

# OWNER'S MANUAL AND SAFETY INSTRUCTIONS

## DEEP TISSUE MASSAGE GUN

ITEM: 97091



**SAVE THIS MANUAL: KEEP THIS MANUAL FOR SAFETY WARNINGS, PRECAUTIONS, ASSEMBLY, OPERATING, INSPECTION, MAINTENANCE AND CLEANING PROCEDURES. WRITE THE PRODUCT'S SERIAL NUMBER ON THE BACK OF THE MANUAL NEAR THE ASSEMBLY DIAGRAM (OR MONTH AND YEAR OF PURCHASE IF PRODUCT HAS NO NUMBER)**

**FOR QUESTIONS PLEASE CALL OUR CUSTOMER SUPPORT: 909.628.0880 MON-FRI 9AM TO 3PM PST**

## WARNING

### GENERAL SAFETY WARNINGS

Read all safety warnings and instructions. Failure to follow the warnings and instructions may result in electric shock, fire and/or serious injury. Save all warnings and instructions for future reference.

The warnings, precautions, and instructions discussed in this instruction manual cannot cover all possible conditions and situations that may occur. It must be understood by the operator that common sense and caution are factors which cannot be built into this product, but must be supplied by the operator. Read carefully and understand all **ASSEMBLY AND OPERATION INSTRUCTIONS** before operating. Failure to follow the safety rules and other basic safety precautions may result in serious personal injury.

- Due to the force of the massage gun, it is recommended to wear clothing on the massage areas. Do not let the massage gun make direct contact with skin.
- Massage time on the same muscle parts should not exceed one minute at a time.
- **DO NOT** use the massage gun on your head nor bony parts of your body.
- Choose a suitable massage head for the specific need. Frequent high speed pressure on the same area can cause abrasions. Stop use immediately if you are experiencing any pain during massage.
- Keep fingers and hair away from the massager to prevent injury.
- Use only the original charger to charge the massage gun.
- Keep the massager away from heat sources.
- Keep the massage gun away from water.
- Inspect the massage gun for any damage prior to each use. If damaged, do not use the massager.
- Modifications of this massage gun is strictly prohibited.
- Do not use high-speed gear or equipment without medical approval if: 1. Pregnancy, Diabetes and Complications, Neuropathy or optic lesions, using a pacemaker, have had a recent surgery, have epilepsy or migraines, disc herniation, relaxation or spondylitis.  
2. Recent joint replacement or IUD or physical health problems. Weak adults and children need to use this machine carefully, accompanied by a healthy adult.
- Consult a physician if you are unsure about using the massage gun.




## SAVE THESE WARNINGS

## BATTERY INFORMATION

1. This machine is equipped with a 12.5V / 2550mAH lithium-ion batteries. This massager can work up to 20 hours continuously.
2. The battery should be charged 4-5 hours when the battery runs out. Never leave the massager charging unattended.
3. Be sure to turn the power off on the massager when not in use.

**WARNING:** The battery in the massager is not removable. Contact your local dealer for a battery replacement.

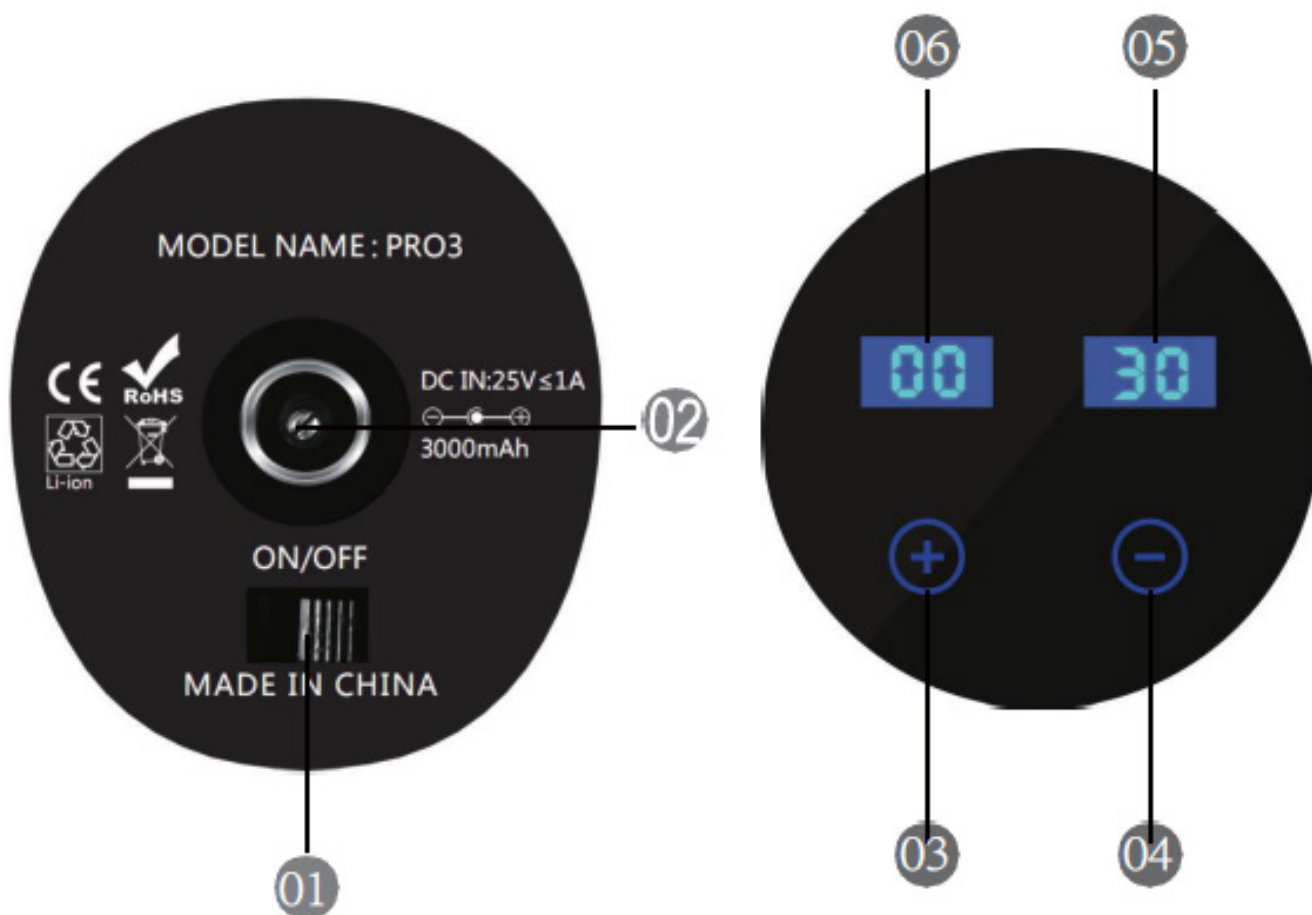

No	Name	No	Name
01	Switch	02	Charge Port
03	Level Plus	04	Level Minus
05	Electricity Indicator	06	Level Indicator ( 1-6)

# FUNCTION INFORMATION

**CHARGING PORT:** When using for the first time, make sure battery is fully charged. Use the original charger, and make sure the machine is not in powered on when charging.

**POWER SWITCH:** Turn the power switch left to power on the massager. Turn the power switch to the right to turn the power off.

**LEVEL PLUS:** Click to increase the massaging power.

**LEVEL MINUS:** Click to decrease the massaging power.

**ELECTRICITY INDICATOR:** Shows the remaining the battery life.

**LEVEL INDICATOR:** **LEVEL 1:** 1500r/min      **LEVEL 2:** 1900r/min      **LEVEL 3:** 2300r/min

**LEVEL 4:** 2700r/min      **LEVEL 5:** 3100r/min      **LEVEL 6:** 3500r/min

**CAUTION:** Levels 5 - 6 are very powerful. Use caution when selecting these two power levels.

## MASSAGER HEADS

Photo	Name	Function
	60MM Globosity head	Suitable for massaging large muscle groups such as neck and thigh of spinal arm
	50MM Globosity head	Suitable for massaging large muscles such as the spine and thighs
	Y head	Suitable for massage of the neck and Achilles tendon
	Small flat head	Suitable for the relaxation and shaping of various muscle parts
	A head	Suitable for various muscle parts such as meridian joint palm and foot base, etc..
	D head	Suitable for the relaxation and shaping of various muscle parts

# MAINTENANCE

1. Wipe the surface of the massage gun with a slightly damp towel. Dry the massager with a soft, dry cloth.
2. When stored for extended periods of time, recharge the battery completely. Keep the massager stored in its original packaging.
3. The massager must be turned off for at least a half an hour after an hour's use. Let the massager cool down for at least 30 minutes before using the machine again.

# DISCLAIMER

## PLEASE READ THE FOLLOWING CAREFULLY

THE MANUFACTURER AND/OR DISTRIBUTOR HAS PROVIDED THE PARTS LIST AND ASSEMBLY DIAGRAM IN THIS MANUAL AS A REFERENCE TOOL ONLY. NEITHER THE MANUFACTURER OR DISTRIBUTOR MAKES ANY REPRESENTATION OR WARRANTY OF ANY KIND TO THE BUYER THAT HE OR SHE IS QUALIFIED TO MAKE ANY REPAIRS TO THE PRODUCT, OR THAT HE OR SHE IS QUALIFIED TO REPLACE ANY PARTS OF THE PRODUCT. IN FACT, THE MANUFACTURER AND/OR DISTRIBUTOR EXPRESSLY STATES THAT ALL REPAIRS AND PARTS REPLACEMENTS SHOULD BE UNDERTAKEN BY CERTIFIED AND LICENSED TECHNICIANS, AND NOT BY THE BUYER. THE BUYER ASSUMES ALL RISK AND LIABILITY ARISING OUT OF HIS OR HER REPAIRS TO THE ORIGINAL PRODUCT OR REPLACEMENT PARTS THERETO, OR ARISING OUT OF HIS OR HER INSTALLATION OF REPLACEMENT PARTS THERETO.

**Record Product's Serial Number Here:** \_\_\_\_\_

**Note:** If product has no serial number, record month and year of purchase instead.

**Note:** Some parts are listed and shown for illustration purposes only and are not available individually as replacement parts.



**SAVE THESE INSTRUCTIONS.**

### Questions, problems, missing parts?



Before returning to your retailer, our exceptional customer service is here to help.

Call Us: 909.628.0880

Hours of Operation: 9am - 4pm (Monday - Friday)

**PRODUCT MADE IN CHINA**