

CAPPUCCINO • SHERMAN

# MICROBIOLOGY

A Laboratory Manual

TENTH EDITION



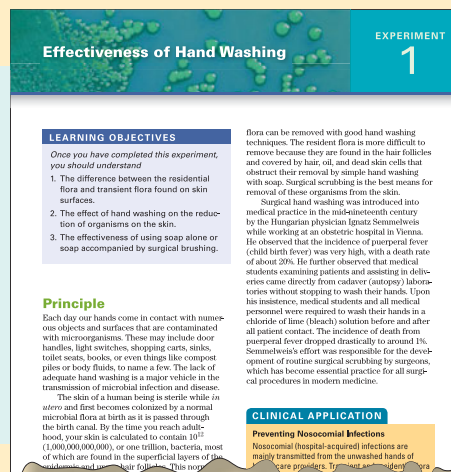
# New Features Make the Micro Lab More Engaging and Applicable to Real Life

Two newly added features in the Tenth Edition of *Microbiology: A Laboratory Manual* speak directly to student interest and engagement: A new experiment on hand washing and a new Clinical Applications feature that appears in each experiment.

NEW!

## Hand-Washing Experiment

A brand-new first experiment emphasizes the importance of effective hand washing and helps introduce microbiology in an engaging way. By learning about the presence of microorganisms on their own hands, this experiment gives students a practical understanding of why aseptic techniques in the lab are so important. A fun, practical lab, this experiment also gets students excited about the course to come.



NEW!

## Clinical Applications

Students often wonder how learning microbiological lab techniques will help them in their future careers. New Clinical Applications added to each experiment make the connection between lab, lecture, and future careers in nursing or other allied health fields.

### CLINICAL APPLICATION

#### Gram Staining: The First Diagnostic

The Gram stain is a diagnostic staining procedure that can be done on body fluids, tissue biopsies, and pure cultures, samples from abscesses when the cause is suspected, and more. Clinically important results are obtained much more rapidly from staining than from culturing the specimen. The results of the Gram stain will aid a clinical lab in determining which tests may be required for identification of the bacterial strain in question. Once the bacterial shape, and orientation are determined, the appropriate choice of antibiotic therapy can be recommended to the patient.

### CLINICAL APPLICATION

#### Cold-Resistant Killers

The field of food science is highly concerned with the temperature-related growth patterns of bacteria. Refrigeration temperatures below 4.4°C are generally considered safe for the short-term storage of food, since most pathogenic bacteria grow very slowly below that temperature. However, some dangerous bacteria are resistant to cold. *Listeria monocytogenes*, which causes a flu-like illness and can be deadly, is capable of doubling its population every 36 hours, even at 4.2°C, and can still attain slow growth below 2°C. *Listeria's* cold tolerance may be due to adaptive genes, prompting research into novel methods of controlling its growth at low temperatures.

### CLINICAL APPLICATION

#### Multiple Drug Therapy

Antibiotic therapy for drug-resistant bacteria and opportunistic pathogen *Pseudomonas aeruginosa* multiple drugs may be used to take advantage of synergistic effects. Research has shown that the use of gentamicin to degrade gram-negative cell wall synthesis. Combination therapies that allow easier entry of kanamycin, which then allows for synergism also allow use of multiple drugs of each drug, which reduces overall toxicity on the patient.

NEW!

## Photos of Protozoa, Yeast, and Fungi

Parts 6 and 7 of the manual now contain even more photos and micrographs of protozoa, yeast, and fungi to help students visualize these important organisms in the lab.

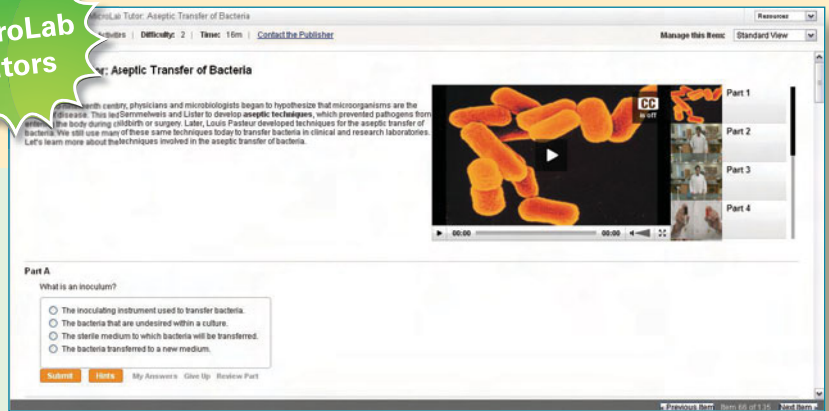


# MasteringMicrobiology® Takes the Microbiology Laboratory Course to the Next Level

**MasteringMicrobiology®** The Mastering online homework, assessment, and tutorial system delivers self-paced tutorials that provide students with individualized coaching set to the instructor's course objectives. MasteringMicrobiology® helps students arrive better prepared for lecture and lab with practice quizzes, coaching activities, tutorials, and more. Instructors can now assign MicroLab Tutors, MicroLab Practicals, Lab Technique Videos, and Lab-specific Quizzes within the MasteringMicrobiology item library. Students can access the MicroLab Tutors, MicroLab Practicals and Lab Technique Videos for review and self-study in the StudyArea of MasteringMicrobiology.


## MicroLab Tutors

**MicroLab Tutors** help instructors and students get the most out of lab time by providing additional opportunities to reinforce lab concepts and techniques. MicroLab Tutors introduce important lab techniques and include each technique's background, purpose, and clinical applications before walking students through the procedure itself. Using a combination of micrographs, video footage in the lab, and, for select topics, 3D molecular animations, MicroLab Tutors help students visualize important lab techniques and processes and make the connection between lecture and lab.




## MicroLab Practical

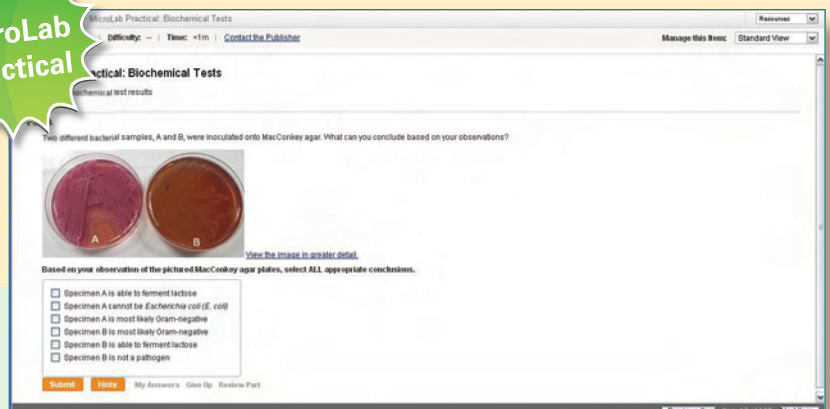**MicroLab Practical** assessments give students extra practice in observing and analyzing test results—putting knowledge into practice. These Lab Practical questions give users an extra opportunity to analyze and interpret important lab tests, procedures, and results in order to get ready for the lab practical!



## Lab Technique Videos

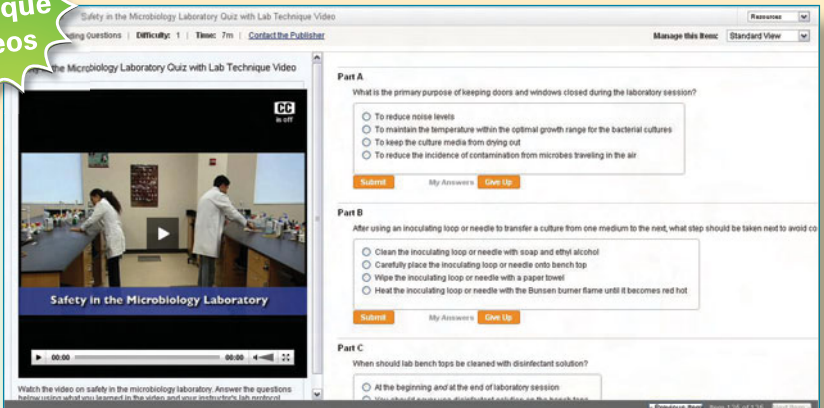**Lab Technique Videos** are 3–5 minute videos, demonstrating specific lab techniques. These videos cover important procedures, such as aseptic technique, Gram staining, and smear preparation. The videos and their associated assessments help students prepare for the wet lab and allow them to review techniques on their own time.



## Lab-Specific Quizzes

**Lab-Specific Quizzes** are newly created for each experiment in the manual. A total of ten multiple-choice and true/false questions per experiment allow instructors to quiz students on microbiological principles, lab techniques, and the theory behind lab procedures. Each quiz question is tied to a Learning Objective for the lab.

## What Instructors are Saying

"These are the types of lab questions I am looking for. Hallelujah!"

► **Pele Rich**  
*North Hennepin Community College*

"I think they are really great!! It puts things in a place where the students can observe firsthand what they should be seeing. I definitely would use them."

► **Tanya Crider**  
*East Mississippi Community College*

# Pearson Custom Library: Lab Manuals Tailored for Your Course

With *The Pearson Custom Library for the Biological Sciences*, instructors can create a lab manual that meets their specific content needs and organization exactly. But that's not all!

Instructors now have the power to **customize their favorite Pearson microbiology lab manuals**. View and select just the labs you need, in the sequence you want, deleting experiments you don't use: Your students pay only for the materials you choose. You can also include your own original content, like your syllabus, course schedule or unique labs. Then preview your custom manual online (with net price and ISBN), and request an evaluation copy to examine before you order. **The Pearson Custom Library is the perfect option to put you in control of your course materials.**

To create the perfect lab manual, visit: [www.pearsoncustomlibrary.com](http://www.pearsoncustomlibrary.com)

## Build Your Manual in Three Easy Steps

In just three easy steps, you can build your ideal customized lab manual that is aligned to your teaching objectives and designed to promote student achievement.

### STEP 1 SEARCH

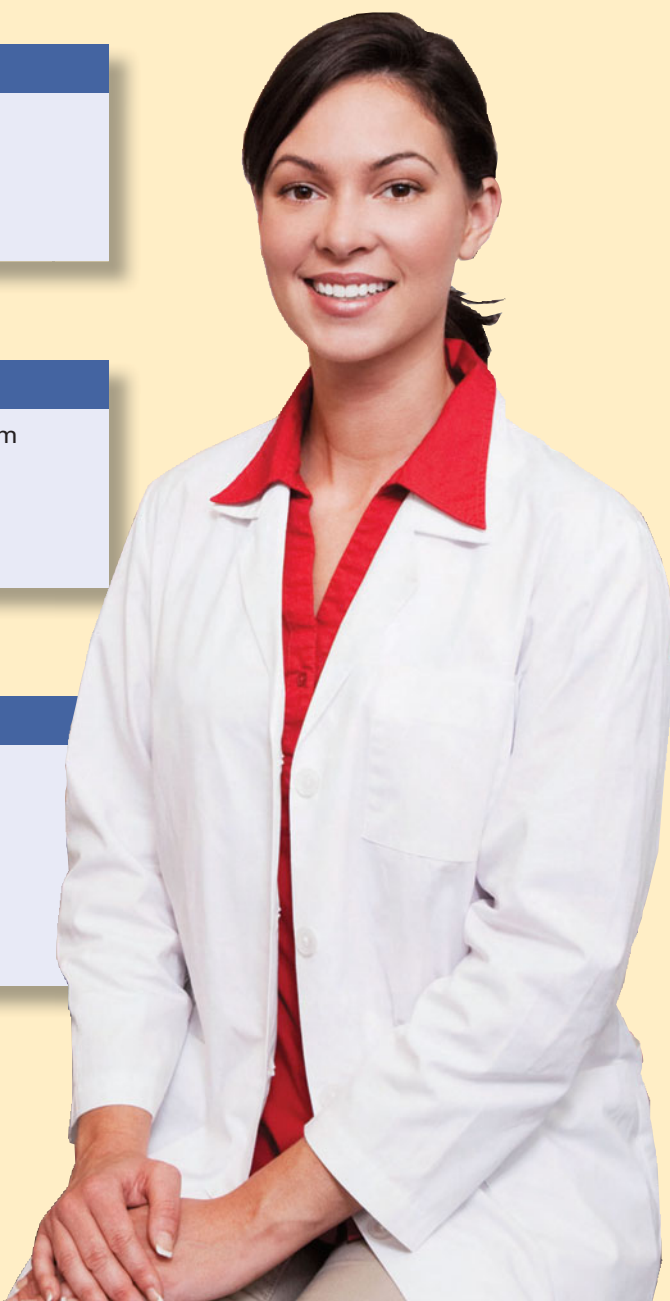
- Search Pearson's vast library of content, including Pearson textbooks, lab manuals, editorially-selected third party readings, and cases.
- Search by author, title, topic, or content type.

### STEP 2 SELECT & PERSONALIZE

- Select only the experiments you use and arrange them in the order you want.
- Upload any original content.
- Choose from a set of diverse and engaging covers.

### STEP 3 PREVIEW & FINALIZE

- Preview Table of Contents
- View a digital review copy instantaneously, or have a bound print copy sent to you.
- Once you decide to adopt, use your automatically generated custom ISBN to place an order at your campus store.





# MICROBIOLOGY

A LABORATORY MANUAL

TENTH EDITION

**James G. Cappuccino**

*SUNY Rockland Community College*

**Natalie Sherman**

*SUNY Rockland Community College*

**PEARSON**

Boston Columbus Indianapolis New York San Francisco Upper Saddle River  
Amsterdam Cape Town Dubai London Madrid Milan Munich Paris Montréal Toronto  
Delhi Mexico City São Paulo Sydney Hong Kong Seoul Singapore Taipei Tokyo

Acquisitions Editor: Kelsey Churchman  
Director of Development: Barbara Yien  
Senior Project Development Editor: Marie Beaugureau  
Editorial Assistant: Ashley Williams  
Managing Editor: Deborah Cogan  
Production Project Manager: Megan Power  
Production Management and Composition: Integra  
Interior Designer: Jerilyn Bockorick  
Cover Designer: Yvo Riezebos  
Image Management: Travis Amos  
Manufacturing Buyer: Stacey Weinberger  
Executive Marketing Manager: Neena Bali  
Text Printer: Webcrafters, Inc.  
Cover Printer: Lehigh-Phoenix Color  
Cover Photo Credit: Andreas Reh / Getty Images

Credits and acknowledgments borrowed from other sources and reproduced, with permission, in this textbook appear on the appropriate page within the text or on p. 531.

Copyright © 2014, 2011, 2008 Pearson Education, Inc. All rights reserved. Manufactured in the United States of America. This publication is protected by Copyright, and permission should be obtained from the publisher prior to any prohibited reproduction, storage in a retrieval system, or transmission in any form or by any means, electronic, mechanical, photocopying, recording, or likewise. To obtain permission(s) to use material from this work, please submit a written request to Pearson Education, Inc., Permissions Department, 1900 E. Lake Ave., Glenview, IL 60025. For information regarding permissions, call (847) 486-2635.

Many of the designations used by manufacturers and sellers to distinguish their products are claimed as trademarks. Where those designations appear in this book, and the publisher was aware of a trademark claim, the designations have been printed in initial caps or all caps.

### Library of Congress

Cappuccino, James G.

Microbiology : a laboratory manual/James G. Cappuccino, Natalie Sherman.—10th ed.  
p.; cm.

Includes bibliographical references and index.

ISBN-13: 978-0-321-84022-6

ISBN-10: 0-321-84022-4

I. Sherman, Natalie. II. Title.

[DNLM: 1. Microbiology—Laboratory Manuals. QW 25]

LC Classification not assigned

579.078—dc23

2012037762

Student edition

ISBN 10: 0-321-84022-4

ISBN 13: 978-0-321-84022-6

Instructor's Review Copy

ISBN 10: 0-321-88451-5

ISBN 13: 978-0-321-88451-0

**PEARSON**

www.pearsonhighered.com

1 2 3 4 5 6 7 8 9 10 11 12—WBC—17 16 15 14 13

# Contents

Preface ix  
Laboratory Safety xiii  
Laboratory Protocol xv

## **PART 1 Basic Laboratory Techniques for Isolation, Cultivation, and Cultural Characterization of Microorganisms 1**

**Introduction 1**  
**Experiment 1:** *Effectiveness of Hand Washing* 7  
**Experiment 2:** *Culture Transfer Techniques* 13  
**Experiment 3:** *Techniques for Isolation of Pure Cultures* 19  
    **Part A:** *Isolation of Discrete Colonies from a Mixed Culture* 19  
    **Part B:** *Isolation of Pure Cultures from a Spread-Plate or Streak-Plate Preparation* 22  
**Experiment 4:** *Cultural Characteristics of Microorganisms* 29

## **PART 2 Microscopy 35**

**Introduction 35**  
**Experiment 5:** *Microscopic Examination of Stained Cell Preparations* 37  
**Experiment 6:** *Microscopic Examination of Living Microorganisms Using a Hanging-Drop Preparation or a Wet Mount* 45  
**Experiment 7:** *The Microscopic Measurement of Microorganisms* 51

## **PART 3 Bacterial Staining 55**

**Introduction 55**  
**Experiment 8:** *Preparation of Bacterial Smears* 59  
**Experiment 9:** *Simple Staining* 65  
**Experiment 10:** *Negative Staining* 71  
**Experiment 11:** *Gram Stain* 75  
**Experiment 12:** *Acid-Fast Stain* 81  
**Experiment 13:** *Differential Staining for Visualization of Bacterial Cell Structures* 87  
    **Part A:** *Spore Stain (Schaeffer-Fulton Method)* 87  
    **Part B:** *Capsule Stain (Anthony Method)* 89

## **PART 4 Cultivation of Microorganisms: Nutritional and Physical Requirements, and Enumeration of Microbial Populations 97**

**Introduction 97**  
**Experiment 14:** *Nutritional Requirements: Media for the Routine Cultivation of Bacteria* 99  
**Experiment 15:** *Use of Differential, Selective, and Enriched Media* 105  
**Experiment 16:** *Physical Factors: Temperature* 115  
**Experiment 17:** *Physical Factors: pH of the Extracellular Environment* 121



**Experiment 18:** *Physical Factors: Atmospheric Oxygen Requirements* 125

**Experiment 19:** *Techniques for the Cultivation of Anaerobic Microorganisms* 131

**Experiment 20:** *Serial Dilution–Agar Plate Procedure to Quantitate Viable Cells* 137

**Experiment 21:** *The Bacterial Growth Curve* 145

## **PART 5 Biochemical Activities of Microorganisms 153**

**Introduction** 153

**Experiment 22:** *Extracellular Enzymatic Activities of Microorganisms* 155

**Experiment 23:** *Carbohydrate Fermentation* 161

**Experiment 24:** *Triple Sugar–Iron Agar Test* 167

**Experiment 25:** *IMViC Test* 171

**Part A:** *Indole Production Test* 172

**Part B:** *Methyl Red Test* 173

**Part C:** *Voges-Proskauer Test* 174

**Part D:** *Citrate Utilization Test* 175

**Experiment 26:** *Hydrogen Sulfide Test* 183

**Experiment 27:** *Urease Test* 187

**Experiment 28:** *Litmus Milk Reactions* 191

**Experiment 29:** *Nitrate Reduction Test* 197

**Experiment 30:** *Catalase Test* 201

**Experiment 31:** *Oxidase Test* 205

**Experiment 32:** *Utilization of Amino Acids* 209

**Part A:** *Decarboxylase Test* 209

**Part B:** *Phenylalanine Deaminase Test* 211

**Experiment 33:** *Genus Identification of Unknown Bacterial Cultures* 215

## **PART 6 The Protozoa 221**

**Introduction** 221

**Experiment 34:** *Free-Living Protozoa* 223

**Experiment 35:** *Parasitic Protozoa* 229

## **PART 7 The Fungi 237**

**Introduction** 237

**Experiment 36:** *Cultivation and Morphology of Molds* 239

**Part A:** *Slide Culture Technique* 239

**Part B:** *Mold Cultivation on Solid Surfaces* 242

**Experiment 37:** *Yeast Morphology, Cultural Characteristics, and Reproduction* 247

**Experiment 38:** *Identification of Unknown Fungi* 255

## **PART 8 The Viruses 261**

**Introduction** 261

**Experiment 39:** *Cultivation and Enumeration of Bacteriophages* 265

**Experiment 40:** *Isolation of Coliphages from Raw Sewage* 271

## **PART 9 Physical and Chemical Agents for the Control of Microbial Growth 277**

**Introduction** 277

**Experiment 41:** *Physical Agents of Control: Moist Heat* 279

**Experiment 42:** *Physical Agents of Control: Environmental Osmotic Pressure* 285

**Experiment 43:** *Physical Agents of Control: Electromagnetic Radiations* 289

**Experiment 44:** *Chemical Agents of Control: Chemotherapeutic Agents* 293

**Part A:** *The Kirby-Bauer Antibiotic Sensitivity Test Procedure* 294

**Part B: Synergistic Effect of Drug Combinations** 296

**Experiment 45: Determination of Penicillin Activity in the Presence and Absence of Penicillinase** 303

**Experiment 46: Chemical Agents of Control: Disinfectants and Antiseptics** 309

**Experiment 55: Bacterial Conjugation** 371

**Experiment 56: Isolation of a Streptomycin-Resistant Mutant** 377

**Experiment 57: The Ames Test: A Bacterial Test System for Chemical Carcinogenicity** 381

## **PART 10 Microbiology of Food** 317

**Introduction** 317

**Experiment 47: Microbiological Analysis of Food Products: Bacterial Count** 319

**Experiment 48: Wine Production** 323

## **PART 11 Microbiology of Water** 327

**Introduction** 327

**Experiment 49: Standard Qualitative Analysis of Water** 329

**Experiment 50: Quantitative Analysis of Water: Membrane Filter Method** 337

## **PART 12 Microbiology of Soil** 343

**Introduction** 343

**Experiment 51: Microbial Populations in Soil: Enumeration** 345

**Experiment 52: Isolation of Antibiotic-Producing Microorganisms and Determination of Antimicrobial Spectrum of Isolates** 351

**Part A: Isolation of Antibiotic-Producing Microorganisms** 352

**Part B: Determination of Antimicrobial Spectrum of Isolates** 353

**Experiment 53: Isolation of Pseudomonas Species by Means of the Enrichment Culture Technique** 357

## **PART 13 Bacterial Genetics** 363

**Introduction** 363

**Experiment 54: Enzyme Induction** 365

## **PART 14 Biotechnology** 387

**Introduction** 387

**Experiment 58: Bacterial Transformation** 389

**Experiment 59: Isolation of Bacterial Plasmids** 397

**Experiment 60: Restriction Analysis and Electrophoretic Separation of Bacteriophage Lambda DNA** 407

## **PART 15 Medical Microbiology** 417

**Introduction** 417

**Experiment 61: Microbial Flora of the Mouth: Determination of Susceptibility to Dental Caries** 419

**Experiment 62: Normal Microbial Flora of the Throat and Skin** 423

**Experiment 63: Identification of Human Staphylococcal Pathogens** 431

**Experiment 64: Identification of Human Streptococcal Pathogens** 439

**Experiment 65: Identification of Streptococcus pneumoniae** 447

**Experiment 66: Identification of Enteric Microorganisms Using Computer-Assisted Multitest Microsystems** 453

**Experiment 67: Isolation and Presumptive Identification of Campylobacter** 463

**Experiment 68: Microbiological Analysis of Urine Specimens** 467

**Experiment 69: Microbiological Analysis of Blood Specimens** 473

**Experiment 70: Species Identification of Unknown Bacterial Cultures** 479

## **PART 16 Immunology 487**

### **Introduction 487**

#### **Experiment 71: Precipitin**

*Reaction: The Ring Test 489*

#### **Experiment 72: Agglutination**

*Reaction: The Febrile Antibody*

*Test 493*

#### **Experiment 73: Immuno- fluorescence 499**

#### **Experiment 74: Enzyme-Linked**

*Immunosorbent Assay 503*

#### **Experiment 75: Sexually**

*Transmitted Diseases: Rapid*

*Immunodiagnostic*

*Procedures 507*

##### **Part A: Rapid Plasma Reagin**

*Test for Syphilis 507*

##### **Part B: Genital Herpes: Isolation**

*and Identification of Herpes*

*Simplex Virus 509*

##### **Part C: Detection of Sexually**

*Transmitted Chlamydial*

*Diseases 510*

## **Appendices**

### **Appendix 1: Scientific**

*Notation 515*

### **Appendix 2: Methods for the**

*Preparation of Dilutions 517*

### **Appendix 3: Microbiological**

*Media 519*

### **Appendix 4: Biochemical Test**

*Reagents 525*

### **Appendix 5: Staining**

*Reagents 528*

### **Appendix 6: Experimental**

*Microorganisms 529*

### **Art & Photo Credits 531**

### **Index 533**

# Dedication



**1932–2001**

*It is with great pride that I dedicate this book to the memory of Natalie Sherman. She was my friend, colleague, and coauthor for 32 years. Her passion for teaching was only exceeded by her ability to teach her students well. They have become the beneficiaries of her unique talent.*

*This page intentionally left blank*

# Preface

Microbiology is a dynamic science. It is constantly evolving as more information is added to the continuum of knowledge, and as microbiological techniques are rapidly modified and refined. The tenth edition of *Microbiology: A Laboratory Manual* continues to provide a blend of traditional methodologies with more contemporary procedures to meet the pedagogical needs of all students studying microbiology. As in previous editions, this laboratory manual provides a wide variety of critically selected and tested experiments suitable for undergraduate students in allied health programs, as well as elementary and advanced general microbiology courses.

## Our Approach

This laboratory manual is designed to guide students in the development of manipulative skills and techniques essential for understanding the biochemical structure and function of a single cell. Its main goal is to encourage students to apply these laboratory skills in the vocational field of applied microbiology and allied health or to further pursue the study of life at the molecular level.

In this manual, comprehensive introductory material is given at the beginning of each major area of study, and specific explanations and detailed directions precede each experiment. This approach augments, enhances, and reinforces course lectures, enabling students to comprehend more readily the concepts and purposes of each experiment. This also provides a review aid if the laboratory and lecture sections are not taught concurrently. The manual should also reduce the time required for explanations at the beginning of each laboratory session and thus make more time available for performing the experiments. Finally, care has been taken to design experimental procedures so that the supplies, equipment, and instrumentation commonly found

in undergraduate institutions will suffice for their successful execution.

## Organization

This manual consists of 75 experiments arranged into 16 parts. The experiments progress from those that are basic and introductory, requiring minimal manipulations, to those that are more complex, requiring more sophisticated skills. The format of each experiment is intended to facilitate presentation of the material by the instructor and to maximize the learning experience. To this end, each experiment is designed as follows:

### Learning Objectives

This introductory section defines the specific principles and/or techniques to be mastered.

### Principle

This is an in-depth discussion of the microbiological concept or technique and the specific experimental procedure.

### Clinical Application

Clinical or medical applications that appear within each experiment help students connect what they are learning in lecture with what they are doing in the lab. For students who intend to have careers as nurses or in other allied health fields, Clinical Applications explain the relevance of each lab technique to their career plans.

### At the Bench

This section signals the beginning of the experiment, and includes the materials, notes of caution, and procedural instructions—all of the things students will need to know at the bench, during the course of the experiment.

## Materials

This comprehensive list helps students and instructors prepare for each laboratory session. Materials appear under one of the following headings:

**Cultures** These are the selected test organisms that have been chosen to demonstrate effectively the experimental principle or technique under study. The choice is also based on their ease of cultivation and maintenance in stock culture. A complete listing of the experimental cultures and prepared slides is presented in Appendix 6.


**Media** These are the specific media and their quantities per designated student group. Appendix 3 lists the composition and method of preparation of all the media used in this manual.

**Reagents** These include biological stains as well as test reagents. The chemical composition and preparation of the reagents are presented in Appendices 4 and 5.

**Equipment** Listed under this heading are the supplies and instrumentation that are needed during the laboratory session. The suggested equipment was selected to minimize expense while reflecting current laboratory technique.

## Procedure

This section provides explicit instructions, augmented by diagrams, that aid in the execution and interpretation of the experiment.

 A caution icon has been placed in experiments that may use **potentially pathogenic materials**. The instructor may wish to perform some of these experiments as demonstrations.

## Lab Report

These tear-out sheets, located at the end of each experiment, facilitate interpretation of data and subsequent review by the instructor. The *Observations and Results* portion of the report provides tables for recording observations and results, and helps the students draw conclusions from and interpret their data. The *Review Questions* aid the instructor in determining the student's ability to understand the experimental concepts and techniques. Questions that call for more critical thinking are indicated by the light bulb icon.



## New to the Tenth Edition

For this tenth edition, the primary aim was to build upon and enrich the student experience. The changes described below are intended to impart the relevance of microbiological lab techniques to student career goals, and to enhance student understanding and performance of each of the microbiological procedures.

### New Clinical Applications

The tenth edition includes new Clinical Applications added to each experiment that help students understand why they're learning individual techniques and connect the lab and lecture to students' future careers as nurses and in other allied health fields. With topics that vary from diagnosing lung infections to selection of effective antibiotics and beyond, clinical applications make the experiments come alive to students.

### New Experiment 1: Effectiveness of Hand Washing

This experiment reinforces the importance of appropriate hand washing and introduces students to microbiology in a fun, easy way. By understanding that microorganisms thrive on their hands, students are also shown the importance of aseptic techniques, both to protect their health and preserve the value of the experiments they are performing.

### New Quiz Questions in *Mastering Microbiology*<sup>®</sup>

Each experiment in the tenth edition is now accompanied by assignable, automatically graded quiz questions in *MasteringMicrobiology*<sup>®</sup>. Each Experiment includes 10 quiz questions correlated to the lab manual Learning Objectives for fast, convenient tracking of course outcomes in the *MasteringMicrobiology* gradebook.

### Additional Protozoa, Fungi, and Yeast Photos Added

Based on reviewer requests, the tenth edition now includes nine new full-color photos and micrographs of protozoa, fungi, and yeast that can be found in Parts 6 and 7 of the manual.

## New “Additional Readings” in the Instructor Guide

Additional reading research articles for each experiment have been added to the Instructor Guide. Articles are relevant to the experiment being performed and are ideal recommendations for further exploration.

## Updates and Revisions

Throughout the manual, updates and revisions have been made to background information, terminology, equipment, and procedural techniques, including the following:

- Added a new procedure and photo for the slide method of the catalase test in Experiment 30.
- Added a new procedure for the filter paper method of the oxidase test in Experiment 31.
- New Observations and Results table in the Lab Report for Experiment 32 covers all experimental organisms.
- Updated photos of *Plasmodium vivax*, *Trypanosoma gambiense*, and *Giardia intestinalis* in Experiment 35.
- Added a new photo showing the minimum inhibitory concentration tube set-up in Experiment 45.
- Updated photo of the development of colonies on a membrane filter plate in Experiment 50.
- Figure 58.3 was altered to agree more precisely with the experimental bacterial transformation procedure.
- Deleted the Agglutination Reaction: Mono-Test for Infectious Mononucleosis from the ninth edition in order to make room for the new hand-washing experiment.

## Instructor Resources

The *Instructor Guide* (ISBN 0-321-86365-8) is a valuable teaching aid for instructors. It has been updated to reflect changes in the main text, and provides the following:

- Laboratory safety protocol for the instructional staff
- Laboratory safety protocol for the technical staff
- New Additional Reading research articles for each experiment

- For each experiment: detailed lists of required materials, procedural points to emphasize, suggestions for optional procedural additions or modifications, helpful tips for preparing or implementing each experiment, and answers to the Review Questions in the student manual
- Appendices with the formulas for the preparation of all media, test reagents, and microbiological stains, as well as the microorganisms required for the performance of each procedure

## Acknowledgments

I wish to express my sincere gratitude to the following instructors for their reviews of the tenth edition. Their comments and direction contributed greatly to the tenth edition.

**Cynthia Anderson**, *Mt. San Antonio College*

**M. Kristie Bompadre**, *Alvernia University*

**Deborah Britt**, *Rhode Island College and Brown University*

**Orianna Carter**, *Ohio University Southern*

**Woo-Suk Chang**, *University of Texas at Arlington*

**Tanya Crider**, *East Mississippi Community College*

**Donald Deters**, *Bowling Green State University*

**Srinivasan Durairaj**, *MacMurray College*

**Sherie Edenborn**, *Chatham University*

**Ralph Fleming**, *Burlington County College*

**Eric Ford**, *East Mississippi Community College*

**Julie Harless**, *Lone Star College–Montgomery*

**Robin R. Hayes**, *Hartnell College*

**Parmeet Jodhka**, *Tarrant County College, Northwest Campus*

**Fred Krick**, *Des Moines Area Community College*

**Jane McLaughlin**, *Orange Coast College*

**Oluwatoyin Osunsanya**, *Muskingum University*

**Adam Parks**, *Shepherd University*

**Diane Pelletier**, *Green River Community College*

**Stacy Pfluger**, *Angelina College*

**Nahed Salama**, *SUNY Rockland Community College*

**George Shiflet**, *Wofford College*

**Amy Sprenkle**, *Salem State University*

**Douglas Stemke**, *University of Indianapolis*

**Chad Welsh**, *Lindenwood University*



I would like to express my sincere gratitude to Dr. Chad Welsh at Lindenwood University and Dr. James Yount at Brevard Community College for reviewing, editing, and writing Clinical Applications for the new edition. Their help made this exciting new feature even stronger. Thanks also to Dr. Welsh for providing the Additional Readings for the Instructor Guide. I am grateful for the work of Dr. Brent Sellinger at the University of Lethbridge. His fine photography and expert experimental setup will be valuable learning tools for students for many years to come.

I wish to extend my appreciation to the staff at Pearson Education, whose expertise and technical skills have guided me over many years. Marie Beaugureau, project editor, provided

invaluable direction in this new edition; her dedication to the highest standards has left its mark on every page of the manual. In addition, I would like to thank Kelsey Churchman, acquisitions editor, Megan Power, production coordinator, and Cheryl Wilson, project manager. I also wish to express my gratitude to the microbiology laboratory technicians at SUNY Rockland Community College—Ms. Joan Grace, who early on performed all the experiments to ensure their success when repeated by the students, and Ms. Roz Wehrman, who followed in Joan's footsteps. Last, but certainly not least, I would like to acknowledge Ms. Allyson Cerreon for her assistance with the hand washing experiment.

*James G. Cappuccino*

## General Rules and Regulations

A rewarding laboratory experience demands strict adherence to prescribed rules for personal and environmental safety. The former reflects concern for your personal safety in terms of avoiding laboratory accidents. The latter requires that you maintain a scrupulously clean laboratory setting to prevent contamination of experimental procedures by microorganisms from exogenous sources.

Because most microbiological laboratory procedures require the use of living organisms, an integral part of all laboratory sessions is the use of aseptic techniques. Although the virulence of microorganisms used in the academic laboratory environment has been greatly diminished because of their long-term maintenance on artificial media, *all microorganisms should be treated as potential pathogens* (organisms capable of producing disease). Thus, microbiology students must develop aseptic techniques (free of contaminating organisms) in the preparation of pure cultures that are essential in the industrial and clinical marketplaces.

The following basic steps should be observed at all times to reduce the ever-present microbial flora of the laboratory environment.

1. Upon entering the laboratory, place coats, books, and other paraphernalia in specified locations—never on bench tops.
2. Keep doors and windows closed during the laboratory session to prevent contamination from air currents.
3. At the beginning and termination of each laboratory session, wipe bench tops with a disinfectant solution provided by the instructor.
4. Do not place contaminated instruments, such as inoculating loops, needles, and pipettes, on bench tops. Loops and needles should be sterilized by incineration, and pipettes should be disposed of in designated receptacles.
5. On completion of the laboratory session, place all cultures and materials in the disposal area as designated by the instructor.
6. Rapid and efficient manipulation of fungal cultures is required to prevent the dissemination of their reproductive spores in the laboratory environment.  
  
To prevent accidental injury and infection of yourself and others, observe the following regulations:
  1. Wash your hands with liquid detergent, rinse with 95% ethyl alcohol, and dry them with paper towels upon entering and prior to leaving the laboratory.
  2. Always use the appropriate safety equipment as determined by your instructor:
    - a. A laboratory coat or apron may be necessary while working in the laboratory. Lab coats protect clothing from contamination or accidental discoloration by staining solutions.
    - b. You may be required to wear gloves while performing the lab exercises. Gloves shield your hands from contamination by microorganisms. They also prevent the hands from coming in direct contact with stains and other reagents.
    - c. Masks and safety goggles may be required to prevent materials from coming in contact with your eyes.
  3. Wear a paper cap or tie back long hair to minimize its exposure to open flames.
  4. Wear closed shoes at all times in the laboratory setting.
  5. Never apply cosmetics or insert contact lenses in the laboratory.
  6. Do not smoke, eat, or drink in the laboratory. These activities are absolutely prohibited.
  7. Carry cultures in a test-tube rack when moving around the laboratory. Likewise, keep cultures in a test-tube rack on the bench tops when not in use. This serves a dual purpose: to prevent accidents and to avoid contamination of yourself and the environment.

8. Never remove media, equipment, or especially, *microbial cultures* from the laboratory. Doing so is absolutely prohibited.
9. Immediately cover spilled cultures or broken culture tubes with paper towels and then saturate them with disinfectant solution. After 15 minutes of reaction time, remove the towels and dispose of them in a manner indicated by the instructor.
10. Report accidental cuts or burns to the instructor immediately.
11. Never pipette by mouth any broth cultures or chemical reagents. Doing so is strictly prohibited. Pipetting is to be carried out with the aid of a mechanical pipetting device only.
12. Do not lick labels. Use only self-stick labels for the identification of experimental cultures.
13. Speak quietly and avoid unnecessary movement around the laboratory to prevent distractions that may cause accidents.

The following specific precautions must be observed when handling body fluids of unknown origin due to the possible transmission of the HIV and hepatitis B viruses in these test specimens.

1. Wear disposable gloves during the manipulation of test materials such as blood, serum, and other body fluids.
2. Immediately wash hands if contact with any of these fluids occurs and also upon removal of the gloves.
3. Wear masks, safety goggles, and laboratory coats if an aerosol might be formed or splattering of these fluids is likely to occur.
4. Decontaminate spilled body fluids with a 1:10 dilution of household bleach, covered with paper toweling, and allowed to react for 10 minutes before removal.
5. Place test specimens and supplies in contact with these fluids into a container of disinfectant prior to autoclaving.

I have read the above laboratory safety rules and regulations and agree to abide by them.

Name: \_\_\_\_\_ Date: \_\_\_\_\_

# Laboratory Protocol

## Student Preparation for Laboratory Sessions

The efficient performance of laboratory exercises mandates that you attend each session fully prepared to execute the required procedures. Read the assigned experimental protocols to effectively plan and organize the related activities. This will allow you to maximize use of laboratory time.

## Preparation of Experimental Materials

**Microscope Slides:** Meticulously clean slides are essential for microscopic work. Commercially precleaned slides should be used for each microscopic slide preparation. However, wipe these slides with dry lens paper to remove dust and finger marks prior to their use. With a glassware marking pencil, label one end of each slide with the abbreviated name of the organism to be viewed.

**Labeling of Culture Vessels:** Generally, microbiological experiments require the use of a number of different test organisms and a variety of culture media. To ensure the successful completion of experiments, organize all experimental cultures and sterile media at the start of each experiment. Label culture vessels with non-water-soluble glassware markers and/or self-stick labels prior to their inoculation. The labeling on each of the experimental vessels should include the name of the test organism, the name of the medium, the dilution of sample (if any), your name or initials, and the date. *Place labeling directly below the cap of the culture tube.* When labeling Petri dish cultures, only the name of the organism(s) should be written on the bottom of the plate, close to its periphery, to prevent obscuring observation of the results. The additional information for the

identification of the culture should be written on the cover of the Petri dish.

## Inoculation Procedures

Aseptic techniques for the transfer or isolation of microorganisms, using the necessary transfer instruments, are described fully in the experiments in Part 1 of the manual. Technical skill will be acquired through repetitive practice.

**Inoculating Loops and Needles:** It is imperative that you incinerate the entire wire to ensure absolute sterilization. The shaft should also be briefly passed through the flame to remove any dust or possible contaminants. To avoid killing the cells and splattering the culture, cool the inoculating wire by tapping the inner surface of the culture tube or the Petri dish cover prior to obtaining the inoculum, or touch the edge of the medium in the plate.

When performing an aseptic transfer of microorganisms, a minute amount of inoculum is required. If an agar culture is used, touch only a single area of growth with the inoculating wire to obtain the inoculum. *Never drag the loop or needle over the entire surface, and take care not to dig into the solid medium.* If a broth medium is used, first tap the bottom of the tube against the palm of your hand to suspend the microorganisms. *Caution: Do not tap the culture vigorously as this may cause spills or excessive foaming of the culture, which may denature the proteins in the medium.*

**Pipettes:** Use only sterile, disposable pipettes or glass pipettes sterilized in a canister. The practice of *pipetting by mouth* has been discontinued to eliminate the possibility of autoinfection by accidentally imbibing the culture or infectious body fluids. Instead, use a mechanical pipetting device to obtain and deliver the material to be inoculated.


## Incubation Procedure

Microorganisms exhibit a wide temperature range for growth. However, for most used in this manual, optimum growth occurs at 37°C over a period of 18 to 24 hours. Unless otherwise indicated in specific exercises, incubate all cultures under the conditions cited above. Place culture tubes in a rack for incubation. Petri dishes may be stacked; however, they *must always be incubated in an inverted position (top down)* to prevent water condensation from dropping onto the surface of the culture medium. This excess moisture could allow the spread of the microorganisms on the surface of the culture medium, producing confluent rather than discrete microbial growth.


## Procedure for Recording Observations and Results

The accurate accumulation of experimental data is essential for the critical interpretation of the observations upon which the final results will be based. To achieve this end, it is imperative that you complete all the preparatory readings that are necessary for your understanding of the basic principles underlying each experiment. Meticulously record all the observed data in the Lab Report of each experiment.

In the experiments that require drawings to illustrate microbial morphology, it will be advantageous to depict shapes, arrangements, and cellular structures enlarged to 5 to 10 times their actual microscopic size, as indicated by the following illustrations. For this purpose a number 2 pencil is preferable. Stippling may be used to depict different aspects of cell structure (e.g., endospores or differences in staining density).

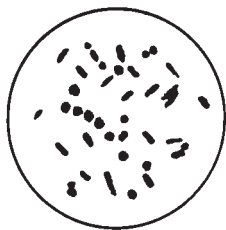

Microscopic drawing

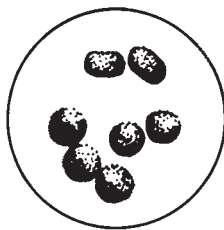

Enlarged drawing

## Review Questions

The review questions are designed to evaluate the student's understanding of the principles and the interpretations of observations in each experiment. Completion of these questions will also serve to reinforce many of the concepts that are discussed in the lectures. At times, this will require the use of ancillary sources such as textbooks, microbiological reviews, or abstracts. The designated critical-thinking questions are designed to stimulate further refinement of cognitive skills.

## Procedure for Termination of Laboratory Sessions

1. Return all equipment, supplies, and chemical reagents to their original locations.
2. Neatly place all capped test tube cultures and closed Petri dishes in a designated collection area in the laboratory for subsequent autoclaving.
3. Place contaminated materials, such as swabs, disposable pipettes, and paper towels, in a biohazard receptacle prior to autoclaving.
4. Carefully place hazardous biochemicals, such as potential carcinogens, into a sealed container and store in a fume hood prior to their disposal according to the institutional policy.
5. Wipe down table tops with recommended disinfectant.
6. Wash hands before leaving the laboratory.

# Basic Laboratory Techniques for Isolation, Cultivation, and Cultural Characterization of Microorganisms

## LEARNING OBJECTIVES

*Once you have completed the experiments in this section, you should be familiar with*

1. The types of laboratory equipment and culture media needed to develop and maintain pure cultures.
2. The types of microbial flora that live on the skin and the effect of hand washing on them.
3. The concept of aseptic technique and the procedures necessary for successful subculturing of microorganisms.
4. Streak-plate and spread-plate inoculation of microorganisms in a mixed microbial population for subsequent pure culture isolation.
5. Cultural and morphological characteristics of microorganisms grown in pure culture.


## Introduction

Microorganisms are ubiquitous. They are found in soil, air, water, food, sewage, and on body surfaces. In short, every area of our environment is replete with them. The microbiologist separates these mixed populations into individual species for study. A culture containing a single unadulterated species of cells is called a **pure culture**. To isolate and study microorganisms in pure culture, the microbiologist requires basic laboratory apparatus and the application of specific techniques, as illustrated in **Figure P1.1**.

## Media

The survival and continued growth of microorganisms depend on an adequate supply of nutrients and a favorable growth environment. For survival, most microbes must use soluble low-molecular-weight substances that are frequently derived from the enzymatic degradation of complex nutrients. A solution containing these

nutrients is a **culture medium**. Basically, all culture media are liquid, semisolid, or solid. A liquid medium lacks a solidifying agent and is called a **broth medium**. A broth medium supplemented with a solidifying agent called **agar** results in a solid or semisolid medium. Agar, an extract of seaweed, is a complex carbohydrate composed mainly of galactose, and is without nutritional value. Agar serves as an excellent solidifying agent because it liquefies at 100°C and solidifies at 40°C. Because of these properties, organisms, especially pathogens, can be cultivated at temperatures of 37.5°C or slightly higher without fear of the medium liquefying. A completely solid medium requires an agar concentration of about 1.5 to 1.8%. A concentration of less than 1% agar results in a **semisolid medium**. A solid medium has the advantage that it presents a hardened surface on which microorganisms can be grown using specialized techniques for the isolation of discrete colonies. Each colony is a cluster of cells that originates from the multiplication of a single cell and

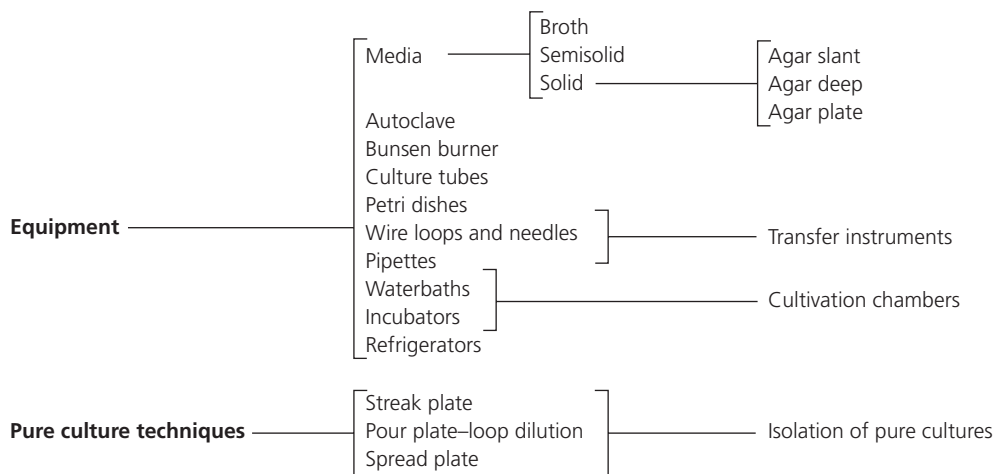


**Figure P1.1** Laboratory apparatus and culture techniques


represents the growth of a single species of microorganism. Such a defined and well-isolated colony is a pure culture. Also, while in the liquefied state, solid media can be placed in test tubes, which are then allowed to cool and harden in a slanted position, producing **agar slants**. These are useful for maintaining pure cultures. Similar tubes that, following preparation, are allowed to harden in the upright position are designated as **agar deep tubes**. Agar deep tubes are used primarily for the study of the gaseous requirements of microorganisms. However, they may be liquefied in a boiling

water bath and poured into Petri dishes, producing **agar plates**, which provide large surface areas for the isolation and study of microorganisms. The various forms of solid media are illustrated in **Figure P1.2**.

In addition to nutritional needs, the environmental factors must also be regulated, including proper pH, temperature, gaseous requirements, and osmotic pressure. A more detailed explanation is presented in Part 4, which deals with cultivation of microorganisms; for now, you should simply bear in mind that numerous types of media are available.

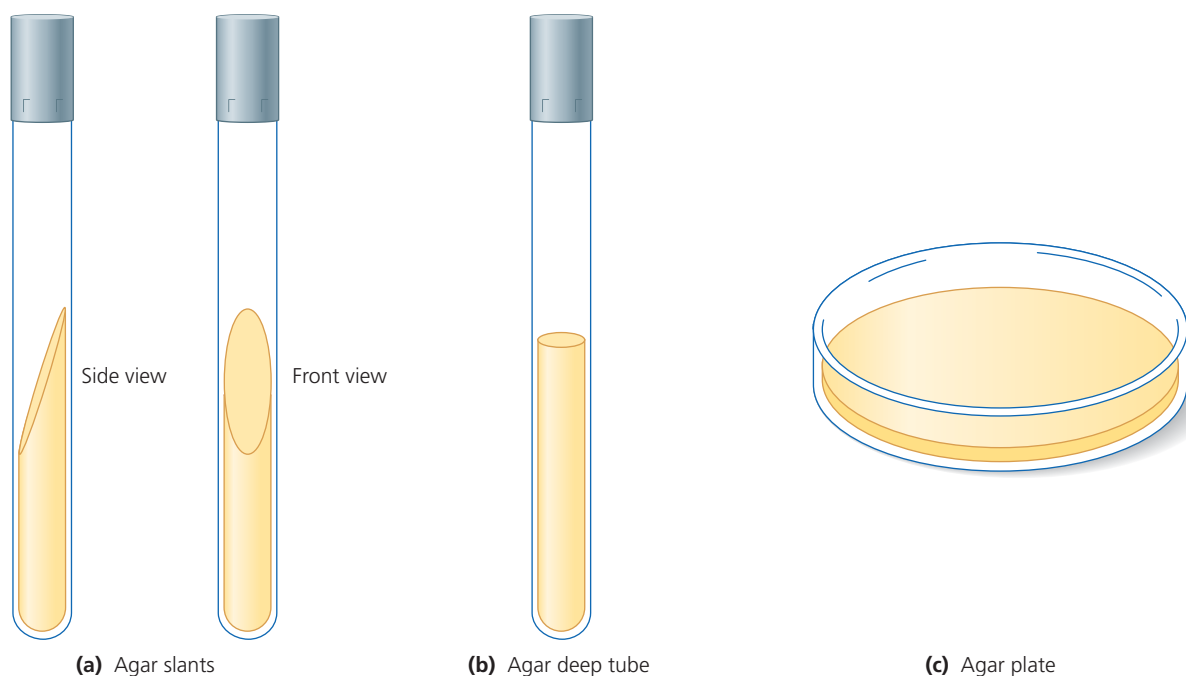

**Figure P1.2** Forms of solid (agar) media

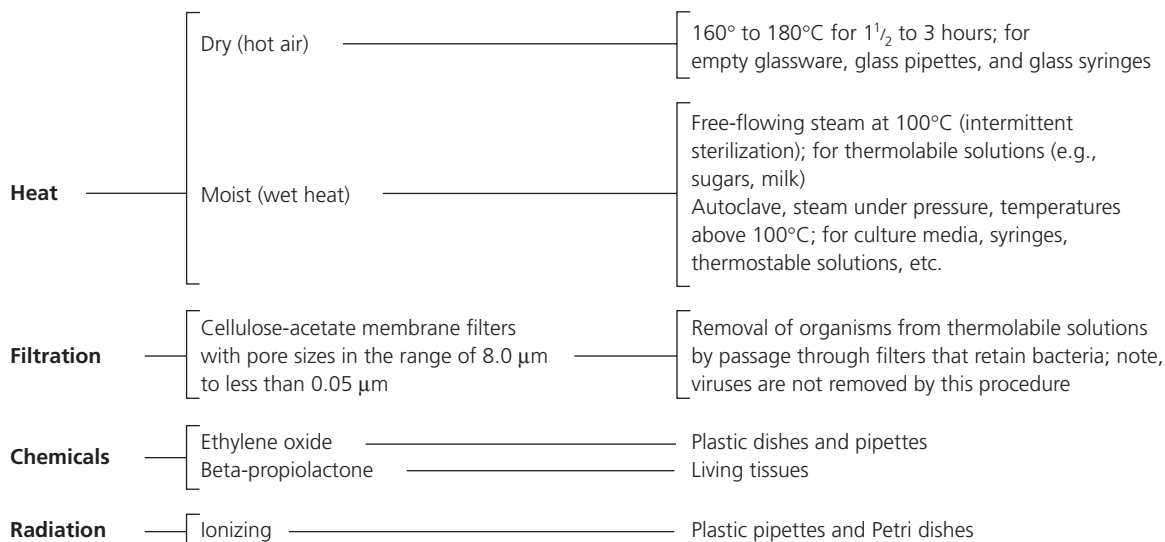

**Figure P1.3 Sterilization techniques**

## Aseptic Technique

Sterility is the hallmark of successful work in the microbiology laboratory, and **sterilization** is the process of rendering a medium or material free of all forms of life. To achieve sterility, it is mandatory that you use sterile equipment and employ **aseptic techniques** when handling bacterial cultures. Although a more detailed discussion is presented in Part 9, which describes the control of microorganisms, **Figure P1.3** is a brief outline of the routine techniques used in the microbiology laboratory.

## Culture Tubes and Petri Dishes

Glass **test tubes** and glass or plastic **Petri dishes** are used to cultivate microorganisms. A suitable nutrient medium in the form of broth or agar may be added to the tubes, while only a solid medium is used in Petri dishes. A sterile environment is maintained in culture tubes by various types of closures. Historically, the first type, a cotton plug, was developed by Schröder and von Dusch in the nineteenth century. Today most laboratories use sleeve-like caps (Morton closures) made of metal, such as stainless steel, or heat-resistant plastics. The advantage of these closures over the cotton plug is that they are labor-saving and, most of all, slip on and off the test tubes easily.

Petri dishes provide a larger surface area for growth and cultivation. They consist of a bottom dish portion that contains the medium and a larger top portion that serves as a loose cover.

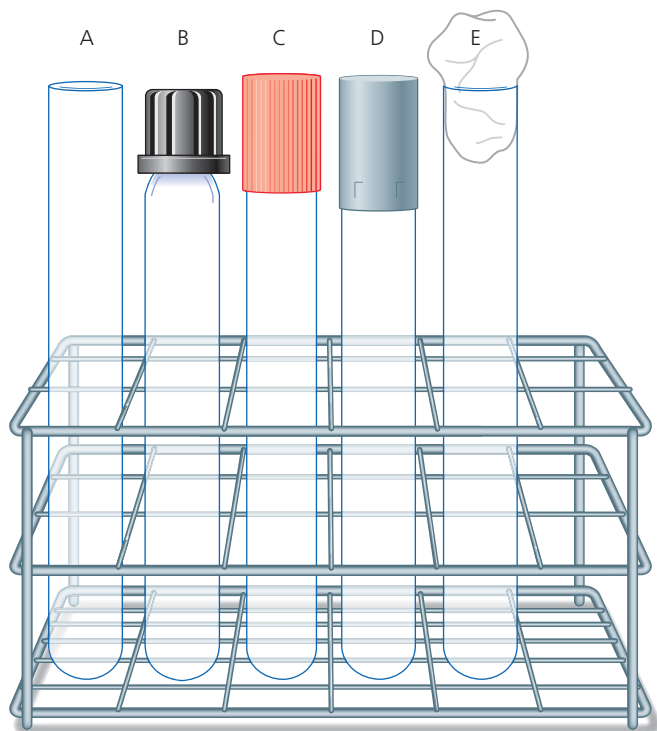
Petri dishes are manufactured in various sizes to meet different experimental requirements. For routine purposes, dishes approximately 15 cm in diameter are used. The sterile agar medium is dispensed to previously sterilized dishes from molten agar deep tubes containing 15 to 20 ml of medium, or from a molten sterile medium prepared in bulk and contained in 250-, 500-, and 1000-ml flasks, depending on the volume of medium required. When cooled to 40°C, the medium will solidify. Remember that *after inoculation, Petri dishes are incubated in an inverted position* (top down) to prevent condensation formed on the cover during solidification from dropping down onto the surface of the hardened agar. **Figure P1.4** illustrates some of the culture vessels used in the laboratory. Built-in ridges on tube closures and Petri dishes provide small gaps necessary for the exchange of air.

## Transfer Instruments

Microorganisms must be transferred from one vessel to another or from stock cultures to various media for maintenance and study. Such a transfer is called **subculturing** and must be carried out under aseptic conditions to prevent possible contamination.


**Wire loops** and needles are made from inert metals such as Nichrome or platinum and are inserted into metal shafts that serve as handles. They are extremely durable instruments and are easily sterilized by incineration in the blue (hottest) portion of the Bunsen burner flame.






- A. Bacteriological tube      D. Metal closure  
 B. Screw cap                      E. Nonabsorbent cotton  
 C. Plastic closure

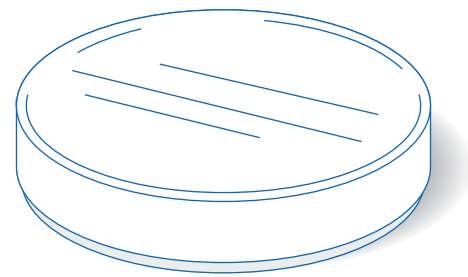(a) Test tube rack with tubes showing various closures



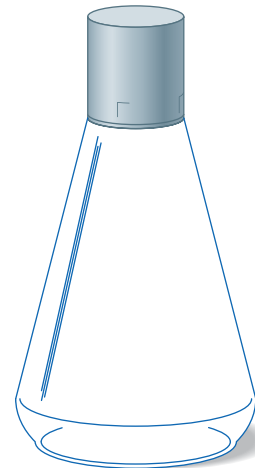(b) Petri dish



(c) DeLong shaker flask with closure

**Figure P1.4 Culture vessels**

A **pipette** is another instrument used for aseptic transfers. Pipettes are similar in function to straws; that is, they draw up liquids. They are made of glass or plastic drawn out to a tip at one end and with a mouthpiece forming the other end. They are calibrated to deliver different volumes depending on requirements. Pipettes may be sterilized in bulk inside canisters, or they may be wrapped individually in brown paper and sterilized in an autoclave or dry-heat oven.


**Figure P1.5** illustrates these transfer instruments. The proper procedure for the use of pipettes will be demonstrated by your instructor.

**⚠ Pipetting by mouth is not permissible!**  
**Pipetting is to be performed with the aid of mechanical pipette aspirators.**

## Cultivation Chambers

The specific temperature requirements for growth are discussed in detail in Part 4. However, a prime requirement for the cultivation of microorganisms is that they be grown at their optimum temperature. An incubator is used to maintain optimum temperature during the necessary growth period. It resembles an oven and is thermostatically controlled so that temperature can be varied depending on the requirements of specific microorganisms. Most incubators use dry heat. Moisture is supplied by placing a beaker of water in the incubator during the growth period. A moist environment retards dehydration of the medium and thereby avoids misleading experimental results.

A thermostatically controlled **shaking waterbath** is another piece of apparatus used to

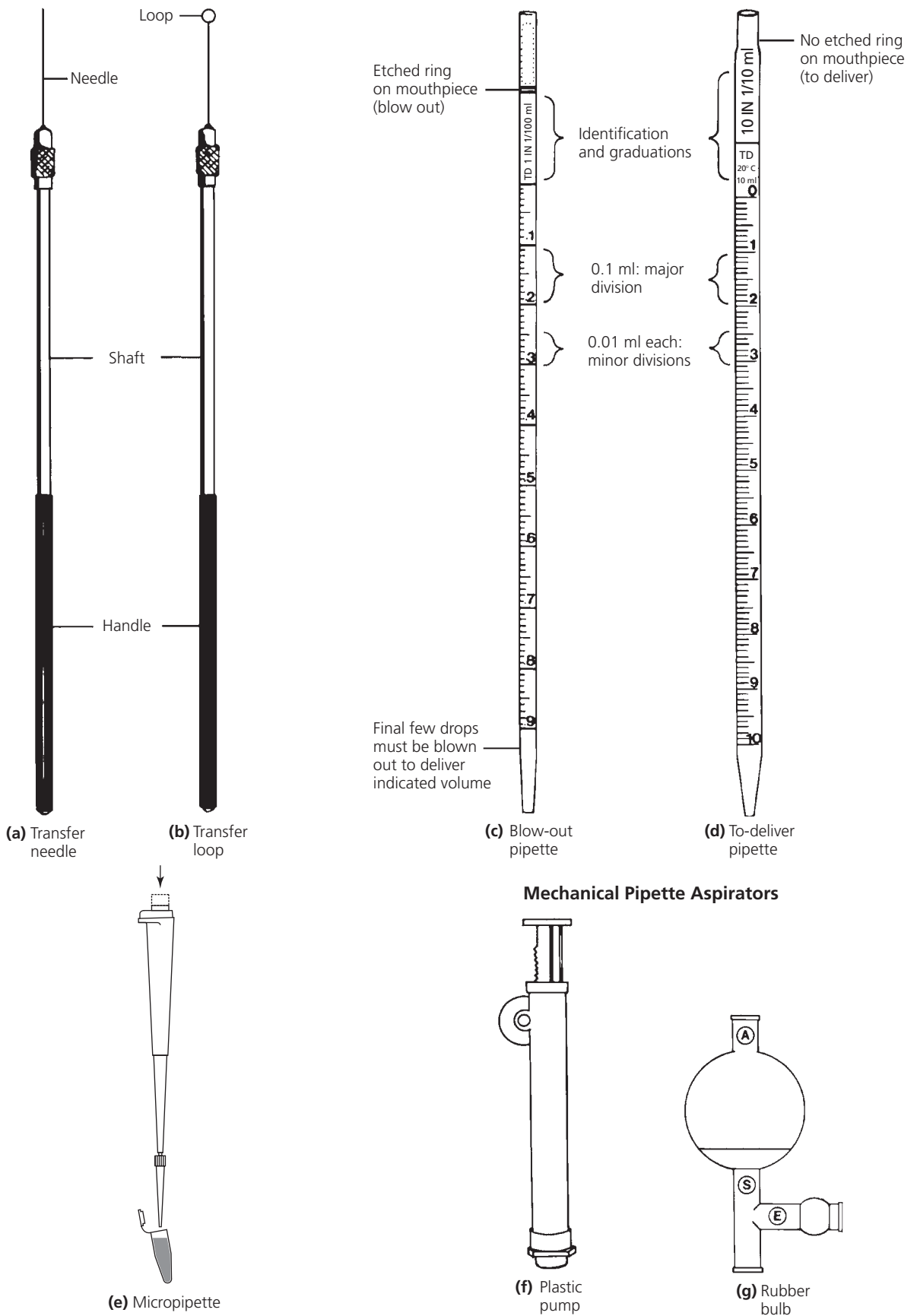

**Figure P1.5** Transfer instruments

cultivate microorganisms. Its advantage is that it provides a rapid and uniform transfer of heat to the culture vessel, and its agitation provides increased aeration, resulting in acceleration of growth. The single disadvantage of this instrument is that it can be used only for cultivation of organisms in a broth medium.

## Refrigerator

A refrigerator is used for a wide variety of purposes such as maintenance and storage of stock cultures between subculturing periods, and storage of sterile media to prevent dehydration. It is also used as a repository for thermolabile solutions, antibiotics, serums, and biochemical reagents.

### LEARNING OBJECTIVES

Once you have completed this experiment, you should understand

1. The difference between the residential flora and transient flora found on skin surfaces.
2. The effect of hand washing on the reduction of organisms on the skin.
3. The effectiveness of using soap alone or soap accompanied by surgical brushing.

## Principle

Each day our hands come in contact with numerous objects and surfaces that are contaminated with microorganisms. These may include door handles, light switches, shopping carts, sinks, toilet seats, books, or even things like compost piles or body fluids, to name a few. The lack of adequate hand washing is a major vehicle in the transmission of microbial infection and disease.

The skin of a human being is sterile while *in utero* and first becomes colonized by a normal microbial flora at birth as it is passed through the birth canal. By the time you reach adulthood, your skin is calculated to contain  $10^{12}$  (1,000,000,000,000), or one trillion, bacteria, most of which are found in the superficial layers of the epidermis and upper hair follicles. This normal flora of microorganisms is called the **resident flora**, the presence of which does not cause negative effects in healthy individuals. In fact, it forms a symbiotic relationship with your skin, which is vital to your health. This beneficial relationship can change in patients who are immunocompromised, or when residential flora accidentally gains entrance to the host via inoculating needles, indwelling catheters, lacerations, and the like. Microorganisms that are less permanent and present for only short periods are termed **transient flora**. This latter

flora can be removed with good hand washing techniques. The resident flora is more difficult to remove because they are found in the hair follicles and covered by hair, oil, and dead skin cells that obstruct their removal by simple hand washing with soap. Surgical scrubbing is the best means for removal of these organisms from the skin.

Surgical hand washing was introduced into medical practice in the mid-nineteenth century by the Hungarian physician Ignatz Semmelweis while working at an obstetric hospital in Vienna. He observed that the incidence of puerperal fever (child birth fever) was very high, with a death rate of about 20%. He further observed that medical students examining patients and assisting in deliveries came directly from cadaver (autopsy) laboratories without stopping to wash their hands. Upon his insistence, medical students and all medical personnel were required to wash their hands in a chloride of lime (bleach) solution before and after all patient contact. The incidence of death from puerperal fever dropped drastically to around 1%. Semmelweis's effort was responsible for the development of routine surgical scrubbing by surgeons, which has become essential practice for all surgical procedures in modern medicine.

### CLINICAL APPLICATION

#### Preventing Nosocomial Infections

Nosocomial (hospital-acquired) infections are mainly transmitted from the unwashed hands of health care providers. Transient and residential flora on health care providers' skin can infect hospital patients whose immune systems are compromised. The cornerstone for the prevention of nosocomial infections is the meticulous hand washing and scrubbing of health care personnel. In the laboratory setting, your normal flora may contaminate patient samples and skew your result, leading to a misdiagnosis. It is important for everyone in the lab to correctly wash their hands before and after handling biological materials.



## Materials

### Media

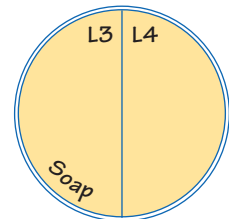
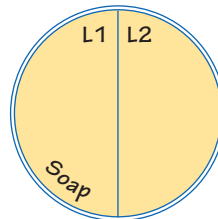
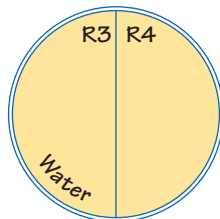
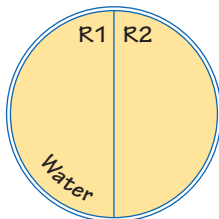
4 nutrient agar plates per student pair

### Equipment

Liquid antibacterial soap, 8 sterile cotton swabs, 2 test tubes of sterile saline, Bunsen burner, glass marking pencil, surgical hand brush, Quebec colony counter, stopwatch.

## Procedure Lab One

- One student will become the washer and the other student the assistant. **The washer must not wash hands before coming to the lab.**
- The assistant will use the glass marking pencil to label the bottoms of the nutrient agar plates. The assistant will mark two plates as “Water” and two plates as “Soap” and draw a line down the middle of each plate to divide each plate in half. For the “Water” plates, label the halves as R1, R2, R3, and R4. For the “Soap” plates, label the halves as L1, L2, L3, and L4. See **Figure 1.1**.
- The assistant will aseptically dip a sterile cotton swab into the first test tube of sterile saline. To do this:
  - First light the Bunsen burner.
  - Uncap the test tube; after removing the cap, keep the cap in your hand with the inner aspect of the cap pointed away from your palm. The cap must never be placed on the laboratory bench because doing so would compromise the aseptic procedure.
  - Flame the neck of the tube by briefly passing it through the flame of the Bunsen burner.
  - Remove the tube from the flame and dip the swab in the tube, soaking it with saline.

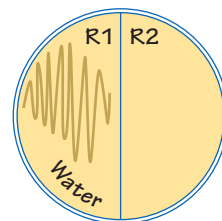


**Figure 1.1** Plate labeling

Avoid touching the sides of the tube with the swab.

The assistant will then rub the moistened cotton swab on the pad of the washer’s **right** thumb.

- The assistant will then aseptically inoculate the half of the nutrient agar plate labeled R1 by streaking the far edge of the plate several times then making a zig zag streak only on the half labeled R1. See **Figure 1.2**. *Caution: Do not gouge the surface of the agar plate.*
- The assistant will turn on the tap on the lab sink, so that the washer can wash the right hand under warm running water, **without soap**, concentrating on the thumb (rubbing the thumb over the right index and middle finger) for one minute. The assistant will turn off the tap. The washer will shake off the excess water from the hand, but not blot dry. The assistant, using a new, dry (not moistened with saline) sterile cotton swab, will obtain a sample from the right thumb pad and inoculate the section of the nutrient agar plate labeled R2 in the same way that R1 was inoculated.
- Repeat step 5 two more times, washing the thumb for 2 minutes and then 3 minutes, respectively. The assistant will use a new, dry sterile cotton swab each time, and will aseptically inoculate R3 and R4, respectively. See Table 1.1.
- The assistant and washer will now move to the left hand. The assistant will aseptically dip the



**Figure 1.2** Plate inoculation

**TABLE 1.1 Inoculation of Nutrient Agar Plates**

WATER—RIGHT THUMB		SOAP—LEFT THUMB	
R1	No wash, damp cotton swab	L1	No wash, damp cotton swab
R2	Wash 1 minute, dry cotton swab	L2	Wash with soap 1 minute, dry cotton swab
R3	Wash 2 minutes, dry cotton swab	L3	Soap and surgical brush 2 minutes, dry cotton swab
R4	Wash 3 minutes, dry cotton swab	L4	Soap and surgical brush 3 minutes, dry cotton swab

sterile cotton swab into the second test tube of sterile saline (following the process from Step 3) and will rub the moistened cotton swab over the pad of the left thumb and aseptically inoculate L1 as shown in Figure 1.2.

8. The assistant will turn on the tap of the lab's sink so that the washer can wet the thumb and index finger of the left hand under warm running water. The assistant will apply one or two drops of liquid soap to the thumb and index finger and the washer will wash for 1 minute by rubbing the thumb over the index finger. Rinse well. Shake off water from the hand but do not blot dry. The assistant will turn off the tap. The assistant will then use a dry, sterile cotton swab to obtain a sample from the washed thumb pad and inoculate L2.
9. Repeat step 8 two more times, not only using soap but also scrubbing the thumb with a surgical brush, for 2 minutes and then 3 minutes, respectively. The washer will obtain the surgical brush and the assistant will add saline to the brush to dampen it, and then add one or two drops of soap to the thumb and also the brush. *Caution: Place the brush bristles up on a dry paper towel between washings.* The assistant will use a new, dry sterile cotton swab each time, and will aseptically inoculate L3 and L4, respectively. Refer back to Table 1.1.
10. Incubate all plates in an inverted position at 37°C for 24 to 48 hours.

## Procedure Lab Two

Examine and record the amount of growth found on each nutrient agar plate. Results may be determined by two methods.

1. **Macroscopically.** Visually observe the presence of growth on the surface of each agar plate in each section. Record your results in your Lab Report as 0 = no growth, 1+ = slight growth, 2+ = moderate growth, 3+ = heavy growth, and 4+ = maximum growth.
2. **Percent Growth Reduction.**
  - a. Count the colonies that appear in each section of the agar plates using a Quebec colony counter. If more than 300 colonies are present, label it as “too numerous to count (TNTC),” if fewer than 30 colonies are present, label it as “too few to count (TFTC).”
  - b. For sections R2, R3, R4 and L2, L3, L4, calculate the percent growth reduction from the first section, using the following equation:

$$\text{Percent reduction} = \frac{[\text{Colonies}(\text{section 1}) - \text{Colonies}(\text{section x})] \div \text{Colonies}(\text{section 1})}{100}$$

X = sections 2, 3, 4 for each hand

*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

- Record the macroscopic observations in the chart below.

Section (Water—Right Thumb)	Time (min)	Growth (0 = none, 1+ = slight, 2+ = moderate, 3+ = heavy, 4+ = maximum)	Section (Soap—Left Thumb)	Time (min)	Growth (0 = none, 1+ = slight, 2+ = moderate, 3+ = heavy, 4+ = maximum)
R1	0		L1	0	
R2	1		L2	1	
R3	2		L3	2	
R4	3		L4	3	

- Record the percent growth reduction in the chart below.

Section (Water— Right Thumb)	Time (min)	Number of Colonies	Percent Reduction	Section (Soap—Left Thumb)	Time (min)	Number of Colonies	Percent Reduction
R1	0		—	L1	0		—
R2	1			L2	1		
R3	2			L3	2		
R4	3			L4	3		


## Review Questions

- Compare the effectiveness of hand washing with water, with soap, and with soap and surgical scrubbing.



2. How does the presence of residential flora influence the infectious process?

3. How does hand washing affect residential versus transient flora?

4.  Why do you think hand washing is necessary when medical and surgical personnel wear gloves during surgery or when examining patients?

### LEARNING OBJECTIVES

Once you have completed this experiment, you should be able to

1. Carry out the technique for aseptic removal and transfer of microorganisms for subculturing.
2. Correctly sterilize inoculating instruments in the flame of a Bunsen burner.
3. Correctly manipulate your fingers to remove and replace the test tube closure.

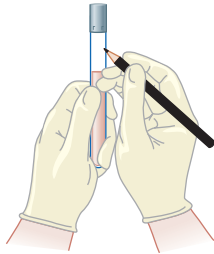
### Principle

Microorganisms are transferred from one medium to another by **subculturing**. This technique is of basic importance and is used routinely in preparing and maintaining stock cultures, as well as in microbiological test procedures.


Microorganisms are always present in the air and on laboratory surfaces, benches, and equipment. They can serve as a source of external contamination and thus interfere with experimental results unless proper aseptic techniques are used during subculturing. Described below are essential steps that you must follow for aseptic transfer of microorganisms. The complete procedure is illustrated in **Figure 2.1**.

1. Label the tube to be inoculated with the name of the organism and your initials.
2. Hold the stock culture tube and the tube to be inoculated in the palm of your hand, secure with your thumb, and separate the two tubes to form a V in your hand.
3. Sterilize an inoculating needle or loop by holding it in the hottest portion of the Bunsen burner flame, until the wire becomes red hot. Then, rapidly pass the upper portion of the handle through the flame. Once flamed, the loop is never put down but is held in the hand and allowed to cool for 10 to 20 seconds.
4. Uncap the tubes by grasping the first cap with your little finger and the second cap with your next finger and lifting the closure upward.  
*Note: Once removed, these caps must be kept in the hand that holds the sterile inoculating loop or needle; thus, the inner aspects of the caps point away from the palm of the hand. They must never be placed on the laboratory bench because doing so would compromise the aseptic procedure.*
5. After removing the closures, flame the necks and mouths of the tubes by briefly passing them through the flame two–three times rapidly. The sterile transfer instrument is further cooled by touching it to the sterile inside wall of the culture tube before removing a small sample of the inoculum.
6. Depending on the culture medium, a loop or needle is used for removal of the inoculum. Loops are commonly used to obtain a sample from a broth culture. Either instrument can be used to obtain the inoculum from an agar slant culture by carefully touching the surface of the solid medium in an area exhibiting growth so as not to gouge the agar. A straight needle is always used when transferring microorganisms to an agar deep tube from both solid and liquid cultures.
  - a. For a slant-to-broth transfer, obtain inoculum from the slant and lightly shake the loop or needle in the broth culture to dislodge the microorganisms.
  - b. For a broth-to-slant transfer, obtain a loopful of broth and place at the base of an agar slant medium. Lightly draw the loop over the hardened surface in a straight or zigzag line, from the base of the agar slant to the top.
  - c. For a slant-to-agar deep transfer, obtain the inoculum from the agar slant. Insert a straight needle to the bottom of the tube in a straight line and rapidly withdraw along the line of insertion. This is called a stab inoculation.
7. Following inoculation, remove the instrument and re flame the necks of the tubes.


## PROCEDURE



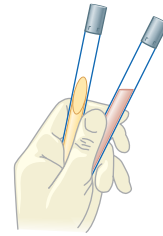
- 1 Label the tube to be inoculated with the name of the organism and your initials.



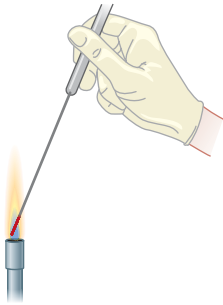
- 2 Place the tubes in the palm of your hand, secure with your thumb, and separate to form a V.



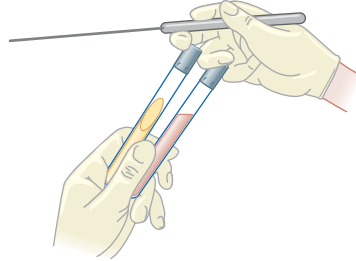
- 3 Flame the needle or loop until the wire is red.



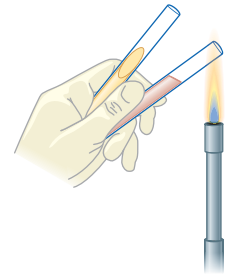
- 4 With the sterile loop or needle in hand, uncapped the tubes.



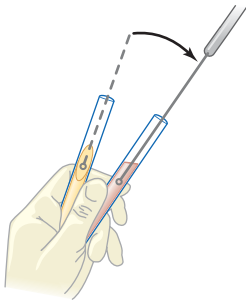
- 5 Flame the necks of the tubes by rapidly passing them through the flame once.



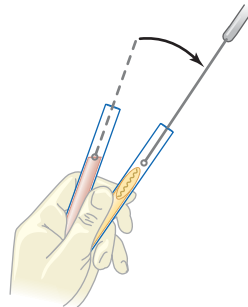
- 6 **Slant-to-broth transfer:** Obtain inoculum from slant and dislodge inoculum in the broth with a slight agitation.



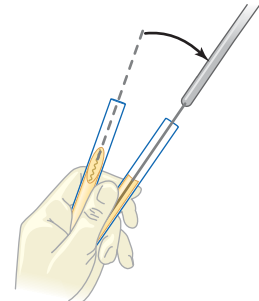
- Broth-to-slant transfer:** Obtain a loopful of broth and place at base of slant. Withdraw the loop in a zigzag motion.



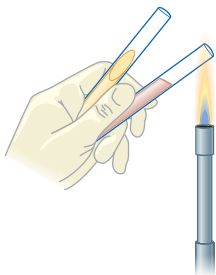
- Slant-to-agar deep transfer:** Obtain inoculum from slant. Insert the needle to the bottom of the tube and withdraw along the line of insertion.



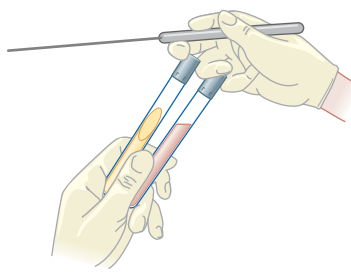- 7 Flame the necks of the tubes by rapidly passing them through the flame once.



- 8 Recap the tubes.

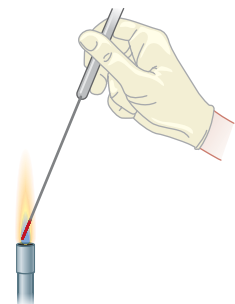

- 9 Re-flame the loop or needle.

Figure 2.1 Subculturing procedure

8. Replace the caps on the same tubes from which they were removed.
9. Reflame the loop or needle to destroy any remaining organisms.

In this experiment you will master the manipulations required for aseptic transfer of microorganisms in broth-to-slant, slant-to-broth, and slant-to-agar deep transfers. The technique for transfer to and from agar plates is discussed in Experiment 3.

## CLINICAL APPLICATION

### Aseptic Inoculation and Transfer

It's mandatory for those working in a microbiology laboratory to learn and perfect the skill of inoculating bacterial specimens on agar plates, in liquid broth, or in semisolid medium, and subsequently be able to subculture the organism from one medium to another. A sterile inoculating needle or loop is the basic instrument of transfer. It is important that you keep in mind that transferring bacterial cultures requires aseptic or sterile techniques at all times, especially if you are working with pathogens. In short, do not contaminate what you are working with and do not contaminate yourself.

## AT THE BENCH



## Materials

### Cultures

24-hour nutrient broth and nutrient agar slant cultures of *Serratia marcescens*.

## Media

Per designated student group: one nutrient broth, one nutrient agar slant, and one nutrient agar deep tube.

## Equipment

Bunsen burner, inoculating loop and needle, and glassware marking pencil.

## Procedure Lab One

1. Label all tubes of sterile media as described in the Laboratory Protocol section on page xv.
2. Following the procedure outlined and illustrated previously (Figure 2.1), perform the following transfers:
  - a. *S. marcescens* broth culture to a nutrient agar slant, nutrient agar deep tube, and nutrient broth.
  - b. *S. marcescens* agar slant culture to a nutrient broth, nutrient agar slant, and nutrient agar deep tube.
3. Incubate all cultures at 25°C for 24 to 48 hours.

## Procedure Lab Two


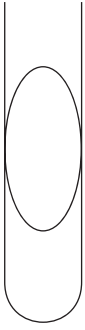

1. Examine all cultures for the appearance of growth, which is indicated by turbidity in the broth culture and the appearance of an orange-red growth on the surface of the slant and along the line of inoculation in the agar deep tube.
2. Record your observations in the chart provided in the Lab Report.

*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

	Nutrient Broth	Nutrient Agar Slant	Nutrient Agar Deep
Growth (+) or (-)	_____	_____	_____
Orange-red pigmentation (+) or (-)	_____	_____	_____
Draw the distribution of growth.			

## Review Questions


1. Explain why the following steps are essential during subculturing:
  - a. Flaming the inoculating instrument *prior to and after* each inoculation.
  - b. Holding the test tube caps in the hand as illustrated in Figure 2.1 on page 14.


c. Cooling the inoculating instrument prior to obtaining the inoculum.

d. Flaming the neck of the tubes immediately after uncapping and before recapping.

2. Describe the purposes of the subculturing procedure.

3. Explain why a straight inoculating needle is used to inoculate an agar deep tube.

4.  There is a lack of orange-red pigmentation in some of the growth on your agar slant labeled *S. marcescens*. Does this necessarily indicate the presence of a contaminant? Explain.

5.  Upon observation of the nutrient agar slant culture, you strongly suspect that the culture is contaminated. Outline the method you would follow to ascertain whether your suspicion is justified.

In nature, microbial populations do not segregate themselves by species but exist with a mixture of many other cell types. In the laboratory, these populations can be separated into **pure cultures**. These cultures contain only one type of organism and are suitable for the study of their cultural, morphological, and biochemical properties.

In this experiment, you will first use one of the techniques designed to produce discrete colonies. Colonies are individual, macroscopically visible masses of microbial growth on a solid medium surface, each representing the multiplication of a single organism. Once you have obtained these discrete colonies, you will make an aseptic transfer onto nutrient agar slants for the isolation of pure cultures.

## PART A Isolation of Discrete Colonies from a Mixed Culture

### LEARNING OBJECTIVE


*Once you have completed this experiment, you should be able to*

1. Perform the streak-plate and/or the spread-plate inoculation procedure to separate the cells of a mixed culture so that discrete colonies can be isolated.

## Principle


The techniques commonly used for isolation of discrete colonies initially require that the number of organisms in the inoculum be reduced. The resulting diminution of the population size ensures that, following inoculation, individual cells will be sufficiently far apart on the surface of the agar medium to separate the different species. The following are techniques that can be used to accomplish this necessary dilution:

1. The **streak-plate** method is a rapid qualitative isolation method. It is essentially a dilution technique that involves spreading a loopful of culture over the surface of an agar plate. Although many types of procedures are performed, the four-way, or quadrant, streak is described. Refer to **Figure 3.1**, which schematically illustrates this technique.
  - a. Place a loopful of culture on the agar surface in Area 1. Flame the loop, cool it by touching an unused part of the agar surface close to the periphery of the plate, and then drag it rapidly several times across the surface of Area 1.
  - b. Reflame and cool the loop, and turn the Petri dish 90°. Then touch the loop to a corner of the culture in Area 1 and drag it several times across the agar in Area 2. The loop should never enter Area 1 again.

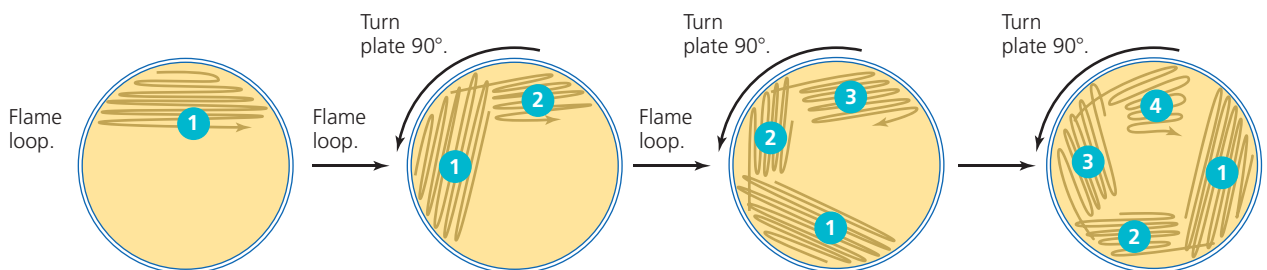


**Figure 3.1** Four-way streak-plate technique



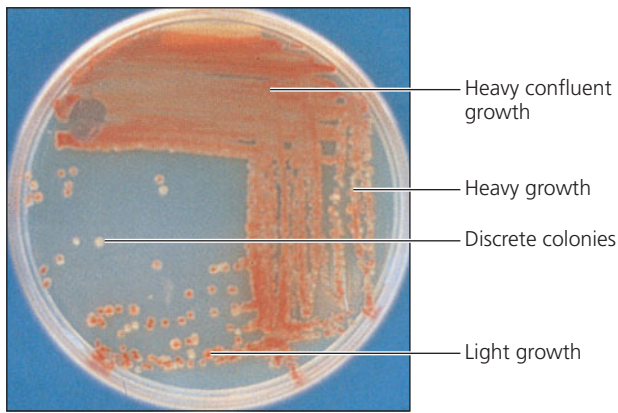

**Figure 3.2** Four-way streak-plate inoculation with *Serratia marcescens*

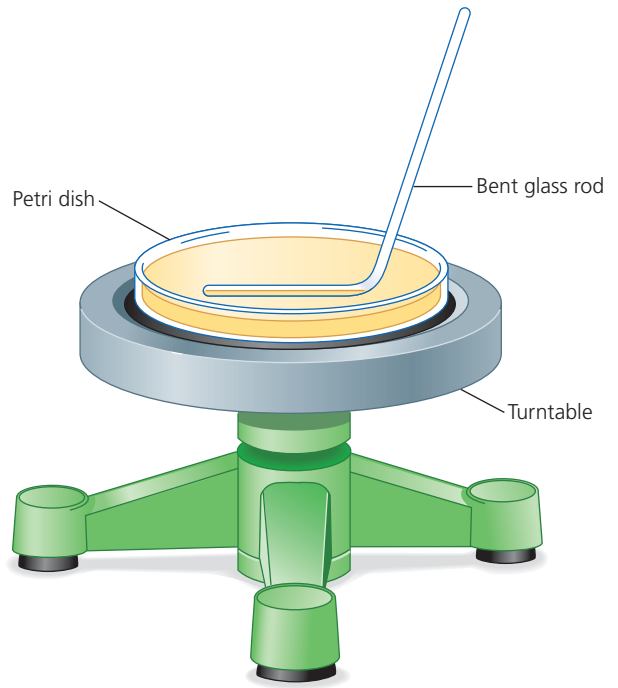

**Figure 3.3** Petri dish turntable

- c. Reflame and cool the loop and again turn the dish  $90^\circ$ . Streak Area 3 in the same manner as Area 2.
  - d. Without reflaming the loop, again turn the dish  $90^\circ$  and then drag the culture from a corner of Area 3 across Area 4, using a wider streak. Don't let the loop touch any of the previously streaked areas. The flaming of the loop at the points indicated is to dilute the culture so that fewer organisms are streaked in each area, resulting in the final desired separation. A photograph of a streak-plate inoculation is shown in **Figure 3.2**.
2. The **spread-plate** technique requires that a previously diluted mixture of microorganisms be used. During inoculation, the cells are spread over the surface of a solid agar medium with a sterile, L-shaped bent glass rod while the Petri dish is spun on a "lazy Susan" turntable (**Figure 3.3**). The step-by-step procedure for this technique is as follows:
    - a. Place the bent glass rod into a beaker and add a sufficient amount of 95% ethyl alcohol to cover the lower, bent portion.
    - b. Place an appropriately labeled nutrient agar plate on the turntable. With a sterile pipette, place one drop of sterile water on the center of the plate, followed by a sterile loopful of *Micrococcus luteus*. Mix gently with the loop and replace the cover.
    - c. Remove the glass rod from the beaker, and pass it through the Bunsen burner flame
      - d. Remove the Petri dish cover and spin the turntable.
      - e. While the turntable is spinning, lightly touch the sterile bent rod to the surface of the agar and move it back and forth. This will spread the culture over the agar surface.
      - f. When the turntable comes to a stop, replace the cover. Immerse the rod in alcohol and reflame.
      - g. In the absence of a turntable, turn the Petri dish manually and spread the culture with the sterile bent glass rod.
  3. The **pour-plate** technique requires a serial dilution of the mixed culture by means of a loop or pipette. The diluted inoculum is then added to a molten agar medium in a Petri dish, mixed, and allowed to solidify. The serial dilution and pour-plate procedures are outlined in Experiment 20.

## CLINICAL APPLICATION

### Isolation of Cultures as a Diagnostic Technique

The isolation of pure cultures is the most important diagnostic tool used in a clinical or research laboratory to uncover the cause of an infection or disease. Before any biochemical or molecular techniques may be used to identify or characterize the causative organism, an individual bacterial colony must be isolated for testing. The isolation of *Staphylococcus aureus* from cultures taken from abscesses or *Streptococcus pyogenes* from a throat culture are two examples of clinical applications of this technique.

## AT THE BENCH



## Materials

### Cultures

24- to 48-hour nutrient broth cultures of a mixture of one part *Serratia marcescens* and three parts *Micrococcus luteus* and a mixture of one part *Escherichia coli* and ten parts *Micrococcus luteus*. For the spread-plate procedure, adjust the cultures to an absorbance (A) of 0.1 at 600 nanometers (nm).

Sources of mixed cultures from the environment could include cultures from a table top, bathroom sink, water fountain, or inside of an incubator. Each student should obtain a mixed culture from one of the environmental sources listed above.

### Media

Three Trypticase™ soy agar plates per designated student group for each inoculation technique to be performed.

### Equipment

Bunsen burner, inoculating loop, turntable, 95% ethyl alcohol, 500-ml beaker, L-shaped bent glass rod, glassware marking pencil, culture tubes containing 1 ml of sterile water, test tube rack, and sterile cotton swabs.

## Procedure Lab One

1. Following the procedures previously described, prepare a spread-plate and/or streak-plate inoculation of each test culture on an appropriately labeled plate.
2. Prepare an environmental mixed culture.
  - a. Dampen a sterile cotton swab with sterile water. Wring out the excess water by pressing the wet swab against the walls of the tube.
  - b. With the moistened cotton swab, obtain your mixed-culture specimen from one of the selected environmental sources listed in the section on cultures.
  - c. Place the contaminated swab back into the tube of sterile water. Mix gently and let stand for 5 minutes.
  - d. Perform spread-plate and/or streak-plate inoculation on an appropriately labeled plate.
3. Incubate all plates in an inverted position for 48 to 72 hours at 25°C.

## Procedure Lab Two

1. Examine all agar plate cultures to identify the distribution of colonies. In the charts provided in Part A of the Lab Report, complete the following:
  - a. Draw the distribution of colonies appearing on each of the agar plate cultures.
  - b. On each of the agar plate cultures, select two discrete colonies that differ in appearance. Using Figure 4.1 on page 30 as a reference, describe each colony as to its  
Form: Circular, irregular, or spreading.  
Elevation: Flat, slightly raised, or markedly raised.  
Pigmentation.  
Size: Pinpoint, small, medium, or large.
2. Retain the mixed-culture plates to perform Part B of this experiment.

## PART B Isolation of Pure Cultures from a Spread-Plate or Streak-Plate Preparation

### LEARNING OBJECTIVE

Once you have completed this experiment, you should be able to

1. Prepare a stock culture of an organism using isolates from mixed cultures prepared on an agar streak plate and/or spread plate.

### Principle

Once discrete, well-separated colonies develop on the surface of a nutrient agar plate culture, each may be picked up with a sterile needle and transferred to separate nutrient agar slants. Each of these new slant cultures represents the growth of a single bacterial species and is designated as a **pure** or **stock culture**.

### CLINICAL APPLICATION

#### Transferring a Colony of Bacteria Daughter Cells

For identification of a bacterial pathogen, a discrete bacterial colony must be transferred from a streak or spread plate to the new testing media. This new culture will consist of daughter cells that are genetic and metabolic clones of the original bacterial cells that were transferred to the plate. This will allow for identification of the unknown bacterial species through its biochemical and molecular characteristics.

## AT THE BENCH



### Materials

#### Cultures

Mixed-culture, nutrient agar streak-plate and/or spread-plate preparations of *S. marcescens* and *M. luteus*, *M. luteus* and *E. coli*, and the environmental specimen plate from Part A.

#### Media

Four Trypticase soy agar slants per designated student group.

#### Equipment

Bunsen burner, inoculating needle, and glassware marking pencil.

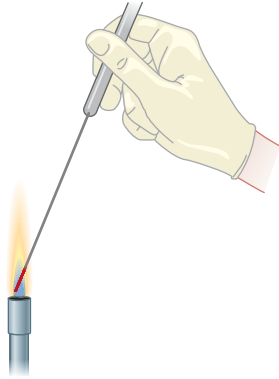
### Procedure Lab One

1. Aseptically transfer, from visibly discrete colonies, the yellow *M. luteus*, the white *E. coli*, the red *S. marcescens*, and a discrete colony from the environmental agar plate specimen to the appropriately labeled agar slants as shown in **Figure 3.4**.
2. Incubate all slants at 37°C for 18 to 24 hours.


### Procedure Lab Two

1. In the chart provided in Part B of the Lab Report, complete the following:
  - a. Draw and indicate the type of growth of each pure-culture isolate, using Figure 4.1 on page 30 as a reference.
  - b. Observe the color of the growth and record its pigmentation.
  - c. Indicate the name of the isolated organisms.


## PROCEDURE



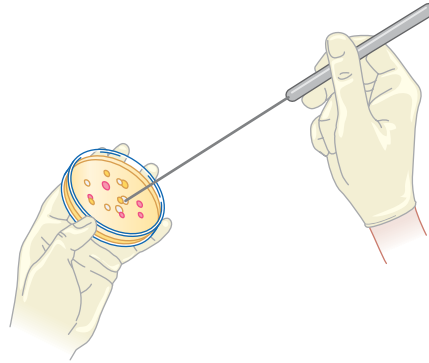
- 1 Flame the straight needle until the entire wire is red.



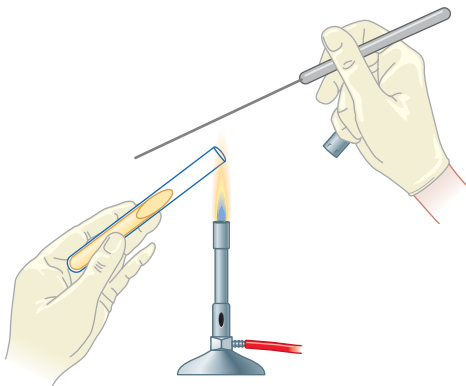
- 2 After isolating a discrete colony on the agar streak plate, touch the straight needle to the surface of the selected colony.



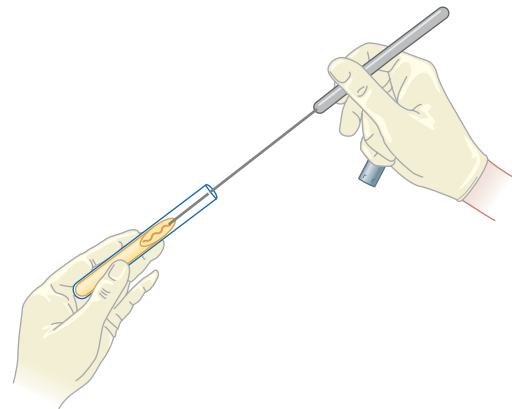
- 3 Uncap the agar slant and pass the neck of the tube rapidly over the Bunsen burner flame.



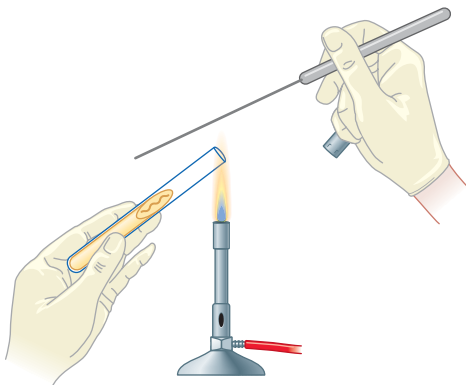- 4 Inoculate the slant by drawing the needle upward in a zigzag motion along the surface of the agar. Do not dig into the agar.



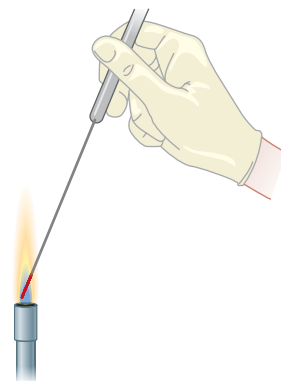- 5 Flame the neck of the tube and recap.



- 6 Flame the inoculating needle.

**Figure 3.4** Procedure for the preparation of a pure culture



*This page intentionally left blank*



Name: \_\_\_\_\_



Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results
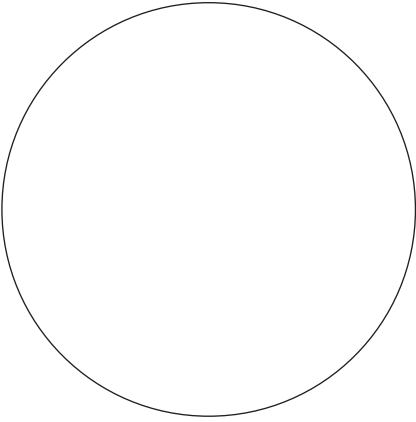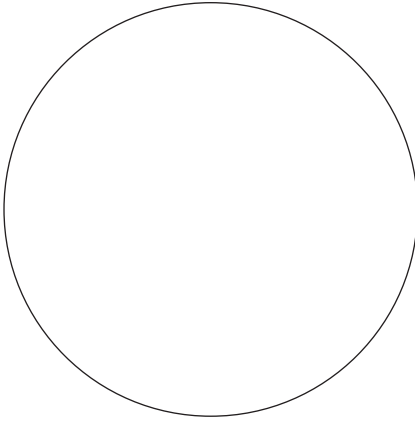
### PART A: Isolation of Discrete Colonies from a Mixed Culture

	SPREAD-PLATE TECHNIQUE			
	<i>S. marcescens</i> and <i>M. luteus</i>		<i>M. luteus</i> and <i>E. coli</i>	
Draw the colonies that appear on each agar plate.				
Colony description:	Isolate 1	Isolate 2	Isolate 1	Isolate 2
Form	_____	_____	_____	_____
Elevation	_____	_____	_____	_____
Pigmentation	_____	_____	_____	_____
Size	_____	_____	_____	_____

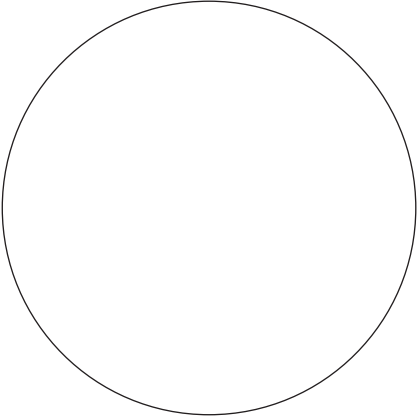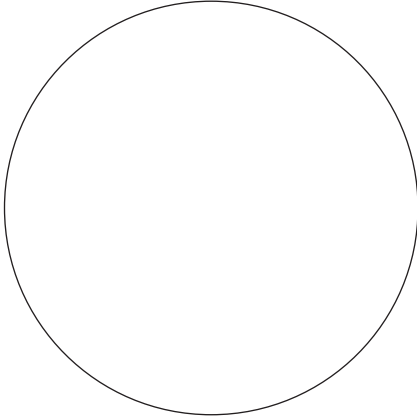	STREAK-PLATE TECHNIQUE			
	<i>S. marcescens</i> and <i>M. luteus</i>		<i>M. luteus</i> and <i>E. coli</i>	
Draw the colonies that appear on each agar plate.				
Colony description:	Isolate 1	Isolate 2	Isolate 1	Isolate 2
Form	_____	_____	_____	_____
Elevation	_____	_____	_____	_____
Pigmentation	_____	_____	_____	_____
Size	_____	_____	_____	_____

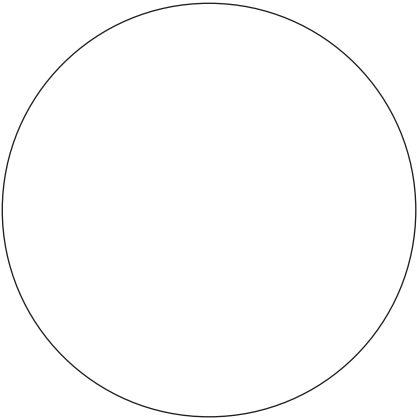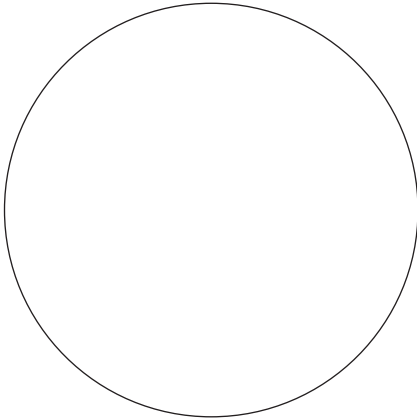	ENVIRONMENTAL SPECIMEN	
	Spread-Plate Technique	Streak-Plate Technique
Draw the colonies that appear on each agar plate.		
Colony description:		
Form	_____	_____
Elevation	_____	_____
Pigmentation	_____	_____
Size	_____	_____







### LEARNING OBJECTIVE

Once you have completed this experiment, you should be able to

1. Determine the cultural characteristics of microorganisms as an aid in identifying and classifying organisms into taxonomic groups.

## Principle

When grown on a variety of media, microorganisms will exhibit differences in the macroscopic appearance of their growth. These differences, called **cultural characteristics**, are used as a basis for separating microorganisms into taxonomic groups. The cultural characteristics for all known microorganisms are contained in *Bergey's Manual of Systematic Bacteriology*. They are determined by culturing the organisms on nutrient agar slants and plates, in nutrient broth, and in nutrient gelatin. The patterns of growth to be considered in each of these media are described below, and some are illustrated in **Figure 4.1**.

### Nutrient Agar Slants


These have a single straight line of inoculation on the surface and are evaluated in the following manner:

1. **Abundance of growth:** The amount of growth is designated as none, slight, moderate, or large.
2. **Pigmentation:** Chromogenic microorganisms may produce intracellular pigments that are responsible for the coloration of the organisms as seen in surface colonies. Other organisms produce extracellular soluble pigments that are excreted into the medium and that also produce a color. Most organisms, however, are nonchromogenic and will appear white to gray.
3. **Optical characteristics:** Optical characteristics may be evaluated on the basis of the amount of light transmitted through the growth. These characteristics are described as **opaque** (no light transmission), **translucent** (partial transmission), or **transparent** (full transmission).
4. **Form:** The appearance of the single-line streak of growth on the agar surface is designated as
  - a. **Filiform:** Continuous, threadlike growth with smooth edges.
  - b. **Echinulate:** Continuous, threadlike growth with irregular edges.
  - c. **Beaded:** Nonconfluent to semiconfluent colonies.
  - d. **Effuse:** Thin, spreading growth.
  - e. **Arborescent:** Treelike growth.
  - f. **Rhizoid:** Rootlike growth.
5. **Consistency:**
  - a. **Dry:** Free from moisture.
  - b. **Buttery:** Moist and shiny.
  - c. **Mucoid:** Slimy and glistening.

### Nutrient Agar Plates

These demonstrate well-isolated colonies and are evaluated in the following manner:

1. **Size:** Pinpoint, small, moderate, or large.
2. **Pigmentation:** Color of colony.
3. **Form:** The shape of the colony is described as follows:
  - a. **Circular:** Unbroken, peripheral edge.
  - b. **Irregular:** Indented, peripheral edge.
  - c. **Rhizoid:** Rootlike, spreading growth.
4. **Margin:** The appearance of the outer edge of the colony is described as follows:
  - a. **Entire:** Sharply defined, even.
  - b. **Lobate:** Marked indentations.
  - c. **Undulate:** Wavy indentations.

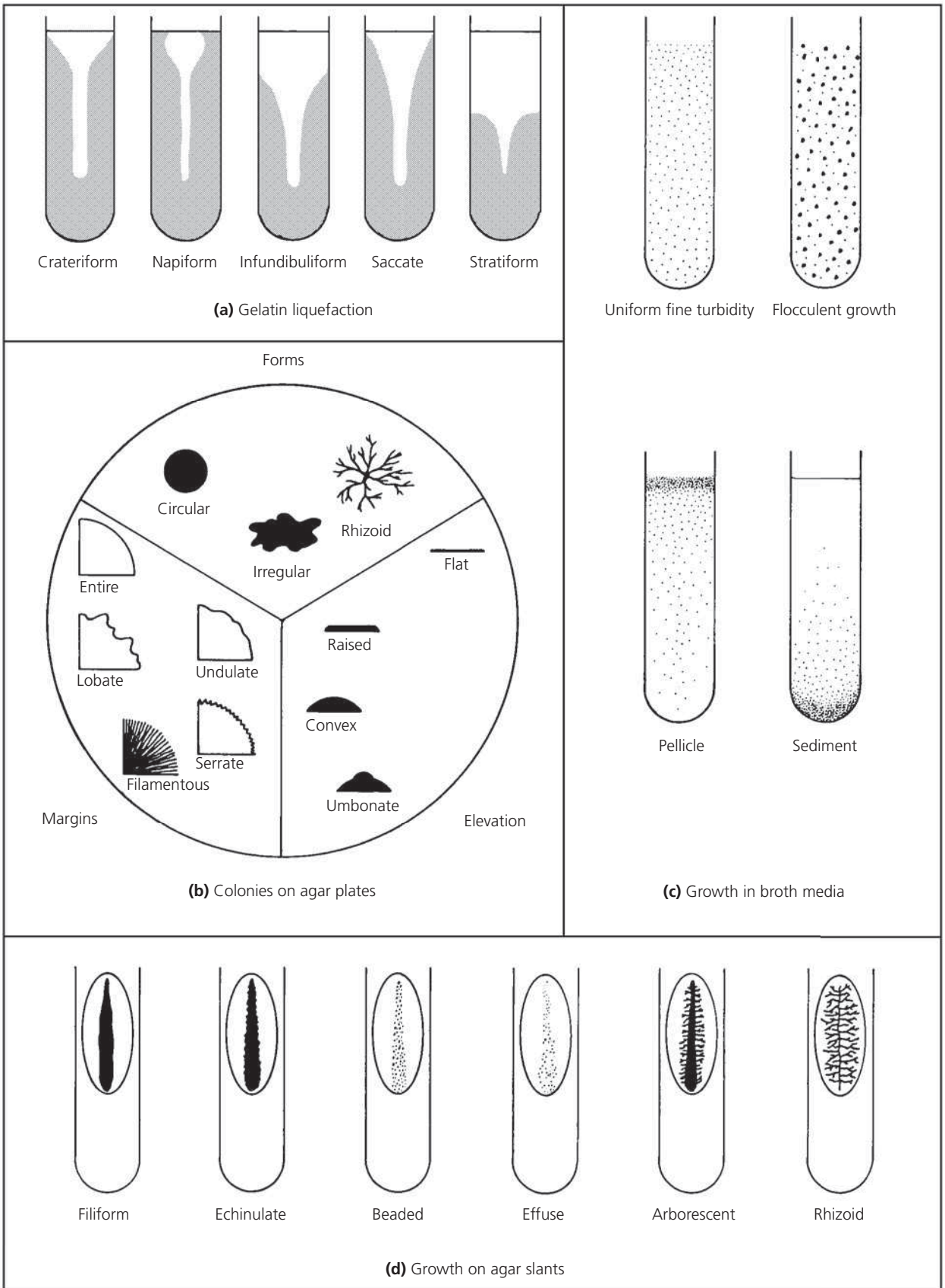

**Figure 4.1** Cultural characteristics of bacteria

- d. **Serrate:** Toothlike appearance.
  - e. **Filamentous:** Threadlike, spreading edge.
5. **Elevation:** The degree to which colony growth is raised on the agar surface is described as follows:
- a. **Flat:** Elevation not discernible.
  - b. **Raised:** Slightly elevated.
  - c. **Convex:** Dome-shaped elevation.
  - d. **Umbonate:** Raised, with elevated convex central region.

## Nutrient Broth Cultures

These are evaluated as to the distribution and appearance of the growth as follows:

1. **Uniform fine turbidity:** Finely dispersed growth throughout.
2. **Flocculent:** Flaky aggregates dispersed throughout.
3. **Pellicle:** Thick, padlike growth on surface.
4. **Sediment:** Concentration of growth at the bottom of broth culture may be granular, flaky, or flocculant.

## Nutrient Gelatin

This solid medium may be liquefied by the enzymatic action of gelatinase. Liquefaction occurs in a variety of patterns:

1. **Crateriform:** Liquefied surface area is saucer-shaped.
2. **Napiform:** Bulbous-shaped liquefaction at surface.
3. **Infundibuliform:** Funnel-shaped.
4. **Saccate:** Elongated, tubular.
5. **Stratiform:** Complete liquefaction of the upper half of the medium.

## CLINICAL APPLICATION

### Examining Colony Growth Characteristics to Aid Identification

Bacterial species each have a characteristic pattern of colony growth in a liquid culture or on a solid medium. While not truly a diagnostic tool, recognition of these patterns of characteristics will aid in a clinical lab setting by helping to minimize the list of potential bacterial species to test for.

## AT THE BENCH



## Materials

### Cultures

24-hour nutrient broth cultures of *Pseudomonas aeruginosa*, *Bacillus cereus*, *Micrococcus luteus*, and *Escherichia coli*. 72- to 96-hour Trypticase soy broth culture of *Mycobacterium smegmatis*.

### Media

Per designated student group: five each of nutrient agar slants, nutrient agar plates, nutrient broth tubes, and nutrient gelatin tubes.

### Equipment

Bunsen burner, inoculating loop and needle, and glassware marking pencil.

## Procedure Lab One

1. Using aseptic technique, inoculate each of the appropriately labeled media listed below in the following manner:
  - a. Nutrient agar slants: With a sterile needle, make a single-line streak of each of the cultures provided, starting at the butt and drawing the needle up the center of the slanted agar surface.
  - b. Nutrient agar plates: With a sterile loop, prepare a streak-plate inoculation of each of the cultures for the isolation of discrete colonies.
  - c. Nutrient broth cultures: Using a sterile loop, inoculate each organism into a tube of nutrient broth. Shake the loop a few times to dislodge the inoculum.
  - d. Nutrient gelatin: Using a sterile needle, prepare a stab inoculation of each of the cultures provided.
2. Incubate all cultures at 37°C for 24 to 48 hours.

## Procedure Lab Two

1. Before beginning observation of all the cultures, place the gelatin cultures in a refrigerator for 30 minutes or in a beaker of crushed ice for a few minutes. The gelatin culture will be the last to be observed.

2. Refer to Figure 4.1 on page 30 and the descriptions presented in the introductory section of Experiment 4 while making the following observations:
  - a. Nutrient agar slants: Observe each of the nutrient agar slant cultures for the amount, pigmentation, form, and consistency of the growth. Record your observations in the chart provided in the Lab Report.
  - b. Nutrient agar plates: Observe a single, well-isolated colony on each of the nutrient agar plate cultures and identify its size, elevation, margin, form, and pigmentation.

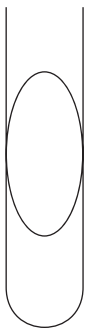
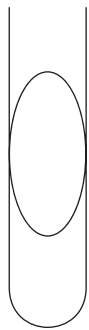



Record your observations in the chart provided in the Lab Report.
  - c. Nutrient broth cultures: Observe each of the nutrient broth cultures for the appearance of growth (flocculation, turbidity, sediment, or pellicle). Record your observations in the chart provided in the Lab Report.
  - d. Nutrient gelatin: Remove gelatin cultures from the refrigerator or beaker of crushed ice, and observe whether liquefaction of the medium has developed and whether the organism has produced gelatinase. Record your observations in the chart provided in the Lab Report.

Name: \_\_\_\_\_

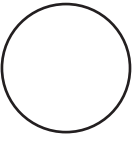
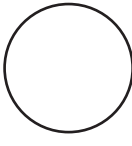
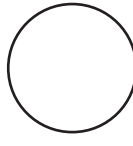
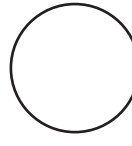
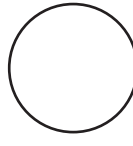
Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

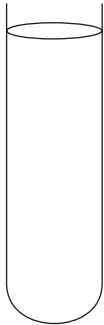



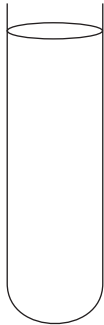
### Nutrient Agar Slants

	NUTRIENT AGAR SLANT CULTURES				
	<i>M. luteus</i>	<i>P. aeruginosa</i>	<i>M. smegmatis</i>	<i>E. coli</i>	<i>B. cereus</i>
Draw the distribution of growth on the slant surface.					
Amount of growth	_____	_____	_____	_____	_____
Pigmentation	_____	_____	_____	_____	_____
Form	_____	_____	_____	_____	_____
Consistency	_____	_____	_____	_____	_____






### Nutrient Agar Plates

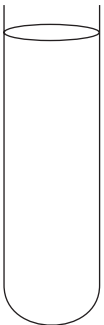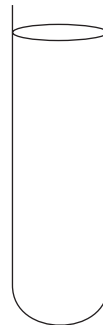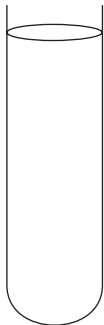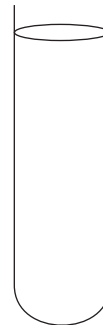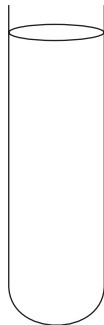	NUTRIENT AGAR PLATES				
	<i>M. luteus</i>	<i>P. aeruginosa</i>	<i>M. smegmatis</i>	<i>E. coli</i>	<i>B. cereus</i>
Draw distribution of colonies.					
Size	_____	_____	_____	_____	_____
Elevation	_____	_____	_____	_____	_____
Margin	_____	_____	_____	_____	_____
Form	_____	_____	_____	_____	_____
Pigmentation	_____	_____	_____	_____	_____

## Nutrient Broth Cultures

	NUTRIENT BROTH CULTURES				
	<i>M. luteus</i>	<i>P. aeruginosa</i>	<i>M. smegmatis</i>	<i>E. coli</i>	<i>B. cereus</i>
Draw the distribution of growth.					
Appearance of growth	_____	_____	_____	_____	_____

## Nutrient Gelatin

	NUTRIENT GELATIN CULTURES				
	<i>M. luteus</i>	<i>P. aeruginosa</i>	<i>M. smegmatis</i>	<i>E. coli</i>	<i>B. cereus</i>
Draw liquefaction patterns.					
Liquefaction (+) or (-)	_____	_____	_____	_____	_____
Type of liquefaction	_____	_____	_____	_____	_____

# Microscopy

## LEARNING OBJECTIVES

Once you have completed the experiments in this section, you should be

1. Familiar with the history and diversity of microscopic instruments.
2. Able to understand the components, use, and care of the brightfield microscope.
3. Able to correctly use the microscope for observation and measurement of microorganisms.

## Introduction

Microbiology, the branch of science that has so vastly extended and expanded our knowledge of the living world, owes its existence to Antoni van Leeuwenhoek. In 1673, with the aid of a crude microscope consisting of a biconcave lens enclosed in two metal plates, Leeuwenhoek introduced the world to the existence of microbial forms of life. Over the years, microscopes have evolved from the simple, single-lens instrument of Leeuwenhoek, with a magnification of  $300\times$ , to the present-day electron microscopes capable of magnifications greater than  $250,000\times$ .

Microscopes are designated as either light microscopes or electron microscopes. The former use visible light or ultraviolet rays to illuminate specimens. They include brightfield, darkfield, phase-contrast, and fluorescent instruments. Fluorescent microscopes use ultraviolet radiations whose wavelengths are shorter than those of visible light and are not directly perceptible to the human eye. Electron microscopes use electron beams (instead of light rays) and magnets (instead of lenses) to observe submicroscopic particles.

## Essential Features of Various Microscopes

**Brightfield Microscope** This instrument contains two-lens systems for magnifying specimens: the ocular lens in the eyepiece and the objective lens located in the nosepiece. The specimen is illuminated by a beam of tungsten light focused on it by a substage lens called a condenser; the result is a specimen that appears dark against a bright background. A major limitation of this system is the absence of contrast between the specimen and the surrounding medium, which makes it difficult to observe living cells. Therefore, most brightfield observations are performed on nonviable, stained preparations.

**Darkfield Microscope** This is similar to the ordinary light microscope; however, the condenser system is modified so that the specimen is not illuminated directly. The condenser directs the light obliquely so that the light is deflected or scattered from the specimen, which then appears bright against a dark background. Living specimens may be observed more readily with darkfield than with brightfield microscopy.



**Phase-Contrast Microscope** Observation of microorganisms in an unstained state is possible with this microscope. Its optics include special objectives and a condenser that make visible cellular components that differ only slightly in their refractive indexes. As light is transmitted through a specimen with a refractive index different from that of the surrounding medium, a portion of the light is refracted (bent) due to slight variations in density and thickness of the cellular components. The special optics convert the difference between transmitted light and refracted rays, resulting in a significant variation in the intensity of light and thereby producing a discernible image of the structure under study. The image appears dark against a light background.

**Fluorescent Microscope** This microscope is used most frequently to visualize specimens that are chemically tagged with a fluorescent dye. The source of illumination is an ultraviolet (UV) light obtained from a high-pressure mercury lamp or hydrogen quartz lamp. The ocular lens is fitted with a filter that permits the longer ultraviolet wavelengths to pass, while the shorter wavelengths are blocked or eliminated. Ultraviolet radiations are absorbed by the fluorescent label, and the energy is re-emitted in the form of a different wavelength in the visible light range. The fluorescent dyes absorb at wavelengths between 230 and 350 nanometers (nm) and emit orange, yellow, or greenish light. This microscope is used primarily for the detection of antigen-antibody reactions. Antibodies are conjugated with a fluorescent dye that becomes excited in the presence of ultraviolet

light, and the fluorescent portion of the dye becomes visible against a black background.

**Electron Microscope** This instrument provides a revolutionary method of microscopy, with magnifications up to 1 million $\times$ . This permits visualization of submicroscopic cellular particles as well as viral agents. In the electron microscope, the specimen is illuminated by a beam of electrons rather than light, and the focusing is carried out by electromagnets instead of a set of optics. These components are sealed in a tube in which a complete vacuum is established. Transmission electron microscopes require specimens that are prepared as thin filaments, fixed and dehydrated for the electron beam to pass freely through them. As the electrons pass through the specimen, images are formed by directing the electrons onto photographic film, thus making internal cellular structures visible. Scanning electron microscopes are used for visualizing surface characteristics rather than intracellular structures. A narrow beam of electrons scans back and forth, producing a three-dimensional image as the electrons are reflected off the specimen's surface.

While scientists have a variety of optical instruments with which to perform routine laboratory procedures and sophisticated research, the compound brightfield microscope is the “workhorse” and is commonly found in all biological laboratories. Although you should be familiar with the basic principles of microscopy, you probably have not been exposed to this diverse array of complex and expensive equipment. Therefore, only the compound brightfield microscope will be discussed in depth and used to examine specimens.

### LEARNING OBJECTIVES

*Once you have completed this experiment, you should be familiar with the*

1. Theoretical principles of brightfield microscopy.
2. Component parts of the compound microscope.
3. Use and care of the compound microscope.
4. Practical use of the compound microscope for visualization of cellular morphology from stained slide preparations.

## Principle

Microbiology is a science that studies living organisms that are too small to be seen with the naked eye. Needless to say, such a study must involve the use of a good compound microscope. Although there are many types and variations, they all fundamentally consist of a two-lens system, a variable but controllable light source, and mechanical adjustable parts for determining focal length between the lenses and specimen (**Figure 5.1**).

## Components of the Microscope

**Stage** A fixed platform with an opening in the center allows the passage of light from an illuminating source below to the lens system above the stage. This platform provides a surface for the placement of a slide with its specimen over the central opening. In addition to the fixed stage, most microscopes have a **mechanical stage** that can be moved vertically or horizontally by means of adjustment controls. Less sophisticated microscopes have clips on the fixed stage, and the slide must be positioned manually over the central opening.

**Illumination** The light source is positioned in the base of the instrument. Some microscopes are equipped with a built-in light source to provide direct illumination. Others are provided with a reversible mirror that has one side flat and the other concave. An external light source, such as a lamp, is placed in front of the mirror to direct the light upward into the lens system. The flat side of the mirror is used for artificial light, and the concave side for sunlight.


**Abbé Condenser** This component is found directly under the stage and contains two sets of lenses that collect and concentrate light as it passes upward from the light source into the lens systems. The condenser is equipped with an **iris diaphragm**, a shutter controlled by a lever that is used to regulate the amount of light entering the lens system.

**Body Tube** Above the stage and attached to the arm of the microscope is the body tube. This structure houses the lens system that magnifies the specimen. The upper end of the tube contains the **ocular** or **eyepiece** lens. The lower portion consists of a movable **nosepiece** containing the **objective lenses**. Rotation of the nosepiece positions objectives above the stage opening. The body tube may be raised or lowered with the aid of **coarse-adjustment** and **fine-adjustment knobs** that are located above or below the stage, depending on the type and make of the instrument.

## Theoretical Principles of Microscopy

To use the microscope efficiently and with minimal frustration, you should understand the basic principles of microscopy: magnification, resolution, numerical aperture, illumination, and focusing.

**Magnification** Enlargement, or magnification, of a specimen is the function of a two-lens system; the **ocular lens** is found in the eyepiece, and the **objective lens** is situated in a

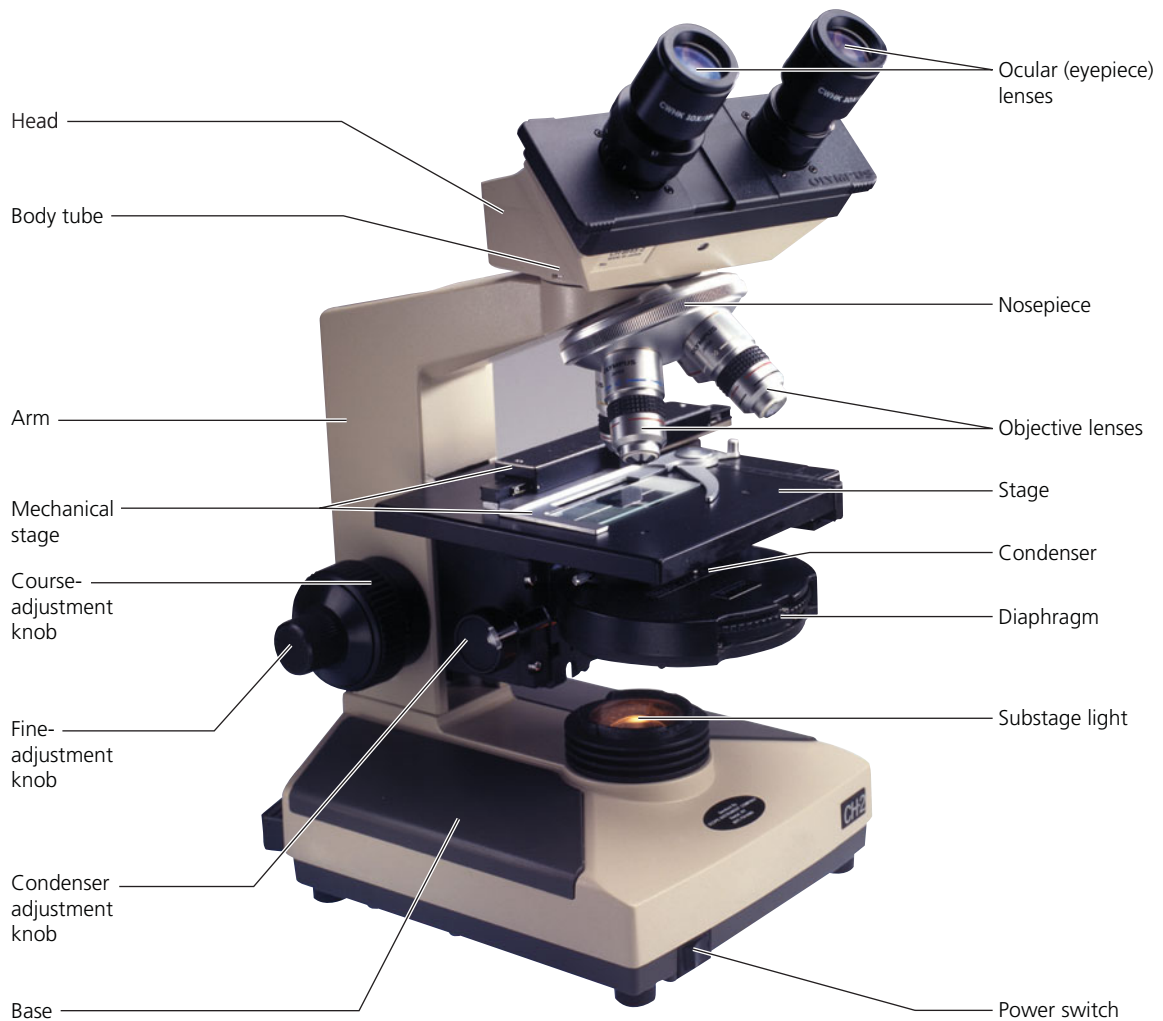

**Figure 5.1** A compound microscope

revolving nosepiece. These lenses are separated by the **body tube**. The objective lens is nearer the specimen and magnifies it, producing the **real image** that is projected up into the focal plane and then magnified by the ocular lens to produce the final image.

The most commonly used microscopes are equipped with a revolving nosepiece containing four objective lenses, each possessing a different degree of magnification. When these are combined with the magnification of the ocular lens, the total or overall linear magnification of the specimen is obtained. This is shown in Table 5.1.

**Resolving Power or Resolution** Although magnification is important, you must be aware that unlimited enlargement is not possible by merely increasing the magnifying power of the lenses or by using additional lenses, because lenses are limited by a property called **resolving power**. By definition, resolving power is how far

apart two adjacent objects must be before a given lens shows them as discrete entities. When a lens cannot discriminate, that is, when the two objects appear as one, it has lost resolution. Increased magnification will not rectify the loss and will, in fact, blur the object. The resolving power of a lens is dependent on the wavelength of light used and the **numerical aperture**, which is a characteristic of each lens and imprinted on each objective. The numerical aperture is defined as a function of the diameter of the objective lens in relation to its focal length. It is doubled by use of the substage condenser, which illuminates the object with rays of light that pass through the specimen obliquely as well as directly. Thus, resolving power is expressed mathematically as follows:

$$\text{resolving power} = \frac{\text{wavelength of light}}{2 \times \text{numerical aperture}}$$

**TABLE 5.1 Overall Linear Magnification**

MAGNIFICATION		TOTAL MAGNIFICATION
OBJECTIVE LENSES	OCULAR LENS	OBJECTIVE MULTIPLIED BY OCULAR
Scanning 4×	10×	40×
Low-power 10×	10×	100×
High-power 40×	10×	400×
Oil-immersion 100×	10×	1000×

Based on this formula, the shorter the wavelength, the greater the resolving power of the lens. Thus, for the same numerical aperture, short wavelengths of the electromagnetic spectrum are better suited for higher resolution than are longer wavelengths.


However, as with magnification, resolving power also has limits. You might rationalize that merely decreasing the wavelength will automatically increase the resolving power of a lens. Such is not the case, because the visible portion of the electromagnetic spectrum is very narrow and borders on the very short wavelengths found in the ultraviolet portion of the spectrum.

The relationship between wavelength and numerical aperture is valid only for increased resolving power when light rays are parallel. Therefore, the resolving power is also dependent on another factor, the **refractive index**. This is the bending power of light passing through air from the glass slide to the objective lens. The refractive index of air is lower than that of glass; as light rays pass from the glass slide into the air, they are bent or refracted so that they do not pass into the objective lens. This would cause a loss of light, which would reduce the numerical aperture and diminish the resolving power of the objective lens. Loss of refracted light can be compensated for by interposing mineral oil, which has the same refractive index as glass, between the slide and the objective lens. In this way, decreased light refraction occurs and more light rays enter directly into the objective lens, producing a vivid image with high resolution (**Figure 5.2**).

**Illumination** Effective illumination is required for efficient magnification and resolving power. Since the intensity of daylight is an uncontrolled variable, artificial light from a tungsten lamp is the most commonly used light source in microscopy. The light is passed through the condenser located beneath the stage. The condenser contains

two lenses that are necessary to produce a maximum numerical aperture. The height of the condenser can be adjusted with the **condenser knob**. Always keep the condenser close to the stage, especially when using the oil-immersion objective.

Between the light source and the condenser is the iris diaphragm, which can be opened and closed by means of a lever, thereby regulating the amount of light entering the condenser. Excessive illumination may actually obscure the specimen because of lack of contrast. The amount of light entering the microscope differs with each objective lens used. A rule of thumb is that *as the magnification of the lens increases,*

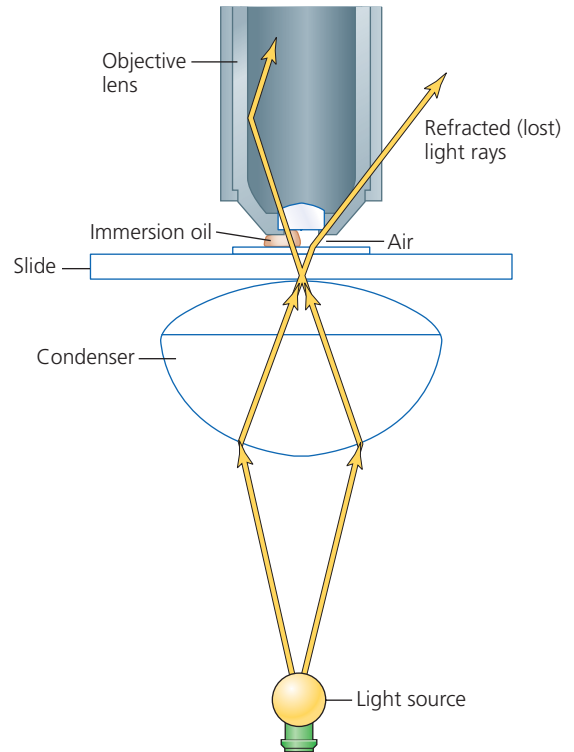**Figure 5.2** Refractive index in air and in mineral oil

the distance between the objective lens and slide, called working distance, decreases, whereas the numerical aperture of the objective lens increases (Figure 5.3).

## Use and Care of the Microscope

You will be responsible for the proper care and use of microscopes. Since microscopes are expensive, you must observe the following regulations and procedures.

The instruments are housed in special cabinets and must be moved by users to their laboratory benches. The correct and only acceptable way to do this is to grip the microscope arm firmly with the right hand and the base with the left hand, and lift the instrument from the cabinet shelf. Carry it close to the body and gently place it on the laboratory bench. This will prevent collision with furniture or coworkers and will protect the instrument against damage.

Once the microscope is placed on the laboratory bench, observe the following rules:

1. Remove all unnecessary materials (such as books, papers, purses, and hats) from the laboratory bench.
2. Uncoil the microscope's electric cord and plug it into an electrical outlet.
3. Clean all lens systems; the smallest bit of dust, oil, lint, or eyelash will decrease the efficiency of the microscope. The ocular, scanning, low-power, and high-power lenses may be cleaned by wiping several times with acceptable lens tissue. Never use paper toweling or cloth on a lens surface. If the oil-immersion lens is gummy or tacky, a piece of lens paper moistened with xylol is used to wipe it clean. The xylol is immediately removed with a tissue moistened with 95% alcohol, and the lens is wiped dry with lens paper. *Note: This xylol cleansing procedure should be performed only by the instructor and only if necessary; consistent use of xylol may loosen the lens.*

The following routine procedures must be followed to ensure correct and efficient use of the microscope.

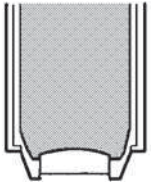
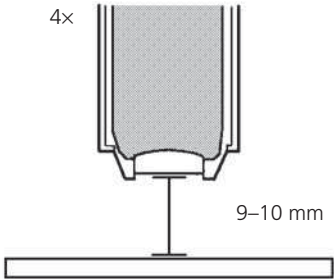
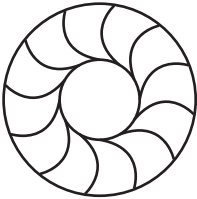
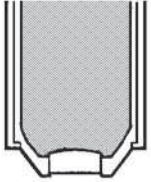
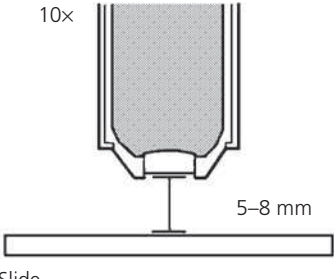
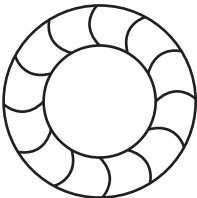
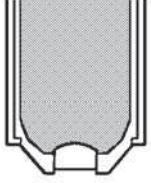
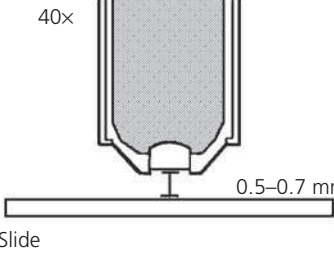
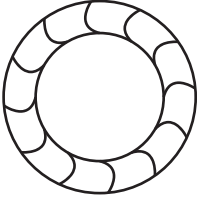
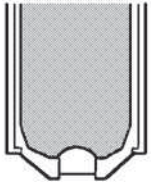
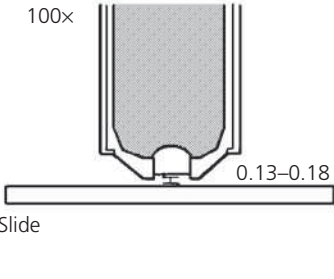
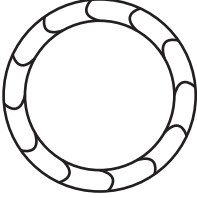
1. Place the microscope slide with the specimen within the stage clips on the fixed stage. Move the slide to center the specimen over the opening in the stage directly over the light source.
2. Raise the microscope stage up as far as it will go. Rotate the scanning lens or low-power lens into position. Lower the body tube with the

coarse-adjustment knob to its lowest position. *Note: Never lower the body tube while looking through the ocular lens.*

3. While looking through the ocular lens, use the fine-adjustment knob, rotating it back and forth slightly, to bring the specimen into sharp focus.
4. Adjust the substage condenser to achieve optimal focus.
5. Routinely adjust the light source by means of the light-source transformer setting, and/or the iris diaphragm, for optimum illumination for each new slide and for each change in magnification.
6. Most microscopes are **parfocal**, which means that when one lens is in focus, other lenses will also have the same focal length and can be rotated into position without further major adjustment. In practice, however, usually a half-turn of the fine-adjustment knob in either direction is necessary for sharp focus.
7. Once you have brought the specimen into sharp focus with a low-powered lens, preparation may be made for visualizing the specimen under oil immersion. Place a drop of oil on the slide directly over the area to be viewed. Rotate the nosepiece until the oil-immersion objective locks into position. *Note: Care should be taken not to allow the high-power objective to touch the drop of oil.* The slide is observed from the side as the objective is rotated slowly into position. This will ensure that the objective will be properly immersed in the oil. The fine-adjustment knob is readjusted to bring the image into sharp focus.
8. During microscopic examination of microbial organisms, it is always necessary to observe several areas of the preparation. This is accomplished by scanning the slide without the application of additional immersion oil. *Note: This will require continuous, very fine adjustments by the slow, back-and-forth rotation of the fine-adjustment knob only.*

On completion of the laboratory exercise, return the microscope to its cabinet in its original condition. The following steps are recommended:

1. Clean all lenses with dry, clean lens paper. *Note: Use xylol to remove oil from the stage only.*
2. Place the low-power objective in position and lower the body tube completely.
3. Center the mechanical stage.

Objective	Working Distance	Diaphragm Opening
Scanning 4x 	4x  9–10 mm Slide	 Reduced
Low power 10x 	10x  5–8 mm Slide	 Not fully opened
High power 40x 	40x  0.5–0.7 mm Slide	 Not fully opened
Oil immersion 100x 	100x  0.13–0.18 mm Slide	 Fully opened

**Figure 5.3** Relationship between working distance, objective, and diaphragm opening

4. Coil the electric cord around the body tube and the stage.
5. Carry the microscope to its position in its cabinet in the manner previously described.

## CLINICAL APPLICATION

### Using Microscopic Examination in the Diagnosis of Tuberculosis

The visualization of stained bacterial cells using a compound light microscope can be the first step in diagnosing microbial infections. For example, a rapid diagnosis for tuberculosis can be made by identifying the unique characteristics of *Mycobacterium tuberculosis* in a stained sample of patient sputum.

## AT THE BENCH



## Materials

### Slides

Commercially prepared slides of *Staphylococcus aureus*, *Bacillus subtilis*, *Aquaspirillum itersonii*, *Saccharomyces cerevisiae*, and a human blood smear.

## Equipment

Compound microscope, lens paper, and immersion oil.

## Procedure

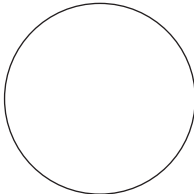
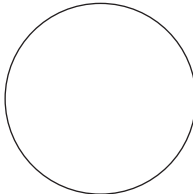
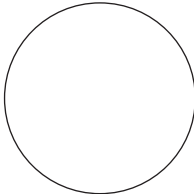
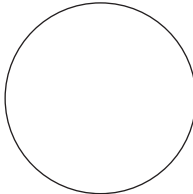
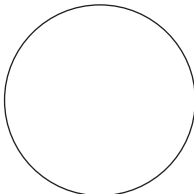
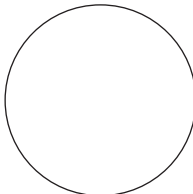
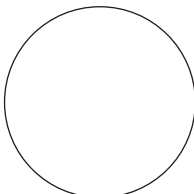
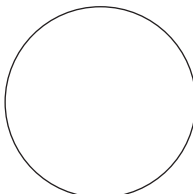
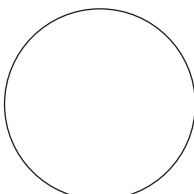
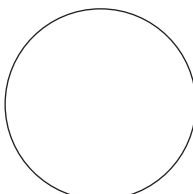
1. Review the parts of the microscope, making sure you know the names and understand the function of each of these components.
2. Review instructions for the use of the microscope, giving special attention to the use of the oil-immersion objective.
3. Examine the prepared slides, noting the shapes and the relative sizes of the cells under the high-power (also called high-dry, because it is the highest power that does not use oil) and oil-immersion objectives.
4. Record your observations in the Lab Report.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_


## Observations and Results

Draw several cells from a typical microscopic field as viewed under each magnification, and give the total magnification for each objective.

	High Power	Oil Immersion
<i>S. aureus</i>		
Magnification	_____	_____
<i>B. subtilis</i>		
Magnification	_____	_____
<i>A. itersonii</i>		
Magnification	_____	_____
<i>S. cerevisiae</i>		
Magnification	_____	_____
Blood smear		
Magnification	_____	_____



## Review Questions

1. Explain why the body tube of the microscope should not be lowered while you are looking through the ocular lens.
2. For what purpose would you adjust each of the following microscope components during a microscopy exercise?
  - a. Iris diaphragm:
  - b. Coarse-adjustment knob:
  - c. Fine-adjustment knob:
  - d. Condenser:
  - e. Mechanical stage control:
3.  As a beginning student in the microbiology laboratory, you experience some difficulties in using the oil-immersion lens. Describe the steps you would take to correct the following problems:
  - a. Inability to bring the specimen into sharp focus.
  - b. Insufficient light while viewing the specimen.
  - c. Artifacts in the microscopic field.

# Microscopic Examination of Living Microorganisms Using a Hanging-Drop Preparation or a Wet Mount

## LEARNING OBJECTIVES

Once you have completed this experiment, you should know how to

1. Microscopically examine living microorganisms.
2. Make a hanging-drop preparation or wet mount to view living microorganisms.

## Principle


Bacteria, because of their small size and a refractive index that closely approximates that of water, do not lend themselves readily to microscopic examination in a living, unstained state. Examination of

living microorganisms is useful, however, to do the following:

1. Observe cell activities such as motility and binary fission.
2. Observe the natural sizes and shapes of the cells, considering that **heat fixation** (the rapid passage of the smear over the Bunsen burner flame) and exposure to chemicals during staining cause some degree of distortion.

In this experiment you will use individual cultures of *Pseudomonas aeruginosa*, *Bacillus cereus*, *Staphylococcus aureus*, and *Proteus vulgaris* for a hanging-drop preparation or a wet mount. Hay infusion or pond water may be substituted or used in addition to the above organisms.

**Figure 6.1** illustrates several organisms commonly found in pond water and hay infusions.

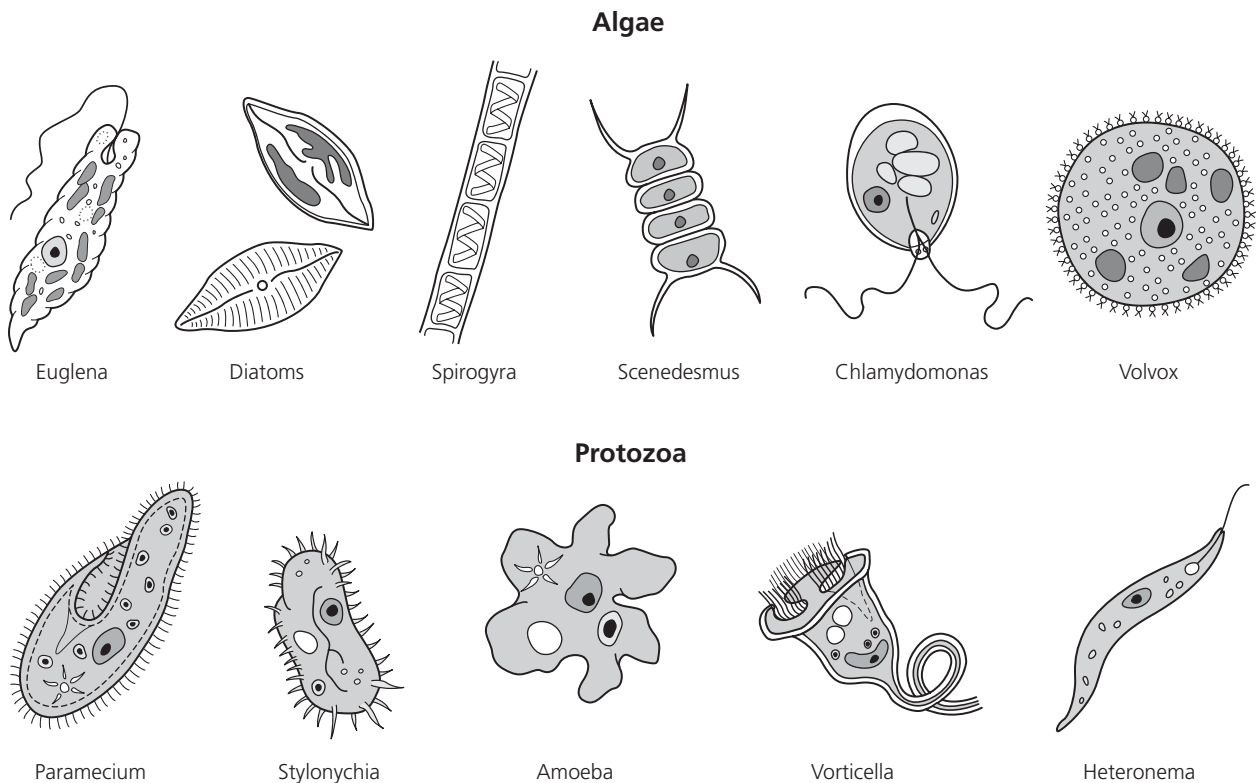

**Figure 6.1** Algae and protozoa commonly found in natural infusions and pond water (drawings are not to scale)

You will observe the preparation(s) microscopically for differences in the sizes and shapes of the cells, as well as for motility, a self-directed movement. It is essential to differentiate between actual motility and **Brownian movement**, a vibratory movement of the cells due to their bombardment by water molecules in the suspension. Hanging-drop preparations and wet mounts make the movement of microorganisms easier to see because they slow down the movement of water molecules.

## CLINICAL APPLICATION

### Observation of Living Bacteria and the Diagnosis of Syphilis

Some microorganisms are difficult or unable to be stained. One such bacteria is *Treponema pallidum*, the causative agent for syphilis. Special stains must be used to stain this bacterium; however, it can be viewed unstained and alive using a darkfield microscope. Under those conditions, its characteristic shape and motility can be observed, leading to a diagnosis of syphilis.

## AT THE BENCH



## Materials

### Cultures

24-hour broth cultures of *P. aeruginosa*, *B. cereus*, *S. aureus*, and *P. vulgaris*; and/or hay infusion broth cultures or pond water. (See Appendix 3 for the preparation of hay infusion broth.)

### Equipment

Bunsen burner, inoculating loop, depression slides, glass slides, coverslips, microscope, petroleum jelly, and cotton swabs.

## Procedure

### Hanging-Drop Preparation

Perform the following steps for each culture provided in this experiment. Steps 1–4 are illustrated in **Figure 6.2**.

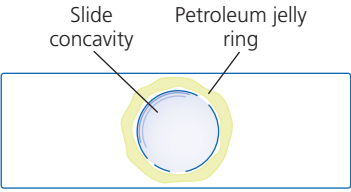
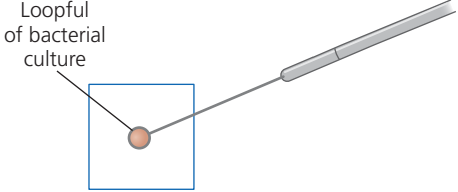
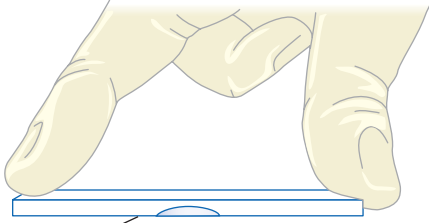
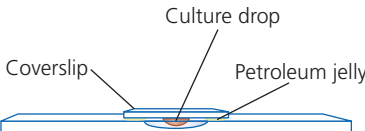
1. With a cotton swab, apply a ring of petroleum jelly around the concavity of the depression slide.
2. Using aseptic technique, place a loopful of the culture in the center of a clean coverslip.
3. Place the depression slide, with the concave surface facing down, over the coverslip so that the depression covers the drop of culture. Press the slide gently to form a seal between the slide and the coverslip.
4. Quickly turn the slide right side up so that the drop continues to adhere to the inner surface of the coverslip.
5. For microscopic examination, first focus on the drop culture under the low-power objective (10 $\times$ ) and reduce the light source by adjusting the Abbé condenser. Repeat using the high-power objective (40 $\times$ ).
6. Examine the hanging-drop preparation and record your observations in the Lab Report.

### Wet Mount

A wet mount may be substituted for the hanging-drop preparation using a similar procedure:

1. With a cotton swab apply a thin layer of petroleum jelly along the edge of the four sides of a coverslip.
2. Using aseptic technique, place a loopful of the culture in the center of a clean coverslip.
3. Place a clean glass slide over the coverslip and press the slide gently to form a seal between the slide and the coverslip.
4. Follow Steps 4 and 5 in the hanging-drop procedure.
5. Examine the wet-mount preparation and record your observations in the Lab Report.

## PROCEDURE

 <p>Slide concavity Petroleum jelly ring</p> <p>Depression slide</p> <p>1 Spread a ring of petroleum jelly around the concavity of the depression slide.</p>	 <p>Loopful of bacterial culture</p> <p>Coverslip</p> <p>2 Place a loopful of the bacterial culture in the center of the coverslip.</p>
 <p>Depression slide Coverslip Culture drop</p> <p>3 Lower the depression slide, with the concavity facing down, onto the coverslip. Press gently to form a seal.</p>	 <p>Culture drop Coverslip Petroleum jelly</p> <p>Hanging-drop preparation</p> <p>4 Turn the hanging-drop preparation over so that the culture drop adheres to the coverslip.</p>

**Figure 6.2** Hanging-drop preparation

*This page intentionally left blank*

Name: \_\_\_\_\_





Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

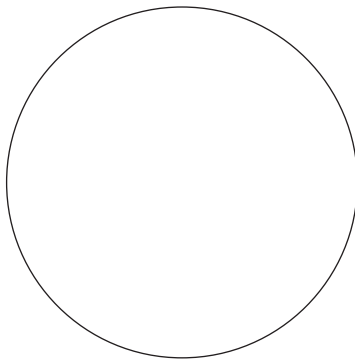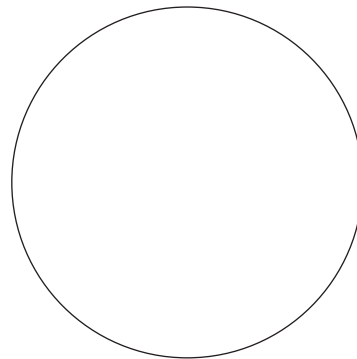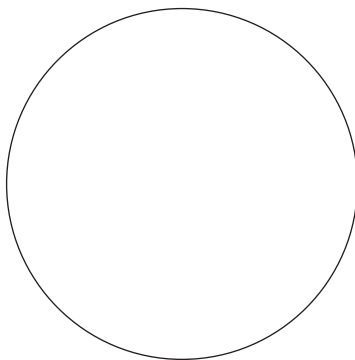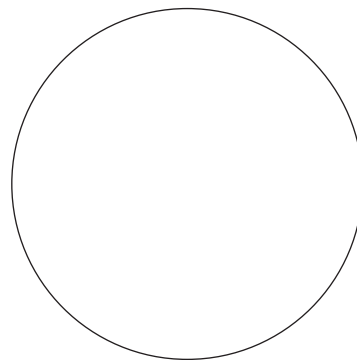
1. Examine the hanging-drop or wet-mount preparation to determine shape and motility of the different bacteria present. Record your results in the chart below.

Organisms	Shape	True Motility or Brownian Movement?
<i>S. aureus</i>		
<i>P. aeruginosa</i>		
<i>B. cereus</i>		
<i>P. vulgaris</i>		


2. Draw a representative field of each of the above organisms.

*S. aureus**P. aeruginosa**B. cereus**P. vulgaris*

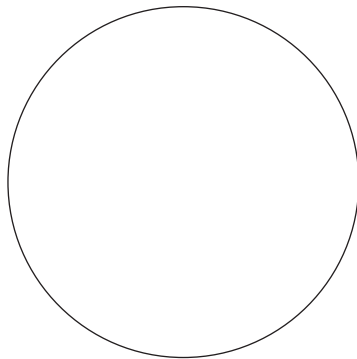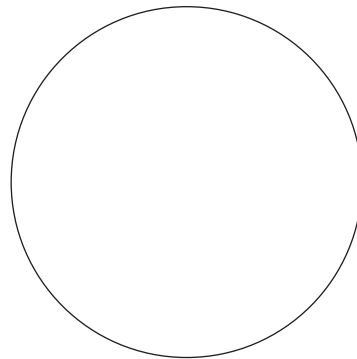3. Draw representative fields of pond water and hay infusion if you used them. Try to identify some of the organisms that you see by referring to Figure 6.1. Note the shape and type of movement in the chart below.





Pond water



Hay infusion

	Pond Water			Hay Infusion		
Shape						
True motility or Brownian movement?						
Organism						

## Review Questions

- Why are living, unstained bacterial preparations more difficult to observe microscopically than stained preparations?
- 
 What is the major advantage of using living cell preparations (hanging-drop or wet mount) rather than stained preparations?
- How do you distinguish between true motility and Brownian movement?
- 
 During the microscopic observation of a drop of stagnant pond water, what criteria would you use to distinguish viable organisms from nonviable suspended debris?

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be able to

1. Calibrate an ocular micrometer.
2. Perform an experimental procedure in the measurement of microorganisms.

## Principle

Determination of microbial size is not as simple as you might assume. Before an accurate measurement of cells can be made, the diameter of the microscopic field must be established by means of optic devices, namely, an ocular micrometer and a stage micrometer. The **ocular micrometer** (Figure 7.1a), which is placed on a circular shelf inside the eyepiece, is a glass disc with graduations etched on its surface. The distance between these graduations will vary depending on the objective being used, which determines the size of the field. This distance is determined by using a **stage micrometer** (Figure 7.1b), a special glass

slide with etched graduations that are 0.01 mm, or 10 micrometers ( $\mu\text{m}$ ), apart.

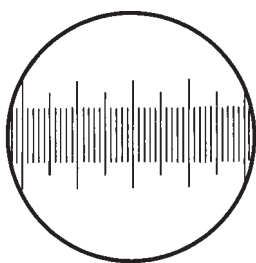
The calibration procedure for the ocular micrometer requires that the graduations on both micrometers be superimposed on each other (Figure 7.1c). This is accomplished by rotating the ocular lens. A determination is then made of the number of ocular divisions per known distance on the stage micrometer. Finally, the calibration factor for one ocular division is calculated as follows:

$$\text{One division on ocular micrometer in mm} = \frac{\text{Known distance between two lines on stage micrometer}}{\text{Number of divisions on ocular micrometer}}$$


**Example:** If 13 ocular divisions coincide with two stage divisions ( $2 \times 0.01 \text{ mm} = \text{known distance of } 0.02 \text{ mm}$ ), then:

$$\begin{aligned} \text{one ocular division} &= \frac{0.02 \text{ mm}}{13} \\ &= 0.00154 \text{ mm, or } 1.54 \mu\text{m} \end{aligned}$$


Once the ocular micrometer is calibrated, the size of a microorganism can easily be determined, first by counting the number of spaces occupied by the organism (Figure 7.1d) and second by



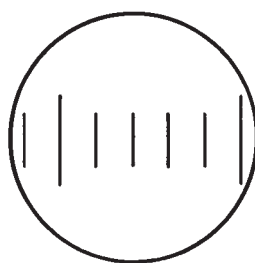
(a) Ocular micrometer: Graduations must be calibrated for each objective on the microscope.



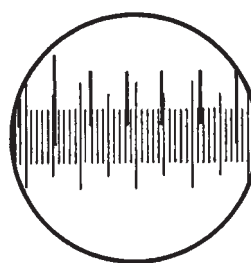(b) Stage micrometer: Graduations are spaced 0.01 mm (10  $\mu\text{m}$ ) apart.



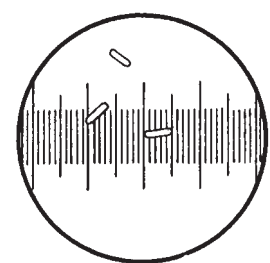(c) Superimposed ocular micrometer over stage micrometer: Thirteen graduations on the ocular micrometer coincide with two graduations on the stage micrometer.



(d) Measurement of microorganisms with the calibrated micrometer: Each cell occupies five divisions on the ocular micrometer.

Figure 7.1 Calibration and use of the ocular micrometer



multiplying this number by the calculated calibration factor for one ocular division.

**Example:** If an organism occupies five spaces on the ocular micrometer, then:

$$\begin{aligned} \text{length of} &= \text{number of} && \text{calibration} \\ \text{organism} &= \text{ocular divisions} && \times \text{factor for one} \\ & \text{occupied} && \text{ocular division} \\ &= 5 \times 1.54 \mu\text{m} = 7.70 \mu\text{m} \end{aligned}$$

In this experiment, you will calibrate an ocular micrometer for the oil-immersion objective and determine the sizes of microorganisms such as bacteria, yeast, and protozoa.

## CLINICAL APPLICATION

### Ocular Micrometers and the Diagnosis of Fungal Infections

Ocular micrometers can measure not only prokaryotic cells, but also eukaryotic cells. Technicians can use ocular micrometers to determine the size of fungal cells found in biopsy tissue in order to differentiate small (1–5  $\mu\text{m}$ ), medium (5–20  $\mu\text{m}$ ), and large (20–200  $\mu\text{m}$ ). The small causal agent of histoplasmosis, *Histoplasma capsulatum*, can easily be distinguished from the large *Coccidioides immitis*, the causal agent of “valley fever.”

## AT THE BENCH



## Materials

### Slides

Prepared slides of yeast cells, protozoa, and bacterial cocci and bacilli.

### Equipment

Ocular micrometer, stage micrometer, microscope, immersion oil, and lens paper.

## Procedure

1. With the assistance of a laboratory instructor, carefully place the ocular micrometer into the eyepiece.
2. Place the stage micrometer on the microscope stage and center it over the illumination source.
3. With the stage micrometer in clear focus under the low-power objective, slowly rotate the eyepiece to superimpose the ocular micrometer graduations over those of the stage micrometer.
4. Add a drop of immersion oil to the stage micrometer, bring the oil-immersion objective into position, and focus, if necessary, with the *fine-adjustment knob only*.
5. Move the mechanical stage so that a line on the stage micrometer coincides with a line on the ocular micrometer at one end. Find another line on the ocular micrometer that coincides with a line on the stage micrometer. Determine the distance on the stage micrometer (number of divisions  $\times$  0.01 mm) and the corresponding number of divisions on the ocular micrometer.
6. Determine the value of the calibration factor for the oil-immersion objective.
7. Remove the stage micrometer from the stage.
8. To determine the size of the cocci on the prepared slides under the oil-immersion objective, do the following:
  - a. Calculate the number of ocular divisions occupied by each of three separate cocci. Record the data in the readings columns of the first observation chart in the Lab Report.
  - b. Determine and record the average of the three measurements.
  - c. Determine the size by multiplying the average by your calculated calibration factor, and record this value.
9. Determine and record the size of the other microorganisms by observing the remaining prepared slides under oil immersion as outlined in Step 8. Since these organisms are not round, both length and width measurements are required.
10. Check that your observations are complete in the Lab Report.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

1. Calibration of ocular micrometer for the oil-immersion objective:

Distance on stage micrometer \_\_\_\_\_

Divisions on ocular micrometer \_\_\_\_\_


Calibration factor for one ocular division \_\_\_\_\_

2. Use the charts below to record your observations and measurements of the microorganisms studied.

	WIDTH OF MICROORGANISMS IN MICROMETERS ( $\mu\text{m}$ )					
	Number of ocular divisions			Average	Calibration Factor	= Size
	<i>Readings</i>					
Organism	1	2	3	Average $\times$ Calibration Factor = Size		
Cocci						
Bacilli						
Yeast						
Protozoa						

	LENGTH OF MICROORGANISMS IN MICROMETERS ( $\mu\text{m}$ )					
	Number of ocular divisions			Average	Calibration Factor	= Size
	<i>Readings</i>					
Organism	1	2	3	Average $\times$ Calibration Factor = Size		
Bacilli						
Yeast						
Protozoa						

## Review Questions

1. Can the same calibration factor be used to determine the size of a micro-organism under all objectives? Explain.
2. If one stage micrometer division contains 12 ocular divisions, then the distance between two lines on the ocular micrometer is \_\_\_\_\_. Show your calculations.
3. If 1 mm is equal to  $1/25,400$  inch, convert the size of a bacterium measuring  $3 \mu\text{m} \times 1.5 \mu\text{m}$  to inches. Show your calculations.
4.  A comparative study of microbial size, using a stage micrometer, requires that you measure *Bacillus subtilis* in both the stained and unstained states.
  - a. Would you expect the measurements to be comparable in both preparations?
  - b. If a variation in size is expected, how would you account for this difference?

# Bacterial Staining

## LEARNING OBJECTIVES

Once you have completed the experiments in this section, you should be familiar with

1. The chemical and theoretical basis of biological staining.
2. Manipulative techniques of smear preparation.
3. Procedures for simple staining and negative staining.
4. The method for performing differential staining procedures, such as the Gram, acid-fast, capsule, and spore stains.

## Introduction

Visualization of microorganisms in the living state is quite difficult, not only because they are minute, but also because they are transparent and practically colorless when suspended in an aqueous medium. To study their properties and to divide microorganisms into specific groups for diagnostic purposes, biological stains and staining procedures in conjunction with light microscopy have become major tools in microbiology.

Chemically, a stain (dye) may be defined as an organic compound containing a benzene ring plus a chromophore and an auxochrome group (**Figure P3.1**).


The stain picric acid may be used to illustrate this definition (**Figure P3.2**).

The ability of a stain to bind to macromolecular cellular components such as proteins or nucleic


acids depends on the electrical charge found on the chromogen portion, as well as on the cellular component to be stained.

**Acidic stains** are anionic, which means that, on ionization of the stain, the chromogen portion exhibits a negative charge and therefore has a strong affinity for the positive constituents of the cell. Proteins, positively charged cellular components, will readily bind to and accept the color of the negatively charged, anionic chromogen of an acidic stain. Structurally, picric acid is an example of an acidic stain that produces an anionic chromogen, as illustrated in **Figure P3.3**.


**Basic stains** are cationic, because on ionization the chromogen portion exhibits a positive charge and therefore has a strong affinity for the negative constituents of the cell. Nucleic acids, negatively charged cellular components, will

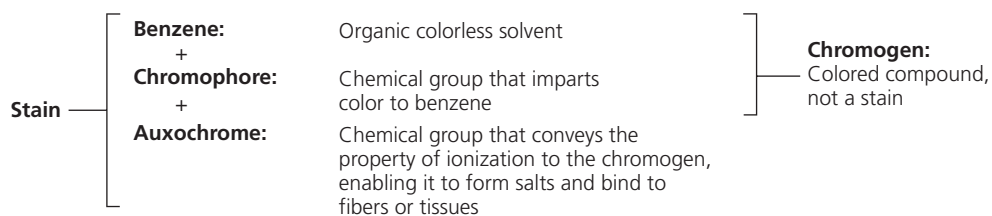


**Figure P3.1** Chemical composition of a stain

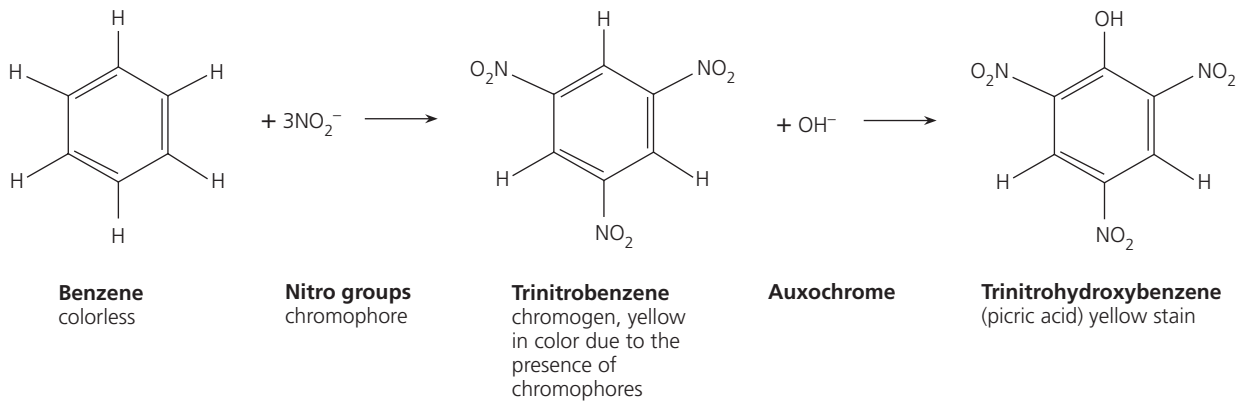

**Figure P3.2** Chemical formation of picric acid

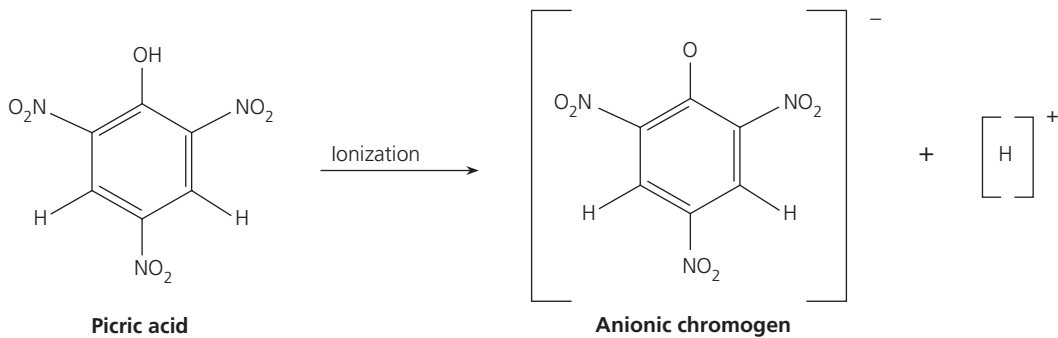

**Figure P3.3** Picric acid: an acidic stain

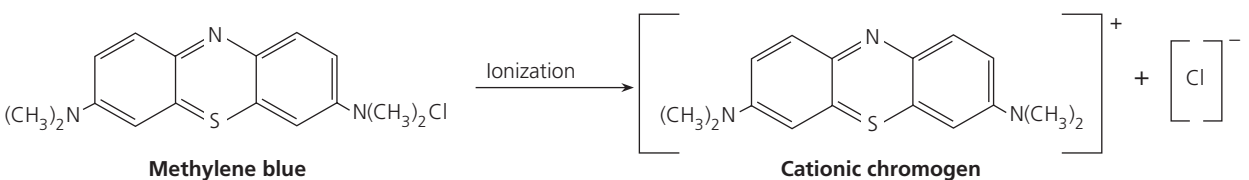

**Figure P3.4** Methylene blue: a basic stain


readily bind to and accept the color of the positively charged, cationic chromogen of a basic stain. Structurally, methylene blue is a basic stain that produces a cationic chromogen, as illustrated in **Figure P3.4**.

**Figure P3.5** is a summary of acidic and basic stains.


Basic stains are more commonly used for bacterial staining. The presence of a negative

charge on the bacterial surface acts to repel most acidic stains and thus prevent their penetration into the cell.

Numerous staining techniques are available for visualization, differentiation, and separation of bacteria in terms of morphological characteristics and cellular structures. A summary of commonly used procedures and their purposes is outlined in **Figure P3.6**.

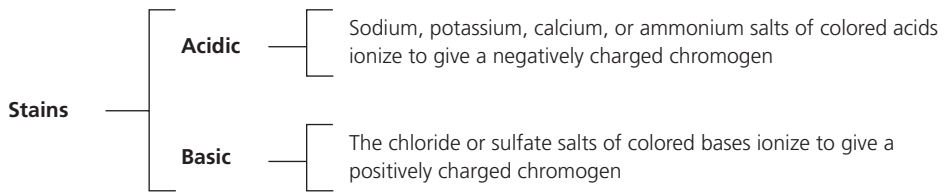

**Figure P3.5 Acidic and basic stains**

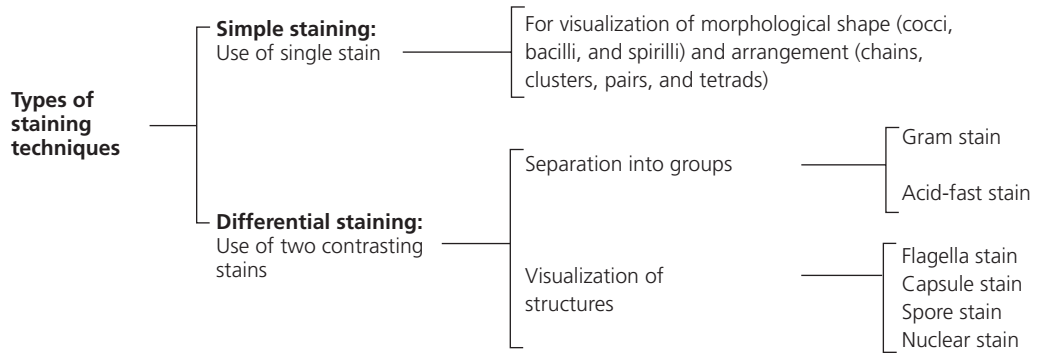

**Figure P3.6 Staining techniques**

*This page intentionally left blank*

# Preparation of Bacterial Smears

## LEARNING OBJECTIVE

Once you have completed this experiment, you should be able to

1. Prepare bacterial smears for the microscopic visualization of bacteria.


## Principle

Bacterial smears must be prepared prior to the execution of any of the staining techniques listed in Figure P3.6 on page 57. Although not difficult, the preparation requires adequate care. Meticulously follow the rules listed below.

- 1. Preparation of the glass microscope slide:** Clean slides are essential for the preparation of microbial smears. Grease or oil from the fingers on slides must be removed by washing the slides with soap and water or scouring powders such as Bon Ami<sup>®</sup>, followed by a water rinse and a rinse of 95% alcohol. After cleaning, dry the slides and place them on laboratory towels until ready for use. *Note: Remember to hold the clean slides by their edges.*
- 2. Labeling of slides:** Proper labeling of the slide is essential. The initials of the organism can be written on either end of the slide with a glassware marking pencil on the surface on which the smear is to be made. Care should be taken that the label does not come into contact with staining reagents.
- 3. Preparation of smear:** It is crucial to avoid thick, dense smears. A thick or dense smear occurs when too much of the culture is used in its preparation, which concentrates a large number of cells on the slide. This type of preparation diminishes the amount of light that can pass through and makes it difficult to visualize the morphology of single cells.
  - a. Broth cultures:** Resuspend the culture by tapping the tube with your finger. Depending on the size of the loop, one or two loopfuls should be applied to the center of the slide with a sterile inoculating loop and spread evenly over an area about the size of a dime. Set the smears on the laboratory table and allow to air-dry.
  - b. Cultures from solid medium:** Organisms cultured in a solid medium produce thick, dense surface growth and are not amenable to direct transfer to the glass slide. These cultures must be diluted by placing one or two loopfuls of water on the center of the slide in which the cells will be emulsified. Transfer of the cells requires the use of a sterile inoculating loop or a needle, if preferred. Only the tip of the loop or needle should touch the culture to prevent the transfer of too many cells. Suspension is accomplished by spreading the cells in a circular motion in the drop of water with the loop or needle. This helps to avoid cell clumping. The finished smear should occupy an area about the size of a nickel and should appear as a translucent, or semitransparent, confluent whitish film (**Figure 8.1**). At this point the smear should be allowed to dry completely. *Note: Do not blow on slide or wave it in the air.*
- 4. Heat fixation:** Unless fixed on the glass slide, the bacterial smear will wash away during the staining procedure. This is avoided by heat fixation, during which the bacterial proteins are coagulated and fixed to

*Note: Smears require only a small amount of the bacterial culture. A good smear is one that, when dried, appears as a thin whitish layer or film. The print of your textbook should be legible through the smear. Different techniques are used depending on whether the smear is made from a broth or solid-medium culture.*



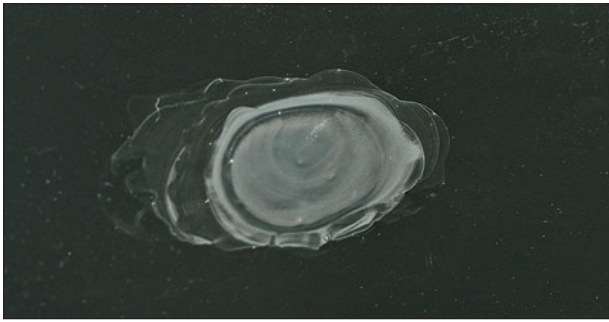

**Figure 8.1** A bacterial smear following fixation

the glass surface. Heat fixation is performed by the rapid passage of the air-dried smear two or three times over the flame of the Bunsen burner.

The preparation of a bacterial smear is illustrated in **Figure 8.2**.

## CLINICAL APPLICATION

### Proper Slide Preparation

Before any staining or visualization of a bacterial sample can take place, a proper smear must be prepared. A smear that is too thick may give a false result due to retention of dye that should have been rinsed away or because the thickness may prevent dye penetration. A smear that is too thin may have too few cells, increasing the time and energy to find the bacteria under magnification. Inconclusive results due to improperly prepared slides may have an impact on patient treatment and outcomes. Good smears are those which allow newsprint to be read through the smear.

## AT THE BENCH



## Materials

### Cultures

24-hour nutrient agar slant culture of *Bacillus cereus* and a 24-hour nutrient broth culture of *Staphylococcus aureus*.

## Equipment

Glass microscope slides, Bunsen burner, inoculating loop and needle, and glassware marking pencil.

## Procedure

### Smears from a Broth Medium

Label three clean slides with the initials of the organism, and number them 1, 2, and 3. Resuspend the sedimented cells in the broth culture by tapping the culture tube with your finger. The next four steps of this procedure are illustrated in Figure 8.2a and c:

1. With a sterile loop, place one loopful of culture on Slide 1, two loopfuls on Slide 2, and three loopfuls on Slide 3, respectively.
2. With a circular movement of the loop, spread the cell suspension into an area approximately the size of a dime.
3. Allow the slide to air-dry completely.
4. Heat fix the preparation. *Note: Pass the air-dried slide through the outer portion of the Bunsen flame to prevent overheating, which can distort the morphology through plasmolysis of the cell wall.*

Examine each slide for the confluent, whitish film or haze and record your results in the Lab Report.

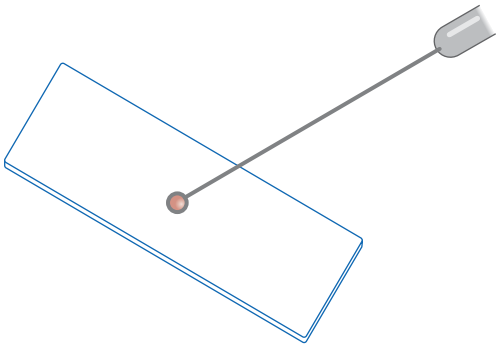
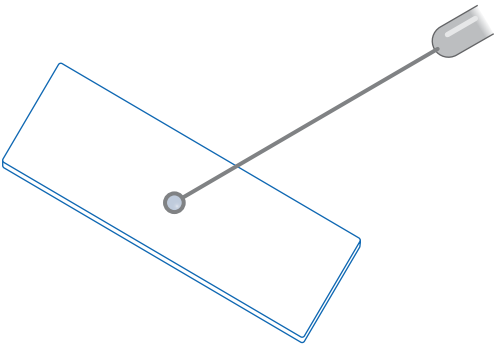
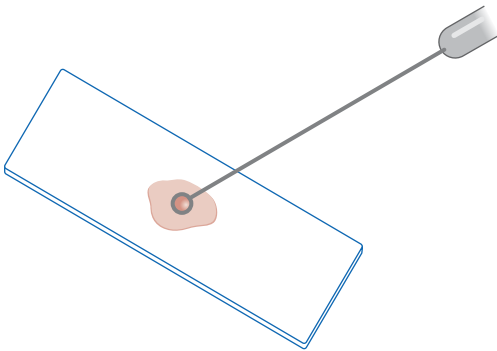
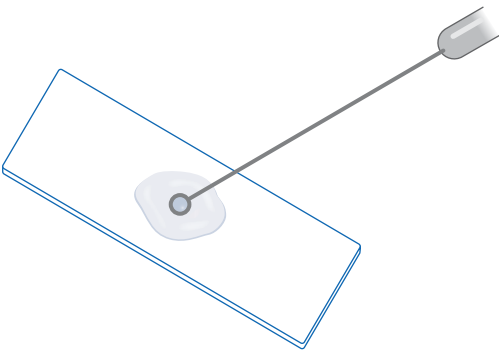
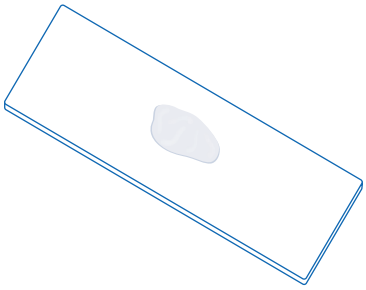
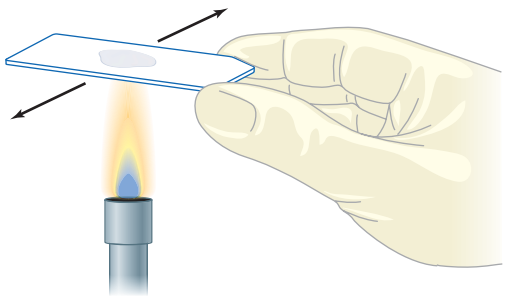
### Smears from a Solid Medium

Label four clean slides with the initials of the organism. Label Slides 1 and 2 with an L for loop, and Slides 3 and 4 with an N for needle. The next four steps of this procedure are illustrated in Figure 8.2b and c:

1. Using a loop, place one to two loops of water on each slide.
2. With a sterile loop, touch the entire loop to the culture and emulsify the cells in water on Slide 1. Then, with a sterile loop, just touch the tip of the loop to the culture and emulsify it in the water on Slide 2. Repeat Steps 1 and 2 using a sterile inoculating needle on Slides 3 and 4.
3. Allow all slides to air-dry completely.
4. Heat fix the preparation.

Examine each slide for the confluent, whitish film or haze and record your results in the Lab Report.

## PROCEDURE

(a) From broth medium	(b) From solid medium
 <p><b>1</b> Place one to two loopfuls of the cell suspension on the clean slide.</p>	 <p><b>1</b> Place one to two loopfuls of water on the center of the slide.</p>
 <p><b>2</b> With a circular movement of the loop, spread the suspension into a thin area approximately the size of a dime.</p>	 <p><b>2</b> Transfer a small amount of the bacterial inoculum from the slant culture into the drop of water. Spread both into a thin area approximately the size of a nickel.</p>
(c) Fixation for solid and broth media	
 <p><b>3</b> Allow the smear to air-dry.</p>	 <p><b>4</b> While holding the slide at one end, quickly pass the smear over the flame of the Bunsen burner two to three times.</p>

**Figure 8.2** Bacterial smear preparation

*This page intentionally left blank*

Name: \_\_\_\_\_


Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

Select the preparation that you think is best and ask your instructor to comment on your choice. Remember that printed material should be legible through a good smear. Indicate by slide number the consistency of smears from both broth and solid cultures that you considered best.

Broth culture \_\_\_\_\_ Solid culture: Loop \_\_\_\_\_ Needle \_\_\_\_\_

## Review Questions

1. Why are thick or dense smears less likely to provide a good smear preparation for microscopic evaluation?
2. Why is it essential that smears be air-dried? Why can't they be gently heated over a flame to speed up the drying process?
3. Why should you be careful not to overheat the smear during the heat-fixing process?
4.  Why do you think the presence of grease or dirt on a glass slide will result in a poor smear preparation? Cite two or three reasons.

*This page intentionally left blank*

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be able to

1. Perform a simple staining procedure.
2. Compare the morphological shapes and arrangements of bacterial cells.

## Principle

In **simple staining**, the bacterial smear is stained with a single reagent, which produces a distinctive contrast between the organism and its background. Basic stains with a positively charged chromogen are preferred because bacterial nucleic acids and certain cell wall components carry a negative charge that strongly attracts and binds to the cationic chromogen. The purpose of simple staining is to elucidate the morphology and arrangement of bacterial cells (**Figure 9.1**). The most commonly used basic stains are methylene blue, crystal violet, and carbol fuchsin.

## CLINICAL APPLICATION

### Quick and Simple Stain

Simple stains are relatively quick and useful methods of testing for the presence of, determining the shape of, or determining the numbers of bacteria present in a sample. Generally involving a single staining step, simple staining methods are not considered differential or diagnostic and will have limited uses. However, this is a quick procedure for determining whether a clinical sample has the presence of a foreign bacterial pathogen.

## AT THE BENCH

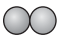






## Materials



### Cultures

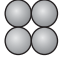24-hour nutrient agar slant cultures of *Escherichia coli* and *Bacillus cereus* and a 24-hour nutrient broth culture of *Staphylococcus aureus*. Alternatively, use the smears prepared in Experiment 8.




**Cocci** are spherical in shape.

- |                                                                                     |                    |                       |
|-------------------------------------------------------------------------------------|--------------------|-----------------------|
|    | (a) Diplococcus    | Diplo = pair          |
|    | (b) Streptococcus  | Strepto = chain       |
|  | (c) Staphylococcus | Staphylo = cluster    |
|  | (d) Tetrad         | Tetrad = packet of 4  |
|  | (e) Sarcina        | Sarcina = packet of 8 |

**Bacilli** are rod-shaped.

- |                                                                                     |                     |                 |
|-------------------------------------------------------------------------------------|---------------------|-----------------|
|  | (a) Diplobacillus   | Diplo = pair    |
|  | (b) Streptobacillus | Strepto = chain |

**Spiral bacteria** are rigid or flexible.

- |                                                                                     |                                           |
|-------------------------------------------------------------------------------------|-------------------------------------------|
|  | (a) Vibrios are curved rods.              |
|  | (b) Spirilla are helical and rigid.       |
|  | (c) Spirochetes are helical and flexible. |

**Figure 9.1** Bacterial shapes and arrangements

## Reagents

Methylene blue, crystal violet, and carbol fuchsin.

## Equipment

Bunsen burner, inoculating loop, staining tray, microscope, lens paper, bibulous (highly absorbent) paper, and glass slides.

## Procedure

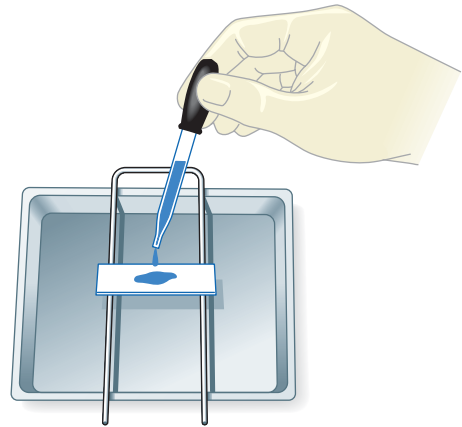
1. Prepare separate bacterial smears of the organisms following the procedure described in Experiment 8. *Note: All smears must be heat fixed prior to staining.*

## Simple Staining


The following steps are illustrated in **Figure 9.2**.

1. Place a slide on the staining tray and flood the smear with one of the indicated stains, using the appropriate exposure time for each: carbol fuchsin, 15 to 30 seconds; crystal violet, 20 to 60 seconds; methylene blue (shown in Figure 9.2), 1 to 2 minutes.
2. Gently wash the smear with tap water to remove excess stain. During this step, hold the slide parallel to the stream of water; in this way you can reduce the loss of organisms from the preparation.
3. Using bibulous paper, blot dry but *do not* wipe the slide.
4. Repeat this procedure with the remaining two organisms, using a different stain for each.
5. Examine all stained slides under oil immersion.
6. In the chart provided in the Lab Report, complete the following:
  - a. Draw a representative field for each organism. Refer to page xvi for proper drawing procedure.
  - b. Describe the morphology of the organisms with reference to their shapes (bacilli, cocci, spirilla) and arrangements (chains, clusters, pairs). Refer to the photographs in **Figure 9.3**.


## PROCEDURE



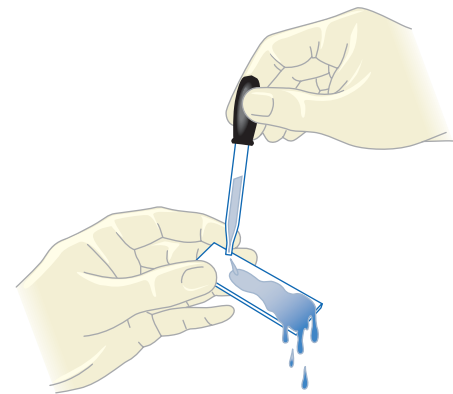- 1 Place slide on the staining tray and flood the smear with methylene blue. Allow 1 to 2 minutes of exposure to the stain.



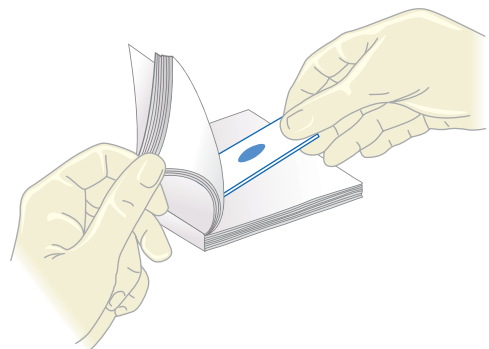- 2 Gently wash the smear with tap water.



- 3 Blot the slide dry with bibulous paper.


**Figure 9.2** Simple staining procedure



(a) Bacilli and diplobacilli (rod-shaped) bacteria



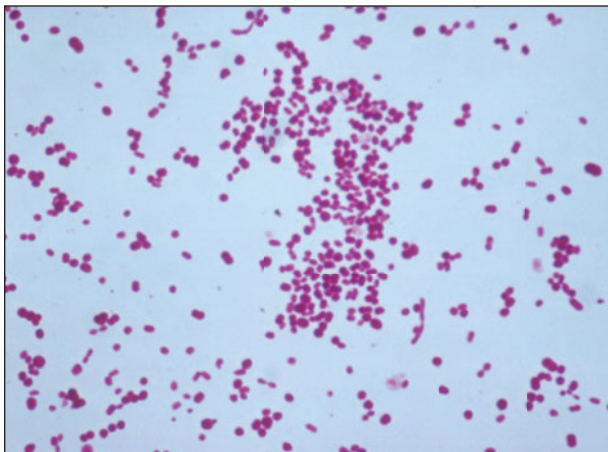(b) Spirilla (spiral-shaped) bacteria



(c) Cocci (spherical-shaped) bacteria: *Staphylococcus*

**Figure 9.3** Micrographs showing bacteria morphology

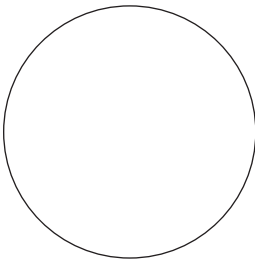
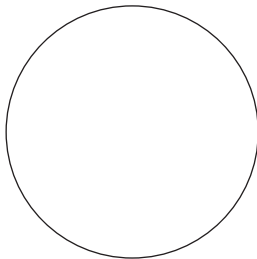
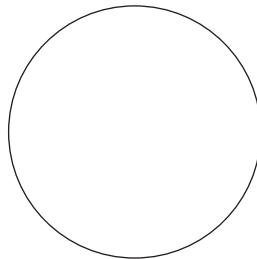


*This page intentionally left blank*

Name: \_\_\_\_\_



Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

	Methylene Blue	Crystal Violet	Carbol Fuchsin
Draw a representative field.			
Organism	_____	_____	_____
Cell morphology:			
Shape	_____	_____	_____
Arrangement	_____	_____	_____
Cell color	_____	_____	_____

## Review Questions

1. Why are basic dyes more effective for bacterial staining than acidic dyes?
  
  
  
  
  
  
2. Can simple staining techniques be used to identify more than the morphological characteristics of microorganisms? Explain.

3.  During the performance of the simple staining procedure, you failed to heat fix your *E. coli* smear preparation. Upon microscopic examination, how would you expect this slide to differ from the correctly prepared slides?
4.  During a coffee break, your friend spills coffee on your lab coat and the fabric is discolored. Is this a true biological stain or simply a compound capable of imparting color? Explain your rationale.

### LEARNING OBJECTIVES

Once you have completed this experiment, you should be able to

1. Perform a negative staining procedure.
2. Understand the benefit obtained from visualizing unstained microorganisms.

## Principle

**Negative staining** requires the use of an acidic stain such as India ink or nigrosin. The acidic stain, with its negatively charged chromogen, will not penetrate the cells because of the negative charge on the surface of bacteria. Therefore, the unstained cells are easily discernible against the colored background.

The practical application of negative staining is twofold. First, since heat fixation is not required and the cells are not subjected to the distorting effects of chemicals and heat, their natural size and shape can be seen. Second, it is possible to observe bacteria that are difficult to stain, such as some spirilla. Because heat fixation is not done during the staining process, keep in mind that the organisms are not killed and *slides should be handled with care*. **Figure 10.1** shows a negative stain of bacilli.



**Figure 10.1** Negative staining: *Bacilli* (1000 ×)

### CLINICAL APPLICATION

#### Detecting Encapsulated Invaders

The principle application of negative staining is to determine if an organism possesses a capsule (a gelatinous outer layer that makes the microorganism more virulent), although it can also be used to demonstrate spore formation. The technique is frequently used in the identification of fungi such as *Cryptococcus neoformans*, an important infectious agent found in bird dropping that is linked to meningial and lung infections in humans.

### AT THE BENCH



## Materials

### Cultures

24-hour agar slant cultures of *Micrococcus luteus*, *Bacillus cereus*, and *Aquaspirillum itersonii*.

### Reagent

Nigrosin.

### Equipment

Bunsen burner, inoculating loop, staining tray, glass slides, lens paper, and microscope.


## Procedure

Steps 1–4 are illustrated in **Figure 10.2**.

1. Place a small drop of nigrosin close to one end of a clean slide.
2. Using aseptic technique, place a loopful of inoculum from the *M. luteus* culture in the drop of nigrosin and mix.

- Place a slide against the drop of suspended organisms at a 45° angle and allow the drop to spread along the edge of the applied slide.
- Push the slide away from the drop of suspended organisms to form a thin smear. Air-dry. *Note: Do not heat fix the slide.*
- Repeat Steps 1–4 for slide preparations of *B. cereus* and *A. itersonii*.
- Examine the slides under oil immersion, and record your observations in the Lab Report.

## PROCEDURE

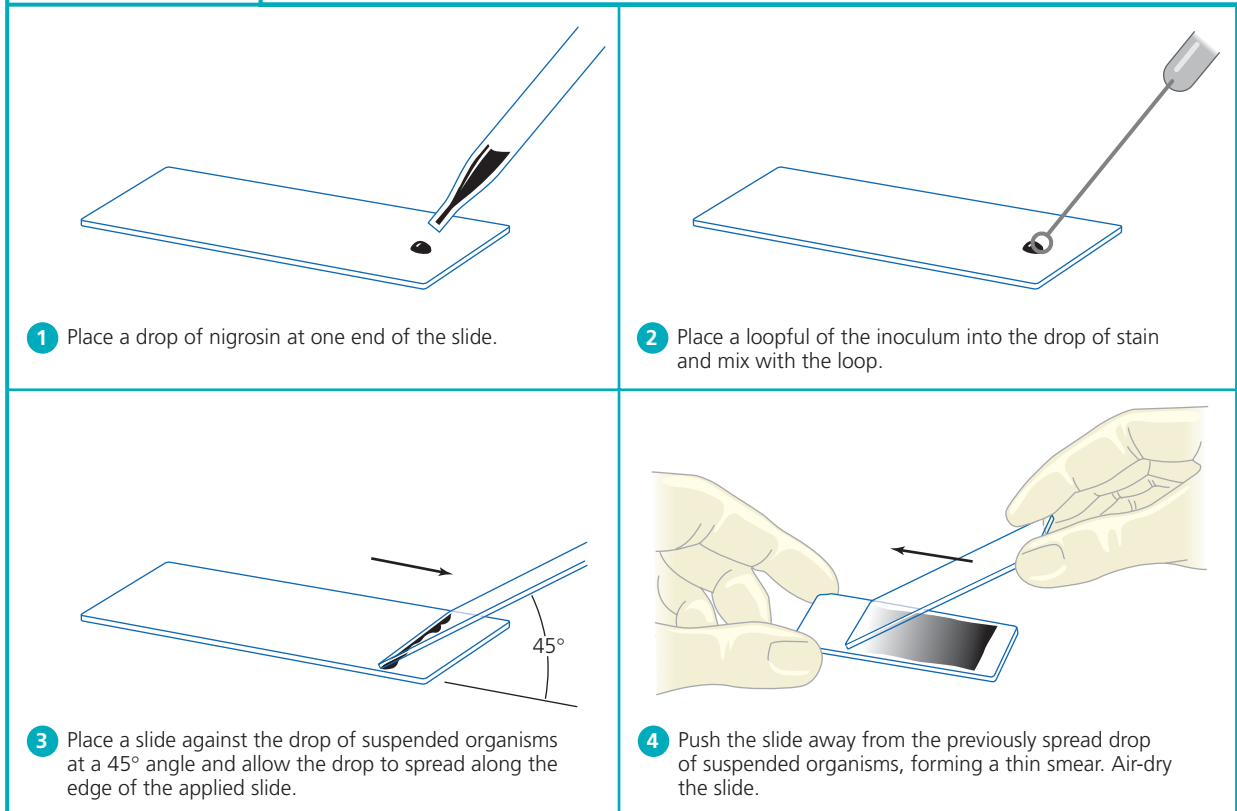




**Figure 10.2** Negative staining procedure

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

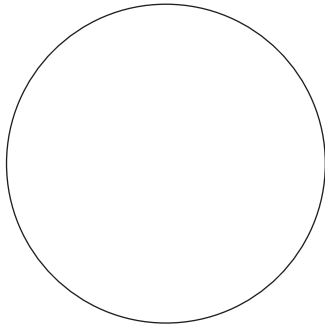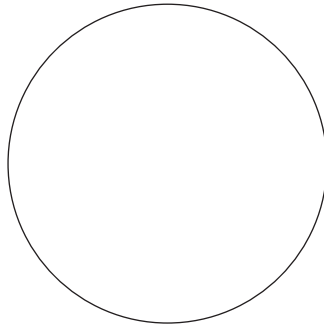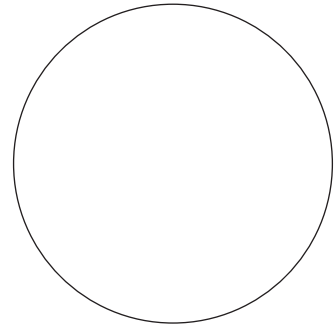1. Draw representative fields of your microscopic observations.

*M. luteus**B. cereus**A. itersonii*

2. Describe the microscopic appearance of the different bacteria using the chart below.

Organism	<i>M. luteus</i>	<i>B. cereus</i>	<i>A. itersonii</i>
Shape			
Arrangement			
Magnification			

## Review Questions

1. Why can't methylene blue be used in place of nigrosin for negative staining? Explain.
2. What are the practical advantages of negative staining?
3. Why doesn't nigrosin penetrate bacterial cells?

## LEARNING OBJECTIVES

Once you have completed this experiment, you should understand

1. The chemical and theoretical basis for differential staining procedures.
2. The chemical basis for the Gram stain.
3. The procedure for differentiating between two principal groups of bacteria: gram-positive and gram-negative.

## Principle

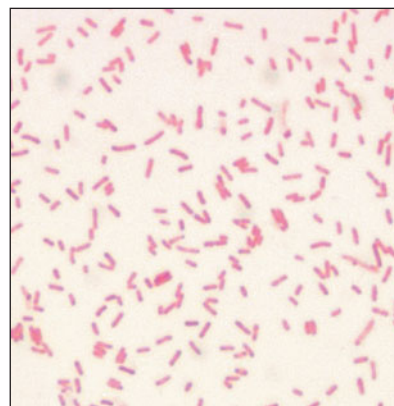
**Differential staining** requires the use of at least four chemical reagents that are applied sequentially to a heat-fixed smear. The first reagent is called the **primary stain**. Its function is to impart its color to all cells. The second stain is a mordant used to intensify the color of the primary stain. In order to establish a color contrast, the third reagent used is the decolorizing agent. Based on the chemical composition of cellular components, the **decolorizing agent** may or may not remove the primary stain from the entire cell or only from certain cell structures. The final reagent, the counterstain, has a contrasting color to that of the primary stain. Following decolorization, if the primary stain is not washed out, the **counterstain** cannot be absorbed and the cell or its components will retain the color of the primary stain. If the primary stain is removed, the decolorized cellular components will accept and assume the contrasting color of the counterstain. In this way, cell types or their structures can be distinguished from each other on the basis of the stain that is retained.

The most important differential stain used in bacteriology is the **Gram stain**, named after Dr. Hans Christian Gram. It divides bacterial cells into two major groups, gram-positive and gram-negative, which makes it an essential tool for classification and differentiation of microorganisms. **Figure 11.1** shows gram-positive and gram-negative stained cells. The Gram stain reaction is based on the difference in the chemical composition of bacterial cell walls.

Gram-positive cells have a thick peptidoglycan layer, whereas the peptidoglycan layer in gram-negative cells is much thinner and surrounded by outer lipid-containing layers. Peptidoglycan is mainly a polysaccharide composed of two chemical subunits found only in the bacterial cell wall. These subunits are *N*-acetylglucosamine and *N*-acetylmuramic acid. With some organisms, as the adjacent layers of peptidoglycan are formed, they are cross-linked by short chains of peptides by means of a transpeptidase enzyme, resulting in the shape and rigidity of the cell wall. In the case of gram-negative bacteria and several of the gram-positive such as the *Bacillus*, the cross-linking of the peptidoglycan layer is direct because the




(a) Gram-positive stain of streptococci



(b) Gram-negative stain of *E. coli*

**Figure 11.1** Gram-stained cells



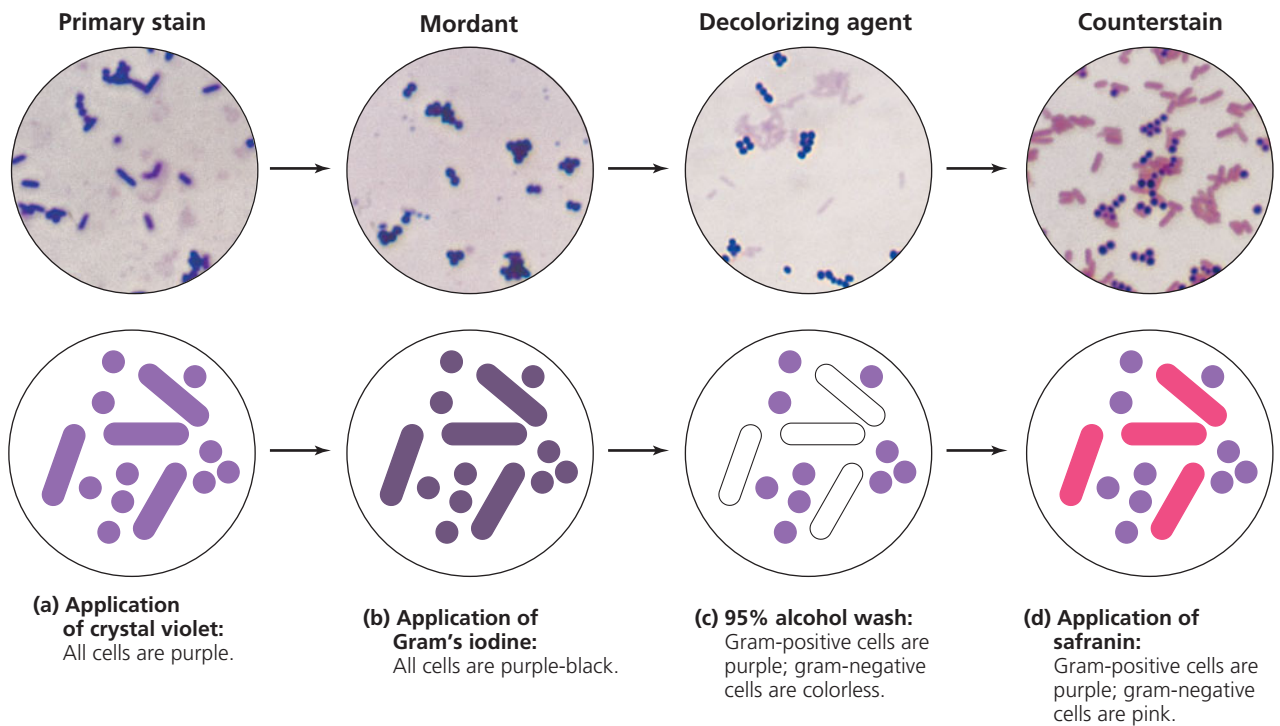

**Figure 11.2** Microscopic observation of cells following steps in the Gram staining procedure

bacteria do not have short peptide tails. Early experiments have shown that a gram-positive cell denuded of its cell wall by the action of lysozyme or penicillin will stain gram-negative.

The Gram stain uses four different reagents. Descriptions of these reagents and their mechanisms of action follow. **Figure 11.2** shows the microscopic appearance of cells at each step of the Gram staining procedure.

## Primary Stain

**Crystal Violet (Hucker's)** This violet stain is used first and stains all cells purple.

## Mordant

**Gram's Iodine** This reagent serves not only as a killing agent but also as a mordant, a substance that increases the cells' affinity for a stain. It does this by binding to the primary stain, thus forming an insoluble complex. The resultant crystal-violet-iodine (CV-I) complex serves to intensify the color of the stain. At this point, all cells will appear purple-black.

## Decolorizing Agent

**Ethyl Alcohol, 95%** This reagent serves a dual function as a protein-dehydrating agent and as a lipid solvent. Its action is determined by two factors, the concentration of lipids and the thickness of the peptidoglycan layer in bacterial cell walls. In gram-negative

cells, the alcohol increases the porosity of the cell wall by dissolving the lipids in the outer layers. Thus, the CV-I complex can be more easily removed from the thinner and less highly cross-linked peptidoglycan layer. Therefore, the washing-out effect of the alcohol facilitates the release of the unbound CV-I complex, leaving the cells colorless or unstained. The much thicker peptidoglycan layer in gram-positive cells is responsible for the more stringent retention of the CV-I complex, as the pores are made smaller due to the dehydrating effect of the alcohol. Thus, the tightly bound primary stain complex is difficult to remove, and the cells remain purple. *Note: Be careful not to over-decolorize the smear with alcohol.*

## Counterstain

**Safranin** This is the final reagent, used to stain pink those cells that have been previously decolorized. Since only gram-negative cells undergo decolorization, they may now absorb the counterstain. Gram-positive cells retain the purple color of the primary stain.

The preparation of adequately stained smears requires that you bear in mind the following precautions:

1. The most critical phase of the procedure is the decolorization step, which is based on the ease with which the CV-I complex is released from the cell. Remember that over-decolorization will result in loss of the primary stain, causing

gram-positive organisms to appear gram-negative. Under-decolorization, however, will not completely remove the CV-I complex, causing gram-negative organisms to appear gram-positive. Strict adherence to all instructions will help remedy part of the difficulty, but individual experience and practice are the keys to correct decolorization.

2. It is imperative that, between applications of the reagents, slides be thoroughly washed under running water or water applied with an eyedropper. This removes excess reagent and prepares the slide for application of the subsequent reagent.
3. The best Gram-stained preparations are made with fresh cultures, that is, not older than 24 hours. As cultures age, especially in the case of gram-positive cells, the organisms tend to lose their ability to retain the primary stain and may appear to be **gram-variable**; that is, some cells will appear purple, while others will appear pink.

## CLINICAL APPLICATION

### Gram Staining: The First Diagnostic Test

The Gram stain is a diagnostic staining procedure that can be done on body fluids, tissue biopsies, throat cultures, samples from abscesses when infection is suspected, and more. Clinically important results are obtained much more rapidly from staining than from culturing the specimen. The results of the Gram stain will aid a clinical lab in determining which additional tests may be required for identification of the bacterial strain in question. Once the bacterial gram type, shape, and orientation are determined, it expedites the appropriate choice of antibiotic needed to treat the patient.

## AT THE BENCH



## Materials

### Cultures

24-hour nutrient agar slant cultures of *Escherichia coli*, *Staphylococcus aureus*, and *Bacillus cereus*.

### Reagents

Crystal violet, Gram's iodine, 95% ethyl alcohol, and safranin.

## Equipment

Bunsen burner, inoculating loop or needle, staining tray, glass slides, bibulous paper, lens paper, and microscope.

## Procedure

### Smear Preparation

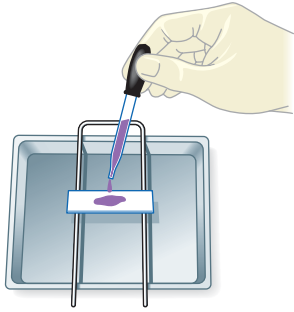
1. Obtain four clean glass slides.
2. Using aseptic technique, prepare a smear of each of the three organisms and on the remaining slide prepare a smear consisting of a mixture of *S. aureus* and *E. coli*. Do this by placing a drop of water on the slide, and then transferring each organism separately to the drop of water with a sterile, cooled loop. Mix and spread both organisms by means of a circular motion of the inoculating loop.
3. Allow smears to air-dry and then heat fix in the usual manner.

### Gram Staining


The following steps are shown in **Figure 11.3**.

1. *Gently* flood smears with crystal violet and let stand for 1 minute.
2. *Gently* wash with tap water.
3. *Gently* flood smears with the Gram's iodine mordant and let stand for 1 minute.
4. *Gently* wash with tap water.
5. Decolorize with 95% ethyl alcohol. *Note: Do not over-decolorize.* Add reagent drop by drop until the alcohol runs almost clear, showing only a blue tinge.
6. *Gently* wash with tap water.
7. Counterstain with safranin for 45 seconds.
8. *Gently* wash with tap water.
9. Blot dry with bibulous paper and examine under oil immersion.
10. As you observe each slide under oil immersion, complete the chart provided in the Lab Report.
  - a. Draw a representative microscopic field.
  - b. Describe the cells according to their morphology and arrangement.
  - c. Describe the color of the stained cells.
  - d. Classify the organism as to the Gram reaction: gram-positive or gram-negative.


## PROCEDURE



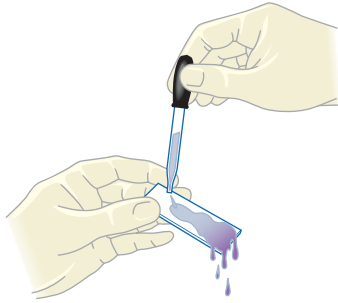
- 1 Gently stain with crystal violet for 1 minute.



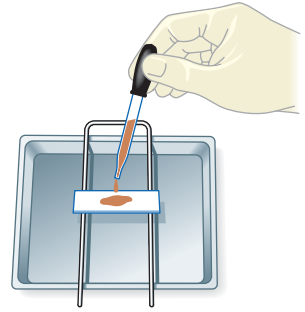
- 2 Gently wash off the stain with tap water.



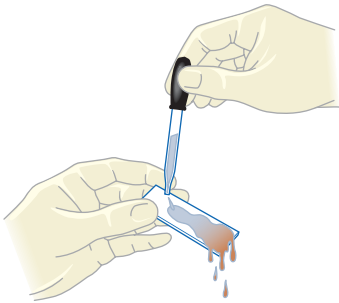
- 3 Gently apply Gram's iodine for 1 minute.



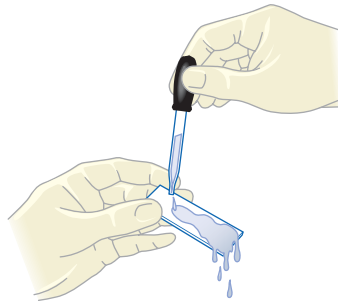
- 4 Gently wash off the Gram's iodine with tap water.



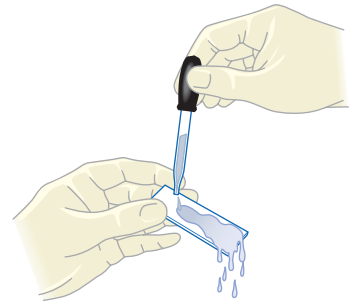
- 5 Add 95% alcohol drop by drop until the alcohol runs almost clear.



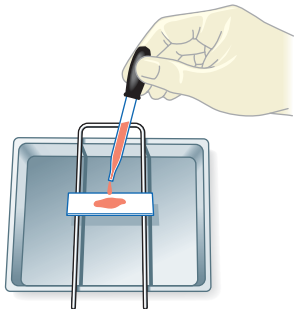
- 6 Gently wash off the 95% alcohol with tap water.



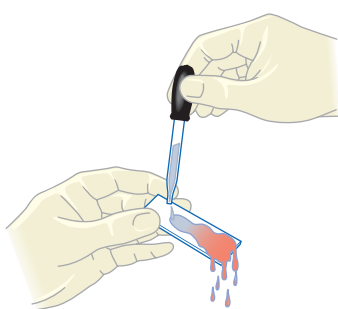- 7 Counterstain with safranin for 45 seconds.



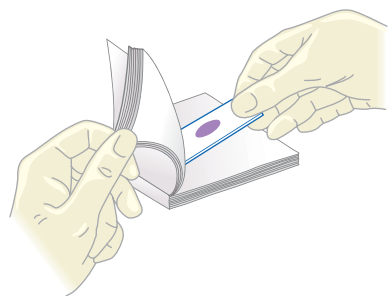- 8 Gently wash off the safranin with tap water.



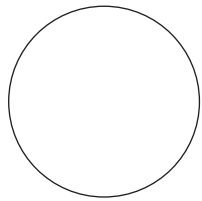
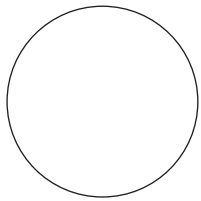
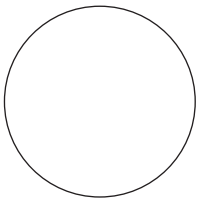
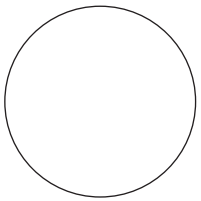
- 9 Blot dry with bibulous paper.

**Figure 11.3** Gram staining procedure

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results


	<i>E. coli</i>	<i>B. cereus</i>	<i>S. aureus</i>	<i>Mixture</i>
Draw a representative field.				
Cell morphology:				
Shape	_____	_____	_____	_____
Arrangement	_____	_____	_____	_____
Cell color	_____	_____	_____	_____
Gram reaction	_____	_____	_____	_____

## Review Questions

1. What are the advantages of differential staining procedures over the simple staining technique?
  
  
  
  
  
  
  
  
  
  
2. Cite the purpose of each of the following reagents in a differential staining procedure.
  - a. Primary stain:
  
  
  
  
  
  
  
  
  
  
  - b. Mordant:

c. Decolorizing agent:

d. Counterstain:

3. Why is it essential that the primary stain and the counterstain be of contrasting colors?
  
  
  
  
  
  
  
  
  
  
4. Which is the most crucial step in the performance of the Gram staining procedures? Explain.
  
  
  
  
  
  
  
  
  
  
5.  Because of a snowstorm, your regular laboratory session was canceled and the Gram staining procedure was performed on cultures incubated for a longer period of time. Examination of the stained *B. cereus* slides revealed a great deal of color variability, ranging from an intense blue to shades of pink. Account for this result.

## LEARNING OBJECTIVES

Once you have completed this experiment, you should understand

1. The chemical basis of the acid-fast stain.
2. The procedure for differentiating bacteria into acid-fast and non-acid-fast groups.

## Principle

While the majority of bacterial organisms are stainable by either simple or Gram staining procedures, a few genera, particularly the members of the genus *Mycobacterium*, are visualized much better by the **acid-fast** method. Since *M. tuberculosis* and *M. leprae* represent bacteria that are pathogenic to humans, the stain is of diagnostic value in identifying these organisms.

The characteristic difference between mycobacteria and other microorganisms is the presence of a thick, waxy (lipoidal) wall that makes penetration by stains extremely difficult. The mycobacteria tend to clump together, and it is difficult to identify individual cells in stained preparations if this clumping effect occurs. To avoid or minimize this phenomenon requires careful preparation of the smear. Place a small drop of water on the slide, suspend the culture in the water, and mix the suspension thoroughly to dislodge and disperse some of the cells. Once the stain has penetrated, however, it cannot be readily removed even with the vigorous use of acid-alcohol as a decolorizing agent (unlike the 95% ethyl alcohol used in the Gram stain). Because of this property, these organisms are called acid-fast, while all other microorganisms, which are easily decolorized by acid-alcohol, are non-acid-fast.

The acid-fast stain uses three different reagents.

## Primary Stain

**Carbol Fuchsin** Unlike cells that are easily stained by ordinary aqueous stains, most species of mycobacteria are not stainable with common

dyes such as methylene blue and crystal violet. Carbol fuchsin, a dark red stain in 5% phenol that is soluble in the lipoidal materials that constitute most of the mycobacterial cell wall, does penetrate these bacteria and is retained. Penetration is further enhanced by the application of heat, which drives the carbol fuchsin through the lipoidal wall and into the cytoplasm. This application of heat is used in the **Ziehl-Neelsen method**. The **Kinyoun method**, a modification of the Ziehl-Neelsen method, circumvents the use of heat by addition of a wetting agent (Tergitol®) to this stain, which reduces surface tension between the cell wall of the mycobacteria and the stain. Following application of the primary stain, all cells will appear red.

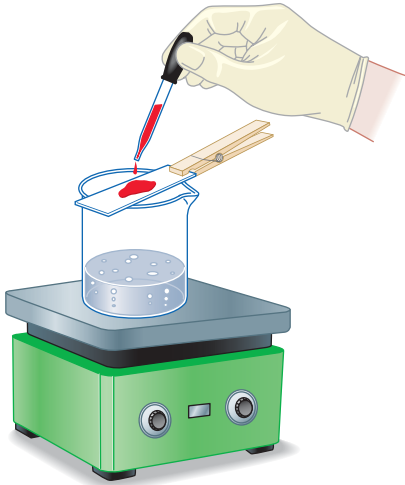
## Decolorizing Agent

**Acid-Alcohol (3% HCl + 95% Ethanol)** Prior to decolorization, the smear is cooled, which allows the waxy cell substances to harden. On application of acid-alcohol, acid-fast cells will be resistant to decolorization since the primary stain is more soluble in the cellular waxes than in the decolorizing agent. In this event, the primary stain is retained and the mycobacteria will stay red. This is not the case with non-acid-fast organisms, which lack cellular waxes. The primary stain is more easily removed during decolorization, leaving these cells colorless or unstained.


## Counterstain

**Methylene Blue** This is used as the final reagent to stain previously decolorized cells. As only non-acid-fast cells undergo decolorization, they may now absorb the counterstain and take on its blue color, while acid-fast cells retain the red of the primary stain.

## PROCEDURE




or

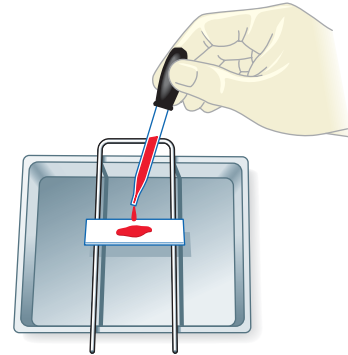


**1a Heat method:** Apply carbol fuchsin and steam over a beaker of boiling water that is placed on a hot plate for 5 minutes. Do not allow the stain to evaporate.

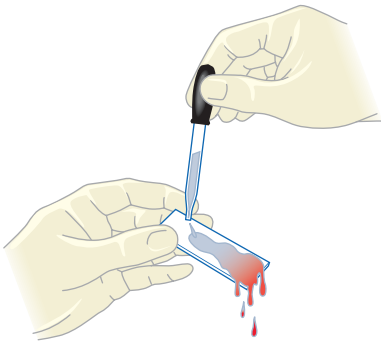
**1b Heatless method:** Apply carbol fuchsin with Tergitol® for 5 to 10 minutes.



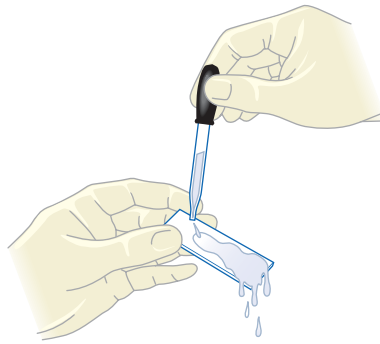
**2** Cool and wash off stain with tap water.



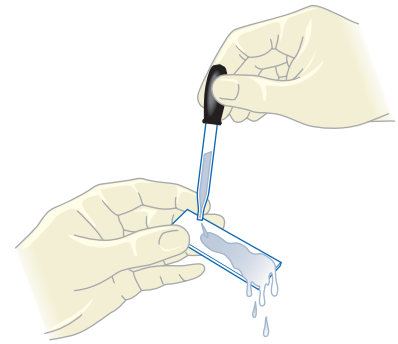
**3** Add acid-alcohol drop by drop until the alcohol runs almost clear.



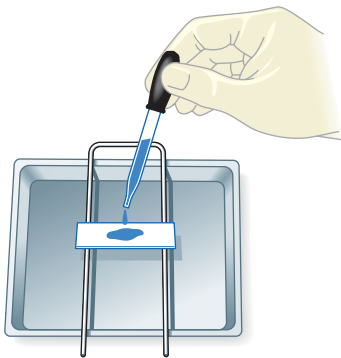
**4** Wash off the acid-alcohol with tap water.



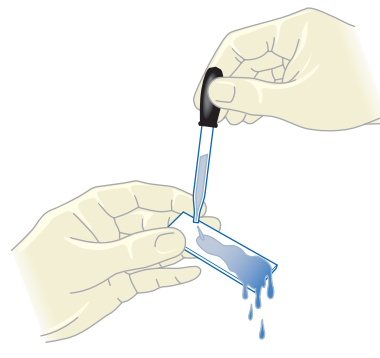**5** Counterstain with methylene blue for 2 minutes.



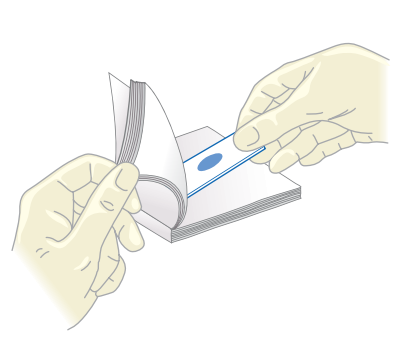**6** Wash off the methylene blue with tap water.



**7** Blot the slide dry with bibulous paper.

**Figure 12.1** Acid-fast staining procedure

## CLINICAL APPLICATION

### Diagnosing Leprosy and Lung Infections

The cell walls of bacteria belonging to the genera *Mycobacterium* or *Nocardia* contain mycolic acid and are resistant to penetration by water-soluble stains such as the Gram stain, which can lead to a “false” gram-positive result. The medical importance of the acid-fast stain is for the diagnosis of the *Mycobacterium* species, which cause tuberculosis, leprosy, and other infections. The genus *Nocardia*, which is the causative agent for lung infections, will also be identified by the acid-fast staining method.

## AT THE BENCH



## Materials

### Cultures

72- to 96-hour Trypticase soy broth culture of *Mycobacterium smegmatis* and 18- to 24-hour culture of *Staphylococcus aureus*.

### Reagents

Carbol fuchsin, acid-alcohol, and methylene blue.

### Equipment

Bunsen burner, hot plate, 250-ml beaker, inoculating loop, glass slides, bibulous paper, lens paper, staining tray, and microscope.

## Procedure

### Smear Preparation

1. Obtain three clean glass slides.
2. Using aseptic technique, prepare a bacterial smear of each organism plus a third mixed smear of *M. smegmatis* and *S. aureus*.
3. Allow smears to air-dry and then heat fix in the usual manner.

### Acid-Fast Staining


Steps 1–7 are pictured in **Figure 12.1**.

1. **a.** Flood smears with carbol fuchsin and place over a beaker of water on a warm hot

plate, allowing the preparation to steam for 5 minutes. *Note: Do not allow stain to evaporate; replenish stain as needed. Also, prevent stain from boiling by adjusting the hot-plate temperature.*

- b.** For a heatless method, flood the smear with carbol fuchsin containing Tergitol® for 5 to 10 minutes.
2. Wash with tap water. Heated slides must be cooled prior to washing.
3. Decolorize with acid-alcohol, adding the reagent drop by drop until the alcohol runs almost clear with a slight red tinge.
4. Wash with tap water.
5. Counterstain with methylene blue for 2 minutes.
6. Wash smear with tap water.
7. Blot dry with bibulous paper and examine under oil immersion.
8. In the chart provided in the Lab Report, complete the following:
  - a. Draw a representative microscopic field for each preparation.
  - b. Describe the cells according to their shapes and arrangements.
  - c. Describe the color of the stained cells.
  - d. Classify the organisms as to reaction: acid-fast or non-acid-fast.

Refer to **Figure 12.2** for a photograph of an acid-fast stain.

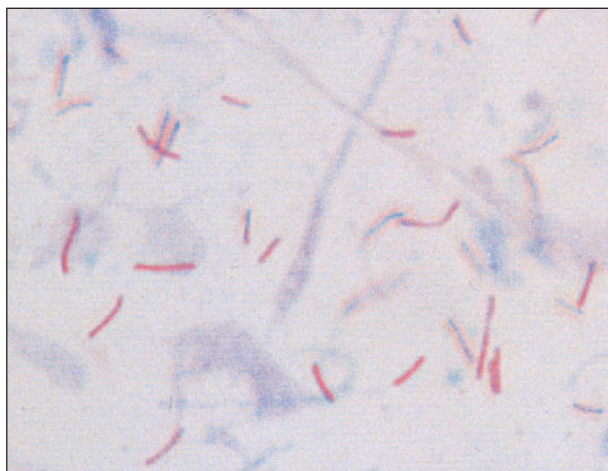

**Figure 12.2** Acid-fast stain of mycobacteria






*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

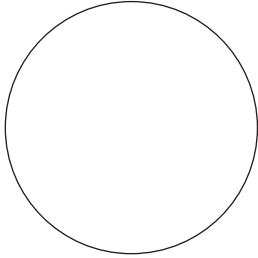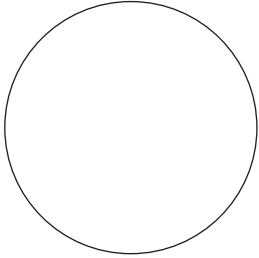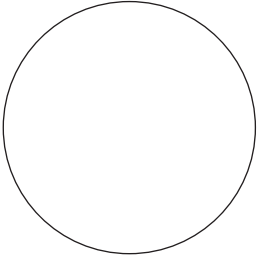
	<i>M. smegmatis</i>	<i>S. aureus</i>	Mixture
Draw a representative field.			
Cell morphology:			
Shape	_____	_____	_____
Arrangement	_____	_____	_____
Cell color	_____	_____	_____
Acid-fast reaction	_____	_____	_____

## Review Questions

1. Why must heat or a surface-active agent be used with application of the primary stain during acid-fast staining?
2. Why is acid-alcohol rather than ethyl alcohol used as a decolorizing agent?

3. What is the specific diagnostic value of this staining procedure?

4.  Why is the application of heat or a surface-active agent not required during the application of the counterstain in acid-fast staining?

5.  A child presents symptoms suggestive of tuberculosis, namely a respiratory infection with a productive cough. Microscopic examination of the child's sputum reveals no acid-fast rods. However, examination of gastric washings reveals the presence of both acid-fast and non-acid-fast bacilli. Do you think the child has active tuberculosis? Explain.

# Differential Staining for Visualization of Bacterial Cell Structures

EXPERIMENT

13

## LEARNING OBJECTIVES

Once you have completed this experiment, you should understand

1. The chemical basis for the spore and capsule stains.
2. The procedure for differentiation between the bacterial spore and vegetative cell forms.
3. The procedure to distinguish capsular material from the bacterial cell.

## PART A Spore Stain (Schaeffer-Fulton Method)

### Principle

Members of the anaerobic genera *Clostridium* and *Desulfotomaculum* and the aerobic genus *Bacillus* are examples of organisms that have the capacity to exist either as metabolically active **vegetative cells** or as highly resistant, metabolically inactive cell types called **spores**. When environmental conditions become unfavorable for continuing vegetative cellular activities, particularly with the exhaustion of a nutritional carbon source, these cells have the capacity to undergo **sporogenesis** and give rise to a new intracellular structure called the **endospore**, which is surrounded by impervious layers called spore coats. As conditions continue to worsen, the endospore is released from the degenerating vegetative cell and becomes an independent cell called a **free spore**. Because of the chemical composition of spore layers,

the spore is resistant to the damaging effects of excessive heat, freezing, radiation, desiccation, and chemical agents, as well as to the commonly employed microbiological stains. With the return of favorable environmental conditions, the free spore may revert to a metabolically active and less resistant vegetative cell through **germination** (see **Figure 13.1**). It should be emphasized that sporogenesis and germination are not means of reproduction but merely mechanisms that ensure cell survival under all environmental conditions.

In practice, the spore stain uses two different reagents.

### Primary Stain


**Malachite Green** Unlike most vegetative cell types that stain by common procedures, the free spore, because of its impervious coats, will not accept the primary stain easily. For further penetration, the application of heat is required. After the primary stain is applied and the smear is heated, both the vegetative cell and spore will appear green.

### Decolorizing Agent

**Water** Once the spore accepts the malachite green, it cannot be decolorized by tap water, which removes only the excess primary stain. The spore remains green. On the other hand, the stain does not demonstrate a strong affinity for vegetative cell components; the water removes it, and these cells will be colorless.

### Counterstain

**Safranin** This contrasting red stain is used as the second reagent to color the decolorized vegetative cells, which will absorb the counterstain and

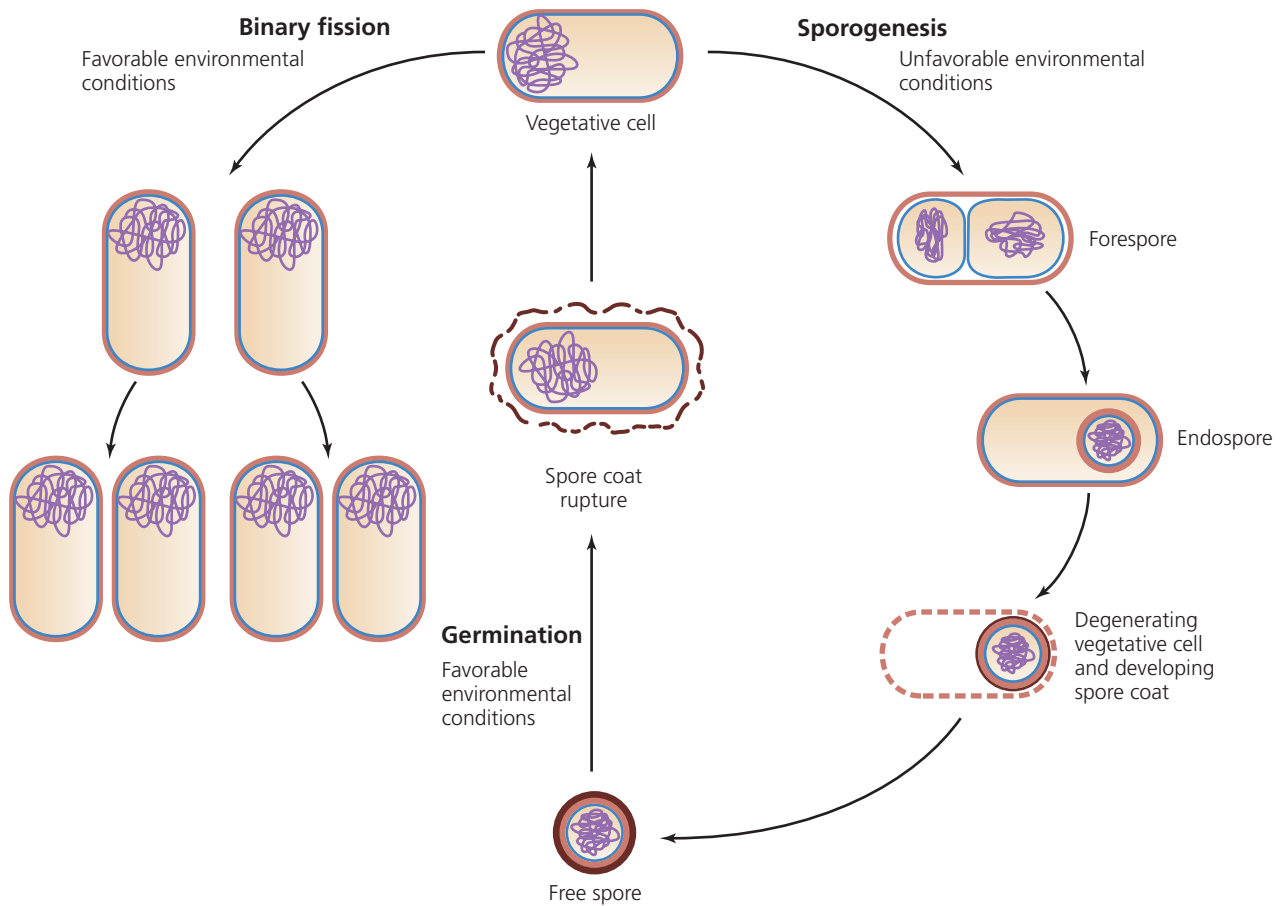


**Figure 13.1** Life cycle of a spore-forming bacterium

appear red. The spores retain the green of the primary stain. A micrograph of spore-stained cells appears in **Figure 13.2**.

## CLINICAL APPLICATION

### Identification of Dangerous Spore-Forming Bacteria

Some spore-forming bacteria can have extremely negative health effects, such as *Bacillus anthracis*, which causes anthrax, and certain *Clostridia* bacteria, which are the causative agents for tetanus, gas gangrene, food poisoning, and pseudomembranous colitis. Differential stains can stain endospores inside bacterial cells as well as free spores to identify these pathogenic bacteria.

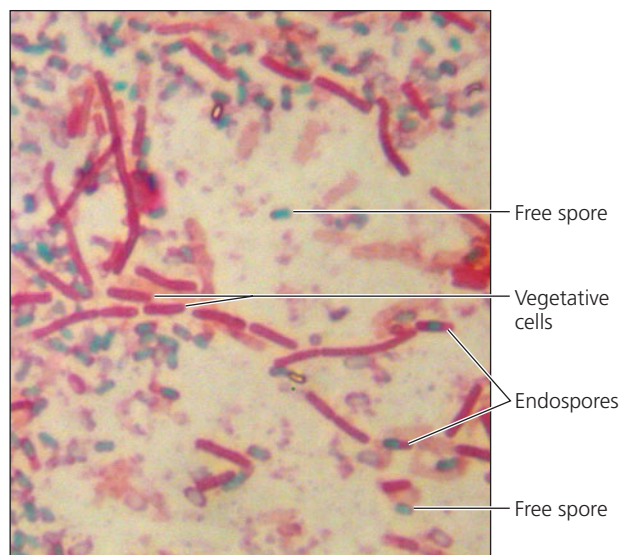

**Figure 13.2** Spore stain showing free spores and vegetative bacilli



## Materials

### Cultures

48- to 72-hour nutrient agar slant culture of *Bacillus cereus* and thioglycollate culture of *Clostridium sporogenes*.

### Reagents

Malachite green and safranin.

### Equipment

Bunsen burner, hot plate, staining tray, inoculating loop, glass slides, bibulous paper, lens paper, and microscope.

## Procedure

### Smear Preparation

1. Obtain two clean glass slides.
2. Make individual smears in the usual manner using aseptic technique.
3. Allow smear to air-dry, and heat fix in the usual manner.

### Spore Staining

Steps 1–5 are illustrated in **Figure 13.3**.

1. Flood smears with malachite green and place on top of a beaker of water sitting on a warm hot plate, allowing the preparation to steam for 2 to 3 minutes. *Note: Do not allow stain to evaporate; replenish stain as needed.* Prevent the stain from boiling by adjusting the hot plate temperature.
2. Remove slides from hot plate, cool, and wash under running tap water.
3. Counterstain with safranin for 30 seconds.
4. Wash with tap water.
5. Blot dry with bibulous paper and examine under oil immersion.
6. In the chart provided in the Lab Report, complete the following:
  - a. Draw a representative microscopic field of each preparation.

- b. Describe the location of the endospore within the vegetative cell as central, subterminal, or terminal on each preparation.
- c. Indicate the color of the spore and vegetative cell on each preparation.

## PART B Capsule Stain (Anthony Method)

### Principle

A **capsule** is a gelatinous outer layer that is secreted by the cell and that surrounds and adheres to the cell wall. It is not common to all organisms. Cells that have a heavy capsule are generally virulent and capable of producing disease, since the structure protects bacteria against the normal phagocytic activities of host cells. Chemically, the capsular material is composed mainly of complex polysaccharides such as levans, dextrans, and celluloses.

**Capsule staining** is more difficult than other types of differential staining procedures because the capsular materials are water-soluble and may be dislodged and removed with vigorous washing. Smears should not be heated because the resultant cell shrinkage may create a clear zone around the organism that is an artifact that can be mistaken for the capsule.

The capsule stain uses two reagents.

### Primary Stain


**Crystal Violet (1% aqueous)** A violet stain is applied to a non-heat-fixed smear. At this point, the cell and the capsular material will take on the dark color.

### Decolorizing Agent


**Copper Sulfate (20%)** Because the capsule is nonionic, unlike the bacterial cell, the primary stain adheres to the capsule but does not bind to it. In the capsule staining method, copper sulfate is used as a decolorizing agent rather than water. The copper sulfate washes the purple primary stain out of the capsular material without removing the stain bound to the cell wall. At the same time, the decolorized capsule absorbs the copper sulfate, and the capsule will now appear blue in contrast to the deep purple color of the cell.

**Figure 13.4** shows the presence of a capsule as a clear zone surrounding the darker-stained cell.

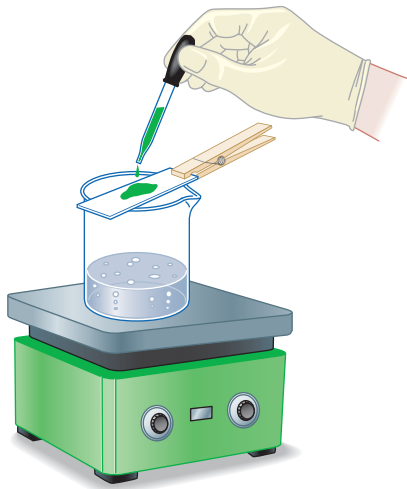
## PROCEDURE



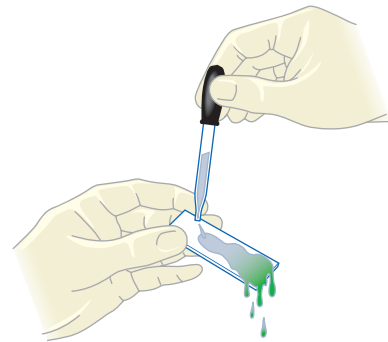
1 Flood smears with malachite green and steam over a beaker of water placed on a hot plate.



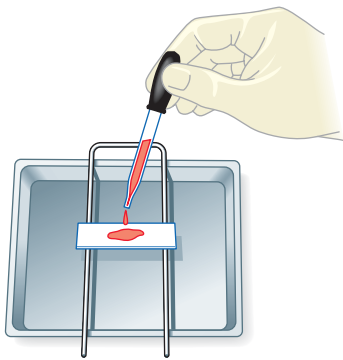
2 Cool and wash off stain with tap water. The water also serves as the decolorizing agent.



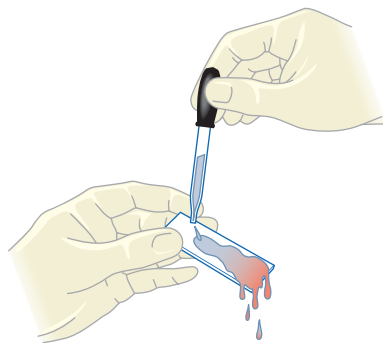3 Counterstain with safranin for 30 seconds.




4 Wash off the safranin with tap water.

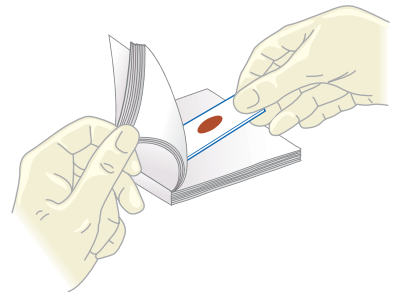

5 Blot the slide dry with bibulous paper.

**Figure 13.3 Spore-staining procedure**

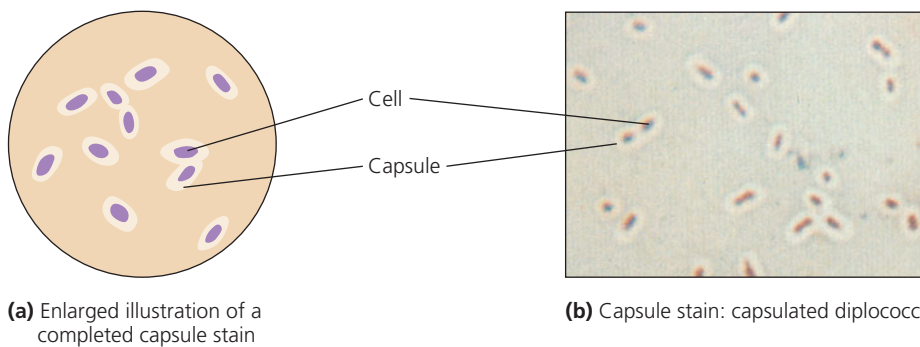

**Figure 13.4 Capsule stain**

## CLINICAL APPLICATION

### Encapsulated Bacterial Pneumonia

The virulence of an organism is increased by the presence of a capsule, since the capsule protects the organism from phagocytosis by white blood cells and inhibits antibody or complement fixation. The water-soluble polysaccharide and/or the polypeptide composition of the bacterial capsule makes staining this feature difficult. Gram-negative bacteria that form capsules include *Haemophilus influenzae* and *Klebsiella pneumoniae*. Gram-positive bacteria that form capsules include *Bacillus anthracis* and *Streptococcus pneumoniae*. If a bacterial infection is not being cleared or responding to antibiotic therapy as expected, staining of isolated organisms to determine the presence of a capsule may be warranted.

## AT THE BENCH



## Materials

### Cultures

48-hour-old skimmed milk cultures of *Alcaligenes viscolactis*, *Leuconostoc mesenteroides*, and *Enterobacter aerogenes*.

### Reagents

1% crystal violet and 20% copper sulfate ( $\text{CuSO}_4 \cdot 5\text{H}_2\text{O}$ ).

## Equipment

Bunsen burner, inoculating loop or needle, staining tray, bibulous paper, lens paper, glass slides, and microscope.

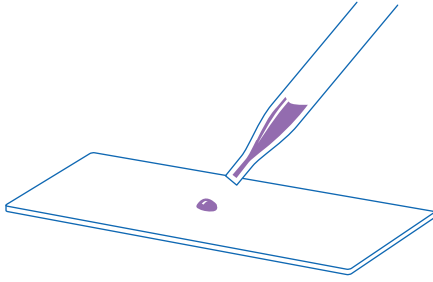
## Procedure

Steps 1–5 are pictured in **Figure 13.5**.


1. Obtain one clean glass slide. Place several drops of crystal violet stain on the slide.
2. Using aseptic technique, add three loopfuls of a culture to the stain and *gently* mix with the inoculating loop.
3. With a clean glass slide, spread the mixture over the entire surface of the slide to create a very thin smear. Let stand for 5 to 7 minutes. Allow smears to air-dry. *Note: Do not heat fix.*
4. Wash smears with 20% copper sulfate solution.
5. *Gently* blot dry and examine under oil immersion.
6. Repeat Steps 1–5 for each of the remaining test cultures.
7. In the chart provided in the Lab Report, complete the following:
  - a. Draw a representative microscopic field of each preparation.
  - b. Record the comparative size of the capsule; that is, small, moderate, or large.
  - c. Indicate the color of the capsule and of the cell on each preparation.




## PROCEDURE



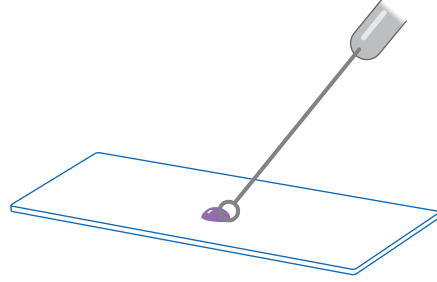
1 Place several drops of crystal violet stain on a clean glass slide.



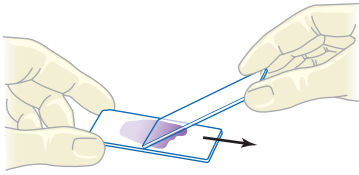
2 Aseptically transfer 3 loopfuls of culture to the stain and gently mix with the loop.



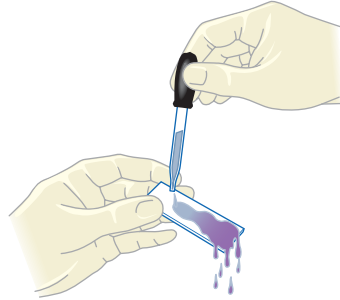3 With a clean glass slide, spread mixture to form a thin smear. Air-dry.



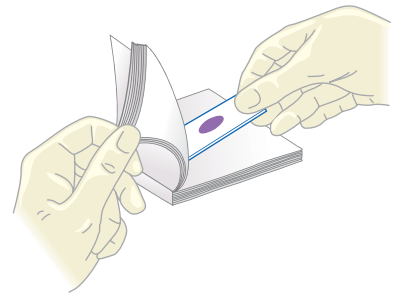4 Wash smear with 20% copper sulfate solution.



5 Gently blot dry with bibulous paper.



**Figure 13.5** Capsule staining procedure

Name: \_\_\_\_\_




Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results
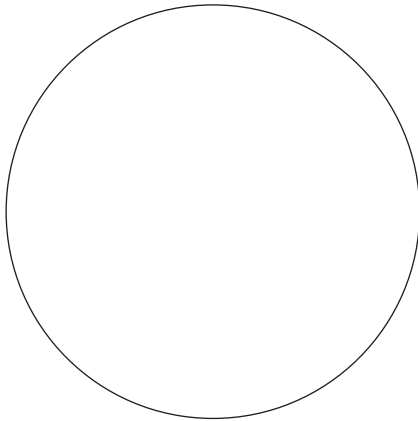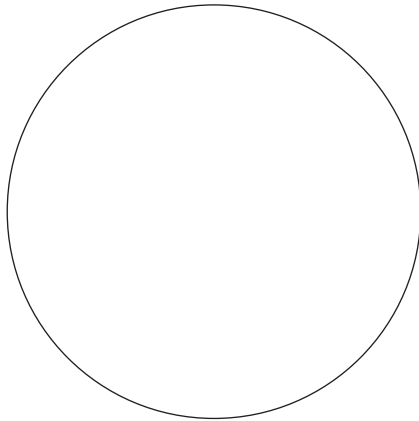
### PART A: Spore Stain

	<i>C. sporogenes</i>	<i>B. cereus</i>
Draw a representative field.		
Color of spores	_____	_____
Color of vegetative cells	_____	_____
Location of endospore	_____	_____

### PART B: Capsule Stain

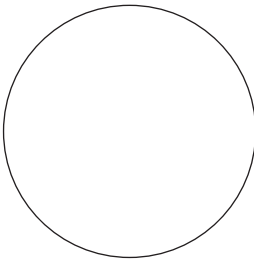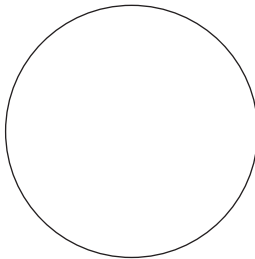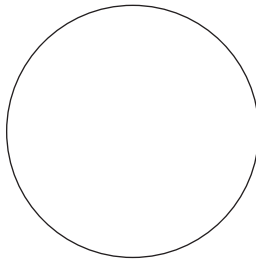	<i>A. viscolactis</i>	<i>L. mesenteroides</i>	<i>E. aerogenes</i>
Draw a representative field.			
Capsule size	_____	_____	_____
Color of capsule	_____	_____	_____
Color of cell	_____	_____	_____



c. You did not apply heat during the application of the primary stain.

4. Explain the medical significance of a capsule.

5. Explain the function of copper sulfate in this procedure.

*This page intentionally left blank*

# Cultivation of Microorganisms: Nutritional and Physical Requirements, and Enumeration of Microbial Populations

## LEARNING OBJECTIVES

Once you have completed the experiments in this section, you should be familiar with

1. The basic nutritional and environmental requirements for the cellular activities of all forms of life.
2. The principles associated with the use of routine and special-purpose media for microbial cultivation.
3. The diversified physical factors essential for microbial cultivation.
4. Specialized techniques for the cultivation of anaerobic microorganisms.
5. The serial dilution–agar plate technique for enumeration of viable microorganisms.
6. The growth dynamics of bacterial populations.

## Introduction

As do all other living organisms, microorganisms require certain basic nutrients and physical factors for the sustenance of life. However, their particular requirements vary greatly. Understanding these needs is necessary for successful cultivation of microorganisms in the laboratory.

## Nutritional Needs

Nutritional needs of microbial cells are supplied in the laboratory through a variety of media. The following list illustrates the nutritional diversity that exists among microbes.

1. **Carbon:** This is the most essential and central atom common to all cellular structures and functions. Among microbial cells, two carbon-dependent types are noted:
  - a. **Autotrophs:** These organisms can be cultivated in a medium consisting solely of inorganic compounds; specifically, they use inorganic carbon in the form of carbon dioxide.
  - b. **Heterotrophs:** These organisms cannot be cultivated in a medium consisting solely of inorganic compounds; they must be supplied with organic nutrients, primarily glucose.
2. **Nitrogen:** This is also an essential atom in many cellular macromolecules, particularly proteins and nucleic acids. Proteins serve as the structural molecules forming the so-called fabric of the cell and as functional molecules, enzymes, that are responsible for the metabolic activities of the cell. Nucleic acids include DNA, the genetic basis of cell life, and RNA, which plays an active role in protein synthesis within the cell. Some microbes use atmospheric nitrogen, others rely on inorganic compounds such as ammonium or nitrate salts, and still others require nitrogen-containing organic compounds such as amino acids.
3. **Nonmetallic elements:** Two major nonmetallic ions are used for cellular nutrition:
  - a. **Sulfur** is integral to some amino acids and is therefore a component of proteins.

Sources include organic compounds such as sulfur-containing amino acids, inorganic compounds such as sulfates, and elementary sulfur.

- b. Phosphorus** is necessary for the formation of the nucleic acids DNA and RNA and also for synthesis of the high-energy organic compound adenosine triphosphate (ATP). Phosphorus is supplied in the form of phosphate salts for use by all microbial cells.
- 4. Metallic elements:**  $\text{Ca}^{2+}$ ,  $\text{Zn}^{2+}$ ,  $\text{Na}^+$ ,  $\text{K}^+$ ,  $\text{Cu}^{2+}$ ,  $\text{Mn}^{2+}$ ,  $\text{Mg}^{2+}$ ,  $\text{Fe}^{2+}$ , and  $\text{Fe}^{3+}$  are some of the metallic ions necessary for continued efficient performance of varied cellular activities. Some of these activities are osmoregulation, regulation of enzyme activity, and electron transport during biooxidation. Remember that these ions are micronutrients and are required in trace concentrations only. Inorganic salts supply these materials.
  - 5. Vitamins:** These organic substances contribute to cellular growth and are essential in minute concentrations for cell activities. They are also sources of coenzymes, which are required for the formation of active enzyme systems. Some microbes require vitamins to be supplied in a preformed state for normal metabolic activities. Some possess extensive vitamin-synthesizing pathways, whereas others can synthesize only a limited number from other compounds present in the medium.
  - 6. Water:** All cells require distilled water in the medium so that the low-molecular-weight nutrients can cross the cell membrane.
  - 7. Energy:** Active transport, biosynthesis, and biodegradation of macromolecules are the metabolic activities of cellular life. These activities can be sustained only if there is a constant availability of energy within the cell. Two bioenergetic types of microorganisms exist:
    - a. Phototrophs** use radiant energy as their sole energy source.
    - b. Chemotrophs** depend on oxidation of chemical compounds as their energy source. Some microbes use organic molecules such as glucose; others utilize inorganic compounds such as  $\text{H}_2\text{S}$  or  $\text{NaNO}_2$ .

## Physical Factors

Three of the most important physical factors that influence the growth and survival of cells are temperature, pH, and the gaseous environment. An understanding of the roles they play in cell metabolism is essential.

- 1. Temperature** influences the rate of chemical reactions through its action on cellular enzymes. Bacteria, as a group of organisms, exist over a wide range of temperatures. However, individual species can exist only within a narrower spectrum of temperatures. Low temperatures slow down or inhibit enzyme activity, thereby slowing down or inhibiting cell metabolism and, consequently, cell growth. High temperatures cause coagulation and thus irreversibly denature thermolabile enzymes. Although enzymes differ in their degree of heat sensitivity, generally temperatures in the range of  $70^\circ\text{C}$  destroy most essential enzymes and cause cell death.
- 2. The pH of the extracellular environment** greatly affects cells' enzymatic activities. Most commonly, the optimum pH for cell metabolism is in the neutral range of 7. An increase in the hydrogen ion concentration resulting in an acidic pH (below 7) or a decrease in the hydrogen ion concentration resulting in an alkaline pH (above 7) is often detrimental. Either increase or decrease will slow down the rate of chemical reactions because of the destruction of cellular enzymes, thereby affecting the rate of growth and, ultimately, survival.
- 3. The gaseous requirement** in most cells is atmospheric oxygen, which is necessary for the biooxidative process of respiration. Atmospheric oxygen plays a vital role in ATP formation and the availability of energy in a utilizable form for cell activities. Other cell types, however, lack the enzyme systems for respiration in the presence of oxygen and therefore must use an anaerobic form of respiration or fermentation.

The following exercises will demonstrate the diversity of nutritional and environmental requirements among microorganisms.

## LEARNING OBJECTIVES

Once you have completed this experiment, you should know how to evaluate

1. The abilities of several types of media to support the growth of different bacterial species.
2. The nutritional needs of the bacteria under study.

## Principle

To satisfy the diverse nutritional needs of bacteria, bacteriologists employ two major categories of media for routine cultivation.

## Chemically Defined Media

These are composed of known quantities of chemically pure, specific organic and/or inorganic compounds. Their use requires knowledge of the organism's specific nutritional needs. The following two chemically defined media are used in this exercise:

1. **Inorganic synthetic broth:** This completely inorganic medium is prepared by incorporating the following salts per 1000 ml of water:

Sodium chloride (NaCl)	5.0 g
Magnesium sulfate (MgSO <sub>4</sub> )	0.2 g
Ammonium dihydrogen phosphate (NH <sub>4</sub> H <sub>2</sub> PO <sub>4</sub> )	1.0 g
Dipotassium hydrogen phosphate (K <sub>2</sub> HPO <sub>4</sub> )	1.0 g
Atmospheric carbon dioxide (CO <sub>2</sub> )	

2. **Glucose salts broth:** This medium is composed of salts incorporated into the inorganic synthetic broth medium plus **glucose**, 5 g per liter, which serves as the sole organic carbon source.

## Complex Media

The exact chemical composition of these media is not known. They are made of extracts of plant and animal tissue and are variable in their chemical composition. Most contain abundant amino acids, sugars, vitamins, and minerals; however, the quantities of these constituents are not known. They are capable of supporting the growth of most heterotrophs. The following two complex media are used in this exercise.

1. **Nutrient broth:** This basic complex medium is prepared by incorporating the following ingredients per 1000 ml of distilled water:

Peptone	5.0 g
Beef extract	3.0 g

**Peptone**, a semidigested protein, is primarily a nitrogen source. The **beef extract**, a beef derivative, is a source of organic carbon, nitrogen, vitamins, and inorganic salts.


2. **Yeast extract broth:** This is composed of the basic artificial medium ingredients used in the nutrient broth plus **yeast extract**, 5 g per liter, which is a rich source of vitamin B and provides additional organic nitrogen and carbon compounds.

The yeast extract broth is an example of an **enriched medium** and is used for the cultivation of **fastidious** microorganisms, organisms that have highly elaborate and specific nutritional needs. These bacteria do not grow or grow poorly on a basic artificial medium and require the addition of one or more growth-supporting substances, enrichments such as additional plant or animal extracts, vitamins, or blood.

## Measuring Turbidity

In this experiment you will evaluate (1) the abilities of media to support the growth of different species of bacteria, and (2) the nutritional needs of the bacteria. You will observe the amount of growth, measured by turbidity, present in each culture following incubation. To evaluate more accurately the amount



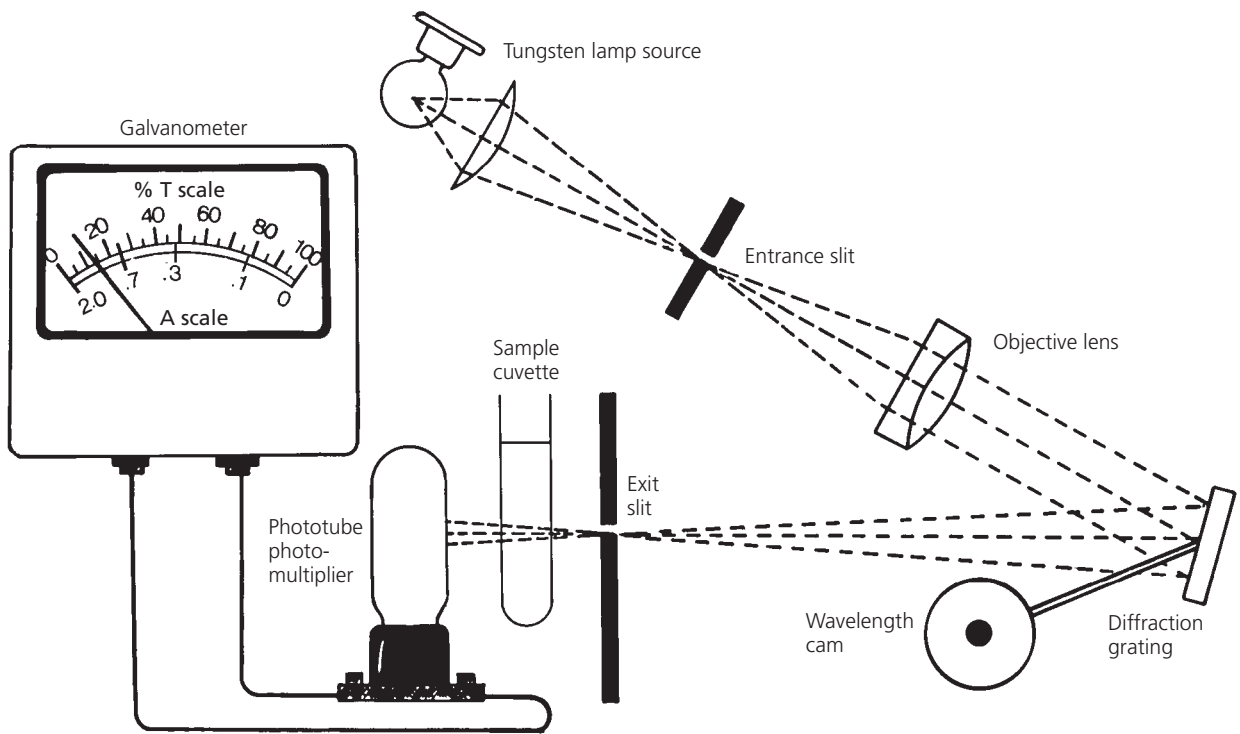

**Figure 14.1** Schematic diagram of a spectrophotometer

of growth, a Bausch & Lomb Spectronic 20 spectrophotometer will be used.

This instrument measures the amount of light transmitted (T) or absorbed (A). It transmits a beam of light at a single wavelength (monochromatic light) through a liquid culture. The cells suspended in the culture interrupt the passage of light, and the amount of light energy transmitted through the suspension is measured on a photoelectric cell and converted into electrical energy. The electrical energy is then recorded on a galvanometer as 0% to 100% T. A schematic representation of a spectrophotometer is shown in **Figure 14.1**.

In practice, the density of a cell suspension is expressed as absorbance (A) rather than percent T, since A is directly proportional to the concentration of cells, whereas percent T is inversely proportional to the concentration of suspended cells. Therefore, as the turbidity of a culture increases, the A increases and percent T decreases, indicating growth of the cell population in the culture. For example, in comparing three cultures with A readings of 0.10 (percent T = 78), 0.30 (percent T = 49), and 0.50 (percent T = 30), the A reading of 0.50 would be indicative of the most abundant growth, and the 0.10 reading would be indicative of the least amount of growth.

**Figure 14.2** shows the Bausch & Lomb Spectronic 20 spectrophotometer.

## CLINICAL APPLICATION

### The Purpose of Specialized Media

The successful cultivation of bacteria requires the use of culture media containing the nutritional and biochemical requirements capable of supporting growth. There is no single medium that can support the growth of all microorganisms. This challenge has been met by the development of a variety of specialized media. For example, the streptococci require media supplemented with blood in order to determine certain properties that are necessary for isolation and species identification. Another example is the thioglycollate medium, which contains thioglycolic acid that removes oxygen from the medium to encourage the growth of certain anaerobic bacteria.

## AT THE BENCH



## Materials Cultures

Saline suspension of 24-hour Trypticase soy broth cultures, adjusted to 0.05 absorbance at a wavelength of 600 nm, of *Escherichia coli*, *Alcaligenes faecalis*, and *Streptococcus mitis*.



**Figure 14.2** The Bausch & Lomb Spectronic 20 spectrophotometer

## Media

Per designated student group: three test tubes (13 × 100 mm) of each type of broth: inorganic synthetic broth, glucose salts broth, nutrient broth, and yeast extract broth.

## Equipment

Bunsen burner, sterile 1-ml serological pipettes, mechanical pipetting device, glassware marking pencil, test tube rack, and Bausch & Lomb Spectronic 20 spectrophotometer.

## Procedure Lab One

1. Using a sterile 1-ml pipette, add 0.1 ml of the *E. coli* culture to one test tube of each of the appropriately labeled media.
2. Repeat Step 1 for inoculation with *A. faecalis* and *S. mitis*.
3. Incubate the test cultures for 24 to 48 hours at 37°C.

## Procedure Lab Two

Follow the instructions below and refer to Figure 14.2 for the use of the spectrophotometer to obtain the absorbance readings of all your cultures.


1. Turn the instrument on 10 to 15 minutes prior to use.
2. Set wavelength at 600 nm.
3. Set percent transmittance to 0% (A to 2) by turning the knob on the left.

4. Read the four yeast extract broth cultures as follows:
  - a. Wipe clean the provided test tube of sterile yeast broth that will serve as the blank for the yeast broth culture readings. Fingerprints on the test tube will obscure the light path of the spectrophotometer.
  - b. Insert the yeast extract broth blank into the tube holder, close the cover, and set the A to 0 (percent T = 100) by turning the knob on the right.
  - c. Shake lightly or tap one of the tubes of yeast extract broth culture to resuspend the bacteria, wipe the test tube clean, and allow it to sit for several seconds for the equilibration of the bacterial suspension.
  - d. Remove the yeast extract broth blank from the tube holder.
  - e. Insert a yeast extract broth culture into the tube holder, close the cover, and read and record the optical density reading in the chart provided in the Lab Report.
  - f. Remove the yeast extract broth culture from the tube holder.
  - g. Reset the spectrophotometer to an A of 2 with the tube holder empty and to an A of 0 with the yeast extract broth blank.
  - h. Repeat Steps c through g to read and record the absorbance of the remaining yeast extract broth cultures.
5. Repeat Step 4 (a–h) to read and record the absorbance of the nutrient broth cultures. Use the provided nutrient broth blank to set the spectrophotometer to an A of 0.

6. Repeat Step 4 (a–h) to read and record the absorbance of the glucose salts broth cultures. Use the provided glucose salts broth blank to set the spectrophotometer to an A of 0.
7. Repeat Step 4 (a–h) to read and record the absorbance of the inorganic synthetic broth cultures. Use the provided inorganic synthetic broth blank to set the spectrophotometer to an A of 0.
8. At the end of the experiment, return all cultures to the area designated for their disposal.
9. Complete the Lab Report.



## Review Questions

1. Explain the advantages of using A readings rather than percent T as a means of estimating microbial growth.
2. Explain the reason for the use of different medium blanks in adjusting the spectrophotometer prior to obtaining A readings.
3. Why are complex media preferable to chemically defined media for routine cultivation of microorganisms?
4. Would you expect a heterotrophic organism to grow in an inorganic synthetic medium? Explain.
5.  A soil isolate is found to grow poorly in a basic artificial medium. You suspect that a vitamin supplement is required.
  - a. What supplement would you use to enrich the medium to support and maintain the growth of the organisms? Explain.
  - b. Outline the procedure you would follow to determine the specific vitamins required by the organism to produce a more abundant growth.

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be familiar with

1. The use and function of specialized media for the selection and differentiation of microorganisms.
2. How an enriched medium like blood agar can also function as both a selective and differential medium.

## Principle

Numerous special-purpose media are available for functions including the following:

1. Isolation of bacterial types from a mixed population of organisms.
2. Differentiation among closely related groups of bacteria on the basis of macroscopic appearance of the colonies and biochemical reactions within the medium.
3. Enumeration of bacteria in sanitary microbiology, such as in water and sewage, and also in food and dairy products.
4. Assay of naturally occurring substances such as antibiotics, vitamins, and products of industrial fermentation.
5. Characterization and identification of bacteria by their abilities to produce chemical changes in different media.

In addition to nutrients necessary for the growth of all bacteria, special-purpose media contain both nutrients and chemical compounds important for specific metabolic pathways in different types of bacteria. In this exercise, three types of media will be studied and evaluated.

## Selective Media

These media are used to select (isolate) specific groups of bacteria. They incorporate chemical substances that inhibit the growth of one type of

bacteria while permitting growth of another, thus facilitating bacterial isolation.

1. **Phenylethyl alcohol agar:** This medium is used for the isolation of most gram-positive organisms. The phenylethyl alcohol is partially inhibitory to gram-negative organisms, which may form visible colonies whose size and number are much smaller than on other media.
2. **Crystal violet agar:** This medium is selective for most gram-negative microorganisms. Crystal violet dye exerts an inhibitory effect on most gram-positive organisms.
3. **7.5% sodium chloride agar:** This medium is inhibitory to most organisms other than halophilic (salt-loving) microorganisms. It is most useful in the detection of members of the genus *Staphylococcus*.


**Figure 15.1** is a photo illustrating the selective effect of phenylethyl alcohol agar, which inhibits the gram-negative organism *E. coli* and selects for the gram-positive organism *S. aureus*.

## Differential/Selective Media

These media can distinguish among morphologically and biochemically related groups of organisms. They incorporate chemical compounds that, following inoculation and incubation, produce a characteristic change in the appearance of bacterial growth and/or the medium surrounding the colonies, which permits differentiation.

Sometimes differential and selective characteristics are combined in a single medium. MacConkey agar is a good example of this because it contains bile salts and crystal violet, which inhibit gram-positive organisms and allow gram-negatives to grow. In addition, it contains the substrate lactose and the pH indicator neutral red, which differentiates the red lactose-fermenting colonies from the translucent nonfermenting colonies. The following media are examples of this type of media:

1. **Mannitol salt agar:** This medium contains a high salt concentration, 7.5% NaCl, which is inhibitory to the growth of most, but not all,

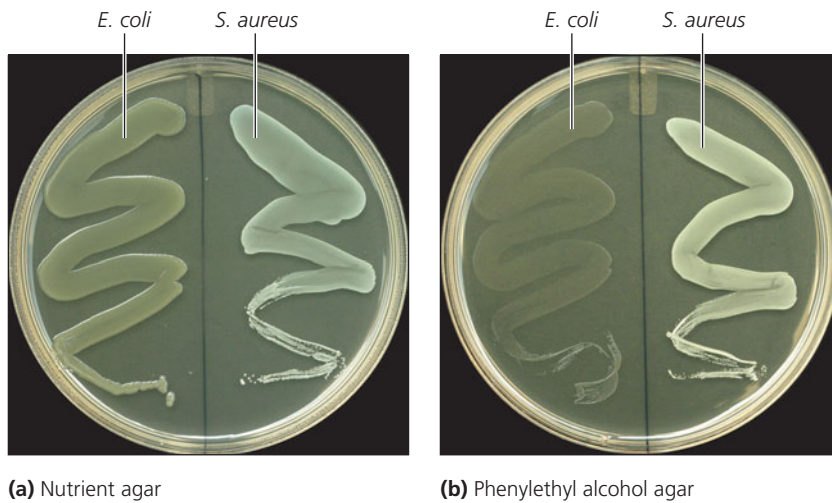

**Figure 15.1** Selective effect of phenylethyl alcohol agar reduces the growth of *E. coli* and selects for *S. aureus*

bacteria other than the staphylococci. The medium also performs a differential function: It contains the carbohydrate mannitol, which some staphylococci are capable of fermenting, and phenol red, a pH indicator for detecting acid produced by mannitol-fermenting staphylococci. These staphylococci exhibit a yellow zone surrounding their growth; staphylococci that do not ferment mannitol will not produce a change in coloration.

2. **MacConkey agar:** The inhibitory action of crystal violet on the growth of gram-positive organisms allows the isolation of gram-negative bacteria. Incorporation of the carbohydrate lactose, bile salts, and the pH indicator neutral red permits differentiation of enteric bacteria on the basis of their ability to ferment lactose. On this basis, enteric bacteria are separated into two groups:

- a. **Coliform bacilli** produce acid as a result of lactose fermentation. The bacteria exhibit a red coloration on their surface. *Escherichia coli* produce greater quantities of acid from lactose than other coliform species. When this occurs, the medium surrounding the growth also becomes pink because of the action of the acid that precipitates the bile salts, followed by absorption of the neutral red.

- b. **Dysentery, typhoid, and paratyphoid bacilli** are not lactose fermenters and therefore do not produce acid. The colonies appear tan and frequently transparent.


3. **Eosin–methylene blue agar (Levine):** Lactose and the dyes eosin and methylene blue permit differentiation between enteric lactose fermenters and nonfermenters as well as identification of the colon bacillus, *E. coli*. The *E. coli* colonies are blue-black with a metallic green sheen caused by the large quantity of acid that is produced and that precipitates the dyes onto the growth's surface. Other coliform bacteria, such as *Enterobacter aerogenes*, produce thick, mucoid, pink colonies on this medium. Enteric bacteria that do not ferment lactose produce colorless colonies, which, because of their transparency, appear to take on the purple color of the medium. This medium is also partially inhibitory to the growth of gram-positive organisms, and thus gram-negative growth is more abundant.

A photographic representation of the effects of selective/differential media is presented in **Figure 15.2**.

## Enriched Media

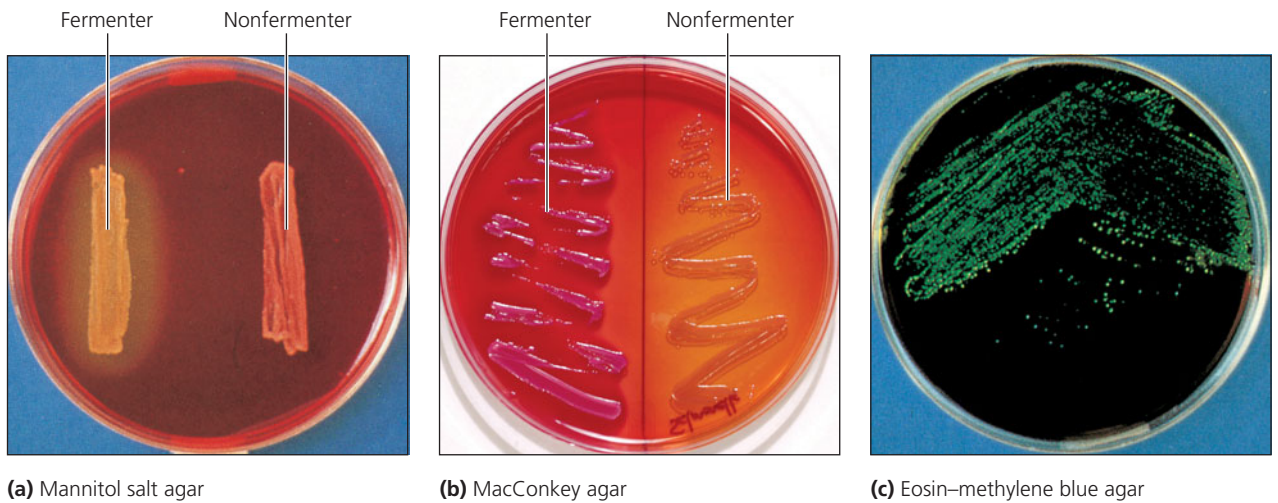Enriched media are media that have been supplemented with highly nutritious materials, such as blood, serum, or yeast extract, for the purpose of cultivating fastidious organisms.

For example, in **blood agar**, the blood incorporated into the medium is an enrichment ingredient for the cultivation of fastidious organisms such as the *Streptococcus* spp. The blood also permits demonstration of the hemolytic properties of some microorganisms, particularly the streptococci,



(a) Mannitol salt agar (b) MacConkey agar (c) Eosin–methylene blue agar


**Figure 15.2** Effects of selective/differential media

whose hemolytic activities are classified as follows:

1. **Gamma hemolysis:** No lysis of red blood cells results in no significant change in the appearance of the medium surrounding the colonies.
2. **Alpha hemolysis:** Incomplete lysis of red blood cells, with reduction of hemoglobin to methemoglobin, results in a greenish halo around the bacterial growth.
3. **Beta hemolysis:** Lysis of red blood cells with complete destruction and use of hemoglobin by the organism results in a clear zone surrounding the colonies. This hemolysis is produced by two types of beta hemolysins,

namely **streptolysin O**, an antigenic, oxygen-labile enzyme, and **streptolysin S**, a nonantigenic, oxygen-stable lysin. The hemolytic reaction is enhanced when blood agar plates are streaked and simultaneously stabbed to show subsurface hemolysis by streptolysin O in an environment with reduced oxygen tension. Based on the hemolytic patterns on blood agar, the pathogenic beta-hemolytic streptococci may be differentiated from other streptococci.

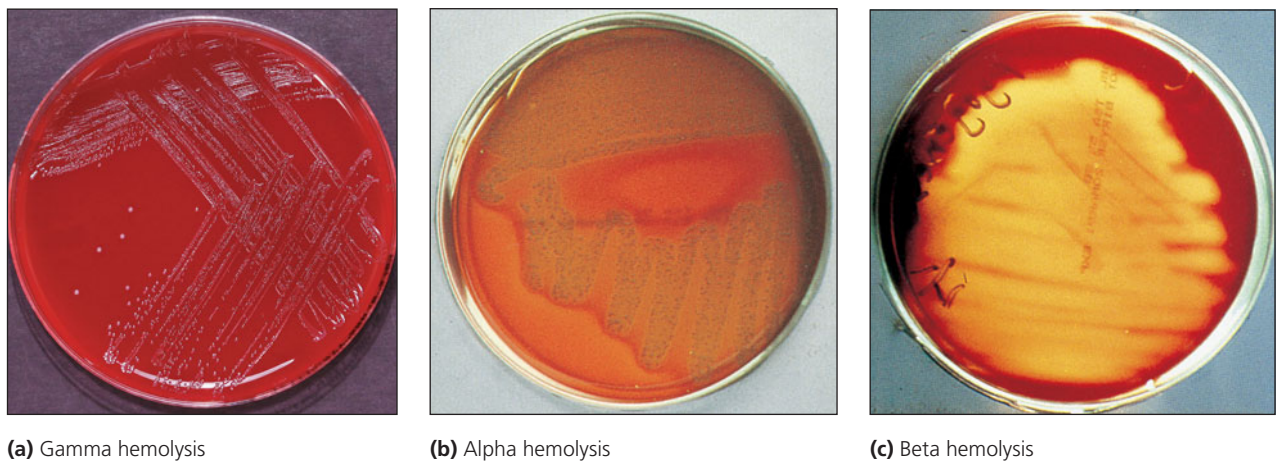**Figure 15.3** shows the different types of hemolysis exhibited by different species of the genus *Streptococcus* on blood agar.



(a) Gamma hemolysis (b) Alpha hemolysis (c) Beta hemolysis

**Figure 15.3** Types of hemolysis exhibited on a blood agar plate



## CLINICAL APPLICATION

### First Steps in Infected Wound Diagnosis

Wounds that have become infected may be swabbed or surgically processed to remove tissue. Once stained samples have revealed infectious agents, cultures are typically made on (1) blood agar for isolation of staphylococci and streptococci bacteria, (2) MacConkey agar for gram-negative rods, and (3) enriched media that can support aerobes or anaerobes, such as thioglycollate broth. Additional media may be used, depending on what was observed microscopically, including Sabouraud dextrose agar for fungi and Löwenstein-Jensen medium for acid-fast rods. Once the microbes are isolated, further tests (which you will learn soon!) would likely be needed for complete identification.

## AT THE BENCH



## Materials

### Cultures

24- to 48-hour Trypticase soy broth cultures of *Enterobacter aerogenes*, *Escherichia coli*, *Streptococcus* var. Lancefield Group E, *Streptococcus mitis*, *Enterococcus faecalis*, *Staphylococcus aureus*, *Staphylococcus epidermidis*, and *Salmonella typhimurium*.

## Media

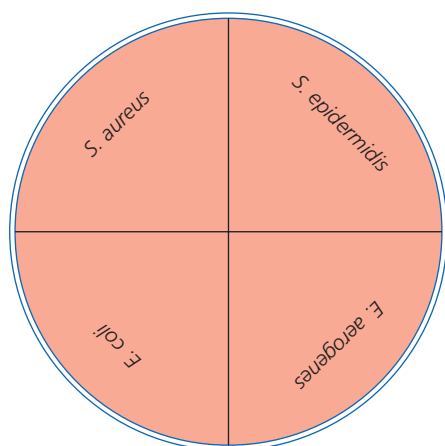
Per designated student group: one each of phenylethyl alcohol agar, crystal violet agar, 7.5% sodium chloride agar, mannitol salt agar, MacConkey agar, eosin–methylene blue agar, and blood agar.

## Equipment


Bunsen burner, inoculating loop, and glassware marking pencil.

## Procedure Lab One

- Using the bacterial organisms listed in Step 2, prepare and inoculate each of the plates in the following manner:
  - Appropriately label the cover of each plate as indicated in the Laboratory Protocol section on page xv.
  - Divide each of the Petri dishes into the required number of sections (one section for each different organism) by marking the *bottom of the dish*. Label each section with the name of the organism to be inoculated, as illustrated in **Figure 15.4a**.
  - Using aseptic technique, inoculate all plates, except the blood agar plate, with the designated organisms by making a single line of inoculation of each organism in its appropriate section (**Figure 15.4b**). Be sure to close the Petri dish and flame the inoculating needle between inoculations of the different organisms.



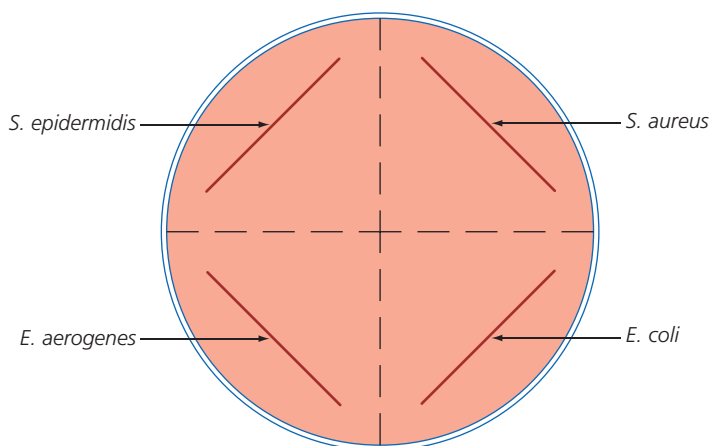(a) Bottom view of plate (Step 1b)



(b) Surface view of plate (Step 1c)

**Figure 15.4** Mannitol salt agar plate preparation and inoculation procedure



*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

Type of Medium	Medium	Bacterial Species	Amount of Growth	Appearance of Growth	Appearance of Medium
Selective	Phenylethyl alcohol agar	<i>E. coli</i>			
		<i>S. aureus</i>			
		<i>E. faecalis</i>			
	Crystal violet agar	<i>E. coli</i>			
		<i>S. aureus</i>			
		<i>E. faecalis</i>			
	7.5% sodium chloride agar	<i>E. coli</i>			
		<i>S. aureus</i>			
		<i>S. epidermidis</i>			
Differential/ Selective	Mannitol salt agar	<i>E. coli</i>			
		<i>E. aerogenes</i>			
		<i>S. aureus</i>			
		<i>S. epidermidis</i>			
	MacConkey agar	<i>E. coli</i>			
		<i>E. aerogenes</i>			
		<i>S. typhimurium</i>			
		<i>S. aureus</i>			
	Eosin–methylene blue agar	<i>E. coli</i>			
		<i>E. aerogenes</i>			
		<i>S. typhimurium</i>			
		<i>S. aureus</i>			
Enriched	Blood agar	<i>S. mitis</i>			
		<i>E. faecalis</i>			
		<i>Streptococcus</i> var. Lancefield Group E			

Indicate the specific selective and/or differential purpose of each of the following media:

a. Phenylethyl alcohol agar:

b. Crystal violet agar:

c. 7.5% sodium chloride agar:

d. Mannitol salt agar:

e. MacConkey agar:




d. High salt concentration in the mannitol salt agar medium:

e. Lactose in the MacConkey agar medium:

f. Phenylethyl alcohol in the phenylethyl alcohol agar medium:

2. Why are crystal violet agar and 7.5% sodium chloride agar considered selective media?

3.  A patient exhibits a boil on his neck. You, as a microbiology technician, are asked to identify the causative organism and determine whether it is pathogenic. Describe the procedure that you would follow to make this determination.

## LEARNING OBJECTIVES

Once you have completed this experiment, you should know

1. The diverse growth temperature requirements of bacteria.
2. How to determine whether the optimum growth temperature is also the ideal temperature for enzyme-regulated cell activities such as pigment production and carbohydrate fermentation.


## Principle

Microbial growth is directly dependent on how temperature affects cellular enzymes. With increasing temperatures, enzyme activity increases until the three-dimensional configuration of these molecules is lost because of denaturation of their protein structure. On the other hand, as the temperature is lowered toward the freezing point,

enzyme inactivation occurs and cellular metabolism gradually diminishes. At 0°C, biochemical reactions cease in most cells.


Bacteria, as a group of living organisms, are capable of growth within an overall temperature range of minus 5°C to 80°C. Each species, however, requires a narrower range that is determined by the heat sensitivity of its enzyme systems. Specific temperature ranges consist of the following **cardinal (significant) temperature points** (Figure 16.1):

1. **Minimum growth temperature:** The lowest temperature at which growth will occur. Below this temperature, enzyme activity is inhibited and the cells are metabolically inactive so that growth is negligible or absent.
2. **Maximum growth temperature:** The highest temperature at which growth will occur. Above this temperature, most cell enzymes are destroyed and the organism dies.
3. **Optimum growth temperature:** The temperature at which the rate of reproduction is most rapid; however, it is not necessarily optimum or ideal for all enzymatic activities of the cell.

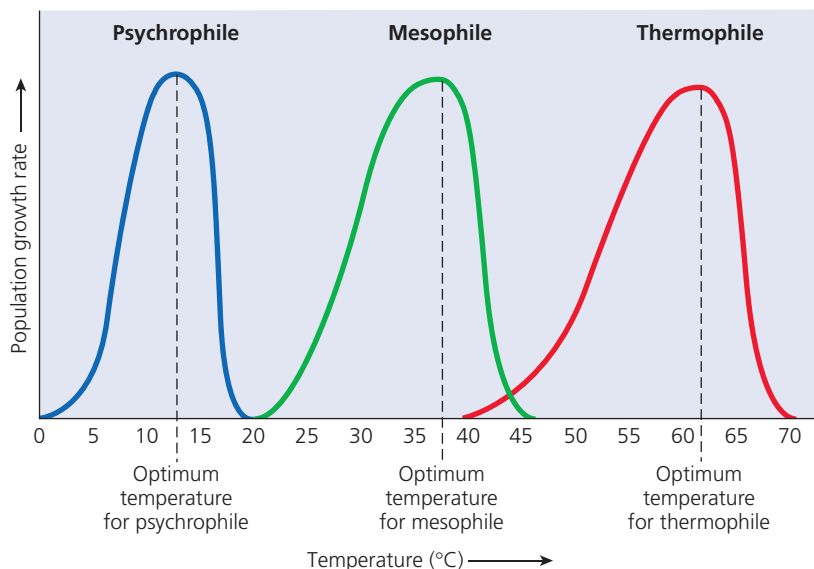

**Figure 16.1** The effect of temperature on the growth of microorganisms



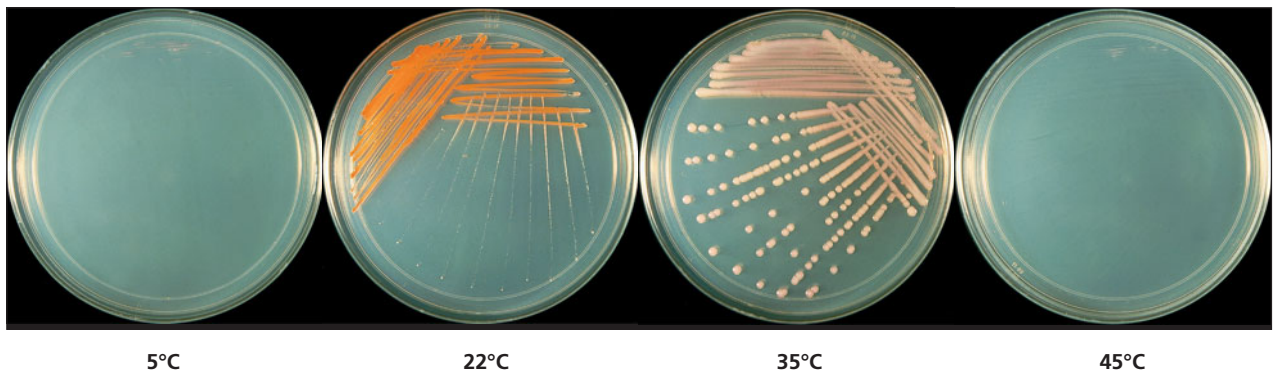

**Figure 16.2** Effect of temperature on bacterial growth and pigmentation

**Figure 16.2** shows the effects of temperature on bacterial growth and pigment production.

All bacteria can be classified into one of three major groups, depending on their temperature requirements:

1. **Psychrophiles:** Bacterial species that will grow within a temperature range of  $-5^{\circ}\text{C}$  to  $20^{\circ}\text{C}$ . The distinguishing characteristic of all psychrophiles is that they will grow between  $0^{\circ}$  and  $5^{\circ}\text{C}$ .
2. **Mesophiles:** Bacterial species that will grow within a temperature range of  $20^{\circ}\text{C}$  to  $45^{\circ}\text{C}$ . The distinguishing characteristics of all mesophiles are their ability to grow at human body temperature ( $37^{\circ}\text{C}$ ) and their inability to grow at temperatures above  $45^{\circ}\text{C}$ . Included among the mesophiles are two distinct groups:
  - a. Those whose optimum growth temperature is in the range of  $20^{\circ}\text{C}$  to  $30^{\circ}\text{C}$  are plant saprophytes.
  - b. Those whose optimum growth temperature is in the range of  $35^{\circ}\text{C}$  to  $40^{\circ}\text{C}$  are organisms that prefer to grow in the bodies of warm-blooded hosts.
3. **Thermophiles:** Bacterial species that will grow at  $35^{\circ}\text{C}$  and above. Two groups of thermophiles exist:
  - a. **Facultative thermophiles:** Organisms that will grow at  $37^{\circ}\text{C}$ , with an optimum growth temperature of  $45^{\circ}\text{C}$  to  $60^{\circ}\text{C}$ .
  - b. **Obligate thermophiles:** Organisms that will grow only at temperatures above  $50^{\circ}\text{C}$ , with optimum growth temperatures above  $60^{\circ}\text{C}$ .

The ideal temperature for specific enzymatic activities may not coincide with the optimum

growth temperature for a given organism. To understand this concept, you will investigate pigment production and carbohydrate fermentation by selected organisms at a variety of incubation temperatures.

1. The production of an endogenous red or magenta pigment by *Serratia marcescens* is determined by the presence of an orange to deep red coloration on the surface of the colonial growth.
2. Carbohydrate fermentation by *Saccharomyces cerevisiae* is indicated by the presence of gas, one of the end products of this fermentative process. Detection of this accumulated gas may be noted as an air pocket, of varying size, in an inverted inner vial (Durham tube) within the culture tube. Refer to Experiment 23 for a more extensive discussion of carbohydrate fermentation.

## CLINICAL APPLICATION

### Cold-Resistant Killers

The field of food science is highly concerned with the temperature-related growth patterns of bacteria. Refrigeration temperatures below  $4.4^{\circ}\text{C}$  are generally considered safe for the short-term storage of food, since most pathogenic bacteria grow very slowly below that temperature. However, some dangerous bacteria are resistant to cold. *Listeria monocytogenes*, which causes a flu-like illness and can be deadly, is capable of doubling its population every 36 hours, even at  $4.2^{\circ}\text{C}$ , and can still attain slow growth below  $2^{\circ}\text{C}$ . *Listeria*'s cold tolerance may be due to adaptive genes, prompting research into novel methods of controlling its growth at low temperatures.



## Materials

### Cultures

24- to 48-hour nutrient broth cultures of *Escherichia coli*, *Bacillus stearothermophilus*, *Pseudomonas savastanoi*, *Serratia marcescens*, and Sabouraud broth culture of *Saccharomyces cerevisiae*.

### Media

Per designated student group: four Trypticase soy agar plates and four Sabouraud broth tubes containing inverted Durham tubes.

### Equipment

Bunsen burner, inoculating loop, refrigerator set at 4°C, two incubators set at 37°C and 60°C, sterile Pasteur pipette, test tube rack, and glassware marking pencil.

## Procedure Lab One

- Score the underside of all plates into four quadrants with a glassware marker. Label each section with the name of the test organism to be inoculated. When labeling the cover of each plate, include the temperature of incubation (4°C, 20°C, 37°C, or 60°C).
- Aseptically inoculate each of the plates with *E. coli*, *B. stearothermophilus*, *P. savastanoi*, and *S. marcescens* by means of a single line of inoculation of each organism in its appropriately labeled section.
- Appropriately label the four Sabouraud broth tubes, including the temperatures of incubation as indicated above.
- Gently shake the *S. cerevisiae* culture to suspend the organisms. Using a sterile Pasteur pipette, aseptically add one drop of the culture into each of the four tubes of broth media.
- Incubate all plates in an inverted position and the broth cultures at each of the four experimental temperatures (4°C, 20°C, 37°C, or 60°C) for 24 to 48 hours.

## Procedure Lab Two

- In the chart provided in the Lab Report, complete the following:
  - Observe all the cultures for the presence of growth. Record your observations: (1+) for scant growth; (2+) for moderate growth; (3+) for abundant growth; and (–) for the absence of growth. Evaluate the amount of growth in the *S. cerevisiae* cultures by noting the degree of developed turbidity.
  - Observe the *S. marcescens* growth on all the plate cultures for the presence or absence of orange to deep red pigmentation. Record the presence of pigment on a scale of 1+ to 3+, and enter (–) for the absence of pigmentation.
  - Observe the *S. cerevisiae* cultures for the presence of a gas pocket in the Durham tube, which is indicative of carbohydrate fermentation. Record your observations using the following designations: (1+) for a minimal amount of gas; (2+) for a moderate amount of gas; (3+) for a large amount of gas; and (–) for the absence of gas.
  - Record and classify the cultures as psychrophiles, mesophiles, facultative thermophiles, or obligate thermophiles.

*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

	<i>Serratia marcescens</i>		<i>Pseudomonas savastanoi</i>	<i>Escherichia coli</i>	<i>Bacillus stearothermophilus</i>	<i>Saccharomyces cerevisiae</i>	
Temperature	Pigment	Growth	Growth	Growth	Growth	Growth	Gas
4°C (refrig.)							
20°C (room temp.)							
37°C (body temp.)							
60°C							
Classification							


Based on your observations of the *S. marcescens* and *S. cerevisiae* cultures, is the optimum growth temperature the ideal temperature for all cell activities? Explain.


## Review Questions

- In the following chart, indicate the types of organisms that would grow preferentially in or on various environments, and indicate the optimum temperature for their growth.

Environment	Type of Organism	Optimum Temperature
Ocean bottom near shore		
Ocean bottom near hot vent		
Hot sulfur spring		
Compost pile (middle)		
High mountain lake		
Center of an abscess		
Antarctic ice		

2. Explain the effects of temperatures above the maximum and below the minimum growth temperatures on cellular enzymes.

3.  If an organism grew at 20°C, explain how you would determine experimentally whether the organism was a psychrophile or a mesophile.

4.  Is it possible for thermophilic organisms to induce infections in warm-blooded animals? Explain.

### LEARNING OBJECTIVE

Once you have completed this experiment, you should be familiar with

1. The pH requirements of microorganisms.

## Principle

Growth and survival of microorganisms are greatly influenced by the pH of the environment, and all bacteria and other microorganisms differ as to their requirements. Based on their pH optima, microorganisms may be classified as acidophiles, neutrophiles, or alkalophiles (Figure 17.1). Each species has the ability to grow within a specific pH range, which may be broad or limited, with the most rapid growth occurring within a narrow optimum range. These specific pH needs reflect the organisms' adaptations to their natural environment. For example, enteric bacteria are capable of survival within a broad pH range, which is characteristic of their natural habitat, the digestive system. Bacterial blood parasites, on the other

hand, can tolerate only a narrow range; the pH of the circulatory system remains fairly constant at approximately 7.4.

Despite this diversity and the fact that certain organisms can grow at extremes of the pH scale, generalities can be made. The specific range for bacteria is between 4 and 9, with the optimum being 6.5 to 7.5. Fungi (molds and yeasts) prefer an acidic environment, with optimum activities at a pH of 4 to 6.

Because a neutral or nearly neutral environment is generally advantageous to the growth of microorganisms, the pH of the laboratory medium is frequently adjusted to approximately 7. Metabolic activities of the microorganism will result in the production of wastes, such as acids from carbohydrate degradation and alkali from protein breakdown, and these will cause shifts in pH that can be detrimental to growth.

To retard this shift, chemical substances that act as **buffers** are frequently incorporated when the medium is prepared. A commonly used **buffering system** involves the addition of equimolar concentrations of  $K_2HPO_4$ , a salt of a weak base, and  $KH_2PO_4$ , a salt of a weak acid. In a medium that has become acidic, the  $K_2HPO_4$  absorbs


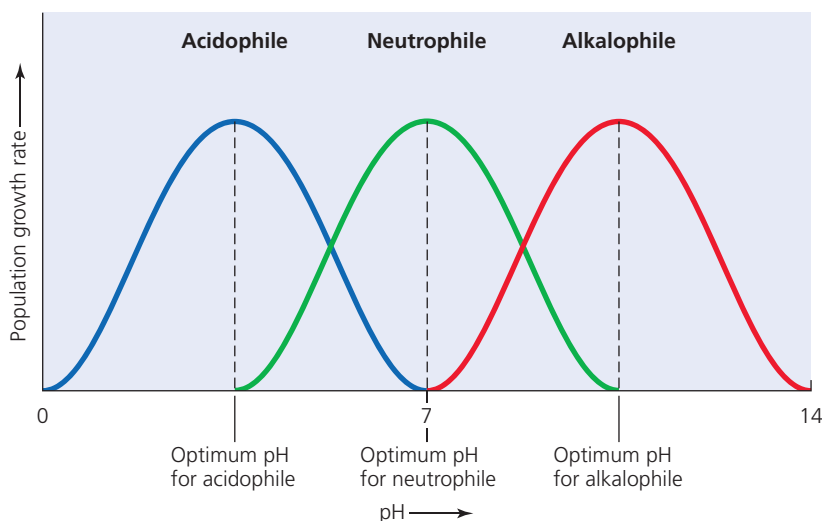



Figure 17.1 The effect of pH on the growth of microorganisms

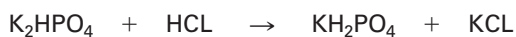excess  $\text{H}^+$  to form a weakly acidic salt and a potassium salt with the anion of the strong acid.




<b>Salt of a weak base</b>	<b>Strong acid</b>	<b>Salt of a weak acid</b>	<b>Potassium chloride salt</b>
----------------------------	--------------------	----------------------------	--------------------------------

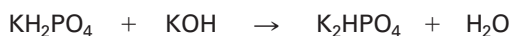In a medium that has become alkaline,  $\text{KH}_2\text{PO}_4$  releases  $\text{H}^+$ , which combines with the excess  $\text{OH}^-$  to form water, and the remaining anionic portion of the weakly acidic salt combines with the cation of the alkali.



<b>Salt of a weak acid</b>	<b>Strong base</b>	<b>Salt of a weak base</b>	<b>Water</b>
----------------------------	--------------------	----------------------------	--------------

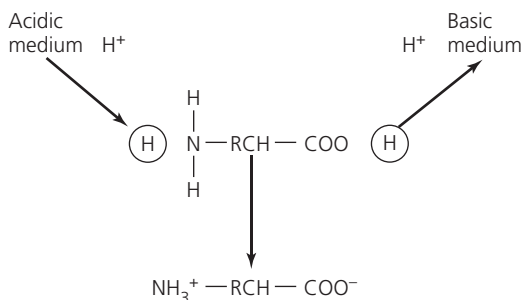Most media contain amino acids, peptones, and proteins, which can act as natural buffers because of their amphoteric nature. For example, amino acids are zwitterions, molecules in which the amino group and the carboxyl group ionize to form dipolar ions. These behave in the following manner:



## CLINICAL APPLICATION

### pH as a Defense Against Infection

Most bacteria grow best at a pH between 6.5 and 7.5, and fungi show optimal growth between a pH of 4 and 6. Many microorganisms are not able to cause stomach infections because the pH of the stomach is 2.0, resembling that of hydrochloric acid. In this way, the acid of the stomach acts as a defense against infection. By the same token, the pH of the skin varies between 4 and 7, with lower ranges (around 5) being the most common, helping prevent many infections of the skin.

## AT THE BENCH



## Materials

### Cultures

Saline suspensions of 24-hour nutrient broth cultures, adjusted to an absorbance (A) of 0.05 at a wavelength of 600 nm, of *Alcaligenes faecalis*, *Escherichia coli*, and *Saccharomyces cerevisiae*.

### Media

Per designated student group: 12 Trypticase soy broth (TSB) tubes, three at each of the following pH designations: 3, 6, 7, and 9. The pH has been adjusted with 1N sodium hydroxide or 1N hydrochloric acid.

### Equipment

Bunsen burner, sterile 1-ml pipettes, mechanical pipetting device, Bausch & Lomb Spectronic 20 spectrophotometer, test tube rack, and glassware marking pencil.

## Procedure Lab One

- Using a sterile pipette, inoculate a series of the appropriately labeled TSB tubes of media, pH values of 3, 6, 7, and 9, with *E. coli* by adding 0.1 ml of the saline culture to each.
- Repeat Step 1 for the inoculation of *A. faecalis* and *S. cerevisiae*, using a new sterile pipette each time.
- Incubate the *A. faecalis* and *E. coli* cultures for 24 to 48 hours at 37°C and the *S. cerevisiae* cultures for 48 to 72 hours at 25°C.

## Procedure Lab Two

- Using the spectrophotometer as described in Experiment 14, determine the absorbance of all cultures. Record the readings in the chart provided in the Lab Report.
- In the second chart provided in the Lab Report, summarize your results as to the overall range and optimum pH of each organism studied.






3. Why are proteins and amino acids considered to be natural buffers?

4. Explain why microorganisms differ in their pH requirements.

5. Will all microorganisms grow optimally at a neutral pH? Explain.

6.  You are instructed to grow *E. coli* in a chemically defined medium containing glucose and  $\text{NH}_4\text{Cl}$  as the carbon and nitrogen sources and also in nutrient broth that contains beef extract and peptone. Both media are adjusted to a pH of 7. With turbidity as an index for the amount of growth in each of the cultures, the following spectrophotometric readings are obtained following incubation:

Time (hours)	ABSORBANCE READINGS	
	Chemically Defined Medium	Nutrient Broth Medium
6	0.100	0.100
12	0.300	0.500
18	0.275	0.900
24	0.125	1.500

Based on the above data, explain why *E. coli* ceased growing in the chemically defined medium but continued to grow in the nutrient broth.

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be familiar with

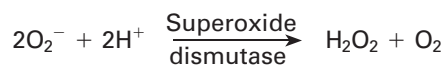
1. The diverse atmospheric oxygen requirements of microorganisms.

## Principle

Microorganisms exhibit great diversity in their ability to use free oxygen ( $O_2$ ) for cellular respiration. These variations in  $O_2$  requirements reflect the differences in biooxidative enzyme systems present in the various species. Microorganisms can be classified into one of five major groups according to their  $O_2$  needs:

1. **Aerobes** require the presence of atmospheric oxygen for growth. Their enzyme system necessitates use of  $O_2$  as the final hydrogen (electron) acceptor in the complete oxidative degradation of high-energy molecules such as glucose.
2. **Microaerophiles** require limited amounts of atmospheric oxygen for growth. Oxygen in excess of the required amount appears to block the activities of their oxidative enzymes and results in death.
3. **Obligate anaerobes** require the absence of free oxygen for growth because their oxidative enzyme system requires the presence of molecules other than  $O_2$  to act as the final hydrogen (electron) acceptor. In these organisms, as in aerobes, the presence of atmospheric oxygen results in the formation of toxic metabolic end products, such as superoxide,  $O_2$ , a free radical of oxygen. However, these organisms lack the enzymes superoxide dismutase
4. **Aerotolerant anaerobes** are fermentative organisms, and therefore they do not use  $O_2$  as a final electron acceptor. Unlike the obligate anaerobes, they produce catalase and/or superoxide dismutase, and thus they are not killed by the presence of  $O_2$ . Hence, these organisms are anaerobes that are termed aerotolerant.
5. **Facultative anaerobes** can grow in the presence or absence of free oxygen. They preferentially use oxygen for aerobic respiration. However, in an oxygen-poor environment, cellular respiration may occur anaerobically, utilizing such compounds as nitrates ( $NO_3^-$ ) or sulfates ( $SO_4^{2-}$ ) as final hydrogen acceptors, or via a fermentative pathway (refer to Experiment 23).

and catalase, whose function is to degrade the superoxide to water and oxygen as follows:



In the absence of these enzymes, small amounts of atmospheric oxygen are lethal, and these organisms are justifiably called obligate anaerobes.

4. **Aerotolerant anaerobes** are fermentative organisms, and therefore they do not use  $O_2$  as a final electron acceptor. Unlike the obligate anaerobes, they produce catalase and/or superoxide dismutase, and thus they are not killed by the presence of  $O_2$ . Hence, these organisms are anaerobes that are termed aerotolerant.
5. **Facultative anaerobes** can grow in the presence or absence of free oxygen. They preferentially use oxygen for aerobic respiration. However, in an oxygen-poor environment, cellular respiration may occur anaerobically, utilizing such compounds as nitrates ( $NO_3^-$ ) or sulfates ( $SO_4^{2-}$ ) as final hydrogen acceptors, or via a fermentative pathway (refer to Experiment 23).

The oxygen needs of microorganisms can be determined by noting their growth distributions following a **shake-tube inoculation**. This procedure requires introduction of the inoculum into a melted agar medium, shaking of the test tube to disperse the microorganisms throughout the agar, and rapid solidification of the medium to ensure that the cells remain dispersed. Following incubation, the growth distribution indicates the organisms' oxygen requirements. Aerobes exhibit surface growth, whereas anaerobic growth is limited to the bottom of the deep tube. Facultative

anaerobes, because of their indifference to the presence or absence of oxygen, exhibit growth throughout the medium. Microaerophiles grow in a zone slightly below the surface. **Figure 18.1** illustrates the shake-tube inoculation procedure and the distribution of growth following an appropriate incubation period.

## CLINICAL APPLICATION

### Differentiating Aerobes and Anaerobes

Samples suspected of containing anaerobes need to be handled carefully and transported promptly to a lab, where they are typically inoculated onto anaerobic blood agar plates and anaerobic broth, as well as onto MacConkey agar and an aerobic blood plate. Growth on aerobic or anaerobic agars will determine oxygen requirements, while comparable growth on both aerobic and anaerobic media suggests a facultative anaerobe.

## AT THE BENCH



## Materials

### Cultures

24- to 48-hour nutrient broth cultures of *Staphylococcus aureus*, *Corynebacterium xerosis*, and *Enterococcus faecalis*; 48- to 72-hour Sabouraud broth cultures of *Saccharomyces cerevisiae* and *Aspergillus niger*; and a 48-hour thioglycollate broth culture of *Clostridium sporogenes*.

### Media

Six brain heart infusion agar deep tubes per designated student group.

## Equipment

Bunsen burner, waterbath, iced waterbath, thermometer, sterile Pasteur pipettes, test tube rack, and glassware marking pencil.

## Procedure Lab One

1. Liquefy the sterile brain heart infusion agar by boiling in a waterbath at 100°C.
2. Cool molten agar to 45°C; check temperature with a thermometer inserted into the waterbath.

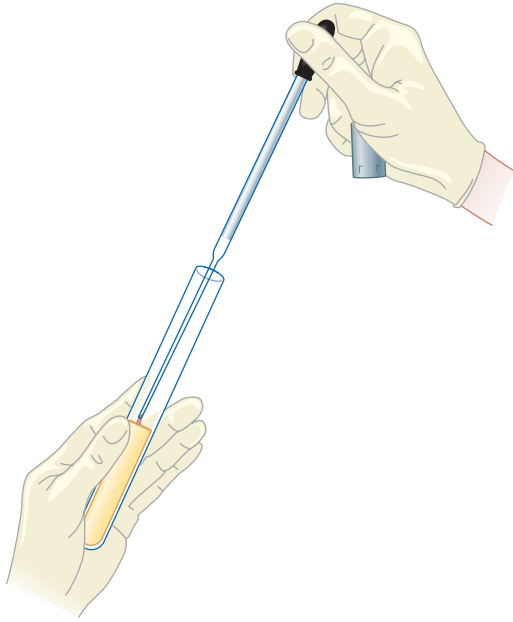
## Determining Oxygen Requirements

1. Using aseptic technique, inoculate each experimental organism by introducing two drops of the culture from a sterile Pasteur pipette into the appropriately labeled tubes of molten agar.
2. Vigorously rotate the freshly inoculated molten infusion agar between the palms of the hands to distribute the organisms.
3. Place inoculated test tubes in an upright position in the iced waterbath to solidify the medium rapidly.
4. Incubate the *S. aureus*, *C. xerosis*, *E. faecalis*, and *C. sporogenes* cultures for 24 to 48 hours at 37°C and the *A. niger* and *S. cerevisiae* cultures for 48 to 72 hours at 25°C.


## Procedure Lab Two

1. Observe each of the experimental cultures for the distribution of growth in each tube.
2. Record your observations and your determination of the oxygen requirements for each of the experimental species in the chart provided in the Lab Report.

## PROCEDURE



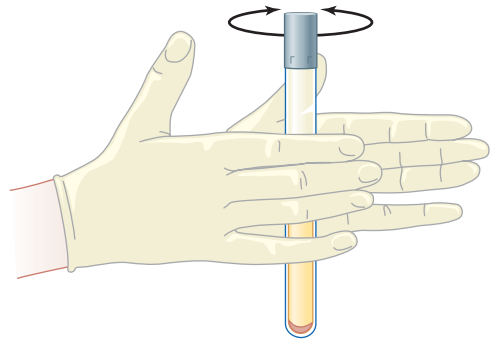**1** Transfer two drops of inoculum from the test culture into a melted agar deep tube.




**2** Disperse the organisms throughout the molten agar medium by rapidly rotating the tube between the palms of your hands.

**3** Cool rapidly by immersion in an iced waterbath.

**4** Incubate at 37°C.

### Distribution of Growth

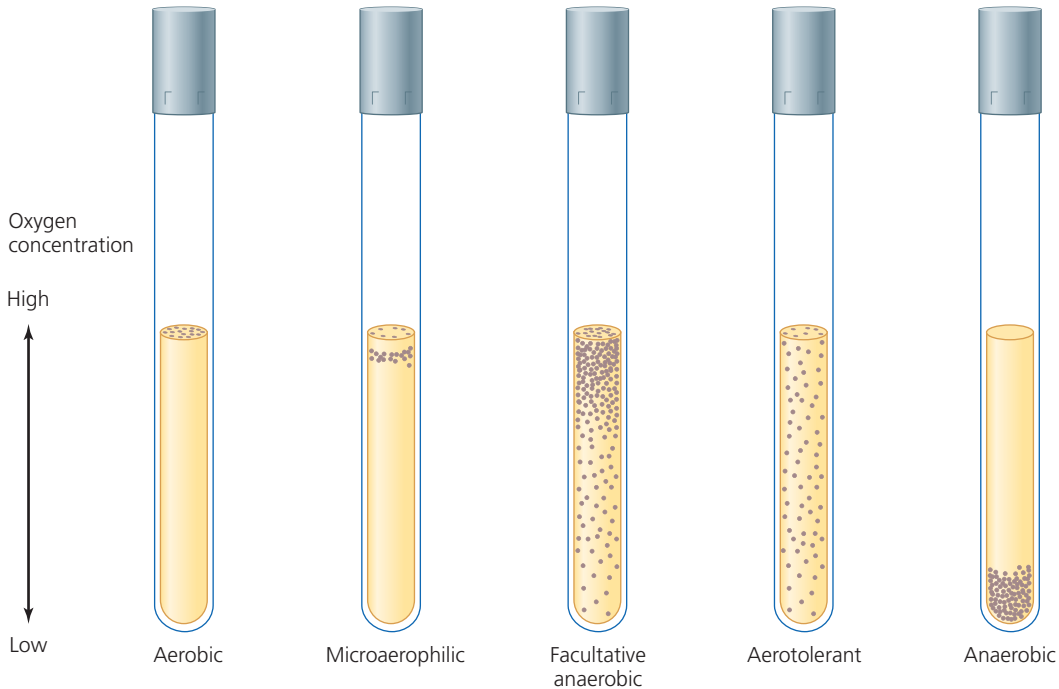

**Figure 18.1** Procedure for determination of oxygen requirements

*This page intentionally left blank*

Name: \_\_\_\_\_




Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

Species	Distribution of Growth	Classification According to Oxygen Requirement
<i>S. aureus</i>		
<i>C. xerosis</i>		
<i>E. faecalis</i>		
<i>A. niger</i>		
<i>S. cerevisiae</i>		
<i>C. sporogenes</i>		

## Review Questions

1. Why is it necessary to place the inoculated molten agar cultures in an iced waterbath for their rapid solidification?
2. As indicated by its oxygen requirements, which group of microorganisms has the most extensive bioenergetic enzyme system? Explain.
3. Account for the inability of aerobes to grow in the absence of O<sub>2</sub>.

4. Account for the subsurface growth of microaerophiles in a shake-tube culture.
  
5.  Consider the culture type in which growth was distributed throughout the entire medium and explain why the growth was more abundant toward the surface of the medium in some cultures, whereas other cultures showed an equal distribution of growth throughout the tubes.
  
6.  Account for the fact that the *C. sporogenes* culture showed a separation within the medium or an elevation of the medium from the bottom of the test tube.
  
7.  Your instructor asks you to explain why the *Streptococcus* species that are catalase negative are capable of growth in the presence of oxygen. How would you respond?

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be familiar with


1. The methods for cultivation of anaerobic organisms.

## Principle

Microorganisms differ in their abilities to use oxygen for cellular respiration. **Respiration** involves the oxidation of substrates for energy necessary to life. A substrate is **oxidized** when it loses a hydrogen ion and its electron ( $H^+e^-$ ). Since the  $H^+e^-$  cannot remain free in the cell, it must immediately be picked up by an electron acceptor, which becomes reduced. Therefore reduction means gaining the  $H^+e^-$ . These are termed **oxidation-reduction (redox)** reactions. Some microorganisms have enzyme systems in which oxygen can serve as an electron acceptor, thereby being reduced to water. These cells have high oxidation-reduction potentials; others have low potentials and must use other substances as electron acceptors.

The enzymatic differences in microorganisms are explained more fully in the section dealing with metabolism (see Part 5). This discussion is limited to cultivation of the strict anaerobes, which cannot be cultivated in the presence of atmospheric oxygen (Figure 19.1). The procedure is somewhat more difficult because it involves sophisticated equipment and media enriched with substances that lower the redox potential. Figure 19.2 shows some of the methods available for anaerobic cultivation.

The following experiment uses fluid thioglycollate medium and the GasPak™ anaerobic system.

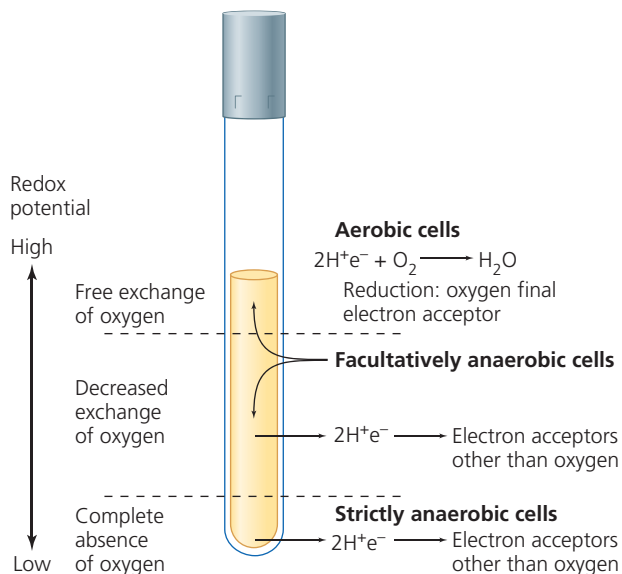


**Figure 19.1** Illustration of redox potentials in an agar deep tube

## CLINICAL APPLICATION

### Oxygen as a Treatment?

The causative agent of gas gangrene, *Clostridium perfringens*, is an anaerobic bacterium that thrives in wounds deprived of circulation and oxygen and can cause limb loss and death. Treatment may involve amputation or surgical removal of infected tissue. Doctors may also prescribe therapy using enriched oxygen delivered to the patient in a hyperbaric chamber. This allows the blood to carry more oxygen to the wounds, slowing the growth of anaerobic microbes. Patients typically undergo five 90-minute sessions lying in a chamber pressurized to 2.5 atmospheres, possibly alleviating the need for surgery.



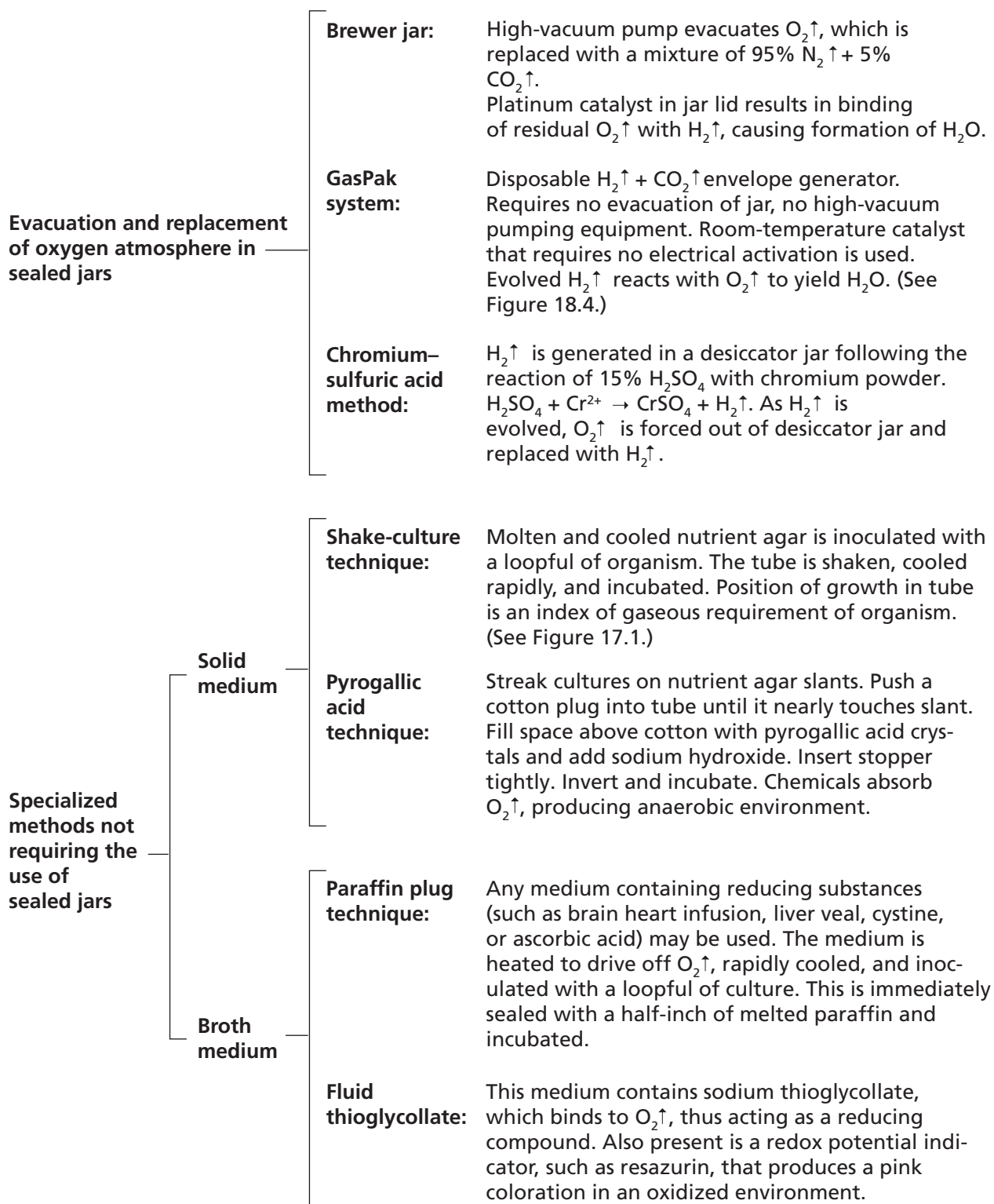

**Figure 19.2** Methods for the cultivation of anaerobic microorganisms



## Materials

### Cultures

24- to 48-hour nutrient broth cultures of *Bacillus cereus*, *Escherichia coli*, and *Micrococcus luteus*; and 48-hour thioglycollate broth culture of *Clostridium sporogenes*.

### Media

Per designated student group: four screw-cap tubes of fluid thioglycollate medium and four nutrient agar plates.

### Equipment

Bunsen burner, inoculating loop, GasPak anaerobic system, test tube rack, and glassware marking pencil.

## Procedure Lab One

### Fluid Thioglycollate Medium

1. For the performance of this procedure, the fluid thioglycollate medium must be fresh. Freshness is indicated by the absence of a pink color in the upper one-third of the medium. If this coloration is present, loosen the screw caps and place the tubes in a boiling water bath for 10 minutes to drive off the

dissolved  $O_2$  from the medium. Cool the tubes to  $45^\circ C$  before inoculation.


2. Aseptically inoculate the appropriately labeled tubes of thioglycollate with their respective test organisms by means of loop inoculations to the depths of the media.
3. Incubate the cultures for 24 to 48 hours at  $37^\circ C$ .

The appearance of the growth of organisms according to their gaseous requirements in thioglycollate medium is shown in **Figure 19.3**.

### GasPak Anaerobic Technique


The GasPak system, shown in **Figure 19.4**, is a contemporary method for the exclusion of oxygen from a sealed jar used for incubation of anaerobic cultures in a nonreducing medium. This system uses a GasPak generator that consists of a foil package that generates hydrogen and carbon dioxide upon the addition of water. A palladium catalyst in the lid of the jar combines the evolved hydrogen with residual oxygen to form water, thereby creating a carbon dioxide environment within the jar that is conducive for anaerobic growth. The establishment of anaerobic conditions is verified by the color change of a methylene blue indicator strip in the jar. This blue indicator becomes colorless in the absence of oxygen.

1. With a glassware marking pencil, divide the bottom of each nutrient agar plate into two sections.

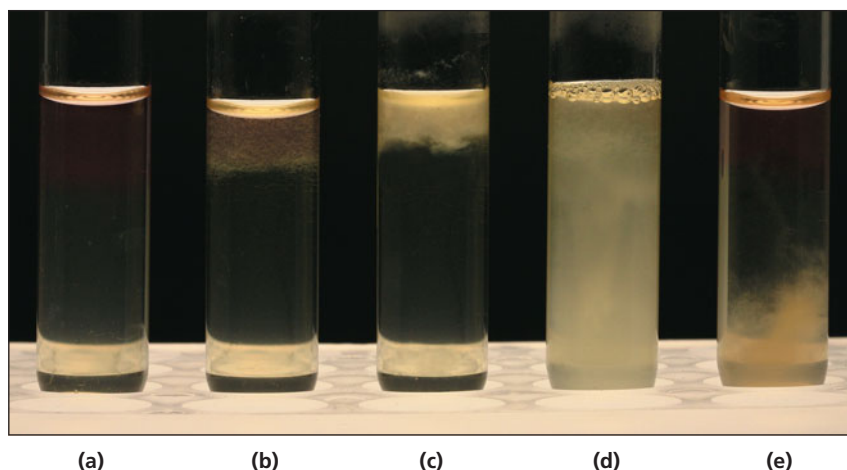

**Figure 19.3** Bacterial growth patterns in thioglycollate broth tubes.

(a) Uninoculated control. (b,c) Uniform growth indicates facultative anaerobic bacteria. (d) Bubbles indicate gas-producing bacteria. (e) Bottom growth indicates anaerobic bacteria.

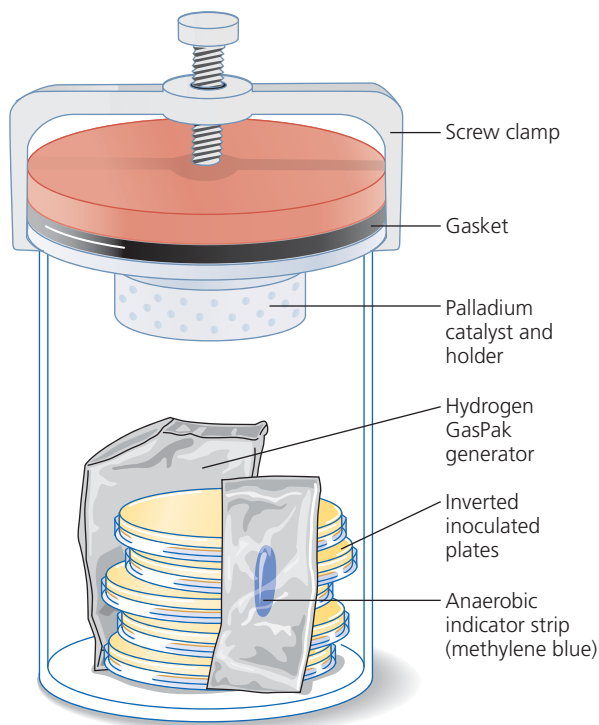

**Figure 19.4 GasPak system**




2. Label each section on two plates with the name of the organism to be inoculated.
3. Repeat Step 2 to prepare a duplicate set of cultures.
4. Using aseptic technique, make a single-line streak inoculation of each test organism in its respectively labeled section on both sets of plates.

5. Tear off the corner of the hydrogen and carbon dioxide gas generator and insert this inside the GasPak jar.
6. Place one set of plate cultures in an inverted position inside the GasPak chamber.
7. Expose the anaerobic indicator strip and place it inside the anaerobic jar so that the wick is visible from the outside.
8. With a pipette, add the required 10 ml of water to the gas generator and quickly seal the chamber with its lid.
9. Place the sealed jar in an incubator at 37°C for 24 to 48 hours. After several hours of incubation, observe the indicator strip for a color change to colorless, which is indicative of anaerobic conditions.
10. Incubate the duplicate set of plates in an inverted position for 24 to 48 hours at 37°C under aerobic conditions.

## Procedure Lab Two

1. Observe the fluid thioglycollate cultures, GasPak system, and aerobically incubated plate cultures for the presence of growth. Record your results in the chart provided in the Lab Report.
2. Based on your observation, record the oxygen requirement classification of each test organism as anaerobe, facultative anaerobe, or aerobe.



3.  Heroin addicts have a high incidence of *Clostridium tetani* infections. Discuss the reasons for the development of this type of infection in these IV drug users.
4.  While you are working in your garden, a tine of the pitchfork accidentally produces a deep puncture wound in the calf of your leg. Discuss the type of infectious process you would be primarily concerned about and why.
5.  The physician who treats your puncture wound opts to insert a drain before applying the dressing. What is the rationale for the insertion of the drain?

# Serial Dilution–Agar Plate Procedure to Quantitate Viable Cells

## EXPERIMENT 20

### LEARNING OBJECTIVES

Once you have completed this experiment, you should understand

1. The diverse methods used to determine the number of cells in a bacterial culture.
2. How to determine quantitatively the number of viable cells in a bacterial culture.

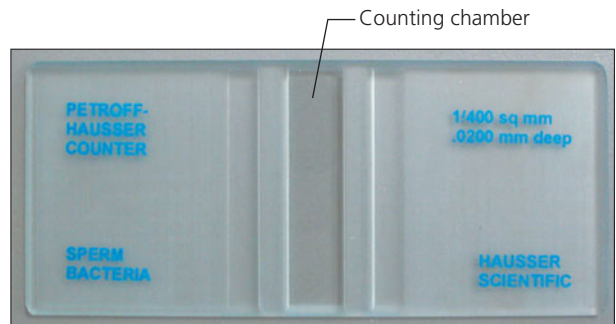
## Principle

Studies involving the analysis of materials such as food, water, milk, and, in some cases, air require quantitative enumeration of microorganisms in the substances. Many methods have been devised to accomplish this, including direct microscopic counts, use of an electronic cell counter such as the Coulter Counter<sup>®</sup>, chemical methods for estimating cell mass or cellular constituents, turbidimetric measurements for increases in cell mass, and the serial dilution–agar plate method.


## Direct Microscopic Counts

Direct microscopic counts require the use of a specialized slide called the **Petroff-Hausser counting chamber**, in which an aliquot of a bacterial suspension is counted and the total number of cells is determined mathematically. The Petroff-Hausser counting chamber is a thick glass microscope slide with a chamber 0.02 mm (1/50 mm) deep in the center. The chamber contains an etched grid and has improved Neubauer rulings (1/400 square mm). The slide and the counting chamber are illustrated in **Figure 20.1**.

The rulings cover  $9 \text{ mm}^2$ . The boundary lines (Neubauer rulings) are the center lines of the groups of three. The center square millimeter is ruled into groups of 16 small squares, and each group is separated by triple lines, the middle one of which is the boundary. The ruled surface is 0.02 mm below the cover glass, which makes



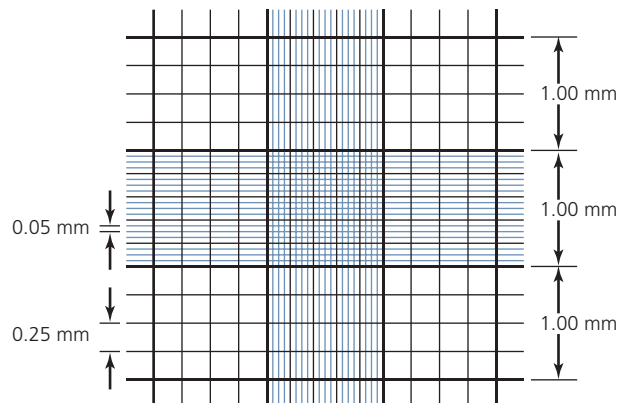(a) Petroff-Hausser counting chamber



(b) Petroff-Hausser counting chamber grid


**Figure 20.1** The Petroff-Hausser chamber

the volume over a square millimeter  $0.02 \text{ mm}^3$  (cubic mm). All cells are counted in this square millimeter.

The number of bacterial cells counted is calculated as follows:

$$\text{number of cells per mm} = \text{number of cells counted} \times \text{dilution} \times 50,000$$

[The factor of 50,000 is used in order to determine the cell count for 1 ml:  $1 \text{ ml} = 1000 \text{ mm}^3 = (50 \text{ times the chamber depth of } 0.02 \text{ mm}) \times 1000$ .] Although rapid, a direct count has the disadvantages that both living and dead cells are counted and that it is not sensitive to populations of fewer than 1 million cells.

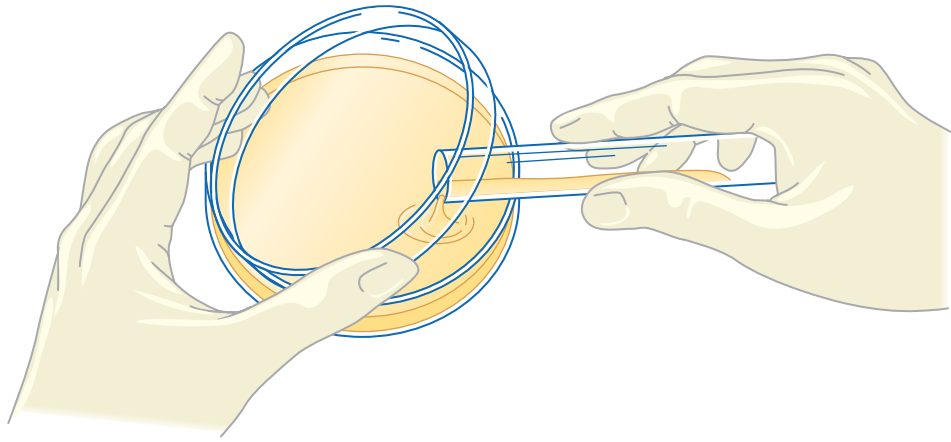

**Figure 20.2** Pour-plate technique

**Breed smears** are used mainly to quantitate bacterial cells in milk. Using stained smears confined to a 1-square-millimeter ruled area of the slide, the total population is determined mathematically. This method also fails to discriminate between viable and dead cells.

### Electronic Cell Counters

The **Coulter Counter** is an example of an instrument capable of rapidly counting the number of cells suspended in a conducting fluid that passes through a minute orifice through which an electric current is flowing. Cells, which are nonconductors, increase the electrical resistance of the conducting fluid, and the resistance is electronically recorded, enumerating the number of organisms flowing through the orifice. In addition to its inability to distinguish between living and dead cells, the apparatus is also unable to differentiate inert particulate matter from cellular material.

### Chemical Methods

While not considered means of direct quantitative analysis, chemical methods may be used to indirectly measure increases both in protein concentration and in DNA production. In addition, cell mass can be estimated by dry weight determination of a specific aliquot of the culture. Measurement of certain metabolic parameters may also be used to quantitate bacterial populations. The amount of oxygen consumed (oxygen uptake) is directly proportional to the increasing number of vigorously growing aerobic cells, and the rate of carbon dioxide production is related to increased growth of anaerobic organisms.

### Spectrophotometric Analysis

Increased turbidity in a culture is another index of growth. With turbidimetric instruments, the amount of transmitted light decreases as the cell population increases, and the decrease in radiant energy is converted to electrical energy and indicated on a galvanometer. This method is rapid but limited because sensitivity is restricted to microbial suspensions of 10 million cells or greater.

### Serial Dilution–Agar Plate Analysis

While all these methods may be used to enumerate the number of cells in a bacterial culture, the major disadvantage common to all is that the total count includes dead as well as living cells. Sanitary and medical microbiology at times require determination of viable cells. To accomplish this, the serial dilution–agar plate technique is used. Briefly, this method involves serial dilution of a bacterial suspension in sterile water blanks, which serve as a diluent of known volume. Once diluted, the suspensions are placed on suitable nutrient media. The **pour-plate technique**, illustrated in **Figure 20.2**, is the procedure usually employed. Molten agar, cooled to 45°C, is poured into a Petri dish containing a specified amount of the diluted sample. Following addition of the molten-then-cooled agar, the cover is replaced, and the plate is gently rotated in a circular motion to achieve uniform distribution of microorganisms. This procedure is repeated for all dilutions to be plated. Dilutions should be plated in duplicate for greater accuracy, incubated overnight, and counted on a **Quebec colony counter** either by hand or by an electronically modified version of this instrument. **Figure 20.3** is an illustration of this apparatus for counting colonies.



**Figure 20.3** Quebec colony counter for the enumeration of bacterial colonies

Plates suitable for counting must contain not fewer than 30 nor more than 300 colonies. See **Figure 20.4**. The total count of the suspension is obtained by multiplying the number of cells per plate by the dilution factor, which is the reciprocal of the dilution.


Advantages of the serial dilution–agar plate technique are as follows:

1. Only viable cells are counted.
2. It allows isolation of discrete colonies that can be subcultured into pure cultures, which may then be easily studied and identified.

Disadvantages of this method are as follows:

1. Overnight incubation is necessary before colonies develop on the agar surface.
2. More glassware is used in this procedure.
3. The need for greater manipulation may result in erroneous counts due to errors in dilution or plating.

The following experiment uses the pour-plate technique for plating serially diluted culture samples.

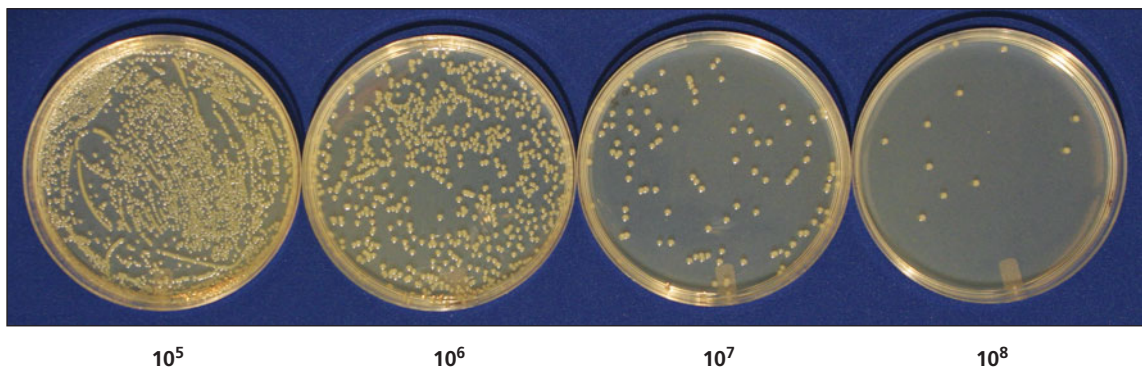

**Figure 20.4** Agar plating method for viable cell counts using dilutions  $1 \times 10^5$ ,  $1 \times 10^6$ ,  $1 \times 10^7$ , and  $1 \times 10^8$

## CLINICAL APPLICATION

### The Multiple Uses of Cell Counts

Determining how many cells are present in a sample is of particular importance in the food and dairy industries, which monitor the number and types of bacteria in their products. Elevated bacteria counts can indicate a sick animal, inadequate sanitation, or improper storage. Viable cell counts are also used in water treatment facilities as well as in wineries and breweries, where the number of yeast cells is monitored. In medical laboratories, sometimes the number of cells and growth rates are used to determine antimicrobial sensitivity as well as the course of infection.

## AT THE BENCH



## Materials

### Culture

24- to 48-hour nutrient broth culture of *Escherichia coli*.

### Media


Per designated student group: six 20-ml nutrient agar deep tubes and seven sterile 9-ml water blanks.

### Equipment

Hot plate, waterbath, thermometer, test tube rack, Bunsen burner, sterile 1-ml serological pipettes, mechanical pipetting device, sterile Petri dishes, Quebec colony counter, manual hand counter, disinfectant solution in a 500-ml beaker, glassware marking pencil, turntable, bent glass rod, and beaker with 95% alcohol.



## PROCEDURE

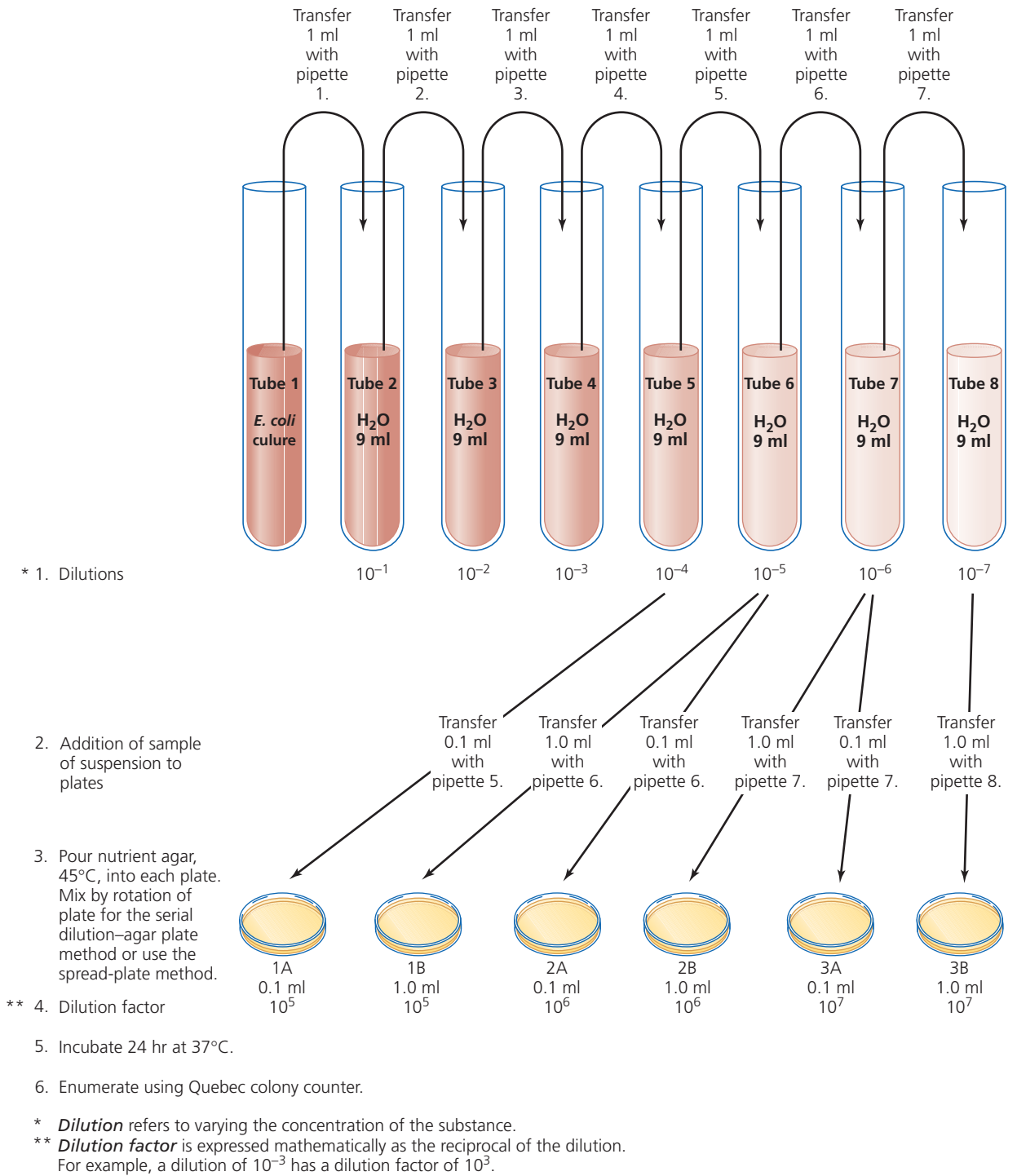


**Figure 20.5** Serial dilution-agar plate procedure

## Procedure Lab One

This procedure is illustrated in **Figure 20.5**. A photograph of the dilutions is shown in **Figure 20.6**.

1. Liquefy six agar deep tubes in an autoclave or by boiling. Cool the molten agar tubes and maintain in a water bath at 45°C.
2. Label the *E. coli* culture tube with the number 1 and the seven 9-ml water blanks as numbers 2 through 8. Place the labeled tubes in a test tube rack. Label the Petri dishes 1A, 1B, 2A, 2B, 3A, and 3B.
3. Mix the *E. coli* culture (Tube 1) by rolling the tube between the palms of your hands to ensure even dispersal of cells in the culture.
4. With a sterile pipette, aseptically transfer 1 ml from the bacterial suspension, Tube 1, to water blank Tube 2. Discard the pipette in the beaker of disinfectant. The culture has been diluted 10 times to  $10^{-1}$ .
5. Mix Tube 2 and, with a fresh pipette, transfer 1 ml to Tube 3. Discard the pipette. The culture has been diluted 100 times to  $10^{-2}$ .
6. Mix Tube 3 and, with a fresh pipette, transfer 1 ml to Tube 4. Discard the pipette. The culture has been diluted 1000 times to  $10^{-3}$ .
7. Mix Tube 4 and, with a fresh pipette, transfer 1 ml to Tube 5. Discard the pipette. The culture has been diluted 10,000 times to  $10^{-4}$ .
8. Mix Tube 5 and, with a fresh pipette, transfer 0.1 ml of this suspension to Plate 1A. Return the pipette to Tube 5 and transfer 1 ml to Tube 6. Discard the pipette. The culture has been diluted 100,000 times to  $10^{-5}$ .
9. Mix Tube 6 and, with a fresh pipette, transfer 1 ml of this suspension to Plate 1B. Return the pipette to Tube 6 and transfer 0.1 ml to Plate 2A. Return the pipette to Tube 6 and transfer 1 ml to Tube 7. Discard the pipette. The culture has been diluted 1,000,000 times to  $10^{-6}$ .
10. Mix Tube 7 and, with a fresh pipette, transfer 1 ml of this suspension to Plate 2B. Return the pipette to Tube 7 and transfer 0.1 ml to Plate 3A. Return the pipette to Tube 7 and transfer 1 ml to Tube 8. Discard the pipette. The culture has been diluted 10,000,000 times to  $10^{-7}$ .
11. Mix Tube 8 and, with a fresh pipette, transfer 1 ml of this suspension to Plate 3B. Discard the pipette. The dilution procedure is now complete.
12. Check the temperature of the molten agar medium to be sure the temperature is 45°C. Remove a tube from the waterbath and wipe the outside surface dry with a paper towel. Using the pour-plate technique, pour the agar into Plate 1A as shown in Figure 20.2 and rotate the plate gently to ensure uniform distribution of the cells in the medium.
13. Repeat Step 12 for the addition of molten nutrient agar to Plates 1B, 2A, 2B, 3A, and 3B.
14. Once the agar has solidified, incubate the plates in an inverted position for 24 hours at 37°C.

*Note: If desired, the spread-plate technique (Experiment 3) may be substituted for the agar pour-plate method described in this experiment. In this case, the dilutions may be placed on the surface of the hardened agar with a sterile pipette and distributed over the surface by means of a bent glass rod and turntable as illustrated in Figure 3.3. Following incubation, cell counts may be made as described in Experiment 20. The resultant cell counts should be the same with either system. The main difference is that there will be no subsurface colonies in the spread-plate method.*

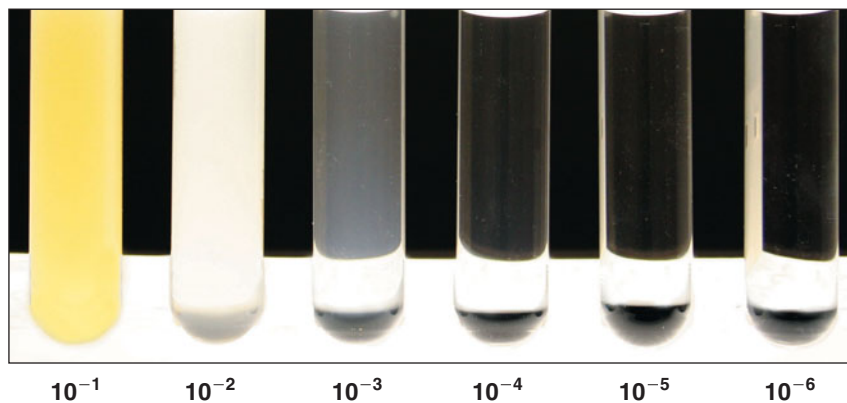

**Figure 20.6** Serial dilution of bacterial culture for quantitation of viable cell numbers

## Procedure Lab Two

1. Using a Quebec colony counter and a mechanical hand counter, observe all colonies on plates. Statistically valid plate counts are only obtained from bacterial cell dilutions that yield between 30 and 300 colonies. Plates with more than 300 colonies cannot be counted and are designated as **too numerous to count—TNTC**; plates with fewer than 30 colonies are designated as **too few to count—TFTC**. Count only plates containing between 30 and 300 colonies. Remember to count all subsurface as well as surface colonies.
2. The number of organisms per ml of original culture is calculated by multiplying the number of colonies counted by the dilution factor:


$$\frac{\text{number of cells per ml}}{\text{number of colonies} \times \text{dilution factor}}$$

### *Examples:*

- a. Colonies per plate = 50  
Dilution factor =  $1 : 1 \times 10^6$  (1 : 1,000,000)  
Volume of dilution added to plate = 1 ml  
 $50 \times 1,000,000 = 50,000,000$  or  
 $(5 \times 10^7)$  CFUs/ml  
(colony-forming units)

- b. Colonies per plate = 50  
Dilution factor =  $1 : 1 \times 10^5$  (1 : 1,00,000)  
Volume of dilution added to plate = 0.1 ml  
 $50 \times 100,000 = 5,000,000$  ( $5 \times 10^6$ )  
cells/0.1 ml  
 $5,000,000 \times 10 = 50,000,000$   
 $(5 \times 10^7)$  CFUs/ml
3. Record your observations and calculated bacterial counts per ml of sample in the Lab Report.
  4. Since the dilutions plated are replicates of each other, determine the average of the duplicate bacterial counts per ml of sample and record in the chart provided in the Lab Report.



3. What are the advantages and disadvantages of the serial dilution–agar plate procedure?
  
  
  
  
  
  
  
  
  
  
4. If 0.1 ml of a  $1 \times 10^{-6}$  dilution plate contains 56 colonies, calculate the number of cells per ml of the original culture.
  
  
  
  
  
  
  
  
  
  
5. How would you record your observation of a plate containing 305 colonies? A plate with 15 colonies?
  
  
  
  
  
  
  
  
  
  
6. Explain the chemical methods for measuring cell growth.
  
  
  
  
  
  
  
  
  
  
7.  Your instructor asks you to determine the number of organisms in a water sample. Observation of your dilution plates reveals the presence of spreading colonial forms on some of the culture plates. What is the rationale for the elimination of these plate counts from your experimental data?

### LEARNING OBJECTIVES

Once you have completed this experiment, you should be able to


1. Understand the population growth dynamics of bacterial cultures.
2. Plot a bacterial growth curve.
3. Determine the generation time of a bacterial culture from the bacterial growth curve.

## Principle

**Bacterial population growth studies** require inoculation of viable cells into a sterile broth medium and incubation of the culture under optimum temperature, pH, and gaseous conditions. Under these conditions, the cells will reproduce rapidly and the dynamics of the microbial growth can be charted in a population growth curve, which is constructed by plotting the increase in cell numbers versus time of incubation. The curve can be used to delineate stages of the growth cycle. It also facilitates measurement of cell numbers and the rate of growth of a particular organism under standardized conditions as expressed by its **generation time**, the time required for a microbial population to double.

The stages of a typical growth curve (Figure 21.1) are as follows:


1. **Lag phase:** During this stage the cells are adjusting to their new environment. Cellular metabolism is accelerated, resulting in rapid biosynthesis of cellular macromolecules, primarily enzymes, in preparation for the next phase of the cycle. Although the cells are increasing in size, there is no cell division and therefore no increase in numbers.
2. **Logarithmic (log) phase:** Under optimum nutritional and physical conditions, the physiologically robust cells reproduce at a uniform

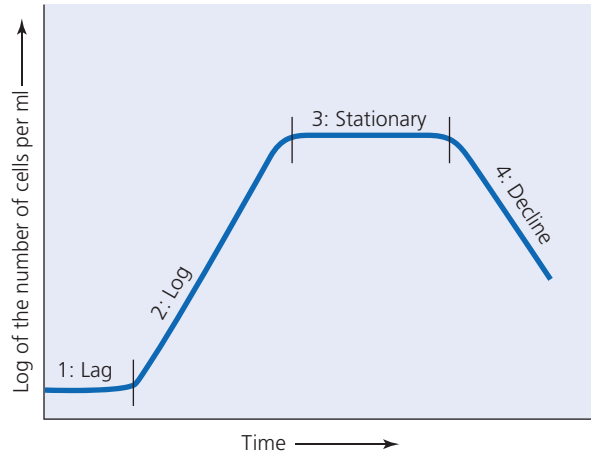

**Figure 21.1** Population growth curve

and rapid rate by binary fission. Thus, there is a rapid exponential increase in population, which doubles regularly until a maximum number of cells is reached. The time required for the population to double is the generation time. The length of the log phase varies, depending on the organisms and the composition of the medium. The average may be estimated to last 6 to 12 hours.

3. **Stationary phase:** During this stage, the number of cells undergoing division is equal to the number of cells that are dying. Therefore there is no further increase in cell number, and the population is maintained at its maximum level for a period of time. The primary factors responsible for this phase are the depletion of some essential metabolites and the accumulation of toxic acidic or alkaline end products in the medium.
4. **Decline, or death, phase:** Because of the continuing depletion of nutrients and buildup of metabolic wastes, the microorganisms die at a rapid and uniform rate. The decrease in population closely parallels its increase during the log phase. Theoretically, the entire population should die during a time interval equal to that of the log phase. This does not occur, however, since a small number of highly resistant organisms persist for an indeterminate length of time.

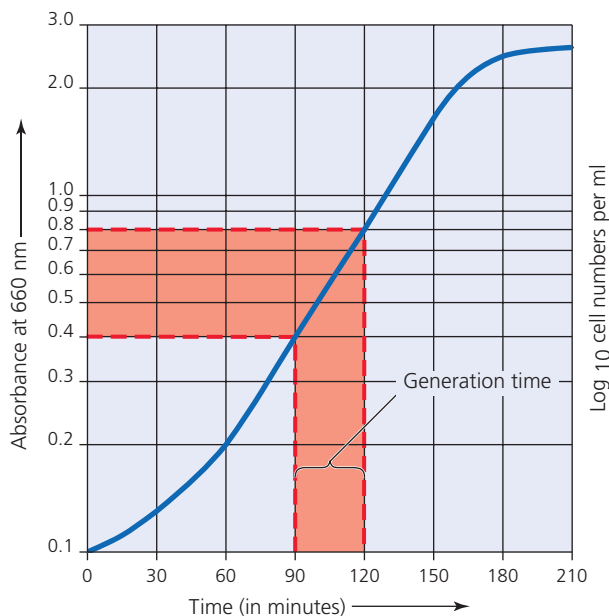

**Figure 21.2** Indirect method of determining generation time

Construction of a complete bacterial growth curve requires that aliquots of a 24-hour shake-flask culture be measured for population size at intervals during the incubation period. Such a procedure does not lend itself to a regular laboratory session. Therefore, this experiment follows a modified procedure designed to demonstrate only the lag and log phases. The curve will be plotted on semilog paper by using two values for the measurement of growth. The direct method requires enumeration of viable cells in serially diluted samples of the test culture taken at 30-minute intervals as described in Experiment 20. The indirect method uses spectrophotometric measurement of the developing turbidity at the same 30-minute intervals, as an index of increasing cellular mass.

You will determine generation time with indirect and direct methods by using data on the growth curve. Indirect determination is made by simple extrapolation from the log phase as illustrated in **Figure 21.2**. Select two points on the absorbance scale, such as 0.2 and 0.4, that represent a doubling of turbidity. Using a ruler, extrapolate by drawing a line between each of the selected absorbances on the ordinate (y-axis) and the plotted line of the growth curve. Then draw perpendicular lines from these endpoints on the plotted line of the growth curve to their respective time intervals on the abscissa (x-axis). With this information, determine the generation time (GT) as follows:

$$GT = t_{(A,0.4)} - t_{(A,0.2)}$$

$$GT = 90 \text{ minutes} - 60 \text{ minutes} = 30 \text{ minutes}$$

The generation time may be calculated directly using the log of cell numbers scale on a growth curve. The following example uses information from a hypothetical growth curve to calculate the generation time directly.

$C_0$  = number of cells at time zero

$C_E$  = number of cells at end of a specified time ( $t$ )

$N$  = number of generations (doublings)

To describe logarithmic growth, the following equation is used:

$$N = (\log C_E - \log C_0) / \log 2$$

Using this formula, the logarithmic tables to the base 10, and the following supplied information, we may now solve for the generation time:

$$C_E = 52,000,000 \text{ cells} \quad \log C_E = 7.7218$$

$$C_0 = 25,000 \text{ cells} \quad \log C_0 = 4.4048$$

$$\log 2 = 0.301$$

$$N = (7.7218 - 4.4048) / 0.301 = 11 \text{ generations}$$

$$\text{generation time (GT)} = \frac{\text{the specified time (} t \text{)}}{\text{number of generations (} N \text{)}}$$

$$t = 180 \text{ minutes}$$

$$GT = 180 / 11 = 16 \text{ minutes}$$

## CLINICAL APPLICATION

### Using Growth Curves to Determine Antimicrobial Resistance

In medical laboratories, growth curves are being mathematically modeled to quickly determine antimicrobial susceptibility. By monitoring turbidity in a series of wells, each containing a test bacterium and a dilution of an antimicrobial agent, the entire growth curve of the bacterium can be determined from early measurements, greatly speeding up the testing process for drugs as well as assessing newly resistant bacteria.

## AT THE BENCH



## Materials Cultures

5- to 10-hour (log phase) brain heart infusion broth culture of *Escherichia coli* with A of 0.08–0.10 at 600 nm.

## Media


Per designated student group: 100 ml of brain heart infusion in a 250-ml Erlenmeyer flask; eighteen 99-ml sterile water blanks; and four 100-ml bottles of nutrient agar.

## Equipment

37°C waterbath shaker incubator, Bausch & Lomb Spectronic 20 spectrophotometer, 13 × 100-mm cuvettes, Quebec colony counter, 24 sterile Petri dishes, 1-ml and 10-ml sterile pipettes, mechanical pipetting device, glassware marking pencil, 1000-ml beaker, and Bunsen burner.

## Procedure Lab One

1. Separate the eighteen 99-ml sterile water blanks into six sets of three water blanks each. Label each set as to time of inoculation ( $t_0$ ,  $t_{30}$ ,  $t_{60}$ ,  $t_{90}$ ,  $t_{120}$ ,  $t_{150}$ ) and the dilution to be effected in each water blank ( $10^{-2}$ ,  $10^{-4}$ ,  $10^{-6}$ ).
2. Label six sets of four Petri dishes as to time of inoculation and dilution to be plated ( $10^{-4}$ ,  $10^{-5}$ ,  $10^{-6}$ ,  $10^{-7}$ ).
3. Liquefy the four bottles of nutrient agar in an autoclave. Cool and maintain at 45°C.
4. With a sterile pipette, add approximately 5 ml of the log phase *E. coli* culture to the flask containing 100 ml of brain heart infusion broth. The approximate initial A ( $t_0$ ) should be 0.08 to 0.1 at 600 nm. Refer to Experiment 13 for proper use of the spectrophotometer.
5. After the  $t_0$  A has been determined, shake the culture flask and aseptically transfer 1 ml to the 99-ml water blank labeled  $t_0$   $10^{-2}$  and continue to dilute serially to  $10^{-4}$  and  $10^{-6}$ . *Note: A new pipette must be used for each subsequent dilution.*
6. Place the culture flask in a waterbath shaker set at 120 rpm at 37°C, and time for the required 30-minute intervals.
7. Shake the  $t_0$  dilution bottle as illustrated in **Figure 21.3**. Plate the  $t_0$  dilutions on the appropriately labeled  $t_0$  plates as shown in **Figure 21.4**. Aseptically pour 15 ml of the molten agar into each plate and mix by gentle rotation.

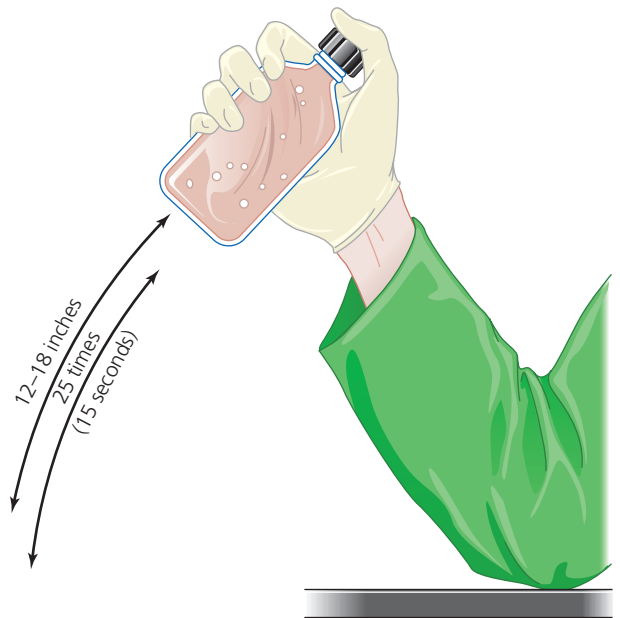

**Figure 21.3** Method for mixing sample in a dilution bottle


8. Thereafter, at each 30-minute interval, shake and aseptically transfer a 5-ml aliquot of the culture to a cuvette and determine its absorbance. Also, aseptically transfer a 1-ml aliquot of the culture into the  $10^{-2}$  water blank of the set labeled with the appropriate time, complete the serial dilution, and plate in the respectively labeled Petri dishes. *Note: A new pipette must be used for each subsequent dilution.*
9. When the pour-plate cultures harden, incubate them in an inverted position for 24 hours at 37°C.

## Procedure Lab Two

1. Perform cell counts on all plates as described in Experiment 20. Cell counts are often referred to as colony-forming units (CFUs) because each single cell in the plate becomes visible as a colony, which can then be counted.
2. Record your results in the Lab Report.



## PROCEDURE

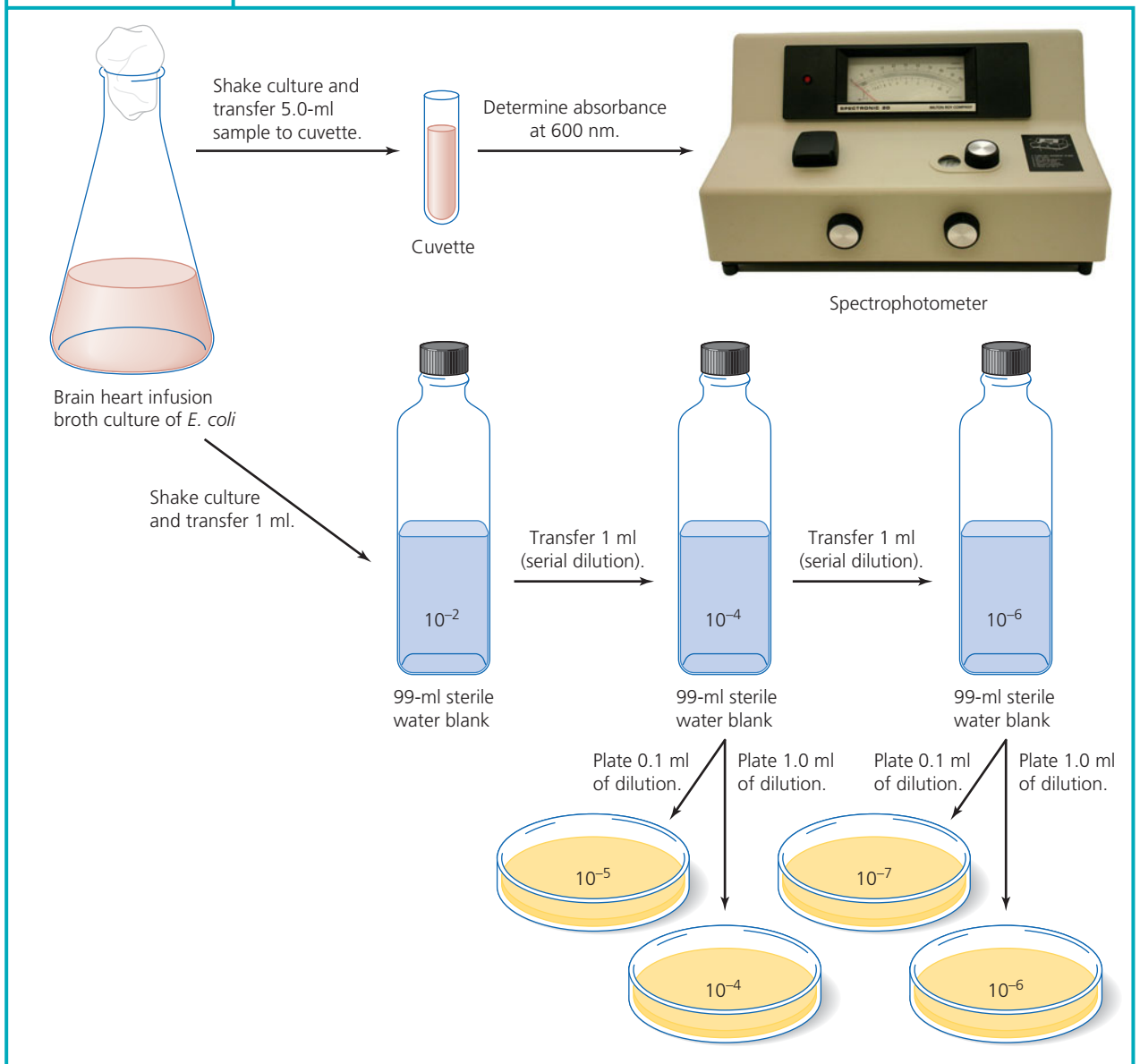

**Figure 21.4** Spectrophotometric and dilution-plateing procedure for use in bacterial growth curves

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results



- Record the absorbances and corresponding cell counts in the chart below.

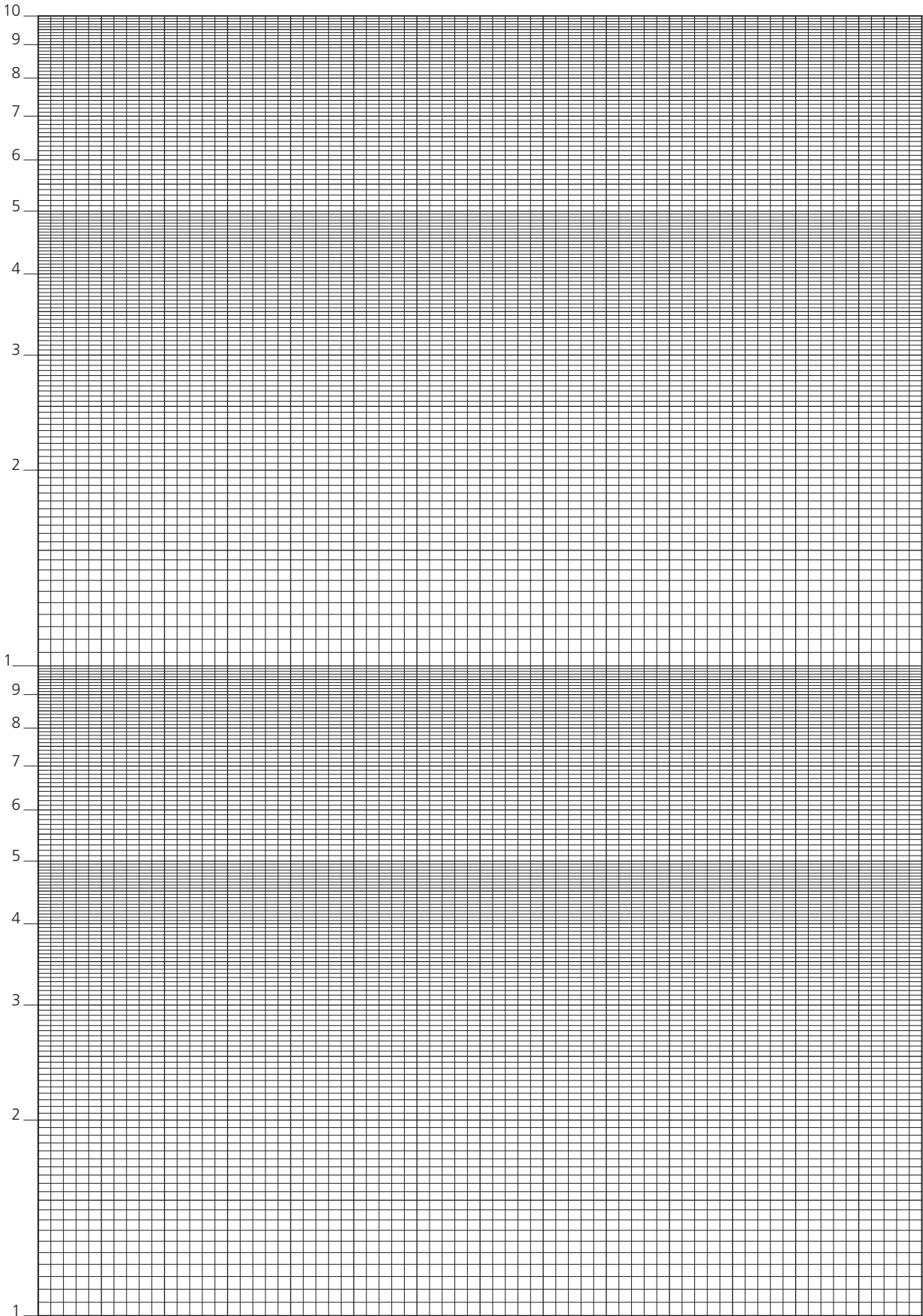
Incubation Time (minutes)	Absorbance at 600 nm	Plate Counts (CFU/ml)	Log of CFU/ml
0			
30			
60			
90			
120			
150			

- On the semilog paper provided on pages 151 and 152:
  - Plot a curve relating the absorbances on the ordinate versus incubation time on the abscissa as shown in Figure 21.2.
  - Plot a population curve with the log of the viable cells/ml on the ordinate and the incubation time on the abscissa. On both graphs, use a ruler to draw the best line connecting the plotted points. The straight-line portion of the curve represents the log phase.
- Calculate the generation time for this culture by the direct method (using the mathematical formula) and by the indirect method (extrapolating from the A scale on the plotted curve). Show calculations, and record the generation time.
  - Direct method:


**b.** Indirect method:

## Review Questions

1. Does the term *growth* convey the same meaning when applied to bacteria and to multicellular organisms? Explain.
  
  
  
  
  
  
  
  
  
  
2. Why do variations in generation time exist:
  - a. Among different species of microorganisms?
  
  
  
  
  
  
  
  
  
  
  - b. Within a single microbial species?
  
  
  
  
  
  
  
  
  
  
3.  The generation time and growth rate of an organism grown in the laboratory can be easily determined by constructing a typical growth curve.
  - a. Would you expect the growth rate of the infectious organisms found in an abscess that developed from a wound to mimic the growth curve obtained in the laboratory? Explain.
  
  
  
  
  
  
  
  
  
  
  - b. Would you expect antibiotic therapy to be effective without any other concurrent treatment of the abscess?
  
  
  
  
  
  
  
  
  
  
4.  Is generation time a useful parameter to indicate the types of media best suited to support the growth of a specific organism? Explain.



Absorbance versus Incubation Time

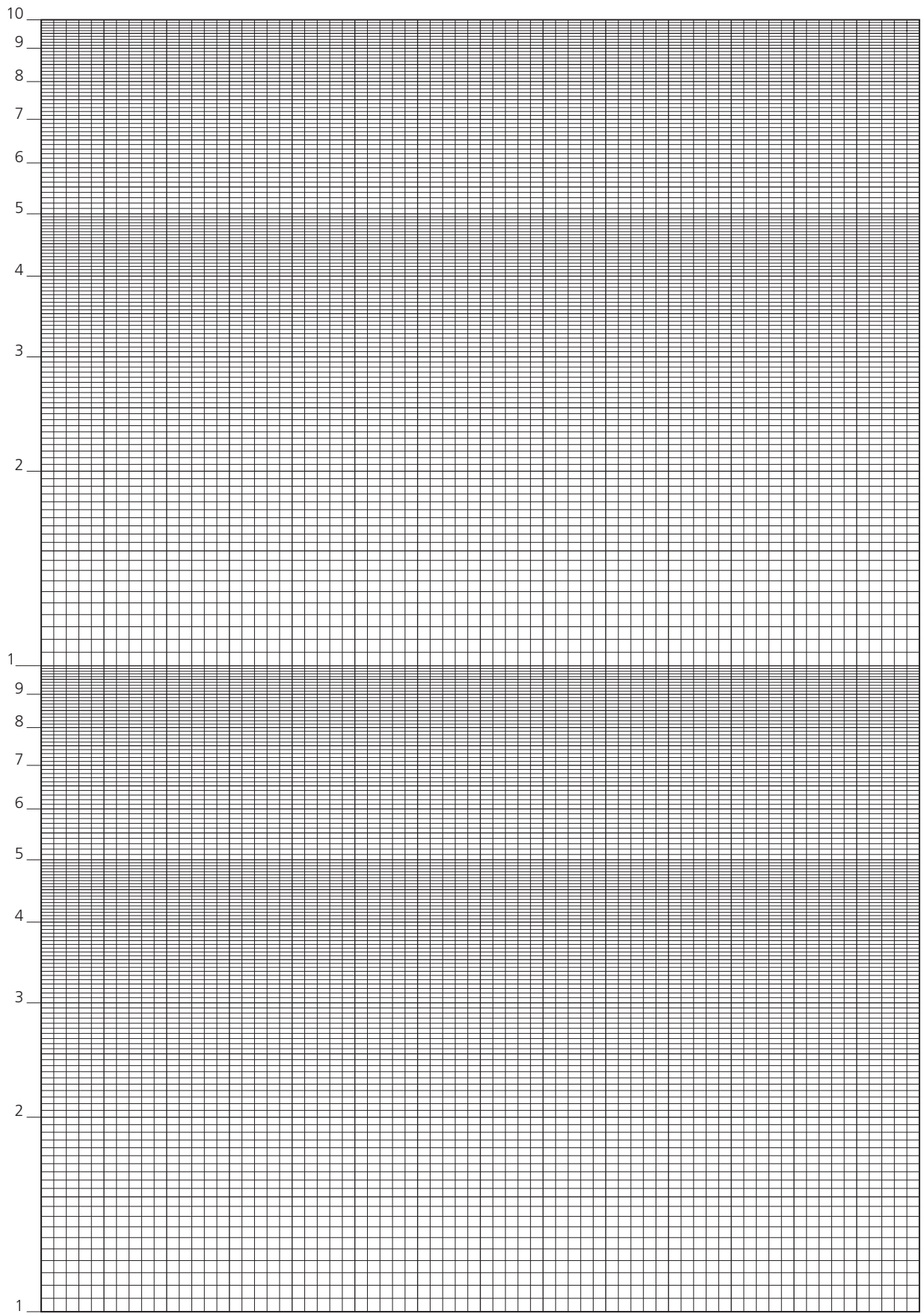

Colony-Forming Units (CFUs) versus Incubation Time

# Biochemical Activities of Microorganisms

## LEARNING OBJECTIVES

*When you have completed the experiments in this section, you should understand*

1. The nature and activities of exoenzymes and endoenzymes.
2. Experimental procedures for differentiation of enteric microorganisms.
3. Biochemical test procedures for identification of microorganisms.

## Introduction

Microorganisms must be separated and identified for a wide variety of reasons, such as

1. Determination of pathogens responsible for infectious diseases.
2. Selection and isolation of strains of fermentative microorganisms necessary for the industrial production of alcohols, solvents, vitamins, organic acids, antibiotics, and industrial enzymes.
3. Isolation and development of suitable microbial strains necessary for the manufacture and the enhancement of quality and flavor in certain food materials such as yogurt, cheeses, and other milk products.
4. Comparison of biochemical activities for taxonomic purposes.

To accomplish these tasks, the microbiologist is assisted by the fact that, just as human beings possess a characteristic and specific set of fingerprints, microorganisms all have their own identifying biochemical characteristics. These so-called biochemical fingerprints are the properties controlled by the cells' enzymatic activity, and they are responsible for bioenergetics, biosynthesis, and biodegradation.

The sum of all these chemical reactions is defined as **cellular metabolism**, and the


biochemical transformations that occur both outside and inside the cell are governed by biological catalysts called **enzymes**.

## Extracellular Enzymes (Exoenzymes)

Exoenzymes act on substances outside of the cell. Most high-molecular-weight substances are not able to pass through cell membranes, and therefore these raw materials—foodstuffs such as polysaccharides, lipids, and proteins—must be degraded to low-molecular-weight materials—nutrients—before they can be transported into the cell. Because of the reactions involved, exoenzymes are mainly **hydrolytic enzymes** that reduce high-molecular-weight materials into their building blocks by introducing water into the molecule. This liberates smaller molecules, which may then be transported into the cell and assimilated.

## Intracellular Enzymes (Endoenzymes)

Endoenzymes function inside the cell and are mainly responsible for synthesis of new protoplasmic requirements and production of cellular energy from assimilated materials. The ability of cells to act on nutritional substrates permeating cell membranes indicates the presence of many endoenzymes capable of transforming the chemically specific substrates into essential materials.

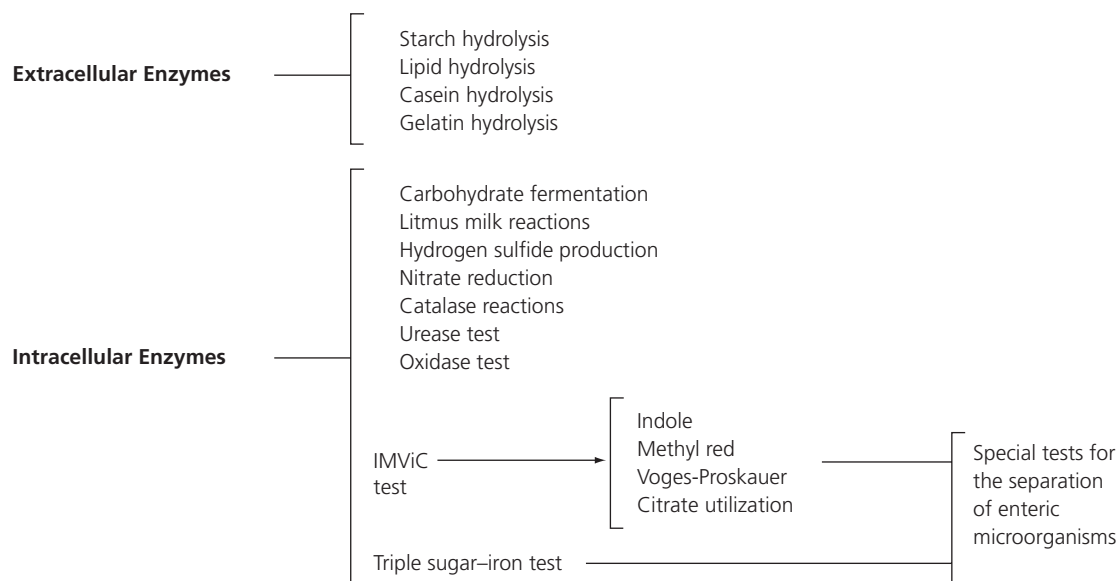

**Figure P5.1** Biochemical activities of microorganisms

This transformation is necessary for cellular survival and function, and it is the basis of cellular metabolism. As a result of these metabolic processes, metabolic products are formed and excreted by the cell into the environment. Assay of these end products not only aids in identification of specific enzyme systems but also serves to identify, separate, and classify microorganisms. **Figure P5.1** represents a simplified schema of experimental procedures used to acquaint students with the intracellular and extracellular enzymatic activities of microorganisms.

The experiments you will carry out in this section can be performed in either of two ways. A short version uses a limited number of organisms to illustrate the possible end product(s) that may result from enzyme action on a substrate. The organisms for this version are designated in the individual exercises.

The alternative, or long, version involves the use of 13 microorganisms. This version provides a complete overview of the biochemical fingerprints of the organisms and supplies the format for their

separation and identification. These organisms were chosen to serve as a basis for identification of an unknown microorganism in Experiment 33. If this alternative version is selected, the following organisms are recommended for use:

- Escherichia coli*
- Enterobacter aerogenes*
- Klebsiella pneumoniae*
- Shigella dysenteriae*
- Salmonella typhimurium*
- Proteus vulgaris*
- Pseudomonas aeruginosa*
- Alcaligenes faecalis*
- Micrococcus luteus*
- Lactococcus lactis*
- Staphylococcus aureus*
- Bacillus cereus*
- Corynebacterium xerosis*

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be able to

1. Understand the function of microbial extracellular enzymes.
2. Determine the ability of microorganisms to excrete hydrolytic extracellular enzymes capable of degrading the polysaccharide starch, the lipid tributyrin, and the proteins casein and gelatin.

## Principle

Because of their large sizes, high-molecular-weight nutrients such as polysaccharides, lipids, and proteins are not capable of permeating the cell membrane. These macromolecules must first be hydrolyzed by specific extracellular enzymes into their respective basic building blocks. These low-molecular-weight substances can then be transported into the cells and used for the synthesis of protoplasmic requirements and energy production. The following procedures are designed to investigate the exoenzymatic activities of different microorganisms.

### Starch Hydrolysis

**Starch** is a high-molecular-weight, branching polymer composed of **glucose** molecules linked together by **glycosidic bonds**. The degradation of this macromolecule first requires the presence of the extracellular enzyme **amylase** for its hydrolysis into shorter polysaccharides, namely **dextrins**, and ultimately into **maltose** molecules. The final hydrolysis of this disaccharide, which is catalyzed by **maltase**, yields low-molecular-weight, soluble **glucose** molecules that can be transported into the cell and used for energy production through the process of glycolysis.

In this experimental procedure, starch agar is used to demonstrate the hydrolytic activities




**Figure 22.1** Starch agar plate. Starch hydrolysis on left; no starch hydrolysis on right.

of these exoenzymes. The medium is composed of nutrient agar supplemented with starch, which serves as the polysaccharide substrate. The detection of the hydrolytic activity following the growth period is made by performing the starch test to determine the presence or absence of starch in the medium. Starch in the presence of iodine will impart a blue-black color to the medium, indicating the absence of starch-splitting enzymes and representing a negative result. If the starch has been hydrolyzed, a clear zone of hydrolysis will surround the growth of the organism. This is a positive result. Positive and negative results are shown in **Figure 22.1**.

### Lipid Hydrolysis

**Lipids** are high-molecular-weight compounds possessing large amounts of energy. The degradation of lipids such as **triglycerides** is accomplished by extracellular hydrolyzing enzymes, called **lipases** (esterases), that cleave the **ester bonds** in this molecule by the addition of water to form the building blocks **glycerol** (an alcohol) and **fatty acids**. **Figure 22.2** shows this reaction. Once assimilated into the cell, these basic components can be further metabolized through aerobic respiration to produce cellular energy, adenosine triphosphate (ATP). The



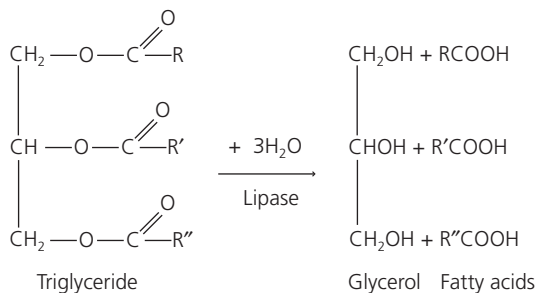

**Figure 22.2 Lipid hydrolysis**

components may also enter other metabolic pathways for the synthesis of other cellular protoplasmic requirements.

In this experimental procedure, tributyrin agar is used to demonstrate the hydrolytic activities of the exoenzyme lipase. The medium is composed of nutrient agar supplemented with the triglyceride tributyrin as the lipid substrate. Tributyrin forms an emulsion when dispersed in the agar, producing an opaque medium that is necessary for observing exoenzymatic activity.

Following inoculation and incubation of the agar plate cultures, organisms excreting lipase will show a zone of **lipolysis**, which is demonstrated by a clear area surrounding the bacterial growth. This loss of opacity is the result of the hydrolytic reaction yielding soluble glycerol and fatty acids and represents a positive reaction for lipid hydrolysis. In the absence of lipolytic enzymes, the medium retains its opacity. This is a negative reaction. Positive and negative results are shown in **Figure 22.3**.


## Casein Hydrolysis

**Casein**, the major milk protein, is a macromolecule composed of **amino acid** subunits linked together by **peptide bonds** (CO—NH). Before their assimilation into the cell, proteins must undergo step-by-step degradation into **peptones**, **polypeptides**, **dipeptides**, and ultimately into their building blocks, **amino acids**. This process is called **peptonization**, or **proteolysis**, and it is mediated by extracellular enzymes called **proteases**. The function of these proteases is to cleave the peptide bond CO—NH by introducing water into the molecule. The reaction then liberates the amino acids, as illustrated in **Figure 22.4**.

The low-molecular-weight soluble amino acids can now be transported through the cell membrane into the intracellular amino acid pool for use in the synthesis of structural and functional cellular proteins.



**Figure 22.3 Tributyrin agar plate.** Lipid hydrolysis on left; no lipid hydrolysis on right.

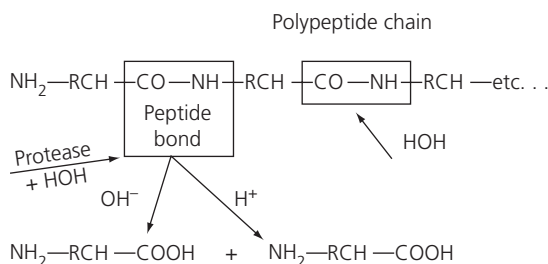

**Figure 22.4 Protein hydrolysis**

In this experimental procedure, milk agar is used to demonstrate the hydrolytic activity of these exoenzymes. The medium is composed of nutrient agar supplemented with milk that contains the protein substrate casein. Similar to other proteins, milk protein is a colloidal suspension that gives the medium its color and opacity because it deflects light rays rather than transmitting them.

Following inoculation and incubation of the agar plate cultures, organisms secreting proteases will exhibit a zone of **proteolysis**, which is demonstrated by a clear area surrounding the bacterial growth. This loss of opacity is the result of a hydrolytic reaction yielding soluble, noncolloidal amino acids, and it represents a positive reaction. In the absence of protease activity, the medium surrounding the growth of the organism remains opaque, which is a negative reaction.

## Gelatin Hydrolysis

Although the value of gelatin as a nutritional source is questionable (it is an incomplete protein, lacking the essential amino acid tryptophan), its value in identifying bacterial species is well established. Gelatin is a protein produced by hydrolysis

of collagen, a major component of connective tissue and tendons in humans and other animals. Below temperatures of 25°C, gelatin will maintain its gel properties and exist as a solid; at temperatures above 25°C, gelatin is liquid. **Figure 22.5** shows gelatin hydrolysis.

Liquefaction is accomplished by some microorganisms capable of producing a proteolytic extracellular enzyme called **gelatinase**, which acts to hydrolyze this protein to **amino acids**. Once this degradation occurs, even very low temperatures of 4°C will not restore the gel characteristic.

In this experimental procedure, you will use nutrient gelatin deep tubes to demonstrate the hydrolytic activity of gelatinase. The medium consists of nutrient broth supplemented with 12% gelatin. This high gelatin concentration results in a stiff medium and also serves as the substrate for the activity of gelatinase.

Following inoculation and incubation for 48 hours, the cultures are placed in a refrigerator at 4°C for 30 minutes. Cultures that remain liquefied produce gelatinase and demonstrate *rapid* gelatin hydrolysis. Re-incubate all solidified cultures for an additional 5 days. Refrigerate for 30 minutes and observe for liquefaction. Cultures that remain liquefied are indicative of *slow* gelatin hydrolysis.

## CLINICAL APPLICATION

### Pathogens and Extracellular Enzymes

Bacteria use enzymes to alter their environments and to gain new sources of nutrients. When known bacterial pathogens are causing symptoms or damage not normally associated with that species, laboratories may test for newly acquired extracellular enzymes. Most known pathogens have been characterized by their abilities to digest proteins (fibronectin and collagen) as well as lipids and starches (glycolipids and glycoproteins).

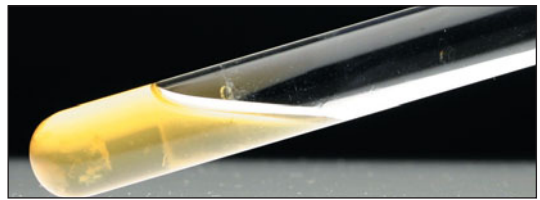
## AT THE BENCH




## Materials

### Cultures

24- to 48-hour trypticase soy broth cultures of *Escherichia coli*, *Bacillus cereus*, *Pseudomonas aeruginosa*, and *Staphylococcus aureus* for the



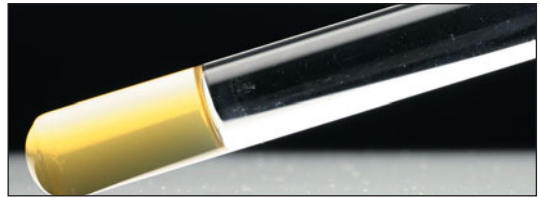(a) Positive for gelatin liquefaction



(b) Negative for gelatin liquefaction

**Figure 22.5** Nutrient gelatin hydrolysis

short version. 24- to 48-hour brain heart infusion broth cultures of the 13 organisms listed on page 154 for the long version.

## Media

Short version: Two plates each of starch agar, tributyrin agar, and milk agar, and three nutrient gelatin deep tubes per designated student group. Long version: Four plates each of starch agar, tributyrin agar, and milk agar, and 14 nutrient gelatin deep tubes per designated student group.

## Reagent

Gram's iodine solution.

## Equipment

Bunsen burner, inoculating loop and needle, glassware marking pencil, test tube rack, and refrigerator.

## Procedure Lab One

1. Prepare the starch agar, tributyrin agar, and milk agar plates for inoculation as follows:
  - a. Short procedure: Using two plates per medium, divide the bottom of each Petri dish into two sections. Label the sections as *E. coli*, *B. cereus*, *P. aeruginosa*, and *S. aureus*, respectively.
  - b. Long procedure: Repeat Step 1a, dividing three plate bottoms into three sections and one plate bottom into four sections for each of the required media, to accommodate the 13 test organisms.

2. Using aseptic technique, make a single-line streak inoculation of each test organism on the agar surface of its appropriately labeled section on the agar plates.
3. Using aseptic technique, inoculate each experimental organism in its appropriately labeled gelatin deep tube by means of a stab inoculation.
4. Incubate all plates in an inverted position for 24 to 48 hours at 37°C. Incubate the gelatin deep tube cultures for 48 hours. Re-incubate all negative cultures for an additional 5 days.

## Procedure Lab Two

### Starch Hydrolysis

1. Flood the starch agar plate cultures with Gram's iodine solution, allow the iodine to remain in contact with the medium for 30 seconds, and pour off the excess.
2. Examine the cultures for the presence or absence of a blue-black color surrounding the growth of each test organism. Record your results in the chart provided in the Lab Report.
3. Based on your observations, determine and record which organisms were capable of hydrolyzing the starch.

### Lipid Hydrolysis

1. Examine the tributyrin agar plate cultures for the presence or absence of a clear area, or zone of lipolysis, surrounding the growth of

each of the organisms. Record your results in the chart provided in the Lab Report.

2. Based on your observations, determine and record which organisms were capable of hydrolyzing the lipid.

### Casein Hydrolysis

1. Examine the milk agar plate cultures for the presence or absence of a clear area, or zone of proteolysis, surrounding the growth of each of the bacterial test organisms. Record your results in the chart provided in the Lab Report.
2. Based on your observations, determine and record which of the organisms were capable of hydrolyzing the milk protein casein.

### Gelatin Hydrolysis

1. Place all gelatin deep tube cultures into a refrigerator at 4°C for 30 minutes.
2. Examine all the cultures to determine whether the medium is solid or liquid. Record your results in the chart provided in the Lab Report.
3. Based on your observations following the 2-day and 7-day incubation periods, determine and record in the Lab Report (a) which organisms were capable of hydrolyzing gelatin and (b) the rate of hydrolysis.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

### Starch and Lipid Hydrolysis

Bacterial Species	STARCH HYDROLYSIS		LIPID HYDROLYSIS	
	Appearance of Medium	Result (+) or (-)	Appearance of Medium	Result (+) or (-)
<i>E. coli</i>				
<i>E. aerogenes</i>				
<i>K. pneumoniae</i>				
<i>S. dysenteriae</i>				
<i>S. typhimurium</i>				
<i>P. vulgaris</i>				
<i>P. aeruginosa</i>				
<i>A. faecalis</i>				
<i>M. luteus</i>				
<i>L. lactis</i>				
<i>S. aureus</i>				
<i>B. cereus</i>				
<i>C. xerosis</i>				

## Casein and Gelatin Hydrolysis

Bacterial Species	CASEIN HYDROLYSIS		GELATIN HYDROLYSIS		
	Appearance of Medium	Result (+) or (-)	Liquefaction (+) or (-)		Rate of Hydrolysis (Slow or Rapid)
			2 days	7 days	
<i>E. coli</i>					
<i>E. aerogenes</i>					
<i>K. pneumoniae</i>					
<i>S. dysenteriae</i>					
<i>S. typhimurium</i>					
<i>P. vulgaris</i>					
<i>P. aeruginosa</i>					
<i>A. faecalis</i>					
<i>M. luteus</i>					
<i>L. lactis</i>					
<i>S. aureus</i>					
<i>B. cereus</i>					
<i>C. xerosis</i>					

## Review Questions

1. Why is the catalytic activity of enzymes essential to ensure and regulate cellular metabolism?
2. Why are microorganisms able to cause dairy products, such as milk, to sour or curdle?
3. Give a reason why it is necessary for polysaccharides, such as starch or cellulose, to be digested outside of the cell even though disaccharides, such as lactose or sucrose, are digestible inside the cell.

## LEARNING OBJECTIVES

Once you have completed this experiment, you should

1. Understand the difference between cellular respiration and fermentation.
2. Be able to determine the ability of microorganisms to degrade and ferment carbohydrates with the production of acid and gas.


## Principle

Most microorganisms obtain their energy through a series of orderly and integrated enzymatic reactions leading to the biooxidation of a substrate, frequently a carbohydrate. The major pathways by which this is accomplished are shown in **Figure 23.1**.

Organisms use carbohydrates differently depending on their enzyme complement. Some organisms are capable of fermenting sugars such as glucose anaerobically, while others use the aerobic pathway. Still others, facultative anaerobes, are enzymatically competent to use both aerobic and anaerobic pathways, and some organisms lack the ability to oxidize glucose by either. In this exercise the fermentative pathways are of prime concern.

In fermentation, substrates such as carbohydrates and alcohols undergo anaerobic dissimilation and produce an organic acid (for example, lactic, formic, or acetic acid) that may be accompanied by gases such as hydrogen or carbon dioxide. Facultative anaerobes are usually the so-called fermenters of carbohydrates. Fermentation is best described by considering the degradation of glucose by way of the **Embden-Meyerhof pathway**, also known as the **glycolytic pathway**, illustrated in **Figure 23.2**.

As the diagram shows, one mole of glucose is converted into two moles of pyruvic acid, which is the major intermediate compound produced by glucose degradation. Subsequent metabolism of

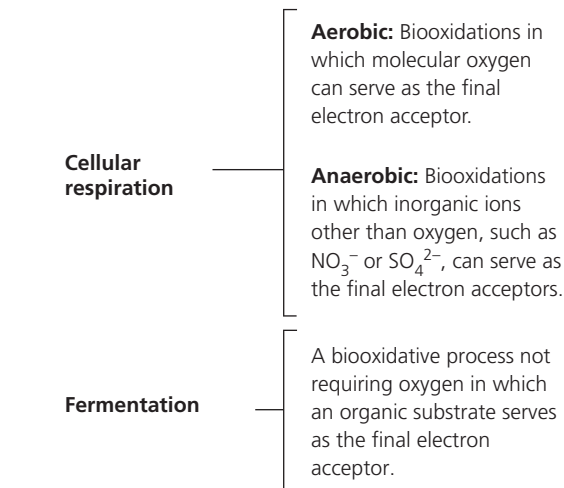


**Figure 23.1** Biooxidative pathways

pyruvate is not the same for all organisms, and a variety of end products result that define their different fermentative capabilities. This can be seen in **Figure 23.3**.


Fermentative degradation under anaerobic conditions is carried out in a fermentation broth tube containing a Durham tube, an inverted inner vial for the detection of gas production as illustrated in **Figure 23.4**. A typical carbohydrate fermentation medium contains

1. Nutrient broth ingredients for the support of the growth of all organisms.
2. A specific carbohydrate that serves as the substrate for determining the organism's fermentative capabilities.
3. The pH indicator phenol red, which is red at a neutral pH (7) and changes to yellow at a slightly acidic pH of 6.8, indicating that slight amounts of acid will cause a color change.

The critical nature of the fermentation reaction and the activity of the indicator make it imperative that all cultures should be observed within 48 hours. Extended incubation may mask acid-producing reactions by production of alkali because of enzymatic action on substrates other than the carbohydrate.

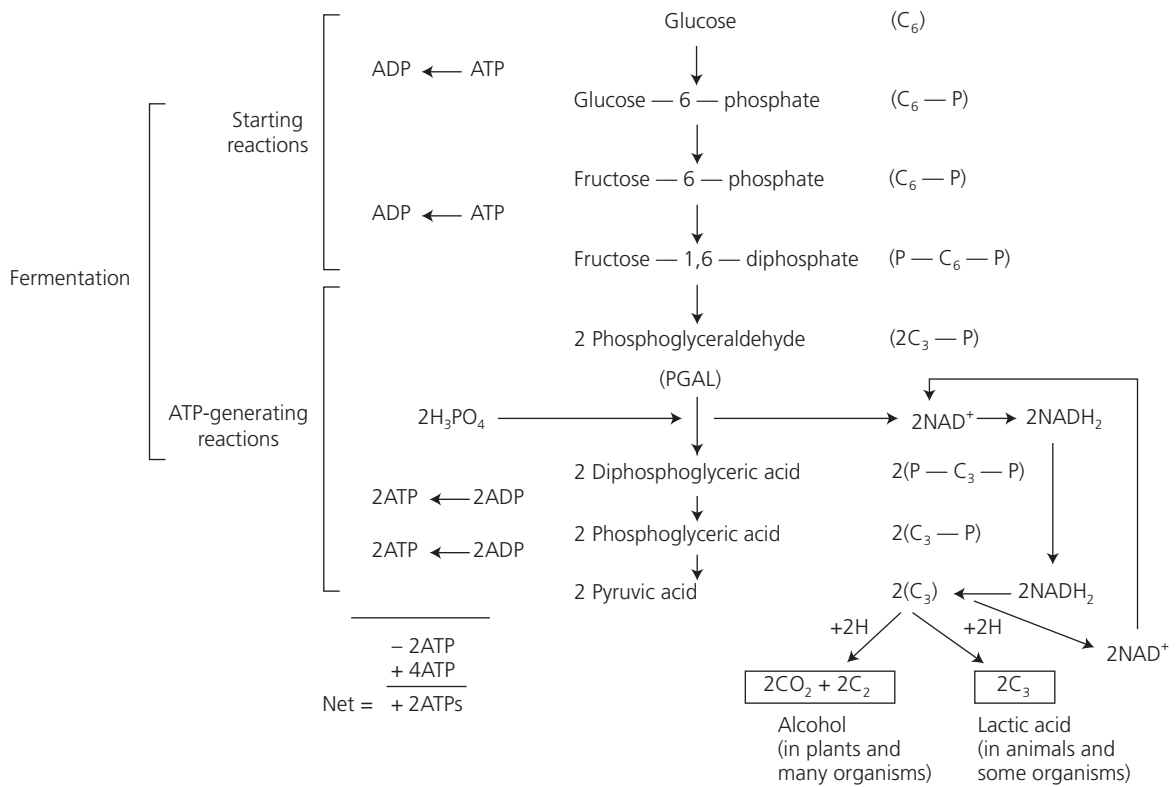


**Figure 23.2** The Embden-Meyerhof pathway

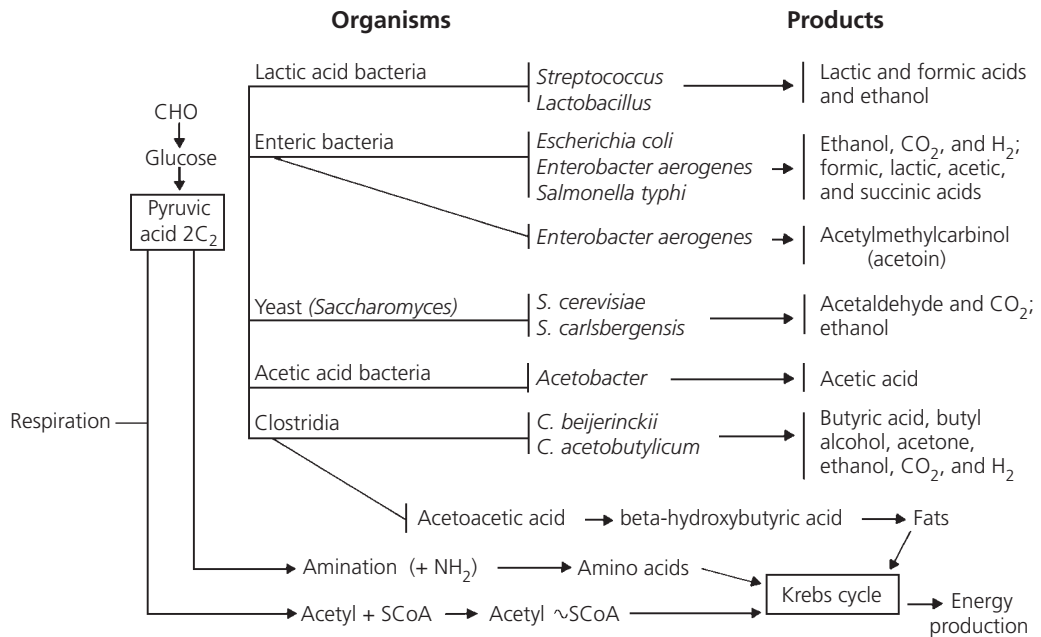

**Figure 23.3** Variations in the use of pyruvic acid

Following incubation, carbohydrates that have been fermented with the production of acidic wastes will cause the phenol red (Figure 23.5a) to turn yellow, thereby indicating a positive reaction (Figures 23.5b and c). In

some cases, acid production is accompanied by the evolution of a gas ( $CO_2$ ) that will be visible as a bubble in the inverted tube (Figure 23.5b). Cultures that are not capable of fermenting a carbohydrate substrate will not change the

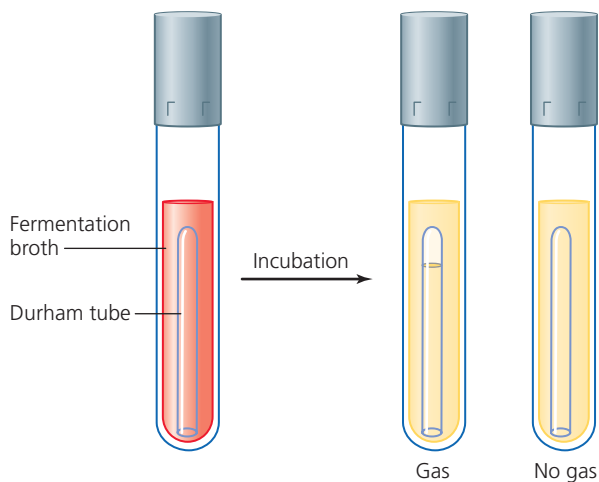

**Figure 23.4** Detection of gas production


indicator, and the tubes will appear red; there will not be a concomitant evolution of gas. This is a negative reaction (**Figure 23.5d**).

The lack of carbohydrate fermentation by some organisms should not be construed as absence of growth. The organisms use other nutrients in the medium as energy sources. Among these nutrients are peptones present in nutrient broth. Peptones can be degraded by microbial enzymes to amino acids that are in turn enzymatically converted by oxidative deamination to ketoamino acids. These are then metabolized through the Krebs cycle for energy production. These reactions liberate ammonia, which accumulates in the medium, forming ammonium hydroxide ( $\text{NH}_4\text{OH}$ ) and producing an alkaline environment. When this occurs, the phenol red turns to a deep red in the now basic medium. This alternative pathway of aerobic respiration is illustrated in **Figure 23.6**.

## CLINICAL APPLICATION

### Using Fermentation Products to Identify Bacteria

The fermentation of carbohydrates assists in the identification of some bacteria by determining what nutrients they are using and what products they produce. The pattern of sugars fermented may be unique to a particular genus, species, or strain. Lactose fermentation is one test that distinguishes between enteric and non-enteric bacteria. Dextrose fermentation allows for the differentiation between the oxidase (+) *Vibrio* and *Pseudomonads* species in patients suffering from septicemia after eating contaminated fish.

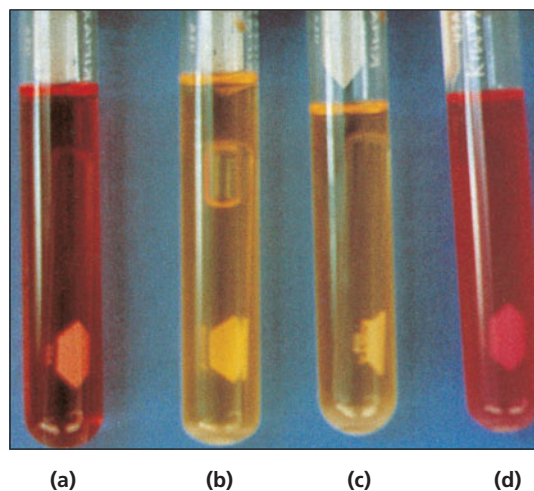

**Figure 23.5** Carbohydrate fermentation test. (a) Uninoculated, (b) acid and gas, (c) acid, and (d) negative.

## AT THE BENCH



## Materials

### Cultures

24- to 48-hour Trypticase soy broth cultures of *Escherichia coli*, *Alcaligenes faecalis*, *Salmonella typhimurium*, and *Staphylococcus aureus* for the short version. 24- to 48-hour Trypticase soy broth cultures of the 13 organisms listed on page 154 for the long version.

### Media

Per designated student group: phenol red lactose, dextrose (glucose), and sucrose broths: 5 of each for the short version, 14 of each for the long version.


### Equipment

Bunsen burner, inoculating loop, and glassware marking pencil.

## Procedure Lab One

- Using aseptic technique, inoculate each experimental organism into its appropriately labeled medium by means of loop inoculation. *Note:* Take care during this step not to shake the fermentation tube; shaking the tube may accidentally force a bubble of air into the inverted



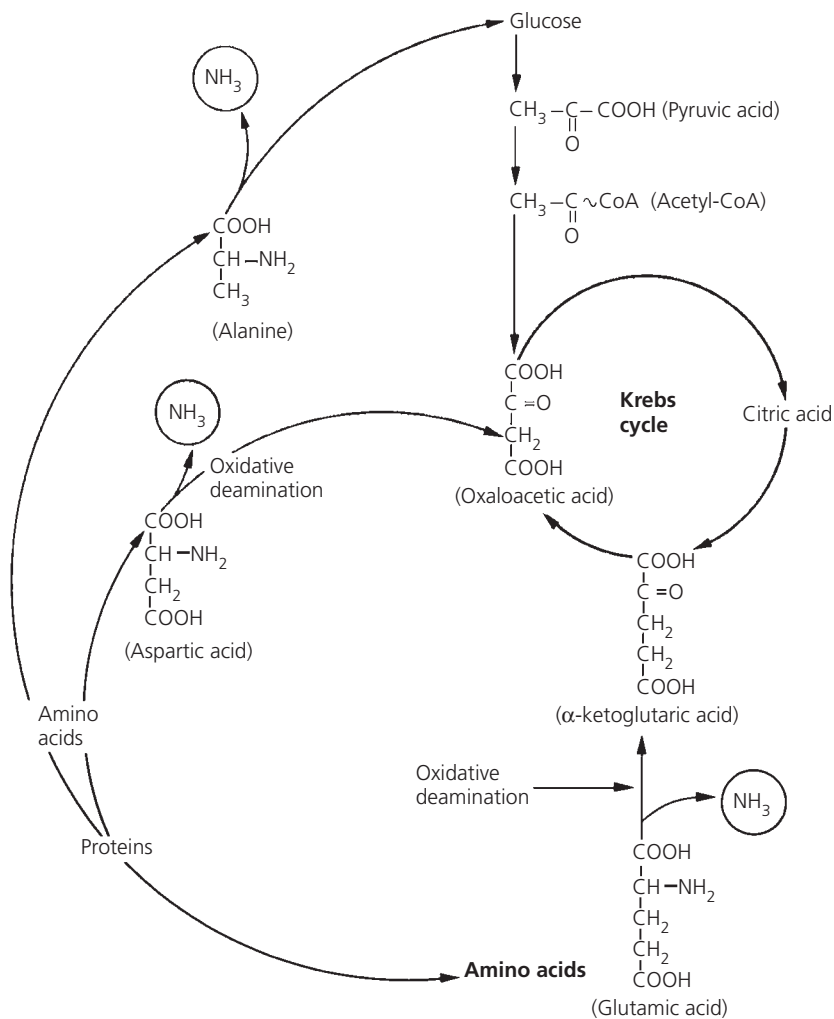

**Figure 23.6** Proteins as energy sources for microbes

gas vial, displacing the medium and possibly rendering a false-positive result. The last tube will serve as a control.

2. Incubate all tubes for 24 hours at 37°C.

2. Based on your observations, determine and record whether or not each organism was capable of fermenting the carbohydrate substrate with the production of acid or acid and gas.

## Procedure Lab Two

1. Examine all carbohydrate broth cultures for color and the presence or absence of a gas bubble. Record your results in the chart provided in the Lab Report.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

Bacterial Species	Lactose Observation (color of medium, bubble in fermentation tube)	Result (A), (A/G), or (–)	Dextrose Observation (color of medium, bubble in fermentation tube)	Result (A), (A/G), or (–)	Sucrose Observation (color of medium, bubble in fermentation tube)	Result (A), (A/G), or (–)
<i>E. coli</i>						
<i>E. aerogenes</i>						
<i>K. pneumoniae</i>						
<i>S. dysenteriae</i>						
<i>S. typhimurium</i>						
<i>P. vulgaris</i>						
<i>P. aeruginosa</i>						
<i>A. faecalis</i>						
<i>M. luteus</i>						
<i>L. lactis</i>						
<i>S. aureus</i>						
<i>B. cereus</i>						
<i>C. xerosis</i>						
Control						

## Review Questions

1. Distinguish between respiration and fermentation.



# Triple Sugar–Iron Agar Test

## LEARNING OBJECTIVES

*Once you have completed this experiment, you should understand a rapid screening procedure that will*

1. Differentiate among members of the Enterobacteriaceae.
2. Distinguish between the Enterobacteriaceae and other groups of intestinal bacilli.

## Principle

The **triple sugar–iron (TSI) agar test** is designed to differentiate among the different groups or genera of the Enterobacteriaceae, which are all gram-negative bacilli capable of fermenting glucose with the production of acid, and to distinguish the Enterobacteriaceae from other gram-negative intestinal bacilli. This differentiation is made on the basis of differences in carbohydrate fermentation patterns and hydrogen sulfide production by the various groups of intestinal organisms.

To facilitate observation of carbohydrate utilization patterns, the TSI agar slants contain lactose and sucrose in 1% concentrations and glucose (dextrose) in a concentration of 0.1%, which permits detection of the utilization of this substrate only. The acid-base indicator phenol red is also incorporated to detect carbohydrate fermentation that is indicated by a change in color of the medium from orange-red to yellow in the presence of acids. The slant is inoculated by means of a stab-and-streak procedure. This requires the insertion of a sterile, straight needle from the base of the slant into the butt. Upon withdrawal of the needle, the slanted surface of the medium is streaked. Following incubation, you will determine the fermentative activities of the organisms as described below.

1. **Alkaline slant (red) and acid butt (yellow) with or without gas production (breaks in the agar butt).** Only glucose fermentation has occurred. The organisms preferentially degrade glucose first. Since this

substrate is present in minimal concentration, the small amount of acid produced on the slant surface is oxidized rapidly. The peptones in the medium are also used in the production of alkali. In the butt the acid reaction is maintained because of reduced oxygen tension and slower growth of the organisms.

2. **Acid slant (yellow) and acid butt (yellow) with or without gas production.**

Lactose and/or sucrose fermentation has occurred. Since these substances are present in higher concentrations, they serve as substrates for continued fermentative activities with maintenance of an acid reaction in both slant and butt.


3. **Alkaline slant (red) and alkaline butt (red) or no change (orange-red) butt.**

No carbohydrate fermentation has occurred. Instead, peptones are catabolized under anaerobic and/or aerobic conditions, resulting in an alkaline pH due to production of ammonia. If only aerobic degradation of peptones occurs, the alkaline reaction is evidenced only on the slant surface. If there is aerobic and anaerobic utilization of peptone, the alkaline reaction is present on the slant and the butt.

For you to obtain accurate results, it is absolutely essential to observe the cultures within 18 to 24 hours following incubation. Doing so will ensure that the carbohydrate substrates have not been depleted and that degradation of peptones yielding alkaline end products has not taken place.

The TSI agar medium also contains sodium thiosulfate, a substrate for hydrogen sulfide ( $H_2S$ ) production, and ferrous sulfate for detection of this colorless end product. Following incubation, only cultures of organisms capable of producing  $H_2S$  will show an extensive blackening in the butt because of the precipitation of the insoluble ferrous sulfide. (Refer to Experiment 26 for a more detailed biochemical explanation of  $H_2S$  production.)

**Figure 24.1** is a schema for the differentiation of intestinal bacilli on the basis of the TSI agar reactions.

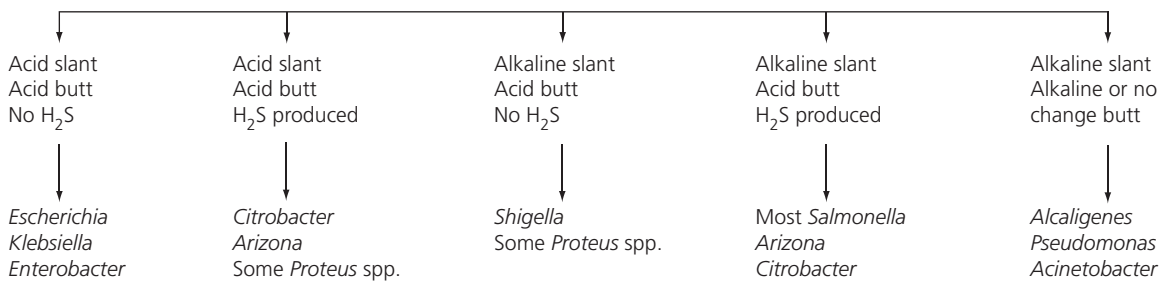

**Figure 24.1** TSI reactions for differentiation of enteric microorganisms

## CLINICAL APPLICATION

### Differentiating Between *Proteus* Species

The TSI test can differentiate enteric organisms based on their abilities to reduce sulfur and ferment carbohydrates. It can be used to separate the three species of *Proteus*, *P. vulgaris*, *P. mirabilis*, and *P. penneri*, all of which are human opportunistic pathogens. *P. mirabilis* causes urinary tract infections and is sensitive to treatment with ampicillin and cephalosporins. *P. vulgaris*, a less common cause of urinary tract infections, is not sensitive to these antibiotics and is found as a nosocomial infectious agent amongst immunocompromised patients.

## AT THE BENCH



## Materials

### Cultures

24-hour Trypticase soy broth cultures of *Pseudomonas aeruginosa*, *Escherichia coli*, *Salmonella typhimurium*, *Shigella dysenteriae*, *Proteus vulgaris*, and *Alcaligenes faecalis* for the short version. 24-hour Trypticase soy broth cultures of the 13 organisms listed on page 154 for the long version.

### Media

Per designated student group: triple sugar–iron agar slants: 7 for the short version, 14 for the long version.

### Equipment

Bunsen burner, inoculating needle, test tube rack, and glassware marking pencil.

## Procedure Lab One


- Using aseptic technique, inoculate each experimental organism into its appropriately labeled tube by means of a stab-and-streak

inoculation. *Note: Do not fully tighten screw cap.* The last tube will serve as a control.

- Incubate for 18 to 24 hours at 37°C.

## Procedure Lab Two

- Examine the color of both the butt and slant of all agar slant cultures (**Figure 24.2**). Based on your observations, determine the type of reaction that has taken place (acid, alkaline, or none) and the carbohydrate that has been fermented (dextrose, lactose, and/or sucrose, all, or none) in each culture. Record your observations and results in the chart provided in the Lab Report.
- Examine all cultures for the presence or absence of blackening within the medium. Based on your observations, determine whether or not each organism was capable of H<sub>2</sub>S production. Record your observations and results in the chart provided in the Lab Report.

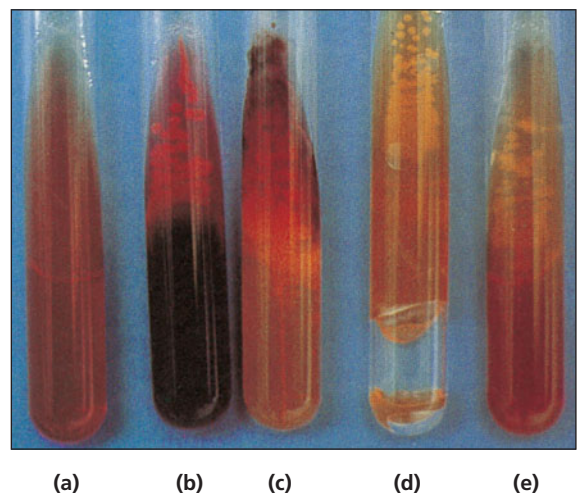

**Figure 24.2** Reactions in triple sugar–iron agar. (a) Uninoculated; (b) alkaline slant/acid butt, H<sub>2</sub>S; (c) alkaline slant/acid butt; (d) acid slant/acid butt, gas; and (e) acid slant/acid butt.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

Bacterial Species	CARBOHYDRATE FERMENTATION			H <sub>2</sub> S PRODUCTION	
	Butt Color and Reaction	Slant Color and Reaction	Carbohydrate Fermented	Blackening	H <sub>2</sub> S (+) or (-)
<i>E. coli</i>					
<i>E. aerogenes</i>					
<i>K. pneumoniae</i>					
<i>S. dysenteriae</i>					
<i>S. typhimurium</i>					
<i>P. vulgaris</i>					
<i>P. aeruginosa</i>					
<i>A. faecalis</i>					
<i>M. luteus</i>					
<i>L. lactis</i>					
<i>S. aureus</i>					
<i>B. cereus</i>					
<i>C. xerosis</i>					
Control					

## Review Questions

1. What is the purpose of the TSI test?



Identification of enteric bacilli is of prime importance in controlling intestinal infections by preventing contamination of food and water supplies. The groups of bacteria that can be found in the intestinal tract of humans and lower mammals are classified as members of the family **Enterobacteriaceae**. They are short, gram-negative, non-spore-forming bacilli. Included in this family are:


- 1. Pathogens** such as members of the genera *Salmonella* and *Shigella*.
- 2. Occasional pathogens** such as members of the genera *Proteus* and *Klebsiella*.
- 3. Normal intestinal flora** such as members of the genera *Escherichia* and *Enterobacter*, which are saprophytic inhabitants of the intestinal tract.

Differentiation of the principal groups of Enterobacteriaceae can be accomplished on the basis of their biochemical properties and enzymatic reactions in the presence of specific substrates. The **IMViC** series of tests (**indole, methyl red, Voges-Proskauer, and citrate utilization**) can be used.

**Figure 25.11** on page 177 shows the biochemical reactions that occur during the IMViC tests. It is designed to assist you in the execution and interpretation of each test.

The following experiments are designed for either a short or long version. The short version uses selected members of the enteric family. The long procedure makes use of bacterial species that do not belong solely to the Enterobacteriaceae. Nonenteric forms are included to acquaint you with the biochemical activities of other organisms grown in these media and to enable you to use these data for further comparisons of both types of bacteria. Selected organisms to be used in the

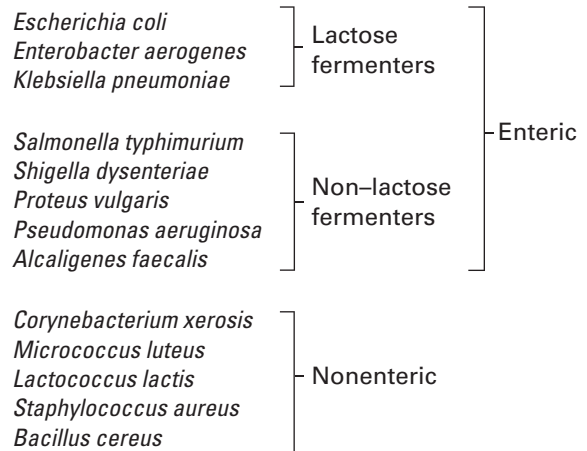long-version procedures are listed below. The enteric organisms are subdivided as lactose fermenters and non-lactose fermenters.



## CLINICAL APPLICATION

### Identification of Enteric Bacteria

The IMViC test is used to identify members of the Enterobacteriaceae, some of which are powerful pathogens such as members of the the genera *Shigella* and *Salmonella*, which cause intestinal infections. Identification of the causative agent may lead to the source of the infection, such as raw food (*Salmonella*) or fecal contamination of food (*Shigella*). This will aid in determining the possible number of individuals that have been exposed and who may require medical attention. This test uses the organisms' biochemical properties and enzymatic reactions on specific substrates as a means of identification.



## PART A Indole Production Test

### LEARNING OBJECTIVE

Once you have completed this test, you should be able to

1. Determine the ability of microorganisms to degrade the amino acid tryptophan.

### Principle

**Tryptophan** is an essential amino acid that can undergo oxidation by way of the enzymatic activities of some bacteria. Conversion of tryptophan into metabolic products is mediated by the enzyme **tryptophanase**. The chemistry of this reaction is illustrated in **Figure 25.1**. This ability to hydrolyze tryptophan with the production of indole is not a characteristic of all microorganisms and therefore serves as a biochemical marker.

In this experiment, SIM agar, which contains the substrate tryptophan, is used. The presence of indole is detectable by adding Kovac's reagent, which produces a cherry red reagent layer. This color is produced by the reagent, which is composed of *p*-dimethylaminobenzaldehyde, butanol, and hydrochloric acid. Indole is extracted from the medium into the reagent layer by the acidified butyl alcohol component and forms a complex with the *p*-dimethylaminobenzaldehyde, yielding the

cherry red color. The chemistry of this reaction is illustrated in **Figure 25.2**.

Cultures producing a red reagent layer following addition of Kovac's reagent are indole-positive. The absence of red coloration demonstrates that the substrate tryptophan was not hydrolyzed and indicates an indole-negative reaction.

### AT THE BENCH



### Materials

#### Cultures

24- to 48-hour Trypticase soy broth cultures of *E. coli*, *P. vulgaris*, and *E. aerogenes* for the short version. 24- to 48-hour Trypticase soy broth cultures of the 13 organisms listed on page 171 for the long version.

#### Media


SIM agar deep tubes per designated student group: 4 for the short version, 14 for the long version.

#### Reagent


Kovac's reagent.

#### Equipment


Bunsen burner, inoculating needle, test tube rack, and glassware marking pencil.

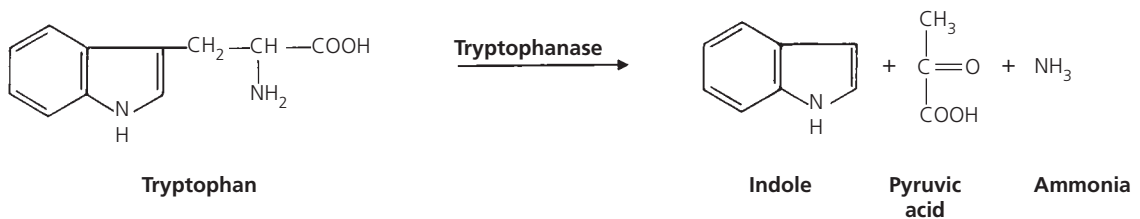

**Figure 25.1** Enzymatic degradation of tryptophan

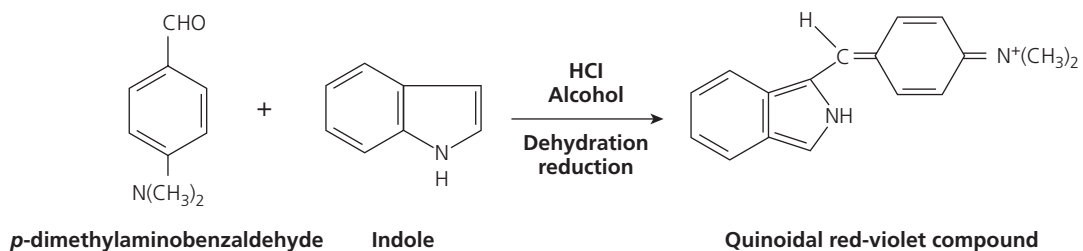

**Figure 25.2** Indole reaction with Kovac's reagent

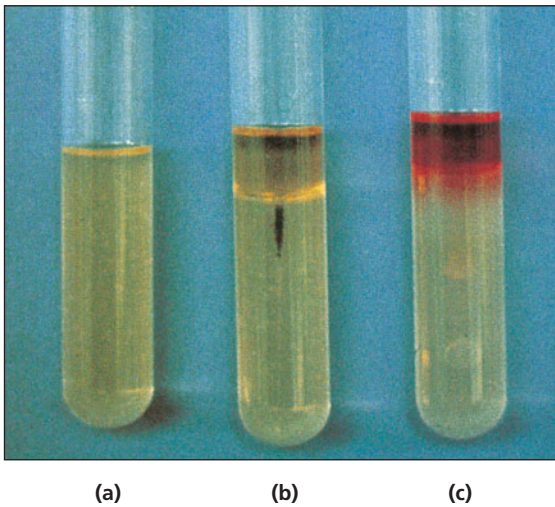

**Figure 25.3** Indole production test. (a) Uninoculated, (b) negative, and (c) positive.

## Procedure Lab One

1. Using aseptic technique, inoculate each experimental organism into its appropriately labeled deep tube by means of a stab inoculation. The last tube will serve as a control.
2. Incubate tubes for 24 to 48 hours at 37°C.

## Procedure Lab Two

1. Add 10 drops of Kovac's reagent to all deep tube cultures and agitate the cultures gently.
2. Examine the color of the reagent layer in each culture (refer to **Figure 25.3**). Record your results in the chart in the Lab Report.
3. Based on your observations, determine and record whether or not each organism was capable of hydrolyzing the tryptophan.

## PART B Methyl Red Test

### LEARNING OBJECTIVES

Once you have completed this experiment, you should be able to

1. Determine the ability of microorganisms to ferment glucose with the production and stabilization of high concentrations of acid end products.
2. Differentiate between all glucose-fermenting enteric organisms, particularly *E. coli* and *E. aerogenes*.

## Principle

The hexose monosaccharide **glucose** is the major substrate utilized by all enteric organisms for energy production. The end products of this process will vary depending on the specific enzymatic pathways present in the bacteria. In this test the pH indicator methyl red detects the presence of large concentrations of acid end products. Although most enteric microorganisms ferment glucose with the production of organic acids, this test is of value in the separation of *E. coli* and *E. aerogenes*.

Both of these organisms initially produce organic acid end products during the early incubation period. The low acidic pH (4) is stabilized and maintained by *E. coli* at the end of incubation. During the later incubation period, *E. aerogenes* enzymatically converts these acids to nonacidic end products such as 2,3-butanediol and acetoin (acetylmethylcarbinol), resulting in an elevated pH of approximately 6. The glucose fermentation reaction generated by *E. coli* is illustrated in **Figure 25.4**.

As shown, the methyl red indicator in the pH range of 4 will turn red, which is indicative of a positive test. At a pH of 6, still indicating the presence of acid but with a lower hydrogen ion concentration, the indicator turns yellow and is a negative test. Production and detection of the nonacidic end products from glucose fermentation by *E. aerogenes* is amplified in Part C of this exercise, the Voges-Proskauer test, which is performed simultaneously with the methyl red test.

## AT THE BENCH




## Materials

### Cultures


24- to 48-hour Trypticase soy broth cultures of *E. coli*, *E. aerogenes*, and *K. pneumoniae* for the short version. 24- to 48-hour Trypticase soy broth cultures of the 13 organisms listed on page 171 for the long version. In Lab Two, aliquots of these experimental cultures must be set aside for the Voges-Proskauer test.

### Media

MR-VP broth per designated student group: 4 for the short version, 14 for the long version.

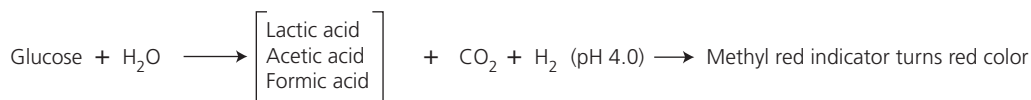

**Figure 25.4** Glucose fermentation reaction with methyl red pH reagent

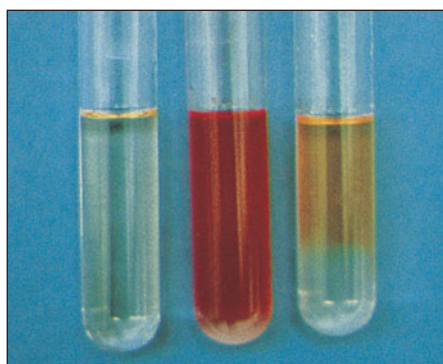

(a) (b) (c)

**Figure 25.5** Methyl red test. (a) Uninoculated, (b) positive, and (c) negative.

## Reagent

Methyl red indicator.

## Equipment

Bunsen burner, inoculating loop, test tubes, and glassware marking pencil.

## Procedure Lab One

1. Using aseptic technique, inoculate each experimental organism into its appropriately labeled tube of medium by means of a loop inoculation. The last tube will serve as a control.
2. Incubate all cultures for 24 to 48 hours at 37°C.

## Procedure Lab Two

1. Transfer approximately one-third of each culture into an empty test tube and set these tubes aside for the Voges-Proskauer test.
2. Add five drops of the methyl red indicator to the remaining aliquot of each culture.
3. Examine the color of all cultures (refer to **Figure 25.5**). Record the results in the chart in the Lab Report.
4. Based on your observations, determine and record whether or not each organism was capable of fermenting glucose with the production and maintenance of a high concentration of acid.

## PART C Voges-Proskauer Test

### LEARNING OBJECTIVE

Once you have completed this experiment, you should be able to


1. Differentiate further among enteric organisms such as *E. coli*, *E. aerogenes*, and *K. pneumoniae*.

## Principle


The Voges-Proskauer test determines the capability of some organisms to produce nonacidic or neutral end products, such as acetylmethylcarbinol, from the organic acids that result from glucose metabolism. This glucose fermentation, which is characteristic of *E. aerogenes*, is illustrated in **Figure 25.6**.

The reagent used in this test, Barritt's reagent, consists of a mixture of alcoholic  $\alpha$ -naphthol and 40% potassium hydroxide solution. Detection of acetylmethylcarbinol requires this end product to be oxidized to a diacetyl compound. This reaction will occur in the presence of the  $\alpha$ -naphthol catalyst and a guanidine group that is present in the peptone of the MR-VP medium. As a result, a pink complex is formed, imparting a rose color to the medium. The chemistry of this reaction is illustrated in **Figure 25.7**.


Development of a deep rose color in the culture 15 minutes following the addition of Barritt's reagent is indicative of the presence of acetylmethylcarbinol and represents a positive result. The absence of rose coloration is a negative result.

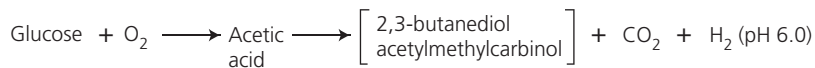

**Figure 25.6** Glucose fermentation by *E. aerogenes*

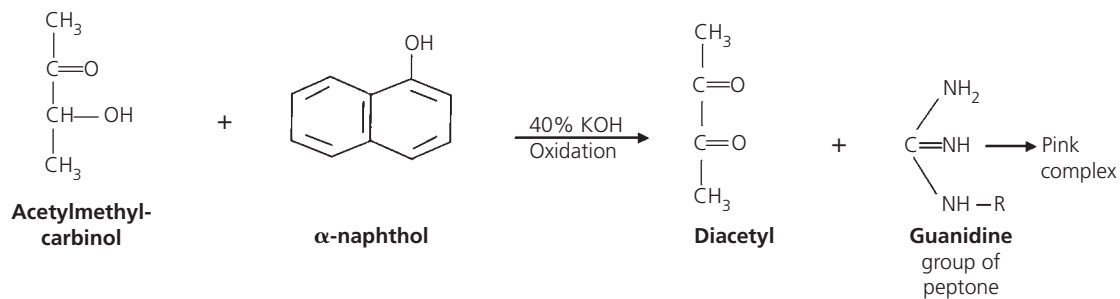

**Figure 25.7** Acetylmethylcarbinol reaction with Barritt's reagent

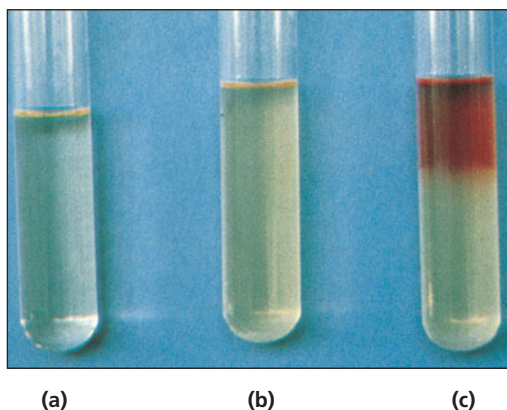

**Figure 25.8** Voges-Proskauer test. (a) Uninoculated, (b) negative, and (c) positive.

## AT THE BENCH



## Materials

### Cultures

24- to 48-hour Trypticase soy broth cultures of *E. coli*, *E. aerogenes*, and *K. pneumoniae* for the short version. 24- to 48-hour Trypticase soy broth cultures of the 13 organisms listed on page 171 for the long version. *Note: Aliquots of these experimental cultures must be set aside from the methyl red test.*

### Reagent

Barritt's reagents A and B.

### Equipment

Bunsen burner, inoculating loop, and glassware marking pencil.

## Procedure Lab One

Refer to the methyl red test in Part B of this exercise.

## Procedure Lab Two


1. To the aliquots of each broth culture separated during the methyl red test, add 10 drops of Barritt's reagent A and shake the cultures. Immediately add 10 drops of Barritt's reagent B and shake. Reshake the cultures every 3 to 4 minutes.
2. Examine the color of the cultures 15 minutes after the addition of Barritt's reagent. Refer to **Figure 25.8**. Record your results in the Lab Report.
3. Based on your observations, determine and record whether or not each organism was capable of fermenting glucose with ultimate production of acetylmethylcarbinol.

## PART D Citrate Utilization Test

### LEARNING OBJECTIVE

*Once you have completed this experiment, you should be able to*

1. Differentiate among enteric organisms on the basis of their ability to ferment citrate as a sole source of carbon.

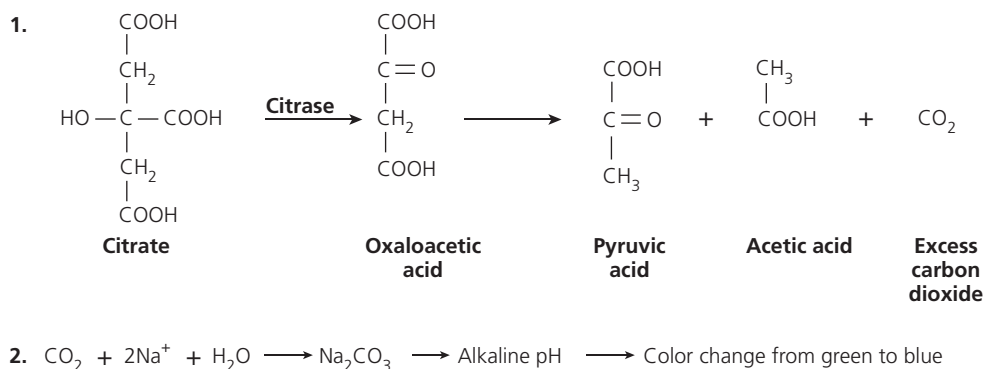

**Figure 25.9** Enzymatic degradation of citrate

## Principle

In the absence of fermentable glucose or lactose, some microorganisms are capable of using **citrate** as a carbon source for their energy. This ability depends on the presence of a **citrate permease** that facilitates the transport of citrate in the cell. Citrate is the first major intermediate in the Krebs cycle and is produced by the condensation of active acetyl with oxaloacetic acid. Citrate is acted on by the enzyme **citrase**, which produces oxaloacetic acid and acetate. These products are then enzymatically converted to pyruvic acid and carbon dioxide. During this reaction the medium becomes alkaline—the carbon dioxide that is generated combines with sodium and water to form sodium carbonate, an alkaline product. The presence of sodium carbonate changes the bromthymol blue indicator incorporated into the medium from green to deep Prussian blue. The chemistry of this reaction is illustrated in **Figure 25.9**.

Following incubation, citrate-positive cultures are identified by the presence of growth on the surface of the slant, which is accompanied by blue coloration. Citrate-negative cultures will show no growth, and the medium will remain green.

## AT THE BENCH



## Materials

### Cultures

24- to 48-hour Trypticase soy broth cultures of *E. coli*, *E. aerogenes*, and *K. pneumoniae* for the short version. 24- to 48-hour Trypticase soy broth cultures of the 13 organisms listed on page 171 for the long version.

## Media


Simmons citrate agar slants per designated student group: 4 for the short version, 14 for the long version.

## Equipment

Bunsen burner, inoculating needle, test tube rack, and glassware marking pencil.

## Procedure Lab One

- Using aseptic technique, inoculate each organism into its appropriately labeled tube by means of streak inoculation. The last tube will serve as a control.
- Incubate all cultures for 24 to 48 hours at 37°C.

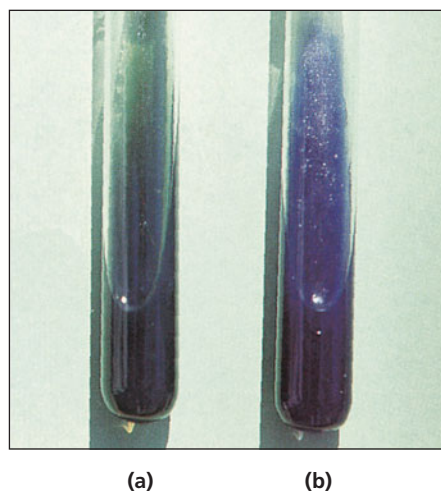


**Figure 25.10** Citrate utilization test. (a) Tube is negative, showing no growth on slant surface. (b) Tube is positive, showing growth on slant surface.

# Procedure Lab Two

1. Examine all agar slant cultures for the presence or absence of growth and coloration of the medium. Refer to **Figure 25.10**. Record your results in the chart in the Lab Report.

2. Based on your observations, determine and record whether or not each organism was capable of using citrate as its sole source of carbon.

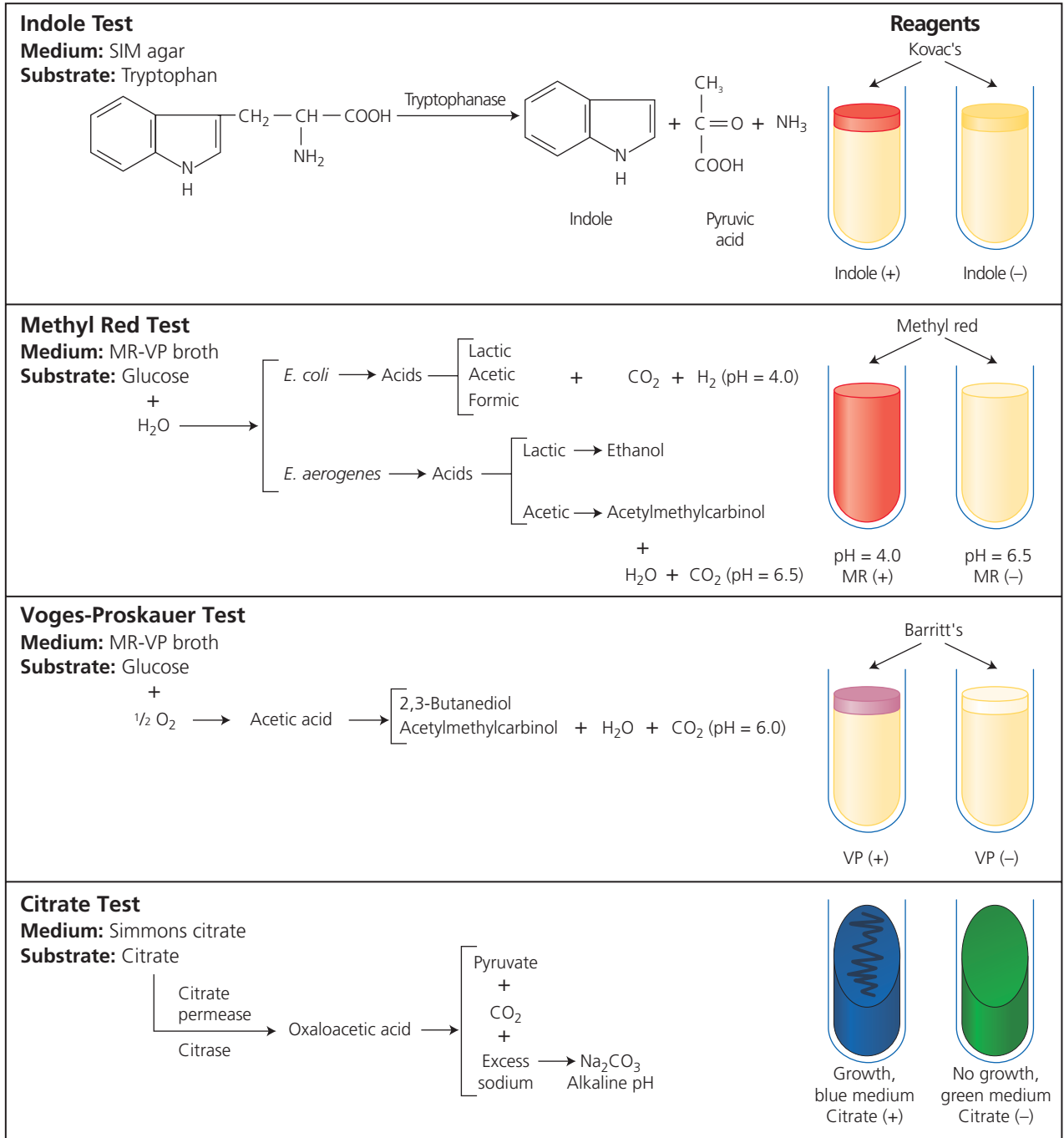

**Figure 25.11** Summary of IMViC reactions

*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

### PART A: Indole Production Test

Bacterial Species	Color of Reagent Layer	Tryptophan Hydrolysis (+) or (-)
<i>E. coli</i>		
<i>E. aerogenes</i>		
<i>K. pneumoniae</i>		
<i>S. dysenteriae</i>		
<i>S. typhimurium</i>		
<i>P. vulgaris</i>		
<i>P. aeruginosa</i>		
<i>A. faecalis</i>		
<i>M. luteus</i>		
<i>L. lactis</i>		
<i>S. aureus</i>		
<i>B. cereus</i>		
<i>C. xerosis</i>		
Control		



## PARTS B, C: Methyl Red Test and Voges-Proskauer Test




Bacterial Species	METHYL RED TEST		VOGES-PROSKAUER TEST	
	Color of Medium	(+) or (-)	Color of Medium	(+) or (-)
<i>E. coli</i>				
<i>E. aerogenes</i>				
<i>K. pneumoniae</i>				
<i>S. dysenteriae</i>				
<i>S. typhimurium</i>				
<i>P. vulgaris</i>				
<i>P. aeruginosa</i>				
<i>A. faecalis</i>				
<i>M. luteus</i>				
<i>L. lactis</i>				
<i>S. aureus</i>				
<i>B. cereus</i>				
<i>C. xerosis</i>				
Control				

## PART D: Citrate Utilization Test

Bacterial Species	Presence or Absence of Growth (+) or (-)	Color of Medium	Citrate Utilization (+) or (-)
<i>E. coli</i>			
<i>E. aerogenes</i>			
<i>K. pneumoniae</i>			
<i>S. dysenteriae</i>			
<i>S. typhimurium</i>			
<i>P. vulgaris</i>			
<i>P. aeruginosa</i>			
<i>A. faecalis</i>			
<i>M. luteus</i>			
<i>L. lactis</i>			
<i>S. aureus</i>			
<i>B. cereus</i>			
<i>C. xerosis</i>			
Control			

## Review Questions

1. Discuss the medical significance of the IMViC series of tests.
2. Explain the chemical mechanism for detecting indole in a bacterial culture.

3. Account for the development of alkalinity in cultures capable of using citrate as their sole carbon source.
4.  In the carbohydrate fermentation test, we found that both *E. coli* and *E. aerogenes* produced the end products acid and gas. Account for the fact that *E. coli* is methyl red–positive and *E. aerogenes* is methyl red–negative.
5.  The end products of tryptophan degradation are indole and pyruvic acid. Why do we test for the presence of indole rather than pyruvic acid as the indicator of tryptophanase activity?
6.  Simmons citrate medium contains primarily inorganic ammonium, potassium, and sodium salts, plus organic citrate. What is the rationale for using a medium with this type of composition for the performance of the citrate utilization test?

## LEARNING OBJECTIVES


Once you have completed this experiment, you will be able to determine

1. The ability of microorganisms to produce hydrogen sulfide from substances such as the sulfur-containing amino acids or inorganic sulfur compounds.
2. Mobility of microorganisms in SIM agar.


## Principle

There are two major fermentative pathways by which some microorganisms are able to produce hydrogen sulfide (H<sub>2</sub>S).

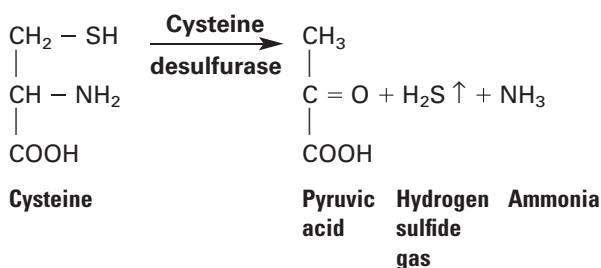**Pathway 1:** Gaseous H<sub>2</sub>S may be produced by the reduction (hydrogenation) of organic sulfur present in the amino acid cysteine, which is a component of peptones contained in the medium. These peptones are degraded by microbial enzymes to amino acids, including the sulfur-containing amino acid cysteine. This amino acid in the presence of a **cysteine desulfurase** loses the sulfur atom, which is then reduced by the addition of hydrogen from water to form bubbles of hydrogen sulfide gas (H<sub>2</sub>S ↑) as illustrated:



**Pathway 2:** Gaseous H<sub>2</sub>S may also be produced by the reduction of inorganic sulfur compounds such as the thiosulfates (S<sub>2</sub>O<sub>3</sub><sup>2-</sup>), sulfates (SO<sub>4</sub><sup>2-</sup>), or sulfites (SO<sub>3</sub><sup>2-</sup>). The medium contains sodium thiosulfate, which certain microorganisms are capable of reducing to sulfite with the liberation of hydrogen sulfide. The sulfur atoms act as hydrogen acceptors during oxidation of the inorganic compound as illustrated in the following:


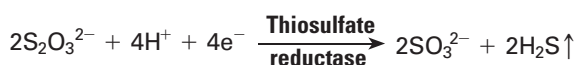


Thiosulfate	Sulfite	Hydrogen Sulfide gas
-------------	---------	----------------------

In this experiment the SIM medium contains peptone and sodium thiosulfate as the sulfur substrates; ferrous sulfate (FeSO<sub>4</sub>), which behaves as the H<sub>2</sub>S indicator; and sufficient agar to make the medium semisolid and thus enhance anaerobic respiration. Regardless of which pathway is used, the hydrogen sulfide gas is colorless and therefore not visible. Ferrous ammonium sulfate in the medium serves as an indicator by combining with the gas, forming an insoluble black ferrous sulfide precipitate that is seen along the line of the stab inoculation and is indicative of H<sub>2</sub>S production. Absence of the precipitate is evidence of a negative reaction. The overall reactions for both pathways and their interpretation are illustrated in **Figure 26.1**.

## Motility

SIM agar may also be used to detect motile organisms. Motility is recognized when culture growth (turbidity) of flagellated organisms is not restricted to the line of inoculation. Growth of nonmotile organisms is confined to the line of inoculation.

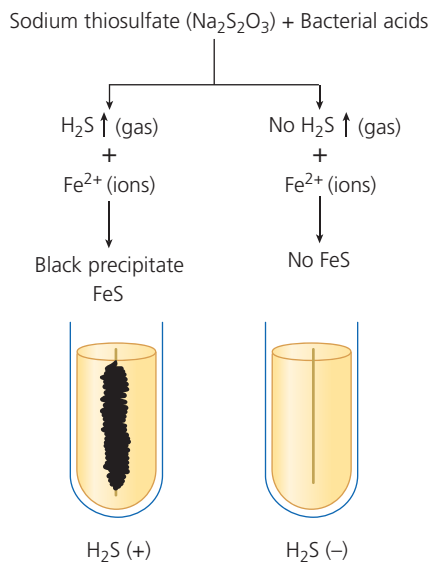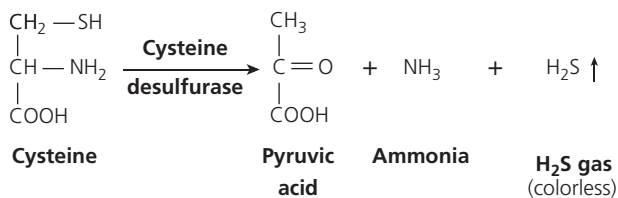

**Figure 26.1** Detection of hydrogen sulfide

## CLINICAL APPLICATION

### Identifying Intestinal Pathogens


While generally considered a self-limiting symptom, diarrhea due to *Proteus* is initially difficult to differentiate from early stages of the more severe bloody diarrhea (dysentery) associated with some *Shigella* or *Salmonella* species. Bacteria belonging to the genera *Salmonella* and *Proteus* enzymatically metabolize inorganic sulfur compounds and sulfur-containing amino acids, producing H<sub>2</sub>S. The hydrogen sulfide test is one way to separate and identify *Shigella dysenteriae*, which does not produce H<sub>2</sub>S, from *Proteus* and *Salmonella*.

## AT THE BENCH

### Materials

#### Cultures

24- to 48-hour Trypticase soy broth cultures of *Enterobacter aerogenes*, *Shigella dysenteriae*, *Proteus vulgaris*, and *Salmonella typhimurium* for

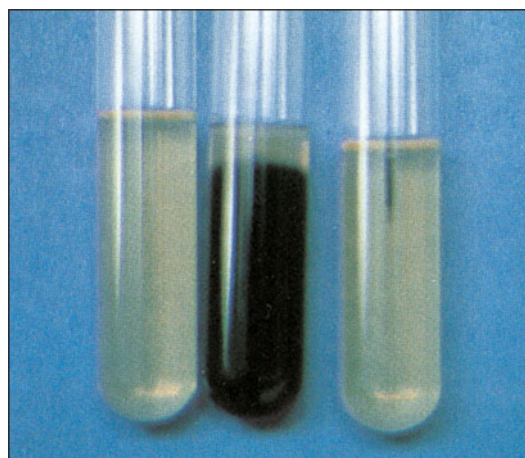

**Figure 26.2** Hydrogen sulfide production test. (a) Negative, (b) positive with motility, and (c) positive with no motility.

the short version. 24- to 48-hour Trypticase soy broth cultures of the 13 organisms listed on page 154 for the long version.

### Media

SIM agar deep tubes per designated student group: 5 for the short version, 14 for the long version.

### Equipment

Bunsen burner, inoculating needle, test tube rack, and glassware marking pencil.

### Procedure Lab One

1. Aseptically inoculate each experimental organism into its appropriately labeled tube by means of stab inoculation. The last tube will serve as a control.
2. Incubate all cultures for 24 to 48 hours at 37°C.

### Procedure Lab Two

1. Examine all SIM cultures for the presence or absence of black coloration along the line of the stab inoculation. Refer to **Figure 26.2**, and record your results in the chart provided in the Lab Report.
2. Based on your observations, determine and record whether or not each organism was capable of producing hydrogen sulfide.
3. Observe all cultures for the presence (+) or absence (-) of motility. Record your results in the chart in the Lab Report.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

Bacterial Species	Color of Medium	H <sub>2</sub> S Production (+) or (-)	Motility (+) or (-)
<i>E. coli</i>			
<i>E. aerogenes</i>			
<i>K. pneumoniae</i>			
<i>S. dysenteriae</i>			
<i>S. typhimurium</i>			
<i>P. vulgaris</i>			
<i>P. aeruginosa</i>			
<i>A. faecalis</i>			
<i>M. luteus</i>			
<i>L. lactis</i>			
<i>S. aureus</i>			
<i>B. cereus</i>			
<i>C. xerosis</i>			
Control			

## Review Questions

1. Distinguish between the types of substrates available to cells for H<sub>2</sub>S production.

2. Explain how SIM medium is used to detect motility.

3. Explain the function of the ferrous ammonium sulfate in SIM agar.



4. Why is *P. vulgaris* H<sub>2</sub>S-positive and *E. aerogenes* H<sub>2</sub>S-negative?



5. A stool specimen of a patient with severe diarrhea was cultured in a series of specialized media for isolation of enteric organisms. The cultures yielded three isolates that were species of *Salmonella*, *Shigella*, and *Escherichia*. Explain why the H<sub>2</sub>S production test would be diagnostically significant.

## LEARNING OBJECTIVE

Once you have completed this experiment, you should be able to


1. Determine the ability of microorganisms to degrade urea by means of the enzyme urease.

## Principle

Urease, which is produced by some microorganisms, is an enzyme that is especially helpful in the identification of *Proteus vulgaris*. Although other organisms may produce urease, their action on the substrate urea tends to be slower than that seen with *Proteus* species. Therefore, this test serves to rapidly distinguish members of this genus from other non-lactose-fermenting enteric microorganisms.

**Urease** is a hydrolytic enzyme that attacks the nitrogen and carbon bond in amide compounds such as urea and forms the alkaline end product ammonia. This chemical reaction is illustrated in **Figure 27.1**.

The presence of urease is detectable when the organisms are grown in a urea broth medium containing the pH indicator phenol red. As the substrate urea is split into its products, the presence of ammonia creates an alkaline environment that causes the phenol red to turn to a deep pink. This is a positive reaction for the presence of urease. Failure of a deep pink color to develop is evidence of a negative reaction.

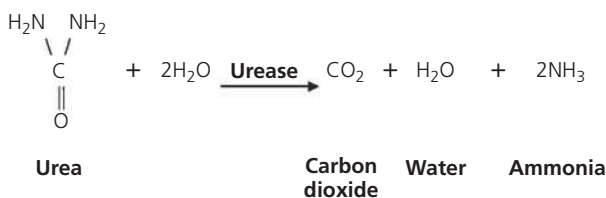

**Figure 27.1** Enzymatic degradation of urea

## CLINICAL APPLICATION

### Pathogens and the Urease Test

The urease test is primarily used to distinguish the small number of urease-positive enterics from other non-lactose-fermenting enteric bacteria. Many enterics can degrade urea but only a few are termed *rapid urease positive organisms*. While part of the normal flora, these commensals have been identified as opportunistic pathogens. Members of the gastrodeudonal commensals are included among this group of organisms.

## AT THE BENCH



## Materials

### Cultures

24- to 48-hour Trypticase soy broth cultures of *Escherichia coli*, *Proteus vulgaris*, *Klebsiella pneumoniae*, and *Salmonella typhimurium* for the short version. 24- to 48-hour Trypticase soy broth cultures of the 13 organisms listed on page 154 for the long version.

### Media

Urea broth per designated student group: 5 for the short version, 14 for the long version.

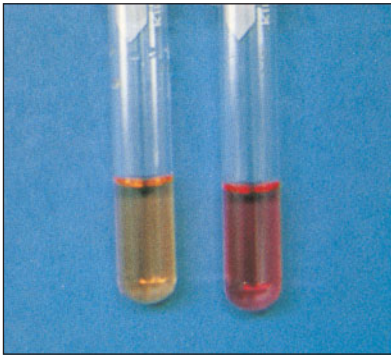
### Equipment

Bunsen burner, inoculating loop, test tube rack, and glassware marking pencil.

## Procedure Lab One

1. Using aseptic technique, inoculate each experimental organism into its appropriately labeled tube by means of loop inoculation. The last tube will serve as a control.
2. Incubate cultures 24 to 48 hours at 37°C.





(a) (b)

**Figure 27.2 Urease test. (a)** Negative and **(b)** positive.

## Procedure Lab Two

1. Examine all urea broth cultures for color (refer to **Figure 27.2**). Record your results in the chart in the Lab Report.
2. Based on your observations, determine and record whether or not each organism was capable of hydrolyzing the substrate urea.

Name: \_\_\_\_\_


Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

Bacterial Species	Color of Medium	Urea Hydrolysis (+) or (-)
<i>E. coli</i>		
<i>E. aerogenes</i>		
<i>K. pneumoniae</i>		
<i>S. dysenteriae</i>		
<i>S. typhimurium</i>		
<i>P. vulgaris</i>		
<i>P. aeruginosa</i>		
<i>A. faecalis</i>		
<i>M. luteus</i>		
<i>L. lactis</i>		
<i>S. aureus</i>		
<i>B. cereus</i>		
<i>C. xerosis</i>		
Control		

## Review Questions

1. Explain the mechanism of urease activity.

2. Explain the function of phenol red in the urea broth medium.
3. Explain how the urease test is useful for identifying members of the genus *Proteus*.
4.  A swollen can of chicken soup is examined by the public health laboratory and found to contain large numbers of gram-negative, H<sub>2</sub>S-positive bacilli. Which biochemical tests would you perform to identify the genus of the contaminant? Justify your test choices.

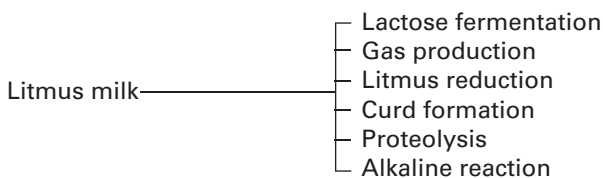
## LEARNING OBJECTIVE

Once you have completed this experiment, you should be able to

1. Differentiate among microorganisms that enzymatically transform different milk substrates into varied metabolic end products.


## Principle

The major milk substrates capable of transformation are the milk sugar lactose and the milk proteins casein, lactalbumin, and lactoglobulin. To distinguish among the metabolic changes produced in milk, a pH indicator, the oxidation-reduction indicator litmus, is incorporated into the medium. Litmus milk now forms an excellent differential medium in which microorganisms can metabolize milk substrates depending on their enzymatic complement. A variety of different biochemical changes result, as follows:



## Lactose Fermentation

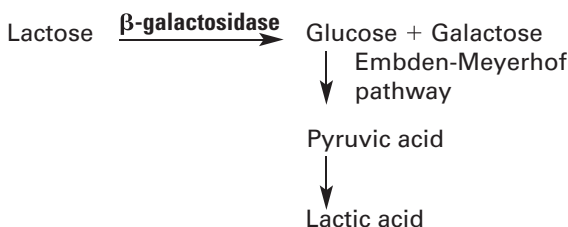Organisms capable of using **lactose** as a carbon source for energy production utilize the inducible enzyme  **$\beta$ -galactosidase** and degrade lactose as follows:




The presence of **lactic acid** is easily detected because litmus is purple at a neutral pH and turns pink when the medium is acidified to an approximate pH of 4.

## Gas Formation

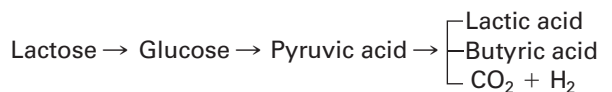The end products of the microbial fermentation of lactose are likely to include the **gases**  $\text{CO}_2 \uparrow + \text{H}_2 \uparrow$ . The presence of gas may be seen as separations of the curd or by the development of tracks or fissures within the curd as gas rises to the surface.

## Litmus Reduction

Fermentation is an anaerobic process involving biooxidations that occur in the absence of molecular oxygen. These oxidations may be visualized as the removal of hydrogen (dehydrogenation) from a substrate. Since hydrogen ions cannot exist in the free state, there must be an immediate and concomitant electron acceptor available to bind these hydrogen ions, or else oxidation-reduction reactions are not possible and cells cannot manufacture energy. In the litmus milk test, **litmus** acts as such an acceptor. While in the oxidized state, the litmus is purple; when it accepts hydrogen from a substrate, it will become reduced and turn white or milk-colored. This oxidation of lactose, which produces lactic acid, butyric acid,  $\text{CO}_2 \uparrow$ , and  $\text{H}_2 \uparrow$ , is as follows:




The excess hydrogen is now accepted by the hydrogen acceptor litmus, which turns white and is said to be reduced.

## Curd Formation

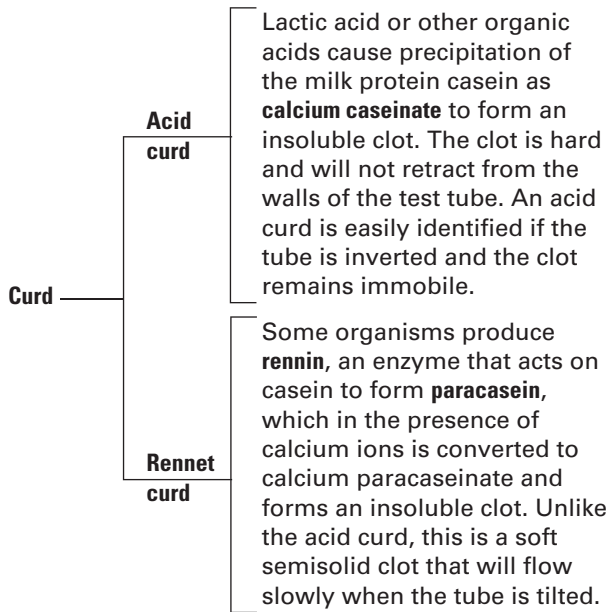The biochemical activities of different microorganisms grown in litmus milk may result in the production of two distinct types of curds (clots). Curds are designated as either acid or rennet,

depending on the biochemical mechanism responsible for their formation.




### Proteolysis (Peptonization)

The inability of some microorganisms to obtain their energy by way of lactose fermentation means they must use other nutritional sources such as proteins for this purpose (see Figure 23.6 on page 164). By means of proteolytic enzymes, these organisms hydrolyze the milk proteins, primarily casein, into their basic building blocks, namely amino acids. This digestion of proteins is accompanied by the evolution of large quantities of ammonia, resulting in an alkaline pH in the medium. The litmus turns deep purple in the upper portion of the tube, while the medium begins to lose body and produces a translucent, brown, wheylike appearance as the protein is hydrolyzed to amino acids.

### Alkaline Reaction

An alkaline reaction is evident when the color of the medium remains unchanged or changes to a deeper blue. This reaction is indicative of the partial degradation of **casein** into **shorter polypeptide chains**, with the simultaneous release of alkaline end products that are responsible for the observable color change.

**Figure 28.1** and **28.2** show the possible litmus milk reactions and their appearance following the appropriate incubation of the cultures.

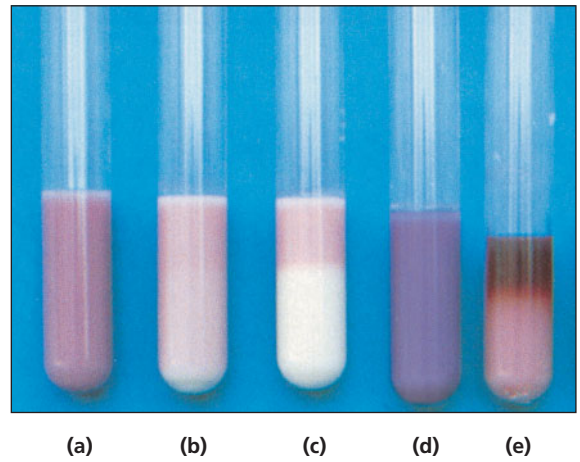

**Figure 28.1** Litmus milk reactions. (a) Uninoculated, (b) acid, (c) acid with reduction and curd, (d) alkaline, and (e) proteolysis.

## CLINICAL APPLICATION

### Differentiating Enterobacteriaceae and *Clostridium*.

The litmus milk test differentiates members of the Enterobacteriaceae from other gram-negative bacilli based on the enterics' ability to reduce litmus. It is also used to differentiate members within the genus *Clostridium*. Watery diarrhea caused by *C. perfringens* (contaminated food) is generally considered self-limiting. But diarrhea caused by *C. difficile* may be associated with antibiotic use that has removed the normal flora of the colon.


## AT THE BENCH



## Materials

### Cultures

24- to 48-hour Trypticase soy broth cultures of *Escherichia coli*, *Alcaligenes faecalis*, *Lactococcus lactis*, and *Pseudomonas aeruginosa* for the short version. 24- to 48-hour Trypticase soy broth cultures of the 13 organisms listed on page 154 for the long version.

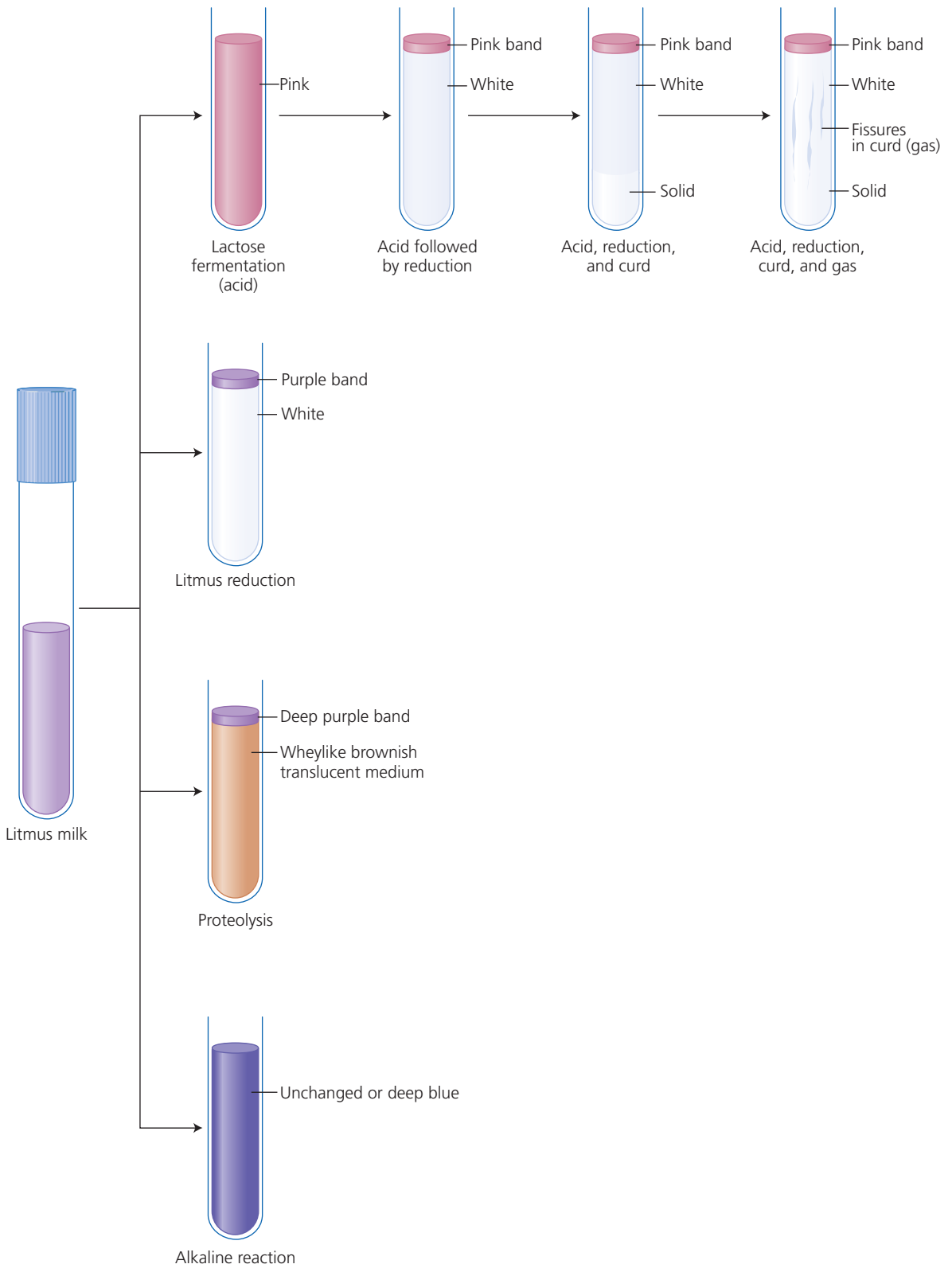

**Figure 28.2** Summary of possible litmus milk reactions

## Media

Litmus milk broth per designated student group: 5 for the short version, 14 for the long version.

## Equipment

Bunsen burner, inoculating loop, test tube rack, and glassware marking pencil.

## Procedure Lab One

1. Using aseptic technique, inoculate each experimental organism into its appropriately labeled tube by means of a loop inoculation. The last tube will serve as a control.
2. Incubate all cultures for 24 to 48 hours at 37°C.

## Procedure Lab Two

1. Examine all the litmus milk cultures for color and consistency of the medium. Record your results in the chart in the Lab Report.
2. Based on your observations, determine and record the type(s) of reaction(s) that have taken place in each culture.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

Bacterial Species	Appearance of Medium	Litmus Milk Reactions
<i>E. coli</i>		
<i>E. aerogenes</i>		
<i>K. pneumoniae</i>		
<i>S. dysenteriae</i>		
<i>S. typhimurium</i>		
<i>P. vulgaris</i>		
<i>P. aeruginosa</i>		
<i>A. faecalis</i>		
<i>M. luteus</i>		
<i>L. lactis</i>		
<i>S. aureus</i>		
<i>B. cereus</i>		
<i>C. xerosis</i>		
Control		


## Review Questions

1. Distinguish between acid and rennet curds.



2. Describe the litmus milk reactions that may occur when proteins are metabolized as an energy source.

3. Explain how the litmus in the litmus milk acts as a redox indicator.

4.  Can a litmus milk culture show a pink band at the top and a brownish translucent layer at the bottom? Explain.

5.  Explain why litmus milk is considered a good differential medium.

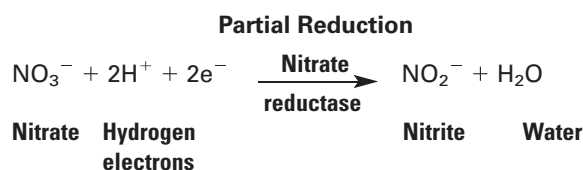
## LEARNING OBJECTIVE

Once you have completed this experiment, you should be able to

1. Determine the ability of some microorganisms to reduce nitrates ( $\text{NO}_3^-$ ) to nitrites ( $\text{NO}_2^-$ ) or beyond the nitrite stage.

## Principle


The reduction of nitrates by some aerobic and facultative anaerobic microorganisms occurs in the absence of molecular oxygen, an anaerobic process. In these organisms anaerobic respiration is an oxidative process whereby the cell uses inorganic substances such as nitrates ( $\text{NO}_3^-$ ) or sulfates ( $\text{SO}_4^{2-}$ ) to supply oxygen that is subsequently utilized as a final hydrogen acceptor during energy formation. The biochemical transformation may be visualized as follows:



Some organisms possess the enzymatic capacity to act further on nitrites to reduce them to ammonia ( $\text{NH}_3^+$ ) or molecular nitrogen ( $\text{N}_2$ ). These reactions may be described as follows:




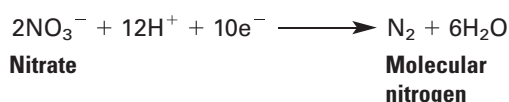or




Nitrate reduction can be determined by cultivating organisms in a nitrate broth medium. The medium is basically a nutrient broth supplemented with 0.1% potassium nitrate ( $\text{KNO}_3$ ) as the nitrate substrate. In addition, the medium is made into a semisolid by the addition of 0.1% agar. The semisolidity impedes the diffusion of oxygen into the medium, thereby favoring the anaerobic requirement necessary for nitrate reduction.

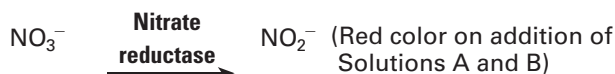
Following incubation of the cultures, an organism's ability to reduce nitrates to nitrites is determined by the addition of two reagents: Solution A, which is sulfanilic acid, followed by Solution B, which is  $\alpha$ -naphthylamine. *Note: This should not be confused with Barritt's reagent.* Following reduction, the addition of Solutions A and B will produce an immediate cherry red color.



Cultures not producing a color change suggest one of two possibilities: (1) nitrates were not reduced by the organism, or (2) the organism possessed such potent **nitrate reductase** enzymes that nitrates were rapidly reduced beyond nitrites to ammonia or even molecular nitrogen. To determine whether or not nitrates were reduced past the nitrite stage, a small amount of zinc powder is added to the basically colorless cultures already containing Solutions A and B. Zinc reduces nitrates to nitrites. The development of red color therefore verifies that nitrates were not reduced to nitrites by the organism. If nitrates were not reduced, a negative nitrate reduction reaction has occurred. If the addition of zinc does not produce a color change, the nitrates in the medium were reduced beyond nitrites to ammonia or nitrogen gas. This is a positive reaction, as shown in **Figure 29.1**. Results of nitrate reduction tests are shown in **Figure 29.2**.

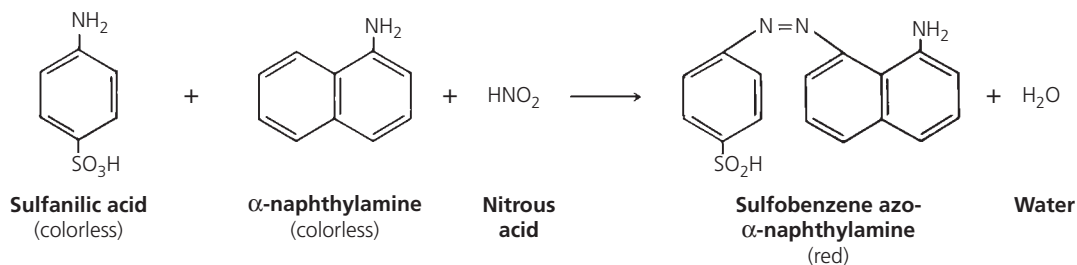

**Figure 29.1** Formation of colored complex indicative of  $\text{NO}_3^-$  reduction

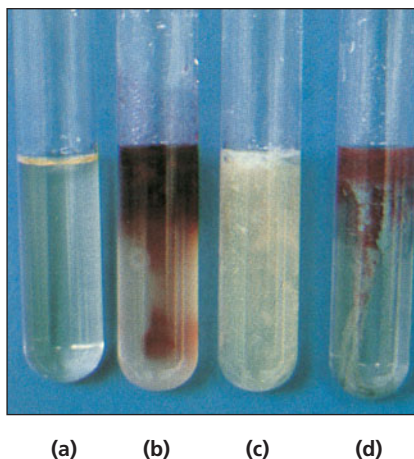

**Figure 29.2** Nitrate reduction tests. (a) Uninoculated, (b) positive with Solutions A + B, (c) positive with Solutions A + B + zinc powder, and (d) negative with Solutions A + B + zinc powder.

## CLINICAL APPLICATION

### Differentiating *Mycobacterium tuberculosis* from Non-tubercle *Mycobacterium*

This test is used to identify intestinal bacteria that are able to reduce nitrates to nitrites. When presented with a patient that exhibits the symptoms of tuberculosis and is positive for tubercles on an x-ray, a sputum sample will be tested for *Mycobacterium*. To distinguish between *Mycobacterium tuberculosis* and other *Mycobacterium* species a nitrate reduction test will be used since *M. tuberculosis* is the only *Mycobacterium* species with this capacity.

## AT THE BENCH



## Materials

### Cultures

24- to 48-hour Trypticase soy broth cultures of *Escherichia coli*, *Alcaligenes faecalis*, and

*Pseudomonas aeruginosa* for the short version. 24- to 48-hour Trypticase soy broth cultures of the 13 organisms listed on page 154 for the long version.

### Media

Trypticase nitrate broth per designated student group: 4 for the short version, 14 for the long version.

### Reagents

Solution A (sulfanilic acid), Solution B ( $\alpha$ -naphthylamine), and zinc powder.

### Equipment

Bunsen burner, inoculating loop, test tube rack, and glassware marking pencil.

## Procedure Lab One

- Using aseptic technique, inoculate each experimental organism into its appropriately labeled tube by means of a loop inoculation. The last tube will serve as a control.
- Incubate all cultures for 24 to 48 hours at 37°C.

## Procedure Lab Two

- Add five drops of Solution A and then five drops of Solution B to all nitrate broth cultures. Observe and record in the Lab Report chart whether or not a red coloration develops in each of the cultures.
- Add a minute quantity of zinc to the cultures in which no red color developed. Observe and record whether or not red coloration develops in each of the cultures.
- On the basis of your observations, determine and record in the Lab Report chart whether or not each organism was capable of nitrate reduction. Identify the end product ( $\text{NO}_2^-$  or  $\text{NH}_3^+/\text{N}_2$ ), if any, that is present.

Name: \_\_\_\_\_


Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

Bacterial Species	Red Coloration with Solutions A and B (+) or (-)	Red Coloration with Zinc (+) or (-)	Nitrate Reductions (+) or (-)	End Products
<i>E. coli</i>				
<i>E. aerogenes</i>				
<i>K. pneumoniae</i>				
<i>S. dysenteriae</i>				
<i>S. typhimurium</i>				
<i>P. vulgaris</i>				
<i>P. aeruginosa</i>				
<i>A. faecalis</i>				
<i>M. luteus</i>				
<i>L. lactis</i>				
<i>S. aureus</i>				
<i>B. cereus</i>				
<i>C. xerosis</i>				
Control				

## Review Questions

1. Explain the function of the 0.1% agar in the nitrate medium.

2. Explain the functions of Solutions A and B.
3. If a culture does not undergo a color change on the addition of Solutions A and B, explain how you would interpret this result.
4. Explain why the development of a red color on the addition of zinc is a negative test.
5.  Discuss the relationship between an organism's ability to reduce nitrate past the nitrite stage and that organism's proteolytic activity.

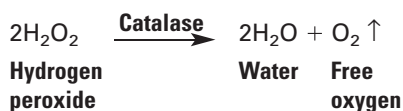
## LEARNING OBJECTIVE

Once you have completed this experiment, you should be able to

1. Determine the ability of some microorganisms to degrade hydrogen peroxide by producing the enzyme catalase.

## Principle

During aerobic respiration, microorganisms produce hydrogen peroxide and, in some cases, an extremely toxic superoxide. Accumulation of these substances will result in death of the organism unless they can be enzymatically degraded. These substances are produced when aerobes, facultative anaerobes, and microaerophiles use the aerobic respiratory pathway, in which oxygen is the final electron acceptor, during degradation of carbohydrates for energy production. Organisms capable of producing **catalase** rapidly degrade hydrogen peroxide as illustrated:



Aerobic organisms that lack catalase can degrade especially toxic superoxides using the enzyme **superoxide dismutase**; the end product of a superoxide dismutase is  $\text{H}_2\text{O}_2$ , but this is less toxic to the bacterial cells than are the superoxides.

The inability of strict anaerobes to synthesize catalase, peroxidase, or superoxide dismutase may explain why oxygen is poisonous to these microorganisms. In the absence of these enzymes, the toxic concentration of  $\text{H}_2\text{O}_2$  cannot be degraded when these organisms are cultivated in the presence of oxygen.

Catalase production can be determined by adding the substrate  $\text{H}_2\text{O}_2$  to an appropriately incubated Trypticase soy agar slant culture. If catalase is present, the chemical reaction mentioned is indicated by

bubbles of free oxygen gas ( $\text{O}_2 \uparrow$ ). This is a positive catalase test; the absence of bubble formation is a negative catalase test. **Figure 30.1** shows the results of the catalase test using (a) the tube method, (b) the plate method, and (c) slide method.

## CLINICAL APPLICATION

Differentiation of *Staphylococci*, *Streptococci*, and Enterobacteriaceae

The catalase test is used for the biochemical differentiation of catalase-positive *Staphylococci* and catalase-negative *Streptococci*, as well as members of the Enterobacteriaceae. With the increasing worry about methicillin-resistant strains of *Staphylococcus* in hospital settings, the catalase test is a quick and easy way to differentiate *S. aureus*, which may be MRSA, from other *Staphylococcus* species that have exhibited lower incidences of methicillin-resistance.

## AT THE BENCH



## Materials

## Cultures

24- to 48-hour Trypticase soy broth cultures of *Staphylococcus aureus*, *Micrococcus luteus*, and *Lactococcus lactis* for the short version. 24- to 48-hour Trypticase soy broth cultures of the 13 organisms listed on page 154 for the long version.

## Media

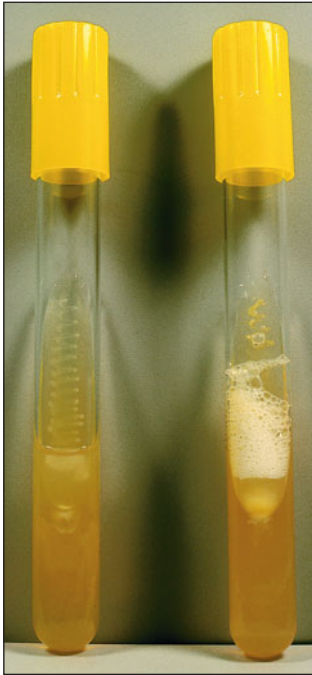
Trypticase soy agar slants per designated student group: 4 for the short version, 14 for the long version.

## Reagent


3% hydrogen peroxide.

## Equipment


*Tube method:* Bunsen burner, inoculating loop, test tube rack, and glassware marking pencil.



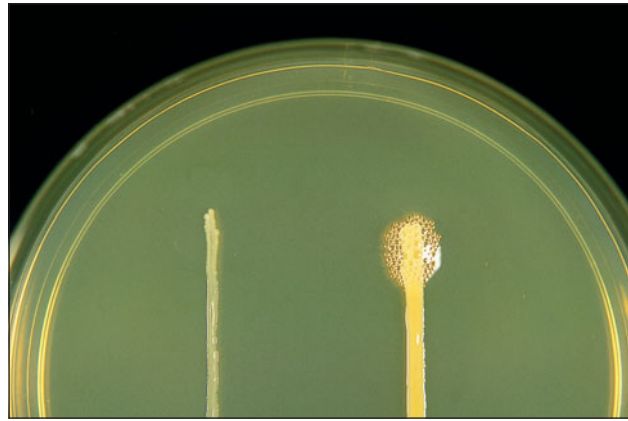(a) Tube method



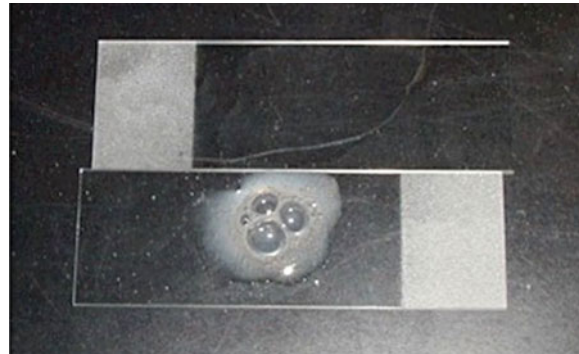(b) Plate method



(c) Slide method

**Figure 30.1 Catalase test.** Negative results are shown on the left and positive results on the right in the (a) tube method and (b) plate method. Negative results are shown on the top and positive results on the bottom in the (c) slide method.

*Slide method:* Bunsen burner, inoculating loop, glassware marking pencil, glass microscope slides (4 for the short version, 14 for the long version), Petri dish and cover.

## Procedure Lab One

### Tube Method

1. Using aseptic technique, inoculate each experimental organism into its appropriately labeled tube by means of a streak inoculation. The last tube will serve as a control.
2. Incubate all cultures for 24 to 48 hours at 37°C.

## Procedure Lab Two

### Tube method

1. Allow three or four drops of the 3% hydrogen peroxide to flow over the entire surface of each slant culture.

2. Examine each culture for the presence or absence of bubbling or foaming. Record your results in the chart in the Lab Report.
3. Based on your observations, determine and record whether or not each organism was capable of catalase activity.

### Slide Method

1. Label slides with the names of the organisms.
2. Using a sterile loop, collect a small sample of the first organism from the culture tube and transfer it to the appropriately labeled slide.
3. Place the slide in the Petri dish.
4. Place one drop of 3% hydrogen peroxide on the sample. Do not mix. Place the cover on the Petri dish to contain any aerosols.
5. Observe for immediate presence of bubble formation. Record your results in the chart in the Lab Report.
6. Repeat Steps 2 through 5 for the remaining test organisms.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

Bacterial Species	PRESENCE OR ABSENCE OF BUBBLING		CATALASE PRODUCTION (+) OR (-)	
	Tube	Slide	Tube	Slide
<i>E. coli</i>				
<i>E. aerogenes</i>				
<i>K. pneumoniae</i>				
<i>S. dysenteriae</i>				
<i>S. typhimurium</i>				
<i>P. vulgaris</i>				
<i>P. aeruginosa</i>				
<i>A. faecalis</i>				
<i>M. luteus</i>				
<i>L. lactis</i>				
<i>S. aureus</i>				
<i>B. cereus</i>				
<i>C. xerosis</i>				
Control				

## Review Questions

1. Explain the toxic effect of O<sub>2</sub> on strict anaerobes.



2. Illustrate the chemical reaction involved in the degradation of hydrogen peroxide in the presence of catalase.

3.  Would catalase be classified as an endoenzyme or an exoenzyme? Explain.

4.  Account for the ability of streptococci to tolerate O<sub>2</sub> in the absence of catalase activity.

## LEARNING OBJECTIVE

Once you have completed this experiment, you should be able to

1. Perform an experimental procedure that is designed to distinguish among groups of bacteria on the basis of cytochrome oxidase activity.

## Principle

Oxidase enzymes play a vital role in the operation of the electron transport system during aerobic respiration. **Cytochrome oxidase** catalyzes the oxidation of a reduced cytochrome by molecular oxygen ( $O_2$ ), resulting in the formation of  $H_2O$  or  $H_2O_2$ . Aerobic bacteria, as well as some facultative anaerobes and microaerophiles, exhibit oxidase activity. The oxidase test aids in differentiation among members of the genera *Neisseria* and *Pseudomonas*, which are oxidase-positive, and Enterobacteriaceae, which are oxidase-negative.

The ability of bacteria to produce cytochrome oxidase can be determined by the addition of the test reagent *p*-aminodimethylaniline oxalate to colonies grown on a plate medium. This light pink reagent serves as an artificial substrate, donating electrons and thereby becoming oxidized to a blackish compound in the presence of the oxidase and free oxygen. Following the addition of the test reagent, the development of pink, then maroon, and finally dark purple coloration on the surface of the colonies is indicative of cytochrome oxidase production and represents a positive test. No color change, or a light pink coloration on the colonies, is indicative of the absence of oxidase activity and is a negative test. The filter paper method may also be used and is described in this experiment.

## CLINICAL APPLICATION

## Test to Distinguish Family Enterobacteriaceae from non-Enterobacteriaceae

The Enterobacteriaceae are cytochrome oxidase-negative, while the genera *Neisseria* and *Pseudomonas* are cytochrome oxidase-positive. The oxidase test is an important tool in the identification of *Neisseria meningitidis*, the causative agent of bacterial meningitis, which has a significant morbidity and mortality rate. In addition, yeast of medical importance such as *Candida* can be separated from *Saccharomyces* and *Torulopsis* by this test.

## AT THE BENCH



## Materials

## Cultures

24- to 48-hour Trypticase soy broth cultures of *Escherichia coli*, *Pseudomonas aeruginosa*, and *Alcaligenes faecalis* for the short version. 24- to 48-hour Trypticase soy broth cultures of the 13 organisms listed on page 154 for the long version.

## Media

Trypticase soy agar plates per designated student group: one for the short version, four for the long version.

## Reagent

*p*-Aminodimethylaniline oxalate (Difco 0329-13-9).

## Equipment

*Plate method:* Bunsen burner, inoculating loop, and glassware marking pencil.

*Filter paper method:* All of the above and filter paper (one for short version, four for long version).

## Procedure Lab One

### Plate Method

1. Prepare the Trypticase soy agar plate(s) for inoculation as follows:
  - a. Short procedure: With a glassware marking pencil, divide the bottom of a Petri dish into three sections and label each section with the name of the test organism to be inoculated.
  - b. Long procedure: Repeat Step 1a, dividing three plates into three sections and one plate into four sections to accommodate the 13 test organisms.
2. Using aseptic technique, make a single-line streak inoculation of each test organism on the agar surface of its appropriate section of the plate(s).
3. Incubate the plate(s) in an inverted position for 24 to 48 hours at 37°C.

## Procedure Lab Two

### Plate Method

1. Add two or three drops of the *p*-aminodimethylaniline oxalate to the surface of the growth of each test organism.
2. Observe the growth for the presence or absence of a color change from pink, to maroon, and finally to purple. Positive test (+), color change in 10–30 seconds; negative test (–), no color change, or light pink color. Refer to **Figure 31.1**. Record the results on the chart in the Lab Report.
3. Based on your observations, determine and record whether or not each organism was capable of producing cytochrome oxidase.



**Figure 31.1** Oxidase test. Negative test, on left, results in no color change, and positive test, on right, results in a color change to purple.

### Filter Paper Method

1. Prepare Petri dishes as described in Lab One Steps 1a and 1b.
2. Place filter paper in Petri dish.
3. With a sterile loop, obtain a heavy loopful of the first test organism and gently smear it on the filter paper.
4. Drop one or two drops of *p*-aminodimethylaniline oxalate reagent on the test organism.
5. Observe the organism for the appearance of a purple color within 30 seconds of contact with the oxidase reagent, indicating a positive test.
6. Repeat Steps 3 to 5 for the remaining test organisms.
7. Record your results in the chart in the Lab Report.

Name: \_\_\_\_\_


Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

Bacterial Species	COLOR OF COLONIES		OXIDASE PRODUCTION (+) OR (-)	
	Plate	Filter paper	Plate	Filter paper
<i>E. coli</i>				
<i>E. aerogenes</i>				
<i>K. pneumoniae</i>				
<i>S. dysenteriae</i>				
<i>S. typhimurium</i>				
<i>P. vulgaris</i>				
<i>P. aeruginosa</i>				
<i>A. faecalis</i>				
<i>M. luteus</i>				
<i>L. lactis</i>				
<i>S. aureus</i>				
<i>B. cereus</i>				
<i>C. xerosis</i>				
Control				

## Review Questions

1. What is the function of cytochrome oxidase?

2. Why are strict aerobes oxidase-positive?
3. The oxidase test is used to differentiate among which groups of bacteria?
4. What is the function of the test reagent in this procedure?
5.  Your instructor asks you to isolate and identify the organisms in an unknown culture. You find that the culture contains two gram-negative bacilli that produce swarming colonies. What biochemical test would you use to identify the bacilli? Justify your answer.

The study of the metabolism of amino acids began in the early part of the twentieth century. Some investigators found that the enteric microorganisms, such as *Proteus* and the so-called Providence species, were able to deaminate a variety of amino acids that provided a vehicle for distinguishing these microorganisms from other members of the large family of the Enterobacteriaceae. It was determined that 11 of the 22 amino acids were deaminated by amino acid oxidases, and it was phenylalanine deaminase that produced the most rapid enzymatic activity. Thus, phenylalanine deaminase became the most widely studied deaminase used to differentiate enteric organisms.

Likewise, some organisms were found to be capable of decarboxylating amino acids, providing a way to differentiate between the enteric genera and species. For instance, lysine decarboxylase is capable of differentiating between *Salmonella* and *Citrobacter*. Ornithine decarboxylase separates *Enterobacter* from *Klebsiella*. Decarboxylase enzymes are numerous, and each is specific for a particular substrate.

It is now evident that decarboxylases and deaminases play a vital role in the utilization of amino acids and the metabolism of nitrogen compounds.

## PART A Decarboxylase Test

### LEARNING OBJECTIVE

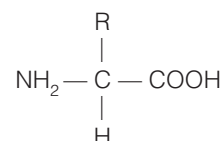
Once you have completed this experiment, you should be able to

1. Identify and differentiate organisms based on their ability to enzymatically degrade amino acid substrates.

## Principle


Every biologically active protein is composed of the 20 essential amino acids. Structurally, amino acids are composed of an alpha carbon ( $-\text{C}-$ ), an amino

group ( $-\text{NH}_2$ ), a carboxyl group ( $-\text{COOH}$ ), and a hydrogen atom ( $-\text{H}$ ). Also attached to the alpha carbon is a side group or an atom designated by an ( $-\text{R}$ ), which differs in each of the amino acids.



**Decarboxylation** is a process whereby some microorganisms that possess decarboxylase enzymes are capable of removing the carboxyl group to yield end products consisting of an **amine** or **diamine** plus **carbon dioxide**. Decarboxylated amino acids play an essential role in cellular metabolism since the amines produced may serve as end products for the synthesis of other molecules required by the cell. Decarboxylase enzymes are designated as adaptive (or induced) enzymes and are produced in the presence of specific amino acid substrates upon which they act. These amino acid substrates must possess at least one chemical group other than an amine ( $-\text{NH}_2$ ) or a carboxyl group ( $-\text{COOH}$ ). In the process of decarboxylation, organisms are cultivated in an acid environment and in the presence of a specific substrate. The decarboxylation end product (amines) results in a shift to a more alkaline pH.

In the clinical or diagnostic microbiology laboratory, three decarboxylase enzymes are used to differentiate members of the Enterobacteriaceae: lysine, ornithine, and arginine. Decarboxylase activity is determined by cultivating the organism in a nutrient medium containing glucose, the specific amino acid substrate, and bromthymol blue (the pH indicator). If decarboxylation occurs, the pH of the medium becomes alkaline despite the fermentation of glucose since the end products (amines or diamines) are alkaline. The function of the glucose in the medium is to ensure good microbial growth and thus more reliable results in the presence of the pH indicator. The presence of each decarboxylase enzyme can be tested for by supplementing decarboxylase broth with the specific

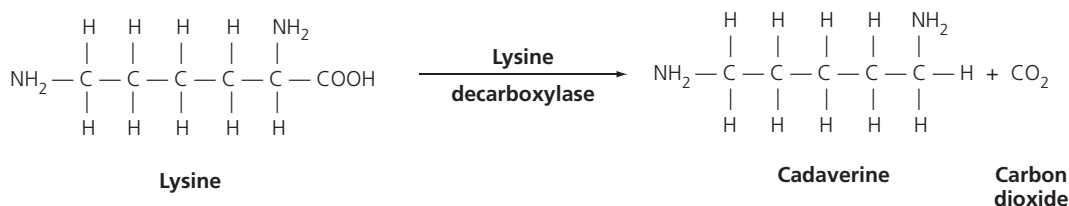

**Figure 32.1** Degradation of lysine

amino acid substrate, namely lysine, arginine, and ornithine. For example, **lysine decarboxylase** degrades L-lysine, forming the diamine end product **cadaverine** plus **carbon dioxide** as illustrated in **Figure 32.1**.

In the experiment that follows, the decarboxylation of L-lysine will be studied. It should be noted that decarboxylation reactions occur under anaerobic conditions that are satisfied by sealing the culture tubes with sterile mineral oil. In the sealed tubes, all of the unbound oxygen is utilized during the organisms' initial growth phase, and the pH of the medium becomes alkaline as carbon dioxide (CO<sub>2</sub>) is produced in the culture tube. A pH indicator such as bromocresol purple is usually incorporated into the medium for the easy detection of pH changes. The production of acid end products will cause the bromocresol purple to change color from purple to yellow, indicating that acid has formed, the medium has been acidified, and the decarboxylase enzymes have been activated. The activated enzyme responds with the production of the alkalizing diamine (cadaverine) and carbon dioxide, which will produce a final color change from yellow back to purple, thereby indicating that L-lysine has been decarboxylated. The development of a turbid purple color verifies a positive test for amino acid decarboxylation. The absence of a purple color indicates a negative result.

## CLINICAL APPLICATION

### Distinguishing between *Enterobacter* Species

The decarboxylase test identifies bacteria based on the production of ammonia from the amino acids lysine, ornithine, and arginine. The decarboxylase test can be used to differentiate the causative agent in many nosocomial infections of immunocompromised patients. The bacterium *Enterobacter aerogenes* is lysine decarboxylase positive while other *Enterobacter* species are negative. It is used primarily to identify bacteria within the Enterobacteriaceae family.

## AT THE BENCH



## Materials

### Cultures

24-hour nutrient broth cultures of *Proteus vulgaris*, *Escherichia coli*, and *Citrobacter freundii* for the short version. 24- to 48-hour nutrient broth cultures of the 13 organisms listed on page 154 for the long version.

### Media


Per designated student group: three tubes of Moeller's decarboxylase broth supplemented with L-lysine (10 gm/l) (labeled LD+), three tubes of Moeller's decarboxylase broth without lysine (labeled LD-).

### Equipment

Bunsen burner, glassware marking pencil, inoculating loop and needle, sterile Pasteur pipettes, rubber bulbs, test tube rack, and sterile mineral oil.

## Procedure Lab One

1. With a glassware marking pencil, label three tubes of the LD+ medium with the name of the organism to be inoculated. Similarly label three tubes of LD- medium. The use of (LD-) control tubes is essential since some bacterial strains are capable of turning substrate-free media positive. *Note: Control tubes should remain yellow after incubation, denoting that only glucose was fermented.* The presence of a positive control tube invalidates the test, and no interpretation is possible.
2. Using aseptic technique, inoculate each experimental organism into its appropriately labeled tube using a loop inoculation.
3. Place a rubber bulb onto a sterile Pasteur pipette and overlay the surface of the inoculated

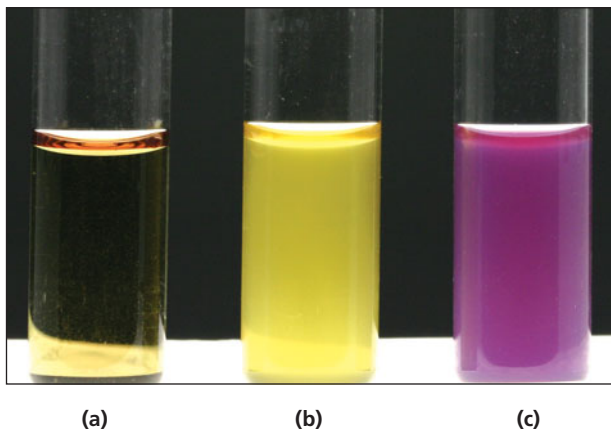

**Figure 32.2** Decarboxylase test. (a) Uninoculated, (b) negative, and (c) positive.

culture tubes with 1 ml of sterile mineral oil. Hold the tubes in a slanted position while adding the mineral oil. *Note: Do not let the tip of the pipette touch the inoculated medium or the sides of the test tube walls.*

- Repeat the above procedure for the remaining test cultures.
- Incubate all tubes at 37°C for 24 to 48 hours.

## Procedure Lab Two

- Examine each culture tube for the presence of a color change. Refer to **Figure 32.2**.
- Based on your observations, determine whether or not each organism was capable of performing decarboxylation of lysine.
- Record your results in the chart in the Lab Report.

## PART B Phenylalanine Deaminase Test

### LEARNING OBJECTIVE

Once you have completed this experiment, you should be able to

- Demonstrate the ability of some organisms to remove the amino group ( $-\text{NH}_2$ ) from amino acids.

## Principle

Microorganisms that contain deaminase enzymes are capable of removing the amino group ( $-\text{NH}_2$ ) from amino acids and other  $\text{NH}_2$ -containing chemical compounds. During this process the amino acid, under the auspices of its specific deaminase, will produce keto acids and ammonia as end products. In the experiment to follow, the amino acid phenylalanine will be deaminated by **phenylalanine deaminase** and converted to the **keto acid phenylpyruvic acid** and **ammonia**. The organisms are cultured on a medium incorporating phenylalanine as the substrate. This chemical reaction is illustrated in **Figure 32.3**.


If the organism possesses phenylalanine deaminase, phenylpyruvic acid will be released into the medium and can be detected by the addition of a 10 to 12% ferric chloride solution to the surface of the medium. If a green color develops, the enzymatic deamination of the substrate has occurred and is indicative of a positive result. The absence of any color change indicates a negative result. The resultant green color produced upon the addition of ferric chloride ( $\text{FeCl}_3$ ) is due to the formation of a keto acid (phenylpyruvic acid). It has been shown that  $\alpha$ - and  $\beta$ -keto acids give a positive color reaction with either alcoholic or aqueous solutions of  $\text{FeCl}_3$ . Phenylpyruvic acid is an  $\alpha$ -keto acid. The results should be read immediately following the addition of the reagent since the color produced fades quickly. When not in use, the ferric chloride reagent should be refrigerated and kept in a dark bottle to avoid exposure to light. The stability of this reagent varies and should be checked weekly with known positive cultures.

## CLINICAL APPLICATION

### Differentiating Intestinal Bacteria

The phenylalanine deaminase test uses the differential medium phenylalanine agar to detect bacteria containing the enzyme phenylalanine deaminase, and is used to differentiate the genera *Proteus*, *Morganella*, and *Providencia* from other gram-negative intestinal bacilli. These genera of enteric and environmental bacteria are known to cause UTIs and gastroenteritis. It is clinically important to distinguish them from other enteric bacteria due to their high level of antibiotic resistance.



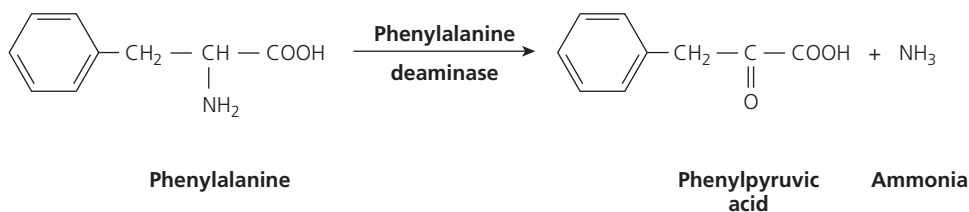

**Figure 32.3** Deamination of phenylalanine

## AT THE BENCH



## Materials

### Cultures

24-hour nutrient broth cultures of *Escherichia coli* and *Proteus vulgaris* for the short version. 24-hour nutrient broth cultures of the 13 organisms listed on page 154 for the long version.

### Media

Two phenylalanine agar slants.

### Reagents

10 to 12% ferric chloride solution.

### Equipment

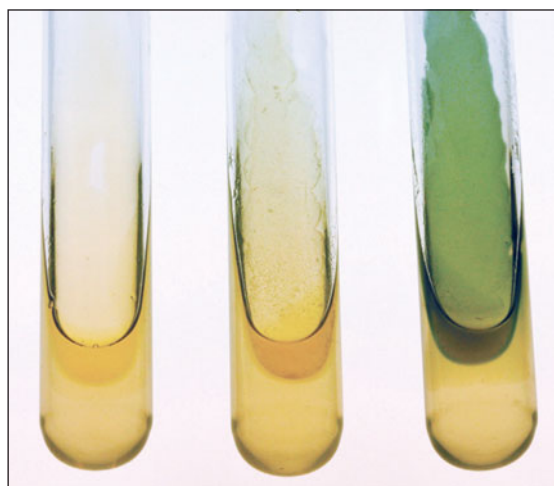
Bunsen burner, glassware marking pencil, Pasteur pipettes, rubber bulbs, test tube racks, and inoculating loop.

## Procedure Lab One

1. Using aseptic technique, inoculate each experimental organism into its appropriately labeled tube using a streak inoculation.
2. Incubate cultures at 37°C for 24 to 48 hours.

## Procedure Lab Two

1. Add 5 to 10 drops of the ferric chloride solution to each agar slant and mix gently. Ferric chloride is a chelating agent and binds to the phenylpyruvic acid to produce a green color on the slant (**Figure 32.4**).



(a) (b) (c)

**Figure 32.4** Phenylalanine deaminase test. (a) Uninoculated, (b) negative, and (c) positive.

2. Based on your observations, determine whether or not each organism was capable of amino acid deamination. *Note: Results should be read immediately following the addition of ferric chloride because the green color fades rapidly.*
3. Record your results in the Lab Report.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_


## Observations and Results

### Parts A, B: Decarboxylase Test and Phenylalanine Deaminase Test

Bacterial Species	DECARBOXYLASE TEST				PHENYLALANINE DEAMINASE TEST	
	Color of Medium		Lysine Decarboxylase (+) or (-)		Color after FeCl <sub>3</sub>	Deamination (+) or (-)
	LD+	LD-	LD+	LD-		
<i>E. coli</i>						
<i>E. aerogenes</i>						
<i>K. pneumoniae</i>						
<i>S. dysenteriae</i>						
<i>S. typhimurium</i>						
<i>P. vulgaris</i>						
<i>P. aeruginosa</i>						
<i>A. faecalis</i>						
<i>M. luteus</i>						
<i>L. lactis</i>						
<i>S. aureus</i>						
<i>B. cereus</i>						
<i>C. xerosis</i>						
<i>C. freundii</i>						
Control						

## Review Questions

1. A negative decarboxylase test is indicated by the production of a yellow color in the medium. Explain the reason for the development of this color.

2. Explain why deaminase activity must be determined immediately following the addition of ferric chloride.
  
3. What is the function of ferric chloride in the detection of deaminase activity?
  
  
  
  
  
  
  
  
  
  
4. Explain why the anaerobic environment is essential for decarboxylation of the substrate to occur.
  
  
  
  
  
  
  
  
  
  
5.  Following a normal delivery, a nurse observes that the urine of the infant has a peculiar odor resembling that of burnt sugar or maple syrup. Subsequent examination by the pediatrician reveals that this child has maple syrup urine disease.
  - a. What is this disease?
  
  
  
  
  
  
  
  
  
  
  - b. How is it treated?

## LEARNING OBJECTIVE

Once you have completed this experiment, you should be able to

1. Use previously studied staining, cultural characteristics, and biochemical procedures for independent genus identification of an unknown bacterial culture.

## Principle

Identification of unknown bacterial cultures is one of the major responsibilities of the microbiologist. Samples of blood, tissue, food, water, and cosmetics are examined daily in laboratories throughout the world for the presence of contaminants. In addition, industrial organizations are constantly screening materials to isolate new antibiotic-producing organisms or organisms that will increase the yield of marketable products such as vitamins, solvents, and enzymes. Once isolated, these unknown organisms must be identified and classified.

The science of classification is called **taxonomy** and deals with the separation of living organisms into interrelated groups. *Bergey's Manual* has been the official, internationally accepted reference for bacterial classification since 1923. The current edition, *Bergey's Manual of Systematic Bacteriology*, arranges related bacteria into 33 groups called sections rather than into the classical taxonomic groupings of phylum, class, order, and family. The interrelationship of the organisms in each section is based on characteristics such as morphology, staining reactions, nutrition, cultural characteristics, physiology, cellular chemistry, and biochemical test results for specific metabolic end products.

At this point you have developed sufficient knowledge of staining methods, isolation techniques, microbial nutrition, biochemical activities, and characteristics of microorganisms to be able to work independently in attempting to identify

the genus of an unknown culture. Characteristics of the major organisms that have been used in experiments thus far are given in Table 33.1. You are to use this table for the identification of the unknown cultures. The observations and results obtained following the experimental procedures are the basis of this identification. *However, you should note that your biochemical results may not be identical to those shown in Table 33.1; they may vary because of variations in bacterial strains (subgroups of a species).* Therefore, it becomes imperative to recall the specific biochemical tests that differentiate among the different genera of the test organisms.

Identification of an unknown culture using a more extensive procedure to differentiate bacterial species is presented in Experiment 70. The rationale for the performance of this exercise later in the semester is twofold. First, you will have acquired expanded knowledge of microbial activities and will be more proficient in laboratory skills. Second, and more important, you will be more cognizant of and more critical in your approach to species identification using dichotomous keys supplemented with *Bergey's Manual*.

## CLINICAL APPLICATION

### Application of Learned Assays to Identify an Unknown Bacterial Pathogen

The role of the clinical laboratory in a hospital is to quickly and efficiently identify the causative agent of a patient's infection. This will entail choosing the correct assays and performing them in the correct order to logically identify the genus and species of the agent.

**TABLE 33.1 Cultural and Biochemical Characteristics of Unknown Organisms**

ORGANISM	GRAM STAIN	AGAR SLANT CULTURAL CHARACTERISTICS	LITMUS MILK REACTION	FERMENTATION														
				LACTOSE	DEXTRROSE	SUCROSE	H <sub>2</sub> S PRODUCTION	NO <sub>3</sub> REDUCTION	INDOLE PRODUCTION	MR REACTION	VP REACTION	CITRATE USE	UREASE ACTIVITY	CATALASE ACTIVITY	OXIDASE ACTIVITY	GELATIN LIQUEFACTION	STARCH HYDROLYSIS	LIPID HYDROLYSIS
<i>Escherichia coli</i>	Rod –	White, moist, glistening growth	Acid, curd ±, gas ±, reduction ±	AG	AG	A±	–	+	+	+	–	–	–	+	–	–	–	–
<i>Enterobacter aerogenes</i>	Rod –	Abundant, thick, white, glistening growth	Acid	AG	AG	AG±	–	+	–	–	+	+	–	+	–	–	–	–
<i>Klebsiella pneumoniae</i>	Rod –	Slimy, white, somewhat translucent, raised growth	Acid, gas, curd ±	AG	AG	AG	–	+	–	±	±	+	+	+	–	–	–	–
<i>Shigella dysenteriae</i>	Rod –	Thin, even, grayish growth	Alkaline	–	A	A±	–	+	±	+	–	–	–	+	–	–	–	–
<i>Salmonella typhimurium</i>	Rod –	Thin, even, grayish growth	Alkaline	–	AG±	A±	+	+	–	+	–	+	–	+	–	–	–	–
<i>Proteus vulgaris</i>	Rod –	Thin, blue-gray, spreading growth	Alkaline	–	AG	AG±	+	+	+	+	–	±	+	+	–	+	–	–
<i>Pseudomonas aeruginosa</i>	Rod –	Abundant, thin, white growth, with medium turning green	Rapid peptonization	–	–	–	–	+	–	–	–	+	–	+	+	+	+	+
<i>Alcaligenes faecalis</i>	Rod* –	Thin, white, spreading, viscous growth	Alkaline	–	–	–	–	–	–	–	–	±	–	+	+	–	–	–
<i>Staphylococcus aureus</i>	Cocci +	Abundant, opaque, golden growth	Acid reduction ±	A	A	A	–	+	–	+	±	–	–	+	–	+	–	+
<i>Lactococcus lactis</i>	Cocci +	Thin, even growth	Acid, rapid reduction with curd	A	A	A	–	–	–	+	–	–	–	–	–	–	–	–
<i>Micrococcus luteus</i>	Cocci +	Soft, smooth, yellow growth	Alkaline	–	–	–	–	±	–	–	–	–	+	+	–	+	+	–
<i>Corynebacterium xerosis</i>	Rod +	Grayish, granular, limited growth	Alkaline	–	A±	A±	–	+	–	–	–	–	–	+	–	–	–	–
<i>Bacillus cereus</i>	Rod +	Abundant, opaque, white waxy growth	Peptonization	–	A	A	–	+	–	–	±	–	–	+	–	+	+	±

Note: AG = Acid and gas; ± = Variable reaction; Rod\* = Coccobacillus



## Materials

### Cultures

Number-coded 24- to 48-hour Trypticase soy agar slant cultures of the 13 bacterial species listed on page 154. You will be provided with one unknown pure culture.

### Media

Two Trypticase soy agar slants, and one each of the following per student: phenol red sucrose broth, phenol red lactose broth, phenol red dextrose broth, SIM agar deep tube, MR-VP broth, tryptic nitrate broth, Simmons citrate agar slant, urea broth, litmus milk, Trypticase soy agar plate, nutrient gelatin deep tube, starch agar plate, and tributyrin agar plate.

### Reagents

Crystal violet; Gram's iodine; 95% ethyl alcohol; safranin; methyl red; 3% hydrogen peroxide; Barritt's reagent, Solutions A and B; Kovac's reagent; zinc powder; and p-aminodimethylaniline oxalate.

### Equipment

Bunsen burner, inoculating loop and needle, staining tray, immersion oil, lens paper, bibulous paper, microscope, and glassware marking pencil.

## Procedure Lab One

1. Perform a Gram stain of the unknown organism. Observe and record in the Lab Report chart the reaction and the morphology and arrangement of the cells.
2. Using aseptic inoculating technique, inoculate two Trypticase soy agar slants by means of a streak inoculation. Following incubation, you will use one slant culture to determine the cultural characteristics of the unknown microorganism. You will use the second as a stock subculture should it be necessary to repeat any of the tests.

3. Exercising care in aseptic technique so as not to contaminate cultures and thereby obtain spurious results, inoculate the media for the following biochemical tests:

Medium	Test
a. Phenol red lactose broth b. Phenol red dextrose broth c. Phenol red sucrose broth	Carbohydrate fermentationc
d. Litmus milk	Litmus milk reactions
e. SIM medium	Indole production H <sub>2</sub> S production
f. Tryptic nitrate broth	Nitrate reduction
g. MR-VP broth	Methyl red test Voges-Proskauer test
h. Simmons citrate agar slant	Citrate utilization
i. Urea broth	Urease activity
j. Trypticase soy agar slant	Catalase activity
k. Starch agar plate	Starch hydrolysis
l. Tributyrin agar plate	Lipid hydrolysis
m. Nutrient gelatin deep tube	Gelatin liquefaction
n. Trypticase soy agar plate	Oxidase test

4. Incubate all cultures for 24 to 72 hours at 37°C.

## Procedure Lab Two

1. Examine a Trypticase soy agar slant culture and determine the cultural characteristics of your unknown organism. Record your results in the Lab Report.
2. Perform biochemical tests on the remaining cultures, making reference to the specific laboratory exercise for each test. Record your observations and results.
3. Based on your results, identify the genus and species of the unknown organism. *Note: Results may vary depending on the strains of each species used and the length of time the organism has been maintained in stock culture. The observed results may not be identical to the expected results. Therefore choose the organism that best fits the results summarized in Table 33.1.*

*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

Description of Unknown's Characteristics		Student
		_____
		Culture no.
		_____
		Organism
		_____
Experimental Procedure	Observations	Results
Gram stain		
Acid-fast stain		
Shape and arrangement		
Cultural characteristics		
Litmus milk reactions		
Carbohydrate fermentations: Lactose		
Dextrose		
Sucrose		
H <sub>2</sub> S production		
Nitrate reduction		
Indole production		
Methyl red test		
Voges-Proskauer test		
Citrate utilization		
Urease activity		
Catalase activity		
Starch hydrolysis		
Lipid hydrolysis		
Gelatin liquefaction		
Oxidase test		



*This page intentionally left blank*

# The Protozoa

## LEARNING OBJECTIVES

Once you have completed the experiments in this section, you should be

1. Familiar with the distinguishing characteristics of protozoans.
2. Able to identify free-living and parasitic protozoans in microscopic views.

## Introduction

The protozoa are a large and diverse group of unicellular, eukaryotic organisms. Most are free-living, but some are parasites. Their major distinguishing characteristics are

1. The absence of a cell wall; some, however, possess a flexible layer, a pellicle, or a rigid shell of inorganic materials outside of the cell membrane.
2. The ability during their entire life cycle or part of it to move by locomotor organelles or by a gliding mechanism.
3. Heterotrophic nutrition whereby the free-living forms ingest particulates such as bacteria, yeast, and algae, while the parasitic forms derive nutrients from the body fluids of their hosts.
4. Primarily asexual means of reproduction, although sexual modes occur in some groups.

Protozoan taxonomy is being continually updated as new technology enables classification based on molecular characteristics. To ease our discussion of protozoans, we follow a more traditional taxonomic scheme, dividing them into four groups based on means of locomotion.

1. **Sarcodina:** Motility results from the streaming of ectoplasm, producing protoplasmic projections called pseudopods (false feet).

Prototypic amoebas include the free-living *Amoeba proteus* and the parasite *Entamoeba histolytica*.

2. **Mastigophora:** Locomotion is effected by one or more whiplike, thin structures called flagella. Free-living members include the genera *Cercomonas*, *Heteronema*, and *Euglena*, which are photosynthetic protists that may be classified as flagellated algae. The parasitic forms include *Trichomonas vaginalis*, *Giardia intestinalis* (formerly called *Giardia lamblia*), and the *Trypanosoma* species.
3. **Ciliophora:** Locomotion is carried out by means of short hairlike projections called cilia, whose synchronous beating propels the organisms. The characteristic example of free-living members of this group is *Paramecium caudatum*, and the parasitic example is *Balantidium coli*.
4. **Sporozoa:** Unlike other members of this phylum, sporozoa do not have locomotor organelles in their mature stage; however, immature forms exhibit some type of movement. All the members of this group are parasites. The most significant members belong to the genus *Plasmodium*, the malarial parasites of animals and humans.

*This page intentionally left blank*

### LEARNING OBJECTIVE

Once you have completed this experiment, you should be familiar with

1. The protozoa of pond water.

## Principle

There are more than 20,000 known species of free-living protozoa. It is not within the scope of this manual to present an in-depth study of this large and diverse population. Therefore in this procedure you will use Table 34.1 and **Figures 34.1** and **Figure 34.2** to become familiar with the general structural characteristics of representative protozoa, and you will identify these in a sample of pond water.

### CLINICAL APPLICATION

#### Wet Mounts for Diagnosis

Wet mount slides, often utilizing stains, are routinely used in the examination of stool samples for infectious protozoans such as *Entamoeba histolytica*. This organism causes amoebic dysentery and has been known to lead to severe liver damage. Although most infections are asymptomatic, carriers can still spread the disease. Diagnosis may require examination of several slide preparations.

### AT THE BENCH



## Materials

### Cultures

Stagnant pond water and prepared slides of amoebas, paramecia, euglenas, and stentors.

### Reagent

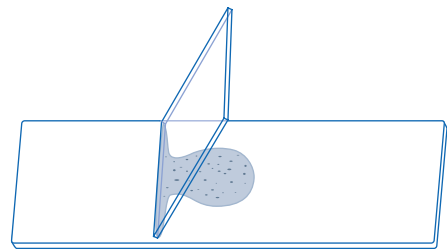
Methyl cellulose.

### Equipment


Microscope, glass slides, coverslips, and Pasteur pipettes.

## Procedure

1. Obtain a drop of pond water from the bottom of the culture and place it in the center of a clean slide.
2. Add a drop of methyl cellulose to the culture to slow down the movement of the protozoa.
3. Apply a coverslip in the following manner to prevent formation of air bubbles:
  - a. Place one edge of the coverslip against the outer edge of the drop of culture.



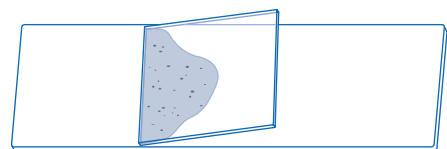
- b. After the drop of culture spreads along the inner aspect of the edge of the coverslip, gently lower the coverslip onto the slide.



4. Examine your slide preparation under scanning, low-power, and high-power objectives with diminished light, and observe for the different protozoa present. Record your results in the Lab Report.


**TABLE 34.1 Structural Characteristics of Free-Living Protozoa**

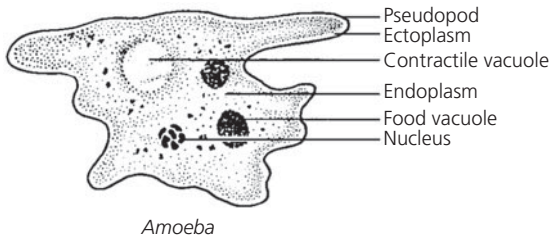Sarcodina




1. **Pseudopods:** Protoplasmic projections that function for locomotion
2. **Ectoplasm:** Outer layer of cytoplasm; clear in appearance
3. **Endoplasm:** Inner cytoplasmic region; granular in appearance
4. **Nucleus:** One present
5. **Food vacuoles:** Contain engulfed food undergoing digestion
6. **Contractile vacuole:** Large, clear circular structure that regulates internal water pressure

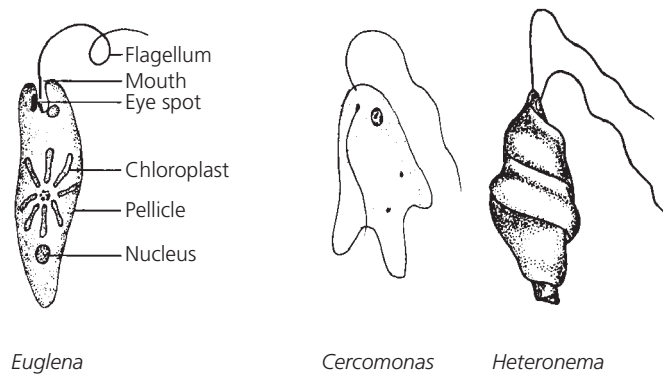Mastigophora




1. **Flagella:** One to several long whiplike structures that function for locomotion
2. **Pellicle:** Elastic layer outside of cell membrane
3. **Mouth:** Present but indistinct
4. **Chloroplast:** Organelles containing chlorophyll; present in photosynthetic forms only
5. **Eye spot:** Light-sensitive pigmented spot
6. **Nucleus:** One present

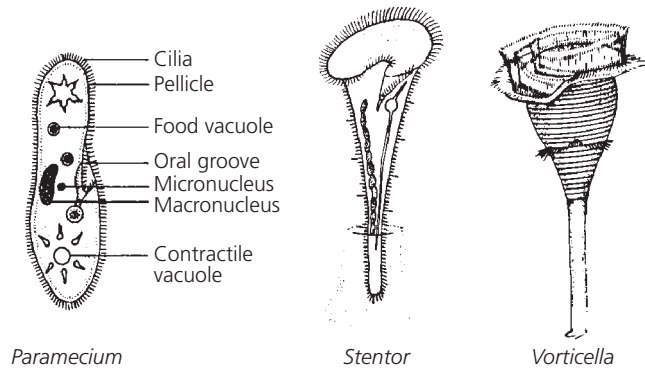Ciliophora



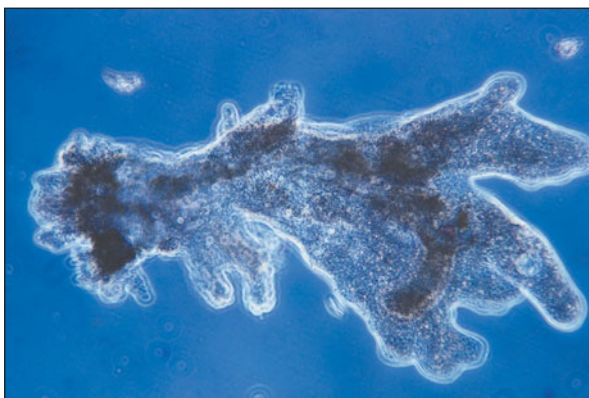1. **Cilia:** Numerous, short, hairlike structures that function for locomotion
2. **Pellicle:** Outermost flexible layer
3. **Contractile vacuole** with radiating canals; regulates osmotic pressure
4. **Oral groove:** Indentation that leads to the mouth and gullet
5. **Food vacuoles:** Sites of digestion of ingested food
6. **Macronucleus:** A large nucleus that functions to control the cell's activities; one to several may be present
7. **Micronucleus:** A small nucleus that functions in conjugation, a mode of sexual reproduction



(a) *Euglena viridis*

(b) *Paramecium caudatum*

**Figure 34.1 Amoeba**



**Figure 34.2 Euglena and Paramecium**

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

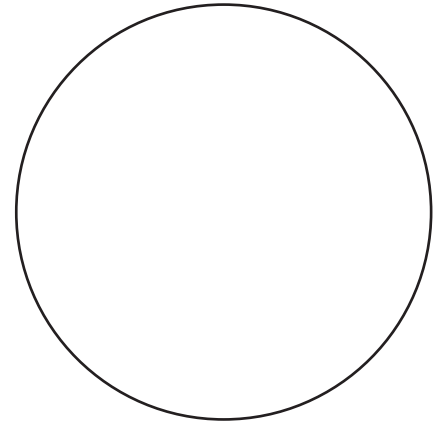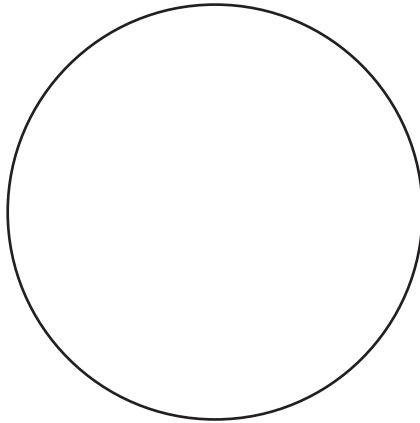1. In the space provided, draw a representative sketch of several of the observed protozoa in stagnant pond water, indicate the magnifications used, and label their structural components. Identify each organism according to its class based on its mode of locomotion and its genus.



Magnification: \_\_\_\_\_

\_\_\_\_\_

Organelles of locomotion: \_\_\_\_\_

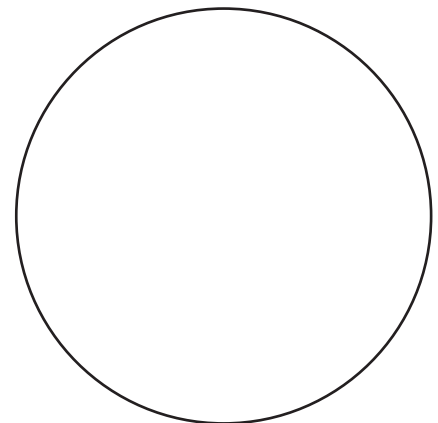
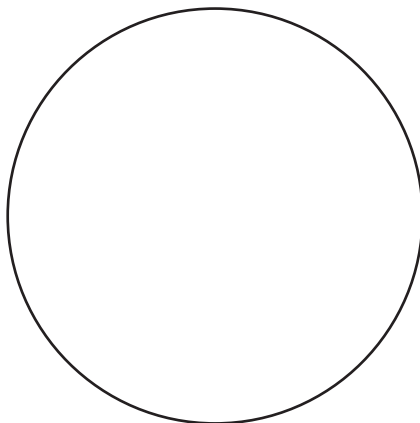
\_\_\_\_\_

Class: \_\_\_\_\_

\_\_\_\_\_

Genus: \_\_\_\_\_

\_\_\_\_\_



Magnification: \_\_\_\_\_

\_\_\_\_\_

Organelles of locomotion: \_\_\_\_\_

\_\_\_\_\_

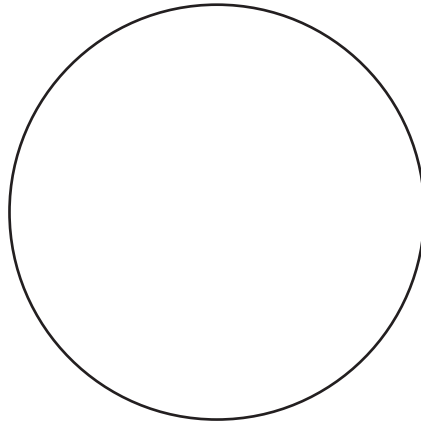
Class: \_\_\_\_\_

\_\_\_\_\_

Genus: \_\_\_\_\_

\_\_\_\_\_

2. Draw representative sketches, indicate magnification, and label the structural components. Identify each organism according to its class based on locomotion and genus.



Amoeba

Magnification:

\_\_\_\_\_

Organelles of locomotion:

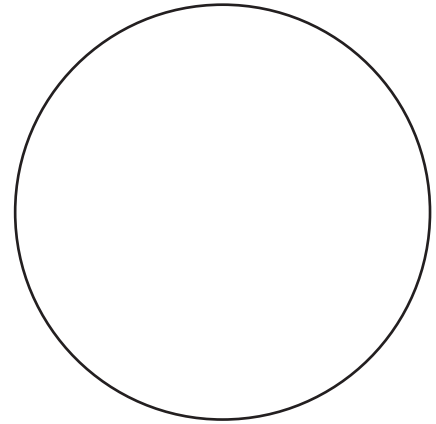
\_\_\_\_\_

Class:

\_\_\_\_\_

Genus:

\_\_\_\_\_




Paramecium

\_\_\_\_\_

\_\_\_\_\_

\_\_\_\_\_

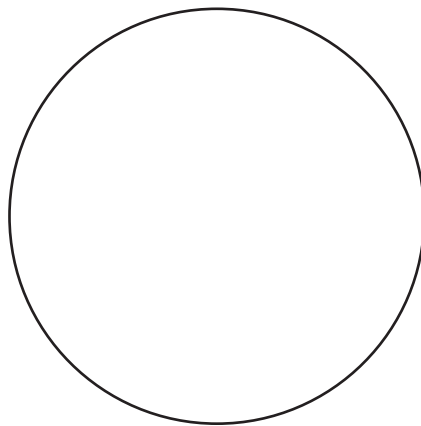\_\_\_\_\_



Euglena

Magnification:

\_\_\_\_\_

Organelles of locomotion:

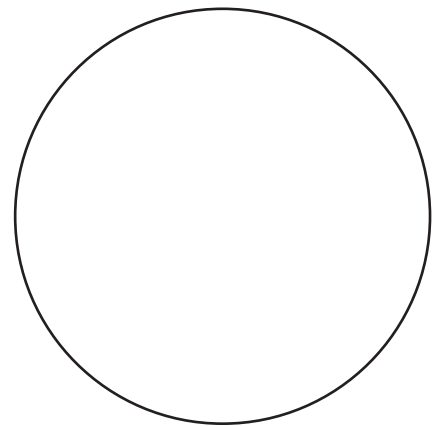
\_\_\_\_\_

Class:

\_\_\_\_\_

Genus:

\_\_\_\_\_



Stentor

\_\_\_\_\_


\_\_\_\_\_

\_\_\_\_\_

\_\_\_\_\_





3.  People with AIDS are vulnerable to toxoplasmosis caused by the protozoan *Toxoplasma gondii*, resulting in infection of lungs, liver, heart, and brain, and often leading to death. About 25% of the world's population is infected, usually without developing symptoms. Why then are people with AIDS so susceptible to this disease?

### LEARNING OBJECTIVE

Once you have completed this experiment, you should be familiar with

1. Parasitic protozoan forms.

## Principle

Unlike the life cycles of the free-living forms, the life cycles of parasitic protozoa vary greatly in complexity. Knowledge of the various developmental stages in these life cycles is essential in the diagnosis, clinical management, and chemotherapy of parasitic infections.

Parasites with the simplest or most direct life cycles not requiring an intermediate host are the following:

1. *Entamoeba histolytica*: A pseudopodian parasite of the class Sarcodina that causes amebic dysentery. Infective, resistant cysts are released from the lumen of the intestine through the feces and are deposited in water, in soil, or on vegetation. On ingestion, the mature quadrinucleated cyst wall disintegrates and the nuclei divide, producing eight active trophozoites (metabolically active cells) that move to the colon, where they establish infection.
2. *Balantidium coli*: The ciliated parasitic protozoan exhibits a life cycle similar to that of *E. histolytica* except that no multiplication occurs within the cyst. This organism resides primarily in the lumen and submucosa of the large intestine. It causes intestinal ulceration and alternating constipation and diarrhea.
3. *Giardia intestinalis*: The intestinal mastigophoric flagellate exhibits a life cycle

comparable to those of the above parasites. This organism is responsible for the induction of abdominal discomfort and severe diarrhea. Diagnosis is made by finding cysts in the formed stool and both cysts and trophozoites in the diarrhetic stool.

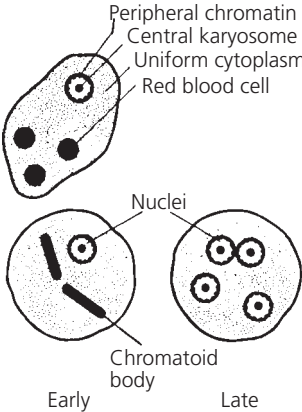
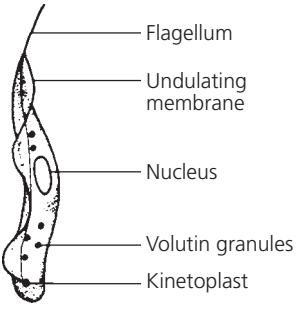
The mastigophoric hemoflagellate responsible for various forms of African sleeping sickness has a more complex life cycle. The *Trypanosoma* must have two hosts to complete its cyclic development: a vertebrate and an invertebrate, blood-sucking insect host. Humans are the definitive hosts harboring the sexually mature forms; the tsetse fly (*Glossina*) and the reduviid bug are the invertebrate hosts in which the developmental forms occur.

Table 35.1 illustrates the morphological characteristics of prototypic members of the parasitic protozoa except the Sporozoa.

Protozoa demonstrating the greatest degree of cyclic complexity are found in the class Sporozoa. They are composed of exclusively obligate parasitic forms, such as members of the genus *Plasmodium*, and are responsible for malaria in both humans and animals. The life cycle requires two hosts, a human being and the female *Anopheles* mosquito. It is significant to note that in this life cycle, the mosquito, and not the human, is the definitive host harboring the sexually mature parasite.

Malaria is initiated when a person is bitten by an infected mosquito, during which time infective sexually mature sporozoites are injected with the insect's saliva. These parasites pass rapidly from the blood into the liver, where they infect the parenchymal cells. This is the **pre-erythrocytic stage**. The parasites develop asexually within the liver cells by a process called **schizogony**, producing **merozoites**. This cycle may be repeated or the merozoites that are released from the ruptured liver cells may now infect red blood cells

**TABLE 35.1** Characteristics of Representative Parasitic Protozoa



CLASS, ORGANISM, AND INFECTION	STRUCTURAL CHARACTERISTICS	LOCOMOTOR ORGANELLES	SITE OF INFECTION	ISOLATION OF PARASITIC FORM
<p>SARCODINA <i>Entamoeba histolytica</i> Infection: Amebic dysentery</p> 	<p><b>Trophozoite:</b> Shape: Variable Nucleus: Discrete nuclear membrane with central karyosome and peripheral chromatin granules Cytoplasm: Clear, red blood cells may be present</p> <p><b>Cyst:</b> Shape: Round to oval with thick wall Nuclei: 1–4 present; mature cyst is quadrinucleated Chromatoid bodies: Sausage-shaped with rounded ends, present in young cysts only</p>	<p>Pseudopods</p> <p>None</p>	<p>Large intestine by ingestion of mature cysts</p> <p>Formed stool</p>	<p>Diarrhetic stool</p> <p>Formed stool</p>
<p>MASTIGOPHORA <i>Trypanosoma gambiense</i> Infection: African sleeping sickness</p> 	<p><b>Trophozoite:</b> Shape: Crescent Nucleus: Large, central, and polymorphic Cytoplasm: Granular</p> <p><b>Cyst:</b> None</p>	<p>Single flagellum along undulating membrane</p>	<p>Peripheral blood-stream by means of tsetse fly vector</p>	<p>Peripheral blood</p>

and initiate the **erythrocytic stage**. During this asexual development, the parasite undergoes a series of morphological changes that are of diagnostic value. These forms are designated as **signet rings, trophozoites, schizonts, segmenters, merozoites, and gametocytes**. The merozoites are capable of reinfesting other blood cells or liver cells. Ingestion of the **microgametocytes** (♂) and **macrogametocytes** (♀) by another mosquito during a blood meal initiates the sexual cycle called **sporogamy**. Male and female gametes give

rise to a zygote in the insect's gut. The zygote is then transformed into an **ookinete** that burrows through the gut wall to form an **oocyst** in which the sexually mature **sporozoites** develop, thereby completing the life cycle.

In this experiment, you will study the parasitic protozoa by using prepared slides and the diagnostic characteristics shown in **Figure 35.1** on page 232 and Table 35.1. The purpose of the experiment is to help you understand life cycles of parasitic protozoa.

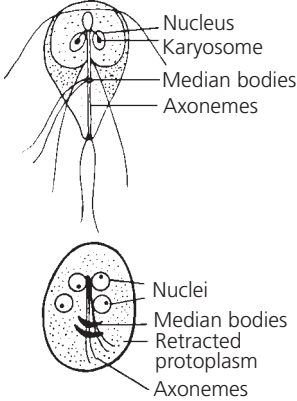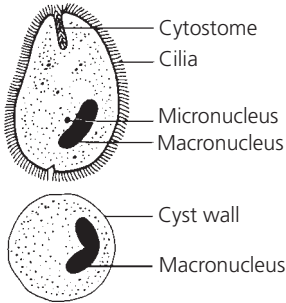**TABLE 35.1 Characteristics of Representative Parasitic Protozoa (continued)**

CLASS, ORGANISM, AND INFECTION	STRUCTURAL CHARACTERISTICS	LOCOMOTOR ORGANELLES	SITE OF INFECTION	ISOLATION OF PARASITIC FORM
<p>MASTIGOPHORA <i>Giardia intestinalis</i> Infection: Dysentery</p> 	<p><b>Trophozoite:</b> Shape: Pear-shaped with concave sucking disc Nuclei: 2 bilaterally located with central karyosome and no peripheral chromatin Cytoplasm: Uniform and clear</p> <p><b>Cyst:</b> Shape: Oval to ellipsoidal Nuclei: 2–4 present and protoplasm retracted from cyst wall Axostyle Parabasal body</p>	<p>4 pairs of flagella</p> <p>4 pairs of flagella within cyst</p>	<p>Small intestine through ingestion of cysts</p>	<p>Diarrhetic stool</p> <p>Formed stool</p>
<p>CILIOPHORA <i>Balantidium coli</i> Infection: Dysentery</p> 	<p><b>Trophozoite:</b> Shape: Oval Nuclei: Kidney-shaped macronucleus and a micronucleus Cytoplasm: Vacuolated</p> <p><b>Cyst:</b> Shape: Round and thick-walled Nuclei: 1 macronucleus and a micronucleus that is not visible</p>	<p>Cilia</p> <p>Cilia within cyst</p>	<p>Large intestine by the ingestion of cysts</p>	<p>Diarrhetic stool</p> <p>Formed stool</p>

## CLINICAL APPLICATION

### Understanding Parasitic Protozoa

Parasitic protozoa can exist extracellularly or intracellularly and possess diverse morphologies. They rapidly reproduce, asexually or sexually, with short generation times. They are highly organ, tissue, or cell specific. Examples are *Plasmodium* species, which colonize red blood cells (malaria); *Trichomonas*, which colonize the urinary tract (vaginal infections); and *Entamoeba*, which colonizes the large intestine (severe diarrhea).


## AT THE BENCH



## Materials

### Prepared Slides

*Entamoeba histolytica* trophozoite and cyst, *Giardia intestinalis* trophozoite and cyst (formerly *G. lamblia*), *Balantidium coli* trophozoite and cyst, *Trypanosoma gambiense*, and *Plasmodium vivax* in human blood smears.

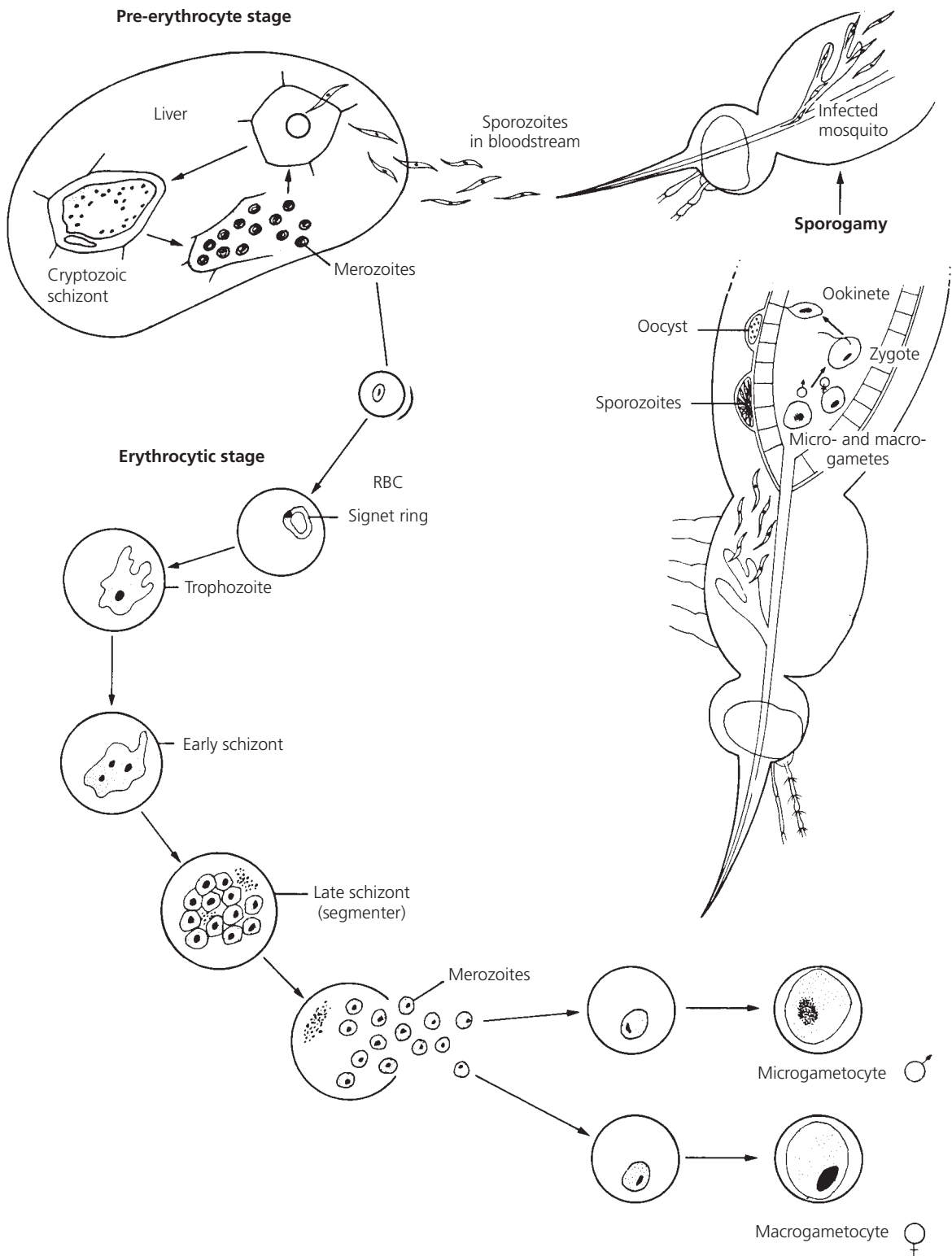

**Figure 35.1** Life cycle of *Plasmodium vivax*

## Equipment

Microscope, immersion oil, and lens paper.

## Procedure

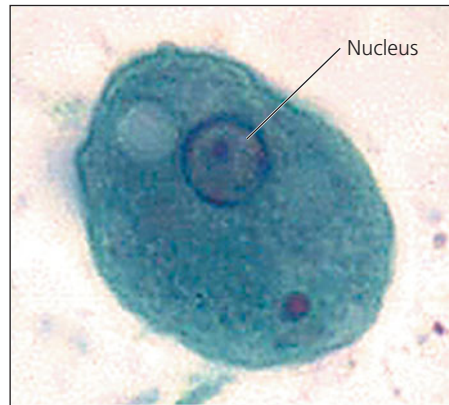
1. Examine all available slides under the oil-immersion objective. Use Table 35.1, Figure 35.1, and the photographs in **Figure 35.2**

through **Figure 35.6** to identify the distinguishing microscopic characteristics of each parasite studied.

2. Record your observations in the Lab Report.




(a) *Entamoeba histolytica* (650 $\times$ )




(b) A trophozoite of *Entamoeba histolytica* (2000 $\times$ )

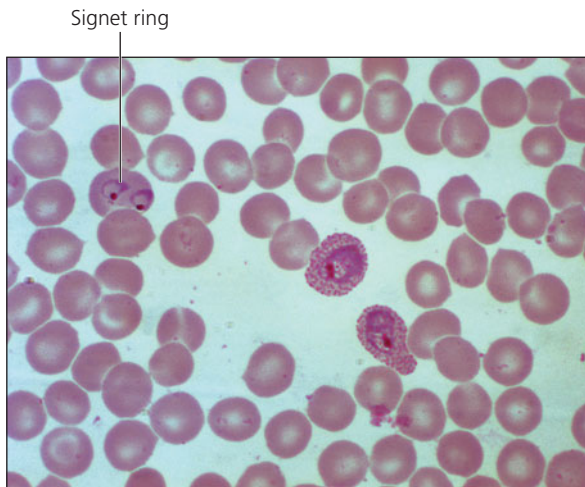**Figure 35.2** *Entamoeba histolytica*. Causative agent of amebic dysentery.




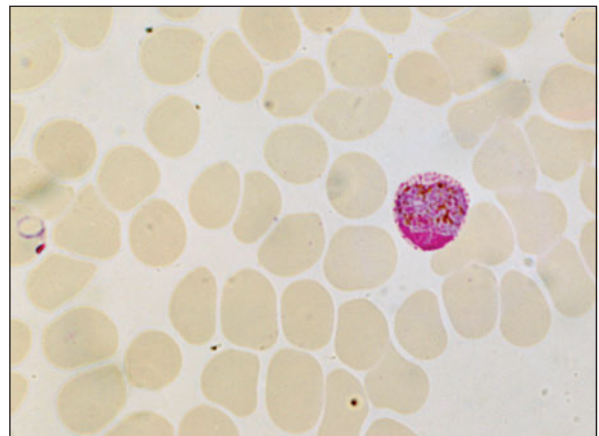(a) Red blood cells infected by the ring stage (signet ring) (850 $\times$ )




(b) Late schizont (segmenter) stage (1100 $\times$ )

**Figure 35.3** *Plasmodium vivax*. Causative agent of malaria.

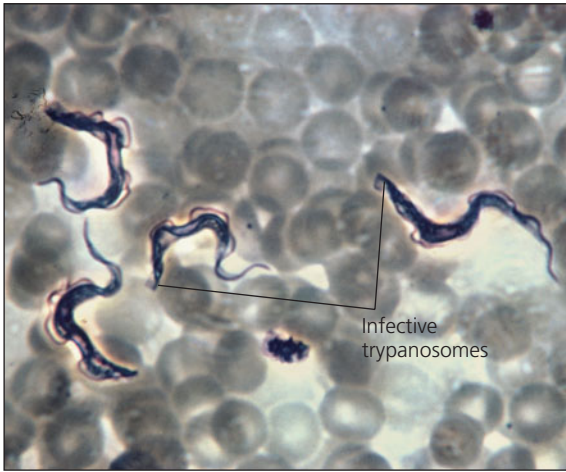

**Figure 35.4** *Trypanosoma gambiense*. Causative agent of African sleeping sickness (1300x).

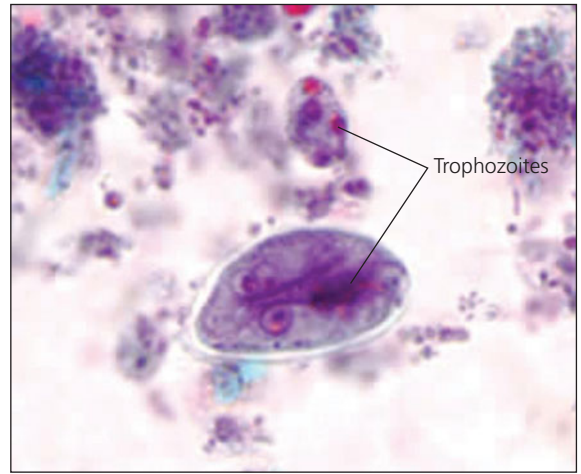

**Figure 35.6** *Giardia intestinalis*. Causative agent of gastrointestinal diarrhea.




**Figure 35.5** A cyst of the ciliated protozoan *Balantidium coli*. The cysts are spherical and lack surface cilia.

Name: \_\_\_\_\_


Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

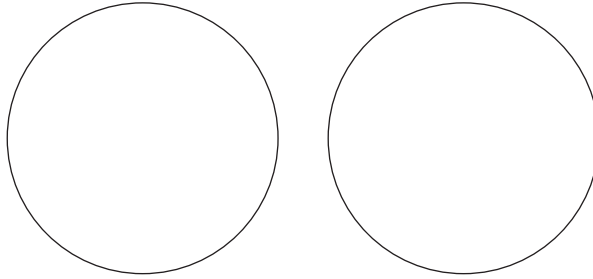Draw representative sketches of the parasitic organisms that you studied, and label the distinguishing structural characteristics you were able to observe.




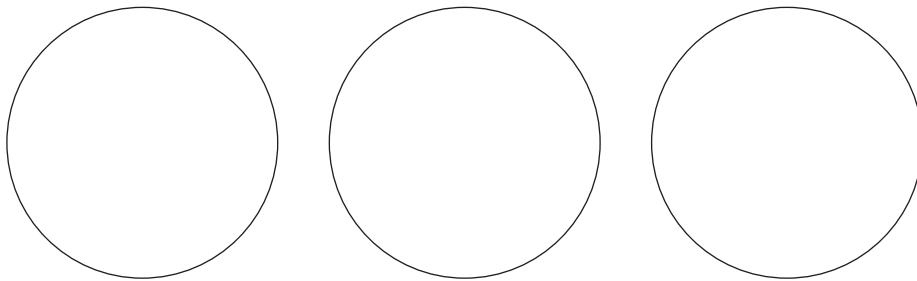*E. histolytica*




*G. intestinalis (G. lamblia)*

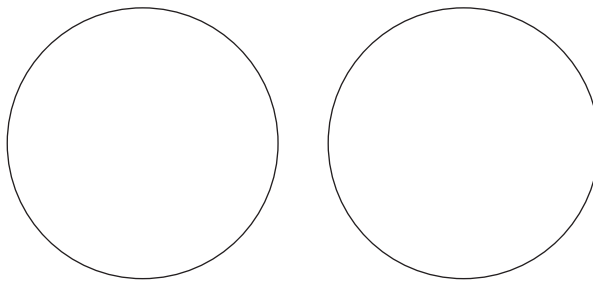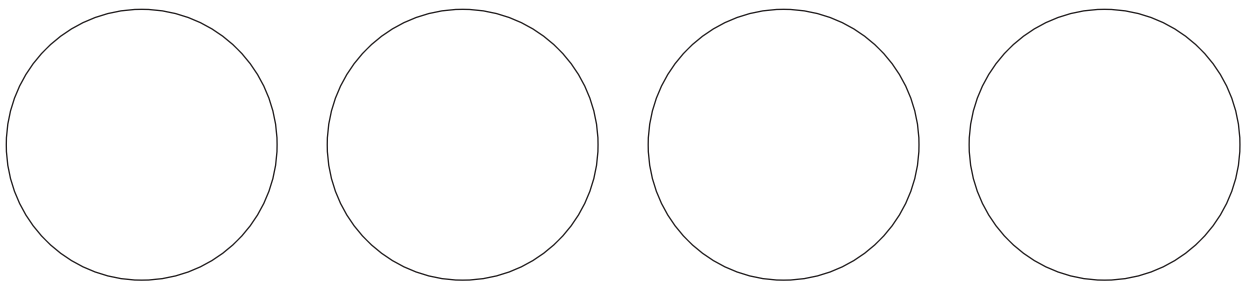*T. gambiense*





*B. coli*



*P. vivax*: Erythrocytic stages



## Review Questions

1. Describe the developmental stages of the malarial parasite during sporogony and schizogony.
2. What role does the invertebrate host play in the life cycle of the trypanosomes? Explain.
3. Distinguish between the pre-erythrocytic and erythrocytic stages in the life cycle of the malarial parasite.
4.  In malarial infections, the sexually mature parasite is found in which host? Is this true for all other protozoan parasitic infections? Explain.
5.  On returning from a trip overseas, an individual with persistent diarrhea is diagnosed as having an *E. histolytica* infection. Fecal examination reveals the presence of blood in the stool, suggesting damage to the intestinal mucosa. Explain why and how the mucosa was compromised by this parasite.

# The Fungi

## LEARNING OBJECTIVES

*When you have completed the experiments in this section, you should be familiar with*

1. The macroscopic and microscopic structures of yeast and molds.
2. The basic mycological culturing and staining procedures.
3. The ability to identify selected common fungal organisms.

## Introduction

The branch of microbiology that deals with the study of fungi (yeasts and molds) is called **mycology**. True fungi are separated into the following four groups on the basis of their sexual modes of reproduction:

1. **Zygomycetes:** Bread and terrestrial molds. Reproductive spores are external and uncovered. Sexual spores are zygospores, and asexual spores are sporangiospores.
2. **Ascomycetes:** Yeasts and molds. Sexual spores, called ascospores, are produced in a saclike structure called an ascus. Conidia are asexual spores produced on a conidiophore.
3. **Basidiomycetes:** Fleshy fungi, toadstools, mushrooms, puffballs, and bracket fungi. Reproductive spores, basidiospores, are separate from specialized stalks called basidia.
4. **Deuteromycetes:** Also called **Fungi Imperfecti** because no sexual reproductive phase has been observed.

The major characteristics of these four groups of fungi are shown in Table P7.1.

Nutritionally, the fungi are heterotrophic, eukaryotic microorganisms that are enzymatically capable of metabolizing a wide variety of organic substrates. Fungi can have beneficial or detrimental effects on humans. Those that inhabit the soil play a vital role in decomposing dead plant and animal tissues, thereby maintaining a fertile soil environment. The fermentative fungi are of industrial importance in producing beer and wine, bakery products, cheeses, industrial enzymes, and antibiotics. The detrimental activities of some fungi include spoilage of foods by rots, mildews, and rusts found on fruit, vegetables, and grains. Some species are capable of producing toxins (for example, aflatoxin) and hallucinogens. A few fungal species are of medical significance because of their capacities to produce diseases in humans. Many of the pathogenic fungi are deuteromycetes and can be divided into two groups based on site of infection. The **superficial mycoses** cause infections of the skin, hair, and nails (for example, ringworm infections). The **systemic mycoses** cause infections of the subcutaneous and deeper tissues such as those of the lungs, genital areas, and nervous system.

**TABLE P7.1 Major Characteristics of the Four Groups of Fungi**

CHARACTERISTICS	GROUP			
	ZYGOMYCETES	ASCOMYCETES	BASIDIOMYCETES	DEUTEROMYCETES
Mycelium	Nonseptate	Septate	Septate	Septate
Asexual spores	Found in sporangium; sporangiospores (nonmotile)	Formed on tip of conidiophore; conidia (nonmotile)	Same as the ascomycetes	Same as the ascomycetes
Sexual spores	Zygospores (motile), found in terrestrial forms; oospores, found in aquatic forms	Ascospores, contained in a saclike structure called the ascus	Basidiospores, carried on the outer surface of a club-shaped cell called the basidium	Fungi Imperfecti—no sexual reproductive phase observed; some members of the ascomycetes and basidiomycetes are Fungi Imperfecti
Common species	Bread molds, mildews, potato blight, <i>Rhizopus</i> species	Cup fungi, ergot, Dutch elm, yeast species	Smuts, rusts, puffballs, toadstools, mushrooms	<i>Aspergillus</i> , <i>Candida</i> , <i>Trichophyton</i> , <i>Cryptococcus</i> , <i>Blastomyces</i> , <i>Histoplasma</i> , <i>Microsporium</i> , and <i>Sporothrix</i>

Molds are the major fungal organisms that can be seen by the naked eye. We have all seen them growing on foods such as bread or citrus fruit as a cottony, fuzzy, black, green, or orange growth, or as a mushroom with a visible cap attached to a stalk, depending on the mold. Examination with a simple hand lens shows that these organisms are composed of an intertwining branching mat called a **mycelium**. The filaments that make up this mycelial mat are called **hyphae**. Most of the mat grows on or in the surface of the nutrient medium so that it can extract nutrients; the mat is therefore called **vegetative mycelium**. Some of the mycelium mat rises upward from the mat and is referred to as **aerial mycelium**. Specialized hyphae are produced from the aerial mycelium and give rise to spores that are the reproductive elements of the mold. **Figures 36.1, 36.2, 36.5, and 36.6** show the reproductive structures of some fungi.

The cultivation, growth, and observation of molds require techniques that differ from those used for bacteria. Mold cultivation requires the use of a selective medium such as Sabouraud agar or potato dextrose agar. These media favor mold growth because their low acidity (pH 4.5 to 5.6) discourages the growth of bacteria, which favor a neutral (pH 7.0) environment. The temperature requirements of molds are also different from those of bacteria, in

that molds grow best at room temperature (25°C). In addition, molds grow at a much slower rate than bacteria do, requiring several days to weeks before visible colonies appear on a solid agar surface. Colony growth is shown in **Figures 36.3 and 36.4**.

## PART A Slide Culture Technique


### LEARNING OBJECTIVES

*When you have completed this experiment you should be*


1. Acquainted with mold cultivation on glass slides.
2. Able to visualize and identify the structural components of molds.

### Principle

Because the structural components of molds are very delicate, even simple handling with an inoculating loop may result in mechanical disruption of their components. The following slide culture

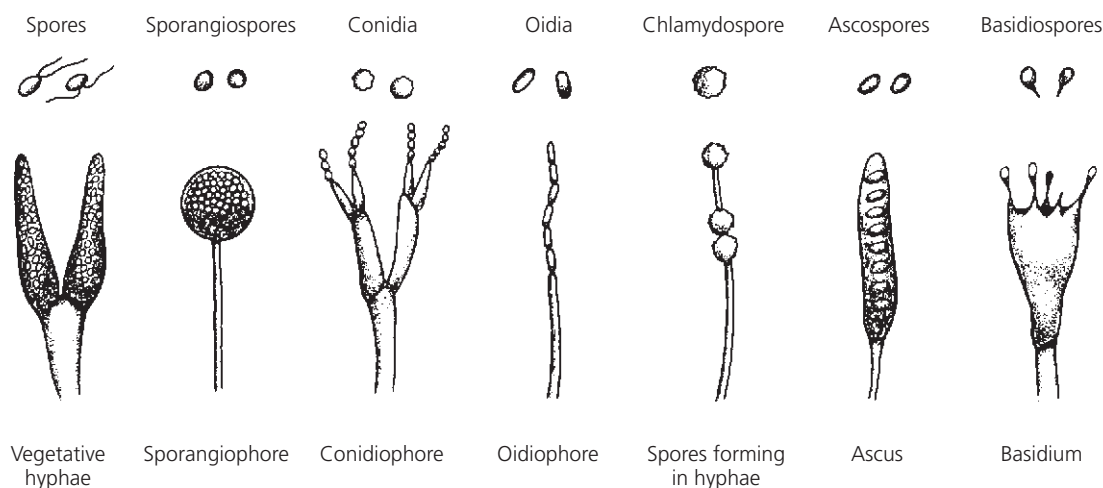


**Figure 36.1** Spore and sporangia types

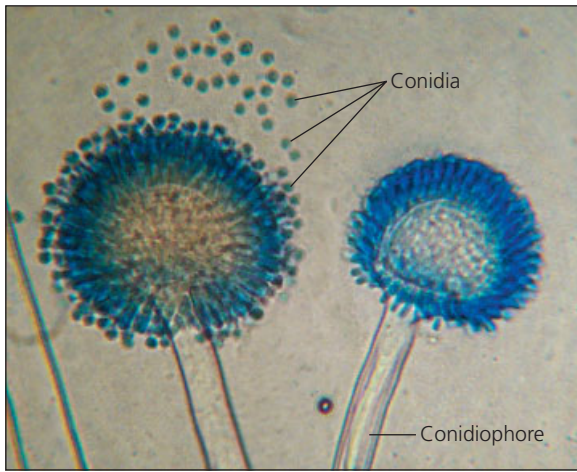

**Figure 36.2** Conidiophore and conidia of mold *Aspergillus niger*




**Figure 36.3** Colony of *Penicillium chrysogenum*

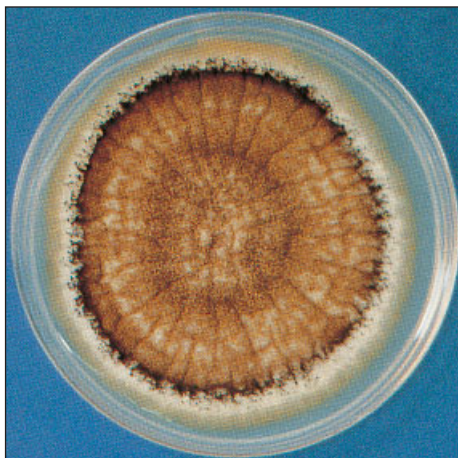

**Figure 36.4** Colony of *Aspergillus niger* on a Sabouraud agar plate



**Figure 36.5** *Mucor mucedo*

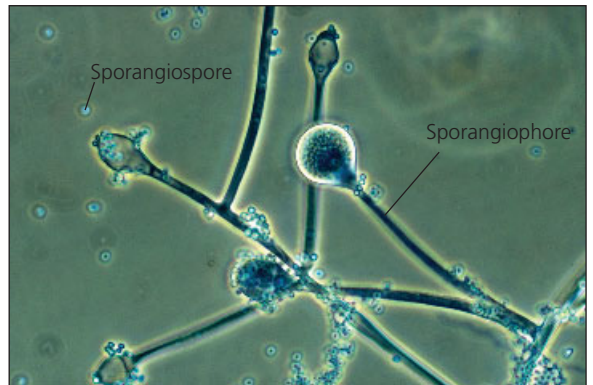

**Figure 36.6** *Rhizopus stolonifer*

technique is used to avoid such disruption. A deep concave slide containing a suitable nutrient medium with an acidic pH, such as Sabouraud agar, is covered by a removable coverslip. Mold spores are deposited in the surface of the agar and incubated in a moist chamber at room temperature. Direct microscopic observation is then possible without fear of disruption or damage to anatomical components. Molds can be identified as to spore type and shape, type of sporangia, and type of mycelium, as shown in Figure 36.1 and Table 38.1 on pages 255–257.

## CLINICAL APPLICATION

### Cultivation of Fungi on Glass Slides

Since sporangia may be damaged during transfer to a glass slide, the slide culture technique prevents the disturbance and damage of the sporangia and other spore structures required for fungi identification. Intact samples can be used to distinguish a fungi like *Aspergillus niger*, which causes the most common fungal infection of the ear, from *Aspergillus flavus*, a fungal pathogen that may result in disseminating infection of the lungs.

## AT THE BENCH



## Materials

### Cultures

7- to 10-day Sabouraud agar cultures of *Penicillium chrysogenum* (formerly called *P. notatum*) and *Aspergillus niger*, *Mucor mucedo*, and *Rhizopus stolonifer*.

### Media

Per student group: one Sabouraud agar deep tube.

### Equipment

Bunsen burner, waterbath, four concave glass slides, four coverslips, petroleum jelly, sterile Pasteur pipettes, toothpicks, four sterile Petri dishes, filter paper, forceps, inoculating loop and needle, four sterile U-shaped bent glass rods, thermometer, dissecting microscope, and beaker with 95% ethyl alcohol.

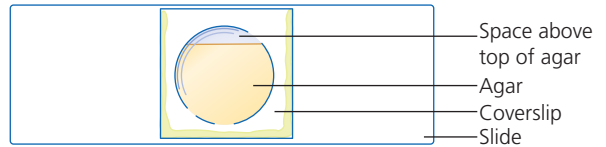
## Procedure Lab One

1. Melt the deep tube of Sabouraud agar in a boiling water bath and cool to 45°C.
2. Place a piece of filter paper in the bottom of each Petri dish, lay a sterile bent glass rod in each dish, and replace the covers.
3. Using forceps, dip the concave slides and coverslips in a beaker of 95% ethyl alcohol, pass through Bunsen burner flame, remove from flame, and hold until all the alcohol has burned off the slides and coverslips.
4. Cool slides and coverslips. Place a slide, concave side up, with a coverslip to one side of the concavity, on the glass rod inside each Petri dish.
5. With a toothpick, add petroleum jelly to three sides surrounding the concavity of each slide. The fourth side will serve as a vent for air.



6. With a sterile Pasteur pipette, add one or two drops of cooled Sabouraud agar to the concavity of each slide.
7. Place a coverslip over the concave portion of each slide so that it is completely sealed.

8. With forceps, stand each slide upright inside its respective Petri dish until the agar solidifies, as illustrated below:



9. When agar is fully hardened, slide coverslips downward with forceps, and with a sterile needle inoculate each prepared slide with the spores from the test cultures.
10. Push the coverslips to their original positions, thereby sealing off the slide.
11. With a Pasteur pipette, moisten the filter paper with sterile water to provide a moist atmosphere. Remoisten filter paper when necessary during the incubation period.
12. Place the slide on the U-shaped bent rod, replace Petri dish cover, and label with the names of the organism and your initials.
13. Incubate the preparations for 7 days at 25°C.

## Procedure Lab Two

1. Examine each mycological slide preparation under the low and high power of a dissecting microscope. Identify the mycelial mat, vegetative and reproductive hyphae, and spores. Use Table 38.1 on pages 255–257 to aid with your identification of mold structures.
2. Record your observations in the Lab Report.

## PART B Mold Cultivation on Solid Surfaces

### LEARNING OBJECTIVES

When you have completed this experiment, you should be

1. Acquainted with the technique of mold cultivation on agar plates.
2. Able to observe and identify colonial characteristics such as growth rate, texture, pigmentation on the surface and reverse side, and folds or ridges on the surface.

## Principle

Cultivating molds on solid surfaces allows you to observe the variations in gross colonial morphology among different genera of molds. These variations in colonial appearance play a major role in the identification of the filamentous fungi. Most microbiologists are familiar with the gross appearance of multicellular fungi, but even to the untrained, the macroscopic differences in colonial growths are obvious and recognizable. For example, most people have seen rotting citrus fruits (lemons and oranges) produce a blue-green velvety growth characteristic of *Penicillium* species. It is also common for stale cheese to show a grayish-white furry growth of *Mucor* species, and the black stalklike appearance of *Rhizopus* molds growing on bread is familiar to many.

In this part of the experiment, you will be able to visualize the gross appearance of the colonial growth of four different molds.

### CLINICAL APPLICATION

#### Isolation of Fungi on Solid Media

Before a fungal species may be identified or studied it must first be isolated. Similar to using an agar plate for isolating a distinct bacterial species, agar plating may be used as a growth medium for the isolation of fungi spores. Once spores have been isolated from individual sporangia, subculturing on solid agar or slides will allow for characterization and genetic studies of the fungus.

### AT THE BENCH



## Materials

### Cultures

7- to 10-day Sabouraud agar cultures of *Aspergillus niger*, *Penicillium chrysogenum* (formerly called *P. notatum*), *Mucor mucedo*, and *Rhizopus stolonifer*.

## Media

Per designated student group: three Sabouraud agar plates and one potato dextrose agar plate.

## Equipment

Bunsen burner, four test tubes containing 2 ml of sterile saline, dissecting microscope, and an inoculating loop.

## Procedure Lab One

1. Label the three Sabouraud agar plates as *Aspergillus niger*, *Penicillium chrysogenum*, and *Mucor mucedo*, and label the fourth plate containing potato dextrose agar as *Rhizopus stolonifer*.
2. Prepare a saline suspension of each mold culture. Label each of the four tubes of saline with the name of the organism. Using a sterile inoculating loop, scrape two loopfuls of mold culture into the corresponding tube of 2 ml of sterile saline and mix well by tapping the tube with your finger.
3. Using aseptic technique, inoculate each of the plates by placing a single loopful of mold suspension in the center of its respective agar plate. *Note: Do not spread the inoculum and do not shake or jostle the plates.*
4. Incubate all plates at room temperature, 25°C, for 7 to 10 days. *Note: Do not invert the plates.*

## Procedure Lab Two

1. Examine each mold plate under the low and high power of a dissecting microscope. Refer to Table 38.1 for your identification of mold structures. *Note: Do not remove Petri dish covers.*
2. Record your observations in the Lab Report.

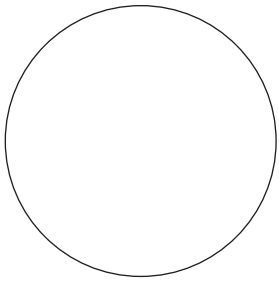
Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_


## Observations and Results

### PART A: Slide Culture Technique


Draw a representative microscopic field under low-power and high-power magnification and label the structural components of each test organism.




Low power

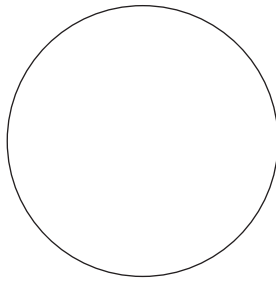*Penicillium chrysogenum (P. notatum)*

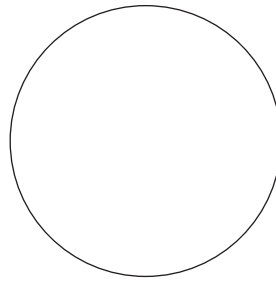
High power




Low power

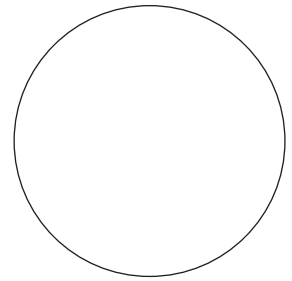*Aspergillus niger*

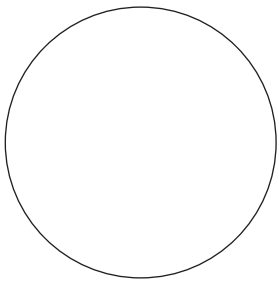
High power




Low power

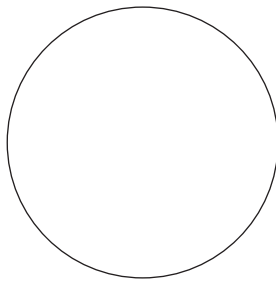*Rhizopus stolonifer*

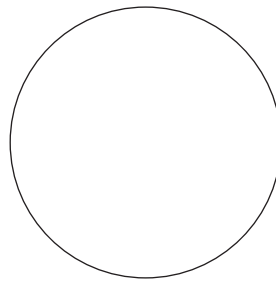High power



Low power

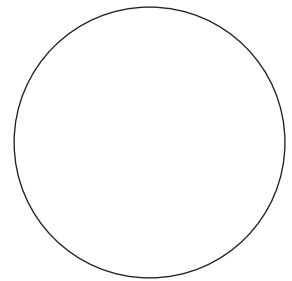
*Mucor mucedo*

High power



## PART B: Mold Cultivation on Solid Surfaces

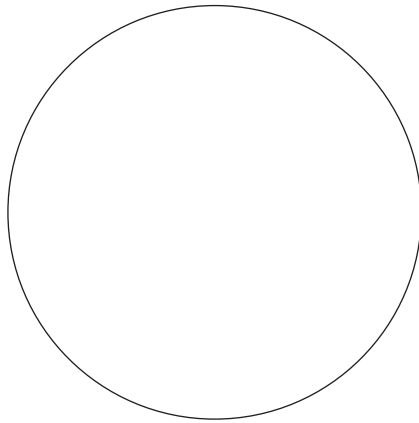Draw sketches of the mold colonies under low power, indicating the extent of growth (diameter in mm), pigmentation, and the presence or absence of aerial hyphae. Refer to Table 38.1 to aid with your identification of mold structures.



*Penicillium chrysogenum (P. notatum)*

Colony diameter (mm):

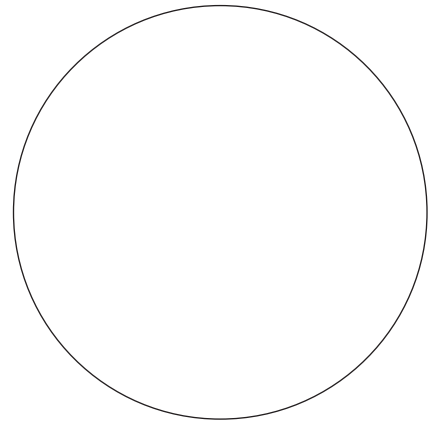
\_\_\_\_\_

Pigmentation:

\_\_\_\_\_

Aerial hyphae (+ or -):

\_\_\_\_\_




*Aspergillus niger*

\_\_\_\_\_

\_\_\_\_\_

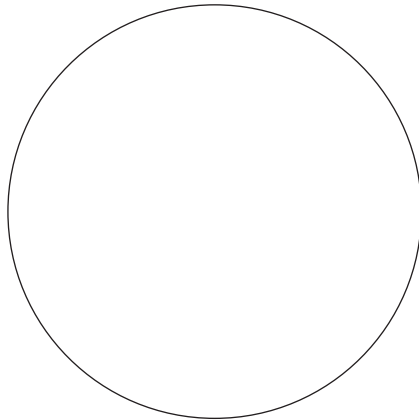\_\_\_\_\_



*Rhizopus stolonifer*

Colony diameter (mm):

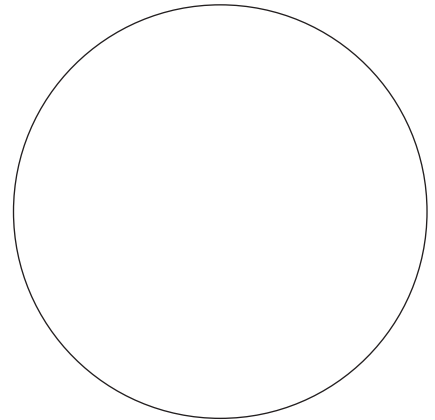
\_\_\_\_\_

Pigmentation:

\_\_\_\_\_

Aerial hyphae (+ or -):

\_\_\_\_\_



*Mucor mucedo*



\_\_\_\_\_

\_\_\_\_\_

\_\_\_\_\_

## Review Questions

1. Cite some beneficial and harmful aspects of molds.
2. What is the advantage of using Sabouraud agar?
3. In the slide culture technique, what is the purpose of the following?
  - a. Moistened filter paper in the Petri dish:
  - b. A U-shaped glass rod in the Petri dish:

4.  What is the advantage of the slide culture technique over that of a simple loop inoculation onto an agar plate (as in Part B)?
5. Why would it be advantageous to observe mold colonies on an agar plate?
6.  Since dimorphism is a property of fungi, how do you account for the fact that molds grow preferentially *in vitro* rather than *in vivo*?

# Yeast Morphology, Cultural Characteristics, and Reproduction

## EXPERIMENT 37

### LEARNING OBJECTIVES

Once you have completed this experiment, you should

1. Know the morphology of different genera of yeast.
2. Understand the growth and fermentative properties of yeast cells.
3. Be familiar with the sexual and asexual modes of reproduction in yeast cells.

### Principle

Yeasts are nonfilamentous unicellular fungi. Yeast cultures resemble bacteria when grown on the surface of artificial laboratory media; however, they are 5 to 10 times larger than bacteria. Yeast colonies are illustrated in **Figure 37.1**. Microscopically, yeast cells may be ellipsoidal, spherical, or in some cases, cylindrical (**Figure 37.2**). Unlike molds, yeast do not have aerial hyphae and supporting sporangia.

Yeast reproduce asexually by **budding** or by **fission**. In budding, an outgrowth from the parent cell (a **bud**) pinches off, producing a daughter

cell (**Figures 37.3a** and **37.4**). Fission occurs in certain species of yeast, such as those in the genus *Schizosaccharomyces*. During fission, the parent cell elongates, its nucleus divides, and it splits evenly into two daughter cells.

Some yeast may also undergo sexual reproduction when two sexual spores conjugate, giving rise to a zygote, or diploid cell. The nucleus of this cell divides by meiosis, producing four new haploid nuclei (sexual spores), called **ascospores**, contained within a structure called the **ascus** (**Figure 37.3b**). When the ascus ruptures, the ascospores are released and conjugate, starting the cycle again.

Yeasts are important for many reasons. *Saccharomyces cerevisiae* is referred to as baker's yeast and is used as the leavening agent in dough. Two major strains of yeast, *Saccharomyces carlsbergensis* and *Saccharomyces cerevisiae*, are used for brewing. The wine industry relies on wild yeast (present on the grape) for the fermentation of grape juice, which is supplemented with *Saccharomyces ellipsoideus* to begin the fermentation. Also, the high vitamin content of yeasts makes them particularly valuable as food supplements. As useful as some yeasts are, there are a few species that can create problems in the food industry or are harmful to humans. Undesired



(a) *Saccharomyces cerevisiae*




(b) *Candida albicans*




(c) *Rhodotorula rubra*


**Figure 37.1** Colonies of yeast cells

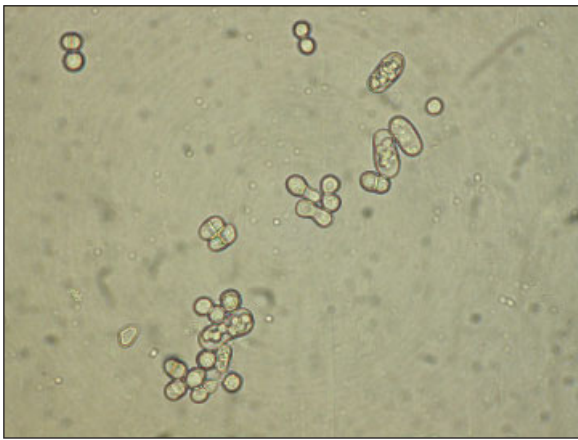


**Figure 37.2** Yeast Cells. *Schizosaccharomyces octosporus*.

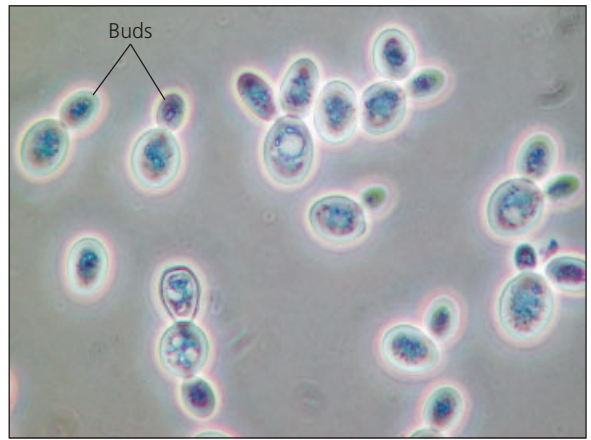

**Figure 37.4** Asexual yeast reproduction by budding. *Saccharomyces cerevisiae*.

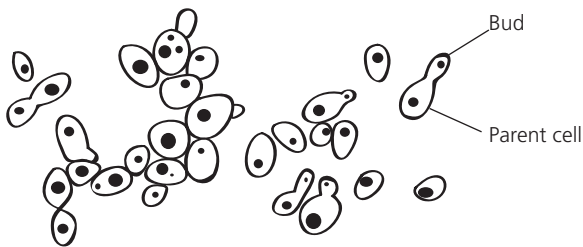

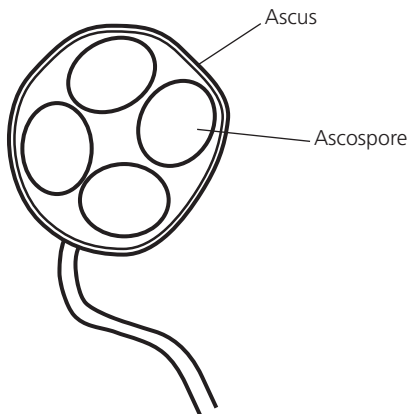(a) Asexual reproductive yeast structures



(b) Sexual reproductive yeast structures

**Figure 37.3** Reproductive structures of yeast

yeast must be excluded from the manufacture of fruit juices, such as grape juice or apple cider, to prevent the fermentation of fruit sugars to alcohol. The contamination of soft cheese by some forms of yeast will destroy the product. Finally, some yeast such as *Candida albicans* are pathogenic and responsible for urinary tract and vaginal infections known as **moniliasis** and infections of the mouth called **thrush**.

The cultural characteristics, the type of reproduction, and the fermentative activities used to identify the different genera of yeast will all be studied in this experiment.

## CLINICAL APPLICATION

### Opportunistic Yeast

Opportunistic yeast may cause serious or life-threatening infections in immunocompromised patients, such as people with cancer or AIDS. The yeast *Candida albicans* can cause relatively minor infections in healthy people, such as thrush or vaginal “yeast infections,” however it can cause a dangerous bloodstream infection called *invasive candidiasis* in those with weakened immune systems. Similarly, the yeast *Cryptococcus neoformans* can cause a pulmonary infection that can lead to meningitis, most often in immunocompromised people.

## AT THE BENCH



## Materials

### Cultures

7-day Sabouraud agar cultures of *Saccharomyces cerevisiae*, *Candida albicans*, *Rhodotorula rubra*, *Selenotila intestinalis*, and *Schizosaccharomyces octosporus*.

## Media

Per designated student group: five tubes each of bromcresol purple glucose broth, bromcresol purple maltose broth, bromcresol purple lactose broth, and bromcresol purple sucrose broth, each containing a Durham tube; two glucose-acetate agar plates; and five test tubes (13 × 100 mm) containing 2 ml of sterile saline.

## Reagents

Water-iodine solution, lactophenol–cotton-blue solution.

## Equipment

Bunsen burner, inoculating loop and needle, 10 glass slides, 10 coverslips, 5 sterile Pasteur pipettes, glassware marking pencil, and microscope.

## Procedure Lab One

### Morphological Characteristics

Prepare a wet mount of each yeast culture in the following manner:

1. Suspend a loopful of yeast culture in a few drops of lactophenol–cotton-blue solution on a microscope slide and cover with a coverslip.
2. Examine all yeast wet-mount slide preparations under low and high power, noting the shape and the presence or absence of budding. Record your observations in the Lab Report.

### Fermentation Studies

1. With a sterile loop, inoculate each experimental organism into appropriately labeled tubes of bromcresol purple glucose, maltose, lactose, and sucrose fermentation broths.
2. Incubate all cultures at 25°C for 4 to 5 days.

### Sexual Reproduction

1. With a glassware marking pencil, divide the bottom of a glucose-acetate agar plate into three sections, and divide another glucose-acetate agar plate in half.

2. Label each section with the name of a test organism.
3. Label each tube of sterile saline with the name of a test organism.
4. With a sterile inoculating loop, suspend a heavy loopful of each test organism into its appropriately labeled tube of saline. Tap the tube with your finger to obtain a uniform cell suspension.
5. With a sterile Pasteur pipette, inoculate one drop of each test organism onto the surface of the appropriately labeled section on an agar plate. *Note: Allow the inoculum to diffuse into the agar for a few minutes. Do not swirl or rotate the plates.*
6. Incubate all plates at 25°C for 7 days. *Note: Visit the laboratory, if possible, during the incubation period and note when sporulation begins.*

## Procedure Lab Two

### Fermentation Studies

1. Examine all fermentation tubes for the presence of growth (turbidity), the presence or absence of acid (change in the color of medium), and the presence or absence of gas (bubble in Durham tube).
2. Record your results in the chart provided in the Lab Report.

## Procedure Lab Three

### Sexual Reproduction

1. Examine the glucose-acetate agar plates for the presence or absence of sporulation.
2. Prepare a water-iodine wet mount using a loopful of culture from each respective section on the glucose-acetate agar plate.
3. Observe the cells using the high-dry objective and record your observations in the Lab Report.

*This page intentionally left blank*

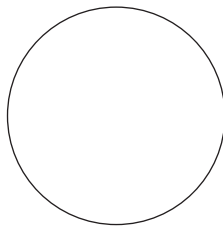
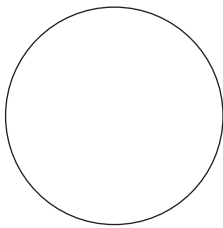
Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_



## Observations and Results

### Morphological Characteristics

Draw a representative field for each organism in the chart below. Note the shape and presence or absence of budding (+ or -).



*Saccharomyces cerevisiae*

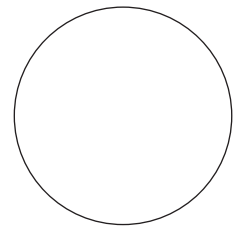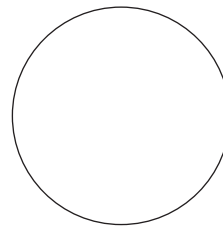



*Candida albicans*

Shape: \_\_\_\_\_



Budding (+ or -): \_\_\_\_\_

Shape: \_\_\_\_\_

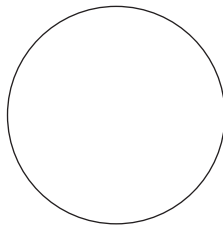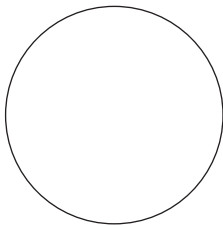Budding (+ or -): \_\_\_\_\_



*Rhodotorula rubra*

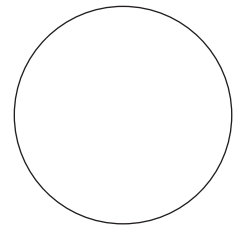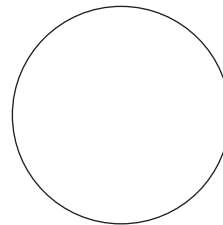



*Selenotila intestinalis*

Shape: \_\_\_\_\_

Budding (+ or -): \_\_\_\_\_

Shape: \_\_\_\_\_

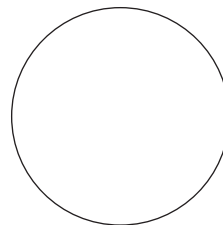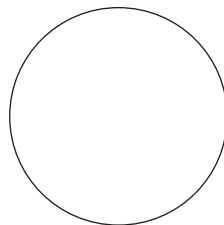Budding (+ or -): \_\_\_\_\_



*Schizosaccharomyces octosporus*

Shape: \_\_\_\_\_

Budding (+ or -): \_\_\_\_\_



## Fermentation Studies

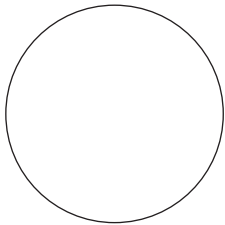
Use a plus (+) or minus (-) in the chart below to record your results.

Organism	GLUCOSE			MALTOSE			LACTOSE			SUCROSE		
	T	A	G	T	A	G	T	A	G	T	A	G
<i>Saccharomyces cerevisiae</i>												
<i>Candida albicans</i>												
<i>Rhodotorula rubra</i>												
<i>Selenotila intestinalis</i>												
<i>Schizosaccharomyces octosporus</i>												


Note: T = turbidity, A = acid, and G = gas

## Sexual Reproduction


In the circles below, draw representative reproductive structures and label the parts.




*Saccharomyces cerevisiae*

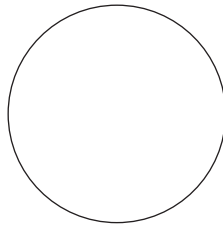


*Candida albicans*

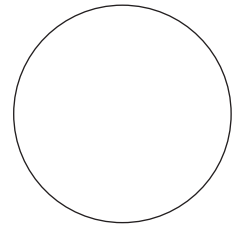

*Schizosaccharomyces octosporus*

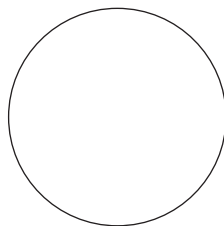

*Rhodotorula rubra*

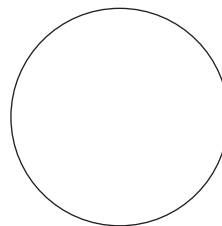


*Selenotila intestinalis*


## Review Questions

1. Indicate the significance of the following structures in the reproductive activities of yeast cells.
  - a. Buds
  - b. Ascus
  - c. Ascospores
2. Why are yeast cells classified as fungi, and how do they differ from other fungi?
3. Why is yeast of industrial importance?

4. Why are yeasts significant from a medical perspective?

5. Why is it necessary to pasteurize fruit juices?

6.  A female patient develops candidiasis (moniliasis) following prolonged antibiotic therapy for a bladder infection caused by *Pseudomonas aeruginosa*. How can you account for the development of this concurrent vaginal infection?

7.  With regard to the fermentation of wine, what kind of wine would be produced if you washed the grapes prior to crushing them?

## LEARNING OBJECTIVE

When you have completed this experiment, you should be able to

1. Identify a fungal unknown based on colonial morphology and microscopic appearance.

## CLINICAL APPLICATION

### Identification of Fungal Infection

When presented with a patient with symptoms that suggest either an intestinal fungal infection or Crohn's Disease complications, isolation and identification of a fungal pathogen is required before the correct treatment may be prescribed. Using culturing techniques, an isolate of the fungi can be identified by mycelium morphology and genetic markers.

## Principle

In this experiment, you will be provided with a number-coded pure culture of a representative fungal organism for cultivation and subsequent identification. Use Table 38.1 to aid in identification of the unknown culture.

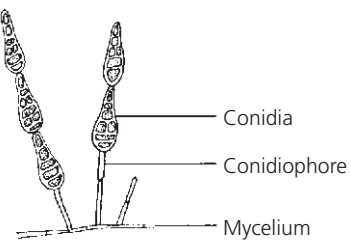
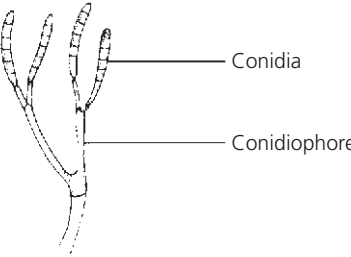
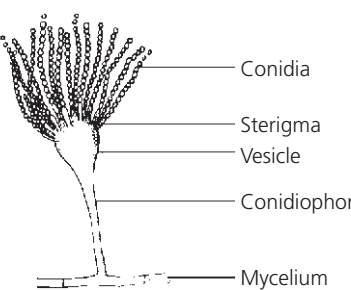
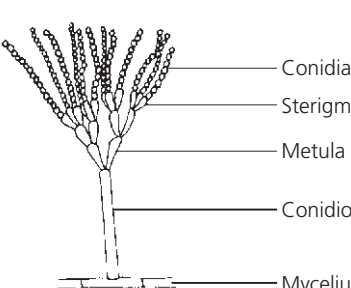
(Text continues on page 258.)

**TABLE 38.1** Identification of Fungi

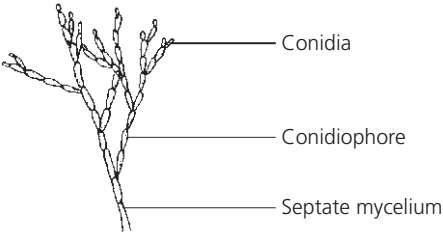
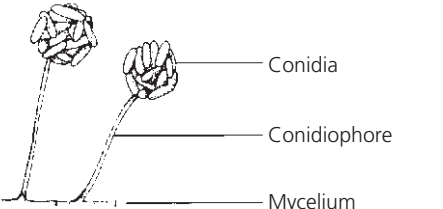
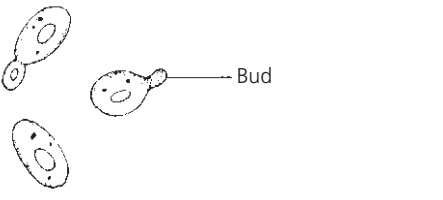
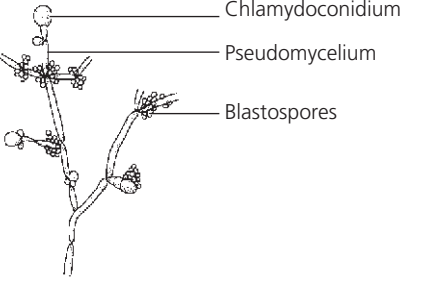
DIAGRAM	COLONIAL MORPHOLOGY	MICROSCOPIC APPEARANCE
<p><b>Molds</b></p> <p>Labels: Sporangium, Columella, Collarete, Sporangiphore, Stolon, Mycelium, Rhizoid</p> <p><i>Rhizopus</i>: Black bread mold; common laboratory contaminant</p>	<p>Rapidly growing white-colored fungus swarms over entire plate; aerial mycelium cottony and fuzzy</p>	<p>Spores are oval, colorless, or brown; nonseptate mycelium gives rise to straight sporangiophores that terminate with black sporangium containing a columella; rootlike hyphae (rhizoids) penetrate the medium</p>
<p>Labels: Sporangium, Columella, Sporangiphore, Mycelium</p> <p><i>Mucor</i>: Food contaminant</p>	<p>Resembles the colonies of <i>Rhizopus</i> except that it lacks rhizoids and collarettes. Sporangiophore arises directly from mycelial mat. Note: Branching sporangiophores may occur with <i>Mucor</i>.</p>	<p>Spores are oval; nonseptate mycelium gives rise to single sporangiophores with globular sporangium containing a columella; there are no rhizoids</p>



**TABLE 38.1 Identification of Fungi (continued)**

DIAGRAM	COLONIAL MORPHOLOGY	MICROSCOPIC APPEARANCE
<p><b>Molds (continued)</b></p>  <p>Conidia Conidiophore Mycelium</p> <p><i>Alternaria:</i> Normally found on plant material; also found in house dust</p>	<p>Grayish-green or black colonies with gray edges rapidly swarming over entire plate; aerial mycelium not very dense, appears grayish to white</p>	<p>Multicelled spores (conidia) are pear-shaped and attached to single conidiophores arising from a septate mycelium</p>
 <p>Conidia Conidiophore</p> <p><i>Fusarium:</i> Found in soil; also likely in eye infections</p>	<p>Woolly, white, fuzzy colonies changing color to pink, purple, or yellow</p>	<p>Spores (conidia) are oval or crescent-shaped and attached to conidiophores arising from a septate mycelium; some spores are single cells, some are multicelled</p>
 <p>Conidia Sterigma Vesicle Conidiophore Mycelium</p> <p><i>Aspergillus:</i> Plant and animal pathogens; some species used industrially</p>	<p>White colonies become greenish-blue, black, or brown as culture matures</p>	<p>Single-celled spores (conidia) in chains developing at the end of the sterigma arising from the terminal bulb of the conidiophore, the vesicle; long conidiophores arise from a septate mycelium</p>
 <p>Conidia Sterigma Metula Conidiophore Mycelium</p> <p><i>Penicillium:</i> Antibiotic-producing citrus fruit contaminant; soil inhabitant</p>	<p>Mature cultures usually greenish or blue-green</p>	<p>Single-celled spores (conidia) in chains develop at the end of the sterigma arising from the metula of the conidiophore; branching conidiophores arise from a septate mycelium</p>

**TABLE 38.1 (continued)**

DIAGRAM	COLONIAL MORPHOLOGY	MICROSCOPIC APPEARANCE
<p><b>Molds (continued)</b></p>  <p>Conidia</p> <p>Conidiophore</p> <p>Septate mycelium</p> <p><i>Cladosporium:</i> Dead and decaying plants</p>	<p>Small, heaped colonies are greenish-black and powdery</p>	<p>Spores (conidia) develop at the end of complex conidiophores arising from a septate mycelium that is usually brownish</p>
 <p>Conidia</p> <p>Conidiophore</p> <p>Mycelium</p> <p><i>Cephalosporium:</i> Antibiotic production</p>	<p>Rapidly growing compact and moist colonies becoming cottony with aerial hyphae that are gray or rose-colored</p>	<p>Single-celled conical or elliptical spores (conidia) held together in clusters at the tips of the conidiophores by a mucoid substance; erect, unbranched conidiophores arise from a septate mycelium</p>
<p><b>Yeast</b></p>  <p>Bud</p> <p><i>Torula:</i> Cheese and food contaminant</p>	<p>Colonies are pink, moist, with unbroken, even edges</p>	<p>Cells are oval, colorless, and reproduce by budding</p>
 <p>Chlamydoconidium</p> <p>Pseudomycelium</p> <p>Blastospores</p> <p><i>Candida:</i> Human pathogen</p>	<p>Colonies are small, round, moist, and colorless, with unbroken, even edges</p>	<p>Yeastlike fungus produces pseudomycelium</p>



## Materials

### Cultures

Number-coded, 7-day Sabouraud broth spore suspensions of *Aspergillus*, *Mucor*, *Penicillium*, *Alternaria*, *Rhizopus*, *Cladosporium*, *Cephalosporium*, *Fusarium*, *Torula*, and *Candida*.

### Media

One Sabouraud agar plate per student.

### Reagent

Lactophenol–cotton-blue solution.

### Equipment

Bunsen burner, dissecting microscope, hand lens, sterile cotton swabs, glass slides, coverslips, inoculating loop, and glassware marking pencil.

## Procedure Lab One

1. With a sterile inoculating loop, inoculate an appropriately labeled Sabouraud agar plate with one of the provided unknown cultures by placing one loopful in the center of the plate.  
*Note: Do not spread culture.*
2. Incubate the plates in a noninverted position for 1 week at 25°C in a moist incubator.

## Procedure Lab Two

1. Observe mold cultures with a hand lens, noting and recording their colonial morphologies.
2. Prepare a wet mount by suspending some of the culture in a few drops of lactophenol–cotton-blue solution. Be gentle to avoid damaging the fungal structures.
3. Examine the preparation under high-power and low-power magnifications with the aid of a dissecting microscope and record your observations in the Lab Report.


Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_


## Observations and Results

1. Draw representative microscopic fields of your culture below.
2. Using Table 38.1 and Figure 4.1 on page 30, identify your unknown fungal organism.
  - a. Color pigmentation
  
  
  
  
  
  
  
  
  
  
  - b. Diameter (mm)
  
  
  
  
  
  
  
  
  
  
  - c. Texture (cottony, smooth, etc.)
  
  
  
  
  
  
  
  
  
  
  - d. Margin (entire, undulating, lobular, etc.)
  
  
  
  
  
  
  
  
  
  
  - e. Aerial hyphae (septate, nonseptate)

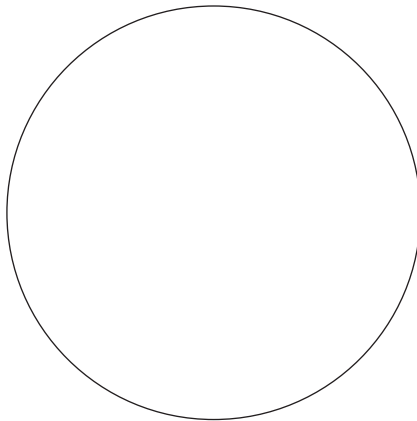Diagram of microscopic appearance:



Low-power magnification

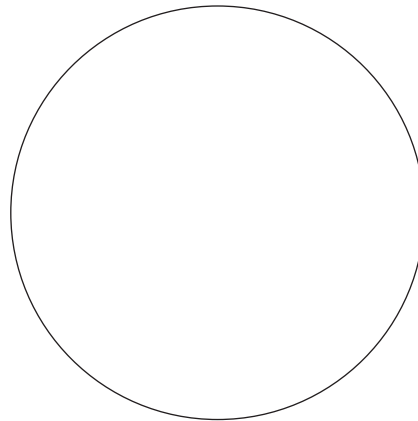

High-power magnification

Number assigned to  
unknown culture: \_\_\_\_\_

Genus of fungal unknown: \_\_\_\_\_



*This page intentionally left blank*

# The Viruses

## LEARNING OBJECTIVES

Once you have completed the experiments in this section, you should know

1. The chemical structures, morphologies, and replicative activities of bacterial viruses (bacteriophages).
2. How to perform a phage dilution procedure for the cultivation and enumeration of bacterial viruses.
3. How to isolate bacteriophages from sewage.

## Introduction

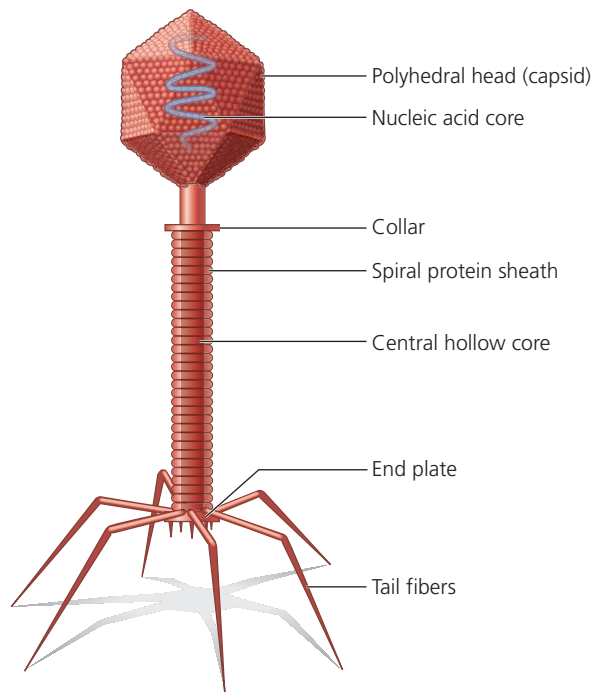
Viruses are noncellular biological entities composed solely of a single type of nucleic acid surrounded by a protein coat called the capsid. Because of their limited and simplistic structures, viruses can be chemically defined as **nucleoproteins**. They are devoid of the sophisticated enzymatic and biosynthetic machinery essential for independent activities of cellular life. This lack of metabolic machinery mandates that they exist as parasites, and they cannot be cultivated outside of a susceptible living cell. Viruses are differentiated from cellular forms of life on the following bases:

1. They are ultramicroscopic and can only be visualized with the electron microscope.
2. They are filterable: They are able to pass through bacteria-retaining filters.
3. They do not increase in size.
4. They must replicate within a susceptible cell.
5. Replication occurs because the viral nucleic acid subverts the synthetic machinery of the host cell (namely, common host cell components and enzyme systems involved in decomposition, synthesis, and bioenergetics) for the purpose of producing new viral components.
6. Viruses are designated either RNA or DNA viruses because they contain one of the nucleic acids but never both.

Much of our knowledge of the mechanism of animal viral infection and replication has been based on our understanding of infection in bacteria by bacterial viruses, called the **bacteriophages**, or **phages**. The bacteriophages were first described in 1915 almost simultaneously by Twort and d'Herelle. The name *bacteriophage*, which in Greek means “to eat bacteria,” was coined by d'Herelle because of the destruction through lysis of the infected cell. Bacteriophages exhibit notable variability in their sizes, shapes, and complexities of structure. The T-even (T2, T4, and T6) phages illustrated in **Figure P8.1** demonstrate the greatest morphological complexity.

Phage replication depends on the ability of the phage particle to infect a suitable bacterial host cell. Infection consists of the following sequential events:

1. **Adsorption:** Tail fibers of the phage particle bind to receptor sites on the host's cell wall.
2. **Penetration (infection):** Spiral protein sheath retracts, and an enzyme, early muramidase, perforates the bacterial cell wall, enabling the phage nucleic acid to pass through the hollow core into the host cell's



The functions of these structural components are as follows:

Component	Function
<b>Capsid</b> (protein coat)	Protection of nucleic acid from destruction by DNases
<b>Nucleic acid core</b>	Phage genome carrying genetic information necessary for replication of new phage particles
<b>Spiral protein sheath</b>	Retracts so that nucleic acid can pass from capsid into host cell's cytoplasm
<b>End plate</b> and <b>tail fibers</b>	Attachment of phage to specific receptor sites on a susceptible host's cell wall


**Figure P8.1 Bacteriophage: Structural components and their functions**

cytoplasm. The empty protein shell remains attached to the cell wall and is called the protein ghost.

3. **Replication:** The phage genome subverts the cell's synthetic machinery, which is then used for the production of new phage components.
4. **Maturation:** During this period, the new phage components are assembled and form complete, mature virulent phage particles.
5. **Release:** Late muramidase (lysozyme) lyses the cell wall, liberating infectious phage

particles that are now capable of infecting new susceptible host cells, thereby starting the cycle over again.

**Virulent phage** particles that infect susceptible host cells always initiate the **lytic cycle** as described above. Other phage particles, called **temperate phages** [that is, **lambda ( $\lambda$ ) phages**], incorporate their nucleic acid into the host's chromosome. Lysis of the host cell does not occur until it is induced by exogenous physical agents such as ultraviolet or ionizing radiation or

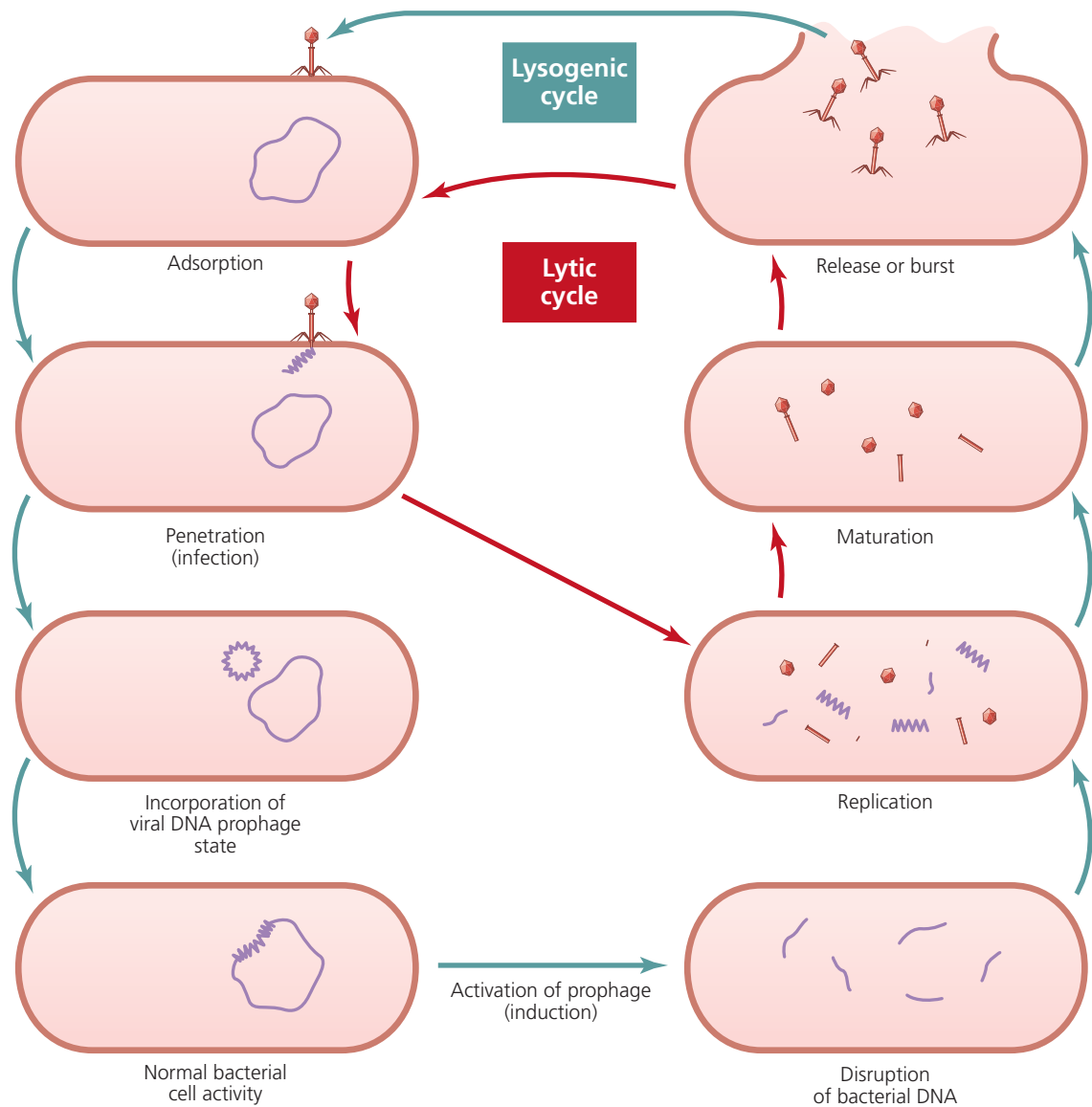

**Figure P8.2** The lytic and lysogenic life cycles of a bacteriophage

chemical mutagenic agents. Bacterial cells containing the incorporated phage nucleic acid, the **prophage**, are called **lysogenic cells**. Lysogenic cells appear and function as normal cells, and they reproduce by fission. When induced by physical or chemical agents, these cells will release a virulent prophage from the host's genome, which then initiates the lytic cycle. **Figure P8.2** illustrates the lytic and lysogenic life cycles of a bacteriophage.

Animal viruses differ structurally from bacteriophages in that they lack the spiral protein sheath, end plate, and tail fibers. Their shapes may be helical or cuboidal (icosahedral, containing 20 triangular facets). Some animal viruses are designated as naked viruses because they are composed solely of nucleocapsids. In others, referred to as enveloped viruses, the nucleocapsid is surrounded by a lipid bilayer that may have glycoproteins associated with it.

The infectious process of the animal virus is very similar to bacteriophage infection. However, there are some notable differences:

1. Adsorption of the virus is to receptor sites that are located on the cell membrane of the host cell instead of the cell wall as in the bacterial host.
2. Viral penetration is accomplished by endocytosis, an energy-requiring, receptor-mediated process in which the entire virus enters the host cell.
3. The uncoating of the animal virus, removal of the capsid, occurs within the host cell; with bacteriophage infection, the phage capsid remains on the outside of the host.
4. The latent period, the time between adsorption and the release of virulent viral particles, is considerably longer—hours to days rather than minutes as in bacteriophage infection.

## LEARNING OBJECTIVE

Once you have completed this experiment, you should be able to


1. Perform techniques for cultivation and enumeration of bacteriophages.

## Principle

This exercise demonstrates the ability of viruses to replicate inside a susceptible host cell. For this purpose, you will be provided with a virulent phage and a susceptible host cell culture. This technique also enables you to enumerate phage particles on the basis of plaque formation in a solid agar medium. **Plaques** are clear areas in an agar medium previously seeded with a diluted phage sample and a host cell culture. Each plaque represents the lysis of a phage-infected bacterial cell.

The procedure requires the use of a double-layered culture technique in which the hard agar serves as a base layer, and a mixture of phage and host cells in a soft agar forms the upper overlay. Susceptible *Escherichia coli* cells multiply rapidly and produce a lawn of confluent growth on the medium. When one phage particle adsorbs to a susceptible cell, penetrates the cell, replicates, and goes on to lyse other host cells, the destroyed cells produce a single plaque in the bacterial lawn (see **Figure 39.1**). Each plaque can be designated as a **plaque-forming unit (PFU)** and used to quantify the number of infective phage particles in the culture.

The number of phage particles contained in the original stock phage culture is determined by

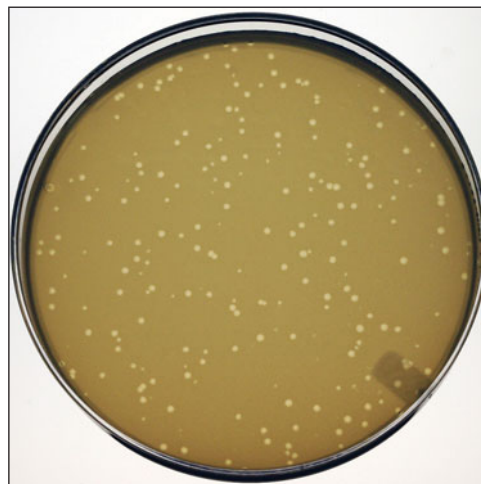

**Figure 39.1** Plaque-forming units (PFUs)

counting the number of plaques formed on the seeded agar plate and multiplying this by the dilution factor. For a valid phage count, the number of plaques per plate should not exceed 300 nor be less than 30.

**Example:** 200 PFUs are counted in a  $10^{-6}$  dilution.

$$(200) \times (10^6) = 200 \times 10^6 \text{ or } 2 \times 10^8 \\ \text{PFUs per ml of stock phage culture}$$

Plates showing greater than 300 PFUs are **too numerous to count (TNTC)**; plates showing fewer than 30 PFUs are **too few to count (TFTC)**.

## CLINICAL APPLICATION

### Identification of Pathogenic Bacteria

Bacterial viruses (bacteriophages) are very common in all natural environments and are directly related to the number of bacteria present. They are most prevalent in soil, intestines of animals, sewage, and seawater. These viral particles have played an important role in the development of all types of viruses. Since many phages are specific about which bacteria they attack, a process called phage typing is used in clinical and diagnostic laboratories for the identification of pathogenic bacteria.

## AT THE BENCH



## Materials

### Cultures

24-hour nutrient broth cultures of *Escherichia coli* B and T2 coliphage.

### Media

Five each of the following per designated student group: tryptone agar plates and tryptone soft agar, 2 ml per tube; and nine tryptone broth tubes, 9 ml per tube.

### Equipment

Bunsen burner, waterbath, thermometer, 1-ml sterile pipettes, sterile Pasteur pipettes, mechanical pipetting device, test tube rack, and glassware marking pencil.

## Procedure Lab One

To perform the dilution procedure as illustrated in **Figure 39.2**, do the following:

1. Label all dilution tubes and media as follows:
  - a. Five tryptone soft agar tubes:  $10^{-5}$ ,  $10^{-6}$ ,  $10^{-7}$ ,  $10^{-8}$ ,  $10^{-9}$ .

- b. Five tryptone hard agar plates:  $10^{-5}$ ,  $10^{-6}$ ,  $10^{-7}$ ,  $10^{-8}$ ,  $10^{-9}$ .

- c. Nine tryptone broth tubes:  $10^{-1}$  through  $10^{-9}$ .

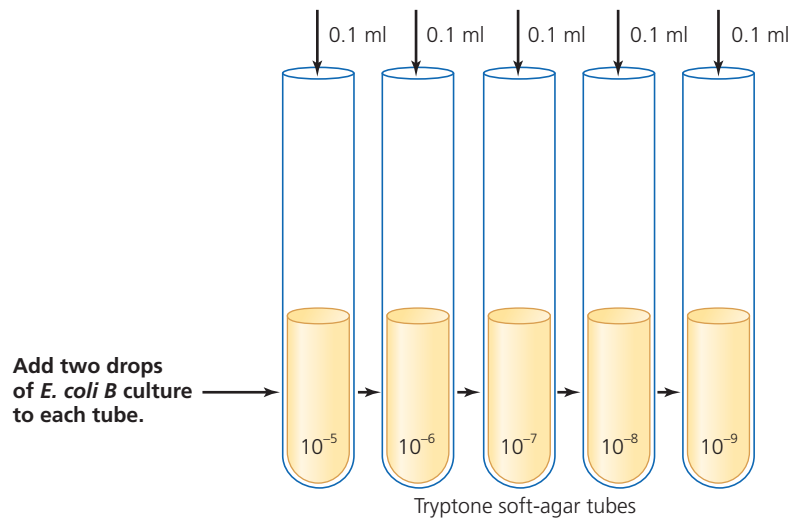
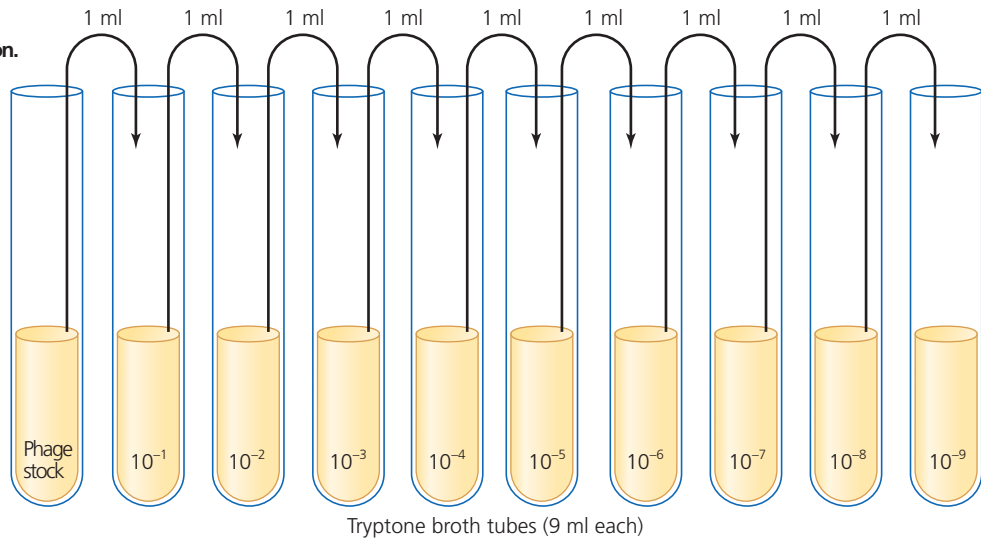
2. Place the five labeled soft tryptone agar tubes into a waterbath. Water should be of a depth just slightly above that of the agar in the tubes. Bring the waterbath to  $100^{\circ}\text{C}$  to melt the agar. Cool and maintain the melted agar at  $45^{\circ}\text{C}$ .
3. With 1-ml pipettes, aseptically perform a 10-fold serial dilution of the provided phage culture using the nine 9-ml tubes of tryptone.
4. To the tryptone soft agar tube labeled  $10^{-5}$ , aseptically add two drops of the *E. coli* B culture with a Pasteur pipette and 0.1 ml of the  $10^{-4}$  tryptone broth phage dilution. Rapidly mix by rotating the tube between the palms of your hands and pour the contents over the hard tryptone agar plate labeled  $10^{-5}$ , thereby forming a double-layered plate culture preparation. Swirl the plate gently and allow to harden.
5. Using separate Pasteur pipettes and 1-ml sterile pipettes, repeat Step 4 for the tryptone broth phage dilution tubes labeled  $10^{-5}$  through  $10^{-8}$  to effect the  $10^{-6}$  through  $10^{-9}$  tryptone soft agar overlays.
6. Following solidification of the soft agar overlay, incubate all plate cultures in an inverted position for 24 hours at  $37^{\circ}\text{C}$ .

## Procedure Lab Two


1. Observe all plates for the presence of plaque-forming units that develop on the bacterial lawn.
2. Count the number of PFUs in the range of 30 to 300 on each plate.
3. Calculate the number of phage particles per ml of the stock phage culture based on your PFU count.
4. Record your results in the chart in the Lab Report.

## PROCEDURE

Perform a 10-fold serial dilution.



Add two drops of *E. coli B* culture to each tube.

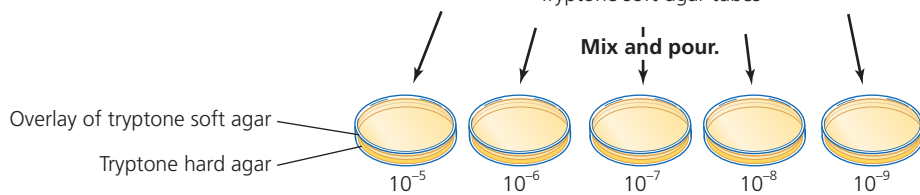

**Figure 39.2** Dilution procedure for cultivation and enumeration of bacteriophages



*This page intentionally left blank*

Name: \_\_\_\_\_


Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

Phage Dilution	Number of PFUs	Calculation: PFUs $\times$ Dilution Factor	PFUs/ml of Stock Phage Culture
$10^{-5}$			
$10^{-6}$			
$10^{-7}$			
$10^{-8}$			
$10^{-9}$			

## Review Questions

1. Discuss the effects of lytic and lysogenic infections on the life cycle of the host cell.
2. Discuss the factors responsible for the transformation of a lysogenic infection to one that is lytic.
3. Distinguish between the replicative and maturation stages of a lytic phage infection.

4. In this experimental procedure, why is it important to use a hard agar with a soft agar overlay technique to demonstrate plaque formation?
  
  
  
  
  
  
  
  
  
  
5. Explain what is meant by plaque-forming units.
  
  
  
  
  
  
  
  
  
  
6. Determine the number of PFUs per ml in a  $10^{-9}$  dilution of a phage culture that shows 204 PFUs in the agar lawn.
  
  
  
  
  
  
  
  
  
  
7.  The release of phage particles from the host bacterium always occurs by lysis of the cell and results in the death of the host. Animal viruses are released by either the lysis of the host cell or exocytosis, a reverse pinocytosis. Regardless of the mechanism of release, most infected cells die, while other viruses may escape the cell without damaging the host cell. Explain.

# Isolation of Coliphages from Raw Sewage

EXPERIMENT

# 40

## LEARNING OBJECTIVE

Once you have completed this experiment, you should be able to

1. Isolate virulent coliphages from sewage.

## Principle

Isolates of bacterial viruses (bacteriophages) can be obtained from a variety of natural sources, including soil, intestinal contents, raw sewage, and some insects such as cockroaches and flies. Their isolation from these environments is not an easy task because the phage particles are usually present in low concentrations. Therefore, isolation requires a series of steps:

1. Collection of the phage-containing sample at its source.
2. Addition of an enriched susceptible host cell culture to the sample to increase the number of phage particles for subsequent isolation.
3. Following incubation, centrifugation of the enriched sample for the removal of gross particles.
4. Filtration of the supernatant liquid through a bacteria-retaining membrane filter.
5. Inoculation of the bacteria-free filtrate onto a lawn of susceptible host cells grown on a soft agar plate medium.
6. Incubation and observation of the culture for the presence of phage particles, which is indicated by plaque formation in the bacterial lawn.

In the following experiment, you will use this procedure, as illustrated in **Figure 40.1**, for the isolation of *Escherichia coli* phage particles from raw sewage. Most bacteriophages that infect *E. coli* (coliphages) are designated by the letter T, indicating types. Seven types have been identified and are labeled T1 through T7. The T-even phages

(T2, T4, and T6) differ from the T-odd phages in that the former vary in size, form, and chemical composition. All of the T phages are capable of infecting the susceptible *E. coli* B host cell.

## CLINICAL APPLICATION

### Phage Therapy

Phage therapy is the therapeutic use of bacteriophages to treat pathogenic bacterial infections. It is mainly used in Russia and the Republic of Georgia and is not universally approved elsewhere. In the West, no phage therapies are authorized for use on humans, although phages for killing food poisoning bacteria (*Listeria*) are now in use. They may also be used as a possible therapy against many strains of drug-resistant bacteria.

## AT THE BENCH



## Materials

### Cultures

Lab One: 5-ml 24-hour broth cultures of *E. coli* B and 45-ml samples of fresh sewage collected in screw-capped bottles. Lab Two: 10-ml 24-hour broth cultures of *E. coli* B.


### Media

Per designated student group: Lab One: One 5-ml tube of bacteriophage nutrient broth, 10 times normal concentration. Lab Two: Five tryptone agar plates and five 3-ml tubes of tryptone soft agar.

### Equipment

Lab One: Sterile 250-ml Erlenmeyer flask and stopper. Lab Two: Sterile membrane filter apparatus, sterile 125-ml Erlenmeyer flask and stopper, 125-ml flask, 1000-ml beaker, centrifuge, Bunsen


## PROCEDURE

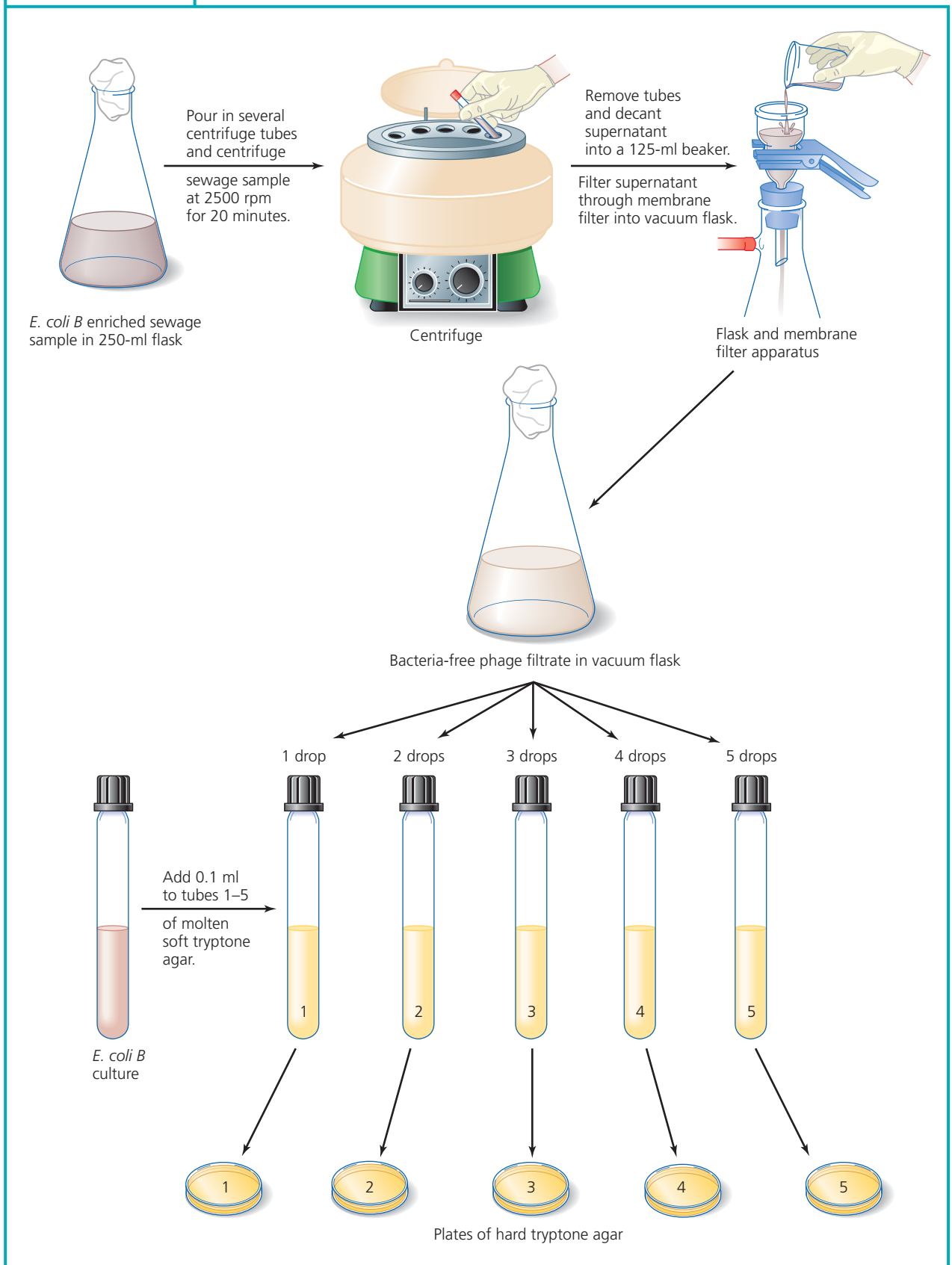

**Figure 40.1** Procedure for isolation of coliphages from raw sewage

burner, forceps, 1-ml sterile disposable pipettes, sterile Pasteur pipette, mechanical pipetting device, test tube rack, and glassware marking pencil.

## Procedure Lab One

 **Use disposable gloves. It is essential to handle raw sewage with extreme caution because it may serve as a vehicle for the transmission of human pathogens.**

### Enrichment of Sewage Sample

1. Aseptically add 5 ml of bacteriophage nutrient broth, 5 ml of the *E. coli* B broth culture, and 45 ml of the raw sewage sample to an appropriately labeled sterile 250-ml Erlenmeyer flask.
2. Incubate the culture for 24 hours at 37°C.

## Procedure Lab Two

### Filtration and Seeding

1. Following incubation, pour the phage-infected culture into a 100-ml centrifuge bottle or several centrifuge tubes and centrifuge at 2500 rpm for 20 minutes.
2. Remove the centrifuge bottle or tubes, being careful not to stir up the sediment, and carefully decant the supernatant into a 125-ml beaker.

3. Pour the supernatant solution through a sterile membrane filter apparatus to collect the bacteria-free, phage-containing filtrate in the vacuum flask below. Refer to Experiment 50 for the procedure in assembling the filter membrane apparatus.
4. Melt the soft tryptone agar by placing the five tubes in a boiling waterbath and cool to 45°C.
5. Label the five tryptone agar plates and the five tryptone agar tubes 1, 2, 3, 4, and 5, respectively.
6. Using a sterile 1-ml pipette, aseptically add 0.1 ml of the *E. coli* B culture to all the molten soft-agar tubes.
7. Using a sterile Pasteur pipette, aseptically add 1, 2, 3, 4, and 5 drops of the filtrate to the respectively labeled molten soft-agar tubes. Mix and pour each tube of soft agar into its appropriately labeled agar plate.
8. Allow agar to harden.
9. Incubate all the plates in an inverted position for 24 hours at 37°C.

## Procedure Lab Three

1. Examine all the culture plates for plaque formation, which is indicative of the presence of coliphages in the culture.
2. Indicate the presence (+) or absence (–) of plaques in each of the cultures in the chart in the Lab Report.

*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

Drops of Phage Filtrate	1	2	3	4	5
Plaque Formation (+) or (-)					


Based on your observations, what is the relationship between the number of plaques observed and the number of drops of filtrate in each culture?

## Review Questions

1. Why is enrichment of the sewage sample necessary for the isolation of phage?
2. How is enrichment of the sewage sample accomplished?



3. How are bacteria-free phage particles obtained?

4.  Why must you exercise caution when handling raw sewage samples?

# Physical and Chemical Agents for the Control of Microbial Growth

## LEARNING OBJECTIVES

Once you have completed the experiments in this section, you should

1. Know the basic methods for inhibiting microbial growth and the modes of antimicrobial action.
2. Be able to demonstrate the effects of physical agents, moist heat, osmotic pressure, and ultraviolet radiation on selected microbial populations.
3. Be able to demonstrate the effects on selected microbial populations of chemical agents used as disinfectants, antiseptics, and antibiotics.

## Introduction

Control of microorganisms is essential in the home, industry, and medical fields to prevent and treat diseases and to inhibit the spoilage of foods and other industrial products. Common methods of control involve chemical and physical agents that adversely affect microbial structures and functions, thereby producing a microbicidal or microbistatic effect. A **microbicidal effect** is one that kills the microbes immediately; a **microbistatic effect** inhibits the reproductive capacities of the cells and maintains the microbial population at a constant size.


## Chemical Methods for Control of Microbial Growth

1. **Antiseptics:** Chemical substances used on living tissue that kill or inhibit the growth of vegetative microbial forms.
2. **Disinfectants:** Chemical substances that kill or inhibit the growth of vegetative microbial forms on nonliving materials.
3. **Chemotherapeutic agents:** Chemical substances that destroy or inhibit the growth of microorganisms in living tissues.

## Physical Methods for Control of Microbial Growth

The modes of action of the different chemical and physical agents of control vary, although they all produce damaging effects to one or more essential cellular structures or molecules in order to cause cell death or inhibition of growth. Sites of damage that can result in malfunction are the cell wall, cell membrane, cytoplasm, enzymes, or nucleic acids. The adverse effects manifest themselves in the following ways.

1. **Cell-wall injury:** This can occur in one of two ways. First, lysis of the cell wall will leave the wall-less cell, called a protoplast, susceptible to osmotic damage, and a hypotonic environment may cause lysis of the vulnerable protoplast. Second, certain agents inhibit cell wall synthesis, which is essential during microbial cell reproduction. Failure to synthesize a missing segment of the cell wall results in an unprotected protoplast.
2. **Cell-membrane damage:** This may be the result of lysis of the membrane, which will cause immediate cell death. Also, the selective nature of the membrane may be affected

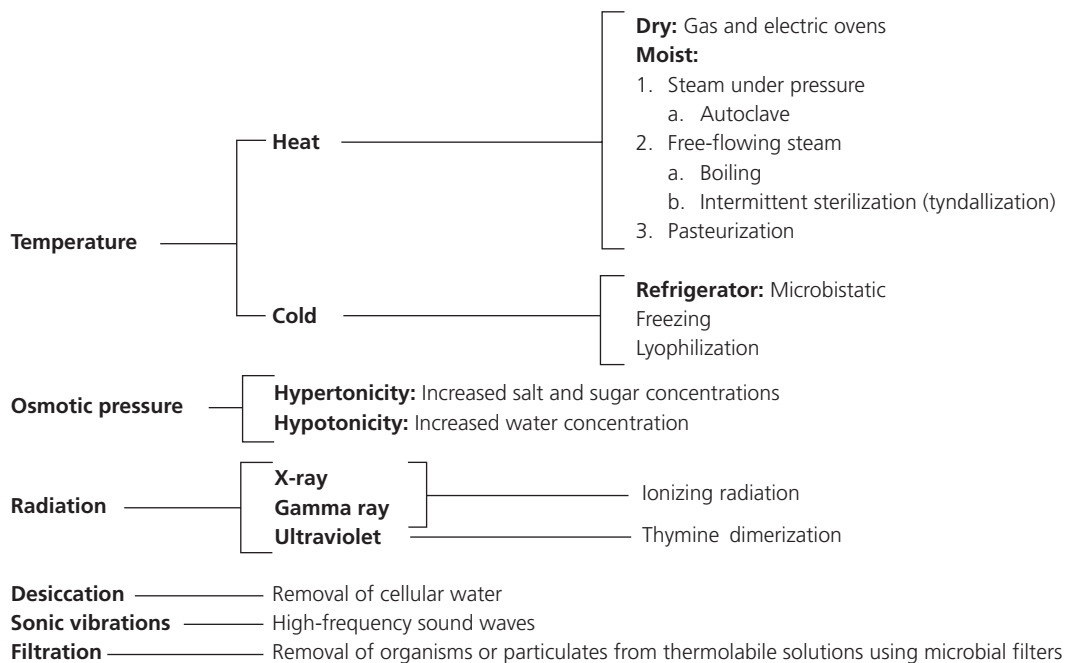

**Figure P9.1** Physical methods used for the control of microbial growth

without causing its complete disruption. As a result, there may be a loss of essential cellular molecules or interference with the uptake of nutrients. In both cases, metabolic processes will be adversely affected.

- Alteration of the colloidal state of cytoplasm:** Certain agents cause denaturing of cytoplasmic proteins. Denaturing processes are responsible for enzyme inactivation and cellular death by irreversibly rupturing the molecular bonds of these proteins and rendering them biologically inactive.
- Inactivation of cellular enzymes:** Enzymes may be inactivated competitively or non-competitively. Noncompetitive inhibition is irreversible and occurs following the application of some physical agent, such as mercuric chloride ( $\text{HgCl}_2$ ), that results in the uncoiling of the protein molecule, rendering it biologically inactive. Competitive inhibition occurs when a natural substrate is forced to compete for the active site on an enzyme surface with

a chemically similar molecular substrate, which can block the enzyme's ability to create end products. Competitive inhibitions are reversible.

- Interference with the structure and function of the DNA molecule:** The DNA molecule is the control center of the cell and may also represent a cellular target area for destruction or inhibition. Some agents have an affinity for DNA and cause breakage or distortion of the molecule, thereby interfering with its replication and role in protein synthesis.

**Figure P9.1** illustrates the acceptable physical methods used for the control of microbial growth.

Awareness of the mode of action of the physical and chemical agents is absolutely essential for their proper selection and application in microbial control. The exercises in this section are designed to acquaint you more fully with several commonly employed agents and their uses.

## LEARNING OBJECTIVE

Once you have completed this experiment, you should understand

1. The susceptibility of microbial species to destruction by the application of moist heat.

## Principle

Temperature has an effect on cellular enzyme systems and therefore a marked influence on the rate of chemical reactions and thus the life and death of microorganisms. Despite the diversity among microorganisms' temperature requirements for growth, extremes in temperature can be used in microbial growth control. Sufficiently low temperatures will inactivate enzymes and produce a static effect. High temperatures destroy cellular enzymes, which become irreversibly denatured.

The application of heat is a common means of destroying microorganisms. Both dry and moist heat are effective. However, moist heat, which (because of the hydrolyzing effect of water and its greater penetrating ability) causes coagulation of proteins, kills cells more rapidly and at lower temperatures than does dry heat. **Sterilization**, the destruction of all forms of life, is accomplished in 15 minutes at 121°C with moist heat (steam) under pressure; dry heat requires a temperature of 160°C to 180°C for 1½ to 3 hours.

Microbes exhibit differences in their resistance to moist heat. As a general rule, bacterial spores require temperatures above 100°C for destruction, whereas most bacterial vegetative cells are killed at temperatures of 60°C to 70°C in 10 minutes. Fungi can be killed at 50°C to 60°C, and fungal spores require 70°C to 80°C for 10 minutes for destruction. Because of this variability, moist heat can either sterilize or disinfect. Common applications include free-flowing steam under pressure (autoclaving), free-flowing steam at 100°C (tyndallization), and the use of lower temperatures (pasteurization).

Free-flowing steam under pressure requires the use of an autoclave, a double-walled metal vessel that allows steam to be pressurized in the outer jacket (see **Figure 41.1**). At a designated pressure, the saturated steam is released into the inner chamber, from which all the air has been evacuated. The steam under pressure in the vacuumed inner chamber is now capable of achieving temperatures in excess of 100°C. The temperature is determined by the pounds of pressure applied per square inch:

Pressure (pounds/inch <sup>2</sup> )	Temperature (°C)
0 (free-flowing steam)	100
10	115
15	121
20	126
25	130

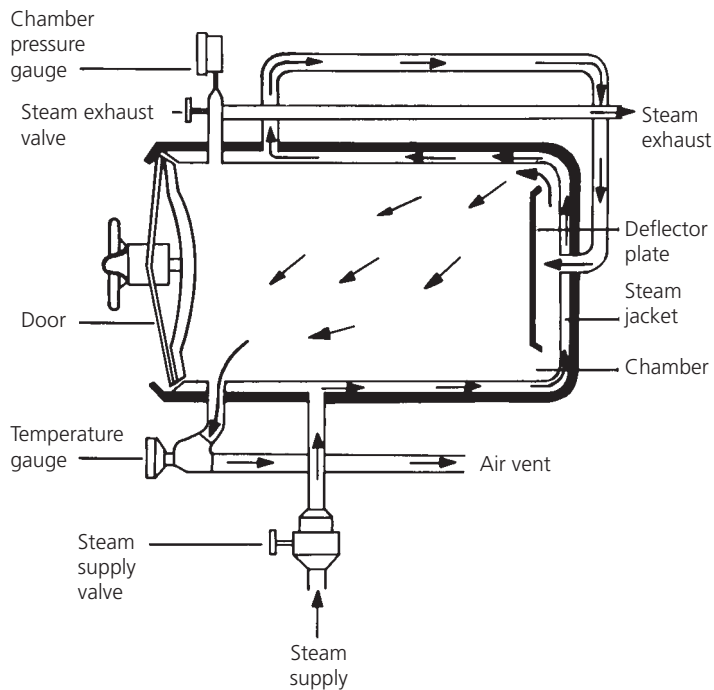
A pressure of 15 pounds/inch<sup>2</sup> achieves a temperature of 121°C and sterilizes in 15 minutes. This is the usual procedure; however, depending on the heat sensitivity of the material to be sterilized, the operating pressure and time conditions can be adjusted.

Application of **free-flowing steam** requires exposure of the contaminated substance to a temperature of 100°C, which is achieved by boiling water. Exposures to boiling water for 30 minutes will result in disinfection only; all vegetative cells will be killed, but not necessarily the more heat-resistant spores.

Another procedure is **tyndallization**, also referred to as intermittent or fractional sterilization. This procedure requires exposure of the material to free-flowing steam at 100°C for 20 minutes on 3 consecutive days with intermittent incubation at 37°C. The steaming kills all vegetative cells. Any spores that may be present germinate during the period of incubation and are



(a) An autoclave



(b) Schematic representation

**Figure 41.1** The autoclave

destroyed during subsequent exposure to a temperature of  $100^{\circ}\text{C}$ . Repeating this procedure for 3 days ensures germination of all spores and their destruction in the vegetative form. Because tyndallization requires so much time, it is used only for sterilization of materials that are composed of thermolabile chemicals and that might be subject to decomposition at higher temperatures.

**Pasteurization** exposes fairly thermolabile products such as milk, wine, and beer for a given period of time to a temperature that is high enough to destroy pathogens and some spoilage-causing microorganisms that may be present, without necessarily destroying all vegetative cells. There are three types of pasteurization: The high-temperature, short-time (HTST) procedure requires a temperature of  $71^{\circ}\text{C}$  for 15 seconds. The low-temperature, long-time (LTLT) method

requires  $63^{\circ}\text{C}$  for 30 minutes, and the ultra high temperature (UHT) approach occurs at  $138^{\circ}\text{C}$  for 2 seconds.

## CLINICAL APPLICATION

### Autoclave Performance Testing

While the original “autoclave” was invented as a pressure cooker for food, modern autoclaves are precision instruments and require maintenance and periodic testing, especially if control of human pathogens is involved. Commonly, a sample of spores of the bacterium *Bacillus stearothermophilus* is sterilized in the chamber with a normal load, and then the sample is allowed to incubate—any growth indicates that the autoclave needs to be serviced.



## Materials

### Cultures


48- to 72-hour nutrient broth cultures (50 ml per 250-ml Erlenmeyer flask) of *Staphylococcus aureus* and *Bacillus cereus*; 72- to 96-hour Sabouraud broth cultures (50 ml per 250-ml Erlenmeyer flask) of *Aspergillus niger* and *Saccharomyces cerevisiae*.

### Media

Per designated student group (pairs or groups of four): five nutrient agar plates, five Sabouraud agar plates, and one 10-ml tube of nutrient broth.

### Equipment

Bunsen burner, 800-ml beaker (waterbath), tripod and wire gauze screen with heat-resistant pad, thermometer, sterile test tubes, glassware marking pencil, and inoculating loop.

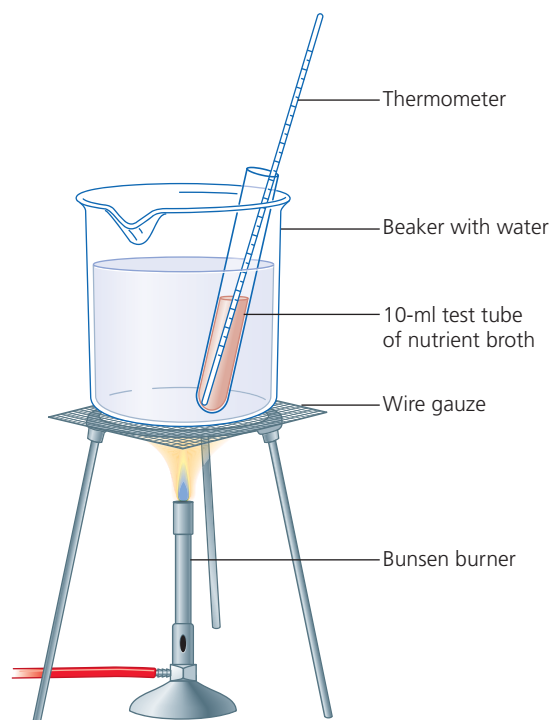

**Figure 41.2** Waterbath for moist heat experiment

## Procedure Lab One

1. Label the covers of each of the nutrient agar and Sabouraud agar plates, indicating the experimental heat temperatures to be used: 25°C (control), 40°C, 60°C, 80°C, and 100°C.
2. Score the underside of all plates with a glassware marking pencil into two sections. On the nutrient agar plates, label one section *S. aureus* and the other *B. cereus*. On the Sabouraud agar plates, label one section *A. niger* and the second *S. cerevisiae*.
3. Using aseptic technique, inoculate the nutrient agar and Sabouraud agar plates labeled 25°C by making a single-line loop inoculation of each test organism in its respective section of the plate.
4. Using a sterile pipette and mechanical pipetter, transfer 10 ml of each culture to four sterile test tubes labeled with the name of the organism and the temperature (40°C, 60°C, 80°C, and 100°C).
5. Set up the waterbath as illustrated in **Figure 41.2**, inserting the thermometer in an uncapped tube of nutrient broth.
6. Slowly heat the water to 40°C; check the thermometer frequently to ensure that it does not exceed the desired temperature. Place the four cultures of the experimental organisms into the beaker and maintain the temperature at 40°C for 10 minutes. Remove the cultures and aseptically inoculate each organism in its appropriate section on the two plates labeled 40°C.
7. Raise the waterbath temperature to 60°C and repeat Step 6 for the inoculation of the two plates labeled 60°C.
8. Raise the waterbath temperature to 80°C and repeat Step 6 for the inoculation of the two plates labeled 80°C.
9. Raise the waterbath temperature to 100°C and repeat Step 6 for the inoculation of the two plates labeled 100°C.
10. Incubate the nutrient agar plate cultures in an inverted position for 24 to 48 hours at 37°C and the Sabouraud agar plate cultures for 4 to 5 days at 25°C in a moist chamber.

## Procedure Lab Two

1. Observe all plates for the amount of growth of the test organisms at each of the temperatures.
2. Record your results in the chart provided in the Lab Report.

*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

- Record your results in the chart as 0 = none; 1+ = slight; 2+ = moderate; 3+ = abundant.

Microbial Species	AMOUNT OF GROWTH				
	25°C	40°C	60°C	80°C	100°C
<i>B. cereus</i>					
<i>S. aureus</i>					
<i>A. niger</i>					
<i>S. cerevisiae</i>					

- List the microbial organisms in order of their increasing heat resistance.

---



---



---



---

## Review Questions


- Account for the microbistatic effect produced by low temperatures as compared to the microbicidal effect produced by high temperatures.
- Cite the advantages of each of the modes of sterilization: tyndallization and autoclaving.



3. Discuss the detrimental effects of control agents on the following: the cytoplasm, the cell wall, nucleic acids, and the cell membrane.

4. Explain why milk is subjected to pasteurization rather than sterilization.

5.  *A. niger* and *B. cereus* cultures used in this experiment contained spores. Why is *B. cereus* more heat resistant?

6.  Account for the fact that aerobic and anaerobic bacterial spore-formers are more heat resistant than the tubercle bacillus, which is also known to tolerate elevated temperatures.

# Physical Agents of Control: Environmental Osmotic Pressure

## LEARNING OBJECTIVE

Once you have completed this experiment, you should understand

1. The possible effect of osmotic pressure environments on microorganisms.

## Principle

**Osmosis** is the net movement of water molecules (solvent) across a semipermeable membrane from a solution of their higher concentration to a solution of their lower concentration. The relative water concentrations of two solutions is determined by their solute concentrations. The **hypertonic solution** possesses a higher osmotic pressure and a higher solute concentration, and therefore it has a lower water concentration; it tends to draw in water. The **hypotonic solution** possesses a lower osmotic pressure and solute concentration, and therefore a higher water concentration; it tends to lose water. If two solutions have equal concentrations of solutes and therefore equal water concentrations, there is no osmosis and the solutions are **isotonic**.

The cell and its environment represent two solutions separated by the semipermeable cell membrane. The cell's cytoplasm contains colloidal and solute particles dispersed in water, as does the cell's environment. The **osmotic pressure of the environment** in relation to that of the cytoplasm of the cell plays a vital role in the life and death of a cell. The ideal environment of animal cells, which are bounded only by the fragile cell membrane, is completely or very nearly isotonic. In a hypertonic, high-pressure environment, all cells lose water by osmosis and become shriveled. This effect is called **plasmolysis**. As water is necessary for the occurrence of many chemical reactions, water loss adversely affects cell metabolism and reproduction. In a hypotonic, low-pressure environment, cells take in water and become swollen. This phenomenon is called **plasmoptysis**. In an environment with sufficiently low osmotic pressure, animal cells undergo **lysis**, which causes their deaths. Microorganisms and other organisms whose cells possess rigid cell walls are not usually susceptible

to lysis in hypotonic environments and usually prefer a slightly hypotonic environment to maintain themselves in a turgid (distended) state. The effects of environmental osmotic pressure on cells unprotected by cell walls are illustrated in **Figure 42.1**.

Osmotic pressure can be used as an antimicrobial agent. Microorganisms are not usually adversely affected by low environmental osmotic pressure because of their small sizes and the presence of rigid cell walls. However, hypertonicity is a commonly used method of inhibiting microbial growth.

Because of their varied habitats, microorganisms are generally well adapted to exist in all types of osmotic pressure environments. Different groups of microorganisms require different degrees of salinity for growth, and as a rule, they can adjust to salt concentrations of 0.5% to 3%. Concentrations of 10% to 15% are inhibitory to the growth of most microbes, except for **halophiles**, which require high salinity concentrations for growth. This sensitivity is the basis of food preservation by the process of salting.

## CLINICAL APPLICATION

### Honey and Wound Healing

Honey has been applied as a wound dressing since at least the time of ancient Egypt, but interest in its antibacterial properties is on the rise today. The exact mechanism of honey's antibacterial actions is not totally understood, and may include anti-inflammatory effects and pH levels inhospitable to bacteria, but its osmotic pressure and tendency to dehydrate bacteria is at least part of its effectiveness.


## AT THE BENCH



## Materials

### Cultures

24- to 48-hour nutrient broth cultures of *Staphylococcus aureus* and *Escherichia coli* and 48-hour *Halobacterium* salt broth culture of *Halobacterium salinarium*.

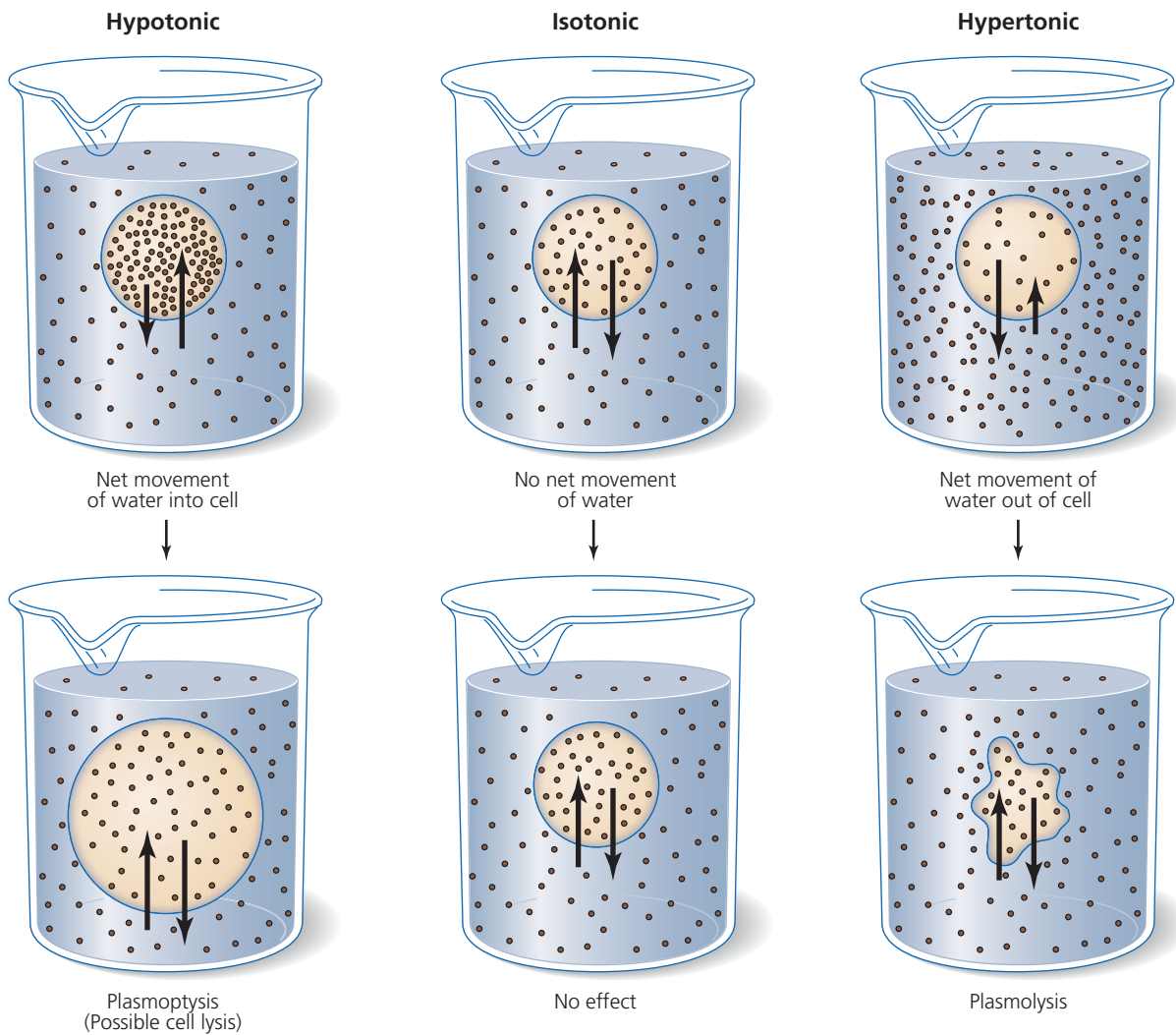

**Figure 42.1** Osmotic environments

## Media

Per student group: one nutrient agar plate of each of the following sodium chloride concentrations: 0.85%, 5%, 10%, 15%, and 25%.

## Equipment

Bunsen burner, inoculating loop, and glassware marking pencil.

## Procedure Lab One

1. Score the underside of each of the five nutrient agar plates into three sections with a glassware marking pencil.
2. Label each of the three sections on each plate with the name of the organism to be inoculated.

3. Using aseptic technique, inoculate each of the agar plates with the three experimental organisms by making a single-line loop inoculation of each organism in its labeled section.
4. Incubate all plates in an inverted position for 4 to 5 days at 25°C.

## Procedure Lab Two

1. Observe each of the nutrient agar plate cultures for the amount of growth on each of the experimental species.
2. Record your observations in the Lab Report.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

- Record your results in the chart as 0 = no growth; 1+ = scant growth; 2+ = moderate growth; 3+ = abundant growth.

Microbial Species	SALT CONCENTRATION OF MEDIUM (%)				
	0.85	5	10	15	25
<i>S. aureus</i>					
<i>E. coli</i>					
<i>H. salinarium</i>					


- For each of the experimental species, indicate (a) the salt concentration range in which growth will occur and (b) the optimal salt concentration.

Microbial Species	Range of Growth	Optimal NaCl Concentration
<i>S. aureus</i>		
<i>E. coli</i>		
<i>H. salinarium</i>		

## Review Questions

- Compare hypertonic, hypotonic, and isotonic solutions and their effects on cells.

2. Explain how hypertonicity can be used as a means of controlling microbial growth.

3.  Early pioneers, traveling to the West, salted their meat products for preservation. Explain the microbiological mechanism by which the salt retarded the spoilage of their meat.

4.  How do halophiles survive in their environment?

## LEARNING OBJECTIVE

Once you have completed this experiment, you should understand

1. The microbicidal effect of ultraviolet (UV) radiation on microorganisms.

## Principle


Certain forms of electromagnetic radiation are capable of producing a lethal effect on cells and therefore can be used for microbial control. Electromagnetic radiations that possess sufficient energy to be microbicidal are the short-wavelength radiations, that is, 300 nm and below. These include UV, gamma rays, and x-rays. The high-wavelength radiations, those above 300 nm, have insufficient energy to destroy cells. The electromagnetic spectrum and its effects on molecules are illustrated in **Figure 43.1**.

**Gamma radiation**, originating from unstable atomic nuclei, and **x-radiation**, originating from outside of the atomic nucleus, are representative of **ionizing forms of radiation**. Both transfer their energy through quanta (photons) to the matter through which they pass, causing excitation and the loss of electrons from molecules in their

paths. This injurious effect is nonspecific in that any molecule in the path of the radiation will undergo ionization. Essential cell molecules can be directly affected through loss of their chemical structures and activity brought about by the ionization. Also, water, the most abundant chemical constituent of cells, commonly undergoes radiation breakdown, with the ultimate production of highly reactive  $H^+$ ,  $OH^-$ , and, in the presence of oxygen,  $HO_2$  free radicals. These may combine with each other, frequently forming hydrogen peroxide ( $H_2O_2$ ), which is highly toxic to cells lacking catalase or other peroxidases, or the highly reactive free radicals may combine with any cellular constituents, again resulting in cell damage.

Because of their high energy content and therefore ability to penetrate matter, x-ray and gamma radiations can be used as means of sterilization, particularly of thermolabile materials. They are not commonly used, however, because of the expense of the equipment and the special facilities necessary for their safe use.

**Ultraviolet light**, which has a lower energy content than ionizing radiations, is capable of producing a lethal effect in cells exposed to the low penetrating wavelengths in the range of 210 nm to 300 nm. Cellular components capable of absorbing ultraviolet light are the nucleic acids; DNA is the primary site of damage. As the pyrimidines especially absorb ultraviolet wavelengths, the

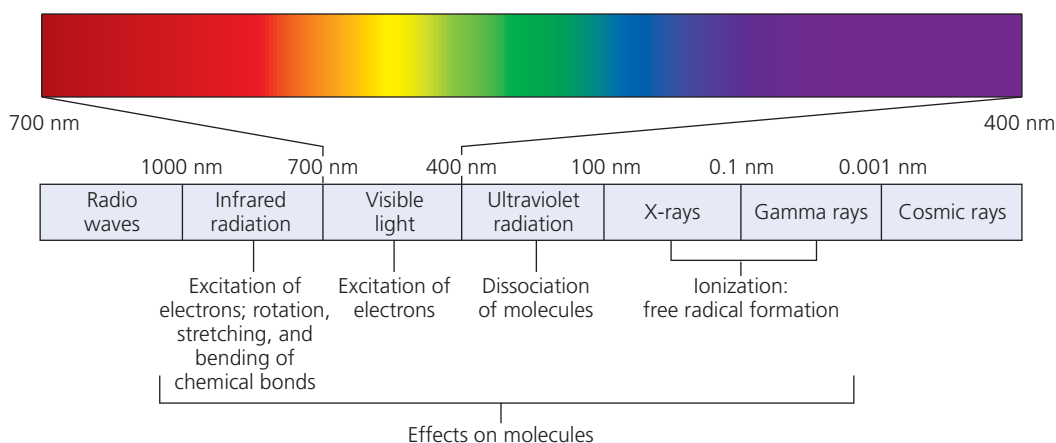

**Figure 43.1** The electromagnetic spectrum and its effects on molecules

major effect of this form of radiation is **thymine dimerization**, which is the covalent bonding of two adjacent thymine molecules on one nucleic acid strand in the DNA molecule. This dimer formation distorts the configuration of the DNA molecule, and the distortion interferes with DNA replication and transcription during protein synthesis.

Some cell types, including some microorganisms, possess enzyme systems for the repair of radiation-induced DNA damage. Two different systems may be operational: (1) The **excision repair system**, which functions in the absence of light; and (2) the **light repair system**, which is made operational by exposure of the irradiated cells to visible light in the wavelength range of 420 nm to 540 nm. The visible light serves to activate an enzyme that splits the dimers and reverses the damage.

Ultraviolet radiation, because of its low penetration ability, cannot be used as a means of sterilization, and its practical application is only for surface or air disinfection.

## CLINICAL APPLICATION

### Radiation Resistant Organisms

As lethal as radiation would seem to be, organisms vary tremendously in their ability to resist harm when exposed. Yeasts and bacterial spores are among the most resistant, while gram-positive bacteria are more resistant than gram-negative ones. Bacterial cocci also tend to be more resistant than bacilli, with viruses the most resistant of all. To effectively neutralize all contaminants, radiation dose and conditions must be carefully controlled.

## AT THE BENCH



## Materials

### Cultures

24- to 48-hour nutrient broth cultures of *Serratia marcescens* and *Bacillus cereus*; sterile saline spore suspension of *Aspergillus niger*.

### Media

Per student group: seven nutrient agar plates.

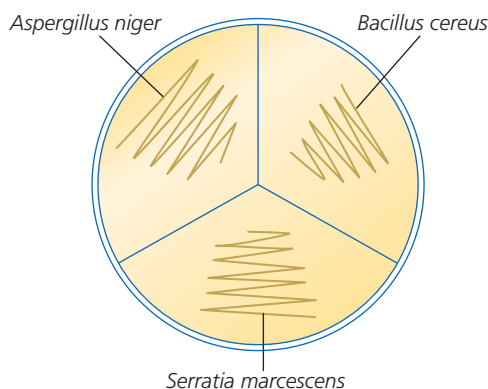
### Equipment

Bunsen burner, inoculating loop, ultraviolet radiation source (254 nm), and glassware marking pencil.

## Procedure Lab One

**⚠️ Wear disposable gloves and do not expose your eyes to the ultraviolet light source.**

1. Divide all nutrient agar plates into three sections by scoring the underside of each plate with a glassware marking pencil.
2. Label each of the sections on each plate with the name of the organism to be inoculated.
3. Using aseptic technique, inoculate all the plates by means of a streak inoculation *specifically* as shown in the following illustration:



4. Label the cover of each inoculated plate with the exposure time to ultraviolet radiation as 0 (control), 15 seconds, 30 seconds, 1 minute, and 3 minutes. Label two plates as 5 minutes; one of these plates will serve as the irradiated, covered control.
5. Irradiate all inoculated plates for the designated period of time by placing them 12 inches below the ultraviolet light source. Make sure first to remove all Petri dish covers except that of the 5-minute irradiated control plate.
6. Incubate all plates in an inverted position for 4 to 5 days at 25°C.

## Procedure Lab Two

1. Observe each of the nutrient agar plate cultures for the amount of growth on each of the microbial species.
2. Record your observations in the chart provided in the Lab Report.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

Record your observations in the chart as 0 = no growth; 1+ = slight growth; 2+ = moderate growth; 3+ = abundant growth.

Microbial Species	TIME OF IRRADIATION						
	Seconds			Minutes			
	0	15	30	1	3	5	5*
<i>B. cereus</i>							
<i>S. marcescens</i>							
<i>A. niger</i>							


\*Irradiated, covered plate.


## Review Questions

1. Discuss the effects of ionizing radiation on cellular constituents.
2. Explain why x-rays can be used for sterilization, whereas ultraviolet rays can be used only for surface disinfection of materials.



3. Explain the mechanism of action of ultraviolet radiation on cells.

4.  Account for the greater susceptibility of *S. marcescens* than that of *B. cereus* to the effects of ultraviolet radiation.

5.  Why is it not essential to shield your hands from ultraviolet light, whereas you must exercise great care to shield your eyes from this type of radiation?

Chemotherapeutic agents are chemical substances used in the treatment of infectious diseases. Their mode of action is to interfere with microbial metabolism, thereby producing a bacteriostatic or bactericidal effect on the microorganisms, without producing a like effect in host cells.

Chemotherapeutic agents act on a number of cellular targets. Their mechanisms of action include inhibition of cell-wall synthesis, inhibition of protein synthesis, inhibition of nucleic acid synthesis, disruption of the cell membrane, and inhibition of folic acid synthesis. These drugs can be separated into two categories:

1. **Antibiotics** are synthesized and secreted by some true bacteria, actinomycetes, and fungi that destroy or inhibit the growth of other


microorganisms. Today, some antibiotics are laboratory synthesized or modified; however, their origins are living cells.

2. **Synthetic drugs** are synthesized in the laboratory.

To determine a therapeutic drug of choice, one must know its mode of action, possible adverse side effects in the host, and the scope of its antimicrobial activity. The specific mechanism of action varies among different drugs, and the short-term or long-term use of many drugs can produce systemic side effects in the host. These vary in severity from mild and temporary upsets to permanent tissue damage (Table 44.1).

**TABLE 44.1 Prototypic Antibiotics**

ANTIBIOTIC	MODE OF ACTION	POSSIBLE SIDE EFFECTS
<b>Penicillin</b>	Prevents transpeptidation of the <i>N</i> -acetylmuramic acids, producing a weakened peptidoglycan structure	Penicillin resistance; sensitivity (allergic reaction)
<b>Streptomycin</b>	Has an affinity for bacterial ribosomes, causing misreading of codons on mRNA, thereby interfering with protein synthesis	May produce damage to auditory nerve, causing deafness
<b>Chloramphenicol</b>	Has an affinity for bacterial ribosomes, preventing peptide bond formation between amino acids during protein synthesis	May cause aplastic anemia, which is fatal because of destruction of RBC-forming and WBC-forming tissues
<b>Tetracyclines</b>	Have an affinity for bacterial ribosomes; prevent hydrogen bonding between the anticodon on the tRNA–amino acid complex and the codon on mRNA during protein synthesis	Permanent discoloration of teeth in young children
<b>Bacitracin</b>	Inhibits cell-wall synthesis	Nephrotoxic if taken internally; used for topical application only
<b>Polymyxin</b>	Destruction of cell membrane	Toxic if taken internally; used for topical application only
<b>Rifampin</b>	Inhibits RNA synthesis	Appearance of orange-red urine, feces, saliva, sweat, and tears
<b>Quinolone</b>	Inhibits DNA synthesis	Affects the development of cartilage

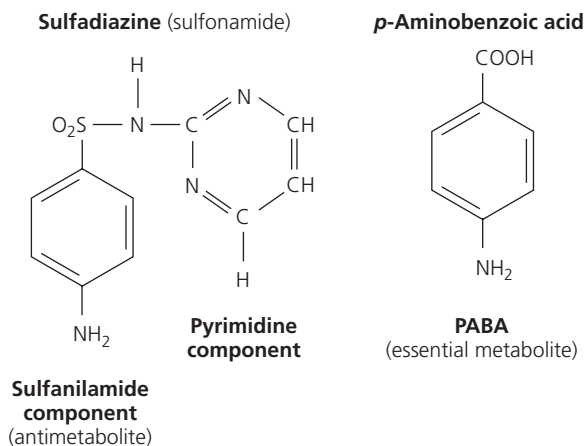

**Figure 44.1** Chemical similarity of sulfanilamide and PABA

## Synthetic Agents

Sulfadiazine (a sulfonamide) produces a static effect on a wide range of microorganisms by a mechanism of action called **competitive inhibition**. The active component of the drug, sulfanilamide, acts as an **antimetabolite** that competes with the **essential metabolite**, p-aminobenzoic acid (PABA), during the synthesis of folic acid in the microbial cell. Folic acid is an essential cellular coenzyme involved in the synthesis of amino acids and purines. Many microorganisms possess enzymatic pathways for folic acid synthesis and can be adversely affected by sulfonamides. Human cells lack these enzymes, and the essential folic acid enters the cells in a preformed state. Therefore, these drugs have no competitive effect on human cells. The similarity between the chemical structure of the antimetabolite sulfanilamide and the structure of the essential metabolite PABA is illustrated in **Figure 44.1**.

## PART A The Kirby-Bauer Antibiotic Sensitivity Test Procedure

### LEARNING OBJECTIVE

*Once you have completed this experiment, you should understand*

1. The Kirby-Bauer procedure for the evaluation of the antimicrobial activity of chemotherapeutic agents.

## Principle


The available chemotherapeutic agents vary in their scope of antimicrobial activity. Some have a limited spectrum of activity, being effective against only one group of microorganisms. Others exhibit broad-spectrum activity against a range of microorganisms. The drug susceptibilities of many pathogenic microorganisms are known, but it is sometimes necessary to test several agents to determine the drug of choice.

A standardized diffusion procedure with filter-paper discs on agar, known as the **Kirby-Bauer method**, is frequently used to determine the drug susceptibility of microorganisms isolated from infectious processes. This method allows the rapid determination of the efficacy of a drug by measuring the diameter of the zone of inhibition that results from diffusion of the agent into the medium surrounding the disc. In this procedure, filter-paper discs of uniform size are impregnated with specified concentrations of different antibiotics and then placed on the surface of an agar plate that has been seeded with the organism to be tested. The medium of choice is Mueller-Hinton agar, with a pH of 7.2 to 7.4, which is poured into plates to a uniform depth of 5 mm and refrigerated after solidification. Prior to use, the plates are transferred to an incubator at 37°C for 10 to 20 minutes to dry off the moisture that develops on the agar surface. The plates are then heavily inoculated with a standardized inoculum by means of a cotton swab to ensure the confluent growth of the organism. The discs are aseptically applied to the surface of the agar plate at well-spaced intervals. Once applied, each disc is gently touched with a sterile applicator stick to ensure its firm contact with the agar surface.

Following incubation, the plates are examined for the presence of growth inhibition, which is indicated by a clear zone surrounding each disc (**Figure 44.2**). The susceptibility of an organism to a drug is assessed by the size of this zone, which is affected by other variables such as

1. The ability and rate of diffusion of the antibiotic into the medium and its interaction with the test organism.
2. The number of organisms inoculated.
3. The growth rate of the organism.

A measurement of the diameter of the zone of inhibition in millimeters is made, and its size is compared to that contained in a standardized

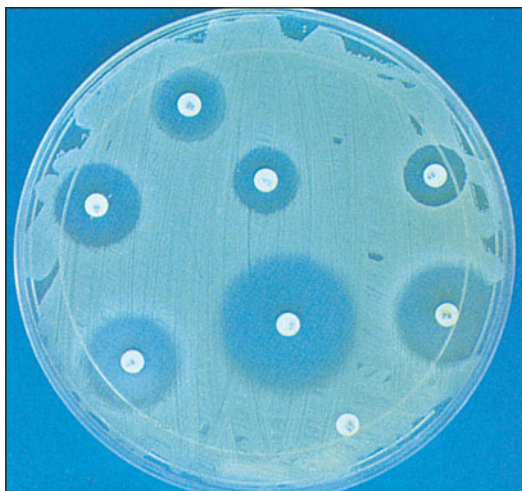

**Figure 44.2** Kirby-Bauer antibiotic sensitivity test

chart, which is shown in Table 44.2 on page 299. Based on this comparison, the test organism is determined to be resistant, intermediate, or susceptible to the antibiotic.

## CLINICAL APPLICATION

### Selection of Effective Antibiotics

Upon isolation of an infectious agent, a chemotherapeutic agent is selected and its effectiveness must be determined. This can be done using the Kirby-Bauer Antibiotic Sensitivity Test. This is the essential tool used in clinical laboratories to select the best agent to treat patients with bacterial infections.

## AT THE BENCH



## Materials

### Cultures

0.85% saline suspensions of *Escherichia coli*, *Staphylococcus aureus*, *Pseudomonas aeruginosa*, *Proteus vulgaris*, *Mycobacterium smegmatis*, *Bacillus cereus*, and *Enterococcus faecalis* adjusted to an absorbance of 0.1 at 600 nm. *Note: For enhanced growth of M. smegmatis, add Tween™ 80 (1 mL/liter of broth medium) and incubate for 3 to 5 days in a shaking waterbath, if available.*

### Media

Per designated student group: seven Mueller-Hinton agar plates.

## Antimicrobial-Sensitivity Discs

Penicillin G, 10  $\mu\text{g}$ ; streptomycin, 10  $\mu\text{g}$ ; tetracycline, 30  $\mu\text{g}$ ; chloramphenicol, 30  $\mu\text{g}$ ; gentamicin, 10  $\mu\text{g}$ ; vancomycin, 30  $\mu\text{g}$ ; and sulfanilamide, 300  $\mu\text{g}$ .

## Equipment

Sensi-Disc™ dispensers or forceps, Bunsen burner, sterile cotton swabs, glassware marking pencil, and millimeter ruler.


## Procedure Lab One

1. Place agar plates right side up in an incubator heated to 37°C for 10 to 20 minutes with the covers adjusted so that the plates are slightly opened.
2. Label the covers of each of the agar plates with the name of the test organism to be inoculated.
3. Using aseptic technique, inoculate all agar plates with their respective test organisms as follows:
  - a. Dip a sterile cotton swab into a well-mixed saline test culture and remove excess inoculum by pressing the saturated swab against the inner wall of the culture tube.
  - b. Using the swab, streak the entire agar surface horizontally, vertically, and around the outer edge of the plate to ensure a heavy growth over the entire surface.
4. Allow all culture plates to dry for about 5 minutes.
5. Using the Sensi-Disc dispenser, apply the antibiotic discs by placing the dispenser over the agar surface and pressing the plunger, depositing the discs simultaneously onto the agar surface (Figure 44.3, Step 1a). Or, if dispensers are not available, distribute the individual discs at equal distances with forceps dipped in alcohol and flamed (Figure 44.3, Step 1b).
6. Gently press each disc down with the wooden end of a cotton swab or with sterile forceps to ensure that the discs adhere to the surface of the agar (Figure 44.3, Step 2). *Note: Do not press the discs into the agar.*
7. Incubate all plate cultures in an inverted position for 24 to 48 hours at 37°C.

## Procedure Lab Two

1. Examine all plate cultures for the presence or absence of a zone of inhibition surrounding each disc.

## PROCEDURE

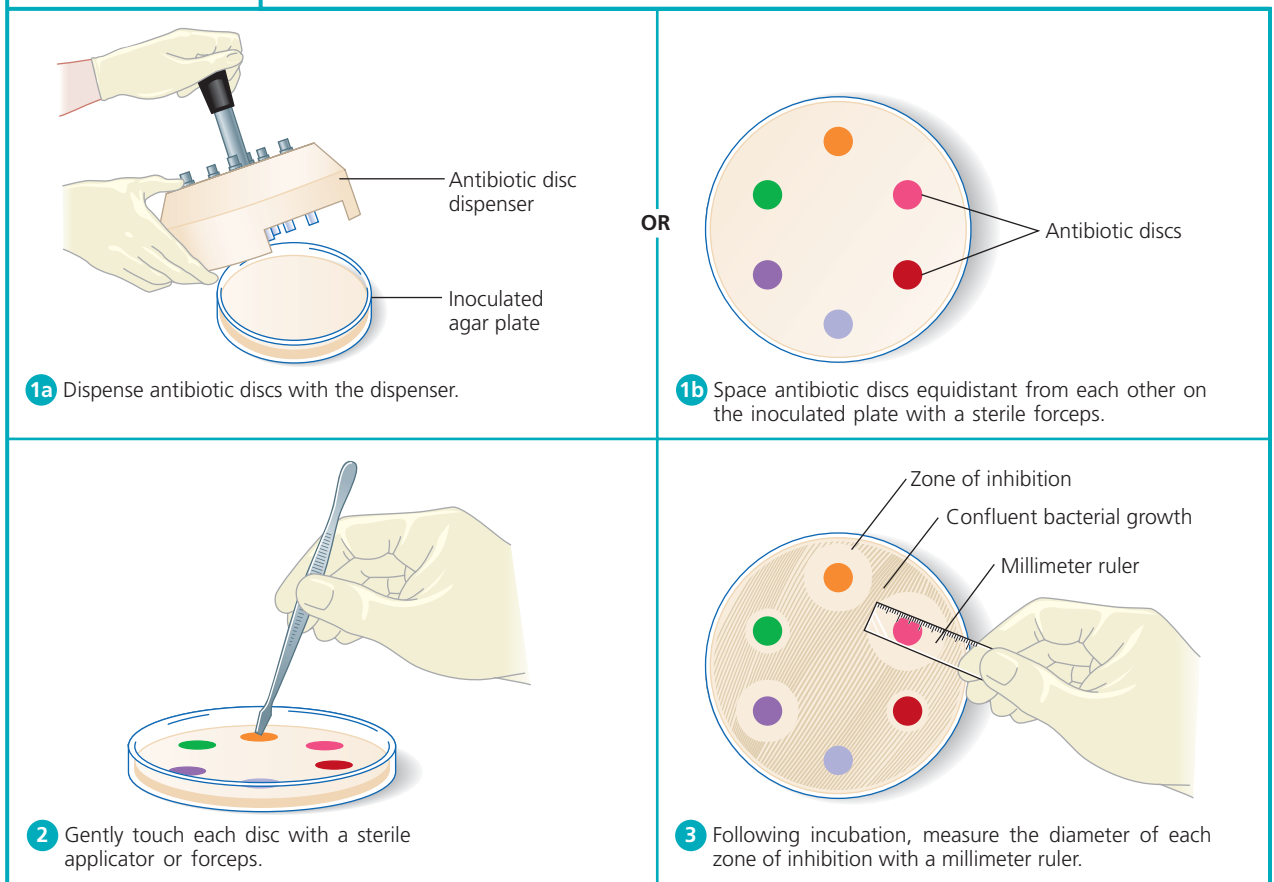

**Figure 44.3** Kirby-Bauer antibiotic sensitivity procedure

- Using a ruler graduated in millimeters, carefully measure each zone of inhibition to the nearest millimeter (Figure 44.3, Step 3). Record your results in the chart provided in the Lab Report.
- Compare your results with Table 44.2 and determine the susceptibility of each test organism to the chemotherapeutic agent. Record your results in the Lab Report.

## PART B Synergistic Effect of Drug Combinations

### LEARNING OBJECTIVE

Once you have completed this experiment, you should be able to

- Perform the disc-agar diffusion technique for determination of synergistic combinations of chemotherapeutic agents.

## Principle

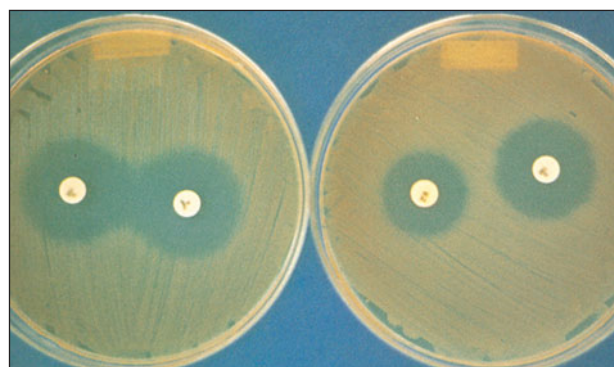
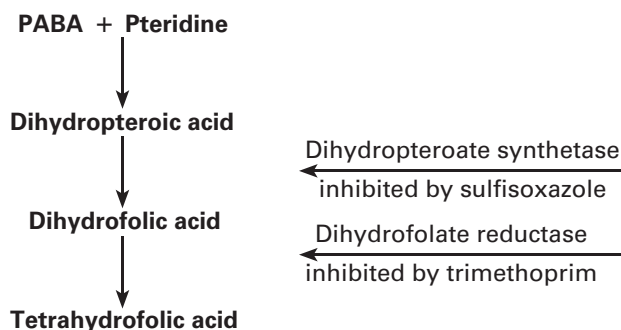
Combination chemotherapy, the use of two or more antimicrobial or antineoplastic agents, is being employed in medical practice with ever-increasing frequency. The rationale for using drug combinations is the expectation that effective combinations might lower the incidence of bacterial resistance, reduce host toxicity of the antimicrobial agents (because of decreased dosage requirements), or enhance the agents' bactericidal activity. Enhanced bactericidal activity is known as **synergism**. Synergistic activity is evident when the sum of the effects of the chemotherapeutic agents used in combination is significantly greater than the sum of their effects when used individually. This result is readily differentiated from an **additive (indifferent) effect**, which is evident when the interaction of two drugs produces a combined effect that is no greater than the sum of their separately measured individual effects.

A variety of *in vitro* methods are available to demonstrate synergistic activity. In this

experiment, a disc–agar diffusion technique will be performed to demonstrate this phenomenon. This technique uses the Kirby-Bauer antibiotic susceptibility test procedure, as described in Part A of this experiment, and requires both Mueller-Hinton agar plates previously seeded with the test organisms and commercially prepared, antimicrobial-impregnated discs. The two discs, representing the drug combination, are placed on the inoculated agar plate and separated by a distance (measured in mm) that is equal to or slightly greater than one-half the sum of their individual zones of inhibition when obtained separately. Following the incubation period, an additive effect is exhibited by the presence of two distinctly separate circles of inhibition. If the drug combination is synergistic, the two inhibitory zones merge to form a “bridge” at their juncture, as illustrated in **Figure 44.4**.

The following drug combinations will be used in this experimental procedure:

- 1. Sulfisoxazole, 150  $\mu\text{g}$ , and trimethoprim, 5  $\mu\text{g}$ .** Both antimicrobial agents are enzyme inhibitors that act sequentially in the metabolic pathway, leading to folic acid synthesis. The antimicrobial effect of each drug is enhanced when the two drugs are used in combination. The pathway thus exemplifies synergism.



(a) Synergistic effect

(b) Additive effect

**Figure 44.4** Synergistic and additive effects of drug combinations

- 2. Trimethoprim, 5  $\mu\text{g}$ , and tetracycline, 30  $\mu\text{g}$ .** The modes of antimicrobial activity of these two chemotherapeutic agents differ; tetracycline acts to interfere with protein synthesis at the ribosomes. Thus, when used in combination, these drugs produce an additive effect.

## CLINICAL APPLICATION

### Multiple Drug Therapy

In antimicrobial therapy for drug-resistant bacteria such as the opportunistic pathogen *Pseudomonas aeruginosa*, multiple drugs may be used to take advantage of synergistic effects. Research has shown that use of ampicillin to degrade gram-negative cell walls allows easier entry of kanamycin, which then inhibits protein synthesis. Combination therapies taking advantage of synergism also allow use of lower doses of each drug, which reduces overall toxic effects on the patient.

## AT THE BENCH



## Materials

### Cultures

0.85% saline suspensions of *Escherichia coli* and *Staphylococcus aureus* adjusted to an absorbance of 0.1 at 600 nm.

### Media

Per designated student group: four Mueller-Hinton agar plates.

### Antimicrobial-Sensitivity Discs

Tetracycline, 30  $\mu\text{g}$ ; trimethoprim, 5  $\mu\text{g}$ ; and sulfisoxazole, 150  $\mu\text{g}$ .

### Equipment

Bunsen burner, forceps, sterile cotton swabs, millimeter ruler, and glassware marking pencil.

## Procedure Lab One

- To inoculate the Mueller-Hinton agar plates, follow Steps 1 through 4 as described under the procedure in Part A of this experiment.

2. Using the millimeter ruler, determine the center of the underside of each plate and mark with a glassware marking pencil.
3. Using the glassware marking pencil, mark the underside of each agar plate culture at both sides from the center mark at the distances specified below:
  - a. *E. coli*-inoculated plate for trimethoprim and sulfisoxazole combination sensitivity: 12.5 mm on each side of center mark.
  - b. *S. aureus*-inoculated plate for trimethoprim and sulfisoxazole combination sensitivity: 14.5 mm on each side of center mark.
  - c. *E. coli*-inoculated plate for trimethoprim and tetracycline combination sensitivity: 14.0 mm on each side of center mark.
  - d. *S. aureus*-inoculated plate for trimethoprim and tetracycline combination: 14.0 mm on each side of center mark.
4. Using sterile forceps, place the antimicrobial discs, in the combinations specified in Step 3, onto the surface of each agar plate culture at the previously marked positions. Gently press each disc down with the sterile forceps to ensure that it adheres to the agar surface.
5. Incubate all plate cultures in an inverted position for 24 to 48 hours at 37°C.

## Procedure Lab Two

1. Examine all agar plate cultures to determine the zone of inhibition patterns exhibited. Distinctly separate zones of inhibition are indicative of an additive effect, whereas a merging of the inhibitory zones is indicative of synergism.
2. Record your observations and results in the chart provided in the Lab Report.

**TABLE 44.2 Zone Diameter Interpretive Standards for Organisms Other Than *Haemophilus* and *Neisseria gonorrhoeae***

ANTIMICROBIAL AGENT	DISC CONTENT	ZONE DIAMETER, NEAREST WHOLE mm		
		RESISTANT	INTERMEDIATE	SUSCEPTIBLE
<b>Ampicillin</b>				
when testing gram-negative bacteria	10 µg	≤ 13	14–16	≥ 17
when testing gram-positive bacteria	10 µg	≤ 28	—	≥ 29
<b>Carbenicillin</b>				
when testing <i>Pseudomonas</i>	100 µg	≤ 13	14–16	≥ 17
when testing other gram-negative organisms	100 µg	≤ 19	20–22	≥ 23
Cefoxitin	30 µg	≤ 14	15–17	≥ 18
Cephalothin	30 µg	≤ 14	16–17	≥ 18
Chloramphenicol	30 µg	≤ 12	13–17	≥ 18
Clindamycin	2 µg	≤ 14	15–20	≥ 21
Erythromycin	15 µg	≤ 13	14–22	≥ 23
Gentamicin	10 µg	≤ 12	13–14	≥ 15
Kanamycin	30 µg	≤ 13	14–17	≥ 18
Methicillin when testing staphylococci	5 µg	≤ 9	10–13	≥ 14
Novobiocin	30 µg	≤ 17	18–21	≥ 22
<b>Penicillin G</b>				
when testing staphylococci	10 units	≤ 28	—	≥ 29
when testing other bacteria	10 units	≤ 14	—	≥ 15
Rifampin	5 µg	≤ 16	17–19	≥ 20
Streptomycin	10 µg	≤ 11	12–14	≥ 15
Tetracycline	30 µg	≤ 14	15–18	≥ 19
Tobramycin	10 µg	≤ 12	13–14	≥ 15
Trimethoprim/sulfamethoxazole	1.25/23.75 µg	≤ 10	11–15	≥ 16
<b>Vancomycin</b>				
when testing enterococci	30 µg	≤ 14	15–16	≥ 17
when testing <i>Staphylococcus</i> spp.	30 µg	—	—	≥ 15
Sulfonamides	250 or 300 µg	≤ 12	—	≥ 17
Trimethoprim	5 µg	≤ 10	—	≥ 16

Source: Clinical and Laboratory Standards Institute. *Performance Standards for Antimicrobial Disk Susceptibility Tests*, Tenth Edition, 2008.



*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

### PART A: Kirby-Bauer Antibiotic Sensitivity Test Procedure

- Record the zone size and the susceptibility of each test organism to the chemotherapeutic agent as resistant (R), intermediate (I), or sensitive (S) in the charts below.

	GRAM-NEGATIVE						ACID-FAST	
	<i>E. coli</i>		<i>P. aeruginosa</i>		<i>P. vulgaris</i>		<i>M. smegmatis</i>	
Chemotherapeutic Agent	Zone Size	Susceptibility	Zone Size	Susceptibility	Zone Size	Susceptibility	Zone Size	Susceptibility
Penicillin								
Streptomycin								
Tetracycline								
Chloramphenicol								
Gentamicin								
Vancomycin								
Sulfanilamide								

	GRAM-POSITIVE					
	<i>S. aureus</i>		<i>E. faecalis</i>		<i>B. cereus</i>	
Chemotherapeutic Agent	Zone Size	Susceptibility	Zone Size	Susceptibility	Zone Size	Susceptibility
Penicillin						
Streptomycin						
Tetracycline						
Chloramphenicol						
Gentamicin						
Vancomycin						
Sulfanilamide						


2. For each of the chemotherapeutic agents, indicate the following:
  - a. The spectrum of its activity as broad or limited.
  - b. The type or types of organisms it is effective against as gram-positive, gram-negative, or acid-fast.

Chemotherapeutic Agent	Spectrum of Activity	Type(s) of Microorganisms
Penicillin		
Streptomycin		
Tetracycline		
Chloramphenicol		
Gentamicin		
Vancomycin		
Sulfanilamide		

## PART B: Synergistic Effect of Drug Combinations

Cultures	Appearance of Zone Inhibition	Synergistic or Additive Effect
<i>E. coli</i> :		
trimethoprim and sulfisoxazole	_____	_____
trimethoprim and tetracycline	_____	_____
<i>S. aureus</i> :		
trimethoprim and sulfisoxazole	_____	_____
trimethoprim and tetracycline	_____	_____

## Review Question

1.  Your experimental results indicate that antibiotics such as tetracycline, streptomycin, and chloramphenicol have a broad spectrum of activity against prokaryotic cells. Why do these antibiotics lack inhibitory activity against eukaryotic cells such as fungi?

# Determination of Penicillin Activity in the Presence and Absence of Penicillinase

EXPERIMENT

# 45

## LEARNING OBJECTIVES


Once you have completed this experiment, you should be able to

1. Employ a broth culture system for the determination of the minimal inhibitory concentration (MIC) of penicillin.
2. Demonstrate the reversal of penicillin inhibition against the test organism in the presence of penicillinase ( $\beta$ -lactamase).

## Principle

In addition to the Kirby-Bauer paper disc–agar diffusion procedure, the broth tube dilution method may be used to determine the susceptibility of an organism to an antibiotic. The latter procedure, in which dilutions of the antibiotic are prepared in the broth medium, also permits the **minimal inhibitory concentration (MIC)** to be determined for the antibiotic under investigation. The **MIC** is the lowest concentration of an antimicrobial agent that inhibits the growth of the test microorganism. Quantitative data of this nature may be used by a clinician to establish effective antimicrobial regimens for the treatment of a bacterial infection in a host. This data is of particular significance when the toxicity of the antibiotic is known to produce major adverse effects in host tissues.

Penicillin is a potent antibiotic produced by the mold *Penicillium chrysogenum* (formerly called *P. notatum*). Sir Alexander Fleming's discovery of penicillin in 1928 provided the world with the first clinically useful antibiotic in the fight to control human infection. The activity of this antibiotic, as illustrated in **Figure 45.1**, is associated with the  $\beta$ -lactam ring within its molecular structure. Shortly after the clinical introduction of benzylpenicillin (penicillin G), pathogenic organisms such as *Staphylococcus aureus* were found to be resistant to this “wonder drug.” Research

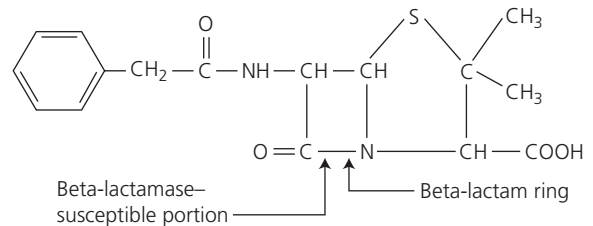


**Figure 45.1** Molecular structure of benzylpenicillin (penicillin G)



**Figure 45.2** Penicillinase activity. Penicillin sensitivity is shown on the left; penicillin resistance is shown on the right.

revealed that some organisms were genetically capable of producing  $\beta$ -lactamase (penicillinase), an enzyme that breaks a bond in the  $\beta$ -lactam ring portion of the molecule. When the integrity of this ring is compromised, the inhibitory activity of the antibiotic is lost. Penicillinase activity is illustrated in **Figure 45.2**.

In this experiment, the MIC of penicillin will be determined against penicillin-sensitive and penicillinase-producing strains of *Staphylococcus*

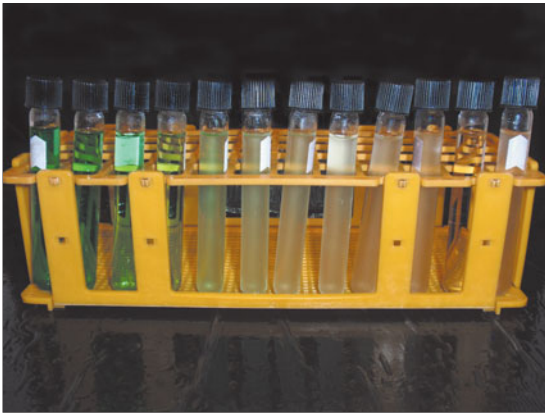

**Figure 45.3** Minimal inhibitory concentration tube set-up

*aureus*. The procedure to be followed involves specific concentrations of the penicillin prepared by means of a twofold serial dilution technique in an enriched broth medium. The tubes containing the antibiotic dilutions are then inoculated with a standardized concentration of the test organism (Figure 45.3). Table 45.1 illustrates the protocol for the antibiotic serial dilution–broth medium setup.

Following incubation, spectrophotometric absorbance readings will be used to determine the presence or absence of growth in the cultures. The culture that shows no growth in the presence of the lowest concentration of penicillin represents the minimal inhibitory concentration of this antibiotic against *S. aureus*.

## CLINICAL APPLICATION

### Wider Capability Seen in $\beta$ -lactamases

Penicillinases are  $\beta$ -lactam ring breakers with specific activity against penicillin, while cephalosporins are generally not affected by them. New gene variants in gram-negative bacteria such as *Klebsiella pneumoniae* and *Neisseria gonorrhoeae* are now producing extended-spectrum  $\beta$ -lactamases (ESBLs), which not only hydrolyze penicillin, but also many cephalosporins and monobactams. These variants have been reported worldwide and now pose significant challenges in infection control.

## AT THE BENCH



## Materials

### Cultures

1:1000 brain heart infusion (BHI) broth dilutions of 24-hour BHI broth cultures of *Staphylococcus aureus* ATCC® 27661™ (penicillin-sensitive strain) and *Staphylococcus aureus* ATCC 27659 (penicillinase-producing strain).

### Media

Per designated student group: 40 ml of brain heart infusion broth in a 100-ml Erlenmeyer flask and 10 ml of sterile aqueous crystalline penicillin G solution (100  $\mu$ g/ml).

**TABLE 45.1** Antibiotic Serial Dilution–Broth Medium Setup

Additions (ml) to:	TUBE NUMBER									
	1	2	3	4	5	6	7	8	9	10
Medium	0	2	2	2	2	2	2	2	2	2
Penicillin	2	2	Serial dilution (See protocol)							0
Test culture	2	2	2	2	2	2	2	2	2	2
Total volume	4	4	4	4	4	4	4	4	4*	4
Penicillin ( $\mu$ g/ml)	50	50	25	12.5	6.25	3.12	1.56	0.78	0.39	0
Control	(–)									(+)

\*After 2 ml discarded

## Equipment

Sterile 13 × 100-mm test tubes, test tube racks, sterile 2-ml and 10-ml pipettes, mechanical pipetting device, Bunsen burner, spectrophotometer, glassware marking pencil, and disinfectant solution in a 500-ml beaker.

## Procedure Lab One

1. Into each of two test tube racks, place a set of 10 sterile 13 × 100-mm test tubes labeled 1 through 10. Label one rack Set I—penicillin-sensitive and the other rack Set II—penicillin-resistant. *Refer to Table 45.1 for Steps 2 through 7.*
2. Using a sterile 10-ml pipette and mechanical pipetting device, add 2 ml of BHI broth to the tubes labeled 2 through 10 in Sets I and II. *Note: Discard the pipette into the beaker of disinfectant.*
3. With a 2-ml sterile pipette, add 2 ml of the penicillin solution to Tubes 1 and 2 in Sets I and II. Discard the pipette. *Note: Mix the contents of the tubes well.*
4. **Set I Serial Dilution:** Using a sterile 2-ml pipette, transfer 2 ml from Tube 2 to Tube 3. Mix well and transfer 2 ml from Tube 3 to Tube 4. Continue this procedure through Tube 9 into beaker. Discard 2 ml from Tube 9. Tube 10

receives no antibiotic and serves as a positive control. Discard the pipette. *Note: Remember to mix the contents of each tube well between transfers.*

5. **Set II Serial Dilution:** Using a sterile 2-ml pipette, repeat Step 4.
6. Using a sterile 2-ml pipette, add 2 ml of the 1:1000 dilution of the *S. aureus* ATCC 27661 (penicillin-sensitive strain) to all tubes in Set I. Discard the pipette.
7. Repeat Step 6 to inoculate all the tubes in Set II with the 1:1000 dilution of *S. aureus* ATCC 27659 (penicillinase-producing strain). Discard the pipette.
8. Incubate both sets of tubes for 12 to 18 hours at 37°C.

## Procedure Lab Two

1. Follow the instructions for the use of the spectrophotometer as outlined in Experiment 13 to determine the absorbance readings for Tubes 2 through 10 in Sets I and II. Use the Number 1 tubes, the negative controls, as your blanks to adjust the spectrophotometer.
2. Record your absorbance readings in the chart in the Lab Report.

*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results



### Absorbance Readings at 600 nm

Tube Number	2	3	4	5	6	7	8	9	10
Penicillin concentration ( $\mu\text{g/ml}$ )	50	25	12.5	6.25	3.12	1.56	0.78	0.39	0
Set I									
Set II									

Set I: Minimal inhibitory concentration: \_\_\_\_\_

Set II: Minimal inhibitory concentration: \_\_\_\_\_

## Review Questions

- 
 Was the ability of some microorganisms to produce  $\beta$ -lactamase present prior to their exposure to the antibiotic penicillin? Explain.
  
- 
 Can the results of an MIC test be used to determine whether an antibiotic is bactericidal or bacteriostatic? If not, set up an experimental procedure to determine whether the effect is bactericidal or bacteriostatic.



*This page intentionally left blank*

Antiseptics and disinfectants are chemical substances used to prevent contamination and infection. Many are available commercially for disinfection and asepsis.

Table 46.1 shows the major groups of antimicrobial agents, their modes and ranges of action, and their practical uses.

*(Text continues on page 312.)*

**TABLE 46.1 Chemical Agents—Disinfectants and Antiseptics**

AGENT	MECHANISM OF ACTION	USE
<b>Phenolic Compounds</b> Phenol	<ol style="list-style-type: none"> <li>1. Germicidal effect caused by alteration of protein structure resulting in protein denaturation.</li> <li>2. Surface-active agent (surfactant) precipitates cellular proteins and disrupts cell membranes.</li> </ol> <p>(Phenol has been replaced by better disinfectants that are less irritating, less toxic to tissues, and better inhibitors of microorganisms.)</p>	<ol style="list-style-type: none"> <li>1. 5% solution: Disinfection.</li> <li>2. 0.5% to 1% solutions: Antiseptic effect and relief of itching as it exerts a local anesthetic effect on sensory nerve endings.</li> </ol>
Cresols	<ol style="list-style-type: none"> <li>1. Similar to phenol.</li> <li>2. Poisonous and must be used externally.</li> <li>3. 50% solution of cresols in vegetable oil, known as Lysol®.</li> </ol>	2% to 5% Lysol solutions used as disinfectants.
Hexachlorophene	<p>Germicidal activity similar to phenol.</p> <p>(This agent is to be used with care, especially on infants, because after absorption it may cause neurotoxic effects.)</p>	<ol style="list-style-type: none"> <li>1. Reduction of pathogenic organisms on skin; added to detergents, soaps, lotions, and creams.</li> <li>2. Effective against gram-positive organisms.</li> <li>3. An antiseptic used topically.</li> </ol>
Resorcinol	<ol style="list-style-type: none"> <li>1. Germicidal activity similar to that of phenol.</li> <li>2. Acts by precipitating cell protein.</li> </ol>	<ol style="list-style-type: none"> <li>1. Antiseptic.</li> <li>2. Keratolytic agent for softening or dissolving keratin in epidermis.</li> </ol>
Hexylresorcinol	Germicidal activity similar to that of phenol.	<ol style="list-style-type: none"> <li>1. Treatment of worm infections.</li> <li>2. Urinary antiseptic.</li> </ol>



**TABLE 46.1 Chemical Agents—Disinfectants and Antiseptics (continued)**

AGENT	MECHANISM OF ACTION	USE
Thymol	<ol style="list-style-type: none"> <li>1. Related to the cresols.</li> <li>2. More effective than phenol.</li> </ol>	<ol style="list-style-type: none"> <li>1. Antifungal activity.</li> <li>2. Treatment of hookworm infections.</li> <li>3. Mouthwashes and gargle solutions.</li> </ol>
<b>Alcohols</b> Ethyl: $\text{CH}_3\text{CH}_2\text{OH}$ Isopropyl: $(\text{CH}_3)_2\text{CHOH}$	<ol style="list-style-type: none"> <li>1. Lipid solvent.</li> <li>2. Denaturation and coagulation of proteins.</li> <li>3. Wetting agent used in tinctures to increase the wetting ability of other chemicals.</li> <li>4. Germicidal activity increases with increasing molecular weight.</li> </ol>	Skin antiseptics: Ethyl—50% to 70%. Isopropyl—60% to 70%.
<b>Halogens</b> Chlorine compounds: Sodium hypochlorite (Dakin's fluid): $\text{NaOCl}$ Chloramine: $\text{CH}_3\text{C}_6\text{H}_4\text{SO}_2\text{NNaCl}$	<ol style="list-style-type: none"> <li>1. Germicidal effect resulting from rapid combination with proteins.</li> <li>2. Chlorine reacts with water to form hypochlorous acid, which is bactericidal.</li> <li>3. Oxidizing agent.</li> <li>4. Noncompetitively inhibits enzymes, especially those dealing with glucose metabolism, by reacting with SH and <math>\text{NH}_2</math> groups on the enzyme molecule.</li> </ol>	<ol style="list-style-type: none"> <li>1. Water purification.</li> <li>2. Sanitation of utensils in dairy and restaurant industries.</li> <li>3. Chloramine, 0.1% to 2% solutions, for wound irrigation and dressings.</li> <li>4. Microbicidal.</li> </ol>
Iodine compounds: Tincture of iodine Povidone-iodine solution (Betadine <sup>®</sup> )	<ol style="list-style-type: none"> <li>1. Mechanism of action is not entirely known, but it is believed that it precipitates proteins.</li> <li>2. Surface-active agent.</li> </ol>	<ol style="list-style-type: none"> <li>1. Tinctures of iodine are used for skin antiseptics.</li> <li>2. Treatment of goiter.</li> <li>3. Effective against spores, fungi, and viruses.</li> </ol>
<b>Heavy Metals</b> Mercury compounds: Inorganic: Mercury bichloride Mercurial ointments	<ol style="list-style-type: none"> <li>1. Mercuric ion brings about precipitation of cellular proteins.</li> <li>2. Noncompetitive inhibition of specific enzymes caused by reaction with sulfhydryl group (SH) on enzymes of bacterial cells.</li> </ol>	<ol style="list-style-type: none"> <li>1. Inorganic mercurials are irritating to tissues, toxic systemically, adversely affected by organic matter, and incapable of acting on spores.</li> <li>2. Mercury compounds are mainly used as disinfectants of laboratory materials.</li> </ol>
Organic mercurials: Mercurochrome (merbromin) Merthiolate (thimerosal) Metaphen (nitromersol) Merbak (acetomerocetol)	<ol style="list-style-type: none"> <li>1. Similar to those of inorganic mercurials, but in proper concentrations are useful antiseptics.</li> <li>2. Much less irritating than inorganic mercurials.</li> </ol>	<ol style="list-style-type: none"> <li>1. Less toxic, less irritating; used mainly for skin asepsis.</li> <li>2. Do not kill spores.</li> </ol>
Silver compounds: Silver nitrate	<ol style="list-style-type: none"> <li>1. Precipitate cellular proteins.</li> <li>2. Interfere with metabolic activities of microbial cells.</li> <li>3. Inorganic salts are germicidal.</li> </ol>	Asepsis of mucous membrane of throat and eyes.

**TABLE 46.1 (continued)**

AGENT	MECHANISM OF ACTION	USE
<b>Surface-Active Agents</b> Wetting agents: Emulsifiers, soaps, and detergents	<ol style="list-style-type: none"> <li>1. Lower surface tension and aid in mechanical removal of bacteria and soil.</li> <li>2. If active portion of the agent carries a negative electric charge, it is called an anionic surface-active agent. If active portion of the agent carries a positive electric charge, it is called a cationic surface-active agent.</li> <li>3. Exert bactericidal activity by interfering with or by depressing metabolic activities of microorganisms.</li> <li>4. Disrupt cell membranes.</li> <li>5. Alter cell permeability.</li> </ol>	Weak action against fungi, acid-fast microorganisms, spores, and viruses.
Cationic agents: Quaternary ammonium compounds Benzalkonium chloride	<ol style="list-style-type: none"> <li>1. Lower surface tension because of keratolytic, detergent, and emulsifying properties.</li> <li>2. Their germicidal activities are reduced by soaps.</li> </ol>	<ol style="list-style-type: none"> <li>1. Bactericidal, fungicidal; inactive against spores and viruses.</li> <li>2. Asepsis of intact skin.</li> <li>3. Disinfectant for operating-room equipment.</li> <li>4. Dairy and restaurant sanitization.</li> </ol>
Anionic agents: Tincture of green soap Sodium tetradecyl sulfate	<ol style="list-style-type: none"> <li>1. Neutral or alkaline salts of high-molecular-weight acids. Common soaps included in this group.</li> <li>2. Exert their maximum activity in an acid medium and are most effective against gram-positive cells.</li> <li>3. Same as all surface-active agents.</li> </ol>	<ol style="list-style-type: none"> <li>1. Cleansing agent.</li> <li>2. Sclerosing agent in treatment of varicose veins and internal hemorrhoids.</li> </ol>
Acids ( $H^+$ ) Alkali ( $OH^-$ )	<ol style="list-style-type: none"> <li>1. Destruction of cell wall and cell membrane.</li> <li>2. Coagulation of proteins.</li> </ol>	Disinfection; however, of little practical value.
<b>Formaldehyde</b> (liquid or gas)	Alkylating agent causes reduction of enzymes.	<ol style="list-style-type: none"> <li>1. Room disinfection.</li> <li>2. Alcoholic solution for instrument disinfection.</li> <li>3. Specimen preservation.</li> </ol>
<b>Ethylene Oxide</b>	Alkylating agent causes reduction of enzymes.	Sterilization of heat-labile material.
<b><math>\beta</math>-Propiolactone</b> (liquid or gas)	Alkylating agent causes reduction of enzymes.	<ol style="list-style-type: none"> <li>1. Sterilization of tissue for grafting.</li> <li>2. Destruction of hepatitis virus.</li> <li>3. Room disinfection.</li> </ol>
<b>Basic Dyes</b> Crystal violet	Affinity for nucleic acids; interfere with reproduction in gram-positive organisms.	<ol style="list-style-type: none"> <li>1. Skin antiseptic.</li> <li>2. Laboratory isolation of gram-negative bacteria.</li> </ol>

The efficiency of all disinfectants and antiseptics is influenced by a variety of factors, including the following:

- 1. Concentration:** The concentration of a chemical substance markedly influences its effect on microorganisms, with higher concentrations producing a more rapid death. Concentration cannot be arbitrarily determined; the toxicity of the chemical to the tissues being treated and the damaging effect on nonliving materials must also be considered.
- 2. Length of exposure:** All microbes are not destroyed within the same exposure time. Sensitive forms are destroyed more rapidly than resistant ones. The longer the exposure to the agent, the greater its antimicrobial activity. The toxicity of the chemical and environmental conditions must be considered in assessing the length of time necessary for disinfection or asepsis.
- 3. Type of microbial population to be destroyed:** Microorganisms vary in their susceptibility to destruction by chemicals. Bacterial spores are the most resistant forms. Capsulated bacteria are more resistant than noncapsulated forms; acid-fast bacteria are more resistant than non-acid-fast; and older, metabolically less-active cells are more resistant than younger cells. Awareness of the types of microorganisms that may be present will influence the choice of agent.
- 4. Environmental conditions:** Conditions under which a disinfectant or antiseptic affects the chemical agent are as follows:
  - a. Temperature:** Cells are killed as the result of a chemical reaction between the agent and cellular component. As increasing temperatures increase the rate of chemical reactions, application of heat during disinfection markedly increases the rate at which the microbial population is destroyed.
  - b. pH:** The pH conditions during disinfection may affect not only the microorganisms but also the compound. Extremes in pH are harmful to many microorganisms and may enhance the antimicrobial action of a chemical. Deviation from a neutral pH may cause ionization of the disinfectant; depending on the chemical agent, this may serve to increase or decrease the chemical's microbicidal action.
  - c. Type of material on which the microorganisms exist:** The destructive power of the compound on cells is due to its

combination with organic cellular molecules. If the material on which the microorganisms are found is primarily organic, such as blood, pus, or tissue fluids, the agent will combine with these extracellular organic molecules, and its antimicrobial activity will be reduced.

Numerous laboratory procedures are available for evaluating the antimicrobial efficiency of disinfectants or antiseptics. They provide a general rather than an absolute measure of the effectiveness of any agent because test conditions frequently differ considerably from those seen during practical use. The agar plate-sensitivity method, a commonly employed procedure, is presented.

## LEARNING OBJECTIVE

*Once you have completed this experiment, you should be able to*

1. Evaluate the effectiveness of antiseptic agents against selected test organisms.


## Principle

This procedure requires the heavy inoculation of an agar plate with the test organism. Sterile, color-coded filter-paper discs are impregnated with a different antiseptic and equally spaced on the inoculated agar plate. Following incubation, the agar plate is examined for zones of inhibition (areas of no microbial growth) surrounding the discs. A zone of inhibition is indicative of microbicidal activity against the organism. Absence of a zone of inhibition indicates that the chemical was ineffective against the test organism. *Note: The size of the zone of inhibition is not indicative of the degree of effectiveness of the chemical agent.* Antiseptic susceptibility is represented in **Figure 46.1**.

## CLINICAL APPLICATION

### MRSA and Disinfection

Methicillin-resistant *Staphylococcus aureus* (MRSA) has achieved notoriety for causing infections that are difficult to treat with conventional antimicrobials, but these strains have also demonstrated resistance to disinfection. One study has shown that resistance to methicillin is directly related to lack of susceptibility to benzalkonium chloride and other disinfectants. It may be that adjusted contact times are necessary to adequately kill these troublesome strains.

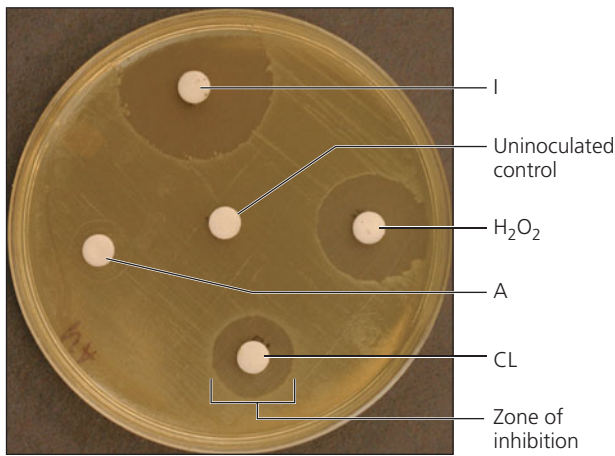

**Figure 46.1 Antiseptic susceptibility test.** Discs are saturated with chlorine bleach (CL), hydrogen peroxide ( $H_2O_2$ ), isopropyl alcohol (A), and tincture of iodine (I).

## AT THE BENCH



## Materials

### Cultures

24- to 48-hour Trypticase soy broth cultures of *Escherichia coli*, *Bacillus cereus*, *Staphylococcus aureus*, and *Mycobacterium smegmatis*, and a 7-day-old Trypticase soy broth culture of *Bacillus cereus*.

### Media

Per designated student group: five Trypticase soy agar plates.

### Antiseptics/Disinfectants

Per designated student group: 10 ml of each of the following dispensed in 25-ml beakers: tincture of iodine, 3% hydrogen peroxide, 70% isopropyl alcohol, and 5% chlorine bleach.

### Equipment

Four different-colored, sterile Sensi-Discs; forceps; sterile cotton swabs; Bunsen burner; and glass-ware marking pencil.

## Procedure Lab One

1. Aseptically inoculate the appropriately labeled agar plates with their respective test organisms by streaking each plate in horizontal and vertical directions and around the edge with a sterile swab.
2. Color-code the Sensi-Discs according to the chemical agents to be used (e.g., red = chlorine bleach).
3. Using forceps dipped in alcohol and flamed, expose five discs of the same color by placing them into the solution of one of the chemical agents. Drain the saturated discs on absorbent paper immediately prior to placing one on each of the inoculated agar plates. Place each disc approximately 2 cm in from the edge of the plate. Gently press the discs down with the forceps so that they adhere to the surface of the agar.
4. Impregnate the remaining discs as described in Step 3. Place one of each of the three remaining colored discs on the surface of each of the five inoculated agar plates equidistant from each other around the periphery of the plate.
5. Incubate all plate cultures in an inverted position for 24 to 48 hours at  $37^\circ\text{C}$ .

## Procedure Lab Two

1. Observe all the plates for the presence of a zone of inhibition surrounding each of the impregnated discs.
2. Record your observations in the chart provided in the Lab Report.

*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

1. Indicate the absence of a zone of inhibition as (0), and the presence of a zone of inhibition as (+).

Bacterial Species	ANTIMICROBIAL AGENT			
	Tincture of Iodine	3% Hydrogen Peroxide	70% Isopropyl Alcohol	5% Chlorine Bleach
<i>E. coli</i> gram-negative				
<i>S. aureus</i> gram-positive				
<i>M. smegmatis</i> acid-fast				
<i>B. cereus</i> spore-former gram-positive				
<i>B. cereus</i> spore-former gram-positive 7-day-old				

2. Indicate which of the antiseptics exhibited microbicidal activity against each of the following groups of microorganisms.

Bacterial Group	Tincture of Iodine	3% Hydrogen Peroxide	70% Isopropyl Alcohol	5% Chlorine Bleach
Gram-negative				
Gram-positive				
Acid-fast				
Spore-former				

3. Which of the experimental chemical compounds appears to have the broadest range of microbicidal activity? The narrowest range of microbicidal activity?



## Review Questions

1. Can the disinfection period (exposure time) be arbitrarily increased? Explain.

2.  A household cleanser is labeled germicidal. Explain what this means to you.

# Microbiology of Food

## LEARNING OBJECTIVES

Once you have completed the experiments in this section, you should be familiar with

1. The endogenous and exogenous organisms that may be found in food products.
2. The analysis of food products as a means of determining their quality from the public health point of view.
3. The microbiological production of wine.

## Introduction

Microbiologists have always been aware that foods, especially milk, have served as important inanimate vectors in the transmission of disease. Foods contain the organic nutrients that provide an excellent medium to support the growth and multiplication of microorganisms under suitable temperatures.

Food and dairy products may be contaminated in a variety of ways and from a variety of sources:

1. **Soil and water:** Food-borne organisms that may be found in soil and water and that may contaminate food are members of the genera *Alcaligenes*, *Bacillus*, *Citrobacter*, *Clostridium*, *Pseudomonas*, *Serratia*, *Proteus*, *Enterobacter*, and *Micrococcus*. The common soil and water molds include *Rhizopus*, *Penicillium*, *Botrytis*, *Fusarium*, and *Trichothecium*.
2. **Food utensils:** The type of microorganisms found on utensils depends on the type of food and the manner in which the utensils were handled.
3. **Enteric microorganisms of humans and animals:** The major members of this group are *Bacteroides*, *Lactobacillus*, *Clostridium*,

*Escherichia*, *Salmonella*, *Proteus*, *Shigella*, *Staphylococcus*, and *Streptococcus*. These organisms find their way into the soil and water, from which they contaminate plants and are carried by wind currents onto utensils or prepared and exposed foods.

4. **Food handlers:** People who handle foods are especially likely to contaminate them because microorganisms on hands and clothing are easily transmitted. A major offending organism is *Staphylococcus*, which is generally found on hands and skin, and in the upper respiratory tract. Food handlers with poor personal hygiene and unsanitary habits are most likely to contaminate foods with enteric organisms.
5. **Animal hides and feeds:** Microorganisms found in water, soil, feed, dust, and fecal debris can be found on animal hides. Infected hides may serve as a source of infection for workers, or the microorganisms may migrate into the musculature of the animal and remain viable following its slaughter.

By enumerating microorganisms in milk and foods, the quality of a particular sample can be determined. Although the microorganisms cannot be identified, the presence of a high number suggests a good possibility that pathogens may be present. Even if a sample contains a low microbial count, it can still transmit infection.

In the laboratory procedures that follow, you will have an opportunity directly and indirectly to enumerate the number of microorganisms present in milk and other food products and to thereby determine the quality of the samples.

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be able to determine

1. The total number of microorganisms present in food products.
2. The number of coliform bacteria in the selected food products.

## Principle

Microorganisms in food may be harmful in some cases, while in others they are beneficial. Certain microorganisms are necessary in preparation of foods such as cheese, pickles, yogurt, and sausage. Other microorganisms, however, are responsible for serious and sometimes fatal food poisoning and spoilage.

## CLINICAL APPLICATION

### Microorganisms in Food

Some of the pathogens that are tested for in food include: *Escherichia coli*, *Listeria monocytogenes*, *Salmonella* species, and fungi of the genus *Aspergillus*.

## AT THE BENCH



## Materials

### Cultures

Samples of fresh vegetables, ground beef, and dried fruit.

### Media

Per designated student group: nine brain heart infusion agar deep tubes, three eosin–methylene blue (EMB) agar plates, three 99-ml sterile water blanks, and three 180-ml sterile water blanks.

## Equipment

Bunsen burner, waterbath, Quebec or electronic colony counter, balance, sterile glassine weighing paper, blender with three sterile jars, sterile Petri dishes, 1-ml pipettes, mechanical pipetting device, inoculation loop, and glassware marking pencil.

## Procedure Lab One


Figure 47.1 illustrates the procedure.

1. Label three sets of three Petri dishes for each of the food samples to be tested and their dilutions ( $10^{-2}$ ,  $10^{-3}$ ,  $10^{-4}$ ). Label the three EMB agar plates with the names of the food.
2. Melt the brain heart infusion agar deep tubes in a waterbath, cool, and maintain at  $45^{\circ}\text{C}$ .
3. Place 20 g of each food sample, weighed on sterile glassine paper, into its labeled blender jar. Add 180 ml of sterile water to each of the blender jars and blend each mixture for 5 minutes. You will have made a 1:10 ( $10^{-1}$ ) dilution of each food sample.
4. Transfer 1 ml of the  $10^{-1}$  ground beef suspension into its labeled 99-ml sterile water blank, thereby effecting a  $10^{-3}$  dilution, and 0.1 ml to the appropriately labeled  $10^{-2}$  Petri dish. Shake the  $10^{-3}$  sample dilution, and using a different pipette, transfer 1 ml to the plate labeled  $10^{-3}$  and 0.1 ml to the plate labeled  $10^{-4}$ . Add a 15-ml aliquot of the molten and cooled agar to each of the three plates. Swirl the plates gently to obtain a uniform distribution, and allow the plates to solidify.
5. Repeat Step 4 for the remaining two  $10^{-1}$  test food sample dilutions.
6. Aseptically prepare a four-way streak plate, as described in Experiment 3, and inoculate each  $10^{-1}$  food sample dilution on its appropriately labeled EMB agar plate.
7. Incubate all plates in an inverted position for 24 to 48 hours at  $37^{\circ}\text{C}$ .

## Procedure Lab Two

1. Following the instructions in the Lab Report, count and record the number of colonies on each plate.

## PROCEDURE

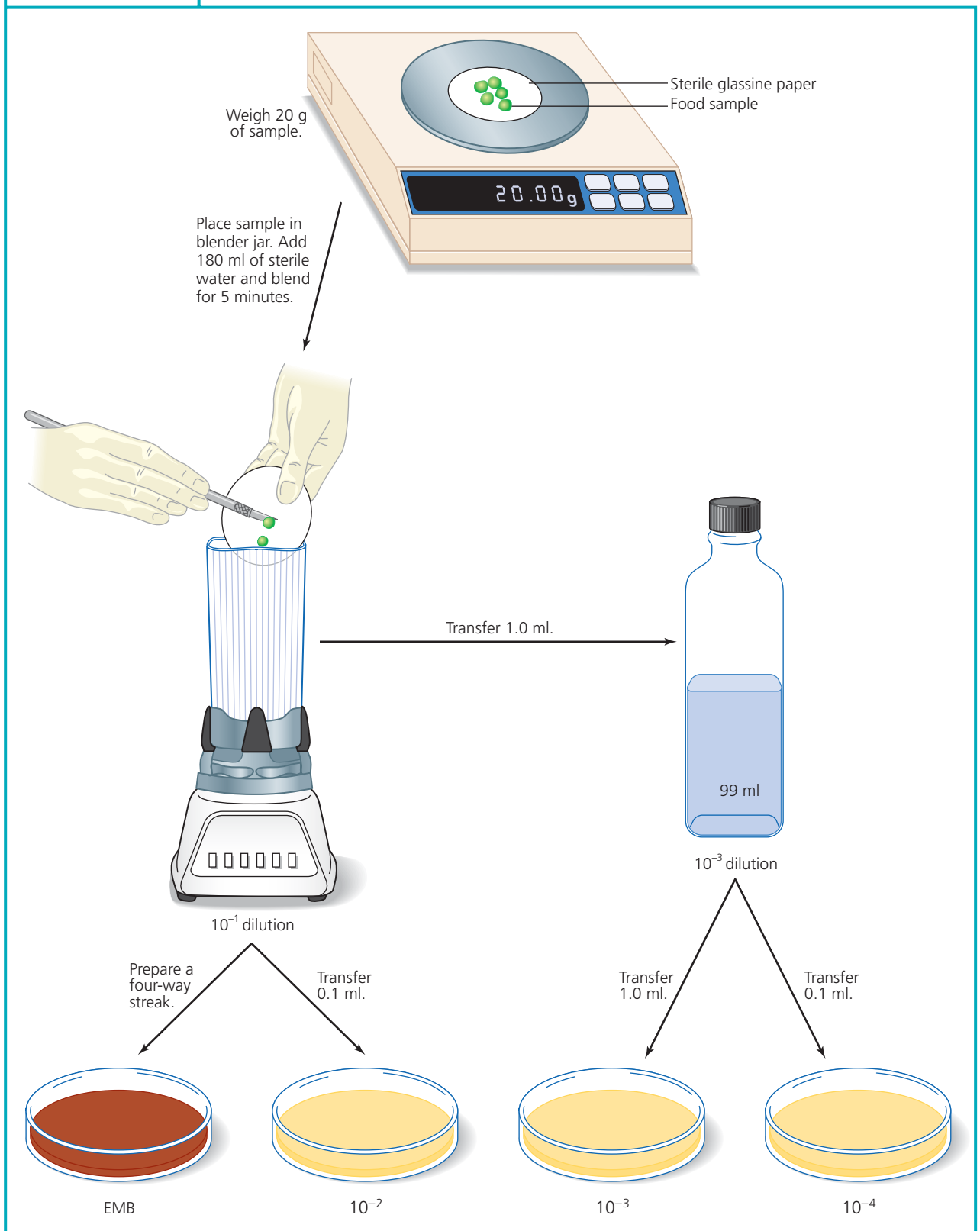

**Figure 47.1** Preparation of a food sample for analysis

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

- Using either the Quebec or electronic colony counter, count the number of colonies on each plate. Count only statistically valid plates that contain between 30 and 300 colonies. Designate plates with fewer than 30 colonies as **too few to count (TFTC)** and plates with more than 300 colonies as **too numerous to count (TNTC)**.
- Determine the number of organisms per milliliter of each food sample on plates not designated as TFTC or TNTC by multiplying the number of colonies counted by the dilution factor.
- Record in the chart below the number of colonies per plate and the number of organisms per milliliter of each food sample.

Type of Food	Dilution	Number of Colonies per Plate	Number of Organisms per ml
Ground beef	$10^{-2}$		
	$10^{-3}$		
	$10^{-4}$		
Fresh vegetables	$10^{-2}$		
	$10^{-3}$		
	$10^{-4}$		
Dried fruits	$10^{-2}$		
	$10^{-3}$		
	$10^{-4}$		


- Examine the eosin–methylene blue agar plate cultures for colonies with a metallic green sheen on their surfaces, which is indicative of *E. coli*. Indicate in the chart below the presence or absence of *E. coli* growth and the possibility of fecal contamination of the food.

Sample	<i>E. coli</i> (+) or (-)	Fecal Contamination (+) or (-)
Ground beef		
Fresh vegetables		
Dried fruit		

## Review Questions

1. Indicate some possible ways in which foods may become contaminated with enteric organisms.

2.  Explain why it is not advisable to thaw and then refreeze food products without having cooked them.

3.  Following a Fourth of July picnic lunch of ham, sour pickles, potato salad, and cream puffs, a group of students was admitted to the hospital with severe gastrointestinal distress. A diagnosis of staphylococcal food poisoning was made. Explain how the staphylococci can multiply in these foods and produce severe abdominal distress.

## LEARNING OBJECTIVES

Once you have completed this experiment, you should understand

1. Wine production by the fermentative activities of yeast cells.


## Principle

Wine is a product of the natural fermentation of the juices of grapes and other fruits, such as peaches, pears, plums, and apples, by the action of yeast cells. This biochemical conversion of juice to wine occurs when the yeast cells enzymatically degrade the fruit sugars fructose and glucose first to acetaldehyde and then to alcohol, as illustrated in **Figure 48.1**.

Grapes containing 20% to 30% sugar concentration will yield wines with an alcohol content of approximately 10% to 15%. Also present in grapes are acids and minerals whose concentrations

are increased in the finished product and that are responsible for the characteristic tastes and bouquets of different wines. For red wine, the crushed grapes must be fermented with their skins to allow extraction of their color into the juice. White wine is produced from the juice of white grapes.

The commercial production of wine is a long and exacting process. First, the grapes are crushed or pressed to express the juice, which is called **must**. Potassium metabisulfite is added to the must to retard the growth of acetic acid bacteria, molds, and wild yeast that are endogenous to grapes in the vineyard. A wine-producing strain of yeast, *Saccharomyces cerevisiae* var. *ellipsoideus*, is used to inoculate the must, which is then incubated for 3 to 5 days under aerobic conditions at 21°C to 32°C. This is followed by an anaerobic incubation period. The wine is then aged for a period of 1 to 5 years in aging tanks or wooden barrels. During this time, the wine is clarified of any turbidity, thereby producing volatile esters that are responsible for characteristic flavors. The clarified product is then filtered, pasteurized at 60°C for 30 minutes, and bottled.

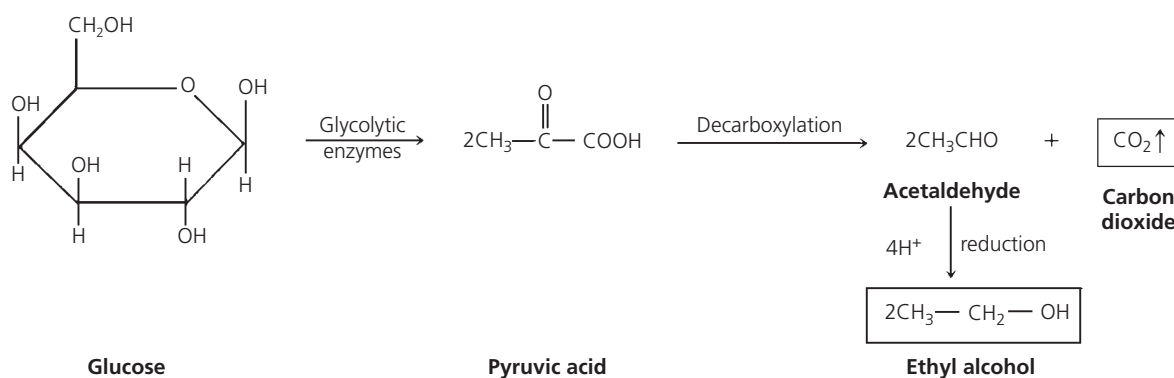

**Figure 48.1** Biochemical pathway for alcohol production



This experiment is a modified method in which white wine is produced from white grape juice. You will examine the fermenting wine at 1-week intervals during the incubation period for:

1. Total acidity (expressed as % tartaric acid): To a 10-ml aliquot of the fermenting wine, add 10 ml of distilled water and 5 drops of 1% phenolphthalein solution. Mix and titrate to the first persistent pink color with 0.1N sodium hydroxide. Calculate total acidity using the following formula:

$$\% \text{ tartaric acid} = \frac{\text{ml alkali} \times \text{normality of alkali} \times 7.5}{\text{weight of sample in g}^*}$$

\*1 ml = 1 g

2. Volatile acidity (expressed as % acetic acid): Following titration, calculate volatile acidity using the following formula:

$$\% \text{ acetic acid} = \frac{\text{ml alkali} \times \text{normality of alkali} \times 6.0}{\text{weight of sample in g}^*}$$

\*1 ml = 1 g

3. Alcohol (expressed as volume %): Optional; can be determined by means of an ebulliometer.
4. Aroma: Fruity, yeastlike, sweet, none.
5. Clarity: Clear, turbid.

## CLINICAL APPLICATION

### Drinking Wine instead of Water for Better Health

For thousands of years mankind has allowed crushed fruits and boiled grains to ferment, creating wine. Wild yeasts and bacteria metabolize and break down the inherent sugars in these liquids, and the fermentation byproduct of alcohol kills all bacteria and protozoa present. Early civilizations drank wine instead of water to protect against diseases. Poorer subjects and young children would drink watered down wine. By replacing water with wine in their daily diet, early civilizations were able to limit their exposure to pathogens.

## AT THE BENCH



## Materials

### Cultures

50 ml of white grape juice broth culture of *Saccharomyces cerevisiae* var. *ellipsoideus* incubated for 48 hours at 25°C.

### Media

Per designated student group: 500 ml of pasteurized Welch's® commercial white grape juice.

### Reagents

1% phenolphthalein solution, 0.1N sodium hydroxide, and sucrose.

### Equipment

1-liter Erlenmeyer flask, one-holed rubber stopper containing a 2-inch glass tube plugged with cotton, pan balance, spatula, glassine paper, 10-ml graduated cylinder, ebulliometer (optional), and burette or pipette for titration.


## Procedure

1. Pour 500 ml of the white grape juice into the 1-liter Erlenmeyer flask. Add 20 g of sucrose and the 50 ml of *S. cerevisiae* grape juice broth culture (10% starter culture). Close the flask with the stopper containing a cotton plugged air vent.
2. After 2 days and 4 days of incubation, add 20 g of sucrose to the fermenting wine.
3. Incubate the fermenting wine for 21 days at 25°C.
4. Using uninoculated white grape juice:
  - a. Perform a titration to determine total acidity and volatile acidity.
  - b. Note aroma and clarity.
  - c. Determine volume % alcohol (optional).
5. Record your results in the chart in the Lab Report.
6. At 7-day intervals, using samples of the fermenting wine, repeat Steps 4a through 4c and record your results in the Lab Report.



3. What are the chemical end products of fermentation?

4. How are white and red wines produced?

5.  Why is wine pasteurized? Would it be preferable to sterilize the wine? Explain.

# Microbiology of Water

## LEARNING OBJECTIVES

Once you have completed the experiments in this section, you should be familiar with

1. The types of microorganisms present in water.
2. The methods to determine the potability of water using standard qualitative and quantitative procedures.

## Introduction

The importance of potable (drinking) water supplies cannot be overemphasized. With increasing industrialization, water sources available for consumption and recreation have been adulterated with industrial as well as animal and human wastes. As a result, water has become a formidable factor in disease transmission. Polluted waters contain vast amounts of organic matter that serve as excellent nutritional sources for the growth and multiplication of microorganisms. The presence of nonpathogenic organisms is not of major concern, but intestinal contaminants of fecal origin are important. These pathogens are responsible for intestinal infections such as **bacillary dysentery, typhoid fever, cholera, and paratyphoid fever.**

The World Health Organization (WHO) estimates that 1.7 million deaths per year result from unsafe water supplies. Most of these are from diarrheal diseases, and 90% of these deaths occur in children living in developing countries where sanitary facilities and potable water are at a minimum. The WHO indicates that about 3.4 million deaths annually are caused by dangerous waterborne enteric bacterial pathogens such as *Shigella dysenteriae*, *Campylobacter jejuni*, *Salmonella typhi*, and *Vibrio cholerae*.

In addition to bacterial infections, unsafe water supplies are responsible for numerous parasitological infections, including helminth diseases such as schistosomiasis and especially guinea worm

(*Dracunculus medinensis*), which infects about 200 million people worldwide each year. Intestinal, hepatic, and pulmonary flukes such as *Fasciolopsis buski*, *Clonorchis sinensis*, and *Paragonimus westermani* are responsible for human infection and are all associated with unsafe water and sanitation. The parasitic protozoa *Entamoeba histolytica*, *Giardia intestinalis* (formerly called *G. lamblia*), and *Balantidium coli* are just a few of the protozoa responsible for major diarrheal disease in humans.

Although waterborne infections occur in the United States, their incidence in comparison to the rest of the world is much lower, and they occur sporadically. This can be attributed to the diligent attention given to our water supplies and sewage disposal systems.

Analysis of water samples on a routine basis would not be possible if each pathogen required detection. Therefore, water is examined to detect *Escherichia coli*, the bacterium that indicates fecal pollution. Since *E. coli* is always present in the human intestine, its presence in water alerts public health officials to the possible presence of other human or animal intestinal pathogens. However, in the tropics and subtropics it is not considered a reliable indicator of fecal pollution because the soil in these regions naturally contains high levels of *E. coli*. Therefore, *E. coli* is present in the water anytime there is surface runoff. Both qualitative and quantitative methods are used to determine the sanitary condition of water.

*This page intentionally left blank*

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be able to

1. Determine the presence of coliform bacteria in a water sample.
2. Obtain an index indicating the possible number of organisms present in the sample under analysis.
3. Confirm the presence of coliform bacteria in a water sample for which the presumptive test was positive.
4. Confirm the presence of coliform bacteria in a water sample, or if necessary, confirm a suspicious or doubtful result from the previous test.

## Principle

The three basic tests to detect coliform bacteria in water are presumptive, confirmed, and completed (Figure 49.1). The tests are performed sequentially on each sample under analysis. They detect the presence of coliform bacteria (indicators of fecal contamination), the gram-negative, non-spore-forming bacilli that ferment lactose with the production of acid and gas that is detectable following a 24-hour incubation period at 37°C.

## The Presumptive Test

The **presumptive test** is specific for detection of coliform bacteria. Measured aliquots of the water to be tested are added to a lactose fermentation broth containing an inverted gas vial. Because these bacteria are capable of using lactose as a carbon source (the other enteric organisms are not), their detection is facilitated by the use of this medium. In this experiment, you will use lactose fermentation broth containing an inverted Durham tube for gas collection.


Tubes of this lactose medium are inoculated with 10-ml, 1-ml, and 0.1-ml aliquots of the water sample. The series consists of at least three groups, each composed of five tubes of the specified medium. The tubes in each group are then inoculated with the designated volume of the water sample, as described under “Procedure: Lab One.” The greater the number of tubes per group, the greater the sensitivity of the test. Development of gas in any of the tubes is *presumptive* evidence of the presence of coliform bacteria in the sample. The presumptive test also enables the microbiologist to obtain some idea of the number of coliform organisms present by means of the **most probable number (MPN) test**. The MPN is estimated by determining the number of tubes in each group that show gas following the incubation period (Table 49.1 on page 334).

## The Confirmed Test

The presence of a positive or doubtful presumptive test immediately suggests that the water sample is nonpotable. Confirmation of these results is necessary because positive presumptive tests may be the result of organisms of noncoliform origin that are not recognized as indicators of fecal pollution.

The **confirmed test** requires that selective and differential media such as eosin-methylene blue (EMB) or Endo agar be streaked from a positive lactose broth tube obtained from the presumptive test. The nature of the differential and selective media was discussed in Experiment 15 but is reviewed briefly here. Eosin-methylene blue contains the dye methylene blue, which inhibits the growth of gram-positive organisms. In the presence of an acid environment, EMB forms a complex that precipitates out onto the coliform colonies, producing dark centers and a green metallic sheen. The reaction is characteristic for *Escherichia coli*, the major indicator of fecal pollution. Endo agar is a nutrient medium containing the dye fuchsin, which is present in the decolorized state. In the presence of acid produced by the coliform bacteria, fuchsin forms a dark pink complex that turns the *E. coli* colonies and the surrounding medium pink.

# PROCEDURE

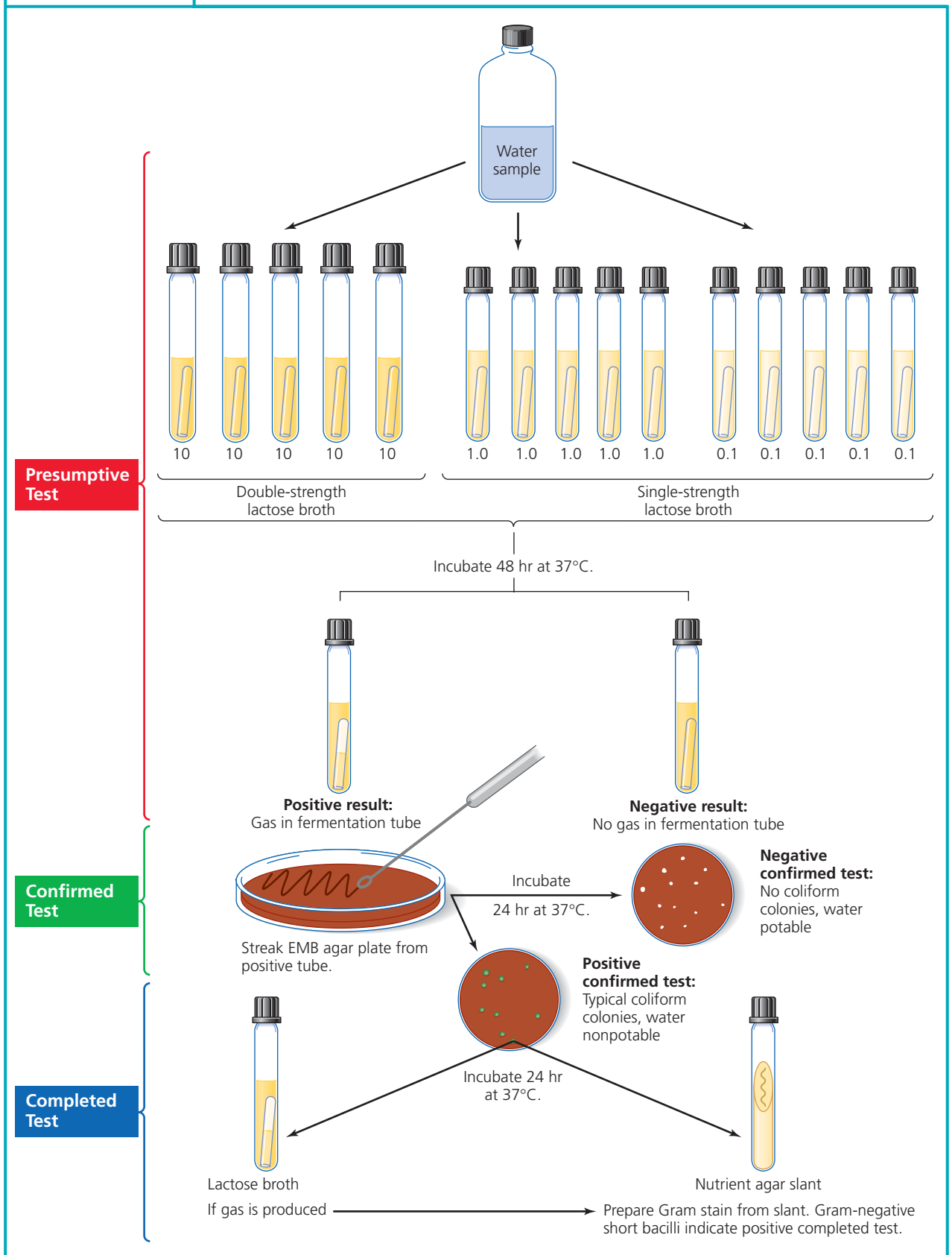

**Figure 49.1** Standard method for bacteriological water analysis

## The Completed Test

The **completed test** is the final analysis of the water sample. It is used to examine the coliform colonies that appeared on the EMB or Endo agar plates used in the confirmed test. An isolated colony is picked up from the confirmatory test plate and inoculated into a tube of lactose broth and streaked on a nutrient agar slant to perform a Gram stain. Following inoculation and incubation, tubes showing acid and gas in the lactose broth and presence of gram-negative bacilli on microscopic examination are further confirmation of the presence of *E. coli*, and they are indicative of a positive completed test.

### CLINICAL APPLICATION

#### Testing for Safe Water

The analysis of water used for human consumption and recreational use is routinely performed for safety. Water sources are regularly tested for the presence of *Escherichia coli* to determine the quality and safety of municipal water supplies. Several testing methods are available for this purpose including: most probable numbers (MPN), ATP testing, membrane filtration, and the use of pour plates.

### AT THE BENCH



## Materials

### Cultures

Lab One: Water samples from sewage plant, pond, and tap. Lab Two: One 24-hour-old positive lactose broth culture from each of the three series from the presumptive test. Lab Three: One 24-hour coliform-positive EMB or Endo agar culture from each of the three series of the confirmed test.

### Media

Lab One (per designated student group): 15 double-strength lactose fermentation broths (LB2X) and 30 single-strength lactose fermentation broths (LB1X). Lab Two (three each per designated student group): eosin–methylene blue agar plates or Endo agar plates. Lab Three (three each per designated student group): nutrient agar slants and lactose fermentation broths.

## Reagents

Lab Three: Crystal violet, Gram's iodine, 95% ethyl alcohol, and safranin.

## Equipment

Lab One: Bunsen burner, 45 test tubes, test tube rack, sterile 10-ml pipettes, sterile 1-ml pipettes, sterile 0.1-ml pipettes, mechanical pipetting device, and glassware marking pencil. Lab Two: Bunsen burner, glassware marking pencil, and inoculating loop. Lab Three: Bunsen burner, staining tray, inoculating loop, lens paper, bibulous paper, microscope, and glassware marking pencil.

## Procedure Lab One

### Presumptive Test

 **Exercise care in handling sewage waste water sample because enteric pathogens may be present.**

1. Set up three separate series consisting of three groups, a total of 15 tubes per series, in a test tube rack; for each tube, label the water source and volume of sample inoculated as illustrated.

Series 1: Sewage water	5 tubes of LB2X-10 ml 5 tubes of LB1X-1 ml 5 tubes of LB1X-0.1 ml
Series 2: Pond water	5 tubes of LB2X-10 ml 5 tubes of LB1X-1 ml 5 tubes of LB1X-0.1 ml
Series 3: Tap water	5 tubes of LB2X-10 ml 5 tubes of LB1X-1 ml 5 tubes of LB1X-0.1 ml

2. Mix sewage plant water sample by shaking thoroughly.
3. Flame bottle and then, using a 10-ml pipette, transfer 10-ml aliquots of water sample to the five tubes labeled LB2X-10 ml.
4. Flame bottle and then, using a 1-ml pipette, transfer 1-ml aliquots of water sample to the five tubes labeled LB1X-1 ml.
5. Flame bottle and then, using a 0.1-ml pipette, transfer 0.1-ml aliquots of water sample to the five tubes labeled LB1X-0.1 ml.





(a) (b) (c) (d) (e)

**Figure 49.2 Possible MPN presumptive test results.** (a) Uninoculated control tube, (b, c) inoculated tubes with no change, (d) inoculated tube with acid production only, and (e) inoculated tube with acid and gas production—the only positive result of the five tubes.

- Repeat Steps 2 through 5 for the tap and pond water samples.
- Incubate all tubes for 48 hours at 37°C.

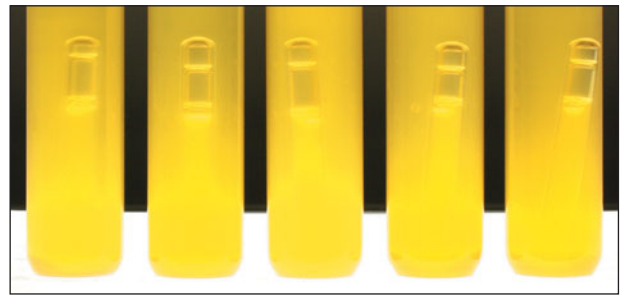
## Procedure Lab Two

### Presumptive Test


- Examine the tubes from your presumptive test after 24 and 48 hours of incubation. Your results are positive if the Durham tube fills 10% or more with gas in 24 hours, doubtful if gas develops in the tube after 48 hours, and negative if there is no gas in the tube after 48 hours. Refer to **Figure 49.2** for a summary of possible MPN presumptive test results. **Figure 49.3** shows actual results from an MPN presumptive test for a water sample. Record your results in the Lab Report.
- Determine the MPN using Table 49.1, and record your results in the Lab Report.

### Confirmed Test


- Label the covers of the three EMB plates or the three Endo agar plates with the source of the water sample (sewage, pond, and tap).
- Using a positive 24-hour lactose broth culture from the sewage water series from the presumptive test, streak the surface of one EMB or one Endo agar plate, as described in Experiment 3, to obtain discrete colonies.
- Repeat Step 2 using the positive lactose broth cultures from the pond and tap water series from the presumptive test to inoculate the remaining plates.



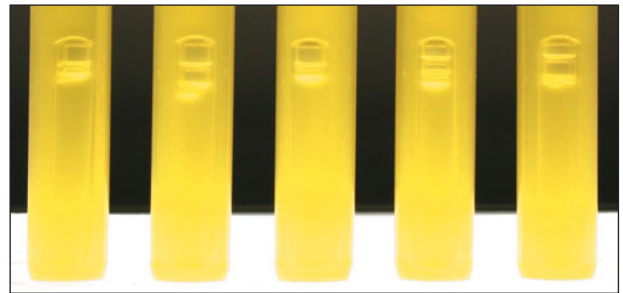(a) Positive results (acid and gas) in five 10-ml tubes



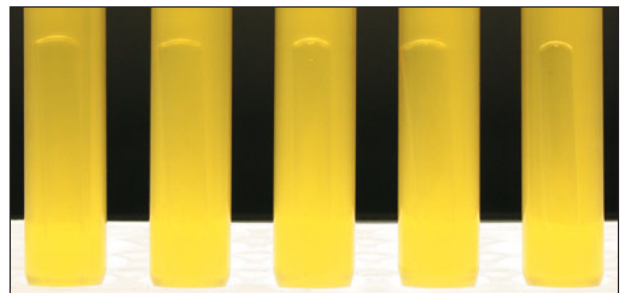(b) Positive results (acid and gas) in five 1-ml tubes



(c) Negative results (acid only) in five 0.1-ml tubes

**Figure 49.3 MPN presumptive test results for a water sample.** The results of this test (5 positive, 5 positive, and 5 negative) indicate 240 coliforms per 100 ml of water (see Table 49.1). This represents a positive presumptive test for the presence of coliforms in the tested water sample.

- Incubate all plate cultures in an inverted position for 24 hours at 37°C.

## Procedure Lab Three

### Confirmed Test

- Examine all the plates from your confirmed test for the presence or absence of *E. coli* colonies (refer to the description of the confirmed test in the experiment introduction, and see Figure 15.2 for an illustration of *E. coli* growth on EMB agar). Record your results in the Lab Report.

2. Based on your results, determine whether each of the samples is potable or nonpotable. The presence of *E. coli* is a positive confirmed test, indicating that the water is nonpotable. The absence of *E. coli* is a negative test, indicating that the water is not contaminated with fecal wastes and is therefore potable. Record your results in the Lab Report.

### Completed Test

1. Label each tube of nutrient agar slants and lactose fermentation broths with the source of its water sample.
2. Inoculate one lactose broth and one nutrient agar slant with a positive isolated *E. coli* colony obtained from each of the experimental water samples during the confirmed test.
3. Incubate all tubes for 24 hours at 37°C.

## Procedure Lab Four

### Completed Test

1. Examine all lactose fermentation broth cultures for the presence or absence of acid and gas. Record your results in the Lab Report.
2. Prepare a Gram stain, using the nutrient agar slant cultures of the organisms that showed a positive result in the lactose fermentation broth (refer to Experiment 11 for the staining procedure).
3. Examine the slides microscopically for the presence of gram-negative short bacilli, which are indicative of *E. coli* and thus nonpotable water. In the Lab Report, record your results for Gram stain reaction and morphology of the cells.

**TABLE 49.1 The MPN Index per 100 ml for Combinations of Positive and Negative Presumptive Test Results When Five 10-ml, Five 1-ml, and Five 0.1-ml Portions of Sample Are Used**

NUMBER OF TUBES WITH POSITIVE RESULTS						NUMBER OF TUBES WITH POSITIVE RESULTS					
FIVE OF 10 ML EACH	FIVE OF 1 ML EACH	FIVE OF 0.1 ML EACH	MPN INDEX PER 100 ML	95% CONFIDENCE LIMITS		FIVE OF 10 ML EACH	FIVE OF 1 ML EACH	FIVE OF 0.1 ML EACH	MPN INDEX PER 100 ML	95% CONFIDENCE LIMITS	
				LOWER	UPPER					LOWER	UPPER
0	0	0	<2	0	6	4	2	1	26	7	67
0	0	1	2	<0.5	7	4	3	0	27	9	78
0	1	0	2	<0.5	7	4	3	1	33	9	78
0	2	0	4	<0.5	11	4	4	0	34	11	93
1	0	0	2	0.1	10	5	0	0	23	7	70
1	0	1	4	0.7	10	5	0	1	31	11	89
1	1	0	4	0.7	12	5	0	2	43	14	100
1	1	1	6	1.8	15	5	1	0	33	10	100
1	2	0	6	1.8	15	5	1	1	46	14	120
2	0	0	5	<0.5	13	5	1	2	63	22	150
2	0	1	7	1	17	5	2	0	49	15	150
2	1	0	7	1	17	5	2	1	70	22	170
2	1	1	9	2	21	5	2	2	94	34	230
2	2	0	9	2	21	5	3	0	79	22	220
2	3	0	12	3	28	5	3	1	110	34	250
3	0	0	8	2	22	5	3	2	140	52	400
3	0	1	11	4	23	5	3	3	180	70	400
3	1	0	11	5	35	5	4	0	130	36	400
3	1	1	14	6	36	5	4	1	170	58	400
3	2	0	14	6	36	5	4	2	220	70	440
3	2	1	17	7	40	5	4	3	280	100	710
3	3	0	17	7	40	5	4	4	350	100	710
4	0	0	13	4	35	5	5	0	240	70	710
4	0	1	17	6	36	5	5	1	350	100	1100
4	1	0	17	6	40	5	5	2	540	150	1700
4	1	1	21	7	42	5	5	3	920	220	2600
4	1	2	26	10	70	5	5	4	1600	400	4600
4	2	0	22	7	50	5	5	5	≥2400	700	---

Sources: pp 9–51, *Standard Methods for the Examination of Water and Wastewater*, 20<sup>th</sup> Edition (1998). M. J. Taras, A. E. Greenberg, R. D. Hoak, and M. C. Rand, eds. American Public Health Association, Washington, D.C. Copyright 1998, American Public Health Association, and *Bacteriological Analytical Manual (BAM)*, 8th Edition, Food and Drug Administration, 1998.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

### Presumptive Test

Using Table 49.1, determine and record the MPN.

**Example:** If gas appeared in all five tubes labeled LB2X-10, in two of the tubes labeled LB1X-1, and in one labeled LB1X-0.1, the series would be read as 5-2-1. From the MPN table, such a reading would indicate approximately 70 microorganisms per 100 ml of water, with a 95% probability that between 22 and 170 microorganisms are present.

Water Sample	GAS															Reading	MPN	95% Probability Range
	LB2X-10					LB1X-1					LB1X-0.1							
	Tube					Tube					Tube							
	1	2	3	4	5	1	2	3	4	5	1	2	3	4	5			
Sewage																		
Pond																		
Tap																		


### Confirmed Test

Water Sample	COLIFORMS		Potable	Nonpotable
	EMB Plate	Endo Agar Plate		
Sewage				
Pond				
Tap				

## Completed Test

Water Source	Lactose Broth A/G (+) or (-)	GRAM STAIN	POTABILITY	
		Reaction/Morphology	Potable	Nonpotable
Sewage				
Pond				
Tap				

## Review Questions

1. What is the rationale for selecting *E. coli* as the indicator of water potability?
2. Why is this procedure qualitative rather than quantitative?
3. Explain why it is of prime importance to analyze water supplies that serve industrialized communities.
4.  Account for the presence of microorganisms in natural bodies of water and sewage systems. What is their function? Explain.

## LEARNING OBJECTIVES

*Once you have completed this experiment, you should be able to*

1. Determine the quality of water samples using the membrane filter method.

## Principle

Bacteria-tight **membrane filters** capable of retaining microorganisms larger than 0.45 micrometer ( $\mu\text{m}$ ) are frequently used for analysis of water. These filters offer several advantages over the conventional, multiple-tube method of water analysis: (1) results are available in a shorter period of time, (2) larger volumes of sample can be processed, and (3) because of the high accuracy of this method, the results are readily reproducible. A disadvantage involves the processing of turbid specimens that contain large quantities of suspended materials; particulate matter clogs the pores and inhibits passage of the specific volume of water.

A water sample is passed through a sterile membrane filter that is housed in a special filter apparatus contained in a suction flask. Following filtration, the filter disc that contains the trapped microorganisms is aseptically transferred to a sterile Petri dish containing an absorbent pad saturated with a selective, differential liquid medium. Following incubation, the colonies present on the filter are counted with the aid of a microscope.

This experiment is used to analyze a series of dilutions of water samples collected upstream and downstream from an outlet of a sewage treatment plant. A total count of coliform bacteria determines the potability of the water sources. Also, the types of fecal pollution, if any, are established by means of a fecal coliform count, indicative of human pollution, and a fecal streptococcal count, indicative of pollution from other animal origins. The ratio of the fecal coliforms to fecal

streptococci per milliliter of sample is interpreted as follows: Between 2 and 4 indicates human and animal pollution;  $>4$  indicates human pollution; and  $<0.7$  indicates poultry and livestock pollution.

## CLINICAL APPLICATION

### Rapid Water Analysis

In the late 1950s the membrane filter method was introduced as an alternative to the most probable number method (MPN). Microbiological analysis of water by the membrane filter procedure is a rapid method that isolates discrete bacteria that are able to be accurately counted, whereas the MPN method only allows for the approximate determination of the number of organisms and does not separate species without further testing.

## AT THE BENCH



## Materials

### Cultures

Water samples collected upstream (labeled U) and downstream (labeled D) from an outlet of a sewage treatment plant.

### Media

Per designated student group for analysis of one water sample: one 20-ml tube of m-Endo broth, one 20-ml tube of m-FC broth, one 20-ml tube of KF broth, four 90-ml sterile water blanks, and one 300-ml flask of sterile water.

### Equipment

Sterile membrane filtration apparatus (i.e., Millipore<sup>®</sup>, Pall<sup>®</sup> Gelman; sterile, plastic, disposable membrane filters), 1-liter suction flask, 15 sterile membrane filters and absorbent pads, 15 sterile 50-mm Petri dishes, 12 10-ml pipettes,

mechanical pipetting device, small beaker of 95% alcohol, membrane forceps, waterproof tape, watertight plastic bags, 44.5°C waterbath, dissecting microscope, and glassware marking pencil.

## Procedure Lab One

The following instructions are for analysis of one of the provided water samples using the Millipore system. Different samples may be assigned to individual groups.



**Use disposable gloves when handling the water samples in this experiment.**

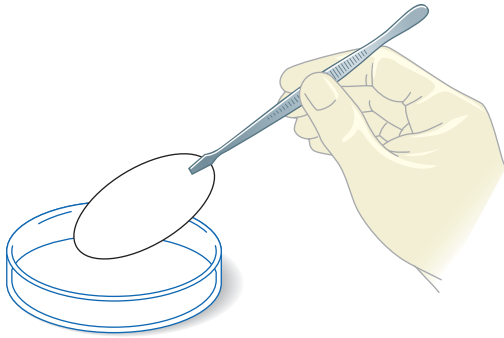
1. Label the four 90-ml water blanks with the source of the water sample and dilution ( $10^{-1}$ ,  $10^{-2}$ ,  $10^{-3}$ , and  $10^{-4}$ ).
2. Using 10-ml pipettes, aseptically perform a 10-fold serial dilution of the assigned undiluted water sample, using the four 90-ml water blanks to effect the  $10^{-1}$ ,  $10^{-2}$ ,  $10^{-3}$ , and  $10^{-4}$  dilutions.
3. Arrange the 15 Petri dishes into three sets of five plates. Label each set as follows:
  - a. For total coliform count (TCC) and dilutions (undiluted,  $10^{-1}$ ,  $10^{-2}$ ,  $10^{-3}$ , and  $10^{-4}$ ).
  - b. For fecal coliform count (FCC) and dilutions as in Step 3a.
  - c. For fecal streptococcal count (FSC) and dilutions as in Step 3a.
3. Aseptically assemble the sterile paper-wrapped membrane filter unit as follows:
  - a. Unwrap and insert the sintered glass filter base into the neck of a 1-liter side-arm suction flask.
  - b. With sterile forceps, place a sterile membrane filter disc, grid side up, on the sintered glass platform.
  - c. Unwrap and carefully place the funnel section of the apparatus on top of the filter disc. Using the filter clamp, secure the funnel to the filter base.
  - d. Attach a rubber hose from the side-arm on the vacuum flask to a vacuum source.
4. Using the highest sample dilution ( $10^{-4}$ ) and a pipette, place 20 ml of the dilution into the funnel and start the vacuum.
  - a. When the entire sample has been filtered, wash the inner surface of the funnel with 20 ml of sterile water.
5. Disconnect the vacuum, unclamp the filter assembly, and with sterile forceps, remove the membrane filter.
6. Place the filter on the medium-saturated pad in the Petri dish labeled TCC,  $10^{-4}$ .
7. Aseptically place a new membrane on the platform, reassemble the filtration apparatus, and repeat Steps 4 through 6 twice, adding the filter discs to the  $10^{-4}$  dilution plates labeled FCC and FSC.
8. Repeat Steps 4a through 7, using 20 ml of the  $10^{-3}$ ,  $10^{-2}$ , and  $10^{-1}$  dilutions and the undiluted samples.
9. Incubate the plates in an inverted position as follows:
  - a. TCC and FSC plates for 24 hours at 37°C.
  - b. FCC plates sealed with waterproof tape and placed in a weighted watertight plastic bag, which is then submerged in a 44.5°C waterbath for 24 hours.

## Membrane Filter Technique


Refer to **Figure 50.1** as you read the instructions below.

1. Using sterile forceps dipped in 95% alcohol and flamed, add a sterile absorbent pad to all Petri dishes.
2. With sterile 10-ml pipettes, aseptically add the following:
  - a. To each pad in the plates labeled TCC, 2 ml of m-Endo broth.
  - b. To each pad in the plates labeled FCC, 2 ml of m-FC broth.


## PROCEDURE



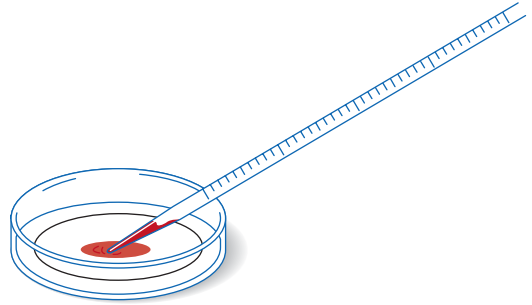
- 1 Aseptically place an absorbent pad in a 50-mm petri dish.



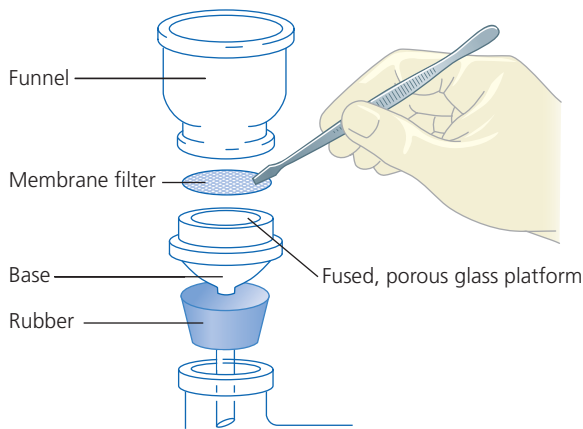
- 2 Saturate the absorbent pad with the specified selective broth medium.



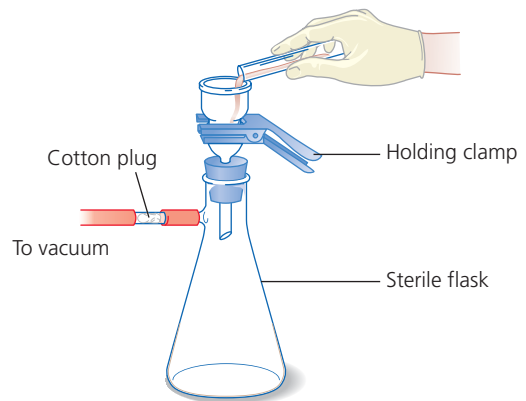
- 3 Assemble the filter apparatus and insert membrane filter.



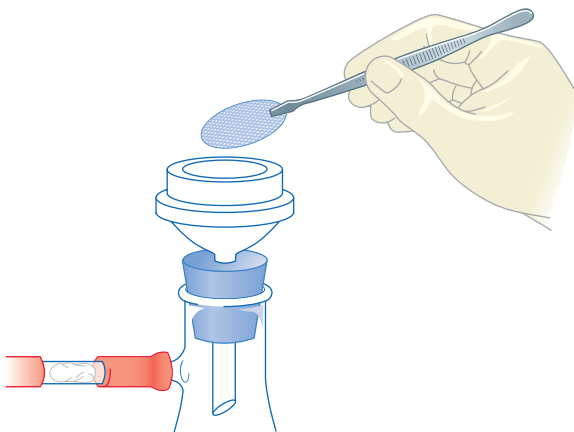- 4 Pour test sample into funnel, filter under vacuum, and rinse with sterile water.




- 5 Aseptically remove filter.

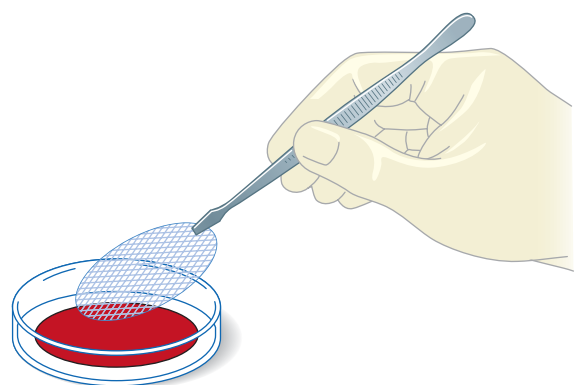

- 6 Place filter in Petri dish on top of medium-saturated pad and incubate.

**Figure 50.1** Membrane filter technique



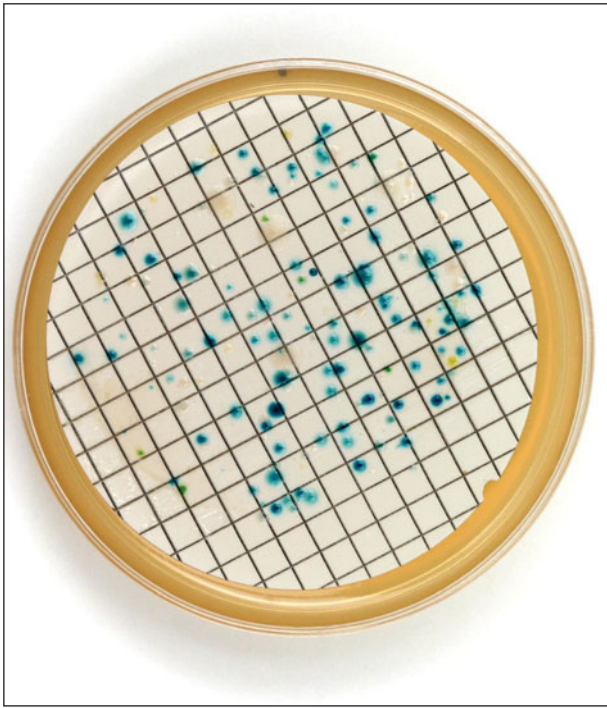

**Figure 50.2** Development of colonies on a membrane filter following incubation

## Procedure Lab Two

- Using sterile membrane forceps, remove the filter discs from the Petri dishes and allow to dry on absorbent paper for 1 hour.
- Using membrane forceps, place each dry filter disc into its Petri dish cover. *Keep the discs within the covers at all times for further observation.*
- Examine all filter discs under a dissecting microscope. Refer to **Figure 50.2**, which shows colonies developing on the membrane filter. Perform colony counts on each set of discs as follows:
  - TCC: Count colonies on m-Endo agar that present a golden metallic sheen (performed on a disc showing 20 to 80 of these colonies).
  - FCC: Count colonies on m-FC agar that are blue (performed on a disc showing 20 to 60 of these colonies).
  - FSC: Count colonies on KF agar that are pink to red (performed on a disc showing 20 to 100 of these colonies).Dilution samples that show fewer colonies than indicated are designated as TFTC, and those showing a greater number of colonies are designated as TNTC.
- For each of the three counts, determine the number of fecal organisms present in 100 ml of the water sample, using the following formula:
$$\frac{\text{colony count} \times \text{dilution factor}}{\text{ml of sample used}} \times 100$$
- Record your results in the Lab Report.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

	UPSTREAM WATER					
	TCC		FCC		FSC	
Dilution	Count	Cells/100 ml	Count	Cells/100 ml	Count	Cells/100 ml
Undiluted						
10 <sup>-1</sup>						
10 <sup>-2</sup>						
10 <sup>-3</sup>						
10 <sup>-4</sup>						

	DOWNSTREAM WATER					
	TCC		FCC		FSC	
Dilution	Count	Cells/100 ml	Count	Cells/100 ml	Count	Cells/100 ml
Undiluted						
10 <sup>-1</sup>						
10 <sup>-2</sup>						
10 <sup>-3</sup>						
10 <sup>-4</sup>						

Determine the fecal coliform to fecal streptococcal (FC:FS) ratio. Record your results in the chart below.

Dilution	UPSTREAM WATER			DOWNSTREAM WATER		
	Cells/ml*		FC:FS Ratio	Cells/ml*		FC:FS Ratio
	FCC	FSC		FCC	FSC	
Undiluted						
1021						
1022						
1023						
1024						

$$*\text{Cells/ml} = \frac{\text{Cells/100 ml}}{100}$$

Based on your FC:FS ratio, indicate the type of fecal pollution, if any, in the two samples:

a. Upstream water sample:

b. Downstream water sample:

## Review Questions

1. What are the advantages of the membrane filter method in the analysis of water samples?
2. What are the disadvantages of the membrane filter method?
3. What is the purpose of determining the FC:FS ratio?
4. Cite some other microbiological applications of the membrane filter technique in environmental studies.

# Microbiology of Soil

## LEARNING OBJECTIVES

Once you have completed the experiments in this section, you should be able to

1. Understand the characteristics and activities of soil microorganisms.
2. Enumerate soil microorganisms.
3. Demonstrate the ability of some soil microorganisms to produce antibiotics.
4. Demonstrate the use of enrichment cultures for the isolation of specific soil microorganisms.

## Introduction

Soil is often thought of as an inert substance by the average layperson. However, contrary to this belief, it serves as a repository for many life forms, including a huge and diverse microbial population. The beneficial activities of these soil inhabitants far outweigh their detrimental effects.


Life on this planet could not be sustained in the absence of microorganisms that inhabit the soil. This flora is essential for degradation of organic matter deposited in the soil, such as dead plant and animal tissues and animal wastes. Hydrolysis of these macromolecules by microbial enzymes supplies and replenishes the soil with basic elemental nutrients. By means of enzymatic transformations, plants assimilate these nutrients into organic compounds essential for their growth and reproduction. In turn, these plants serve as a source of nutrition for animals and man. Thus, many soil microorganisms play a vital role in a number of elemental cycles, such as the nitrogen cycle, the carbon cycle, and the sulfur cycle.

## Nitrogen Cycle

The nitrogen cycle is concerned with the enzymatic conversion of complex nitrogenous compounds in the soil and atmosphere into nitrogen compounds that plants are able to use for the

synthesis of essential macromolecules such as nucleic acids, amino acids, and proteins. The four distinct phases in this cycle are as follows:

1. **Ammonification:** Soil microorganisms sequentially degrade nitrogenous organic compounds derived from dead plants and animals deposited in the soil. The degraded nitrogenous organic compounds are converted to inorganic nitrogen compounds and then to ammonia.
2. **Nitrification:** In this two-step process, (1) ammonia is oxidized to nitrite ions ( $\text{NO}_2^-$ ) by an aerobic species of *Nitrosomonas*, and then (2) nitrites are converted to nitrate ions ( $\text{NO}_3^-$ ) by another aerobic species, *Nitrobacter*. Nitrates are released into the soil and are assimilated as a nutritional source by plants.
3. **Denitrification:** Nitrates ( $\text{NO}_3^-$ ) that are not used by plants are reduced to gaseous nitrogen ( $\text{N}_2\uparrow$ ) and are liberated back into the atmosphere by certain groups of microorganisms.
4. **Nitrogen fixation:** This vital process involves the chemical combination of gaseous nitrogen ( $\text{N}_2\uparrow$ ) with other elements to form fixed nitrogen (nitrogen-containing compounds), which are useful for plant growth. The two types of microorganisms involved in this process are free-living and symbiotic. Free-living microorganisms include *Azotobacter*,

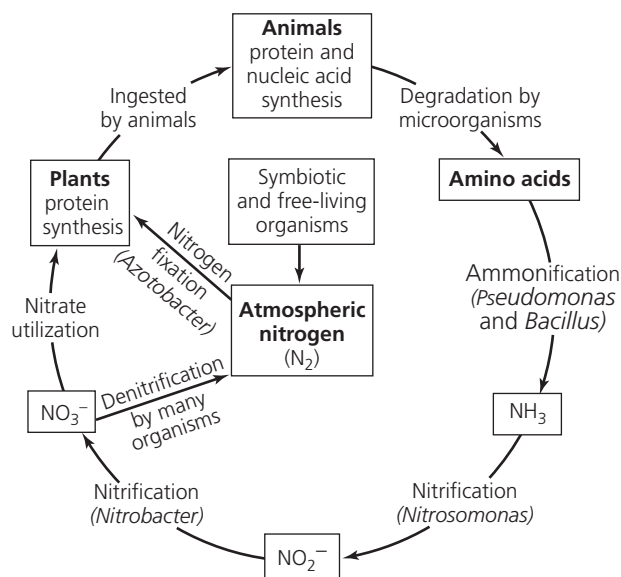

**Figure P12.1** The nitrogen cycle

*Pseudomonas*, *Clostridium*, and *Bacillus* as well as some species of yeast. Symbiotic microorganisms, such as *Rhizobium*, grow in tumorlike nodules in the roots of leguminous plants, and use nutrients in the plant sap to fix gaseous nitrogen as ammonia for its subsequent assimilation into plant proteins. Animals then consume the leguminous plants and convert plant protein to animal protein, completing the process. The nitrogen cycle is shown in **Figure P12.1**.

## Carbon Cycle

Carbon dioxide is the major carbon source for the synthesis of organic compounds. The carbon cycle is basically represented by the following two steps:

1. Oxidation of organic compounds to carbon dioxide with the production of energy and heat by heterotrophs.
2. Fixation of carbon dioxide into organic compounds by green plants and some bacteria, the autotrophic soil flora.

## Sulfur Cycle

Elemental sulfur and proteins cannot be utilized by plants for growth. They must first undergo enzymatic conversions into inorganic

sulfur-containing compounds. The basic steps in the sulfur cycle are

1. Degradation of proteins into hydrogen sulfide (H<sub>2</sub>S) by many heterotrophic microorganisms.
2. Oxidation of H<sub>2</sub>S to sulfur (S) by a number of bacterial genera, such as *Beggiatoa*.
3. Oxidation of sulfur to utilizable sulfate (SO<sub>4</sub><sup>2-</sup>) by several chemoautotrophic genera, such as *Thiobacillus*.

Some soil microorganisms also play a role in the enzymatic transformation of other elements, such as phosphorus, iron, potassium, zinc, manganese, and selenium. These biochemical changes make the minerals available to plants in a soluble form.

Many members of the soil flora, because of their fermentative and synthetic capabilities, play an important role in the synthesis of a variety of industrial products:

1. **Food.** *Penicillium* spp. are used in the production of such cheeses as Camembert, Roquefort, and Brie.
2. **Beverages.** *Saccharomyces* spp. are utilized in the wine, beer, and ale industries.
3. **Vitamins.** *Eremothecium ashbyii* and *Pseudomonas denitrificans*, respectively, synthesize riboflavin (vitamin B<sub>2</sub>) and cobalamin (vitamin B<sub>12</sub>).
4. **Enzymes.** Amylases, pectinases, and proteases are produced by *Aspergillus* spp.
5. **Antibiotics.** *Penicillium* spp. (penicillin), *Streptomyces* spp. (kanamycins and tetracyclines), and *Bacillus* spp. (bacitracin).
6. **Steroids.** *Rhizopus*, *Streptomyces*, and *Curvularia* are microorganisms that are used to carry out specific reactions, bioconversions, to aid in the manufacture of these lipid compounds.
7. **Industrial chemicals.** *Clostridium acetobutylicum* is used in the production of acetone and butanol, and *Aspergillus niger* is used in the synthesis of citric acid.

The major adverse effect of soil organisms is the ability of some species to produce disease in plants and animals. Soil-borne human pathogens include members of the spore-forming bacterial genera *Clostridium* and *Bacillus*, and some fungal genera, such as *Cryptococcus* and *Coccidioides*.

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be

1. Familiar with the microbial soil flora.
2. Able to determine the number of bacteria and fungi present in a soil sample.

## Principle

Soil contains myriads of microorganisms, including bacteria, fungi, protozoa, algae, and viruses. The most prevalent are bacteria, including the moldlike actinomycetes, and fungi:

<b>Simple bacteria</b>	Predominantly members of the orders Pseudomonadales and Eubacteriales
<b>Actinomycetes</b> (moldlike bacteria)	Predominantly members of the genus <i>Streptomyces</i> , characterized by pleomorphism and filamentous structure
<b>Fungi</b>	Predominantly members of the zygomycetes ( <i>Rhizopus</i> , <i>Mucor</i> , and <i>Absidia</i> ) and deuteromycetes ( <i>Penicillium</i> , <i>Aspergillus</i> , <i>Alternaria</i> , <i>Stemphylium</i> , and <i>Cladosporium</i> )

It is essential to bear in mind that the soil environment differs from one location to another and from one period of time to another. Therefore, factors such as moisture, pH, temperature, gaseous oxygen content, and organic and inorganic composition of soil are crucial in determining the specific microbial flora of a particular sample.

Just as the soil differs, microbiological methods used to analyze soil also vary. A single technique cannot be used to count all the different types of microorganisms present in a given soil

sample because no one laboratory cultivation procedure can provide all the physical and nutritional requirements necessary for the growth of a greatly diverse microbial population. In this experiment, only the relative numbers of bacteria, actinomycetes, and fungi are determined. The method used is the serial dilution–agar plate procedure described in Experiment 20. Different media are employed to support the growth of these three types of microorganisms: glycerol yeast agar for the isolation of actinomycetes, Sabouraud agar for the isolation of fungi, and nutrient agar for the isolation of bacteria. The glycerol yeast agar and Sabouraud agar are supplemented with 10  $\mu\text{g}$  of chlortetracycline (Aureomycin) per milliliter of medium to inhibit the growth of bacteria.

## CLINICAL APPLICATION

### Soil Testing

The enumeration of organisms in soil helps to establish the level of soil fertility as well as the types and kinds of pathogens it contains. From a clinical view, many bacterial pathogens originate from a soil environment. Current thought is that the ability of *Bacillus* species (for example, *B. subtilis* and *B. anthracis*) and *Mycobacterium* species (for example, *M. tuberculosis*) to survive in a soil environment—one that contains low nutrients, low moisture, and that necessitates sporulation or slow growth—aids the bacteria in infecting human tissues and surviving the immune response.

## AT THE BENCH



## Materials

### Soil

1-g sample of finely pulverized, rich garden soil in a flask containing 99 ml of sterile water; flask labeled 1:100 dilution ( $10^{-2}$ ).

## Media

Per designated student group: four glycerol yeast agar deep tubes, four Sabouraud agar deep tubes, four nutrient agar deep tubes, and two 99-ml flasks of sterile water.

## Equipment

Bunsen burner, 12 Petri dishes, Quebec colony counter, mechanical hand counter, sterile 1-ml pipettes, mechanical pipetting device, L-shaped bent glass rod, turntable (optional), 95% alcohol in a 500-ml beaker, and glassware marking pencil.

## Procedure Lab One

Refer to **Figure 51.1** as you read the following instructions.

1. Liquefy the glycerol yeast, Sabouraud, and nutrient agar deep tubes in an autoclave or by boiling. Cool the molten agar tubes and maintain in a waterbath at 45°C.
2. Divide the Petri dishes into three groups of four; using a glassware marking pencil, label the groups as nutrient agar, glycerol yeast extract agar, and Sabouraud agar. Then, label each set of Petri dishes as follows:
  - Nutrient agar:  $10^{-4}$ ,  $10^{-5}$ ,  $10^{-6}$ , and  $10^{-7}$  (to be used for enumeration of bacteria).
  - Glycerol yeast extract agar:  $10^{-3}$ ,  $10^{-4}$ ,  $10^{-5}$ , and  $10^{-6}$  (to be used for enumeration of actinomycetes).
  - Sabouraud agar:  $10^{-2}$ ,  $10^{-3}$ ,  $10^{-4}$ , and  $10^{-5}$  (to be used for enumeration of fungi).
3. With a glassware marking pencil, label the soil sample flask as Flask 1, and label the 99-ml sterile water Flasks 2 and 3.
4. Vigorously shake the provided soil sample dilution of 1:100 ( $10^{-2}$ ) approximately 30 times, with your elbow resting on the table.
5. With a sterile 1-ml pipette, transfer 1 ml of the provided soil sample dilution to Flask 2 and shake vigorously as before. The final dilution is 1:10,000 ( $10^{-4}$ ).
6. Using another sterile 1-ml pipette, transfer 1 ml of Dilution 2 to Flask 3 and shake vigorously as before. The final dilution is 1:1,000,000 ( $10^{-6}$ ).

7. Using sterile 1-ml pipettes and aseptic technique, add the proper amount of each dilution into each Petri dish as indicated in a–c and shown in Figure 51.1.

**a. For actinomycetes**—in plates labeled glycerol yeast extract agar:

Transfer 0.1 ml of Dilution 1 into plate to effect a  $10^{-3}$  dilution.

Transfer 1 ml of Dilution 2 into plate to effect a  $10^{-4}$  dilution.

Transfer 0.1 ml of Dilution 2 into plate to effect a  $10^{-5}$  dilution.

Transfer 1 ml of Dilution 3 into plate to effect a  $10^{-6}$  dilution.

**b. For molds**—in plates labeled Sabouraud agar:

Transfer 1 ml of Dilution 1 into plate to effect a  $10^{-2}$  dilution.

Transfer 0.1 ml of Dilution 1 into plate to effect a  $10^{-3}$  dilution.

Transfer 1 ml of Dilution 2 into plate to effect a  $10^{-4}$  dilution.

Transfer 0.1 ml of Dilution 2 into plate to effect a  $10^{-5}$  dilution.

**c. For bacteria**—in plates labeled nutrient agar:

Transfer 1 ml of Dilution 2 into plate to effect a  $10^{-4}$  dilution.

Transfer 0.1 ml of Dilution 2 into plate to effect a  $10^{-5}$  dilution.

Transfer 1 ml of Dilution 3 into plate to effect a  $10^{-6}$  dilution.


Transfer 0.1 ml of Dilution 3 into plate to effect a  $10^{-7}$  dilution.

8. Check the temperature of the molten agar medium to be sure that the temperature is 45°C. Remove the tubes from the waterbath and wipe the outside surface dry with a paper towel. Using the pour-plate technique, pour the liquefied agar into the plates as shown in Figure 20.2 on page 138 and rotate gently to ensure uniform distribution of the cells in the medium.
9. Incubate the plates in an inverted position at 25°C. Perform colony counts on nutrient agar plate cultures in 2 to 3 days and on the remaining agar plate cultures in 4 to 7 days.

## PROCEDURE


### (a) For actinomycetes:

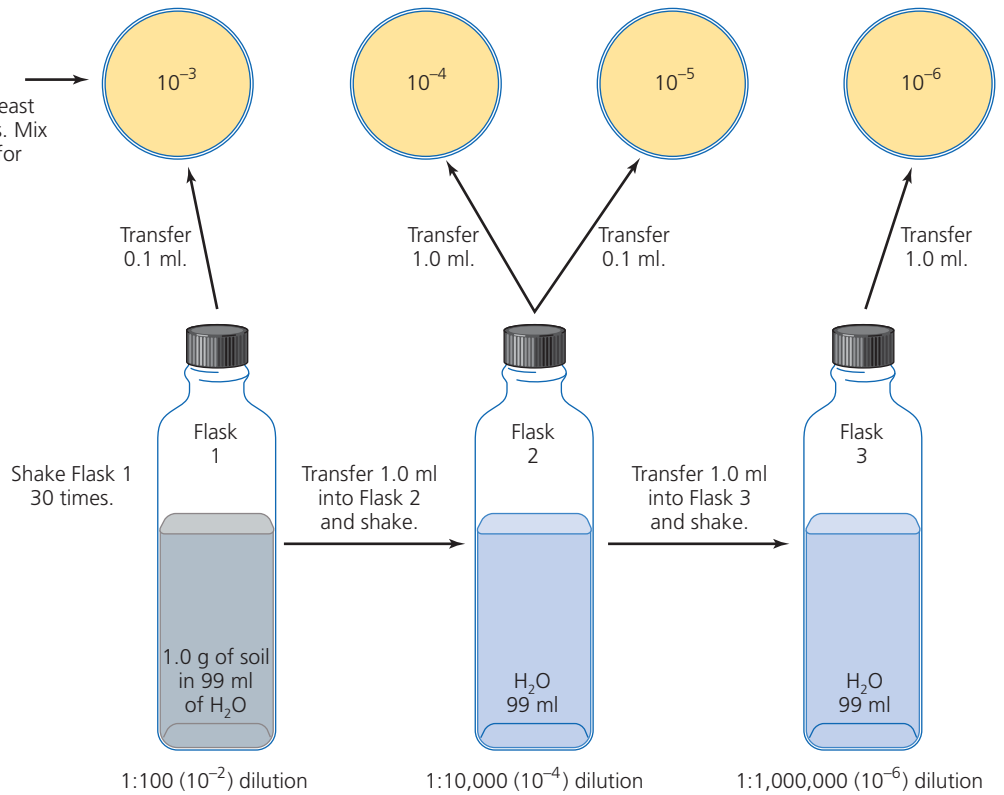Pour 45°C glycerol yeast agar into Petri plates. Mix by rotation of plate for serial dilution–agar plate method.



### (b) For molds:

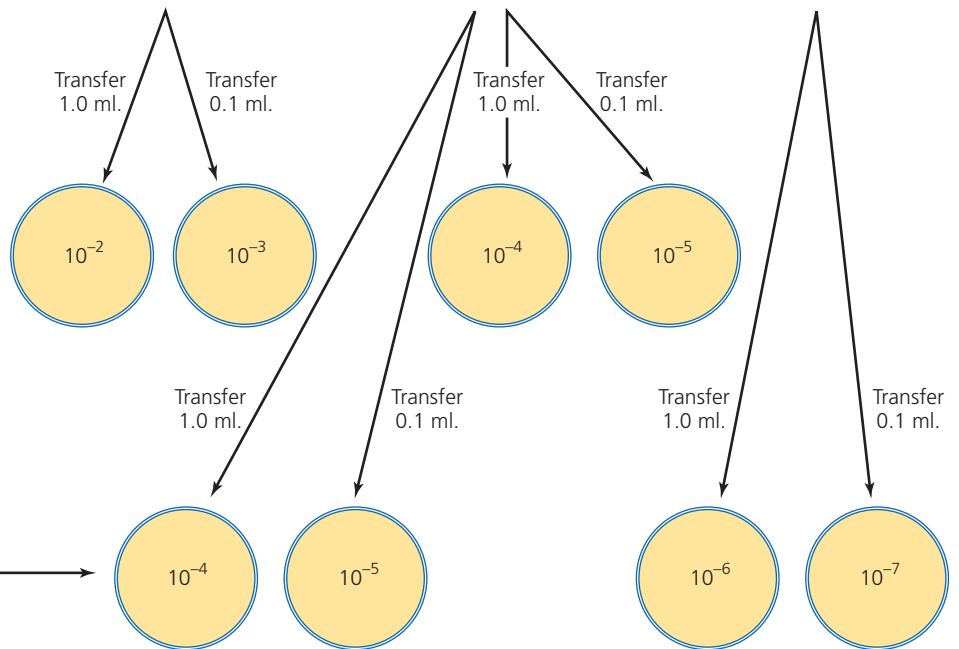Pour 45°C Sabouraud agar into Petri plates. Mix by rotation of plate for serial dilution–agar plate method.



### (c) For bacteria:

Pour 45°C nutrient agar into Petri plates. Mix by rotation of plate for serial dilution–agar plate method.



**Figure 51.1** Procedure for enumeration of soil microorganisms



## Procedure Lab Two

1. Using an electronic colony counter or a Quebec colony counter and a mechanical hand counter, observe all the colonies on each nutrient agar plate 2 to 3 days after incubation begins. Plates with more than 300 colonies cannot be counted and should be designated as **too numerous to count (TNTC)**; plates with fewer than 30 colonies should be designated as **too few to count (TFTC)**. Count only plates with between 30 and 300 colonies.
2. Determine the number of organisms per milliliter of original culture on all plates other

than those designated as TFTC or TNTC by multiplying the number of colonies counted by the dilution factor. Refer to Experiment 20 for examples of the calculation of cell counts.

3. Record your observations and calculated cell counts per gram of sample in the Lab Report chart.

## Procedure Lab Three

1. Repeat Steps 1–3 from Lab Two for the Sabouraud agar and glycerol yeast extract agar plates 4 to 7 days after incubation begins.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_


## Observations and Results

Organism	Dilution	Number of Colonies	Organisms per Gram of Soil
Bacteria	$10^{-4}$		
	$10^{-5}$		
	$10^{-6}$		
	$10^{-7}$		
Actinomycetes	$10^{-3}$		
	$10^{-4}$		
	$10^{-5}$		
	$10^{-6}$		
Molds	$10^{-2}$		
	$10^{-3}$		
	$10^{-4}$		
	$10^{-5}$		

Based on your results, which of the three types of soil organisms was most abundant in your sample? Least abundant?

## Review Questions

1. Would you expect to be able to duplicate your results if a soil sample were taken from the same location at a different time of the year? Explain.

2. In the experiment performed, why wasn't the same medium used for enumeration of all three types of soil organisms?
  
3. Would you expect to be able to isolate an anaerobic organism from any of your cultures? Explain.
  
4. Explain why most microorganisms are present in the upper layers of the soil.
  
5.  Following the nuclear disaster at Chernobyl, the regional microbial flora was destroyed. What impact did this have on higher forms of plant and animal life in this area?

# Isolation of Antibiotic-Producing Microorganisms and Determination of Antimicrobial Spectrum of Isolates

EXPERIMENT

52

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be able to

1. Isolate antibiotic-producing microorganisms.
2. Determine the spectrum of antimicrobial activity of the isolated antibiotic.

## Principle

**Soil** is the major repository of microorganisms that produce **antibiotics** capable of inhibiting the growth of other microorganisms. Clinically useful antibiotics have been isolated from five groups of soil microorganisms—*Streptomyces*, *Amycolatopsis* (including some species formerly classified as *Streptomyces*), *Bacillus*, *Penicillium*, and *Acremonium*—that represent three microbial types, namely, actinomycetes, true bacteria, and molds.

Although soils from all parts of the world are continually screened in industrial laboratories for the isolation of new antibiotic-producing microorganisms, industrial microbiology is directing its energies toward chemical modification of existing antibiotic substances. This is accomplished by adding or replacing chemical side chains, reorganizing intramolecular bonding, or producing mutant microbial strains capable of excreting a more potent form of the antibiotic. The establishment of chemical congeners has been responsible for the circumvention of antibiotic resistance, minimizing adverse side effects in the host and increasing the effective spectrum of a given antibiotic.

In Part A of this experiment, you will use the **crowded-plate** technique for isolation of antibiotic-producing microorganisms from two soil samples, one of which is seeded with

*Streptomyces griseus* to serve as a positive control. Figure 52.1 illustrates the procedure to be followed. In Part B, isolates exhibiting antibiotic activity will be screened against several different microorganisms to establish their effectiveness.

## CLINICAL APPLICATION

### Testing New Antibiotics

Soil is the major reservoir housing microorganisms that produce antibiotics, which are used offensively to reduce competition for available nutrients. The most prolific antibiotic producers are within the phylum *Actinobacteria*. The genus *Streptomyces* are the major producers of currently used antibiotics (for example, neomycin, streptomycin, and tetracyclin) along with the genus *Actinomycetes* (for example, erythromycin). Interestingly, hundreds of new antibiotics are isolated annually using procedures such as the crowded plate technique, but most have a limited spectrum and only a few are found to be clinically acceptable.

## AT THE BENCH



## Materials


### Cultures

For Part B: 24-hour Trypticase soy broth cultures of *Escherichia coli*, *Staphylococcus aureus*, *Mycobacterium smegmatis*, and *Pseudomonas aeruginosa*.

### Soil Suspensions

For Part A: 1:500 dilution of soil sample suspension (0.1 g of soil per 50 ml of tap water) to serve as an unknown; 1:500 dilution of soil sample

## PROCEDURE

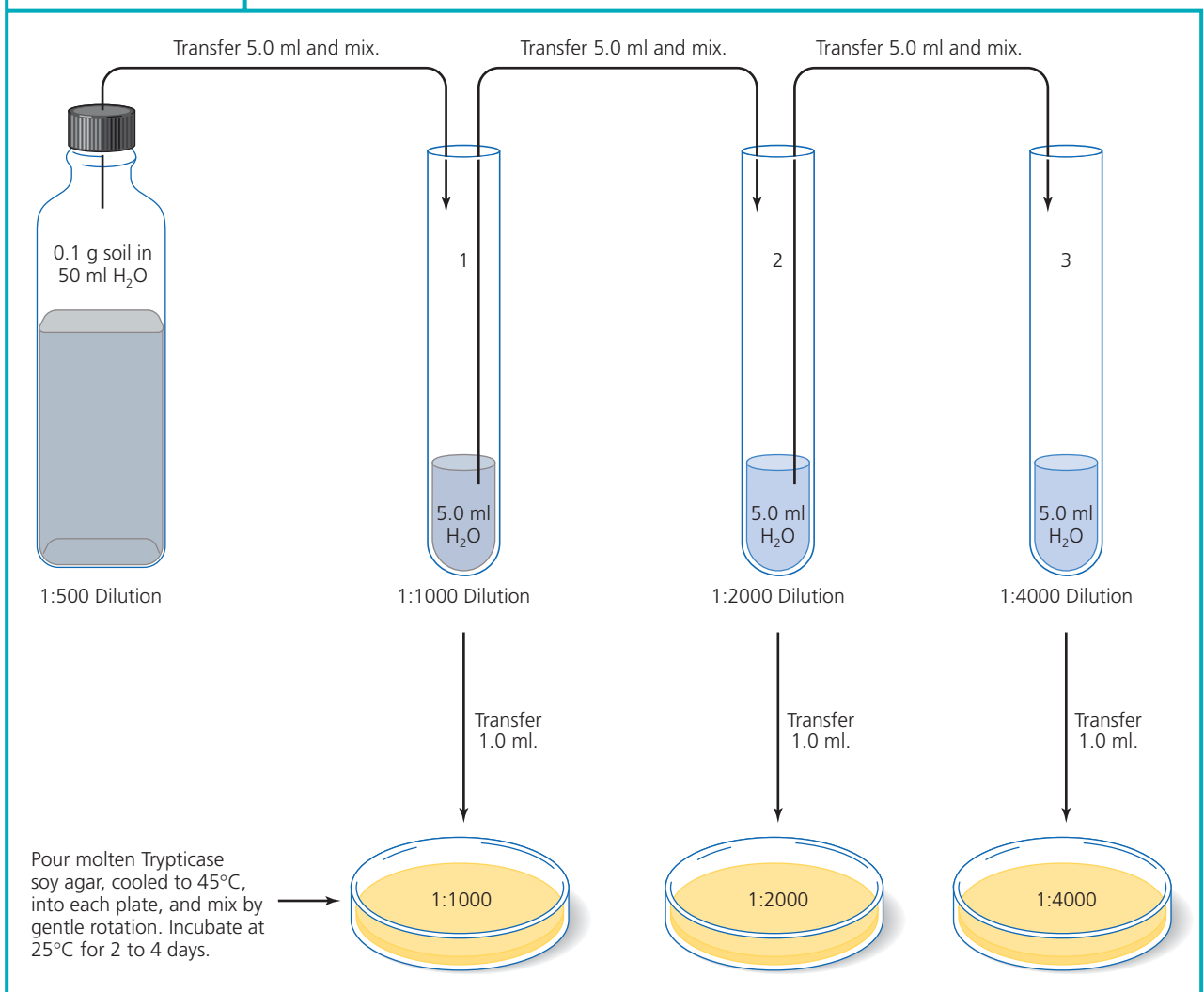

**Figure 52.1** Crowded-plate technique for isolation of antibiotic-producing microorganisms

seeded with *S. griseus* (0.1 g of soil per 50 ml of tap water) to serve as a positive control.

### Media

Per designated student group: Part A: Six 15-ml Trypticase soy agar deep tubes, and two Trypticase soy agar slants. Part B: Two Trypticase soy agar plates.

### Equipment

Part A: 500-ml beaker, test tubes, test tube rack, sterile Petri dishes, inoculating needle, hot plate, thermometer, 1-ml and 5-ml pipettes, mechanical pipetting device, and magnifying hand lens. Part B: Bunsen burner, inoculating loop, and glassware marking pencil.

## PART A Isolation of Antibiotic-Producing Microorganisms

### Procedure Lab One

1. Label two sets of three sterile Petri dishes with the types of soil samples being used and dilutions (1:1000, 1:2000, and 1:4000).
2. Place six Trypticase soy agar deep tubes into a beaker of water and bring to 100°C on a hot plate. Once agar is liquefied, add cool water to the waterbath. Cool to 45°C, checking the temperature with a thermometer.
3. Prepare a serial dilution of the unknown and positive control 1:500 soil samples as follows (refer to **Figure 52.1**):

- Label three test tubes 1, 2, and 3. With a pipette, add 5 ml of tap water to each tube.
  - Shake the provided 1:500 soil sample thoroughly for 5 minutes to effect a uniform soil-water suspension.
  - Using a 5-ml pipette, transfer 5 ml from the 1:500 dilution to Tube 1 and mix. The final dilution is 1:1000.
  - Using another pipette, transfer 5 ml from Tube 1 to Tube 2 and mix. The final dilution is 1:2000.
  - Using another pipette, transfer 5 ml from Tube 2 to Tube 3 and mix. The final dilution is 1:4000.
  - Using separate 1-ml pipettes, transfer 1 ml of the 1:1000, 1:2000, and 1:4000 dilutions to their appropriately labeled Petri dishes.
  - Pour one tube of molten Trypticase soy agar, cooled to 45°C, into each plate and mix by gentle rotation.
  - Allow all plates to solidify.
- Incubate all plates in an inverted position for 2 to 4 days at 25°C.

## Procedure Lab Two

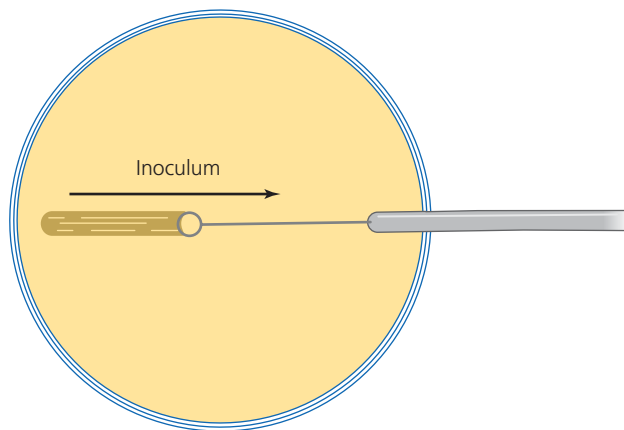
- Examine all crowded-plate dilutions for colonies exhibiting zones of growth inhibition. Use a hand magnifying lens if necessary. Record in the Lab Report the number of colonies showing zones of inhibition.
- Aseptically isolate one colony showing a zone of growth inhibition from each soil culture with an inoculating needle and streak onto Trypticase soy agar slants labeled with the soil sample from which the isolate was obtained.
- Incubate the slants for 2 to 4 days at 25°C. These will serve as stock cultures of antibiotic-producing isolates to be used in Part B.

## PART B Determination of Antimicrobial Spectrum of Isolates

### Procedure Lab One

- Label the Trypticase soy agar plates with the soil sample source of the isolate.
- Using aseptic technique, make a single-line streak inoculation of each isolate on the


surface of an agar plate so as to divide the plate in half as shown:



- Incubate the plates in an inverted position for 3 to 5 days at 25°C.

### Procedure Lab Two

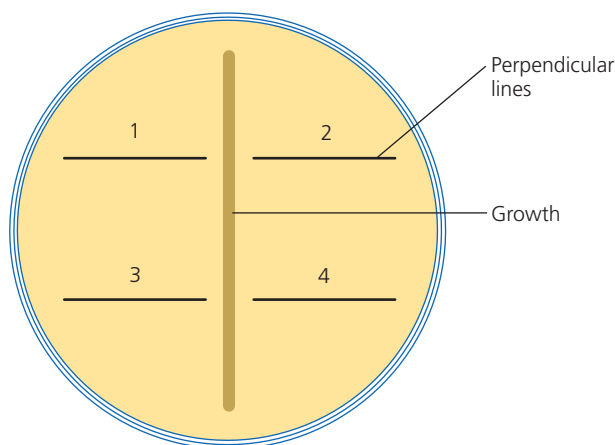- Following incubation, on the bottom of each plate draw four lines perpendicular to the growth of the antibiotic-producing isolate as shown:



- Aseptically make a single-line streak inoculation of each of the four test cultures following the inoculation template on each plate. Start close to, but not touching, the growth of the antibiotic-producing isolate and streak toward the edge of the plate.
- Incubate the plates in an inverted position for 24 hours at 37°C.

### Procedure Lab Three

- Examine all plates for inhibition of test organisms, and record your observations in the Lab Report.

*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

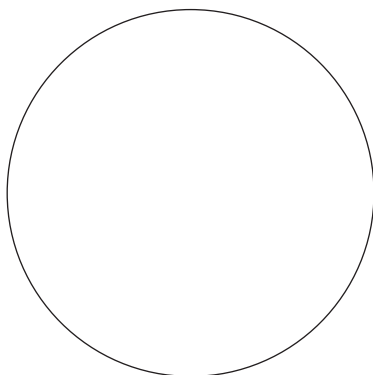
## Observations and Results

### PART A: Isolation of Antibiotic-Producing Microorganisms


Soil Sample	NUMBER OF COLONIES SHOWING INHIBITION ZONE		
	Dilutions		
	1:1000	1:2000	1:4000
Unknown			
Positive control			

### PART B: Determination of Antimicrobial Spectrum of Isolates

1. Draw a representation of the observed antibiotic activity against the test organisms.



Antibiotic-Producing Isolate 1

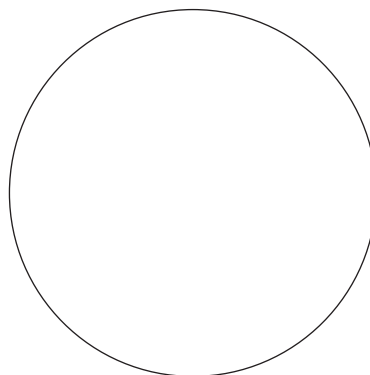


Antibiotic-Producing Isolate 2

2. Based on your observations, record in the chart the presence (+) or absence (–) of antibiotic activity against each of the test organisms and the spectrum of antimicrobial activity (broad or narrow).

Soil Sample	TEST ORGANISMS				Spectrum
	<i>E. coli</i> Gram-Negative	<i>S. aureus</i> Gram-Positive	<i>P. aeruginosa</i> Gram-Negative	<i>M. smegmatis</i> Acid-Fast	
Unknown					
Positive control					



## Review Questions

1. Why is it frequently advantageous to modify antibiotics in industrial laboratories?
2. Is the ability to produce antibiotics limited only to bacterial species? Explain.
3.  Do you feel that sufficient test organisms were used in Part B to determine fully the spectrum of activity of each isolated antibiotic? Explain.

# Isolation of *Pseudomonas* Species by Means of the Enrichment Culture Technique

## LEARNING OBJECTIVES

Once you have completed this experiment, you should

1. Understand the enrichment culture technique for the isolation of a specific microbial cell type.

## Principle


The enrichment culture technique is used for the isolation of a specific type of microorganism from an environment that is replete with different types of microbes. In such an environment, the desired organism may be present only in very small numbers because of the competitive activities of this diverse microbial population. Under these circumstances, the use of conventional enriched media is not suitable for the selection of a specific cell type. These special-purpose media are supplemented with a variety of enriching nutrients capable of supporting the growth of many organisms rather than a single cell type in the test sample. Enrichment broths, on the other hand, are designed to contain a limited number of specific substrates that will preferentially promote the growth of the desired microorganisms.

The enrichment culture technique employs such a specifically designed enrichment broth for the initial inoculation of the test sample. Once growth occurs in the primary culture, it is sequentially transferred into a fresh medium of the same composition until the desired microorganisms are predominant in the culture. These organisms are capable of exponential growth because of their ability to adapt to the medium and to enzymatically use the incorporated substrate(s) as an

energy source. Most of the competitors, however, are incapable of utilizing the substrate(s) and therefore remain in the lag phase of the growth curve. In some instances the organisms to be isolated do not grow more rapidly than their competitors. Instead, they produce a growth inhibitor that greatly suppresses the growth of the competing population. After the serial transfer through the broth medium, the culture is streaked on an agar plate of the same composition as the enrichment broth for the isolation and subsequent identification of the discrete colonies.

The use of the enrichment culture technique has a wide range of applications in clinical, industrial, and environmental microbiology. Enrichment methods may be used to isolate and cultivate specific soil microorganisms for the production of industrial products such as steroids, enzymes, and vitamins. Likewise, a beneficial environmental application may involve the isolation by enrichment of petroleum-utilizing microorganisms such as *Pseudomonas* that would be capable of degrading environmentally destructive oil spills in waterways.

In this experimental procedure, a compost or a rich garden soil sample will be used to isolate *Pseudomonas* species by means of the enrichment culture procedure. Members of the genus *Pseudomonas* can utilize mandelic acid aerobically as their sole carbon and energy source. Therefore, this compound is the most important factor in the enrichment broth, which also contains a number of inorganic salts. The pseudomonads are gram-negative, motile organisms that generally produce a diffusible yellow-green pigment. In addition, they commonly reduce nitrates ( $\text{NO}_3^-$ ) and produce an alkaline or proteolytic reaction in litmus milk. The schema for the experimental procedure to be followed is illustrated in **Figure 53.1**.

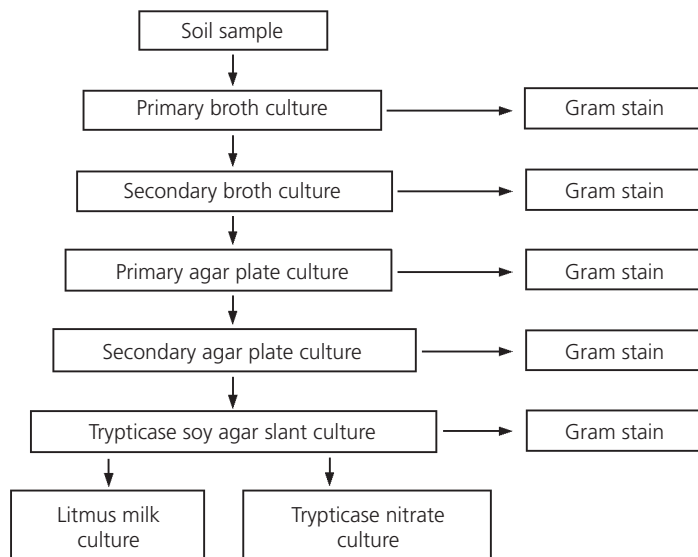

**Figure 53.1** Enrichment culture procedure schema

## CLINICAL APPLICATION

### Medical Use for the Enrichment Culture Technique

Medically, the enrichment culture technique is routinely used for the isolation of intestinal pathogens from fecal samples when these organisms may be present only in low concentrations during the infectious process. With hundreds of different bacterial species composing our intestinal flora, identifying a new bacterial pathogen, such as *Salmonella* or a new strain of *E. coli*, within that population through normal plating techniques may not be possible. By increasing the number of bacteria present in a media that is enriched, thus lowering competition, bacterial species with low numbers may increase their percentage of the population and increase the chances of their identification.

## Media

Per designated student group: Two Erlenmeyer flasks containing 20 ml of basal salts broth supplemented with 2 ml of 2.5% mandelic acid, two agar plates of the same composition as the broth, one Trypticase nitrate broth, one litmus milk, and one Trypticase soy agar slant.

## Reagents

Crystal violet, Gram's iodine, 95% ethanol, safranin, Solution A (sulfanilic acid), Solution B (alpha-naphthylamine), and zinc powder. *Note: Solutions A and B are not Barritt's reagent.*

## Equipment

Sterile 10-ml, 5-ml, and 1-ml pipettes, mechanical pipetting device, microspatula, Bunsen burner, staining tray, glass slides, lens paper, bibulous paper, inoculating loop, and glassware marking pencil.

## AT THE BENCH



## Materials

### Cultures

Rich garden soil or compost sample.

## Procedure Lab One

### Primary Broth Culture Preparation

1. Inoculate an appropriately labeled Erlenmeyer flask containing the enrichment broth by adding an amount of the soil sample equivalent to the size of a pea with a microspatula. Gently swirl the flask to mix the culture.

2. Incubate the primary broth culture for 24 hours at 30°C.

## Procedure Lab Two

### Secondary Broth Culture Preparation

1. Examine the primary culture for presence of growth. If growth is not present, return the flask to the incubator for an additional 24 hours.
2. If growth is present, aseptically transfer 1 ml of the primary culture to an appropriately labeled Erlenmeyer flask containing fresh enrichment medium. Swirl the flask.
3. Incubate the secondary broth culture for 24 hours at 30°C.
4. Prepare and examine a Gram-stained smear from the primary culture. Record your observations of cellular morphology and Gram reaction in the Lab Report.
5. Refrigerate the primary broth culture.

## Procedure Lab Three

### Primary Agar Plate Preparation

1. If growth is present in the secondary broth culture, aseptically perform a four-way streak inoculation on the appropriately labeled agar plate of the enrichment medium (refer to Experiment 3).
2. Incubate the agar plate culture in an inverted position for 24 hours at 30°C.
3. Prepare and examine a Gram-stained smear of the secondary broth culture. Record your observations of cellular morphology and Gram reaction in the Lab Report.
4. Refrigerate the secondary broth culture.

## Procedure Lab Four

### Secondary Agar Plate Preparation

1. Examine the primary plate culture for the presence of discrete colonies. Record your observations of the cultural characteristics of these colonies in the Lab Report. Using a discrete colony:
  - a. Aseptically prepare and examine a Gram-stained smear. Record your observations of cellular morphology and Gram reaction in the Lab Report.

- b. Aseptically perform a four-way streak inoculation on an appropriately labeled agar plate of the enrichment medium.

2. Incubate the secondary agar plate culture in an inverted position for 24 hours at 30°C.
3. Refrigerate the primary agar plate culture.

## Procedure Lab Five

### Pure Culture Isolation

1. Examine the secondary agar plate culture. Record your observations of the cultural characteristics of these colonies in the Lab Report. If the cultural characteristics of discrete colonies appear to be similar:
  - a. Prepare and examine a Gram-stained smear from a discrete colony. Record your observations of cellular morphology and Gram reaction in the Lab Report.
  - b. Pick a discrete colony and aseptically inoculate a Trypticase soy agar slant by means of a streak inoculation.
2. Incubate the agar slant culture for 24 to 48 hours at 30°C.
3. Refrigerate the secondary agar plate culture.

## Procedure Lab Six

### Genus Identification of Isolate

1. Prepare and examine a Gram-stained smear from the Trypticase agar slant culture. Record your observations of cellular morphology and Gram reaction in the Lab Report.
2. Using the Trypticase agar slant culture, aseptically inoculate the appropriately labeled tubes of Trypticase nitrate broth and litmus milk by means of a loop inoculation.
3. Incubate the litmus milk and Trypticase nitrate broth cultures for 24 to 48 hours at 30°C.

## Procedure Lab Seven

1. Observe the litmus milk culture. Determine the type of reaction that has taken place (refer to Experiment 28), and record in the Lab Report.
2. Perform the nitrate reduction test on the Trypticase nitrate broth culture (refer to Experiment 29). Record your results in the Lab Report.

*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

### Gram Reactions and Colony Characteristics

Culture	Gram Stain; Cellular Morphology	Cultural Characteristics
Primary broth culture		
Secondary broth culture		
Primary agar plate culture		
Secondary agar plate culture		
Trypticase soy agar slant culture		

### Litmus Milk Reaction

Record the type of reaction below.

### Nitrate Reduction Test

Record whether or not the organism was capable of nitrate reduction (+ or -) below.



# Bacterial Genetics

## LEARNING OBJECTIVES

Once you have completed the experiments in this section, you should be able to demonstrate the applicability of bacterial test systems in genetic-related studies. The procedures include

1. Enzyme induction.
2. Transfer of genetic material by means of conjugation.
3. Isolation of a streptomycin-resistant mutant.
4. Detection of potential chemical carcinogens.

## Introduction

In recent years, bacteria have proved to be essential organisms in research into the structure and function of DNA, the universal genetic material. Their use is predicated on the following:

1. Their haploid genetic state, which allows the phenotypic, observable expression of a genetic trait in the presence of a single mutant gene.
2. Their rapid rate of growth, which permits observation of transmission of a trait through many generations.
3. The availability of large test populations, which allows isolation of spontaneous mutants and their induction by chemical and physical mutagenic agents.
4. Their low cost of maintenance and propagation, which make it possible to perform a large number of experimental procedures.

In the following experiments, bacterial test systems are used to demonstrate enzyme induction, screening for chemical carcinogens, and the genetic phenomena of mutation and genetic transfer. The last two mechanisms introduce

genetic variability, which is essential for evolutionary survival in asexually reproducing bacterial populations.

**Point mutations** are permanent, sudden qualitative alterations in genetic material that arise as a result of the addition, deletion, or substitution of one or more bases in the region of a single gene. As a result, one or more amino acid substitutions occur during translation, and a protein that may be inactive, reduced in activity, or entirely different is synthesized. **Spontaneous mutations** are the result of the chemical and physical components in the organism's natural environment. The rate at which they occur is extremely low in all organisms. For example, in *Escherichia coli*, the spontaneous mutation rate at a single locus (specific site on the DNA) is estimated to be about  $1 \times 10^{-7}$ , and the possibility of a mutation at any locus in the genome is approximately  $1 \times 10^{-4}$ . **Induced mutations** are genetic changes resulting from the organism's exposure to an artificial physical or chemical mutagen, that is, an agent capable of inducing a mutation. The resultant mutations are of the same type that occur spontaneously;



however, their rate is increased, and in some cases dramatically so.

Transfer of genetic material and its subsequent incorporation into the bacterial genome is also a source of genetic variation in some bacteria. This transfer may occur by means of the following:

**1. Conjugation:** A mating process between “sexually” differentiated bacterial strains

that allows unidirectional transfer of genetic material.

- 2. Transduction:** A bacteriophage-mediated transfer of genetic material from one cell to another.
- 3. Transformation:** A genetic alteration in a cell, resulting from the introduction of free DNA from the environment across the cell membrane.

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be able to

1. Understand the mechanism of the lactose operon.
2. Understand the factors affecting the expression of the  $\beta$ -galactosidase gene.

## Principle

Although bacteria possess a single chromosome, each cell is capable of synthesizing hundreds of different enzymes. Studies have shown that these enzymes are not present within the cells in equal concentrations. Some enzymes, called **constitutive enzymes**, are synthesized at a constant rate regardless of conditions in the cell's environment. Synthesis of other enzymes, called **adaptive enzymes**, occurs only when necessary, and it is subject to regulatory mechanisms that are dependent on the environment. One such mechanism, **induction**, requires the presence of a substrate, the inducer, in the environment to initiate synthesis of its specific enzyme, called an **inducible enzyme**. An extensively studied inducible enzyme in *E. coli* is  **$\beta$ -galactosidase**, which acts on the disaccharide lactose to yield the monosaccharides glucose and galactose. The gene for  $\beta$ -galactosidase is a member of a cluster of genes, called an **operon**, that is involved in the metabolism of lactose. The member genes of the lactose (*lac*) operon function as a unit, all being transcribed only when the inducer, lactose, is present in the surrounding medium. See **Figure 54.1**.

To illustrate  $\beta$ -galactosidase induction, two test strains of *E. coli* will be used: a prototrophic (wild type) strain (lactose-positive) and an auxotrophic (mutant) strain (lactose-negative), which carries a mutation in the gene for  $\beta$ -galactosidase

as well as a mutation in the lactose operon regulatory gene. Both test strains will be grown in the following media:

1. Inorganic synthetic medium lacking an organic carbon and energy source that is required by the heterotrophic *E. coli*.
2. Inorganic synthetic medium plus glucose, which can be utilized by both strains as a carbon and energy source.
3. Inorganic synthetic medium plus lactose, which can be utilized only by the prototrophic strain.

Orthonitrophenyl- $\beta$ -D-galactoside (ONPG), a colorless analog of lactose, can serve as the substrate for the induction of  $\beta$ -galactosidase synthesis. As the inducer, it is hydrolyzed to galactose and a yellow nitrophenolate ion. Following a short incubation period, growth in all the cultures will be determined by spectrophotometry. Induction of  $\beta$ -galactosidase synthesis and activity will be indicated by the appearance of a yellow color in the medium following addition of ONPG, which occurs only in the presence of the nitrophenolate ion. Absence of this macroscopically visible color change indicates that enzyme induction in the lactose-negative strain did not occur.

## CLINICAL APPLICATION

### Enzyme Inducers and Cancer

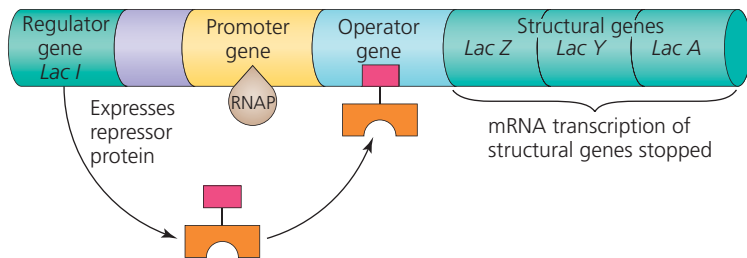
Inducer molecules can include hormones produced by the body as well as toxins and drugs. Both enzyme induction and inhibition are used by the body to control a number of interactions that play a role in many cellular reactions from digestion to cell death. One important type of current research is the deliberate induction of human enzymes that protect against environmental carcinogens. Such intervention may provide advance protection against cell damage.




(a) Lactose operon.

### Genes of lactose operon


- Regulator gene (*Lac I*): produces repressor protein.
- Promoter gene: binding site for RNA polymerase
- Operator gene: binding site for repressor protein
- Structural genes:  
*Lac Z*: codes for  $\beta$ -galactosidase  
*Lac Y*: codes for galactoside permease  
*Lac A*: codes for galactoside transacetylase



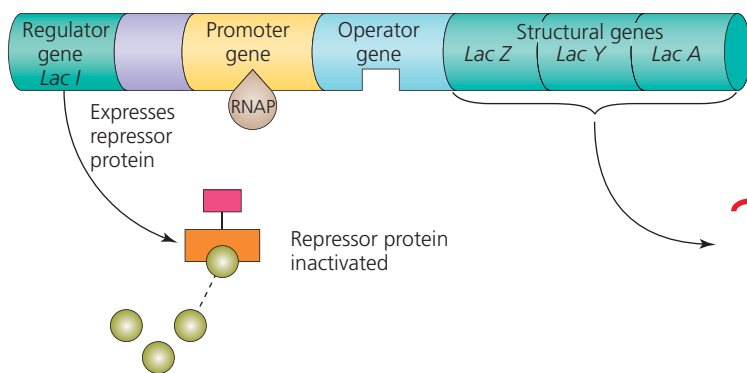(b) No lactose present. The regulator gene (*Lac I*) expresses the repressor protein. Because no lactose is present, the repressor protein binds to the operator gene, blocking the RNA polymerase and stopping mRNA transcription of structural genes.



(c) Lactose present. Lactose acts as an inducer, binding to the repressor protein and inactivating it. The repressor protein cannot bind to the operator gene; therefore mRNA transcription of the structural genes can proceed.

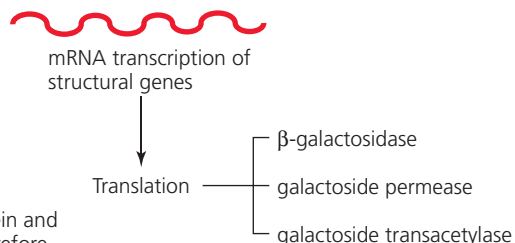

**Figure 54.1 Enzyme induction: The mechanism of operation of the lactose operon**

## AT THE BENCH



## Materials

### Cultures

25-ml inorganic synthetic broth suspensions of 12-hour nutrient agar cultures of a lactose-positive *E. coli* strain (ATCC e 23725) and a lactose-negative *E. coli* strain (ATCC e 23735) adjusted to an absorbance of 0.1 at 600 nm.

### Media

Per designated student group: dropper bottles of sterile 10% glucose, 10% lactose, and water.

## Reagents

Dropper bottles of toluene and orthonitrophenyl- $\beta$ -D-galactoside (ONPG).

## Equipment

1-ml and 5-ml sterile pipettes, mechanical pipetting device, six sterile 13-  $\times$  100-mm test tubes, test tube racks, six sterile 25-ml Erlenmeyer flasks, Bausch & Lomb Spectronic 20 spectro-photometer, shaking waterbath incubator, and glassware marking pencil.

## Procedure

1. Label three sterile test tubes and three sterile 25-ml Erlenmeyer flasks as "Lac<sup>+</sup>" (lactose-positive) and the name of the substrate to be

- added (glucose, lactose, or water). Similarly label three sterile tubes and flasks “Lac<sup>-</sup>” (lactose-negative) for each test organism.
- Using sterile 5-ml pipettes, aseptically transfer 5 ml of the Lac<sup>+</sup> and Lac<sup>-</sup> inorganic synthetic broth cultures to their respectively labeled test tubes.
  - Using a sterile 1-ml pipette, aseptically add 0.5 ml of the glucose and lactose solutions and 0.5 ml of sterile distilled water to the appropriately labeled tubes.
  - Determine the absorbance of all cultures at a wavelength of 600 nm. Record your results in the Lab Report.
  - Aseptically transfer each culture to its appropriately labeled flask. (*Note: If side-arm flasks are available, additions and absorbance readings may be made directly.*)
  - Incubate all flasks for 2 hours in a shaking waterbath at 37°C and 100 strokes per minute.
  - Following incubation, transfer all cultures back to their appropriately labeled test tubes.
  - Determine and record in the Lab Report the absorbance for each culture at a wavelength of 600 nm. Based on your observations, indicate whether growth has occurred in each of the cultures.
  - To each culture, add 5 drops of toluene and shake vigorously (toluene ruptures the cells, releasing intact enzymes).
  - To each culture, add 5 drops of ONPG solution.
  - Incubate all cultures for 40 minutes at 37°C.
  - Following the addition of ONPG, observe the cultures for the presence of yellow coloration indicative of  $\beta$ -galactosidase synthesis and activity. In the Lab Report, record the colors of your cultures and the presence (+) or absence (-) of the  $\beta$ -galactosidase activity.

*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

Cultures	ABSORBANCE AT 600 nm		Growth (+) or (-)	Color of Culture with ONPG	$\beta$ -Galactosidase (+) or (-)
	Prior to Incubation	Following Incubation			
<b>Lac<sup>+</sup> <i>E. coli</i></b> Glucose					
Lactose					
Water					
<b>Lac<sup>-</sup> <i>E. coli</i></b> Glucose					
Lactose					
Water					

Explain the absence of growth in some of the cultures.

## Review Questions

1. Distinguish between constitutive enzymes and inducible enzymes.



## LEARNING OBJECTIVES

Once you have completed this experiment, you should be able to

1. Demonstrate genetic recombination in bacteria by the process of conjugation.


## Principle

Genetic variability is essential for the evolutionary success of all organisms. In diploid eukaryotes, the processes of **crossing over** (exchange of genetic material between homologous chromosomes) and **meiosis** contribute to this variability. In haploid, asexually reproducing prokaryotic organisms, genetic recombination may occur by **conjugation**, **transduction**, and **transformation**. In this experiment, only the process of conjugation is considered.

**Conjugation** is a mating process during which a unidirectional transfer of genetic material occurs at physical contact between two “sexually” differentiated cell types. This differentiation, or existence of different mating strains in some bacteria, is determined by the presence of a **fertility factor**, or **F factor**, within the cell. Cells that lack the F factor are recipients (females) of the genetic material during conjugation and are designated as  $F^-$ . Cells possessing the F factor have the ability to act as genetic donors (males) during mating. If this F factor is extrachromosomal (a **plasmid** or **episome**), the cells are designated as  $F^+$ ; most commonly only the F factor is transferred during conjugation. If this factor becomes incorporated into the bacterial chromosome, there is a transfer of chromosomal genes, although generally not involving the entire chromosome or the F factor. The resulting cells are designated **Hfr**, for **high-frequency recombinants**.

In this experiment, you will prepare a mixed culture representing a cross between an Hfr prototrophic (wild type) strain of *E. coli* that is

streptomycin-sensitive and an  $F^-$  auxotrophic (mutant) *E. coli* strain that requires threonine (thr), leucine (leu), and thiamine (thi) and is streptomycin-resistant (Str-r). Following a short incubation period, isolation of only the threonine and leucine recombinants will be performed by plating the mixed culture on a minimal medium containing streptomycin and thiamine. The streptomycin is incorporated into the medium to inhibit the growth of the wild-type, streptomycin-sensitive (Str-s) parental Hfr cells. The thiamine is required as an essential growth factor for the thiamine-negative ( $thi^-$ ) recombinant cells. Because of its distant location on the chromosome, this marker will not be transferred during the short mating period. A genetic map denotes the time in minutes required for the transfer of a marker (operon) from the donor cell to the recipient cell. **Figure 55.1** is the genetic map showing the site of origin of transfer and location of relevant markers in this experiment.

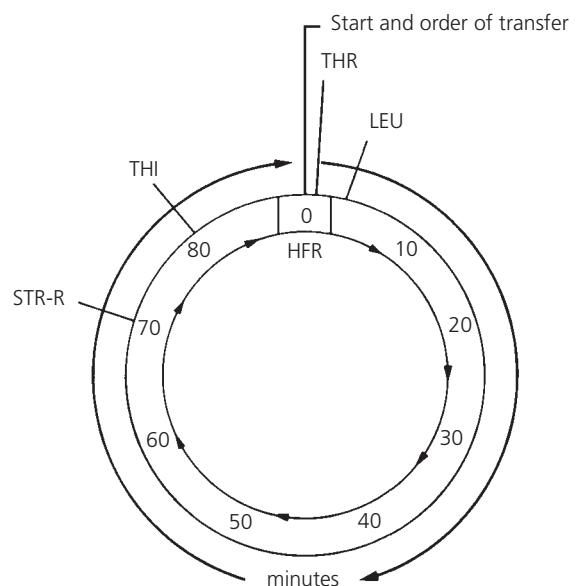

**Figure 55.1** Genetic map of *Escherichia coli*



## CLINICAL APPLICATION

### Antibiotic Resistance

Conjugation is a major cause of the spread of antibiotic resistance and represents a serious problem in antibiotic therapy of immunosuppressed patients. Bacteria that carry several resistant genes are called multi-drug-resistant superbugs. The indiscriminate use of antibiotics both within the health care profession and the illegal sale of drugs without prescriptions are mainly responsible for the increased spread of antibiotic resistance.

## AT THE BENCH



## Materials

### Cultures

12-hour nutrient broth cultures of  $F^-$  *E. coli* strain  $thr^-$ ,  $leu^-$ ,  $thi^-$ , and  $Str-r$  (ATCC e 23724); and Hfr *E. coli* strain  $Str-s$  (ATCC e 23740).

### Media

Per designated student group: three plates of minimal medium plus streptomycin and thiamine.

### Equipment

Bunsen burner, beaker with 95% ethyl alcohol, L-shaped bent glass rod, 1-ml sterile pipettes, mechanical pipetting device, sterile 13- × 100-mm test tube, waterbath shaker, and glassware marking pencil.

## Procedure Lab One

1. With separate sterile 1-ml pipettes, aseptically transfer 1 ml of the  $F^-$  *E. coli* culture and 0.3 ml of the Hfr *E. coli* culture into the sterile 13 × 100-mm test tube.

2. Mix by gently rotating the culture between the palms of your hands.
3. Incubate the culture for 30 minutes at 37°C in a waterbath shaker at the lowest speed setting.
4. Appropriately label two minimal plus streptomycin and thiamine agar plates, to be used for the control plates in Step 5.
5. To prepare control plates of the parental Hfr and  $F^-$  *E. coli* strains, aseptically add 0.1 ml of each *E. coli* strain to its appropriately labeled agar plate.
6. Use the spread-plate technique as shown in **Figure 55.2** and as instructed below.
  - a. Dip the bent glass rod into the beaker of 95% ethyl alcohol.
  - b. Sterilize the glass rod by flaming with a Bunsen burner.
  - c. Remove the glass rod from the Bunsen burner, allow flame to extinguish, and cool the glass rod.
  - d. Spread the inoculum over the agar surface by rotating the plate.
7. Following incubation of the mixed culture, vigorously agitate it to terminate the genetic transfer.
8. Appropriately label a minimal plus streptomycin and thiamine plate. Aseptically add 0.1 ml of the mixed culture. Spread the inoculum over the entire surface with a sterile glass rod.
9. Incubate all plates in an inverted position for 48 hours at 37°C.

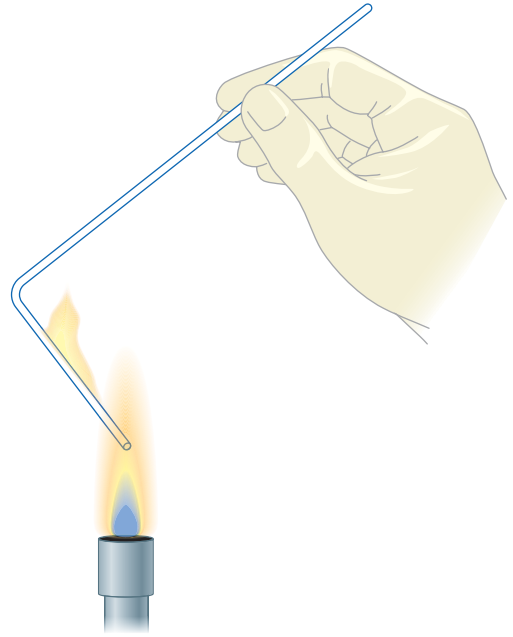
## Procedure Lab Two

1. Observe all plates for growth of colonies.
2. Record your observations in the Lab Report.


## PROCEDURE




**a** Dip the bent glass rod into the beaker of 95% ethyl alcohol.



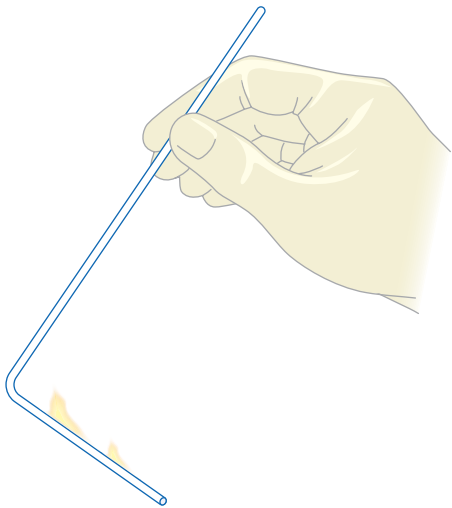**b** Sterilize the glass rod by flaming with a Bunsen burner.



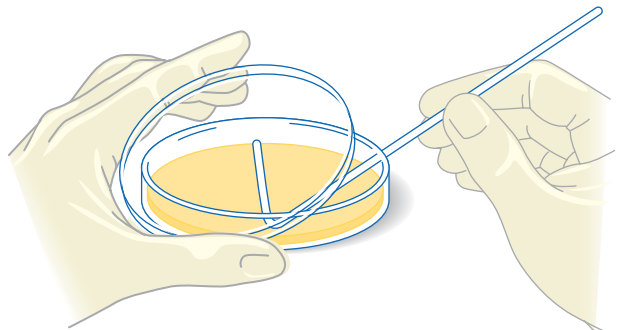**c** Remove from Bunsen burner, allow flame to extinguish, and cool the glass rod.



**d** Spread the inoculum over the agar surface by rotating the plate.

**Figure 55.2** Spread-plate technique

*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

1. Observe all plates for the presence (+) or absence (–) of colonies. Record your results in the chart.

	Hfr <i>E. coli</i> Plate	F <sup>–</sup> <i>E. coli</i> Plate	Mixed-Culture Plate
Growth (+) or (–)			


2. Do you expect any growth to be present on the two parental *E. coli* minimal agar plates? Explain.
3. Did genetic recombination occur? Explain how your observations support your answer.

## Review Questions

1. Explain how genetic variations may be introduced in eukaryotic and prokaryotic cells.

2. Explain the significance of the F factor.

3. Distinguish between F<sup>+</sup> and Hfr bacterial strains.

4.  Explain the importance of the streptomycin marker in the parental *E. coli* strains.

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be able to

1. Isolate a streptomycin-resistant mutant in a prototrophic bacterial population by means of the gradient-plate technique.


## Principle

**Mutation**, a change in the base sequence of a single gene, although infrequent, is one of the sources of genetic variability in cells. In some instances, these changes enable the cell to survive in an otherwise deleterious environment. An example of such a genetic adaptation is the development of **antibiotic resistance** in a small population of microorganisms prior to the advent and large-scale use of these agents. This microbial characteristic of antibiotic resistance is of major clinical importance because the number of drug-resistant microbial strains continues to increase. This occurs due to their extensive use and frequent misuse, which over the years have selected for the drug-resistant strains by their microbicidal effects on the sensitive cell forms. These agents select for the resistant mutant and do not act as inducers of the mutation.

In a drug-resistant organism, the mutated gene enables the cell to circumvent the antimicrobial effect of the drug by a variety of mechanisms, including the following:

1. The production of an enzyme that alters the chemical structure of the antibiotic, as in penicillin resistance.
2. A change in the selective permeability of the cell membrane, as in streptomycin resistance.
3. A decrease in the sensitivity of the organism's enzymes to inhibiting mechanisms, as in the resistance to streptomycin, which interferes with the translation process at the ribosomes.
4. An overproduction of a natural substrate (metabolite) to compete effectively with the drug (antimetabolite), as in the resistance to sulfonamides, which produce their antimicrobial effect by competitive inhibition.

The following procedure is designed to allow you to isolate a streptomycin-resistant mutant from a prototrophic (wild type, streptomycin-sensitive) *Escherichia coli* culture by means of the **gradient-plate technique**. This requires preparation of a double-layered agar plate as illustrated in **Figure 56.1**. The lower, slanted agar-medium layer lacks streptomycin. When poured over the lower slanted layer, the molten agar medium containing the antibiotic will produce a streptomycin concentration gradient in the surface layer. Following a spread-plate inoculation of the prototrophic test

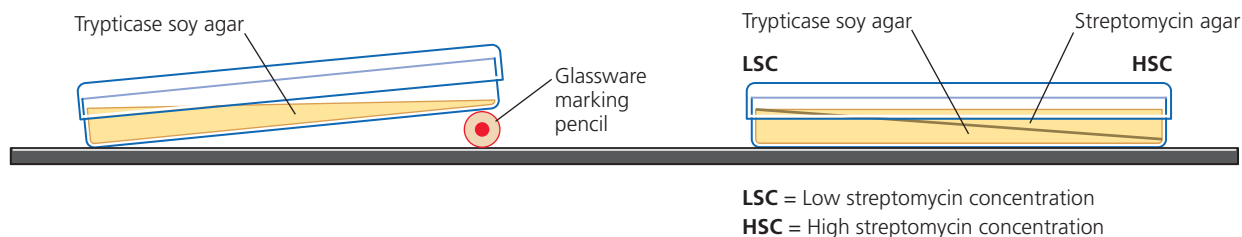

**Figure 56.1** Preparation of a streptomycin gradient plate

culture and incubation, the appearance of colonies in a region of high streptomycin concentration is indicative of streptomycin-resistant mutants.

## CLINICAL APPLICATION

### Searching for Resistance Mutations

Once strains of bacteria that are known to be resistant are isolated, microbiologists attempt to find the gene or genes responsible. A gene called *sasX*, found in methicillin-resistant *Staphylococcus aureus* (MRSA), has almost doubled in frequency over the past decade. It is located in a mobile genetic element that allows its easy transfer to other bacteria. This gene helps the bacterium to more effectively colonize nasal tissues and evade the host's immune system. Targeting this gene may provide a route for highly effective therapies.

## AT THE BENCH



## Materials

### Cultures

24-hour nutrient broth culture of *E. coli*.

### Media

Per designated student group: two 10-ml Trypticase soy agar deep tubes.

### Reagent

Stock streptomycin solution (10 mg per 100 ml of sterile distilled water).

### Equipment

Sterile Petri dish (100 × 15mm), sterile 1-ml pipettes, mechanical pipetting device, inoculating loop, bent glass rod, beaker with 70% ethanol, waterbath, and glassware marking pencil.

## Procedure Lab One

1. In a hot waterbath, melt two Trypticase soy agar tubes. Cool and maintain at 45°C.

2. Place a pencil under one end of a sterile Petri dish, pour in a sufficient amount of the molten agar medium to cover the entire bottom surface, and allow to solidify in the slanted position.
3. Using a sterile 1-ml pipette, add 0.1 ml of the streptomycin solution to a second tube of molten Trypticase soy agar. Mix by rotating the tube between the palms of your hands.
4. Place the dish in a horizontal position, pour in a sufficient amount of the molten agar medium containing streptomycin to cover the gradient agar layer, and allow to solidify.
5. With a sterile 1-ml pipette, add 0.2 ml of the *E. coli* test culture. With an alcohol-dipped and flamed bent glass rod, spread the culture over the entire agar surface as illustrated in Figure 55.2 on page 373.
6. Incubate the appropriately labeled culture in an inverted position for 48 hours at 37°C.

## Procedure Lab Two

1. Observe the plate for the appearance of discrete colonies and indicate their positions in the “Initial Incubation” diagram in the Lab Report.
2. Select one or two isolated colonies present in the middle of the streptomycin concentration gradient. With a sterile inoculating loop, streak the selected colonies toward the high-concentration end of the plate.
3. Incubate the plate in an inverted position for 48 hours at 37°C.

## Procedure Lab Three

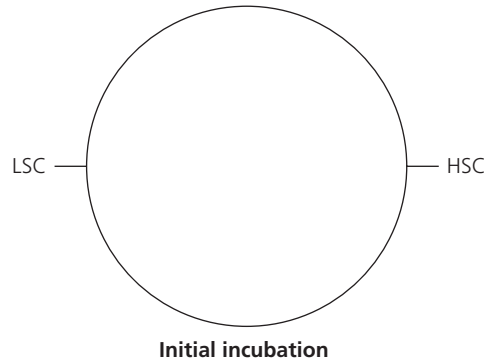
1. Observe the plate for a line of growth from the streaked colonies into the area of high streptomycin concentration. Growth in this area is indicative of streptomycin-resistant mutants.
2. Indicate the observed line(s) of growth in the “Second Incubation” diagram in the Lab Report.

Name: \_\_\_\_\_


Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

1. Indicate the positions of discrete colonies in the diagram below (LSC = low streptomycin concentration; HSC = high streptomycin concentration).



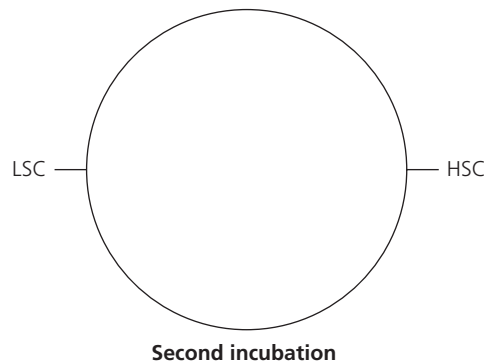2. Indicate the observed line(s) of growth in the diagram below.



## Review Questions

1. What mechanisms are responsible for antibiotic resistance?



2. Why is it necessary to use an antibiotic gradient-plate preparation for isolation of mutants?

3.  Why has there been an increase in drug-resistant bacterial strains in recent years?

4.  Does the streptomycin in the medium cause the mutations? Explain.

# The Ames Test: A Bacterial Test System for Chemical Carcinogenicity

EXPERIMENT

57

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be able to

1. Screen for potential chemical carcinogens using a bacterial test system.

## Principle

Our exposure to a wide variety of chemical compounds has increased markedly over the past decades. Oncological epidemiologists strongly suspect that the intrusion of these chemicals in the form of industrial pollutants, pesticides, food additives, hair dyes, cigarette smoke, and the like may play a significant role in the induction of malignant transformations in humans. From a genetic aspect there is strong evidence linking **carcinogenicity** to **mutagenicity**. Research indicates that approximately 90% of the chemicals proved to be carcinogens are mutagens; they cause cancer by inducing mutations in somatic cells. These mutations are most frequently a result of base substitutions, the substitution of one base for another in the DNA molecule, and frameshift mutations, a shift in the reading frame of a gene because of the addition or deletion of a base.

In view of the rapid advent of new products and new industrial processes with their resultant pollutants, it is essential to determine their potential genetic hazards. Despite the fact that mammalian cell structure and human enzymatic pathways differ from those in bacteria, the chemical nature of DNA is common to all organisms; this permits the use of bacterial test systems for the rapid detection of possible mutagens and therefore possible carcinogens.

The **Ames test** is a simple and inexpensive procedure that uses a bacterial test organism to screen for mutagens. The test organism is a histidine-negative ( $\text{his}^-$ ) and biotin-negative ( $\text{bio}^-$ ) auxotrophic strain of *Salmonella typhimurium*

that will not grow on a medium deficient in histidine unless a back mutation to  $\text{his}^+$  (histidine-positive) has occurred. It is recognized that the mutagenic effect of a chemical is frequently influenced by the enzymatic pathways of an organism, whereby nonmutagens are transformed into mutagens and vice versa when introduced into human systems. In mammals, this toxification or detoxification frequently occurs in the liver. The Ames test generally requires the addition of a liver homogenate, S-9, which serves as a source of activating enzymes, to make this bacterial system more comparable to a mammalian test system.

In the Ames test, by means of the spot method, molten agar containing the test organism, S-9 mix, and a trace of histidine to allow the bacteria to undergo the several cell divisions necessary for mutation to occur is poured on a minimal agar plate. A disc impregnated with the test chemical is then placed in the center of the test plate. Following diffusion of the test compound from the disc, a concentration gradient of the chemical is established. Following incubation, a qualitative indication of the mutagenicity of the test chemical can be determined by noting the number of colonies present on the plate. Each colony represents a  $\text{his}^- \rightarrow \text{his}^+$  revertant. A positive result, indicating mutagenicity, is obtained when an obvious increase in the number of colonies is evident compared to the number of spontaneous revertants on the negative control plate.

In the following procedure, you will perform a modified Ames test; you will not use the S-9 mix to test for the mutagenicity of nitro compounds, which, as in humans, are activated by the bacterial nitroreductases. Four minimal agar plates are inoculated with the *S. typhimurium* test organism. One plate, the negative control, is not exposed to a test chemical. Any colonies developing on this plate are representative of spontaneous  $\text{his}^- \rightarrow \text{his}^+$  mutations. The second plate, the positive control, is exposed to a known nitrocarcinogen, 2-nitrofluorene. The remaining two plates are used to determine the mutagenicity of two commercial hair dyes.

## CLINICAL APPLICATION

### Testing for Cancer-Causing Chemicals

The Ames test is a procedure used for the identification of mutagenic chemical and physical agents. The test was named after Bruce Ames, who invented the test in the 1970s. While the Ames test does not detect all mutagenic chemicals, it is used in the pharmaceutical industry to test drugs prior to use in clinical trials and also in the cosmetic industry to check on the mutagenic potential of makeup. A positive Ames test results in the rejection of the drug or agent for further development and testing.

## AT THE BENCH



## Materials

### Cultures

24-hour nutrient broth cultures of *S. typhimurium*, Strain TA 1538 (ATCC e 29631).

### Media

Per designated student group: four minimal agar plates and four 2-ml top agar tubes.

### Reagents

Sterile biotin-histidine solution, 2-nitrofluorene dissolved in alcohol, and two commercial hair dyes.

### Equipment

1-ml sterile pipettes, mechanical pipetting device, sterile discs, forceps, waterbath, Bunsen burner, and glassware marking pencil.

## Procedure Lab One

**⚠ Wear disposable gloves and a laboratory coat when handling 2-nitrofluorene. For disposal of this chemical, place excess into a sealable container and put it inside a fume hood for subsequent removal according to your institution's policy for disposal of hazardous materials.**


Refer to **Figure 57.1** as you read the following instructions.

1. Label three minimal agar plates with the name of the test chemical to be used. Label the fourth plate as a negative control.
2. Melt four tubes of top agar in a hot waterbath and maintain the molten agar at 45°C.
3. To each molten top agar tube, aseptically add 0.2 ml of the sterile biotin-histidine solution and 0.1 ml of the *S. typhimurium* test culture. Mix by rotating the test tube between the palms of your hands.
4. Aseptically pour the top agar cultures onto the minimal agar plates and allow to solidify.
5. Using sterile forceps, dip each disc into its respective test chemical solution and drain by touching the side of the container.
6. Place the chemical-impregnated discs in the center of the respectively labeled minimal agar plates. Place a sterile disc on the plate labeled negative control. With the sterile forceps, *gently* press down on the discs so that they adhere to the surface of the agar.
7. Incubate all plates in an inverted position for 24 hours at 37°C.

## Procedure Lab Two

1. Count the number of large colonies present on each plate and record on the chart in the Lab Report.
2. Determine and record the number of chemically induced mutations by subtracting the number of colonies on the negative control plate, representative of spontaneous mutations, from the number of colonies on each test plate.
3. Determine and record in the Lab Report the relative mutagenicity of the test compounds on the basis of the number of induced mutations: If below 10, (-); if greater than 10, (1+); if greater than 100, (2+); and if greater than 500, (3+).

## PROCEDURE

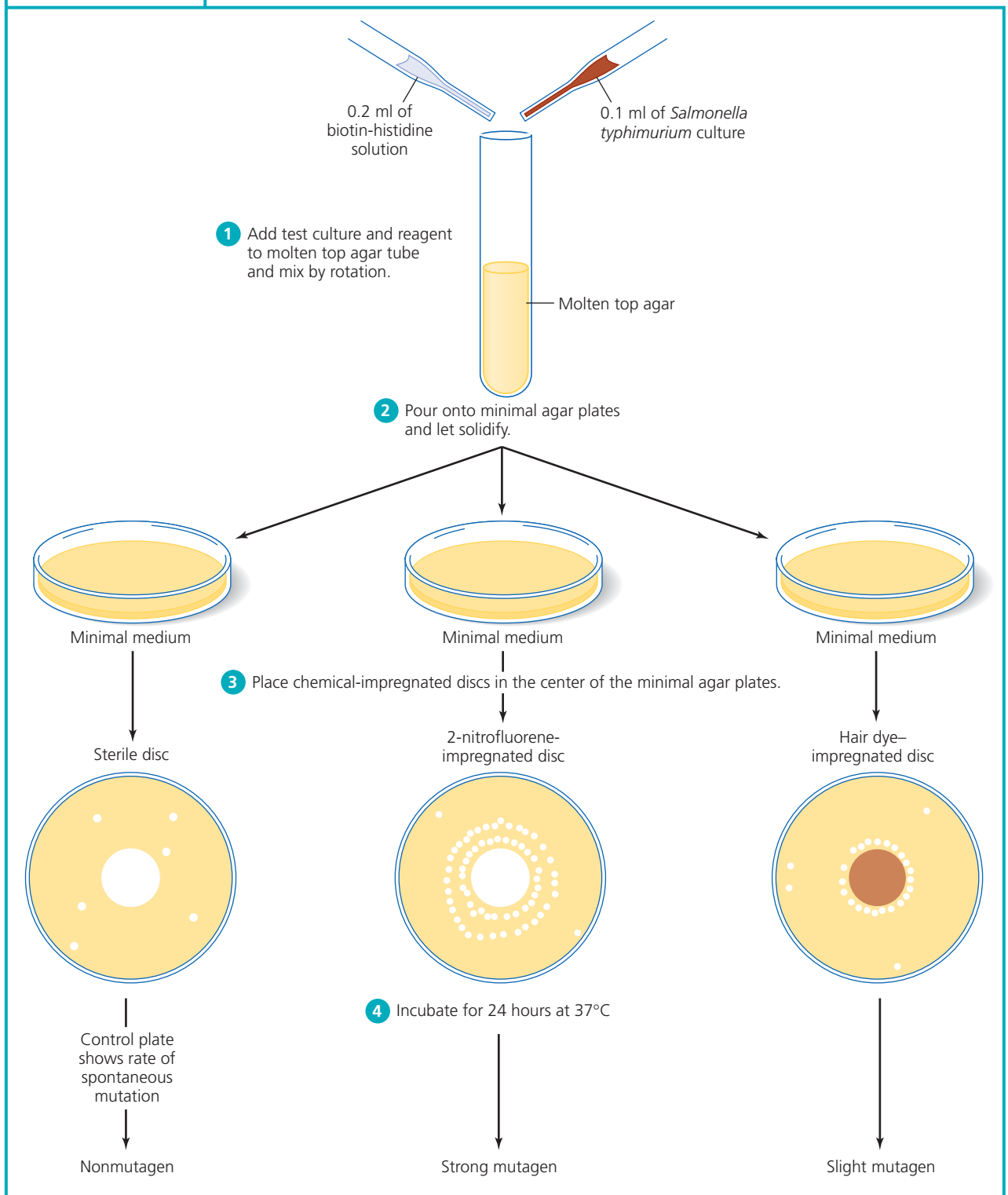

**Figure 57.1** The Ames test

*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_


## Observations and Results


Test Chemical	Number of Colonies	Number of Induced Mutations	Degree of Mutagenicity (-), (1+), (2+), or (3+)*
Negative control			
2-Nitrofluorene			
Hair dye 1			
Hair dye 2			


\* If below 10, (-); if greater than 10, (1+); if greater than 100, (2+); and if greater than 500, (3+).

## Review Questions

1. What is the purpose of the S-9 in the Ames test?
  
  
  
  
  
  
  
  
  
  
2. What is the purpose of the biotin-histidine solution in the Ames test?

3.  What is the relationship between chemical carcinogenicity and mutagenicity?

4.  What are the advantages of using bacterial systems instead of mammalian systems to test for chemical carcinogenicity?

5.  What are the disadvantages of using bacterial systems instead of mammalian systems to test for chemical carcinogenicity?

# Biotechnology

## LEARNING OBJECTIVES

*Once you have completed the experiments in this section, you should be familiar with the*

1. Process of bacterial transformation.
2. Isolation of bacterial plasmids.
3. Use of restriction endonucleases to cut (digest) DNA molecules into small fragments.
4. Use of agarose gel electrophoresis to separate DNA fragments and determine their size.

## Introduction

In the past, microbiologists focused on macromolecular cellular components. In recent years, however, this focus has rapidly shifted, and biological studies can now be easily carried out at the molecular level. **Molecular biology** is the study of life at the molecular level. This rapidly expanding field is based on the ability to manipulate the genes of an organism. As technology has expanded, so too have researchers' abilities to modify an organism's genetic structure by removing the genetic material from one organism and combining it with that of another organism. Commonly, this is referred to as **genetic engineering**, or **recombinant DNA technology**.

Most often, genomic manipulation begins with isolating plasmid DNA and foreign DNA and cutting them with the same restriction enzyme. The foreign DNA inserts itself into the plasmid and the recombinant plasmid is introduced into a bacterium. The bacteria are then grown and screened to select the desired gene. An example of the use of this technology is the exploration of human genetic structure to identify and ameliorate genetic diseases, particularly through prenatal diagnosis.

Scientists are also interested in the genetic engineering of hormones (e.g., pancreatic insulin), which are normally obtained only in their natural state. Genetically engineered organisms also help produce vaccines and antibiotic substances and remediate environmental toxins that threaten human and environmental health. The role of genetic engineering technology in clinical and forensic science is commonplace today. For example, DNA profiling is used to resolve paternity disputes, and it was used to identify the bones of Czar Nicholas, the last Russian czar, and his family. This technology is also used extensively in criminology as a tool to help establish guilt or innocence of individuals involved in criminal cases.

The many applications of recombinant DNA technology have brought about significant advances in medicine, pharmacology, basic research, industry, and agriculture. In medicine, advances in biotechnology now allow insulin-deficient diabetics to use different types of genetically engineered insulin (e.g., rapid acting, moderate acting, or slow acting) to meet diverse medical needs, rather than relying on the slower-acting insulin obtained from



the pancreas of slaughtered animals. Other proteins such as somatostatin, human growth hormone, interferon, interleukin-2—a regulator of the immune system—and blood clotting factor VIII have all been cloned using genetic engineering and are in clinical use. Biotechnology plays a prominent role in agriculture as well. Plants have been made resistant to herbicides such as glyphosate and pathogens such as the European corn borer. These few examples illustrate the amazing achievements made by the use of recombinant

DNA technology. More fascinating biotechnological advances are anticipated in the years ahead.

In the experiments that follow, you will explore some of the basic techniques in biotechnology. You will transform an antibiotic-susceptible cell into an antibiotic-resistant one by means of a plasmid. You will isolate plasmid DNA from bacterial cells and cut it into small fragments using restriction endonucleases. Finally, using agarose gel electrophoresis, you will separate these DNA fragments and determine their size.

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be able to

1. Transform a competent ampicillin-susceptible strain of *Escherichia coli* into one that is ampicillin-resistant by means of a DNA plasmid.
2. Visualize transformed cells using a color marker gene carried in the plasmid.
3. Calculate the efficiency of transformation.




## Principle


Transformation is a process whereby small pieces of host cell genomic DNA are able to enter a recipient cell and become incorporated into a homologous area on the recipient cell's genome. Historically, transformation had its origin in the pioneering experiments of Fred Griffith in London in the late 1920s. Working with *Streptococcus pneumoniae*, he noted that when an encapsulated smooth (S) strain that was lethal for mice was heat-killed and mixed with a living culture of an avirulent nonencapsulated rough (R) strain and then injected into mice, the result was fatal (**Figure 58.1**). Subsequent isolation of the organisms from the tissues of the dead mice revealed that the rough avirulent strain had been converted to a smooth, encapsulated, and lethal strain of *S. pneumoniae*.

This unusual experiment by Griffith, done long before DNA was determined to be the genetic basis for life, was simply termed by him

as a “transformation.” Today, in retrospect, we realize that this experiment proved to be the first indication of gene activity and the first demonstration of genetic recombination in bacteria. Later, Avery, McLeod, and McCarty, research scientists at the then Rockefeller Institute, were able to show that the transforming factor in Griffith's experiment was not a protein, as had been previously suspected, but a little-studied organic chemical called deoxyribonucleic acid (DNA).

During the **transformation** process, the donor cells forcibly lyse, releasing small segments of DNA containing 10 to 20 genes. These small segments have the ability to pass through the cell wall and cell membrane of a **competent cell** (a cell that is able to take up DNA from its environment). During naturally occurring transformations, a double-stranded DNA segment passes through the cell wall and into the cell's cytoplasm, and if there is sufficient sequence similarity, the foreign DNA undergoes homologous recombination with the recipient chromosome. The genome of the recipient cell has now been modified to contain DNA with genetic characteristics of the donor cell. Naturally occurring transformations are of great interest medically because they may serve as a vehicle for genetic exchange among pathogenic organisms. Interestingly, it appears that a larger percentage of pathogenic bacteria, such as *Streptococcus pneumoniae*, *Neisseria gonorrhoeae*, and *Haemophilus influenzae*, is capable of natural transformation than the nonpathogenic bacteria. This raises the intriguing possibility that the exchange of genetic material allows pathogenic cells to acquire the ability to evade a host's bodily defenses.

-  Encapsulated smooth strain (S cell)
-  Nonencapsulated rough strain (R cell)
-  Heat-killed S cell

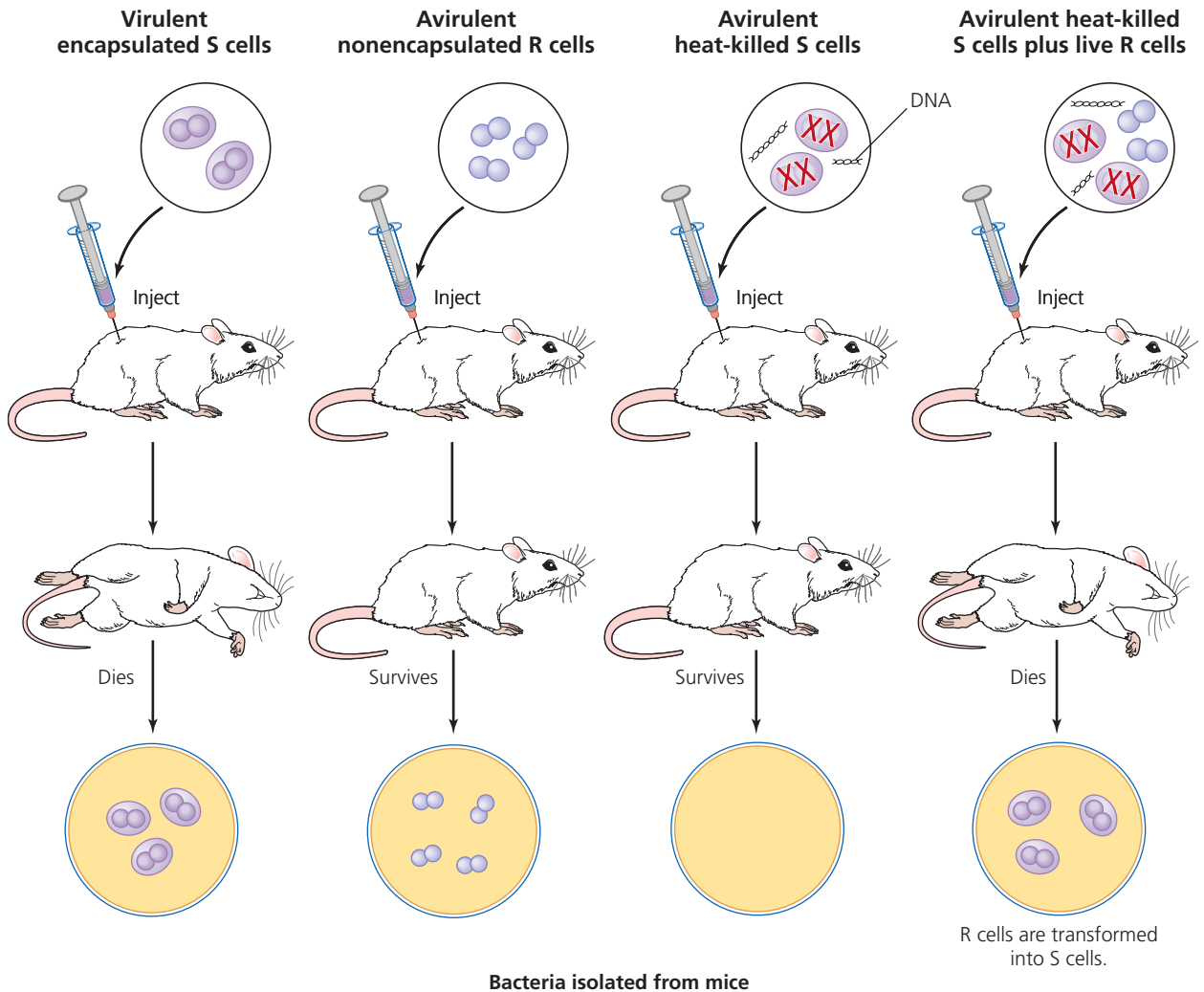

**Figure 58.1 Griffith's experiment on transformation.** An avirulent rough strain of *Streptococcus pneumoniae* becomes lethal when transformed after mixing with the DNA of a heat-killed virulent smooth strain of *S. pneumoniae*.

Not all bacteria are naturally transformable, however, and methods have been developed to produce competency in various types of cells and transform those cells artificially. This process was initiated in the 1970s, when it was shown that treating a recipient cell with a cold calcium chloride ( $\text{CaCl}_2$ ) solution allows the passage of donor DNA into the cell. The porosity of the cell

wall is already almost sufficient to allow the passage of intact DNA; it is the cell membrane that is the true barrier, and its permeability is altered by this drastic treatment with  $\text{CaCl}_2$ , allowing DNA to pass through the membrane and into the cell. With our rapidly advancing knowledge in the field of molecular genetics, it is now possible to artificially induce transformations by the use of **plasmids**.


Plasmids are small, circular pieces of extrachromosomal DNA with a length of 5,000 to 100,000 base pairs (bp), capable of autonomous replication in the bacterial cytoplasm. Another membrane-altering method is **electroporation**. In this method, cells are suspended in a DNA solution and subjected to high-voltage electric impulses that destabilize the cell membrane, resulting in increased permeability and enabling DNA to pass into the cells. **Transduction** is a method of horizontal passage of genetic material from one bacterial cell to another by means of a bacteriophage. **Conjugation** occurs when bacterial DNA is transferred from one cell to another via the formation of a protoplasmic bridge, called a conjugative, or sex, pilus.

In the following experiment, the transformation of *E. coli* is artificial because the cells must be treated with a salt concentration and temperature shocks, an environment that is not natural for these cells. You will use a competent strain of *E. coli* in which the *lac* operon has been deleted, leaving it devoid of a  $\beta$ -galactosidase gene. The procedure directs you to suspend the *E. coli* in cold  $\text{CaCl}_2$  and then introduce an ampicillin-resistant plasmid (pBLU<sup>®</sup>). The plasmid confers ampicillin resistance because it carries the  $\beta$ -lactamase gene  $\text{amp}^r$ , as shown in **Figure 58.2**. This suspension is incubated in ice and then heated. The cold shock and heat shock in the presence of  $\text{CaCl}_2$  alters the permeability of the outer surfaces of the cell and facilitates the passage of the DNA into the cellular cytoplasm.

## CLINICAL APPLICATION

### Transferring Genes Between Bacteria

Bacterial transformation is the exchange of genetic material (DNA) between strains of bacteria by the transfer of a fragment of naked DNA from a donor cell to a recipient cell, followed by recombination in the recipient chromosome. Called “horizontal gene transfer,” this allows for the transfer of genes between bacteria of the same or different species. Human and other genes are routinely transferred into bacteria in order to synthesize products for medical and commercial use, such as human insulin, human growth hormone, and vaccines.

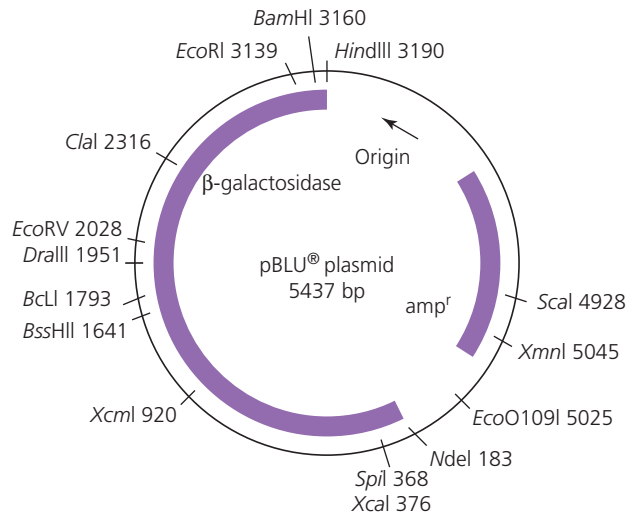

**Figure 58.2** Plasmid map of pBLU

## AT THE BENCH



## Materials

### Cultures

18- to 24-hour Luria-Bertani (LB) agar base streak plate cultures of *Escherichia coli*.

### Media

Per designated student group: two LB agar base plates, three LB agar base plates plus ampicillin (Amp), three LB agar base plates plus ampicillin/X-Gal, and one tube of LB broth. *Note: X-Gal is a colorimetric analog of lactose that is cleaved by  $\beta$ -galactosidase to yield blue-colored colonies.*

### Reagent


50 mM  $\text{CaCl}_2$  solution.

**Plasmid** pBLU is 5437 bp long and has the gene for  $\beta$ -lactamase (penicillinase) and the gene for  $\beta$ -galactosidase.

### Equipment

Sterile plastic 13- × 100-mm test tubes or plastic 1.5-ml centrifuge tubes, adjustable micropipette (0.5 to 100  $\mu\text{l}$ ) with sterile plastic micropipette

## PROCEDURE

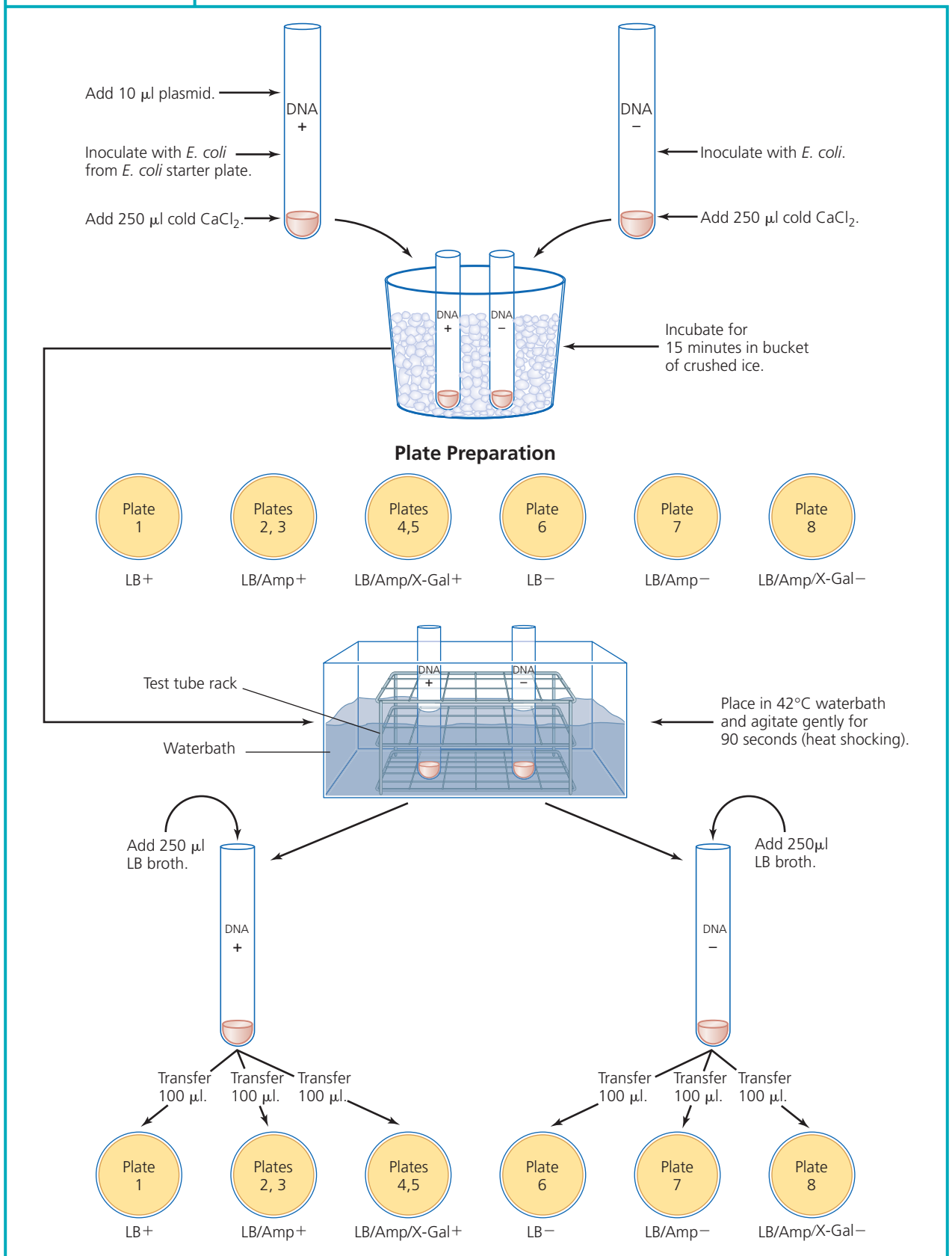

**Figure 58.3** Transformation procedure


tips (10 to 100  $\mu\text{l}$ ) (or 1.0-ml graduated, individually wrapped, disposable plastic transfer pipettes), glass beads (6-mm diameter), glassware marking pencil, disposable plastic inoculating loops (standard wire loops may be used), Bunsen burner, waterbath, 500-ml beaker of crushed ice, 500-ml beaker labeled “waste,” 500-ml beaker containing disinfectant solution, and a Quebec colony counter (or permanent marker); if the spread-plate method is used, bent glass rod, beaker of alcohol, and turntable.

## Procedure Lab One

Refer to **Figure 58.3** as you read the numbered instructions.

*Note: If using a plastic transfer pipette (Figure 58.4), it is essential that you calibrate it to deliver a volume of 100  $\mu\text{l}$  (0.1 ml), required for plasmid transfer. Once calibrated, it should be marked with a permanent glassware marker and retained to be used as a guide in the transformation experiment. Ask your instructor for help if needed.*

- With a glassware marking pencil, label two 13-  $\times$  100-mm test tubes, one as “DNA+” and the other as “DNA-.” The DNA+ tube will receive the plasmid.
- Using a sterile pipette, transfer 250  $\mu\text{l}$  (0.25 ml) of ice-cold  $\text{CaCl}_2$  solution into each tube.
- Place both tubes in a 500-ml beaker of crushed ice.
- Using a sterile inoculation loop, obtain a large mass of cells approximately 5 mm in size (about the size of a pencil eraser) from the *E. coli* starter plate, and inoculate the tube labeled “DNA+.” *Note: Be sure to immerse the loop directly into the  $\text{CaCl}_2$  and shake the loop vigorously to dislodge the inoculum. Discard the plastic loops in the beaker labeled “waste” or sterilize the wire loop by flaming it.*
- Disperse the cells by gently tapping the tube with your finger until a uniform milky-white translucent cell suspension is obtained.
- Repeat Steps 4 and 5 to inoculate the tube marked “DNA-,” using an equal amount of inoculum and a sterile inoculating loop.
- Using a sterile pipette, deliver 10  $\mu\text{l}$  (0.01 ml) of pBLU plasmid into the DNA+ tube. Tap the tube several times with your finger to ensure complete mixing of the plasmid and cell suspension. *Note: Discard the plastic tip or*
- disposable pipette into a beaker containing disinfectant solution.*
- Return the DNA+ tube to the crushed-ice beaker and incubate for 15 minutes. During this time, label your agar plates as described in Step 9.
- Label the 8 LB agar plates as follows:  
Plate 1: LB+  
Plates 2, 3: LB/Amp+  
Plates 4, 5: LB/Amp/X-Gal+  
Plate 6: LB-  
Plate 7: LB/Amp-  
Plate 8: LB/Amp/X-Gal-
- Remove both tubes from ice after 15 minutes; place them in a test tube rack and immediately into a 42°C waterbath with *gentle agitation* for 90 seconds (heat shocking).
- Return both tubes to the crushed-ice beaker for 1 minute.
- With a sterile pipette, add 250  $\mu\text{l}$  (0.25 ml) of LB broth to both the DNA+ and DNA- tubes. Tap the tubes with your finger to achieve uniform cell suspension. (These are the transformation tubes.)
- Incubate both tubes in a test tube rack for 10 minutes. *Note: This is the **recovery period**, when the cells convert their newly modified genotype into a functionally ampicillin-resistant phenotype.*

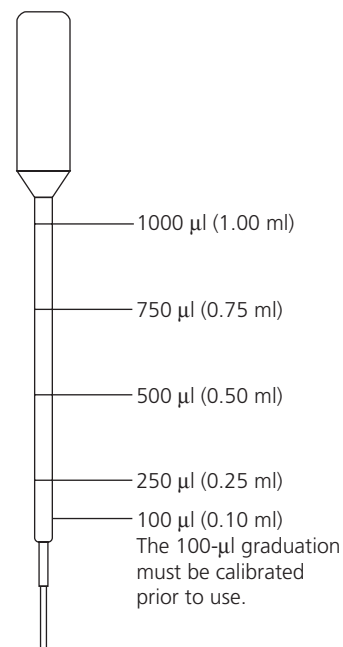

**Figure 58.4** Plastic transfer pipette

14. Using a new plastic micropipette tip or disposable plastic transfer pipette for each inoculation, inoculate 100  $\mu$ l (0.1 ml) of cells from the DNA+ transformation tube onto the surface of LB plates 1–5, and inoculate 100  $\mu$ l (0.1 ml) of cells from the DNA– transformation tube onto plates 6–8.
15. Place six sterile 6-mm glass beads on the surface of each inoculated plate. Replace the cover and spread the cell suspension by gently moving the plate up and down and then side to side a few times. *Note: Do not swirl or rotate the plate.* If the spread-plate technique is used, refer to Experiment 3 (page 20) for the correct procedure.
16. Repeat Steps 14 and 15 for the remaining plates.
17. Allow the plates to set for a few minutes so that the inoculum may be absorbed by the agar.
18. Remove the glass beads from the plate by lifting the cover slightly while holding the plate vertically over the beaker of disinfectant, allowing the beads to leave the plate. *Note: This step may be eliminated if the spread-plate procedure is used.*
19. Incubate all plates at 37°C for 24 to 36 hours or at room temperature for 48 to 72 hours.

## Procedure Lab Two

1. In the Lab Report, predict whether each plate will experience growth or no growth. Use a plus (+) sign for growth and a minus (–) for no growth.
2. Without removing the cover of the Petri plates, observe the colonies through the bottom of each plate.
3. Perform a colony count on each plate using an electronic colony counter if available, or use a permanent marker to mark each colony on the bottom of the plate as it is counted. Plates with more than 300 colonies should be designated as **TNTC (too numerous to count)**; plates with fewer than 30 colonies are designated as **TFTC (too few to count)**. Record your results in the Lab Report.
4. For each plate, did transformation occur? Record your results in the Lab Report.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

Plate Number	Designation	Growth + or -	Transformation Yes or No	Number of Colonies
1	LB+			
2	LB/Amp+			
3	LB/Amp+			
4	LB/Amp/X-Gal+			
5	LB/Amp/X-Gal+			
6	LB-			
7	LB/Amp-			
8	LB/Amp/X-Gal-			

Determine the sum of the colonies on the two experimental plates for LB/AMP+ and the two LB/AMP/X-Gal+ plates. Determine the averages and use these figures to calculate the transformation efficiency, using the following protocol.

Avg. LB/Amp+ colonies \_\_\_\_\_ Avg. LB/Amp/X-Gal+ colonies \_\_\_\_\_

- a.** The concentration of the plasmid used in this experiment was  $0.005 \mu\text{g}/\mu\text{l}$ . Keep in mind that you used  $10 \mu\text{l}$ .

Total mass of plasmid = volume  $\times$  concentration

Total mass of plasmid = \_\_\_\_\_

- b.** Calculate the total volume of the cell suspension.

Total volume = volume of  $\text{CaCl}_2$  solution + volume of plasmid + volume of LB broth

Total volume = \_\_\_\_\_

- c.** Calculate the fraction of cell suspension spread on each plate.

Fraction spread = volume of suspension spread on plate / total volume

Fraction spread = \_\_\_\_\_



- d. Calculate the mass of DNA plasmid in the cell suspension spread on each plate.

Mass of DNA plasmid spread = total mass of plasmid  $\times$  fraction spread

Mass of DNA plasmid spread = \_\_\_\_\_

- e. Calculate the transformation efficiency (the number of colonies per  $\mu\text{g}$  of plasmid DNA).

Transformation efficiency = average of colonies counted/mass of plasmid spread

Express your results using scientific notation, as found in Appendix 1.

TRANSFORMATION EFFICIENCY	
LB/Amp+	LB/Amp/X-Gal+


## Review Questions

- Why do you think that it is necessary to autoclave transformed cells immediately after the termination of the experiment?
- How would you explain to an untrained neighbor the process scientists use to change the genes (genotype) of an organism to give it useful new traits?
- At a neighboring lab bench, students find that there are a few white colonies scattered among a few hundred blue colonies on their LB/Amp/X-Gal+ plates.
  - How would you explain this result?
  - To confirm your hypothesis, you ask to see their LB/Amp/X-Gal- plate. What do you see that supports your hypothesis?

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be able to

1. Isolate plasmids from plasmid-bearing bacteria.
2. Separate plasmids using agarose gel electrophoresis.
3. Compare electrophoretic mobilities of plasmids.

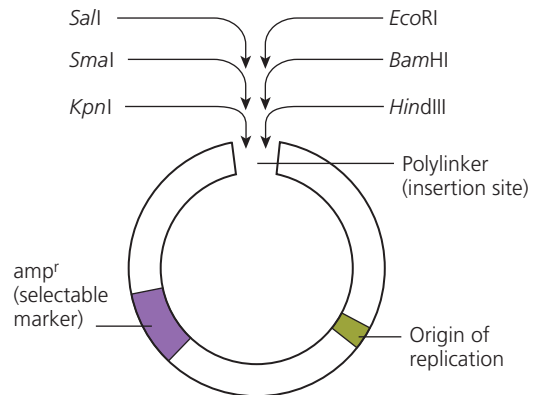

**Figure 59.1** Illustration of a plasmid, showing a selectable marker and a multiple cloning site

## Principle


Isolating pure DNA is a necessary step in studies that incorporate cloning, gene sequencing, gene mapping, or any other recombinant DNA technique. Many microorganisms contain small pieces of circular DNA called **plasmids** that exist separately from the host-cell genome. In studies that use recombinant DNA techniques, plasmid DNA is often preferred over chromosomal DNA because it is smaller and easier to manipulate. Plasmids that are commonly used as cloning vectors possess three important elements: an origin of replication, which allows the plasmid to be replicated independently of the cell's chromosome; a selectable marker, so the presence of the plasmid in the cell can be detected; and a cloning site into which a gene can be inserted. Although a selectable marker is not a required element of a plasmid, it can be useful in order to signal that the plasmid has been incorporated into the host cell. Plasmid genes that code for resistance to antibiotics are able to confer this resistance to the host cell, and a test of the host cell reveals that the transfer has occurred. In fact, antibiotic-resistance genes are among the most commonly used selectable markers.

Because plasmids are circular and capable of self-replication, they are able to serve as vectors for transportation of cloned fragments of DNA

into other cells for genetic engineering purposes. To do this, plasmids must have a multiple cloning site, or **polylinker**, which is a DNA segment with several unique insertion sites for restriction endonucleases located next to each other as shown in **Figure 59.1**.

Gel electrophoresis is a technique used to separate different sizes of DNA fragments from a sample of DNA. Because DNA is negatively charged, when a sample is loaded into a porous agarose gel and subjected to an electric current with a positive charge at one end of the gel and a negative charge at the other, DNA fragments will migrate through the pores in the gel, toward the positively charged end. Different-sized fragments of linear and circular DNA move through the gel at different speeds, thus traveling different distances in the gel over a set time period. Larger, longer pieces snake their way through the pores more slowly, while shorter, smaller pieces move more quickly and travel further toward the positively charged end of the gel.

There are two types of circular DNA: closed and nicked. Closed circular DNA has all of its nucleotides linked with phosphodiester bonds and is supercoiled. Nicked circular DNA has at least one broken phosphodiester linkage. Nicked DNA is sometimes referred to as “relaxed” because

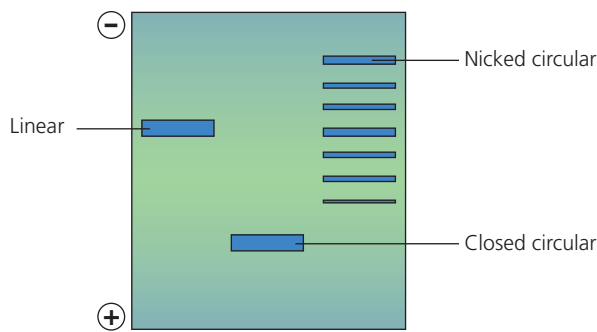

**Figure 59.2** Gel migration pattern for linear, nicked circular, and closed circular DNA

some of the tension present in covalently coiled and twisted DNA has been released. **Figure 59.2** illustrates the relative distance that each type of DNA described previously (linear, closed circular, and nicked circular) travels in an electrophoresed agarose gel.

Isolating plasmids is a multistep process, which involves rupturing a plasmid-bearing bacterium, using a variety of reagents to remove cellular components, and suspending plasmid DNA in an aqueous solution. After a plasmid-bearing organism is cultured, cells are lysed using alkali to release the plasmid DNA. The cellular debris is then precipitated by using a detergent and potassium acetate. Following centrifugation, the pellet that forms is removed, and alcohol is added to the supernatant to precipitate the DNA. The DNA precipitate is resuspended in Tris-EDTA buffer.

In the following experiment, two different strains of plasmid-bearing *Escherichia coli* will be used—*E. coli*-1 and *E. coli*-2. Half of the class will isolate *E. coli*-1 plasmid DNA, and the other half will isolate *E. coli*-2 plasmid DNA.

## CLINICAL APPLICATION

### Plasmids and Genetic Engineering

Plasmids are mostly found in bacteria and are used in recombinant DNA research to transfer genes between cells. Plasmids that confer antibiotic resistance (R Plasmids) have been of special interest because of their medical importance, and also because of their significant role in genetic engineering.

## AT THE BENCH



## Materials

### Cultures

24-hour Luria-Bertani broth, 50  $\mu\text{g/ml}$  of ampicillin, cultures of plasmid-bearing *Escherichia coli* ATCC 39991 (plasmid designation pIVEV) and *Escherichia coli* ATCC 53100 (plasmid designation pDGR-2).

### Reagents

Glucose-Tris-ethylenediaminetetraacetic acid (EDTA) buffer, Tris-EDTA buffer, Tris-acetate-EDTA buffer, 5M potassium acetate (KOAc), sodium hydroxide containing 1% sodium dodecyl-sulfate (NaOH/SDS), 95% ethanol at 0°C, 70% ethanol, molten agarose at 55°C, gel electrophoresis running dye, Carolina Blu stain or 0.025% methylene blue stain, and *Hind*III-cut bacteriophage lambda ( $\lambda$ ) DNA (used as the standard for comparing fragment sizes).

*Note: The formulation for some of these buffers and reagents may be found in Appendix 4. Although it is cheaper to prepare your own solutions, your instructor may have ordered a kit containing premade solutions.*

### Equipment

Microcentrifuge, 2-ml microcentrifuge tubes, digital micropipette—10, 100, and 200  $\mu\text{l}$ , small and large micropipette tips, waterbath, rubber micro test tube racks capable of floating, glassware marking pencil, ice bucket, crushed ice, a light box or overhead projector, millimeter ruler, agarose gel casting tray, staining tray, plastic sandwich size bags, and an electrophoretic apparatus.

## Procedure

### Using a Micropipette

Before the start of the experiment, familiarize yourself with the use of a micropipette, the function of which is to accurately deliver microliter volumes of solution. Not all micropipettes work the same way. Some are designed to deliver a fixed

volume, while others can deliver variable volumes. Your instructor will demonstrate the proper handling and use of these expensive instruments.

Using samples of colored water, practice using a micropipette, attaching different-sized micropipette tips and delivering various sample volumes to digestion tubes.

1. Set the scale on the pipette to the volume you wish to deliver.
2. Place a tip on the micropipette by pushing it firmly onto the pipette.
3. Depress the plunger to the first stop. This is necessary to remove all of the air from the tip.
4. To load the pipette, dip the pipette tip into the solution and release the plunger slowly to draw up the solution.
5. Touch the end of the tip to the side of the tube to remove any excess solution.
6. To deliver the solution, touch the side of the micropipette tip to the inside of the tube receiving the solution to produce a capillary.
7. Depress the plunger to the *first stop* and then continue depressing the plunger to the *second stop* to deliver the full volume of sample, blowing out the last bit in the tip.
8. Continue depressing the plunger while you remove the pipette tip from the tube. *Note: Releasing the plunger before removing the tip of the pipette from the tube will cause fluid to suck back into the tip.*

Before the laboratory session, *E. coli*-1 and *E. coli*-2 were inoculated in their medium and grown overnight. Before the start of the lab, your instructor dispensed 1 ml of culture into a microcentrifuge tube and spun it for 1 minute in a centrifuge. The supernatant was discarded, and the pellet retained. Another 1 ml of culture was added to the tube, and the process was repeated.

## Isolating the Plasmid

Obtain a microcentrifuge tube from your instructor with a retained pellet labeled “EC-1” or “EC-2.” With a glassware marking pencil, label the tube with your group name or number to identify it later. Refer to **Figure 59.3** as you complete the following steps.

1. Add 100  $\mu\text{l}$  of GTE (glucose, Tris, and EDTA) buffer to your tube and resuspend the pelleted cells by tapping with your finger or mixing by vortex. *Note: The EDTA in the buffer chelates the divalent metal ions,  $\text{Ca}^{2+}$  and  $\text{Mg}^{2+}$ , which destabilizes the cell membrane and*

## PROCEDURE

Prior to the lab, your instructor centrifuged 1 ml of *E. coli* culture for 1 minute. The supernatant was discarded and the pellet retained. Another 1 ml of *E. coli* culture was added to the pellet. It was centrifuged for 1 minute. The supernatant was discarded and the pellet retained for this experiment.

**Figure 59.3** Procedure for isolating bacterial plasmid DNA

*inhibits the activity of DNases. The glucose maintains the osmolarity, preventing the buffer from bursting the cell.*

2. Add 200  $\mu$ l of NaOH/SDS solution and mix gently by inversion about four to five times. Incubate the tube in an ice bucket for 5 minutes. *Note: This is a highly alkaline solution that lyses the cell, releasing the cytoplasm into the buffer, and separates the chromosomal DNA into single strands (ssDNA) and complexes with cellular proteins.*
3. Remove the tube from the ice bucket. Then add 500  $\mu$ l of potassium acetate (KOAc) and mix thoroughly by gentle inversion. *Note: The KOAc promotes the precipitation of chromosomal ssDNA and large RNA molecules, which are insoluble in this salt.*
4. Reincubate the tube in the ice bucket for another 5 minutes.
5. Remove the tube from the ice bucket and centrifuge for 5 minutes. Be sure the tubes are balanced in the centrifuge. *Note: In this step, pellets form from all of the cellular debris and organic molecules precipitated in the previous steps.*
6. Carefully decant the supernatant solution into a new microcentrifuge tube. *Note: The plasmid remains in the supernatant solution. The pellet and the tube are discarded.*
7. Add 1 ml of 95% ethanol at 0°C to the supernatant fluid in the new tube. *Note: The ethanol precipitates the plasmid.*
8. Incubate the plasmid in the ice bucket for 15 minutes.
9. Centrifuge the tube for 15 minutes to make the precipitated plasmid form a pellet.
10. Decant and discard the supernatant. *Note: Care must be taken not to shake the tube before or after decanting the supernatant. Do not be concerned if you do not see a pellet. It is there, provided that you were careful during the decanting step.*
11. Add 500  $\mu$ l of cold 70% ethanol to the pellet and gently tap the tube with your finger or rock the tube back and forth. *Note: This step washes the plasmid by removing the excess salt. The plasmid is insoluble in ethanol.*
12. Centrifuge the tube for 5 minutes.
13. Decant and discard the supernatant fluid.
14. Allow the pellet to dry for about 15 minutes until you no longer smell alcohol.


15. Add 100  $\mu$ l of TAE (Tris-acetate-EDTA) buffer to resuspend the pellet.
16. The plasmid may be placed in the freezer until the next lab class, or you may proceed to the electrophoresis step. *Note: If the electrophoresis is to be done during this class period, practice loading and casting the gel, which is described next.*

## Casting the Agarose Gel

*Note: Not all casting trays are the same. Your instructor will indicate which type will be used and whether there are special considerations during the setup.*

Your instructor prepared the 0.8% agarose gel in a 1X TAE buffer solution before class and maintained it at 55°C in a waterbath. It is ready to pour. One or two drops of Carolina Blu stain was added to the agarose buffer solution to give a small tinge of blue to the gel. At this concentration, the pores that form the gel lattice are such that they allow the free migration of the cut DNA fragments between 0.5 and 10 kilobases (kb).

1. Place the casting tray inside the casting tray box on a level surface.
2. Close off the ends of the tray with the rubber dams by tightening the knob on the top of the casting tray box.
3. Place a well-forming comb in the first notch at the end of the casting tray.
4. Pour 60–70 ml of agarose solution that has been cooled to 55°C into the tray. Use a toothpick or applicator stick to move the bubbles to the edge of the gel before it solidifies.
5. Allow the gel to solidify completely. It should be firm to the touch after 20 minutes.
6. Slowly remove the rubber dams and *very gently* remove the well-forming comb by pulling it straight up. *Note: Use extreme care not to damage or tear the wells.*
7. Place the gel on the platform in the electrophoresis box so that the formed wells are properly oriented toward the anode (negative pole with black cord). Because DNA is negatively charged, the cut DNA fragments will migrate to the cathode (positive pole with the red cord). Refer to **Figure 59.4** to see the proper setup of an electrophoretic apparatus.
8. Fill the electrophoresis box with TAE buffer to a level that just covers the gel, about 2 mm. Make sure that all of the wells are filled with the buffer.

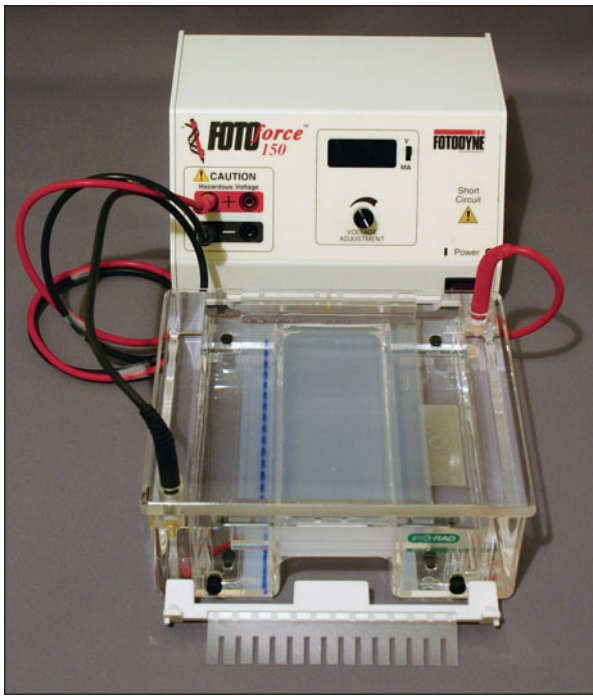

**Figure 59.4** Setup of agarose gel unit for DNA electrophoresis


## Practicing Gel Loading

Before loading your sample into the wells of the agarose gel, practice this challenging technique. Your instructor will demonstrate the proper method for loading the wells. Each student should practice on a gel that has been prepared earlier by your instructor, not on the gel to be used for running the samples.

1. Load the pipette with 22  $\mu\text{l}$  of loading gel.
2. Hold the pipette with both hands and dip the tip slightly through the buffer covering the gel, with the tip barely in the well.
3. Slowly discharge the contents of the pipette.  
*Note: The loading gel contains sucrose, which is heavier than the DNA and will pull the sample into the well.*
4. Practice the technique until you are comfortable with it.

## Electrophoresing of the Plasmids

1. Add 18  $\mu\text{l}$  of plasmid in the 1X TAE buffer to a new microcentrifuge tube. Then add 4  $\mu\text{l}$  of gel electrophoresis running dye to the tube.
2. Add 18  $\mu\text{l}$  of *Hind*III-cut lambda ( $\lambda$ ) DNA and 4  $\mu\text{l}$  of the gel electrophoresis running dye. With your glass marking pencil, label this tube

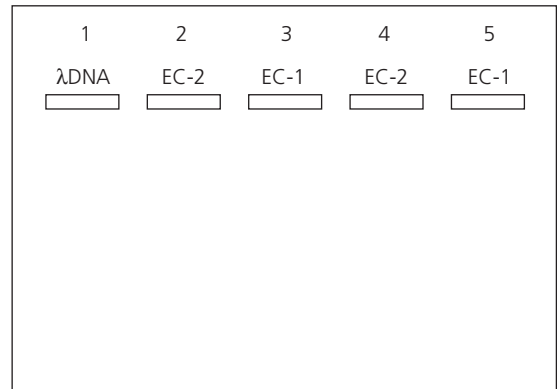

**Figure 59.5** Gel loading scheme

*Hind*III. This DNA has been cut into six linear fragments with the *Hind*III restriction enzyme. The fragments (bands) are various sizes: 23 kb, 9.4 kb, 6.6 kb, 4.4 kb, 2.3 kb, and 2.0 kb.

3. When the wells are ready to be loaded, make a diagram so that you will know the position of your sample in the agarose gel.
4. Fill the wells by designating the EC-1 samples as odd-numbered groups and EC-2 samples as even-numbered groups, as shown in **Figure 59.5**.
5. After the samples are loaded into the wells, place the lid on the electrophoresis gel box. Check that the power switch is turned to the “off” position and then attach the electrical leads (red to red and black to black) from the power supply to the box.
6. Turn the power pack on and adjust the rheostat dial to 110V.
7. Electrophorese the gel for 30 to 40 minutes or until the leading edge of the bromphenol blue dye (the dye in the loading gel) has traveled roughly three-fourths of the distance to the edge of the gel.
8. Turn the rheostat to zero and turn off the power. Disconnect the leads and remove the cover from the gel box.

## Staining the Gel

1. Put on a pair of disposable laboratory gloves.
2. Lift the gel tray out of the electrophoresis box, and slide the gel into a staining tray containing approximately 100 ml of Carolina Blu stain or 0.025% of methylene blue stain.
3. Allow the gel to remain in the stain for 30 to 40 minutes.



**Figure 59.6** Agarose gel stained with methylene blue following the electrophoresis of plasmid DNA

4. Pour off the stain into a waste beaker. Transfer the gel to a staining tray containing 100 ml of distilled water and allow the gel to decolorize (destain) for another 30 minutes. Frequent destaining with fresh distilled water for 2 minutes increases the intensity of the bands. For best results, let the gel destain overnight in a small volume of water. *Note: If the gel is left overnight in a large volume of water, it may destain too much.*
5. Pour off the water, carefully remove the gel from the staining tray, and place it in a plastic sandwich-size bag or wrap it in a piece of clear plastic wrap. *Note: Be careful to keep the gel flat as you place it in the bag or plastic wrap.*
6. The gel can be placed in the refrigerator until the next lab period.

**Figure 59.6** shows an agarose gel stain with methylene blue following electrophoresis of plasmid DNA.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

1. Tape a millimeter ruler to your light box or to the glass on your overhead projector. If you prefer, you may use a millimeter ruler to measure the migration of the plasmid.
2. Align your gel so that the front end of the well is set at the zero point on the ruler.
3. Measure the migration distances from the front of the well to the front edge of the band and record the distance in the following chart.
4. Prepare a standard curve on the semilog paper provided on page 405 by plotting the distance traveled in millimeters on the x-axis versus the size of the fragment of *Hind*III-cut  $\lambda$  DNA in kilobases. Record your results in the following chart.

*Note: 0.8% agarose gel has pore sizes that will allow the free movement of nucleic acids between 0.5 and 10 kb. Therefore, draw the best-fit straight line for all bands except the 23-kb band.*

Migration Distances of Lambda DNA						
Kilobases	23	9.4	6.6	4.4	2.3	2.0
Millimeters						

5. Determine the number of bands in each plasmid and use a ruler to measure the migration distance in centimeters.

Band #	1	2	3	4	5	Total
<i>E. coli-1</i>						
<i>E. coli-2</i>						

6. Draw a diagram of your agarose gel and indicate which bands are linear, closed circular, or nicked circular.

1	2	3	4	5
$\lambda$ DNA	EC-2	EC-1	EC-2	EC-1
<input type="text"/>	<input type="text"/>	<input type="text"/>	<input type="text"/>	<input type="text"/>



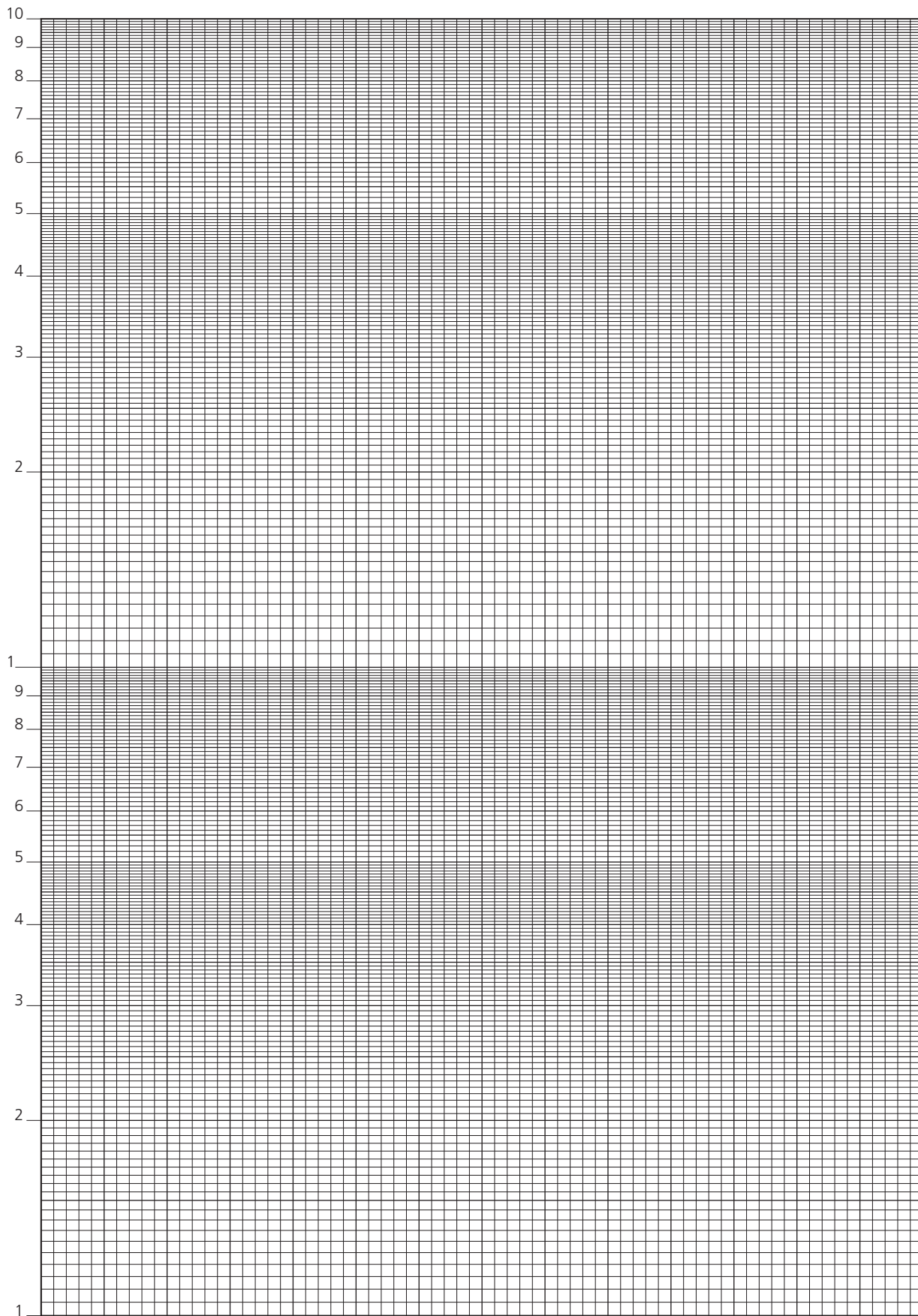
7. Determine the size of linear DNA segments that would migrate the same distance as the various forms of the plasmids, using the standard curve. Record your results below:

Linear DNA	PLASMID DNA	
	<i>E. coli-1</i>	<i>E. coli-2</i>

## Review Questions

1. Why is plasmid DNA preferred for genetic engineering studies?
2. What are selectable markers and why are they important to cloning vectors?
3. What is the rationale for using the following solutions for the isolation of plasmids?
  - a. EDTA:
  - b. Sodium dodecylsulfate (SDS):
  - c. Potassium acetate:

(Questions continue on page 406.)



Migration Distance (mm) versus Fragment Size (kb)

4. Alcohol is considered to be a significant reagent for the isolation of nucleic acids (RNA and DNA). Why is this so?

5. What is nicked circular DNA and why is it termed “relaxed”?

6.  When might you not be able to use a standard curve to determine the size of a plasmid?

7.  When plasmids are isolated from bacterial cells, they may exist in a number of forms.

a. List the different forms that may be found.

b. Which do you think would migrate the fastest and farthest in an electrophoresis experiment and why?

# Restriction Analysis and Electrophoretic Separation of Bacteriophage Lambda DNA

## LEARNING OBJECTIVES


Once you have completed this experiment, you should be able to

1. Cut DNA into small fragments by using restriction endonucleases.
2. Separate DNA fragments through agarose gel electrophoresis.
3. Determine the length of DNA fragments in kilobases.

## Principle

Through their pioneering work, Werner Arber and Hamilton Smith discovered that bacteria possess enzymes that can act as molecular scissors and cut DNA molecules into smaller fragments. These enzymes, called **endonucleases**, are able to differentiate between DNA endogenous to an organism and foreign DNA, such as that of infecting bacteriophages. Endonucleases can cut foreign DNA, rendering it nonfunctional, which spares the bacterium from infection. For scientists, the discovery of restriction endonucleases has been vital in advancing research over the past 40 years because small DNA fragments are much easier to manipulate than longer DNA strands. Now scientists can accurately map a gene's location on a chromosome and describe its base-pair sequence. **Restriction endonucleases** are also being used to develop DNA recombinants for commercial use, detect genetic defects, map restriction sites on plasmids, and create DNA profiles for use in medicine and forensics.

Endonucleases recognize palindromic sequences, four to six base pairs long, on DNA molecules. In everyday usage, a palindrome is a word that is spelled the same way forward and backward. For example, the word "racecar" is a common palindrome. On a double-stranded DNA molecule, a palindrome is a sequence of base pairs that reads the same on one strand 5' to 3' as

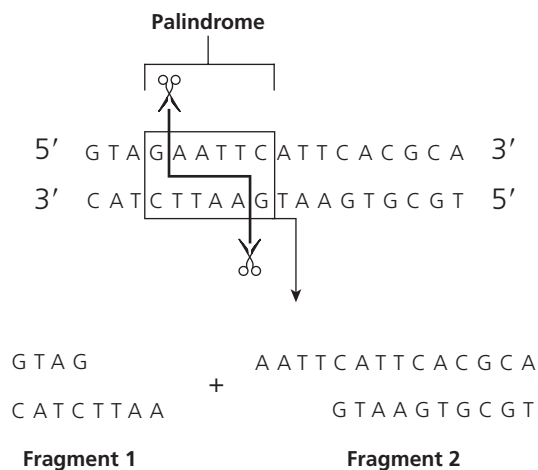


**Figure 60.1** Palindrome for *EcoRI* endonuclease

it does on the other strand 5' to 3'. Keep in mind that reading 5' to 3' on one strand is the opposite direction of reading 5' to 3' on the other. Each endonuclease has its own unique restriction site. **Figure 60.1** shows an example of a palindromic base-pair sequence and the cutting site for the restriction enzyme *EcoRI*. In the figure, *EcoRI* cuts the molecule between guanine and adenine, producing two fragments with staggered ends.

The key property of endonucleases is that they recognize and digest, or cut, one specific sequence of nucleotides on a DNA molecule and cut this same sequence every time. Several endonucleases make staggered cuts in the double-stranded molecule, producing single strands of DNA with cohesive, or sticky, ends that allows them to combine with complementary single-stranded DNA. Other endonucleases cut DNA sequences straight through both strands, producing blunt ends. **Figure 60.2** illustrates the restriction sites of some commonly used endonucleases. The arrows indicate the cutting sites on each strand. The endonucleases that produce sticky, staggered ends are clearly distinguishable in Figure 60.2 from those that produce blunt ends.

DNA fragments cut with the same restriction enzyme can pair with one another. The sticky ends of different strands will join together because of

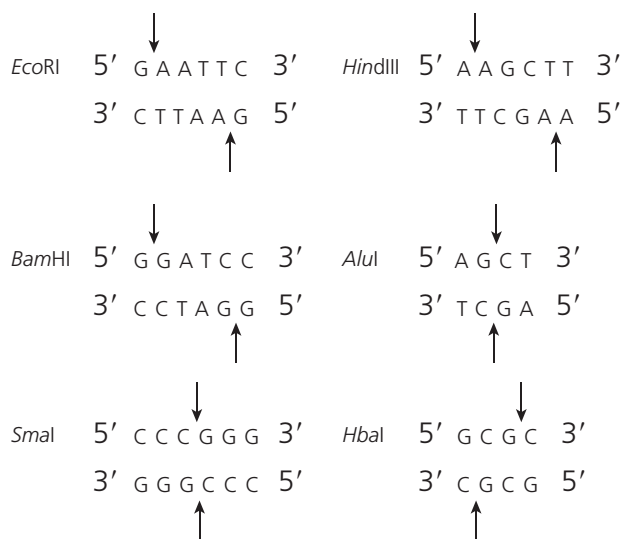

**Figure 60.2** Illustration of restriction sites for common endonucleases that leave blunt and staggered ends

the formation of hydrogen bonds between complementary bases. However, joined fragments lack phosphodiester bonds between guanine and adenine, and nicks form as a result. These nicks are annealed by DNA ligase enzymes. Under optimum environmental conditions (salt concentration, pH, and temperature), restriction endonucleases will cut a strand of DNA into a number of varying-sized fragments. The exact number and size of the fragments depends on the location and number of restriction sites for the enzyme.

Restriction enzymes are named based on the genus and species of bacteria from which they are obtained. The first letter of the genus name is followed by the first two letters of the species name. For example, an endonuclease from *Escherichia coli* is named *Eco*. If a bacterium produces more than one restriction enzyme, each endonuclease is differentiated by Roman numerals. If the enzyme is coded for on a resistance factor, it is further designated with an "R." Thus *EcoRI* is one of several endonucleases produced by *E. coli* and is coded for on a restriction site. Other widely used endonucleases are obtained from *Haemophilus influenzae D* (*HindIII*, which cuts between adenine bases) and *Bacillus amyloliquefaciens H* (*BamHI*, which cuts between two guanine bases).

In the following experiment, you will use endonucleases to cut bacteriophage lambda ( $\lambda$ ) DNA, containing 48,502 base pairs (48.5 kb), into fragments. You will separate the fragments by using agarose gel electrophoresis, and you will determine the size of each.

## CLINICAL APPLICATION

### Enzyme Digestion to Isolate Genes of Interest

It was quickly recognized that restriction enzymes would provide a new tool for the investigation of gene organization, function, and expression. Current medical research is examining how restriction enzymes can be used to insert the genes required for insulin production into diabetic patients.

## AT THE BENCH



## Materials

### DNA Source

Bacteriophage  $\lambda$  (200  $\mu$ l).

### Restriction Endonucleases

*EcoRI*, *HindIII*, and *BamHI*.

### Reagents

Tris-acetate buffer, type-specific buffers for *EcoRI*, *HindIII*, and *BamHI*, electrophoresis loading dye, Carolina Blu or 0.025% methylene blue stain, and 0.8% agarose in 1X TAE buffer. *Note: Formulas for the preparation of type-specific buffers may be found in Appendix 4.*

### Equipment

Plastic 1.5-ml microcentrifuge tubes, microcentrifuge, adjustable micropipettes (0.5 to 10  $\mu$ l), (5 to 10  $\mu$ l), (10 to 100  $\mu$ l), large and small fine-point micropipette tips, waterbath, ice bucket, crushed ice, staining tray, disposable gloves, glassware marking pencil, Bunsen burner or hot plate, 250-ml Erlenmeyer flask, beaker for waste, micro test tube racks, electrophoretic apparatus, millimeter ruler, and a light box or overhead projector.

## Procedure

*Note: In the following procedure, the steps for practicing micropipette use, casting a gel, and loading samples were covered in Experiment 59: Isolation of Bacterial Plasmids. If you completed these practice steps in Experiment 59, you may want to skip them here. If you are uncomfortable with any of the techniques, practice them again.*

**TABLE 60.1** Microliters ( $\mu\text{l}$ ) per Digestion Tube

TUBE	LAMBDA DNA	dH <sub>2</sub> O	10X RESTRICTION BUFFER	ENZYMES			TOTAL VOLUME
				<i>EcoRI</i>	<i>HindIII</i>	<i>BamHI</i>	
B	6	10	2	0	0	2	20
E	6	10	2	2	0	0	20
H	6	10	2	0	2	0	20
L	6	12	2	0	0	0	20

## Using a Micropipette

Before the start of the experiment, familiarize yourself with the use of a micropipette, the function of which is to accurately deliver microliter volumes of solution. Not all micropipettes work the same way. Some are designed to deliver a fixed volume, while others can deliver variable volumes. Your instructor will demonstrate the proper handling and use of these expensive instruments.

Using samples of colored water, practice using a micropipette, attaching different-sized micropipette tips and delivering various sample volumes to digestion tubes.

1. Set the scale on the pipette to the volume you wish to deliver.
2. Place a tip on the micropipette by pushing it firmly onto the pipette.
3. Depress the plunger to the first stop. This is necessary to remove all of the air from the tip.
4. To load the pipette, dip the pipette tip into the solution and release the plunger slowly to draw up the solution.
5. Touch the end of the tip to the side of the tube to remove any excess solution.
6. To deliver the solution, touch the side of the micropipette tip to the inside of the tube receiving the solution to produce a capillary.
7. Depress the plunger to the *first stop* and then continue depressing the plunger to the *second stop* to deliver the full volume of sample, blowing out the last bit in the tip.
8. Continue depressing the plunger while you remove the pipette tip from the tube. *Note: Releasing the plunger before removing the tip of the pipette from the tube will cause fluid to suck back into the tip.*

## Digesting of Lambda ( $\lambda$ ) DNA

 **Be sure to wear gloves, as enzymes on your skin degrade DNA in the experiment.**

1. Obtain a sample of lambda DNA from the instructor.
2. With a glassware marking pencil, label four microcentrifuge tubes with your name or group number followed by an “L” for the uncut DNA, “E” for *EcoRI*, “H” for *HindIII*, and “B” for *BamHI*.
3. Using a new pipette tip for each reagent, add the reagents to the digestion tubes in the following order:
  - a. Lambda DNA
  - b. Deionized or distilled water (dH<sub>2</sub>O)
  - c. Restriction enzyme buffer 10X
  - d. Restriction endonucleases (10 units/ $\mu\text{l}$ )

*Note: The restriction enzyme must be added last to the digestion tubes. Addition of the endonucleases before DNA or buffer may inactivate the endonuclease. Each reagent is added with a new pipette tip to avoid contaminating the digestion tubes.*
4. The addition of the reagents to each tube may be made following the scheme in Table 60.1.
5. Pulse centrifuge or tap your finger on each tube several times to mix the reagents.
6. Place all tubes in a foam rubber test tube rack or a suitable microcentrifuge rack and incubate them in the waterbath at 37°C for 60 minutes.

7. The digestion tubes may be stored in the refrigerator until the next class period. If you are continuing with the experiment now, place the tubes in an ice bucket and proceed to the next step.

## Casting the Agarose Gel

*Note: Not all casting trays are the same. Your instructor will indicate which type will be used and whether there are special considerations during the setup.*

Your instructor prepared the 0.8% agarose gel in a 1X TAE buffer solution before class and maintained it at 55°C in a waterbath. It is ready to pour. One or two drops of Carolina Blu stain was added to the agarose buffer solution to give a small tinge of blue to the gel. At this concentration, the pores that form the gel lattice are such that they allow the free migration of the cut DNA fragments between 0.5 and 10 kb.

Refer to Figure 59.4 on page 401 to see the proper setup of an electrophoretic apparatus.

1. Place the casting tray inside the casting tray box on a level surface.
2. Close off the ends of the tray with the rubber dams by tightening the knob on the top of the casting tray box.
3. Place a well-forming comb in the first notch at the end of the casting tray.
4. Pour 60–70 ml of agarose solution that has been cooled to 55°C into the tray. Use a toothpick or applicator stick to move the bubbles to the edge of the gel before it solidifies.
5. Allow the gel to solidify completely. It should be firm to the touch after 20 minutes.
6. Slowly remove the rubber dams and *very gently* remove the well-forming comb by pulling it straight up. *Note: Use extreme care not to damage or tear the wells.*
7. Place the gel on the platform in the electrophoresis box so that the formed wells are properly oriented toward the anode (negative pole with black cord). Because DNA is negatively charged, the cut DNA fragments will migrate to the cathode (positive pole with the red cord).
8. Fill the electrophoresis box with TAE buffer to a level that just covers the gel, about 2 mm. Make sure that all of the wells are filled with the buffer.


## Practicing Gel Loading

Before loading your sample into the wells of the agarose gel, practice this challenging technique. Your instructor will demonstrate the proper method for loading the wells. Each student should practice on a gel that has been prepared earlier by your instructor, not on the gel to be used for running the samples.

1. Load the pipette with 22  $\mu\text{l}$  of loading gel.
2. Hold the pipette with both hands and dip the tip slightly through the buffer covering the gel, with the tip barely in the well.
3. Slowly discharge the contents of the pipette. *Note: The loading gel contains sucrose, which is heavier than the DNA and will pull the sample into the well.*
4. Practice the technique until you are comfortable with it.

## Loading the DNA Digests into the Wells and Electrophoresing the Samples

1. Remove the digestion tubes from the ice bucket and add 4  $\mu\text{l}$  of 6X loading dye to each tube.
2. Pulse centrifuge or tap your finger on each tube several times so that the contents of the tube move to the bottom.
3. Set the dial on the micropipette to deliver 24  $\mu\text{l}$  (20  $\mu\text{l}$  of restriction digests plus 4  $\mu\text{l}$  of loading dye).
4. Deliver each of the four enzyme digests to separate wells in the agarose gel.
5. Remember the order of your samples and the position of each in the agarose gel. Because the gel cannot be marked, you should draw a diagram of the gel and label the position of your samples as shown in **Figure 60.3**.
6. After the samples are loaded into the wells, place the lid on the electrophoresis gel box. Check that the power switch is turned to the “off” position and then attach the electrical leads (red to red and black to black) from the power supply to the box.
7. Turn the power pack on and adjust the rheostat dial to 110V.
8. Electrophorese the gel for 30 to 40 minutes or until the leading edge of the bromphenol blue dye (the dye in the loading gel) has traveled

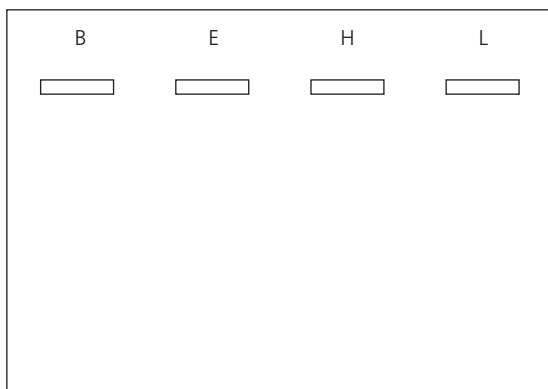

**Figure 60.3** Example of gel loading scheme

roughly three-fourths of the distance to the edge of the gel.

9. Turn the rheostat to zero and turn off the power. Disconnect the leads and remove the cover from the gel box.

## Staining the Gel

1. Put on a pair of disposable laboratory gloves.
2. Lift the gel tray out of the electrophoresis box, and slide the gel into a staining tray containing approximately 100 ml of Carolina Blu stain.
3. Allow the gel to remain in the stain for 30 to 40 minutes.
4. Pour off the stain into a waste beaker. Transfer the gel to a staining tray containing 100 ml of distilled water and allow the gel to decolorize (destain) for another 30 minutes. Frequent destaining with fresh distilled water for 2 minutes increases the intensity of the bands. For best results, let the gel destain overnight in a small volume of water. *Note: If the gel is left overnight in a large volume of water, it may destain too much.*
5. Pour off the water, carefully remove the gel from the staining tray, and place it in a plastic sandwich-size bag or wrap it in a piece of clear plastic wrap. *Note: Be very careful to keep the gel flat as you place it in the bag or plastic wrap.*
6. The gel can be placed in the refrigerator until the next lab period. Refer to Figure 59.6, page 402, for a photo of an electrophoresed and stained gel.



*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

1. Tape a millimeter ruler to the surface of a light box or to the glass surface of an overhead projector.
2. Place the stained gel (inside plastic bag) next to the zero point on the ruler and measure the distances that each fragment (band) migrated. Measure the distance from the front of the well to the front of the band. Record your results in the following chart.

	Migration Distances of Fragments (mm)					
Uncut $\lambda$ DNA						
<i>Bam</i> HI cut $\lambda$ DNA						
<i>Eco</i> RI cut $\lambda$ DNA						
<i>Hind</i> III cut $\lambda$ DNA						

Linear DNA fragments migrate at rates inversely proportional to the  $\log_{10}$  of their molecular weight and base-pair length.

3. A graph (standard curve) can be constructed by plotting known kilobase-pair fragments versus distances migrated from the wells to the front of the fragment. The six kilobase-pair fragment sizes for *Hind*III are well established and can be used to plot a standard curve.
4. Once the fragment sizes are measured and distances traveled are plotted on semilog paper found on page 415, a best-fit straight line can be drawn. The size of each unknown fragment can be determined by drawing a vertical line from the migration distance (mm) on the x-axis up to the point on the curve that intersects that straight line. From there, draw a horizontal line to the fragment size on the y-axis.

5. In the following table, the kilobase lengths of *Hind*III are provided. From the standard curve, use the migration distance you have measured to determine the base-pair lengths for the three restriction enzymes. Record your results in the table.

<i>Hind</i> III*		<i>Bam</i> HI		<i>Eco</i> RI		λ DNA	
Distance (mm)	Actual kb	Distance (mm)	Calculated kb	Distance (mm)	Calculated kb	Distance (mm)	Calculated kb
	27.4*						
	23.1*						
	9.4						
	6.6						
	4.4						
	2.3						
	2.0						

\*Note: Remember that 0.8% agarose allows the free migration of DNA in the range of 0.5 to 10 kb. Therefore, the 27.4- and 23.1-kb fragments will not be detected.

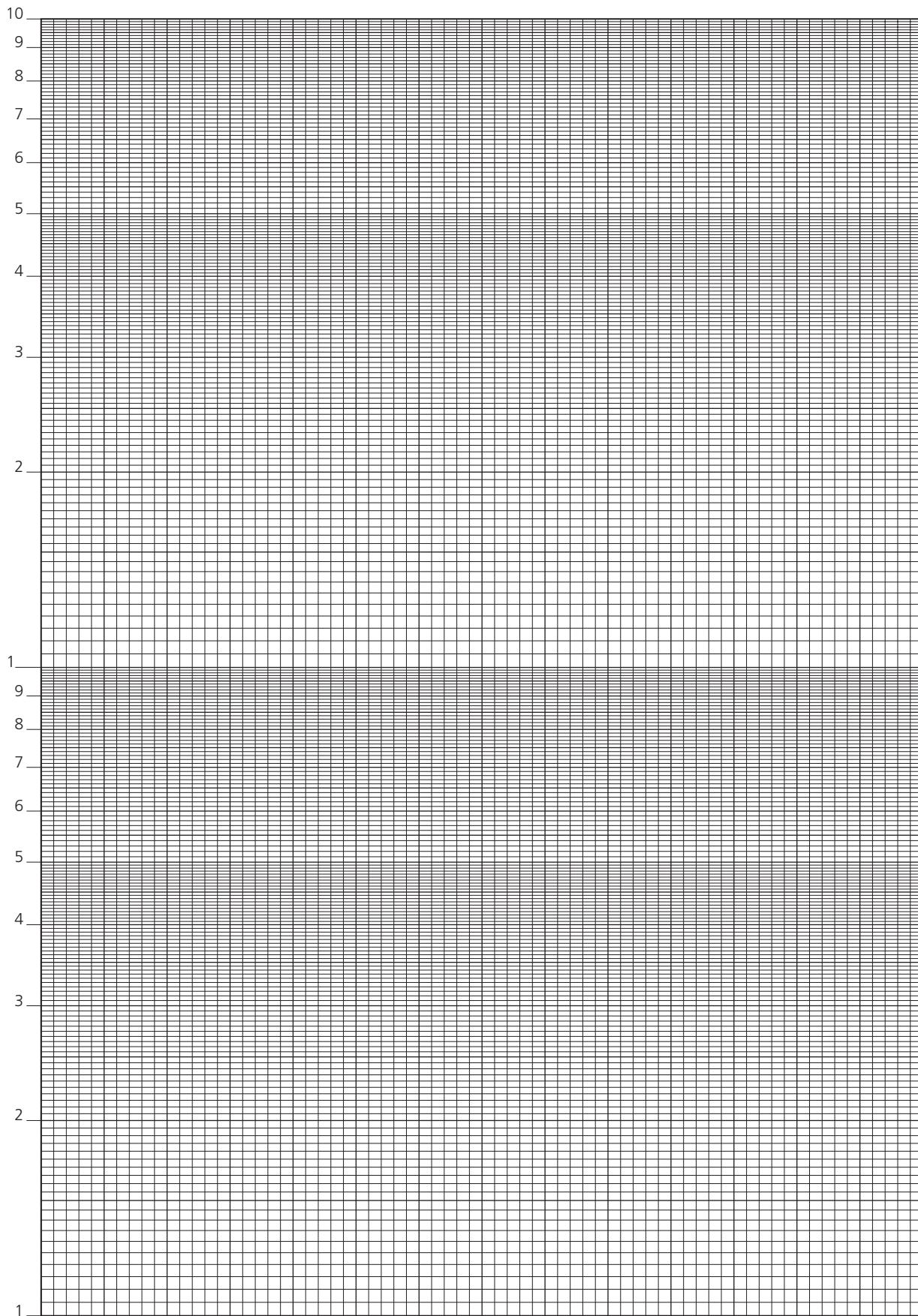
6. Calculate the fragment lengths of *Eco*RI and *Bam*HI from the standard curve and compare them to the actual kilobase lengths listed in the following chart.

<i>Hind</i> III		<i>Bam</i> HI			<i>Eco</i> RI		
Actual kb	Distance (mm)	Actual kb	Distance (mm)	Calculated Length	Actual kb	Distance (mm)	Calculated Length
27.4*		16.8*			24.6*		
23.1*		12.3			21.2*		
9.4		7.2			7.4		
6.6		6.7*			5.8*		
4.3		6.5*			5.6*		
2.3		5.6*			4.9		
2.0		5.5*			3.5		

\*Note: These fragments appear as a single band.



7. Compare and contrast your calculated kilobase pair from the standard curve with the actual kilobase pair for the restriction endonucleases.
- List those that were most accurate.
  - List those that were least accurate.

(See Review Questions on page 416.)



Migration Distance (mm) versus Fragment Size (kb)

## Review Questions

1. Why were the DNA digestions carried out at 37°C?
2. Would any or all of these endonucleases cut the DNA of another bacteriophage or bacterium?
3. Why were specific restriction buffers needed for each of the restriction enzyme digests?
4. What could account for low endonuclease activity?
5. Why are the restriction enzymes added last to the digestion mixtures?
6.  Assume you have one organism with a gene for ampicillin resistance and another organism with a gene for luciferinase. How would you isolate the gene from one organism and connect it with the gene of the other organism?
7.  How would restriction enzymes play a role in having an organism produce a protein that it normally doesn't make?

# Medical Microbiology

## LEARNING OBJECTIVES

Once you have completed the experiments in this section, you should be familiar with the

1. Characteristics and methodology for isolating and identifying selected pathogenic microorganisms.
2. Indigenous microbial flora of selected human anatomical sites.

## Introduction

Although microorganisms are ubiquitous and their benefits to humans have been recorded, a small group of organisms remains a focus of concern: They are the pathogens, whose existence makes medical or clinical microbiology an especially important science.

That living agents are capable of inducing infections (*contagium vivum*) was first put forward by the monk Fracastoro in Verona about 500 years ago. In 1659 Kircher reported the presence of minute motile organisms in the blood of plague victims. Two hundred years after Fracastoro developed his initial concept, the germ theory of disease was formulated by Plenciz based on Leeuwenhoek's revolutionary microscopic observation of microorganisms. Perhaps the most important contributions to microbiology were made by Pasteur, Koch, and Lister during the **Golden Era of Microbiology** from 1870 to 1920. These investigators and their students recorded the observations and discoveries that cemented the cornerstone of medical microbiology. The body of knowledge that has accrued since these early years has made clinical microbiology a

major component of laboratory or diagnostic medicine. The major responsibility of this science is isolating and identifying infectious pathogens to enable physicians to treat patients with infectious disease prudently, intelligently, and rapidly.

Many of the experiments so far described have application in the field of clinical microbiology. Among these are isolation and identification of unknown cultures, the use of selective and differential media, and biochemical tests used to separate and identify various microorganisms. Although studying all of the bacterial pathogens responsible for human illness is not possible here, routine experiments for isolating and identifying some of the most frequently encountered infectious organisms and microorganisms that constitute the indigenous flora of the human body are included. The pathogens chosen are pyogenic cocci, members of the genera *Staphylococcus* and *Streptococcus*, the *Enterobacteriaceae*, and the organisms suspected in formation of dental caries. Experimental procedures designed for the detection and presumptive identification of microorganisms in blood and urine, which are normally sterile body fluids, have also been incorporated

into this section. Organisms that naturally reside in or on body surfaces and constitute the body's **normal flora** are also examined.

The need for the expeditious detection and identification of pathogens has led to the development of rapid testing methods. These are microbiologically and immunologically based and can be performed quickly and without the need for sophisticated and expensive equipment. Some prototypic experiments using these rapid methods are included along with the traditional procedures.



**Many of the organisms that are used, although attenuated by having been subcultured on artificial complex media for many generations, must be viewed as potential pathogens and therefore handled with respect. At this point in your training, your manipulative skills should be sufficiently developed, allowing you to perform aseptically in any medical, hospital, or clinical laboratory setting to prevent infection of yourself and others.**

# Microbial Flora of the Mouth: Determination of Susceptibility to Dental Caries

EXPERIMENT

61

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be

1. Familiar with the organisms responsible for dental caries.
2. Able to perform experiments that demonstrate the host's susceptibility to formation of caries.

## Principle

A variety of microorganisms are known to be involved in the formation of dental caries, including *Lactobacillus acidophilus*, *Streptococcus mutans*, and *Actinomyces odontolyticus*. These organisms in the oral flora produce organic acids, particularly lactic acid, by fermenting carbohydrates that adhere to the surface of the teeth. In the continued presence of lactic acid, dental enamel undergoes decalcification and softening, which result in the formation of tiny perforations called dental caries.

The actual mechanism of action of these organisms is still unclear. However, it has been noted that *S. mutans* excretes an enzyme called **dextranucrase** (glycosyl transferase), which is capable of polymerizing sucroses into a large polymer, dextran, plus the monosaccharide fructose. This polysaccharide clings tenaciously to the teeth and forms dental plaque, in which streptococci reside and ferment fructose with the formation of lactic acid (Figure 61.1).

Similarly, *L. acidophilus* produces lactic acid as an end product of carbohydrate fermentation. Oral lactobacilli are capable of metabolizing glucose found in the mouth, producing organic acids that reduce the oral acid concentration to a pH of less than 5. At this pH, decalcification occurs and dental decay begins.

One of the best microbiological methods for determining susceptibility to dental caries is the **Snyder test**. This test measures the amount of acid produced by the action of the lactobacilli on glucose. The test employs a differential medium, Snyder agar (pH 4.7), which contains glucose and

the pH indicator bromocresol green that gives the medium a green color.

Following incubation, Snyder agar cultures containing lactobacilli from the saliva will show glucose fermentation with the production of acid, which tends to lower the pH to 4.4, the level of acidity at which dental caries form. At this pH the green medium turns yellow. A culture demonstrating a yellow color within 24 to 48 hours is suggestive of the host's susceptibility to the formation of dental caries. A culture that does not change color is indicative of lower susceptibility.

## CLINICAL APPLICATION

### Preventing Dental Caries

Factors that help to control the development of dental caries are proper oral hygiene, consumption of adequate fluoride, and moderation in the consumption of foods that cause decay. Foods likely to lead to decay are sticky, highly processed, and high in fermentable carbohydrates such as breads, muffins and dried fruits. Also, the use of products to control oral pH might be helpful to ensure that bacteria that cause caries will not flourish.

## AT THE BENCH



## Materials

### Cultures

Organisms of the normal oral flora present in saliva.


### Media

Per designated student group: two Snyder test agar deep tubes.

### Equipment

Bunsen burner, ice-water bath, 1-in. square blocks of paraffin, sterile 1-ml pipettes, mechanical pipetting device, sterile test tubes, and glassware marking pencil.



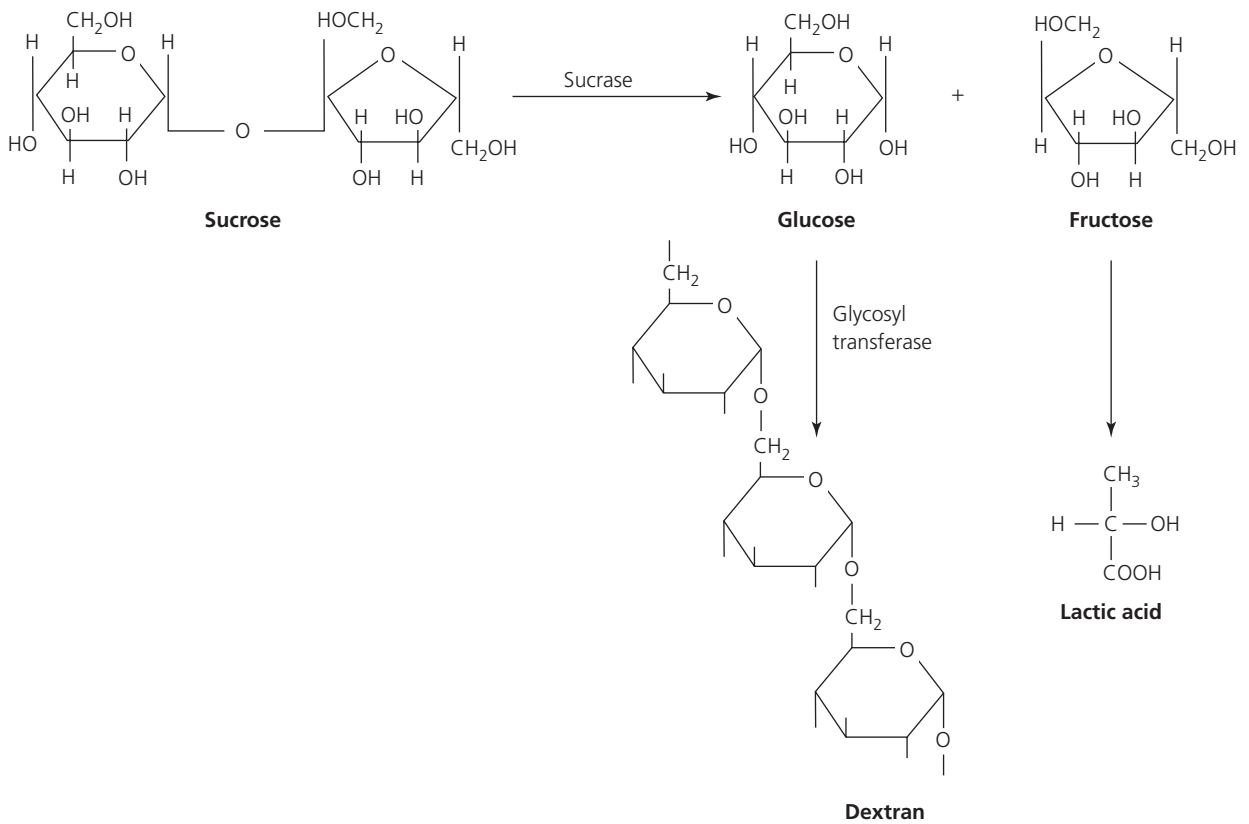

**Figure 61.1** Degradation of sucrose and subsequent conversion of glucose into dextran by *Streptococcus mutans*

## Procedure Lab One


1. Melt two appropriately labeled Snyder agar deep tubes and cool to 45°C.
2. Chew one square of paraffin for 3 minutes *without swallowing the saliva*. As saliva develops, collect it in a sterile test tube.
3. Vigorously shake the collected saliva sample and transfer 0.2 ml of saliva with a sterile pipette into one of the Snyder test medium tubes that have been cooled to 45°C. *Note: Don't let the pipette touch the sides of the tubes or the agar.*
4. Mix the contents of the tube thoroughly by rolling the tube between the palms of your hands or by tapping it with your finger.
5. Rapidly cool the inoculated tube of Snyder agar in an ice-water bath.
6. Repeat Steps 3 through 5 to inoculate the second tube.
7. Incubate both tubes for 72 hours at 37°C. Observe cultures at 24, 48, and 72 hours.

## Procedure Lab Two

1. Examine the Snyder test cultures daily during the 72-hour incubation period for a change

in the color of the culture medium. Use an uninoculated tube of the medium as a control. **Figure 61.2** shows a positive and a negative Snyder test.

2. Record the color of the cultures in the Lab Report.

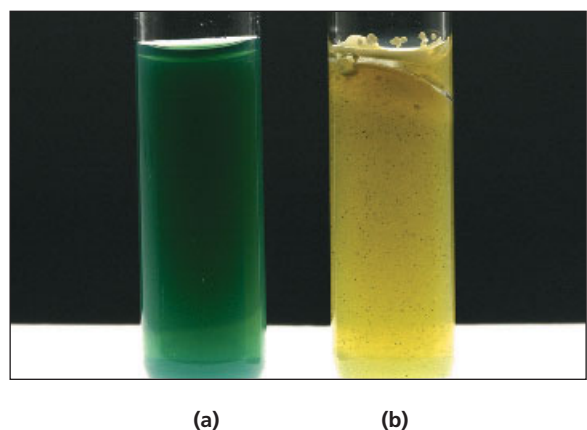

**Figure 61.2** Snyder test. (a) No change in the color indicates a negative result. (b) The color change to yellow indicates a positive result.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

Using Table 61.1 to interpret your observations, record your findings about susceptibility to caries in the chart below.

**TABLE 61.1** Assessment of Susceptibility to Dental Caries

CARIES ACTIVITY	HOURS OF INCUBATION		
	24	48	72
Marked	Positive	...	...
Moderate	Negative	Positive	...
Slight	Negative	Negative	Positive
Negative	Negative	Negative	Negative

*Source:* Courtesy of Difco Laboratories, Inc., Detroit, Michigan.

Positive: Complete color change; green is no longer dominant.

Negative: No color change or a slight color change; medium retains green color throughout.

Tube Number	COLOR OF SNYDER TEST CULTURES			Caries Susceptibility (Yes or No)
	24 hr	48 hr	72 hr	

## Review Questions

- How would you explain the differential nature of the Snyder agar medium as used for the detection of dental caries?

2. How would you explain the mechanism responsible for the formation of dental caries by resident microorganisms?

3. What is the function of the paraffin in this procedure?

4. Based on your results, what is your tendency to form dental caries?

Is this result consistent with your dental history?

5.  Are all members of the resident flora of the mouth capable of initiating dental caries? Explain.

6.  What is the ideal time of day to perform this procedure? Why?

## LEARNING OBJECTIVE

Once you have completed this experiment, you should be able to

1. Identify microorganisms that normally reside in the throat and skin.

## Principle

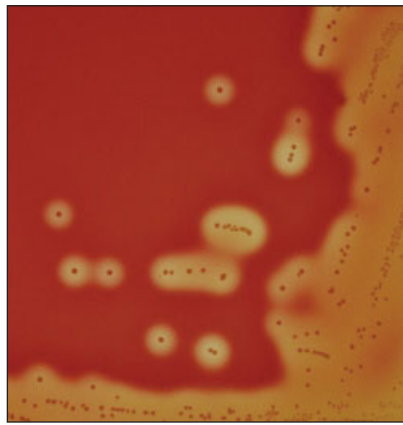
The normal flora are regularly found in specific areas of the body. This specificity is far from arbitrary and depends on environmental factors such as pH, oxygen concentration, amount of moisture present, and types of secretions associated with each anatomical site. Native microbial flora are broadly located as follows:

1. **Skin:** Staphylococci (predominantly *Staphylococcus epidermidis*), streptococci (alpha-hemolytic, nonhemolytic), enterococci, diphtheroid bacilli, yeasts, and fungi.
2. **Eye conjunctiva:** Staphylococci, streptococci, diphtheroids, and neisseriae.
3. **Upper respiratory tract:** Staphylococci, streptococci (alpha-hemolytic, nonhemolytic, and *Streptococcus pneumoniae*), enterococci, diphtheroids, spirochetes, and members of the genera *Moraxella* (formerly called *Branhamella*), *Neisseria*, and *Haemophilus*.
4. **Mouth and teeth:** Anaerobic spirochetes and vibrios, fusiform bacteria, staphylococci, and anaerobic levan-producing and dextran-producing streptococci responsible for dental caries.
5. **Intestinal tract:** In the upper intestine, predominantly lactobacilli and enterococci. In the lower intestine and colon, 96% to 99% is composed of anaerobes such as members of the genera *Bacteroides*, *Lactobacillus*, *Clostridium*, and *Streptococcus*, and 1% to 4% is composed of aerobes, including coliforms, enterococci, and a small number of *Proteus*, *Pseudomonas*, and *Candida* species.
6. **Genitourinary tract:** Staphylococci, streptococci, lactobacilli, gram-negative enteric bacilli, clostridia, spirochetes, yeasts, and protozoa such as *Trichomonas* species.


In this exercise, you will study the resident flora of the throat and skin. Since these sites represent sources of mixed microbial populations, you will perform streak-plate inoculations, as outlined in Experiment 3, to effect their separations. The discrete colonies thus formed can be studied morphologically, biochemically, and microscopically to identify the individual genera of these mixed flora.

The procedure used to identify the native flora of the throat involves the following steps:

1. A *blood agar plate* is inoculated to demonstrate the alpha-hemolytic and beta-hemolytic reactions of some streptococci and staphylococci. Hemolytic reactions on blood agar are shown in **Figure 62.1**. A distinction between these two genera can be made based on their colonial and microscopic appearances. The streptococci typically form pinpoint colonies on blood agar, whereas the staphylococci form larger pinhead colonies that might show a golden coloration. When viewed under a microscope, the streptococcal cells form chains of varying lengths, whereas the staphylococci are arranged in clusters.
2. A *chocolate agar plate* is inoculated to detect *Neisseria* spp. by means of the oxidase test. Members of this genus are recognized when the colonies develop coloration that is pink to dark purple upon addition of *p*-aminodimethylaniline oxalate following incubation. **Figure 62.2** shows colonies growing on chocolate agar from a throat culture.
3. A *Mueller-Hinton tellurite* or *Tinsdale agar plate* is inoculated to demonstrate the presence of diphtheroids, which appear as black, pinpoint colonies on this medium (**Figure 62.3**). This coloration is due to the diffusion of the




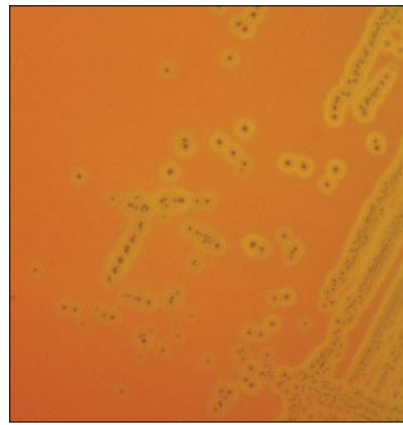(a) Beta hemolysis




(b) Alpha hemolysis

**Figure 62.1** Beta- and alpha-hemolytic reactions on blood agar

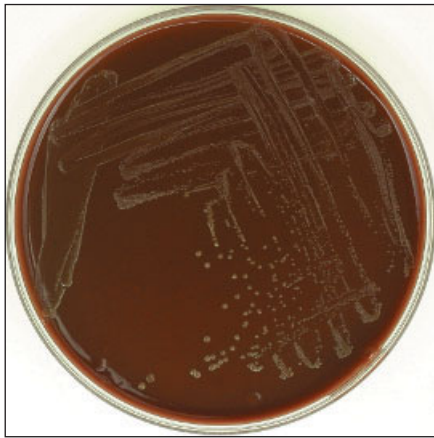

**Figure 62.2** Colony growth on chocolate agar from a throat culture

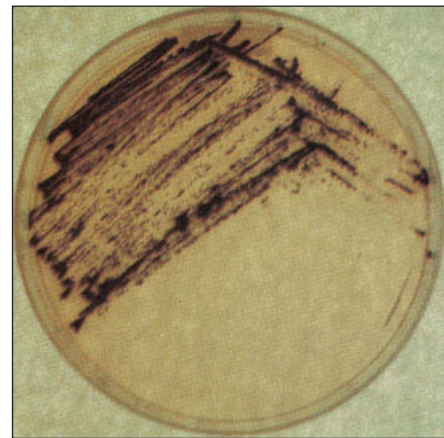

**Figure 62.3** Mueller-Hinton tellurite agar plate. Growth of black, pinpoint colonies indicates the presence of diphtheroids.

tellurite ions into the bacterial cells and their subsequent reduction to tellurium metal, which precipitates inside the cells.


The procedure used to identify the native flora of the skin involves the following steps:

1. A *blood agar plate* inoculated to determine the presence of hemolytic microorganisms, specifically the staphylococci and streptococci: Differentiation between these two genera may be made as previously described.
2. A *mannitol salt agar plate* inoculated for the isolation of the staphylococci: The generally avirulent staphylococcal species can be differentiated from the pathogenic *Staphylococcus aureus* because the latter is able to ferment mannitol, causing yellow coloration of this medium surrounding the growth. **Figure 62.4**
3. A *Sabouraud agar plate* inoculated to detect yeasts and molds: Yeast cells will develop pigmented or nonpigmented colonies that are elevated, moist, and glistening. Mold colonies will appear as fuzzy, powdery growths arising from a mycelial mat in the agar medium. **Figure 62.5** shows yeast colonies and a mold colony.
4. *Chocolate agar plate* inoculated to detect *Neisseria* spp. The presence of *Neisseria* spp. produces pink to purple to black colonies on this medium.
5. *Mueller-Hinton tellurite* or *Tinsdale media* inoculated to detect *Corynebacterium* spp. (diphtheroids). These colonies are black in appearance.

## CLINICAL APPLICATION

### Skin Flora and Acne

The bacterial population on a single human's skin is about  $10^{12}$  organisms. A normal flora of microorganisms colonizes the human skin at birth as it passes through the birth canal, and typically inhabits the superficial layers of the epidermis and upper parts of the hair follicles. They consist mainly of *Staphylococcus epidermidis*, *Micrococcus*, and corynebacteria such as *Propionibacterium*. *Propionibacterium acnes* is normally found in low concentrations, but overgrows in the anaerobic environment of a blocked hair follicle, producing acne.

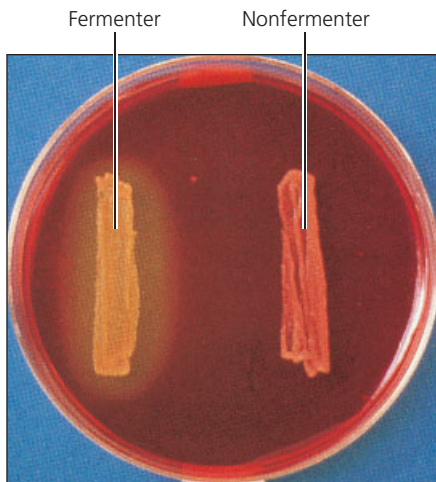

**Figure 62.4** Mannitol salt agar plate showing a fermenter and a nonfermenter organism



(a) Yeast colonies



(b) Mold colony

**Figure 62.5** Sabouraud agar plate. (a) Yeast colonies have an elevated, moist, and glistening appearance. (b) A mold colony shows fuzzy, powdery growth.

## AT THE BENCH



## Materials

### Media

Per designated student group: two blood agar plates, two mannitol salt agar plates, one chocolate agar plate, one Mueller-Hinton tellurite or Tinsdale agar plate, one Sabouraud agar plate, and two 5-ml sterile saline tubes.

### Reagents

Crystal violet, Gram's iodine, safranin, 1% *p*-aminodimethylaniline oxalate, and lactophenol-cotton-blue.

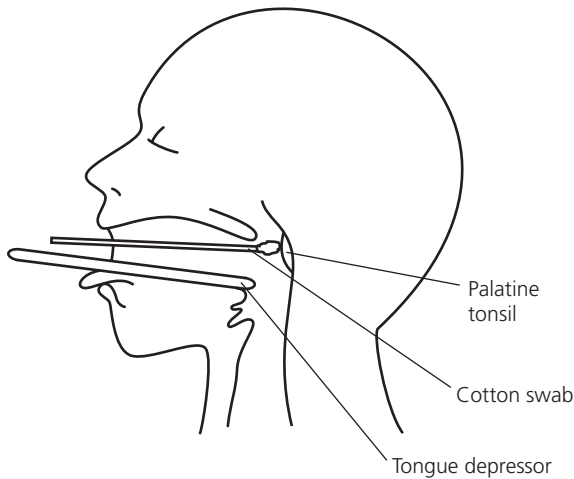
### Equipment

Sterile cotton swabs, tongue depressors, desiccator jar with candle, microscope, glass slides, Bunsen burner, glassware marking pencil, and disposable gloves.

## Procedure Lab One

 You must wear disposable gloves in Steps 1–3.

1. Place a tongue depressor on the extended tongue and with a sterile cotton swab, obtain a specimen from the palatine tonsil by rotating the swab vigorously over its surface without touching the tongue, as illustrated.



2. Inoculate a tube of sterile saline with the swab and mix until you have a uniform suspension.
3. Using a sterile inoculating loop, inoculate one plate each of blood agar, chocolate agar, mannitol salt agar, and Mueller-Hinton tellurite or Tinsdale agar, all previously labeled with the source of the specimen, by means of a four-way streak inoculation as described in Experiment 3.
4. Using a sterile cotton swab moistened in sterile saline, obtain a specimen from the skin by rubbing the swab vigorously against the palm of your hand.
5. Inoculate a tube of sterile saline with the swab and mix the solution.
6. Inoculate one plate each of blood agar, mannitol salt agar, and Sabouraud agar, as described in Step 3.
7. Incubate the inverted chocolate agar plate in a CO<sub>2</sub> incubator, in a CO<sub>2</sub> incubation bag, or in a candle jar. If you use the candle jar, place a lighted candle in a desiccator jar and cover the jar tightly to effect a 5% to 10% CO<sub>2</sub> environment required for the growth of the *Neisseria*. Incubate the jar for 48 hours at 37°C.
8. Incubate the inverted Sabouraud agar plate for 48 hours at 25°C and the remaining plates for 48 hours at 37°C.

## Procedure Lab Two

### Selection and Differentiation of Skin and Throat Isolates

1. Examine the blood agar plate cultures for zones of hemolysis (refer to Figure 62.1 and Experiment 15).
2. Add the *p*-aminodimethylaniline oxalate to the surface of the growth on the chocolate agar plate. Observe for the appearance of a pink-to-purple-to-black color on the surface of any of the colonies (Figure 62.2).
3. Examine the Mueller-Hinton tellurite or Tinsdale agar plate for the presence of black colonies (Figure 62.3).
4. Examine the Sabouraud agar plate for the appearance of moldlike growth (Figure 62.5).
5. Examine the mannitol salt agar plate for the presence of growth that is indicative of staphylococci. Then examine the color of the medium surrounding the growth. A yellow color is indicative of *S. aureus* (refer to Figure 62.4 and Experiment 15).
6. Record your observations in the Lab Report and indicate the types of organisms that may be present in each specimen.

### Staining and Morphological Characteristics of Skin and Throat Isolates

1. Prepare two Gram-stained smears from each of the blood agar cultures, choosing well-isolated colonies that differ in their cultural appearances and demonstrate hemolytic activity. Observe microscopically for the Gram reaction and the size, shape, and arrangement of the cells. Record your observations in the Lab Report and attempt to identify each isolate.
2. Prepare two lactophenol–cotton-blue–stained smears of organisms obtained from discrete colonies that differ in appearance on the Sabouraud agar culture (refer to Experiment 36). Observe microscopically, draw a representative field in the Lab Report, and attempt to identify the fungi by referring to Experiment 38.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

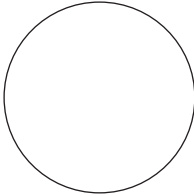
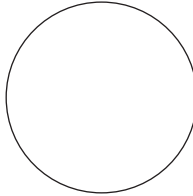
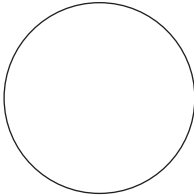
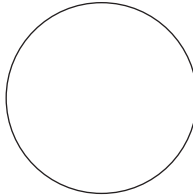
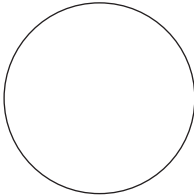
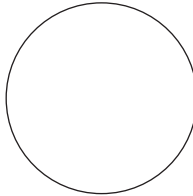
## Observations and Results

### Selection and Differentiation of Skin and Throat Isolates


Cultures	Throat Specimen	Skin Specimen
Blood agar: <i>Staphylococcus</i> spp. <i>Streptococcus</i> spp. Type of hemolysis: alpha beta		
Chocolate agar: <i>Neisseria</i> spp. (+) or (-) pink to purple to black		
Mueller-Hinton tellurite or Tinsdale: <i>Corynebacterium</i> spp. (+) or (-) black colonies		
Sabouraud agar: Fungal colonies (+) or (-)		
Mannitol salt agar: <i>Staphylococcus aureus</i> Other <i>Staphylococcus</i> spp. ( <i>S. epidermidis</i> , <i>S. saprophyticus</i> ) (+) or (-) growth Color of medium		
Types of organisms present		



## Staining and Morphological Characteristics of Skin and Throat Isolates

<b>Skin Specimen</b>	Isolate 1	Isolate 2
Draw a representative field.		
Gram reaction		
Morphology		
Organism		
<b>Throat Specimen</b>	Isolate 1	Isolate 2
Draw a representative field.		
Gram reaction		
Morphology		
Organism		
<b>Sabouraud Agar Colonies Specimen</b>	Isolate 1	Isolate 2
Draw a representative field.		
Morphology		
Organism		

## Review Questions

1.  Why are some microorganisms termed “normal flora,” and of what value are they to the well-being of the host?
  
2. A 6-year-old female is taken to her pediatrician for a checkup. As the doctor takes the child’s history, her mother reports that the child had a severe sore throat several weeks earlier that regressed without treatment. Upon examination the pediatrician notes that the child has a systolic heart murmur consistent with mitral insufficiency and suspects that she has rheumatic fever.
  - a. How was the earlier pharyngitis related to the subsequent development of rheumatic fever?
  
  - b. Rheumatic fever is diagnosed on clinical and serological findings. What test should be done to diagnose rheumatic fever?
  
  - c. How are rheumatic fever patients treated?

**3.** A 35-year-old female underwent serious abdominal surgery involving extensive bowel resection. She was maintained postoperatively on a regimen of intravenous broad-spectrum antibiotics. Three days postoperative she spiked a fever without a clear source. She complains of vaginal discomfort. Blood cultures reveal the presence of an ovoid cell that reproduced by budding.

**a.** Based on this observation, what do you think this organism is?

**b.** Is it part of the normal flora in humans?

**c.** How did the treatment with broad-spectrum antibiotics predispose the patient to infection with this organism?

## LEARNING OBJECTIVES

Once you have completed this experiment, you should understand

1. The medical significance of the staphylococci.
2. Selected laboratory procedures designed to differentiate among the major staphylococcal species.

## Principle

The genus *Staphylococcus* is composed of both pathogenic and nonpathogenic organisms. The three major species are *S. aureus*, *S. saprophyticus*, and *S. epidermidis*. Strains of the last two species are generally avirulent; however, under special circumstances in which a suitable portal of entry is provided, *S. epidermidis* may be the etiological agent for **skin lesions** and **endocarditis**, and *S. saprophyticus* has been implicated in some **urinary tract infections**. **Figure 63.1** is a streak-plate culture of *Staphylococcus aureus*.

Infections are primarily associated with *S. aureus* pathogenic strains that are often responsible for the formation of **abscesses**, localized pus-producing lesions. These lesions most commonly occur in the skin and its associated structures, resulting in **boils**, **carbuncles**, **acne**, and **impetigo**. Infections of internal organs and tissues are not uncommon, however, and include **pneumonia**, **osteomyelitis** (abscesses in bone and bone marrow), **endocarditis** (inflammation of the endocardium), **cystitis** (inflammation of the urinary bladder), **pyelonephritis** (inflammation of the kidneys), **staphylococcal enteritis** due to enterotoxin contamination of foods, and, on occasion, **septicemia**.

Strains of *S. aureus* produce a variety of metabolic end products, some of which may play roles in the organisms' pathogenicity. Included among these are **coagulase**, which causes clot formation; **leukocidin**, which lyses white blood cells; **hemolysins**, which are active against red blood



**Figure 63.1** Streak-plate culture of *Staphylococcus aureus*. Produces colonies that are circular, convex, smooth, and cream-colored to golden yellow in appearance.

cells; and **enterotoxin**, which is responsible for a type of gastroenteritis. Additional metabolites of a nontoxic nature are **DNase**, **lipase**, **gelatinase**, and the fibrinolysin **staphylokinase**.

When there is a possibility of staphylococcal infection, isolation of *S. aureus* is of clinical importance. These virulent strains can be differentiated from other staphylococci and identified by a variety of laboratory tests, some of which are illustrated in Table 63.1.

In this exercise you will distinguish among the staphylococcal species by performing traditional test procedures, a computer-assisted multitest procedure, or a newer rapid latex agglutination test.

## Traditional Procedures

The traditional procedures involve the following steps:

1. **Mannitol salt agar:** This medium is selective for salt-tolerant organisms such as staphylococci. Differentiation among the staphylococci is predicated on their ability to ferment mannitol.

**TABLE 63.1 Laboratory Tests for Differentiation of Staphylococcal Species**

TEST	<i>S. aureus</i>	<i>S. epidermidis</i>	<i>S. saprophyticus</i>
Mannitol salt agar			
Growth	+	+	+
Fermentation	+	–	–
Colonial pigmentation	Generally golden yellow	White	White
Coagulase	+	–	–
DNase	+	–	–
Hemolysis	Generally beta	–	–
Novobiocin sensitivity	Sensitive	Sensitive	Resistant

Following incubation, mannitol-fermenting organisms, typically *S. aureus* strains, exhibit a yellow halo surrounding their growth, and nonfermenting strains do not. (Refer to Figure 62.4.) It should be noted that other salt-tolerant microorganisms, such as the enterococci, are capable of growth on mannitol salt agar. These two genera are easily differentiated by performing a catalase test. *Staphylococcus* will grow in the presence of catalase while the enterococci will not.

- 2. Coagulase test:** Production of coagulase is indicative of an *S. aureus* strain. The enzyme acts within host tissues to convert fibrinogen to fibrin. It is theorized that the fibrin meshwork that is formed by this conversion surrounds the bacterial cells or infected tissues, protecting the organism from nonspecific host resistance mechanisms such as phagocytosis and the antistaphylococcal activity of normal serum. In the coagulase tube test for bound and free coagulase, a suspension of the test organism in citrated plasma is prepared and the inoculated plasma is then periodically examined for fibrin formation, or coagulation. Clot formation within 4 hours is interpreted as a positive result and indicative of a virulent *S. aureus* strain. The absence of coagulation after 24 hours of incubation is a negative result, indicative of an avirulent strain (Figure 63.2).
- 3. Deoxyribonuclease (DNase) test:** Generally, coagulase-positive staphylococci also produce the hydrolytic enzyme DNase; thus this test can be used to reconfirm the identification of *S. aureus*. The test organism is grown on an agar medium containing DNA. Following incubation, DNase activity is determined by the addition of

0.1% toluidine blue to the surface of the agar. DNase-positive cultures capable of DNA hydrolysis will show a rose-pink halo around the area of growth. The absence of this halo is indicative of a negative result and the inability of the organism to produce DNase (Figure 63.3).

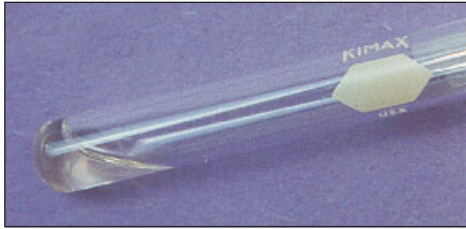
- 4. Novobiocin sensitivity:** This test is used to distinguish *S. epidermidis* from *S. saprophyticus*. The Mueller-Hinton agar plate is heavily seeded with the test organism to produce a confluent growth on the agar surface. After the seeding, a 30- $\mu$ g novobiocin antibiotic disc is applied to the agar surface. Following incubation, the sensitivity of an organism to the antibiotic is determined by the Kirby-Bauer method as described in Experiment 44 and as shown in Figure 63.4.

### STAPH-IDENT System Procedure

A computer-assisted procedure is the API<sup>®</sup> (Analytical Profile Index) STAPH-IDENT<sup>®</sup> system (developed by Analytab Products, Division of Sherwood Medical, Plainview, New York). STAPH-IDENT is a rapid, computer-based micromethod for the separation and identification of the newly proposed 13 species of staphylococci. The system consists of 10 microcupules containing dehydrated substrates for the performance of conventional and chromogenic tests. The addition of a suspension of the test organism serves to hydrate the media and to initiate the biochemical reactions. The identification of the staphylococcal species is made with the aid of the differential charts or the STAPH-IDENT Profile Register that is part of the system (Table 63.2 on page 435), or both.



(a) Positive coagulase test




(b) Negative coagulase test

**Figure 63.2 Coagulase test.** (a) Clot formation indicates a positive result; (b) the absence of coagulation is a negative result.

### Latex Agglutination Procedure

The latex agglutination test as a rapid diagnostic slide test for *Staphylococcus aureus*. The Remel BactiStaph<sup>®</sup> diagnostic kit (Fisher Health Care) uses protein-coated latex particles that are able to detect the clumping factor (bound coagulase and protein A) that causes the *S. aureus* to adhere to the black latex particles, producing a visible agglutination.

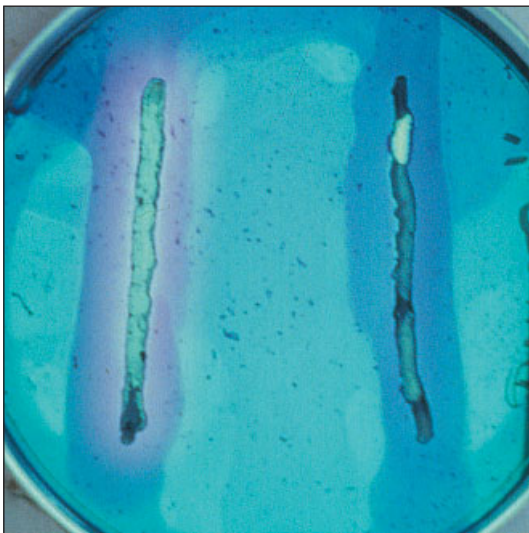

**Figure 63.3 DNase test.** A rose-pink halo around the area of growth on the left side of the plate indicates a positive result, while the absence of a halo on the right is a negative result.

## CLINICAL APPLICATION

### Staphyloxanthin

*Staphylococcus aureus* is one of the most common species of staphylococci to cause human disease, producing many types of skin infections as well as life-threatening diseases like meningitis, osteomyelitis, endocarditis, and toxic shock syndrome. Its pathogenic success is due to its immune-evasive properties, mainly through the production of its yellow pigment staphyloxanthin. This pigment behaves as a virulence factor that helps the organism evade the immune system of the host. Blocking the synthesis of staphyloxanthin may present a unique and vital target for antimicrobials.


## AT THE BENCH



### Materials

#### Cultures

24-hour Trypticase soy agar slant cultures of *Staphylococcus epidermidis*, *Staphylococcus saprophyticus* (ATCC 15305), and *Staphylococcus aureus* (ATCC e 27660). Number-coded, 24-hour blood agar cultures of the above organisms for the STAPH-IDENT system.

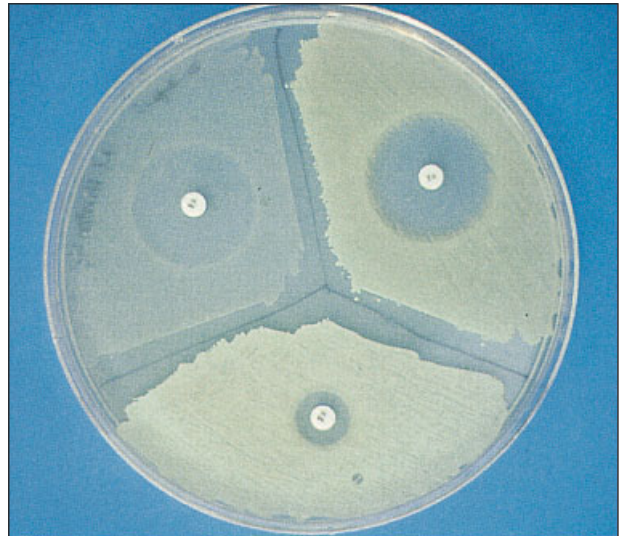

**Figure 63.4 Novobiocin test.** *Staphylococcus aureus* and *Staphylococcus epidermidis* (on top) are sensitive to the antibiotic, while *Staphylococcus saprophyticus* (on bottom) is resistant.

## Media

Per designated student group: three mannitol salt agar plates, one DNA agar plate, three Mueller-Hinton agar plates, and the STAPH-IDENT system.

## Reagents

Citrated human or rabbit plasma, 0.1% toluidine blue, 0.85% saline (pH 5.5–7.0), McFarland barium sulfate standards, and BactiStaph diagnostic kit (latex agglutination test).

## Equipment

Bunsen burner, inoculating loop, 13 × 100-mm test tubes, 15 × 150-mm test tubes, sterile Pasteur pipettes, 1-ml sterile pipettes, mechanical pipetting device, sterile cotton swabs, 30- $\mu$ g novobiocin antibiotic discs, glassware marking pencil, metric ruler, forceps, and beaker with 95% ethyl alcohol.

## Procedure Lab One

### Traditional Procedures

1. Preparation of DNA agar plate culture:
  - a. With a glassware marking pencil, divide the bottom of the plate into three sections. Label each section with the name of the organism to be inoculated.
  - b. Aseptically make a single line of inoculation of each test organism in its respective sector on the agar plate.
2. Preparation of agar plate cultures for novobiocin sensitivity determination:
  - a. Label the three Mueller-Hinton agar plates with the name of the test organism to be inoculated. Inoculate each plate with its respective organism according to the Kirby-Bauer procedure as outlined in Experiment 44.
  - b. Using alcohol-dipped and flamed forceps, aseptically apply a novobiocin antibiotic disc to the surface of each inoculated plate. *Gently* press the discs down with sterile forceps to ensure that they adhere to the agar surface.
3. Preparation of mannitol salt agar plate cultures: Aseptically make a single line of inoculation of each test organism in the center of the appropriately labeled agar plates.
4. Incubation of all plate cultures: Incubate them in an inverted position for 24 to 48 hours at 37°C.

5. Coagulase test procedure:
  - a. Label three 13 × 100-mm test tubes with the name of the organism to be inoculated.
  - b. Aseptically add 0.5 ml of a 1:4 dilution of citrated rabbit or human plasma and 0.1 ml of each test culture to its appropriately labeled test tube.
  - c. Examine the bacterial plasma suspensions for clot formation at 5 minutes, 20 minutes, 1 hour, and 4 hours after inoculation by holding the test tubes in a slanted position. Record your observations and results in the Lab Report.
  - d. At the end of the laboratory session, place all tubes that are coagulase-negative in an incubator for 20 hours at 37°C.

### STAPH-IDENT System Procedure

1. Prepare strip:
  - a. Dispense 5 ml of tap water into incubation tray.
  - b. Place API strip in incubation tray.
2. Prepare inoculum:
  - a. Add 2 ml of 0.85% saline (pH 5.5–7.0) to a sterile 15 × 150-mm test tube.
  - b. Using a sterile swab, pick up a sufficient amount of inoculum to prepare a saline suspension with a final turbidity that is equivalent to a No. 3 McFarland (BaSO<sub>4</sub>) turbidity standard. *Note: Be sure to use suspension within 15 minutes of preparation.*
3. With a sterile Pasteur pipette, add 2 or 3 drops of the inoculum to each microcupule.
4. Place plastic lid on tray and incubate for 5 hours at 37°C.

### Latex Agglutination Procedure

1. Label three of the provided slides (cards) with the name of the organism to be inoculated.
2. Place one drop of *Staphylococcus* latex reagent in the center of the circle on the provided slide.
3. Using an applicator stick or sterile needle, spread one colony of each organism in the reagent of its respective slide.
4. Spread the mixture over the entire circle.
5. Rotate the slide in a circular motion for 60 seconds.
6. Observe all slides for the presence or absence of agglutination. A positive agglutination

reaction usually occurs in 15 seconds and is indicated by a clumping together of the black latex suspension, followed by the loss of the black background. A negative reaction results in little or no agglutination and no loss of the black background within 60 seconds.

- Record your results as positive (+) or as negative (–) in the chart provided in the Lab Report.

## Procedure Lab Two

### Traditional Procedures

- Examine the coagulase-negative tubes, and record your observations in the Lab Report.
- Examine the mannitol salt agar plate. Note and record the following in the Lab Report.
  - Presence (+) or absence (–) of growth of each test organism.
  - Color of the medium surrounding the growth of each test organism.
  - Whether each test organism is a mannitol fermenter (+) or non-mannitol fermenter (–).
- Flood the DNA agar plate with 0.1% toluidine blue. Observe for the delayed development of a

rose-pink coloration surrounding the growth of each test organism. Record your color observation and indicate the presence (+) or absence (–) of DNase activity in the Lab Report.

- With a metric ruler, measure the size of the zone of inhibition, if present, surrounding each of the novobiocin discs on the agar plates. A zone of inhibition of 17 mm or less is indicative of novobiocin resistance, whereas a zone greater than 17 mm indicates that the organism is sensitive to this antibiotic. Record the susceptibility of each test organism to novobiocin as sensitive (S) or resistant (R) in the Lab Report.

### STAPH-IDENT System Procedure

- Interpret your STAPH-IDENT system reactions on the basis of the observed color changes in each of the microcupules described in the chart in the Lab Report. Report your color observations and results as (+) or (–) for each test in the Lab Report.
- Construct a four-digit profile for your unknown organisms using the guidelines provided in the Lab Report.

**TABLE 63.2 API STAPH-IDENT Profile Register**

PROFILE	IDENTIFICATION	PROFILE	IDENTIFICATION
0 040	STAPH CAPITIS	2 000	STAPH SAPROPHYTICUS NOVO R
0 060	STAPH HAEMOLYTICUS		STAPH HOMINIS NOVO S
0 100	STAPH CAPITIS	2 001	STAPH SAPROPHYTICUS
0 140	STAPH CAPITIS	2 040	STAPH SAPROPHYTICUS NOVO R
0 200	STAPH COHNII		STAPH HOMINIS NOVO S
0 240	STAPH CAPITIS	2 041	STAPH SIMULANS
0 300	STAPH CAPITIS	2 061	STAPH SIMULANS
0 340	STAPH CAPITIS	2 141	STAPH SIMULANS
0 440	STAPH HAEMOLYTICUS	2 161	STAPH SIMULANS
0 460	STAPH HAEMOLYTICUS	2 201	STAPH SAPROPHYTICUS
0 600	STAPH COHNII	2 241	STAPH SIMULANS
0 620	STAPH HAEMOLYTICUS	2 261	STAPH SIMULANS
0 640	STAPH HAEMOLYTICUS	2 341	STAPH SIMULANS
0 660	STAPH HAEMOLYTICUS	2 361	STAPH SIMULANS
1 000	STAPH EPIDERMIDIS	2 400	STAPH HOMINIS NOVO S
1 040	STAPH EPIDERMIDIS		STAPH SAPROPHYTICUS NOVO R
1 300	STAPH AUREUS	2 401	STAPH SAPROPHYTICUS
1 540	STAPH HYICUS (An)	2 421	STAPH SIMULANS
1 560	STAPH HYICUS (An)	2 441	STAPH SIMULANS
2 541	STAPH SIMULANS	2 461	STAPH SIMULANS
2 561	STAPH SIMULANS	6 101	STAPH XYLOSUS
		6 121	STAPH XYLOSUS





**TABLE 63.2 API STAPH-IDENT Profile Register (continued)**

PROFILE	IDENTIFICATION		PROFILE	IDENTIFICATION	
2 601	STAPH SAPROPHYTICUS		6 221	STAPH XYLOSUS	
2 611	STAPH SAPROPHYTICUS		6 300	STAPH AUREUS	
2 661	STAPH SIMULANS		6 301	STAPH XYLOSUS	
2 721	STAPH COHNII (SSP1)		6 311	STAPH XYLOSUS	
2 741	STAPH SIMULANS		6 321	STAPH XYLOSUS	
2 761	STAPH SIMULANS		6 340	STAPH AUREUS	COAG+
				STAPH WARNERI	COAG-
3 000	STAPH EPIDERMIDIS		6 400	STAPH WARNERI	
3 040	STAPH EPIDERMIDIS		6 401	STAPH XYLOSUS	XYL + ARA+
3 140	STAPH EPIDERMIDIS			STAPH SAPROPHYTICUS	XYL - ARA-
3 540	STAPH HYICUS (An)		6 421	STAPH XYLOSUS	
3 541	STAPH INTERMEDIUS (An)		6 460	STAPH WARNERI	
3 560	STAPH HYICUS (An)		6 501	STAPH XYLOSUS	
3 601	STAPH SIMULANS	NOVO	6 521	STAPH XYLOSUS	
	STAPH SAPROPHYTICUS	NOVO R	6 600	STAPH WARNERI	
			6 601	STAPH SAPROPHYTICUS	XYL - ARA-
4 060	STAPH HAEMOLYTICUS			STAPH XYLOSUS	XYL + ARA+
4 210	STAPH SCIURI		6 611	STAPH XYLOSUS	
4 310	STAPH SCIURI		6 621	STAPH XYLOSUS	
4 420	STAPH HAEMOLYTICUS		6 700	STAPH AUREUS	
4 440	STAPH HAEMOLYTICUS		6 701	STAPH XYLOSUS	
4 460	STAPH HAEMOLYTICUS		6 721	STAPH XYLOSUS	
4 610	STAPH SCIURI		6 731	STAPH XYLOSUS	
4 620	STAPH HAEMOLYTICUS				
4 660	STAPH HAEMOLYTICUS		7 000	STAPH EPIDERMIDIS	
4 700	STAPH AUREUS	COAG+	7 021	STAPH XYLOSUS	
	STAPH SCIURI	COAG-	7 040	STAPH EPIDERMIDIS	
4 710	STAPH SCIURI		7 141	STAPH INTERMEDIUS (An)	
			7 300	STAPH AUREUS	
5 040	STAPH EPIDERMIDIS		7 321	STAPH XYLOSUS	
5 200	STAPH SCIURI		7 340	STAPH AUREUS	COAG-
5 210	STAPH SCIURI	COAG+	7 401	STAPH XYLOSUS	
5 300	STAPH AUREUS	COAG-	7 421	STAPH XYLOSUS	
	STAPH SCIURI		7 501	STAPH INTERMEDIUS (An)	COAG+
5 310	STAPH SCIURI			STAPH XYLOSUS	COAG-
5 600	STAPH SCIURI		7 521	STAPH XYLOSUS	
5 610	STAPH SCIURI	COAG+	7 541	STAPH INTERMEDIUS (An)	
5 700	STAPH AUREUS	COAG-	7 560	STAPH HYICUS (An)	
	STAPH SCIURI		7 601	STAPH XYLOSUS	
5 710	STAPH SCIURI		7 621	STAPH XYLOSUS	
5 740	STAPH AUREUS		7 631	STAPH XYLOSUS	
			7 700	STAPH AUREUS	
6 001	STAPH XYLOSUS	XYL + ARA+	7 701	STAPH XYLOSUS	
	STAPH SAPROPHYTICUS	XYL - ARA-	7 721	STAPH XYLOSUS	
6 011	STAPH XYLOSUS		7 740	STAPH AUREUS	
6 021	STAPH XYLOSUS				

*Note: The API STAPH-IDENT Profile Register uses "STAPH" as an abbreviation for Staphylococcus. The correct scientific abbreviation is S. aureus, for example.*

*Source: Staph-Ident, Analytab Products, Division of Sherwood Medical, Plainview, New York.*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

### Traditional Procedures

Staphylococcal Species	APPEARANCE OF PLASMA: CLOTTED (+) OR UNCLOTTED (-)					Coagulase (+) or (-)
	5 min	20 min	1 hr	4 hr	24 hr	
<i>S. aureus</i>						
<i>S. epidermidis</i>						
<i>S. saprophyticus</i>						

Procedure	<i>S. aureus</i>	<i>S. epidermidis</i>	<i>S. saprophyticus</i>
Mannitol salt agar:			
Growth	_____	_____	_____
Color of medium	_____	_____	_____
Fermentation	_____	_____	_____
DNA agar:			
Color of medium	_____	_____	_____
DNase activity	_____	_____	_____
Novobiocin sensitivity:			
Growth inhibition in mm	_____	_____	_____
Susceptibility—(R) or (S)	_____	_____	_____

### STAPH-IDENT System Procedure

MICROCUPULE		INTERPRETATION OF REACTIONS		REACTION RESULTS	
No.	Substrate	Positive	Negative	Color	(+) or (-)
1	PHS <i>p</i> -Nitrophenyl-phosphate, disodium salt	Yellow	Clear or straw-colored		
2	URE Urea	Purple to red-orange	Yellow or yellow-orange		
3	GLS <i>p</i> -Nitrophenyl- $\beta$ -D-glucopyranoside	Yellow	Clear or straw-colored		



## STAPH-IDENT System Procedure (continued)

MICROCUPULE		INTERPRETATION OF REACTIONS		REACTION RESULTS	
No.	Substrate	Positive	Negative	Color	(+) or (-)
4	MNE Mannose	Yellow or yellow-orange	Red or orange		
5	MAN Mannitol				
6	TRE Trehalose				
7	SAL Salicin				
8	GLC <i>p</i> -Nitrophenyl- $\beta$ -D-glucuronide	Yellow	Clear or straw-colored		
9	ARG Arginine	Purple to red-orange	Yellow or yellow-orange		
10	NGP 2-Naphthyl- $\beta$ -D-galactopyranoside	Add 1–2 drops of STAPH-IDENT reagent			
		Plum-purple (mauve)	Yellow or colorless		

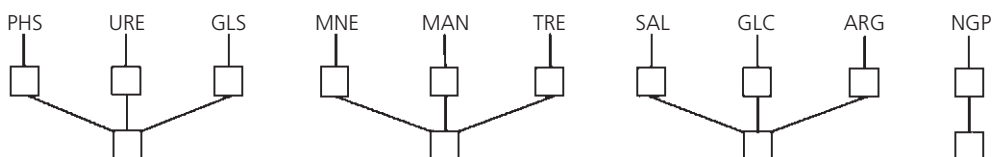
Construct a four-digit profile for your unknown organism as follows: A four-digit profile is derived from the results obtained with STAPH-IDENT. The 10 biochemical tests are divided into four groups, as follows:

PHS	MNE	SAL	NGP
URE	MAN	GLC	
GLS	TRE	ARG	

Only positive reactions are assigned a numerical value. The value depends on the location within the group.

- A value of 1 for the first biochemical in each group (e.g., PHS, MNE)
- A value of 2 for the second biochemical in each group (e.g., URE, MAN)
- A value of 4 for the third biochemical in each group (e.g., GLS, TRE)
- A value of 0 for all negative reactions

A four-digit number is obtained by totaling the values of each of the groups.



Using Table 63.2 on pages 435–436 and your four-digit profile number, identify your organism.

**Unknown organism:** \_\_\_\_\_

## Latex Agglutination Procedure

Record the presence of agglutination as (+), and the absence of agglutination as (–).

	<i>S. aureus</i>	<i>S. epidermidis</i>	<i>S. saprophyticus</i>
Agglutination			
No agglutination			

## LEARNING OBJECTIVES

Once you have completed this experiment, you should understand

1. The medical significance of streptococci.
2. Selected laboratory procedures designed to differentiate streptococci on the basis of their hemolytic activity and biochemical patterns associated with the Lancefield group classifications.

## Principle

Members of the genus *Streptococcus* are perhaps responsible for a greater number of infectious diseases than any other group of microorganisms. Morphologically, they are cocci that divide in a single plane, forming chains. They form circular, translucent to opaque, pinpoint colonies on solid media. All members of this group are gram-positive, and many are nutritionally fastidious, requiring enriched media such as blood for growth.

The streptococci are classified by means of two major methods: (1) their **hemolytic activity**, and (2) the **serologic classification of Lancefield**. The observed hemolytic reactions on blood agar are of the following three types:

1. ( $\alpha$ ) **Alpha hemolysis**, an incomplete form of hemolysis, produces a green zone around the colony.  $\alpha$ -Hemolytic streptococci, the *Streptococcus viridans* species, are usually nonpathogenic opportunists. In some instances, however, they are capable of inducing human infections such as **subacute endocarditis**, which may precipitate valvular damage and heart failure if untreated. *Streptococcus pneumoniae*, the causative agent of **lobar pneumonia**, will be studied in a separate experiment.
2. ( $\beta$ ) **Beta hemolysis**, a complete destruction of red blood cells, exhibits a clear zone of approximately 2 to 4 times the diameter of the colony. The streptococci capable of producing  $\beta$ -hemolysins are most frequently associated with pathogenicity.
3. ( $\gamma$ ) **Gamma hemolysis** is indicative of the absence of any hemolysis around the colony. Most commonly,  $\gamma$ -hemolytic streptococci are avirulent.

These hemolytic reactions are shown in

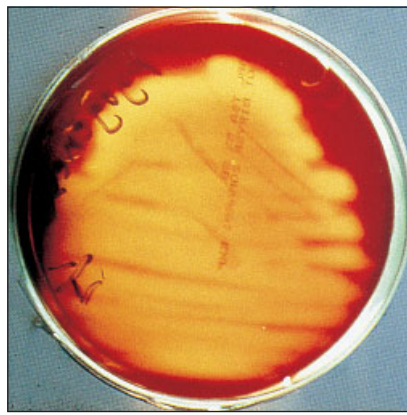
### Figure 64.1.

Lancefield classified the streptococci into 20 **serogroups**, designated **A** through **V**, omitting **I** and **J**, based on the presence of an antigenic group-specific hapten called the **C-substance**. This method of classification generally implicates the members of Groups A, B, C, and D in human infectious processes.

$\beta$ -Hemolytic streptococci belonging to **Group A**, and collectively referred to as *Streptococcus pyogenes*, are the human pathogens of prime importance. Members of this group are the main etiological agents of human respiratory infections such as **tonsillitis**, **bronchopneumonia**, and **scarlet fever**, as well as skin disorders such as **erysipelas** and **cellulitis**. In addition, these organisms are responsible for the development of complicating infections, namely **glomerulonephritis** and **rheumatic fever**, which may surface when primary streptococcal infections either go untreated or are not completely eradicated by antibiotics. The  $\beta$ -hemolytic streptococci found in **Group B** are indigenous to the vaginal mucosa and have been shown to be responsible for **puerperal fever** (childbirth fever), a sometimes-fatal **neonatal meningitis**, and **endocarditis**. Members of **Group C** are also  $\beta$ -hemolytic and have been implicated in **erysipelas**, **puerperal fever**, and **throat infections**. The enterococci formerly classified as Group D streptococci have been reclassified and are now considered a separate genus. The enterococci differ significantly from other members of Group D, such as *S. bovis*, which may be the etiological agent of urinary tract infections. Enterococci such as *Enterococcus faecalis* may cause infections to the lungs, urinary



(a) Alpha hemolysis



(b) Beta hemolysis



(c) Gamma hemolysis

**Figure 64.1** Types of hemolytic reactions on blood agar

tract infections, or bloodstream through an intestinal laceration or poor personal hygiene. The enterococci tend to be antibiotic-resistant, particularly to penicillin and more recently to vancomycin.

The virulence of the streptococci is associated with their ability to produce a wide variety of extracellular metabolites. Included among these are the **hemolysins** ( $\alpha$  and  $\beta$ ), **leukocidins** that destroy phagocytes, and the **erythrogenic toxin** responsible for the rash of scarlet fever. Also of medical significance are three metabolic end products that facilitate the spread of the organisms, thereby initiating secondary sites of streptococcal infection. These metabolites are **hyaluronidase** (the spreading factor), which hydrolyzes the tissue cement hyaluronic acid; **streptokinase**, a

fibrinolysin; and the **nucleases**, ribonuclease and deoxyribonuclease, which destroy viscous tissue debris.

Although the different groups of streptococci have similar colonial morphology and microscopic appearance, they can be separated and identified by the performance of a variety of laboratory tests. Toward this end, you will perform laboratory procedures to differentiate among the medically significant streptococci on the basis of their Lancefield group classification and their hemolytic patterns. Table 64.1 will aid in this separation.


Identification of Group A streptococci involves the following procedures:

1. **Bacitracin sensitivity test:** A filter-paper disc impregnated with 0.04 unit of bacitracin

**TABLE 64.1** Laboratory Differentiation of Streptococci

GROUP:	A	B	C	D	K, H, N	
ORGANISMS:	<i>S. pyogenes</i>	<i>S. agalactiae</i>	<i>S. equi</i>	<i>S. bovis</i> NON-ENTEROCOCCI	<i>S. salivarius</i> <i>S. sanguis</i> <i>S. mitis</i>	<i>E. faecalis</i> ENTEROCOCCI
Hemolysis	$\beta$	$\beta$	$\beta$	$\alpha \rightarrow \gamma$	$\alpha$	$\alpha \rightarrow \gamma$
Bacitracin sensitivity	S	R	R	R	R	R
CAMP test	—	+	—	—	—	—
Bile esculin hydrolysis	—	—	—	+	—	+
6.5% NaCl medium	NG	NG	NG	NG	NG	G
Growth at 10°C	NG	NG	NG	NG	NG	G
Growth at 45°C	NG	NG	NG	NG or G	NG	G

NG = no growth; G = growth; S = sensitive; R = resistant

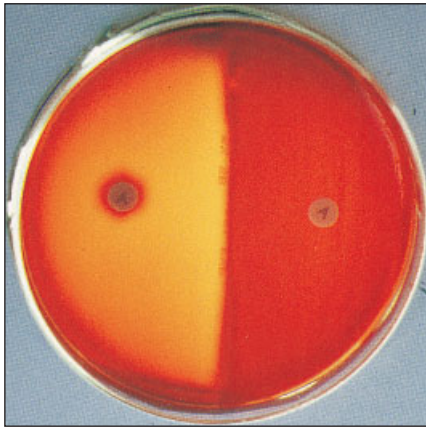

**Figure 64.2 Bacitracin sensitivity test.** Positive for beta-hemolytic Group A streptococci on the left; negative on the right.


is applied to the surface of a blood agar plate previously streaked with the organism to be identified. Following incubation, the appearance of a zone of growth inhibition surrounding the disc is indicative of Group A streptococci. Absence of this zone suggests a non-Group A organism. **Figure 64.2** shows the result of a bacitracin sensitivity test.

2. **Directigen™ test:** A rapid, non-growth-dependent immunological procedure for the detection of the Group A antigen, developed by Becton Dickinson and Company. In this test, a clinical specimen is subjected to reagents designed to extract the Group A antigen, which is then mixed with a reactive and a negative control latex. Agglutination with the reactive latex is indicative of Group A streptococci.


Group B streptococci are identified with the **CAMP test** (named for Christie, Atkins, and Munch-Petersen). Group B streptococci produce a peptide, the CAMP substance, that acts in concert with the  $\beta$ -hemolysins produced by some strains of *Staphylococcus aureus*, causing an increased hemolytic effect. Following inoculation and incubation, the resultant effect appears as an arrow-shaped zone of hemolysis adjacent to the central streak of *S. aureus* growth. The non-Group B streptococci do not produce this reaction. **Figure 64.3** illustrates the CAMP reactions.

Identification of Group D streptococci involves the following:

1. **Bile esculin test:** In the presence of bile, Group D streptococci hydrolyze the glycoside esculin to 6,7-dihydroxy coumarin that reacts

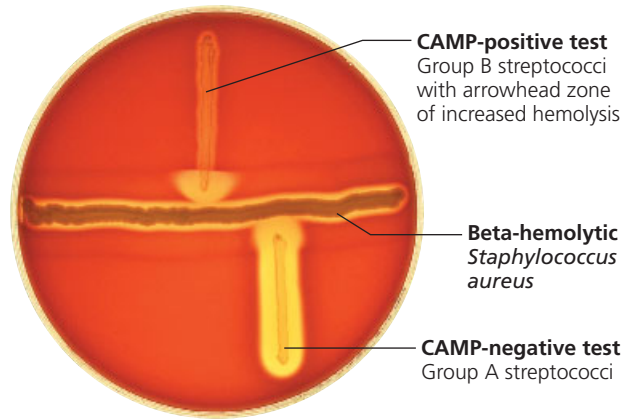

**Figure 64.3 CAMP reactions**

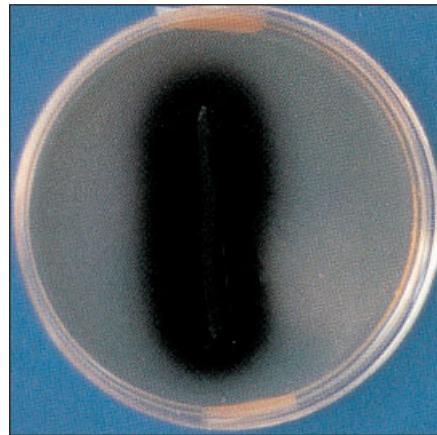

**Figure 64.4 Positive bile esculin test.** A brown-to-black coloration of the medium indicates positive identification of Group D streptococci.

with the iron salts in the medium to produce a brown-to-black coloration of the medium following incubation (**Figure 64.4**). Lack of this dark coloration is indicative of a non-Group D organism.


2. **6.5% sodium chloride broth:** The enterococci can be separated from the non-enterococci by the ability of the former to grow in this medium. This reaction is shown in **Figure 64.5**.

Hemolytic activity is identified with a blood agar medium. The pathogenic streptococci, primarily the  $\beta$ -hemolytic, can be separated from the generally avirulent  $\alpha$ - and  $\gamma$ -hemolytic streptococci by the type of hemolysis produced on blood agar, as previously described.

## CLINICAL APPLICATION
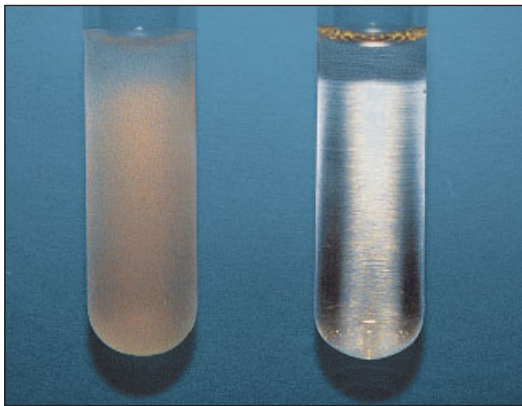
### Streptococci Infections

Medically the streptococci are of significant importance because they are responsible for a wide variety of infections, many of which are pyogenic (pus-producing). *Streptococcus agalactiae* (Lancefield group B) may colonize the vagina as well as the upper respiratory tract of humans, and is the most frequent cause of neonatal pneumonia in the United States. Meanwhile, *Streptococcus pyogenes* (Lancefield group A) causes necrotizing fasciitis, a rare but devastating infection that destroys skin, muscle, and underlying tissue. The CAMP test is used to identify Group A *Streptococcus pyogenes* from Group B *Streptococcus agalactiae*.



(a)

(b)

**Figure 64.5** 65% sodium chloride test. (a) Growth indicates the presence of Group D enterococci. (b) The absence of growth indicates the presence of Group D non-enterococci.

## AT THE BENCH



## Materials

### Cultures

24-hour blood agar slant cultures of *Streptococcus pyogenes* (ATCC 12385), *Enterococcus faecalis*, *Streptococcus bovis*, *Streptococcus agalactiae*, *Streptococcus mitis*, and *Staphylococcus aureus* (ATCC 25923).

## Media

Per designated student group: five blood agar plates, three bile esculin agar plates, and three 6.5% sodium chloride broths.

## Reagents

Directigen Rapid Group A Strep Test (Becton Dickinson and Company), crystal violet, Gram's iodine, ethyl alcohol, safranin, and Taxo™ A discs (0.04 unit of bacitracin).

## Equipment

Bunsen burner, inoculating loop, staining tray, lens paper, bibulous paper, microscope, sterile cotton swabs, glassware marking pencil, sterile 12 × 75-mm test tubes, sterile Pasteur pipettes, sterile applicators, 95% ethyl alcohol in beaker, forceps, and mechanical rotator.

## Procedure Lab One

1. Prepare a Gram-stained preparation of each streptococcal culture and observe under oil immersion. Record in the Lab Report your observations of cell morphology and Gram reaction.
2. Prepare the blood agar plate cultures to identify the type of hemolysis as follows:
  - a. With a glassware marking pencil, divide the bottoms of two blood agar plates to accommodate the five test organisms. Label each section with the name of the culture to be inoculated.
  - b. Using aseptic inoculating technique, make a single line streak of inoculation of each organism in its respective sector on the blood plates.
3. Prepare the blood agar plate cultures for the **bacitracin** test as follows:
  - a. With a glassware marking pencil, label the covers of two blood agar plates with the names of the organisms to be inoculated, *S. pyogenes* and *S. agalactiae*.
  - b. Using a sterile cotton swab, inoculate the agar surface of each plate with its respective test organism by streaking first in a horizontal direction, then vertically to ensure a heavy growth over the entire surface.
  - c. Using alcohol-dipped and flamed forceps, apply a single 0.04-unit bacitracin disc to

- the surface of each plate. Gently touch each disc to ensure its adherence to the agar surface.
4. Prepare a blood agar plate culture for the **CAMP test** as follows:
    - a. Using a sterile inoculating loop, make a single line of inoculation along the center of the plate using the *S. aureus* culture.
    - b. With a sterile loop, inoculate *S. pyogenes* on one side and perpendicular to the central *S. aureus* streak, starting about 5 mm from the central streak and extending toward the periphery of the agar plate.
    - c. On the opposite side of the central streak, but not directly opposite the *S. pyogenes* line of inoculation, repeat Step 4b using *S. agalactiae*.
  5. Prepare the bile esculin agar plate cultures as follows:
    - a. Label the three bile esculin plates with the names of the organisms to be inoculated, *S. bovis*, *S. mitis*, and *E. faecalis*.
    - b. Aseptically inoculate each plate with its test organism by making several lines of inoculation on the agar surface.
  6. Prepare 6.5% sodium chloride broth cultures as follows:
    - a. Label three tubes of 6.5% sodium chloride broth with the names of the organisms to be inoculated, *S. bovis*, *E. faecalis*, and *S. mitis*.
    - b. With a sterile loop, inoculate each tube with its organism.
  7. Conduct the **Directigen test** procedure as follows:
    - a. Label two sterile 12 × 75-mm test tubes as *S. pyogenes* and *S. agalactiae*.
    - b. Add 0.3 ml of Reagent 1 to both test tubes.
    - c. Using a sterile cotton swab, transfer the test organisms into their respectively labeled test tubes. *Note: These samples will emulate the throat swabs obtained in a clinical solution.*
    - d. Add 1 drop of Reagent 2 to each test tube. Mix by rotating the swab against the side of the tube. Allow the swabs to remain in the test tubes for 3 minutes.
    - e. Add 1 drop of Reagent 3 to both tubes and mix.
    - f. Remove swabs after extracting as much liquid as possible by rolling them against the sides of the tubes.
    - g. Place 1 drop of negative antigen control on both circles in Column A of test slide.
    - h. Place 1 drop of positive antigen control on both circles in Column B of test slide.
    - i. Dispense 1 drop of each streptococcal sample on both circles in Columns C and D, respectively.
    - j. Using a new sterile applicator for each specimen, spread each specimen within the confines of both circles in Columns A, B, C, and D.
    - k. Add 1 drop of reactive latex to the top row of circles.
    - l. Add 1 drop of control latex to the bottom row of circles.
    - m. Place the slide on a mechanical rotator for 4 minutes under a moistened humidifying cover.
    - n. Compare the agglutination seen in the upper “reactive latex” circles with the consistency of the latex in the bottom “control latex” circles. Any agglutination in the top circles distinct from any background granules seen in the bottom circles indicates Group A streptococci.
  8. Incubate all tubes and plates in an inverted position for 24 hours at 37°C.

## Procedure Lab Two

1. Examine the two blood agar plates for bacitracin activity. Record in the Lab Report your observations of the presence (+) or absence (–) of a zone of inhibition of any size surrounding the discs.
2. Examine the blood agar plate for the CAMP reaction. Record your observations of the presence (+) or absence (–) of increased arrow-shaped hemolysis.
3. Examine the bile esculin plates for the presence (+) or absence (–) of a brown-black coloration in the medium and record your observations.



4. Observe the 6.5% sodium chloride broth cultures for the presence (+) or absence (–) of growth and record your observations.
5. Examine the two blood agar plates for the presence and type of hemolysis produced by each of the test organisms. Record your observations of the appearance of the medium surrounding the growth and the type of hemolytic reaction that has occurred— $\alpha$ ,  $\beta$ , or  $\gamma$ .
6. Observe the Directigen test slide for the presence (+) or absence (–) of agglutination in the reactive and control latex circles. Based on your observations, indicate the Lancefield group classification of each test organism. Record your results.
7. Based on your observations, classify each test organism according to its Lancefield group.
8. Check that all of your observations have been recorded in the Lab Report.

Name: \_\_\_\_\_



Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

Procedure	<i>S. pyogenes</i>	<i>S. agalactiae</i>	<i>S. bovis</i>	<i>E. faecalis</i>	<i>S. mitis</i>
Gram stain: Morphology	_____	_____	_____	_____	_____
Reaction	_____	_____	_____	_____	_____
Bacitracin test: Zone of inhibition	_____	_____	_____	_____	_____
CAMP test: Increased hemolysis	_____	_____	_____	_____	_____
Bile esculin test: Color of medium	_____	_____	_____	_____	_____
Result: (+) or (-)	_____	_____	_____	_____	_____
6.5% NaCl broth: Growth	_____	_____	_____	_____	_____
Hemolytic activity: Appearance of medium	_____	_____	_____	_____	_____
Type of hemolysis	_____	_____	_____	_____	_____
Directigen test: Agglutination (+) or (-) in:					
Reactive circle	_____	_____	_____	_____	_____
Control circle	_____	_____	_____	_____	_____
Lancefield group	_____	_____	_____	_____	_____
Group classification					

## Review Questions

1. How do the purposes of the bacitracin and CAMP tests differ?

2. What is the mechanism of the bile esculin test?
  
  
  
  
  
  
  
  
  
  
3. Why is it important medically to distinguish between the enterococci and the non-enterococci?
  
  
  
  
  
  
  
  
  
  
4.  Why can some streptococci produce secondary sites of infection?
  
  
  
  
  
  
  
  
  
  
5.  The streptococci are known to be fastidious organisms that require an enriched medium for growth. How would you account for the fact that a medium enriched with blood (blood agar) is the medium of preference for growth of these organisms?

## LEARNING OBJECTIVE

Once you have completed this experiment, you should understand

1. Laboratory procedures to differentiate between *Streptococcus pneumoniae* and other  $\alpha$ -hemolytic streptococci.


## Principle

The pneumococcus *Streptococcus pneumoniae* is the major  $\alpha$ -hemolytic, streptococcal pathogen in humans. It serves as an etiological agent of **lobar pneumonia**, an infection characterized by acute inflammation of the bronchial and alveolar membranes. These organisms are gram-positive cocci, tapered or lancet-shaped at their edges, that occur in pairs or as short, tight chains. The large, thick **capsules** formed *in vivo* are responsible for antiphagocytic activity, which is believed to enhance the organisms' virulence. In addition, the pneumococci produce  **$\alpha$ -hemolysis** on blood agar plates. Because of these properties (short-chain formation,  $\alpha$ -hemolysis, and failure of the capsule to stain on Gram staining), the organisms closely resemble *Streptococcus viridans* species. The *S. pneumoniae* can be differentiated from other  $\alpha$ -hemolytic streptococci on the basis of the following laboratory tests:

Test	<i>S. pneumoniae</i>	<i>S. mitis</i>
Hemolysis	$\alpha$	$\alpha$
Bile solubility	+	-
Optochin sensitivity	+	-
Inulin fermentation	+	-
Quellung reaction	+	-
Mouse virulence	+	-

Brief descriptions of the tests and their mechanisms follow:

1. **Bile solubility test:** In the presence of surface-active agents such as **bile** and **bile salts** (sodium desoxycholate or sodium dodecyl sulfate), the cell wall of the pneumococcus undergoes lysis. Other members of the  $\alpha$ -hemolytic streptococci will not be lysed by these agents and are bile-insoluble. Following incubation, bile-soluble cultures will appear clear, and bile-insoluble cultures will be turbid.
2. **Optochin sensitivity test:** This is a growth inhibition test in which 6-mm filter-paper discs impregnated with 5 mg of **ethylhydrocupreine hydrochloride** (optochin) and called Taxo P discs are applied to the surface of a blood agar plate streaked with the test organisms. The *S. pneumoniae*, being sensitive to this surface-active agent, are lysed with the resultant formation of a zone of inhibition greater than 15 mm surrounding the P disc.

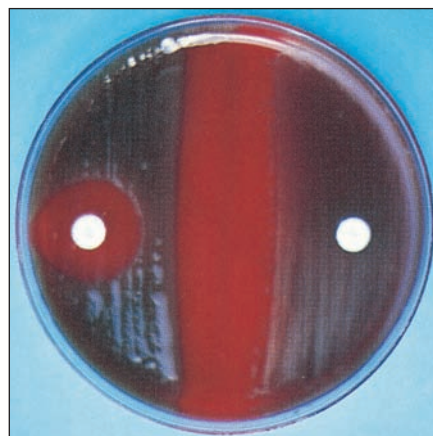

**Figure 65.1** Optochin sensitivity test. The formation of a zone of inhibition greater than 15 mm on the left indicates the presence of alpha-hemolytic *S. pneumoniae*. No zone of inhibition on the right indicates the presence of other alpha-hemolytic streptococcal species.

Nonpneumococcal  $\alpha$ -hemolytic streptococci are resistant to optochin and fail to show a zone of inhibition or produce a zone less than 15 mm. Sensitivity to optochin is illustrated in **Figure 65.1**.

- Inulin fermentation:** The pneumococci are capable of fermenting inulin, while most other  $\alpha$ -hemolytic streptococci are non-inulin fermenters. Following incubation, the **acid** resulting from inulin fermentation will change the color of the culture from red to yellow. Cultures that are not capable of fermenting inulin will not exhibit a color change, which is a negative test result.
- Quellung (Neufeld) reaction:** This **capsular swelling** reaction is a sensitive and accurate method of determining the presence of *S. pneumoniae* in sputum. The reaction of the pneumococcal capsular polysaccharide, a hapten antigen, with an omnivalent capsular antiserum (Abcam, Inc.) produces a microscopically visible swollen capsule surrounding the *S. pneumoniae* organisms.
- Mouse virulence test:** Laboratory white mice are highly susceptible to infection by *S. pneumoniae* and resistant to other streptococcal infections. Intraperitoneal injection of 0.1 ml of pneumococcus-infected sputum will kill the mouse. Examination of the peritoneal fluid by Gram stain and culture will reveal the presence of *S. pneumoniae*.

In the following experiment, you will use hemolytic patterns, bile solubility, the Quellung reaction, the optochin test, and the inulin fermentation test for laboratory differentiation of *S. pneumoniae* from other  $\alpha$ -hemolytic streptococci.

## CLINICAL APPLICATION

### Pneumococcus Infections

*Streptococcus pneumoniae*, formerly called *Diplococcus pneumoniae*, appears as a lancet-shaped diplococcus and is unlike all other cocci. The pneumococcus, as it is called, is the causative agent of lobar pneumonia (lung), otitis media (middle ear), and meningitis (meninges) infections. It is currently the leading invasive bacterial disease in children and the elderly. Presently a vaccine is available for people who are designated as high risk for infection with this organism.

## AT THE BENCH



## Materials

### Cultures

24-hour blood agar slant cultures of *Streptococcus pneumoniae* and *Streptococcus mitis*.

### Media

Per designated student group: one blood agar plate, two phenol red inulin broth tubes, and four 13 × 75-mm tubes containing 1 ml of nutrient broth.

### Reagents

Crystal violet, Gram's iodine, ethyl alcohol, safranin, methylene blue, 10% sodium desoxycholate, commercially available Taxo P discs (5 mg of optochin), and omnivalent pneumococcal antiserum.

## Equipment

Bunsen burner, waterbath, inoculating loop, glass slides, coverslips, sterile cotton swabs, sterile 1-ml serological pipettes, mechanical pipetting device, 95% ethyl alcohol in beaker, forceps, and glassware marking pencil.

## Procedure Lab One

### 1. Bile solubility test

- a. Label two nutrient broth tubes *S. pneumoniae* and two other tubes *S. mitis*.
- b. Aseptically add 2 loopfuls of the test organisms to the appropriately labeled sterile test tubes to effect a heavy suspension.
- c. Aseptically add 0.5 ml of sodium desoxycholate to one tube of each test culture. The remaining two cultures will serve as controls.
- d. Incubate the tubes in a waterbath at 37°C for 1 hour.
- e. After incubation, examine the tubes for the presence or absence of turbidity in each culture. Record your observations of the appearance (clear or turbid) and bile solubility of each test organism in the Lab Report.

### 2. Optochin test

- a. With a glassware marking pencil, divide the bottom of a blood agar plate into two equal sections and label one section *S. pneumoniae* and the other *S. mitis*.
- b. Using a sterile cotton swab, heavily inoculate the surface of each section with its respective test organism in a horizontal and then vertical direction, being careful to stay within the limits of each section.
- c. Using alcohol-dipped and flamed forceps, apply a single Taxo P disc (optochin) to the surface of the agar in each section of the inoculated plate. Touch each disc slightly to ensure its adherence to the agar surface.
- d. Incubate the plate in an inverted position for 24 to 48 hours at 37°C.

### 3. Inulin fermentation test

- a. Label two phenol red inulin broth tubes with the name of each test organism to be inoculated.
- b. Using aseptic technique and loop inoculation, inoculate each experimental organism in its appropriately labeled tube of medium.
- c. Incubate the tube cultures for 24 to 48 hours at 37°C.

### 4. Quellung reaction

- a. Spread a loopful of each test culture on a separate labeled clean glass slide and allow the slides to air-dry.
- b. Place a loopful of the omnivalent capsular antiserum and a loopful of methylene blue on each of two coverslips.
- c. Place the coverslips over the dried bacterial smears. Prepare a Gram-stained preparation of each test organism and observe under oil immersion. Record your observations of cell morphology and Gram reaction in the Lab Report.

## Procedure Lab Two

1. Examine blood agar plates for the presence of hemolysis and optochin activity by measuring the zone of inhibition, if any, surrounding the disc. Record the measurement in the Lab Report and indicate whether each organism is optochin-sensitive (zone of inhibition greater than 15 mm) or optochin-resistant (no zone or less than 15 mm).
2. Observe the inulin fermentation broth cultures containing phenol red and record the color of each culture and whether it is indicative of a positive (+) or negative (-) result in the Lab Report.
3. Examine slides of the Quellung reaction under oil immersion and indicate in the Lab Report the presence (+) or absence (-) of capsular swelling surrounding the blue-stained cells.

*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results


Procedure	<i>S. pneumoniae</i>	<i>S. mitis</i>
Bile solubility test: Appearance of culture	_____	_____
Bile solubility	_____	_____
Gram stain: Morphology	_____	_____
Reaction	_____	_____
Optochin test: Zone of inhibition in mm	_____	_____
Resistant or sensitive	_____	_____
Inulin fermentation: Color of medium	_____	_____
Fermentation (+) or (-)	_____	_____
Quellung reaction: Capsular swelling (+) or (-)	_____	_____


## Review Questions

1. Why is it clinically important to distinguish *S. pneumoniae* from other  $\alpha$ -hemolytic streptococci?



2. How would you separate *S. pneumoniae* from other  $\alpha$ -hemolytic streptococci?

3.  What are secondary pneumonias? Why do they develop most frequently following viral infections?

4.  Why did it require many years of research to develop an effective, long-term pneumococcal vaccine?

# Identification of Enteric Microorganisms Using Computer-Assisted Multitest Microsystems

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be familiar with

1. The members of the family Enterobacteriaceae.
2. Laboratory procedures designed to identify enteric pathogens using commercial multitest microsystems.

## Principle

The Enterobacteriaceae are a significant group of bacteria that are endogenous to the intestinal tract or that may gain access to this site via a host's ingestion of contaminated food and water. The family consists of a number of genera whose members vary in their capacity to produce disease. The *Salmonella* and *Shigella* are considered to be pathogenic. Members of other genera, particularly *Escherichia* and *Enterobacter*, and to a lesser extent *Klebsiella* and *Proteus*, constitute the natural flora of the intestines and are generally considered to be avirulent. Remember, however, that all can produce disease under appropriate conditions.

The Enterobacteriaceae are gram-negative, short rods. They are mesophilic, nonfastidious organisms that multiply in many foods and water sources. They are all non-spore-formers and susceptible to destruction by common physical and chemical agents. They are resistant to destruction by low temperatures and can therefore frequently survive in soil, sewage, water, and many foods for extended periods.


From a medical point of view, the pathogenic Enterobacteriaceae are salmonellae and shigellae. Salmonellae are responsible for enteric fevers, **typhoid**, the milder **paratyphoid**, and **gastroenteritis**. In typhoid, *Salmonella typhi* penetrates the intestinal mucosa and enters the bloodstream, thus infecting organs such as the gallbladder, intestines, liver, kidney, spleen, and heart. Ulceration

of the intestinal wall, caused by the release of the lipopolysaccharide endotoxin into the blood over a long febrile period, and enteric symptoms are common. **Gastroenteritis** is caused by a number of *Salmonella* species. Symptoms associated with this type of food poisoning include abdominal pain, nausea, vomiting, and diarrhea, which develop within 24 hours of ingestion of contaminated food and last for several days.

Several shigellae are responsible for **shigellosis**, a bacillary dysentery that varies in severity. Ulceration of the large intestine, explosive diarrhea, fever, and dehydration occur in the more severe cases.

Isolation and identification of enteric bacteria from feces, urine, blood, and fecally contaminated materials are of major importance in the diagnosis of enteric infections. Although the Enterobacteriaceae are morphologically alike and in many ways metabolically similar, laboratory procedures for the identification of these bacteria are based on differences in biochemical activities (**Figure 66.1**).

In the past, several **multitest systems** have been developed for differentiation and identification of members of the Enterobacteriaceae. They use microtechniques that incorporate a number of media in a single unit. At least six multitest systems are commercially available. The obvious advantages of these units are the need for minimal storage space, the use of less media, the rapidity with which results may be obtained, and the applicability of the results to a computerized system for identification of organisms. There are also certain disadvantages with these systems, including difficulty in obtaining the proper inoculum size since some media require heavy inoculation while others need to be lightly inoculated, the possibility of media carryover from one compartment to another, and the possibility of using inoculum of improper age. Despite these difficulties, when properly correlated with other properties such as Gram stain and colonial morphology on specialized solid media, these systems are acceptable for the identification of Enterobacteriaceae. The most frequently used systems are discussed.

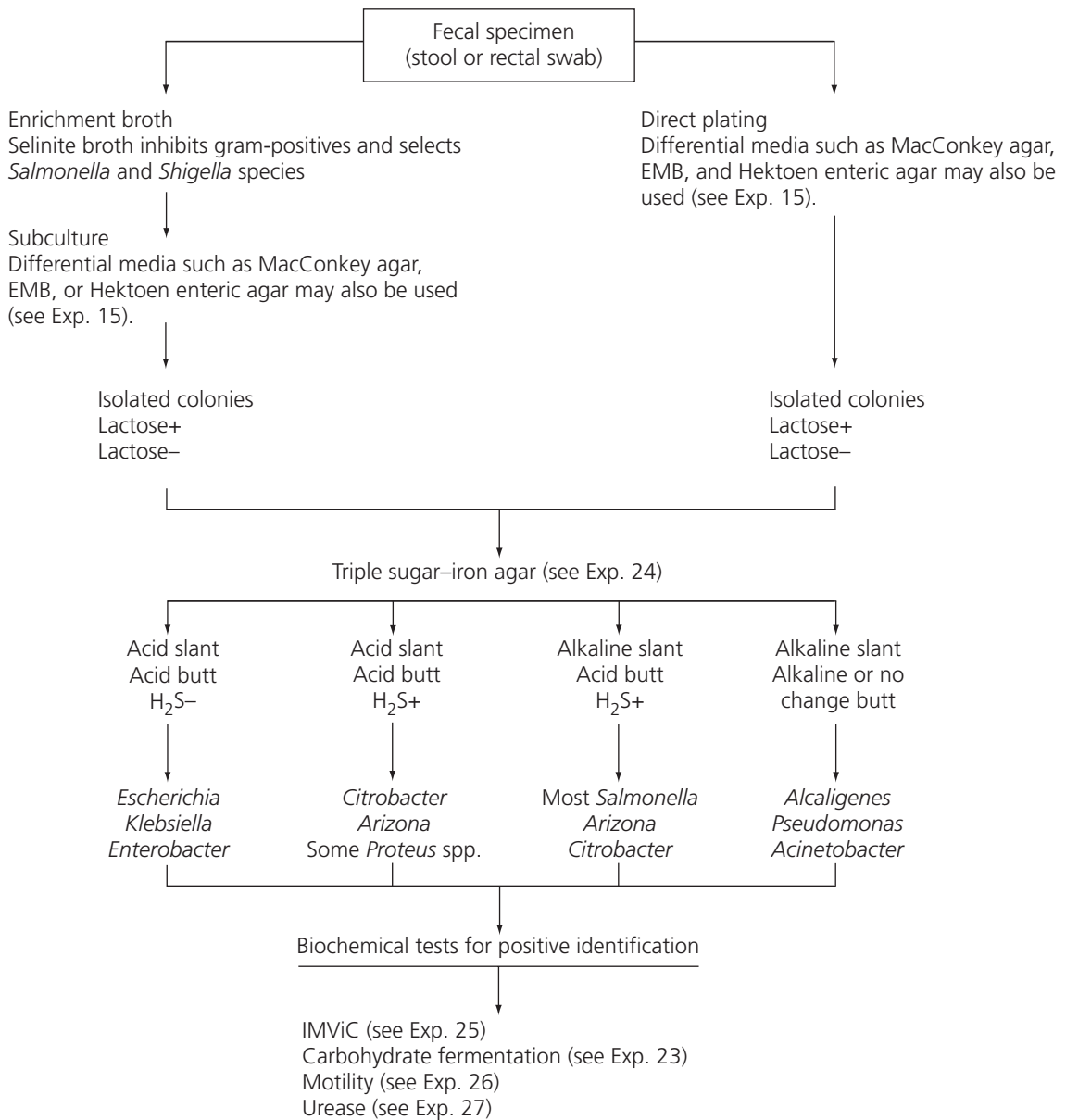


**Figure 66.1** Conventional laboratory procedures for isolation and identification of enteric microorganisms

## Enterotube II Multitest System and ENCISE II

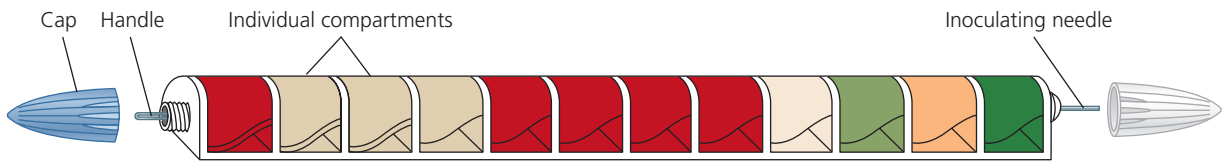The **Enterotube™ II Multitest System** (Roche Diagnostics, Division of Hoffmann-La Roche, Inc.) consists of a single tube containing 12 compartments (**Figure 66.2a**) and a self-enclosed inoculating needle. This needle can touch a single isolated colony and then in one operation be drawn through all 12 compartments, thereby inoculating all of the test media. In this manner, 15 standard

biochemical tests can be performed in one inoculating procedure. Following incubation, the color changes that occur in each of the compartments are interpreted according to the manufacturer's instructions to identify the organisms (**Figure 66.2b**). This method has been further refined to permit identification of the enteric bacteria by means of a computer-assisted system called **ENCISE** (Enterobacteriaceae numerical coding and identification system for Enterotube).



(a) Diagram of Enterotube




(b) Inoculated and uninoculated control

**Figure 66.2** Enterotube II multitest system



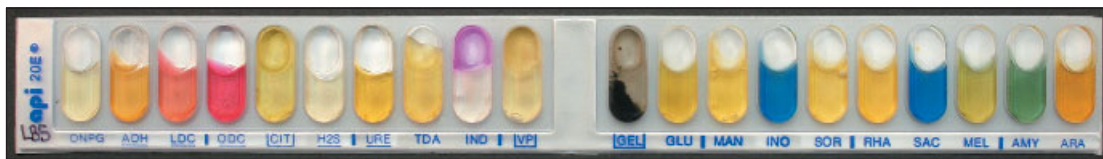(a) Uninoculated control



(b) Inoculated with test organism


**Figure 66.3** The API 20-E system

## API (Analytical Profile Index) System

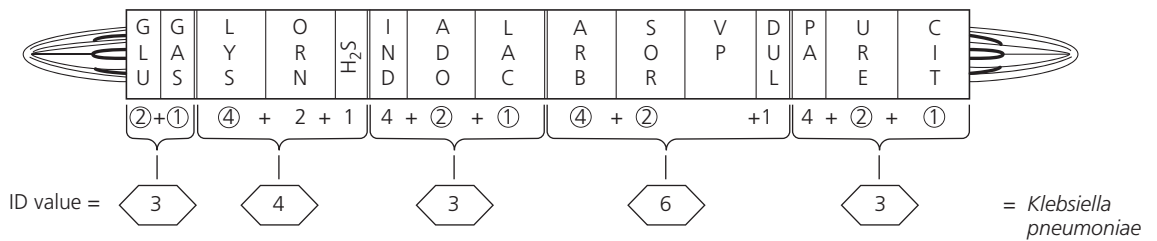The API<sup>®</sup> 20-E employs a plastic strip composed of 20 individual microtubes, each containing a dehydrated medium in the bottom and an upper cupule as shown in **Figure 66.3**. The media become hydrated during inoculation of a suspension of the test organism, and the strip is then incubated in a plastic-covered tray to prevent evaporation. In this manner, 22 biochemical tests are performed. Following incubation, identification of the organism is made by using differential charts supplied by the manufacturer or by means


of a computer-assisted system called **PRS** (Profile Recognition System). PRS includes an API coder, profile register, and selector.

In the following experiment, you will inoculate an Enterotube and an API strip with an unknown enteric organism. Following incubation, you will make your identification by two methods: (1) the traditional method of noting the characteristic color changes and interpreting them according to manufacturers' instructions, and (2) the computer-assisted methods illustrated in **Figure 66.4**.



1. Each positive reaction is indicated by circling the number directly below its compartment.
2. The circled numbers in each bracket are added together, and the sum is placed in the box below.
3. The resultant 5-digit number (ID value) is then located in the computer coding manual to identify the organism.

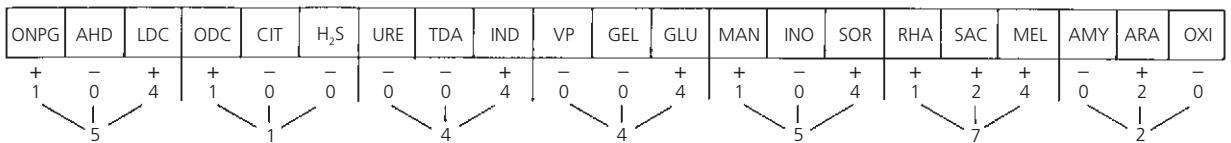(a) The Enterotube® II



7-digit profile number = 5 1 4 4 5 7 2 = *E. coli*

1. The 21 tests are divided into seven groups of three each.
2. A value of 1 is assigned to the first positive test in each group.
3. A value of 2 is assigned to the second positive test in each group.
4. A value of 4 is assigned to the third positive test in each group.
5. A 7-digit number is obtained by totaling the positive values of each of the seven groups of three. This number is located in the analytical profile index to identify the organism.

(b) The API® strip

**Figure 66.4** Computer-assisted techniques for the identification of Enterobacteriaceae

## CLINICAL APPLICATION

### Enterobacteriaceae Infections

The Enterobacteriaceae are a very diverse group of bacteria that commonly inhabit the human colon, but can cause a variety of infections throughout the body. In the hospital environment these often result from colonization of intravascular catheters, leading to bacteremia that can progress rapidly to sepsis and septic shock. Once identification of the infectious agent has been made, treatment with effective antimicrobials may be used.

## AT THE BENCH



## Materials

### Cultures

Number-coded, 24-hour Trypticase soy agar streak plates of *Escherichia coli*, *Salmonella typhimurium*, *Klebsiella pneumoniae*, *Enterobacter aerogenes*, *Shigella dysenteriae*, and *Proteus vulgaris*.

## Media

Per designated student group: one Enterotube II, one API 20-E strip, and one 5-ml tube of 0.85% sterile saline.

## Reagents

Sterile mineral oil, 10% ferric chloride, Kovac's reagent, VP reagent for API system, nitrate reduction reagents, Barritt's reagent (VP test reagent for Enterotube II system), 1.5% hydrogen peroxide, and 1% *p*-aminodimethylaniline oxalate (oxidase reagent).

## Equipment

Bunsen burner, inoculating loop, 5-ml pipette, mechanical pipetting device, sterile Pasteur pipettes, glassware marking pencil, API profile recognition system and differential identification charts, and Enterotube II ENCISE pads and color reaction charts.

## Procedure Lab One

### Enterotube II System

1. Familiarize yourself with the components of the system: screw caps at both ends, medium-containing compartments, self-enclosed inoculating needle, plastic side bar, and blue-taped section.
2. Label the Enterotube II with your name and the number of the unknown culture supplied by the instructor.
3. Remove the screw caps from both ends of the Enterotube II. Using the inoculating needle contained in the Enterotube II, aseptically pick some inoculum from an isolated colony on the provided streak-plate culture.
4. Inoculate the Enterotube II as follows:
  - a. Twist the needle in a rotary motion and withdraw it slowly through all 12 compartments.
  - b. Replace the needle in the tube and with a rotary motion push the needle into the first three compartments (GLU/GAS, LYS, and ORN). The point of the needle should be visible in the H<sub>2</sub>S/IND compartment.
5. Break the needle at the exposed notch by bending, discard the needle remnant, and replace the caps at both ends. The presence of the needle in the three compartments maintains anaerobiosis, which is necessary for glucose fermentation, CO<sub>2</sub> production, and the decarboxylation of lysine and ornithine.
5. Remove the blue tape covering the ADO, LAC, ARB, SOR, VP, DUL/PA, URE, and CIT compartments. Beneath this tape are tiny air vents that provide aerobic conditions in these compartments.
6. Place the clear plastic slide band over the GLU/GAS compartment to contain the wax, which may be spilled by the excessive gas production of some organisms.
7. Incubate the tube on a flat surface for 24 hours at 37°C.

### API 20-E System

1. Familiarize yourself with the components of the system: incubation tray, lid, and the strip with its 20 microtubes.
2. Label the elongated flap on the incubation tray with your name and the number of the unknown culture supplied by the instructor.
3. With a pipette, add approximately 5 ml of tap water to the incubation tray.
4. Using a sterilized loop, touch an isolated colony on the provided streak-plate culture, transfer the inoculum to a 5-ml tube of sterile saline, and mix well to effect a uniform suspension.
5. Remove the API strip from its sterile envelope and place it in the incubation tray.
6. Tilt the incubation tray. Using a sterile Pasteur pipette containing the bacterial saline suspension, fill the tube section of each compartment by placing the tip of the pipette against the side of the cupule. Fill the cupules in the CIT, VP, and GEL microtubes with the bacterial suspension.
7. Using a sterile Pasteur pipette, fill the cupules of the AHD, LDC, ODC, and URE microtubes with sterile mineral oil to provide an anaerobic environment.
8. Cover the inoculated strip with the tray lid and incubate for 18 to 24 hours at 37°C.

## Procedure Lab Two

### Enterotube II System

1. Observe all reactions in the Enterotube II except IND and VP, and interpret your observations using the manufacturer's instructions. Record your observations and results in the Lab Report.
2. Perform the IND and VP tests as follows:
  - a. Place the Enterotube II in a rack with the GLU and VP compartments facing downward.
  - b. With a needle and a syringe, gently pierce the plastic film of the H<sub>2</sub>S/IND compartment and add 2 or 3 drops of Kovac's reagent. Read the results after 1 minute.
  - c. As in Step 2b, add 2 drops of Barritt's reagent to the VP compartment and read the results after 20 minutes.
  - d. Record your IND and VP observations and results in the Lab Report.
3. Based on your results, identify your unknown organism using the manufacturer's color identification charts.
4. Determine and record in the Lab Report the five-digit ID value as described in **Figure 66.4a** on page 456. Identify your unknown organism by referring to the computer coding manual.

### API 20-E System

1. Observe all reactions in the API strip that do not require addition of a test reagent, and interpret your observations using the manufacturer's instructions. Record your observations and results in the Lab Report.
2. Add the required test reagents in the following order: Kovac's reagent to IND, VP reagent to VP (read the result after 15 minutes), ferric chloride to TDA, nitrate reagents to GLU, and oxidase reagent to OXI. Note color changes and interpret your observations according to the manufacturer's instructions. Record your observations and results in the Lab Report.
3. Based on your results, identify your unknown organism using the differential identification chart.
4. Determine and record in the Lab Report the seven-digit profile number as described in **Figure 66.4b** on page 456. Identify your unknown organism by referring to the Profile Recognition System.

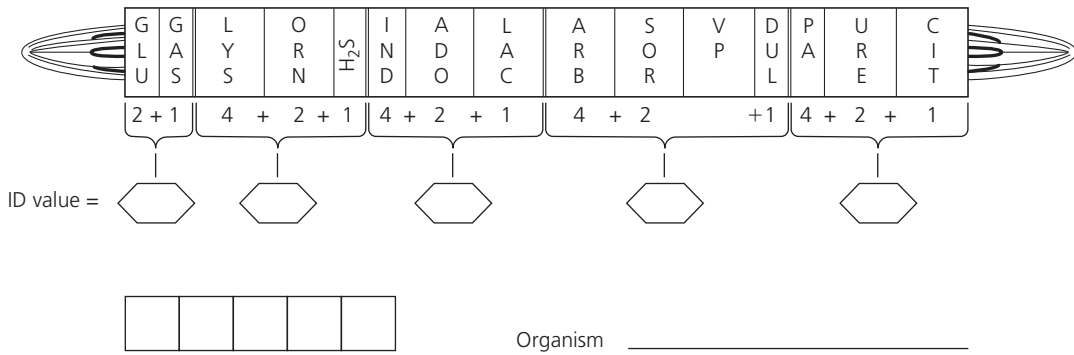
Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_


## Observations and Results

Code	Name	API 20-E		ENTEROTUBE II	
		Appearance (color)	Result (+) or (-)	Appearance (color)	Result (+) or (-)
ONPG	$\beta$ -Galactosidase				
AHD	Arginine dihydrolase				
LDC/LYS	Lysine decarboxylase				
ODC/ORN	Ornithine decarboxylase				
CIT	Citrate				
H <sub>2</sub> S	Hydrogen sulfide				
URE	Urease				
TDA	Tryptophan deaminase				
IND	Indole				
VP	Acetonin				
GEL	Gelatin				
GLU	Glucose				
MAN	Mannitol				
INO	Inositol				
SOR	Sorbitol				
RHA	Rhamnose				
SAC	Sucrose				
MEL	Melibiose				
AMY	Amygdalin				
ARA/ARB	Arabinose				
OXI	Oxidase				
ADO	Adonitol fermentation				
GAS	Gas production				
PHE/PA	Phenylalanine				
LAC	Lactose				
DUL	Dulcitol				
<b>Organism</b>					





**Determination of Enterotube II five-digit identification number**

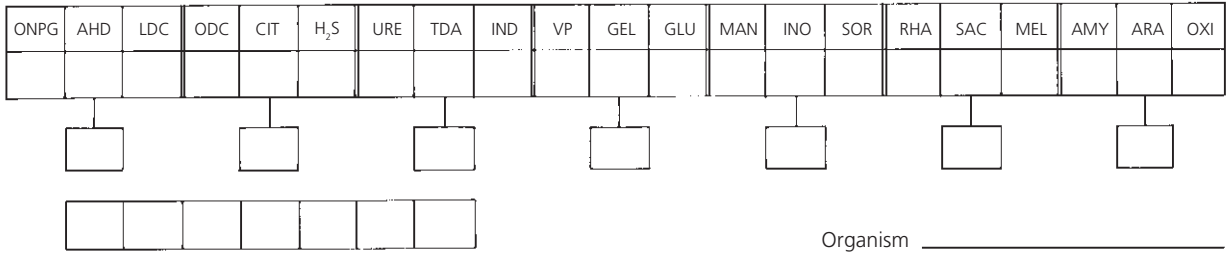

**Determination of API 20-E seven-digit profile number**

**Review Questions**


1. What are the advantages of multitest systems?

Disadvantages?

2. What Enterobacteriaceae are of medical significance?

List and describe the infections caused by these organisms.

3. Would similar results be obtained by use of the computer-assisted method and the traditional color-change method?

4.  What is the clinical justification for the use of a rapid test procedure such as the Enterotube II System for the identification of enteric microorganisms?

*This page intentionally left blank*

## LEARNING OBJECTIVE

Once you have completed this experiment, you should understand

1. The laboratory procedures required for the isolation, cultivation, and presumptive identification of the genus *Campylobacter*.

## Principle

Clinicians are aware of the medical significance of *Campylobacter* strains as the etiological agents of enteric infections. The incidence of enteritis caused by *Campylobacter jejuni* equals or exceeds that of salmonellosis or shigellosis. The clinical syndrome, although varying in severity, is generally characterized by acute gastroenteritis accompanied by the rapid onset of fever, headache, muscular pain, malaise, nausea, and vomiting. Twenty-four hours following this acute phase, diarrhea develops that may be mucoid, bloody, bile-stained, and watery. The precise epidemiology of the infection is not clear; however, contact with animals, waterborne organisms, and fecal-oral transmission remain suspect.

The organisms (*campylo*, curved; *bacter*, rod) were formerly called vibrios because of their curved and spiral morphology. In the early 1980s they were reclassified in the genus *Campylobacter*. They are gram-negative and curved or spiral, with a single flagellum located at one or both poles of the cell. In pure culture, two types of colonies have been recognized and designated as Types I and II. The more commonly observed Type I colonies are large, flat, and spread with uneven margins. They are nonhemolytic, watery, and grayish. Type II colonies are also nonhemolytic, but they are smaller (1 to 2 mm), with unbroken edges. They are convex and glistening.

Initially, the isolation of *Campylobacter* organisms from fecal specimens was difficult because of their microaerophilic nature and their 42°C optimal growth temperature. Furthermore, in

the absence of selective media, their growth was masked by the overgrowth of other enteric organisms, and they were often overlooked on primary isolation. This situation has been rectified with the development of selective media that are designed specifically for isolating *Campylobacter* species and that inhibit the growth of other enteric organisms. These media are nutritionally enriched and supplemented with 5% to 10% sheep or horse blood. In addition they contain three to five antimicrobial agents, depending on the medium. For example, cephalosporins, one of the antimicrobial agents present in the Campy-BAP medium, is selective for *C. jejuni* and inhibits the species *C. intestinalis*, which is rarely responsible for enteric infections.

The most essential requirement for cultivating campylobacteria is a microaerophilic incubation atmosphere. High concentrations of oxygen are toxic to these organisms, and an atmosphere of 3% to 10% carbon dioxide and 5% to 10% oxygen is optimal for their growth. The incubation temperature for *C. jejuni* is 42°C. At this temperature the organism grows optimally, while growth of *C. intestinalis* is inhibited.

In the experiment to follow, a simulated fecal specimen (a culture containing an attenuated strain of *C. jejuni* and other enteric organisms) is used. You will attempt to isolate the *Campylobacter* organisms by using the following two procedures:

1. A conventional method uses MacConkey agar directly, circumventing enrichment procedures, using a mixed simulated fecal population as the test culture.
2. A special method employs Campy-BAP agar and the CampyPak<sup>®</sup> and GasPak<sup>®</sup> jar, which are illustrated in **Figure 67.1**.

Presumptive identification is made on the basis of colonial morphology and the microscopic appearance of the organisms obtained from a typical isolated colony. You may perform the catalase and oxidase tests as described in Experiments 30 and 31 for further presumptive identification. In the case of *C. jejuni*, both tests should be positive.



**Figure 67.1** CampyPak and GasPak jar

## CLINICAL APPLICATION

### Traveler's Diarrhea

*Campylobacter* is the most common cause of bacterial diarrheal diseases worldwide. It is also referred to as "Traveler's Diarrhea." Poultry and poultry products have been associated with *Campylobacter* infections. Other foods have also been implicated in its transmission. *Campylobacter jejuni* and *Campylobacter coli* are the two most clinically significant and may be isolated from the intestinal tract of poultry. They are slow-growing organisms and are identified by biochemical, immunological, and molecular techniques.

## AT THE BENCH



## Materials

### Cultures

Mixed saline suspensions of *Campylobacter jejuni* cultured on a sheep blood-enriched

medium, *Salmonella typhimurium*, and *Escherichia coli*.

### Media

Per designated student group: one Campy-BAP agar plate and one MacConkey agar plate.

### Reagents

Crystal violet, Gram's iodine, 95% ethyl alcohol, and 0.8% carbol fuchsin.

### Equipment

Bunsen burner, inoculating loop, glassware marking pencil, CampyPak and GasPak jars, and 10-ml pipettes.

## Procedure Lab One

1. Aseptically perform a four-way streak inoculation as described in Experiment 3 for the isolation of discrete colonies on both appropriately labeled agar plates.
2. Place the inoculated Campy-BAP agar plate in the GasPak jar in an inverted position. Following the manufacturer's instructions, open the CampyPak envelope and place it in the jar. With a pipette, add 10 ml of water to each envelope and immediately seal the jar to establish a microaerophilic environment.
3. Incubate the jar for 48 hours at 42°C.
4. Incubate the MacConkey agar plate culture in an inverted position for 48 hours at 37°C.

## Procedure Lab Two



1. Observe both plate cultures for the presence of discrete colonies. Record your observations in the chart provided in the Lab Report.
2. Prepare a Gram stain, using 0.8% carbol fuchsin as the counterstain, of a representative colony agar plate culture. Observe microscopically and record in the Lab Report the microscopic morphology and Gram reaction of each preparation.
3. Based on your observations, identify your isolates and record in the Lab Report.
4. Optional: Perform the catalase and oxidase tests on the representative isolates.

Name: \_\_\_\_\_




Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

- In the chart below, diagram the appearance of representative colonies on both plates and describe their colonial characteristics. Also, note and record the color of the medium surrounding the representative colonies on the MacConkey plate. (Refer to Experiment 15 for an explanation of the selective and differential nature of MacConkey agar.)

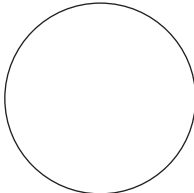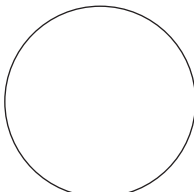Plate Culture	Diagram of Colonies	Colonial Characteristics	Color of Medium
Campy-BAP agar			
MacConkey agar			

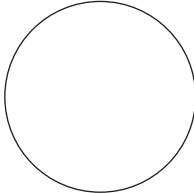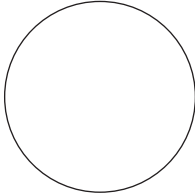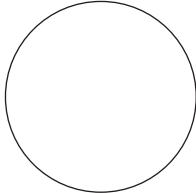- Record your observations of the Gram reactions in the chart below.

Gram Stain Preparation	Campy-BAP Plate Isolate	MACCONKEY AGAR PLATE	
		Isolate 1	Isolate 2
Draw a representative field.			
Microscopic morphology			
Gram reaction			

3. Based on your observations, identify your isolates:

Campy-BAP agar culture isolate: \_\_\_\_\_

MacConkey agar culture Isolate 1: \_\_\_\_\_

MacConkey agar culture Isolate 2: \_\_\_\_\_

## Review Questions

1. How would you describe the clinical syndrome induced by *C. jejuni*?

2. What are the purposes of the antimicrobial agents present in the selective media used for the isolation of *Campylobacter*?

3. How may *C. jejuni* be separated from *C. intestinalis*?

4.  Why might members of *Campylobacter* not be isolated from a stool specimen in a diagnostic laboratory?

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be familiar with

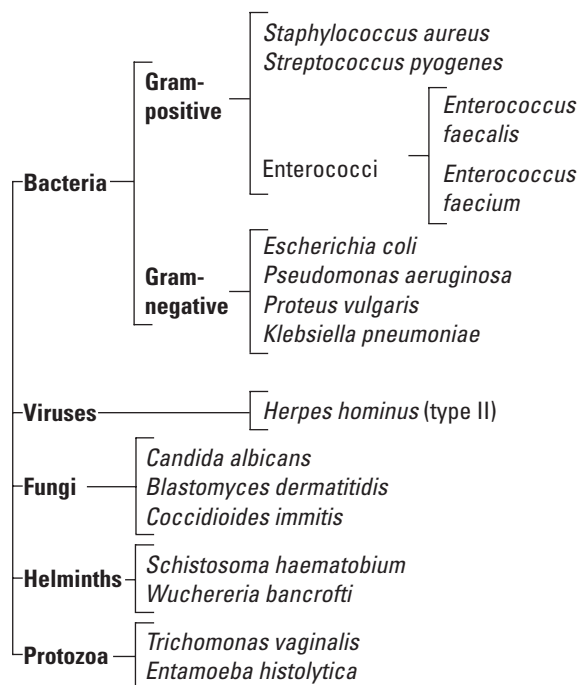
1. The organisms responsible for infections of the genitourinary tract.
2. Laboratory methods for detection of bacteriuria and identification of microorganisms associated with the urinary tract.

## Principle

The anatomical structure of the mammalian urinary system is such that the external genitalia and the lower aspects of the urethra are normally contaminated with a diverse population of microorganisms. The tissues and organs that compose the remainder of the urinary system, the bladder, ureters, and kidneys, are sterile, and therefore urine that passes through these structures is also sterile. When pathogens gain access to this system, they can establish infection. Some etiological agents of urinary tract diseases are illustrated on this page.

Urinary tract infections may be limited to a single tissue or organ, or they may spread upward and involve the entire system. Infections such as **cystitis** involve the bladder but may spread through the ureters to the kidneys. Infections limited to the ureters and kidneys are called **pyelitis**. **Glomerulonephritis** is an inflammation that results in the destruction of renal corpuscles; **pyelonephritis** results in the destruction of renal tubules. Organisms other than bacteria may also act as etiological agents of urogenital infections. *Trichomonas vaginalis*, a pathogenic flagellated protozoan, is commonly found in the vagina, and under appropriate conditions, it is responsible for a severe inflammatory **vaginitis**. *Candida albicans*, a pathogenic yeast, is normally found in low numbers in the intestines. Under suitable conditions, such as the use of antibacterial antibiotics, which disrupt the normal intestinal flora and allow *Candida* to proliferate, it can enter the urogenital


systems, where it gives rise to vaginal infections. *Schistosoma haematobium* is a pathogenic fluke, a helminth, responsible for severe bladder infections.



The initial step in diagnosis of a possible urinary tract infection is laboratory examination of a urine specimen. The sample must be collected midstream in a sterile container following adequate cleansing of the external genitalia. It is imperative to culture the freshly voided, unrefrigerated urine sample immediately to avoid growth of normal indigenous organisms, which may overtake the growth of the more slowly growing pathogens. In this event the infectious organism might be overlooked, resulting in an erroneous diagnosis.

Clinical evaluation of the specimen requires a quantitative determination of the microorganisms per ml of urine. Urine in which the bacterial count per ml exceeds 100,000 ( $10^5$ ) denotes significant **bacteriuria** and is indicative of a urinary tract infection. Urine in which counts range from 0 to 1000 per ml are generally normal.



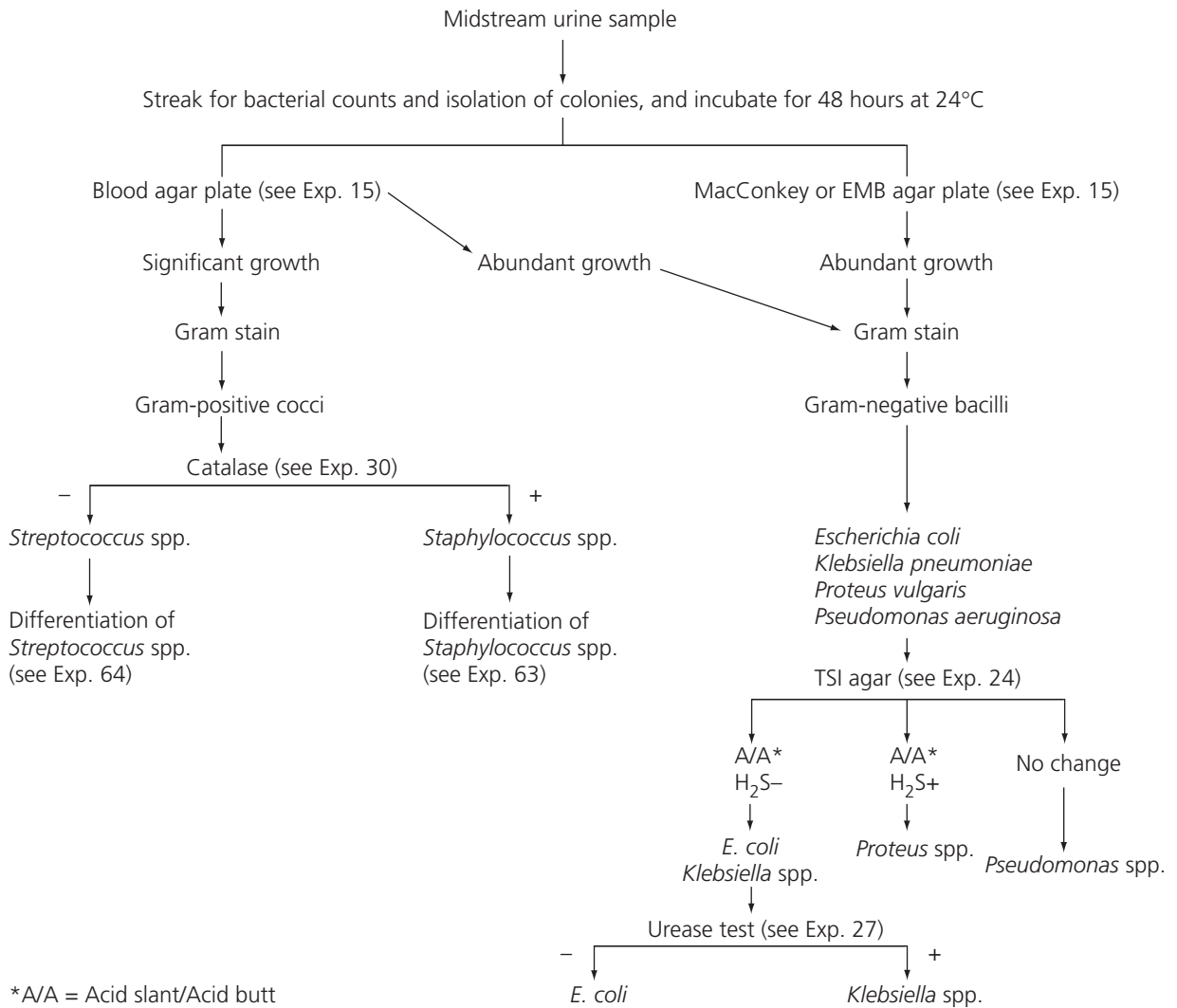

**Figure 68.1** Laboratory procedures for the isolation and identification of urinary tract pathogens

In the conventional method, a urine sample is streaked over the surface of an agar medium with a special loop calibrated to deliver a known volume. Following incubation, the number of isolated colonies present on the plate is determined and multiplied by a factor that converts the volume of urine to 1 ml. The final calculation is then equal to the number of organisms per ml of sample.

**Example:** Twenty-five colonies were present on a plate inoculated with a loop calibrated to deliver 0.01 ml of a urine specimen.

$$\begin{array}{rcccl} \text{number of} & \times & \text{factor that} & = & \text{organisms} \\ \text{colonies} & & \text{converts} & & \text{per ml} \\ & & 0.01 \text{ ml to } 1 \text{ ml} & & \\ 25 & \times & 100 & = & 2500 \text{ organisms} \\ & & & & \text{per ml} \end{array}$$

If the specimen is turbid, dilution is necessary prior to culturing. In this case, conventional

10-fold dilutions are prepared in physiological saline to effect a final dilution of 1:1000 (see Experiment 20). Each of the dilutions ( $10^{-1}$ ,  $10^{-2}$ , and  $10^{-3}$ ) is then streaked on the surface of a suitable agar plate medium for isolation of colonies. Following incubation, the number of microorganisms per ml of sample is determined by the following formula:

$$\begin{array}{l} \text{organisms per ml} = \text{number of colonies} \\ \quad \times \text{factor that converts the} \\ \quad \quad \quad \text{volume of urine to } 1 \text{ ml} \\ \quad \times \text{dilution factor} \end{array}$$

**Example:** Twenty-five colonies were counted on a  $10^{-2}$  dilution plate inoculated with a loop calibrated to deliver 0.01 ml of urine.

Calculation:

$$25 \times 100 \times 100 = 250,000 \text{ organisms per ml}$$

On determination of bacteriuria, identification of the infectious organism can be accomplished by the laboratory procedures outlined in **Figure 68.1**.

A newer, less conventional, and less time-consuming method uses a diagnostic urine-culture tube, Bacturcult<sup>®</sup>, devised by Wampole Laboratories (**Figure 68.2**). Bacturcult is a sterile, disposable plastic tube coated on the interior with a special medium that allows detection of the bacteriuria and a presumptive class identification of urinary bacteria.

Following incubation of the Bacturcult urine culture, bacteriuria can be detected with a bacterial count. This is performed by placing the counting strip around the Bacturcult tube over an area of even colony distribution and counting the number of colonies within the circle. The average number of colonies counted is interpreted in Table 68.1.


For the presumptive identification of bacteria, the medium contains two substrates, lactose and urea, and the pH indicator phenol red. Depending on the organism's enzymatic action on these substrates, differentiation of urinary bacteria into three groups following incubation is possible based on observable color changes that occur in the culture:

**Group I:** *E. coli* and *Enterococcus*—yellow.

**Group II:** *Klebsiella*, *Staphylococcus*, and *Streptococcus*—rose to orange.

**Group III:** *Proteus* and *Pseudomonas*—purplish-red.

Mixed cultures do not always produce clear-cut color changes, however. Therefore, if additional testing is required, the discrete colonies that develop on the medium can be used as the source for subculturing into other media.

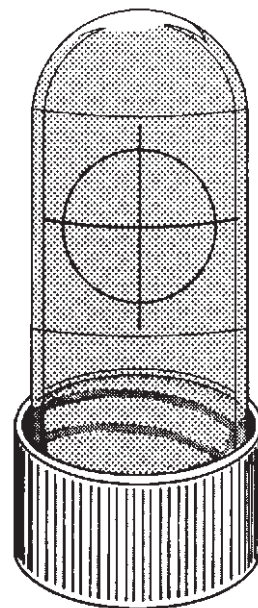

**Figure 68.2** Bacturcult culture tube

In this experiment, seeded saline cultures will be used to simulate urine specimens. This is done to minimize the risk of using a potentially infectious body fluid, urine, as the test sample. The conventional procedure performed with the calibrated loop will be used to determine the number of cells in the specimens. The Bacturcult tube will be used for enumeration and presumptive group identification. If your instructor desires to emulate more closely a clinical evaluation of urine, then a mixed seeded culture must be used. Representative colonies isolated from the blood agar streak-plate culture for detection of bacteriuria can then be identified following the schema in Figure 68.1.

**TABLE 68.1** Bacturcult: Interpretation of Colony Counts

AVERAGE NUMBER OF COLONIES WITHIN CIRCLE	APPROXIMATE NUMBER OF BACTERIA PER ML	DIAGNOSTIC SIGNIFICANCE
<25	<25,000	Negative bacteriuria
25 to 50	25,000 to 100,000	Suspicious*
<50	<100,000	Positive bacteriuria

Source: Wampole Laboratories Division, Carter-Wallace, Inc., Cranbury, NJ 08512. Reprinted with permission.

\*Additional testing recommended.

## CLINICAL APPLICATION

### The Oldest Clinical Test

Urinary tract infections are among the most frequently occurring problems in clinical medicine. Urine is composed of 95% water with the remainder consisting mainly of urea, uric acid, ammonia, hormones, sloughed squamous cells, proteins, salts, and minerals. Urinalysis is performed for the diagnosis of metabolic or systemic diseases that affect kidney function, for disorders of the kidney and urinary tract, screening for drug abuse, and monitoring patients with diabetes. Urinalysis is considered to be the oldest clinical test, with physical examination of urine for diagnosis having been performed as long as 6000 years ago. Hippocrates, in the 4th century BCE, first realized that urine was a filtrate from the kidneys.

## AT THE BENCH



## Materials

### Cultures

Six saline cultures, each seeded with one of the following 24-hour cultures: *Enterococcus faecalis*, *Staphylococcus aureus*, *Proteus vulgaris*, *Escherichia coli*, *Pseudomonas aeruginosa*, and *Klebsiella pneumoniae*. Optional: Saline culture seeded with a gram-positive and a gram-negative organism.

### Media

Per designated student group: three blood agar plates, three sterile 9-ml tubes of saline, and six Bacturcult culture tubes.

### Equipment

Bunsen burner, calibrated 0.01-ml platinum loop, glassware marking pencil, sterile 1-ml pipettes, and mechanical pipetting device.

## Procedure Lab One

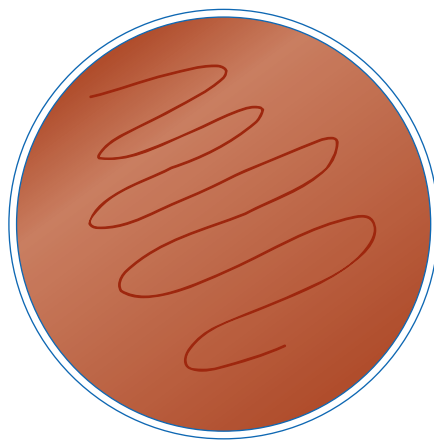
### Bacturcult

1. Label each Bacturcult tube with the name of the bacterial organism present in the urine sample.

2. Fill each tube almost to the top with urine.
3. Immediately pour the urine out of each tube, allowing all the fluid to drain for several seconds. Replace the screw cap securely.
4. Immediately prior to incubation, loosen the cap on each tube by turning the screw cap counterclockwise for one-half turn.
5. Incubate the tubes with the caps down for 24 hours at 37°C.

## Calibrated Loop for Bacterial Counts

1. Label the three 9-ml sterile saline tubes and the three blood agar plates  $10^{-1}$ ,  $10^{-2}$ , and  $10^{-3}$ , respectively.
2. Using the three 9-ml saline blanks, aseptically prepare a 10-fold dilution of the urine sample to effect  $10^{-1}$ ,  $10^{-2}$ , and  $10^{-3}$  dilutions.
3. With a calibrated loop, aseptically add 0.01 ml of the  $10^{-1}$  urine dilution to the appropriately labeled blood agar plate and streak for isolation of colonies as illustrated.



4. Repeat Step 3 to inoculate the remaining urine sample dilutions.
5. Incubate all plates in an inverted position for 24 hours at 37°C.

## Procedure Lab Two

1. Determine the number of colonies in each of the Bacturcult urine cultures (refer to Lab Report for further instructions).
2. Record your results in the Lab Report.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

### Bacturcult Procedure

- Determine the number of colonies in each of the Bacturcult urine cultures as follows:
  - Place the counting strip around the tube over an area of even colony distribution and count the number of colonies within the circle.
  - Repeat the count in another area of the tube.
  - Average the two counts.
  - Record in the Lab Report the average number of colonies counted within the circle.
- Based on your colony count, determine and record in the Lab Report the approximate number of bacteria per ml of each sample and its diagnostic significance as negative bacteriuria, suspicious, or positive bacteriuria.
- Observe and record in the Lab Report the color of the medium in each of the urine cultures and the presumptive bacterial group.



Urine Culture	Number of Colonies	Number of Bacteria per ml	Diagnostic Significance	Color of Medium	Presumptive Group
<i>E. faecalis</i>					
<i>S. aureus</i>					
<i>K. pneumoniae</i>					
<i>P. vulgaris</i>					
<i>P. aeruginosa</i>					
<i>E. coli</i>					

### Calibrated Loop Procedure

Determine the number of colonies on each blood agar culture plate and calculate the number of organisms per ml of the urine. Record your results in the Lab Report.

Urine Sample Dilution	Number of Colonies	Organisms per ml of Sample	Bacteriuria (+) or (-)
$10^{-1}$			
$10^{-2}$			
$10^{-3}$			

## Review Questions

1. What types of urinary infections may be caused by different microorganisms?
2. How is a clinical diagnosis of a bacteriuria established?
3. If five colonies were counted on a  $10^{-3}$  dilution plate streaked with 0.01 ml of urine, what was the number of organisms per ml of the original specimen, and is this count indicative of bacteriuria? Explain.
4.  How accurate is a laboratory analysis of a 24-hour, unrefrigerated, non-midstream urine sample? Explain.
5.  A male patient is diagnosed as having a urinary tract infection. A urine culture is ordered by his physician. She requests that a voided specimen be used rather than a catheterized sample. Why does she make this request?

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be familiar with

1. The microorganisms most frequently associated with septicemia.
2. Laboratory methods for the isolation and presumptive identification of the etiological agents of septicemia.

## Principle

Blood is normally a sterile body fluid. This sterility may be breached, however, when microorganisms gain entry into the bloodstream during the course of an infectious process. The transient occurrence of bacteria in the blood is designated as **bacteremia** and implies the presence of nonmultiplying organisms in this body fluid.

Bacteremias may be encountered in the course of some bacterial infections such as pneumonia, meningitis, typhoid fever, and urinary tract infections. A bacteremia of this nature does not present a life-threatening situation because the bacteria are present in low numbers and the activity of the host's innate (nonspecific) immune system is generally capable of preventing further systemic invasion of tissues. A more dangerous and clinically alarming syndrome is **septicemia**, a condition characterized by the rapid multiplication of microorganisms, with the possible elaboration of their toxins into the bloodstream. The clinical picture frequently present in septicemia is that of septic shock, which is recognized by a severe febrile episode with chills, prostration, and a drop in blood pressure.

A large and diverse microbial population has been implicated in septicemia. The major offenders include the following:


1. Gram-negative bacteria, because of their endotoxic properties, are the most frequently

encountered etiological agents of the serious complications of septicemia. Among these agents are *Haemophilus influenzae*, *Neisseria meningitidis*, *Serratia marcescens*, *Escherichia coli*, *Pseudomonas aeruginosa*, and *Salmonella* spp. Less frequently implicated are *Francisella tularensis* and members of the genera *Campylobacter* and *Brucella*.

2. Gram-positive bacteria that generally do not produce the presenting signs of septic shock include primarily members of the genera *Streptococcus* and *Staphylococcus*.
3. *Candida albicans* is the major fungal invader of the bloodstream.

In the clinical setting, to facilitate the rapid initiation of effective chemotherapy, a culture of the suspect blood sample is required for the isolation and identification of the offending organisms. A blood sample is drawn and cultured in an appropriate medium under both aerobic and anaerobic conditions. Over a period of 3 to 7 days, the cultures are observed for turbidity and Gram-stained smears are prepared to ascertain the presence of microorganisms in the blood. Upon detection of microbial growth in the cultures, transfers onto a variety of specialized agar media are made for the identification of the infectious agent. The schema for this protocol is shown in **Figure 69.1**.

Two methods are outlined in this exercise. Either method or both methods may be used for the isolation and presumptive identification of the microorganisms in the experimental culture. Both procedures use a simulated blood specimen: a prepared culture containing blood previously seeded with selected microorganisms. The traditional method is a modification of the schema shown in Figure 69.1. This procedure requires the preparation of Gram-stained smears for the morphological study of the organisms and the inoculation of selected agar media for their isolation and preliminary identification. The alternative method uses the commercially available **BBL Septi-Chek™ System**, a single unit composed of the Septi-Chek culture bottle and the Septi-Chek slide

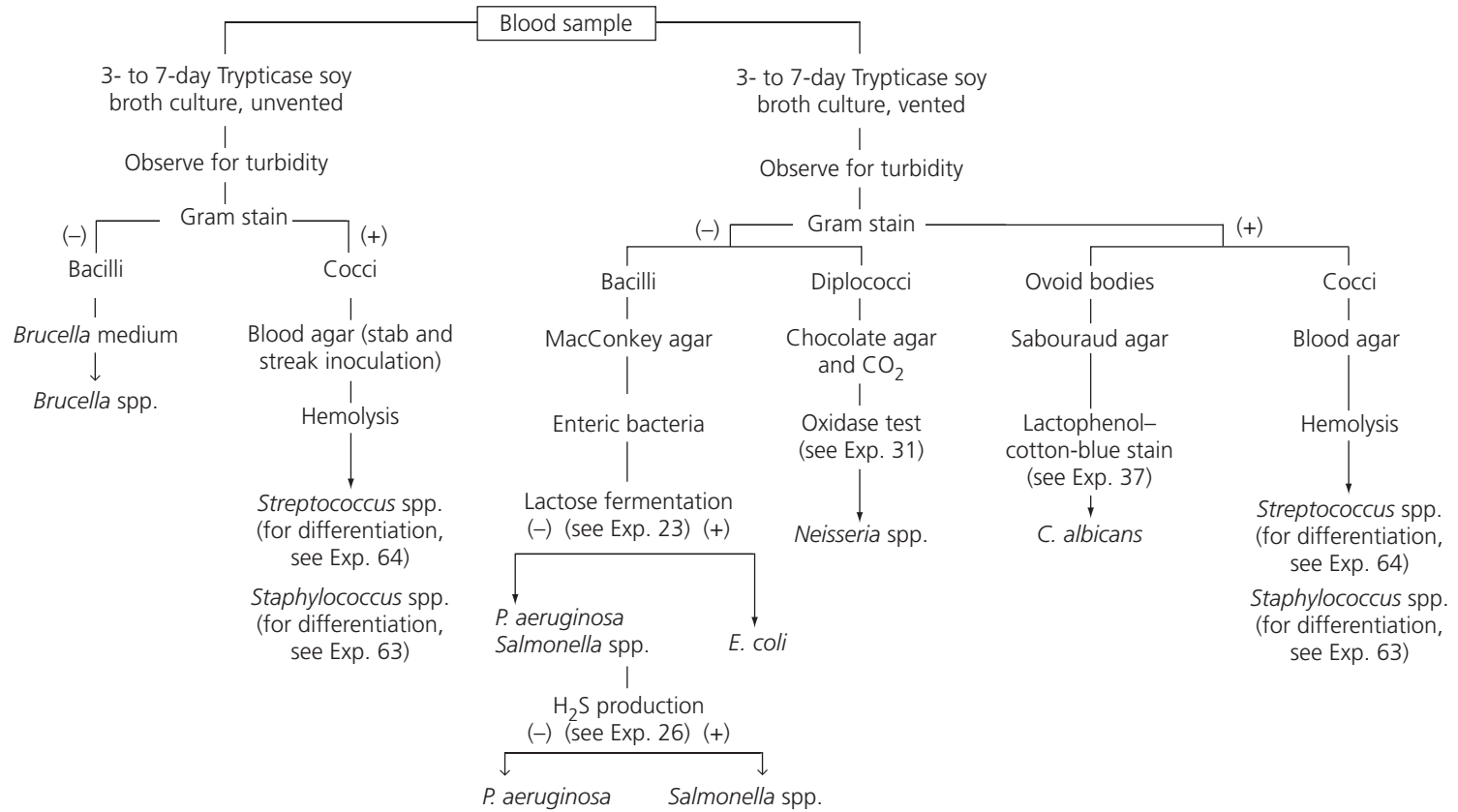

**Figure 69.1** Schema for the isolation and identification of the etiological agents of septicemia

as illustrated in **Figure 69.2**. The culture-bottle component permits the qualitative determination of the presence of microorganisms in the blood sample, and the slide component is designed for the simultaneous subculturing of the organisms onto a plastic slide containing three differential media (chocolate, MacConkey, and malt agar). Differential growth on these media provides preliminary information as to the nature of the infectious agent and isolated colonies for further study.

## CLINICAL APPLICATION

### Drawing Blood for Cultures

Normally, drawing blood for hematological analysis simply requires cleansing of the skin with an alcohol pad, but those draws intended for microbiological testing require a different protocol. The area of the draw is cleaned thoroughly with alcohol followed by a disinfectant such as Chloraprep® One-Step. All palpation after cleansing of the skin must be done with sterile gloves, and the phlebotomist wears a face shield. While it is nearly impossible to eliminate all bacteria from the skin, such techniques attempt to minimize contamination of drawn blood from contact with the skin flora that could produce false positive blood cultures.

## AT THE BENCH



## Materials

### Culture

48- to 72-hour simulated blood culture prepared as follows: 10 ml of citrated blood, obtained from a blood bank, or 10 ml of saline seeded with 2 drops each of *Escherichia coli*, *Neisseria perflava*, and *Saccharomyces cerevisiae*, each adjusted to an absorbance of 0.1 at 600 nm, in 90 ml of Trypticase soy broth containing 0.05% of sodium polyantholesulfonate (SPS) used to prevent clotting of the blood sample.

### Media

Per designated student group: one each of blood agar plate, MacConkey agar plate, chocolate agar plate, Sabouraud agar plate, and Septi-Chek System.



**Figure 69.2** Septi-Chek System

## Reagents

Crystal violet, Gram's iodine, 95% ethyl alcohol, safranin, lactophenol-cotton-blue stain, 70% isopropyl alcohol, and 1% *p*-aminodimethylaniline oxalate.

## Equipment

Sterile 20-gauge, 1½-inch needles; sterile 1-ml and 10-ml syringes; Bunsen burner; staining tray; inoculating loop; glass microscope slide; lens paper; bibulous paper; microscope; glassware marking pencil; and disposable gloves.

## Procedure Lab One

**⚠ Use gloves throughout the procedure.**

1. Swab the rubber stopper of the blood-culture bottle with 70% isopropyl alcohol and allow to air-dry.
2. Using a sterile needle and 1-ml syringe, aseptically remove 0.5 ml of the blood culture by penetrating the rubber stopper.

**⚠ Dispose of the needle and syringe, as a single unit, into the provided puncture-proof receptacle.**

3. To prepare a smear, place a small drop of the culture on a clean glass slide and spread evenly with an inoculating loop.



4. Place 1 drop of culture in one corner of the appropriately labeled blood agar plate and prepare a four-way streak inoculation as described in Experiment 3.
5. Repeat Step 4 to inoculate the MacConkey, chocolate, and Sabouraud agar plates.
6. Incubate the agar plate cultures in an inverted position for 24 to 48 hours as follows: Sabouraud agar culture at 25°C, chocolate agar culture in a 10% CO<sub>2</sub> atmosphere at 37°C, and the remaining cultures at 37°C.
7. Follow the **Septi-Chek System procedure** as follows:
  - a. Remove the protective top of the screw cap of the culture bottle, disinfect the rubber stopper with 70% isopropyl alcohol, and allow to air-dry.
  - b. Using the 10-ml syringe, aseptically transfer 10 ml of the experimental culture to the appropriately labeled Septi-Chek culture bottle.
  - c. Aseptically vent the bottle for aerobic incubation.
  - d. Replace the protective top of the screw cap on the bottle.
  - e. Gently invert the bottle two or three times to disperse the blood evenly throughout the medium.
  - f. Incubate the culture for 4 to 6 hours at 37°C.
  - g. Attach the Septi-Chek slide according to the manufacturer's instructions.
  - h. Tilt the combined system to a horizontal position and hold until the liquid medium enters the slide chamber and floods the agar surfaces. While maintaining this position, rotate the entire system one complete turn to ensure that all agar surfaces have come in contact with the liquid medium. Return the system to an upright position.
  - i. Incubate the system in an upright position at 37°C.
  - j. Check the culture bottle daily for turbidity and the slide for visible colony formation.

## Procedure Lab Two

1. Examine the blood agar plate culture for the presence (+) or absence (–) of hemolytic activity. (Refer to Figure 62.1.) If hemolysis is present, determine the type observed. Record your observations in the Lab Report.
2. For the performance of the oxidase test, add *p*-aminodimethylaniline oxalate to the surface of the growth on the chocolate agar plate. The presence of pink-to-purple colonies is indicative of *Neisseria* spp. (Refer to Figure 62.2.) Record your observations and the oxidase test results in the Lab Report.
3. Examine the MacConkey agar plate culture for determination of lactose fermentation. Lactose fermenters exhibit a pink-to-red halo in the medium, a red coloration on the surface of their growth, or both a halo and red coloration. (Refer to Figure 15.2b.) Record your observations and indicate the presence or absence of lactose fermenters in the Lab Report.
4. Examine the Sabouraud agar plate culture for the presence of growth. Prepare a lactophenol–cotton-blue–stained smear from an isolated colony (see Experiment 37). Examine the smear microscopically for the presence of large ovoid bodies indicative of the yeast cells. Record your morphological observations in the Lab Report.
5. Observe the Septi-Chek slide system for the presence of growth on the three agar surfaces. If growth is present on:
  - a. Medium 1 (MacConkey agar), examine for fermentative patterns as described in Step 3 and record your observations in the Lab Report.
  - b. Medium 2 (chocolate agar), perform the oxidase test as described in Step 2 and record your observations in the Lab Report.
  - c. Medium 3 (malt agar), prepare and examine microscopically a lactophenol–cotton-blue–stained smear as described in Step 4. Record your observations in the Lab Report.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_


## Observations and Results

Culture	Traditional Procedure	Septi-Chek System
<b>Blood agar</b> Hemolysis: (+) or (-) Type of hemolysis	_____ _____	_____ _____
<b>Chocolate agar</b> Color of colonies Oxidase test: (+) or (-)	_____ _____	_____ _____
<b>MacConkey agar</b> Color of colonies Color of medium Lactose fermentation: (+) or (-)	_____ _____ _____	_____ _____ _____
<b>Sabouraud or malt agar</b> Cell morphology	_____	_____
<b>Presumptive identification of organisms present</b>	_____	_____


## Review Questions

1. Differentiate between septicemia and bacteremia, and explain the medical significance of each.

2. Why are blood samples cultured in both vented and unvented systems?

3.  A 15-year-old boy is admitted to the hospital and presents the following symptoms: chills, fever, increased pulse rate, and a drop in blood pressure. The patient indicates that these symptoms have occurred intermittently. The physician suspects a bacteremia and orders a series of three blood cultures over a 24-hour period. Explain the following:
- a. Why did the physician order more than one blood culture?

b. Why does blood culture medium contain an anticoagulant?

4.  Prior to the introduction of antibiotic therapy, what was the prognosis for patients with septicemia? What significant factors played roles in recovery in the absence of antibiotics?

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be able to

1. Identify an unknown bacterial species by the use of dichotomous keys and *Bergey's Manual of Systematic Bacteriology*.

## Principle

At this point in the course, you have developed the manipulative skills and the cognitive microbiological knowledge to identify microorganisms beyond their genus classification to the level of their species identification. Therefore, in this experiment, you will use dichotomous keys, *Bergey's Manual of Systematic Bacteriology*, and information accrued from previously performed laboratory procedures to help identify the species of an unknown culture.

In Experiment 33, “Genus Identification of Unknown Bacterial Cultures,” you were required to use a variety of biochemical tests to successfully accomplish the experimental purpose. Your review of the required procedures and ensuing results should indicate that only a few of these tests were actually necessary, in most instances, for the identification of the unknown culture. Similarly, species identification can be accomplished by using a limited number of carefully selected laboratory procedures. Notice that what appears to be a spurious result in some cases, one that departs from the expected norm for a particular species, may be attributable to strain differences within the given species. These nonconforming results may be verified by the use of *Bergey's Manual* to ascertain the existence of variable biochemical test results for the particular species being studied.

In this experimental procedure, you will receive a mixed culture containing a gram-positive and a gram-negative organism. The protocol will require (1) Gram staining, (2) streak plating for

observation of colonial characteristics, (3) use of selective media for the preparation of pure cultures, (4) the performance of appropriate biochemical tests as indicated in the dichotomous keys outlined in **Figure 70.1** on pages 480–481 and **Figure 70.2** on pages 482–483, and (5) information in *Bergey's Manual*.

## CLINICAL APPLICATION

### New Molecular Techniques for Rapid Species Identification

Once bacteria from blood or other tissues has been cultured, the organisms must be positively identified. While biochemical and serological tests are the norm for such identification, a recently developed technique of mass spectrometry using matrix-assisted laser desorption/ionization (MALDI) offers a quicker (less than one hour after detection in blood) way to identify organisms. This technique releases key molecules from the organisms in question and, using analysis of the size-to-charge ratios of the molecules and specialized computer software, provides accurate identification of infectious organisms and may provide a future alternative or addition to both biochemical and genomic identification schemes.

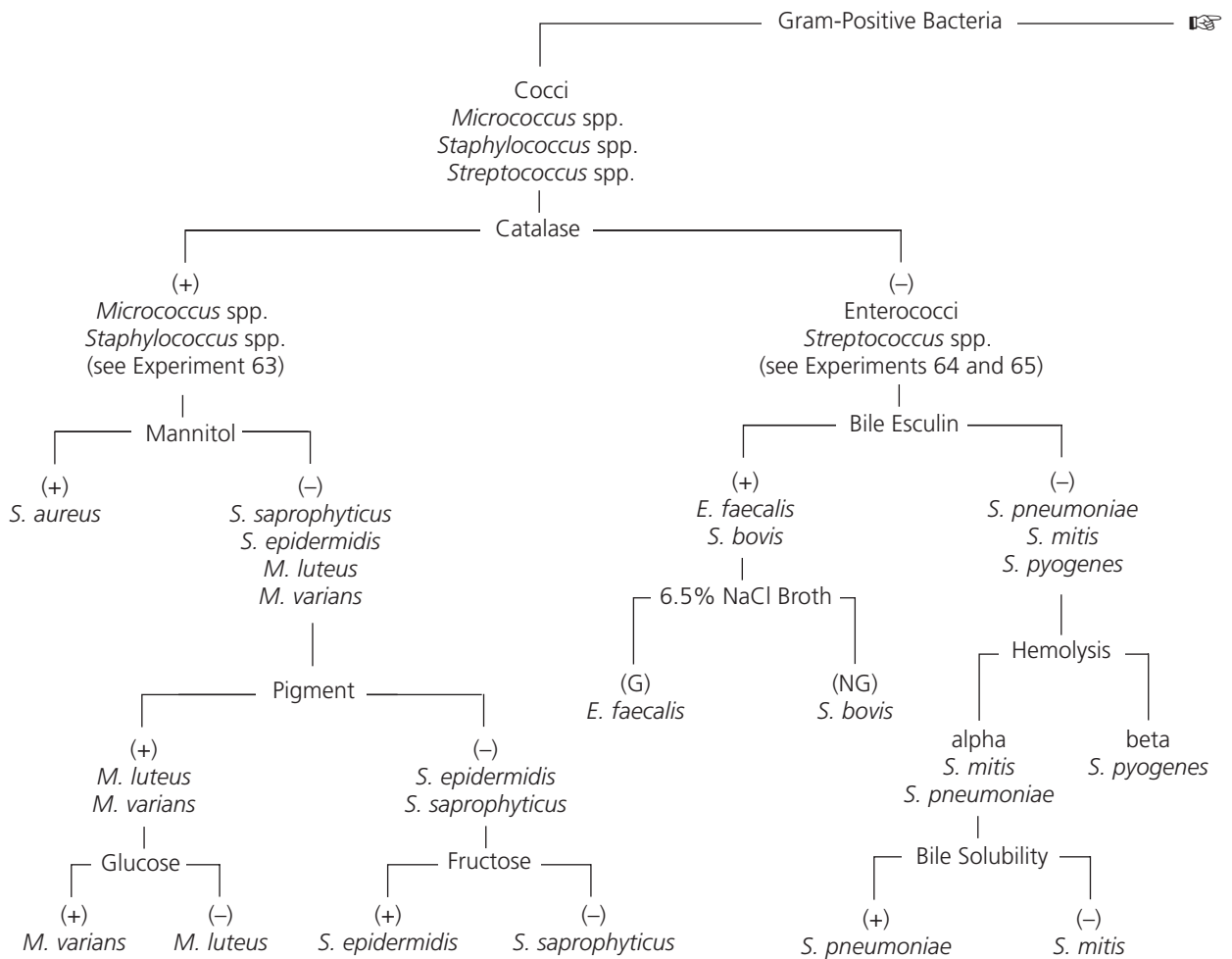
## AT THE BENCH



## Materials

### Cultures

Per student: number-coded, 24- to 48-hour mixed Trypticase soy broth cultures each containing a gram-positive and a gram-negative organism selected from the species listed in Figures 70.1 and 70.2.



NG=No growth; G=Growth; A/G=Acid and gas; A=Acid only

**Figure 70.1** Schema for the identification of gram-positive bacteria

## Media

Per student: one Trypticase soy agar plate, two Trypticase soy agar slants, one Trypticase soy broth, one phenylethyl alcohol agar plate, and one MacConkey agar plate.

Required media for the biochemical tests listed in Figures 70.1 and 70.2 should be available on your request.

## Reagents

Crystal violet, Gram's iodine, 95% ethyl alcohol, safranin, and required reagents for the interpretation of the biochemical reactions listed in Figures 70.1 and 70.2.

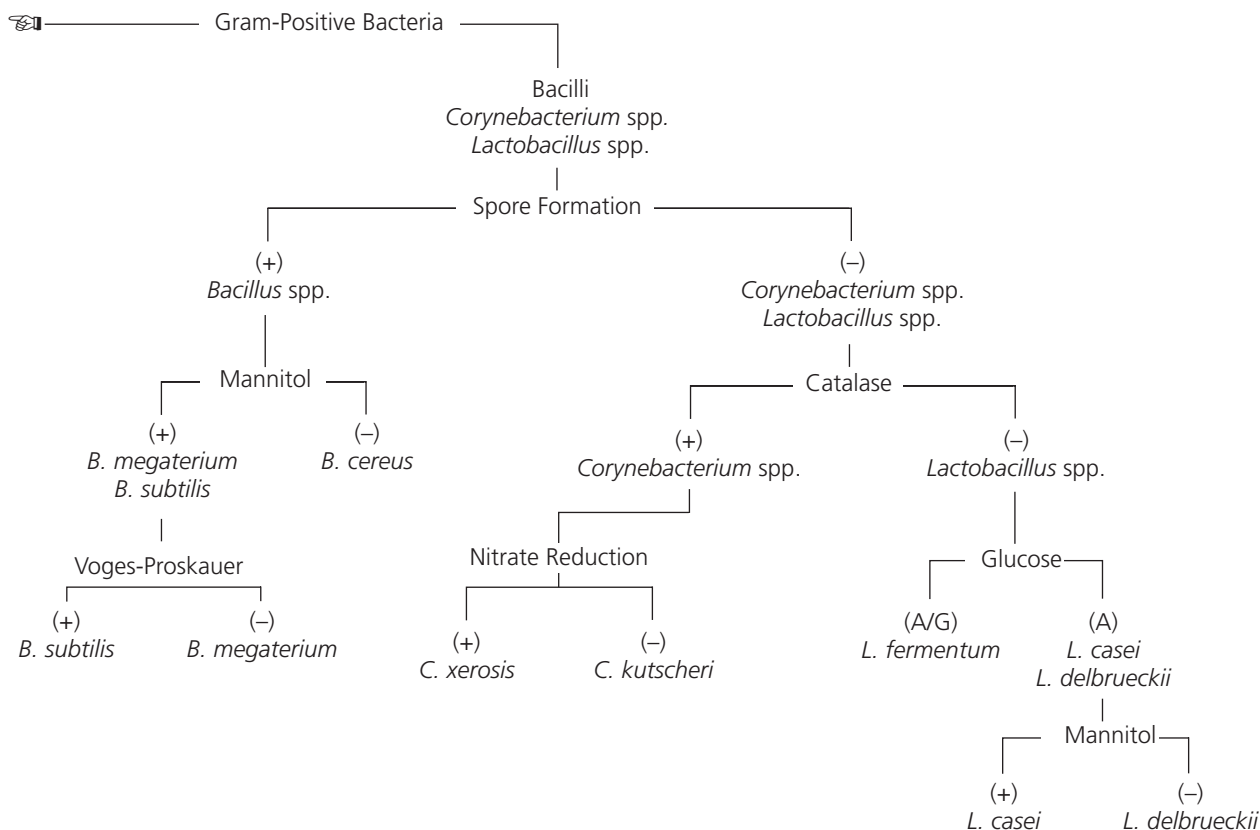
## Equipment

Bunsen burner, inoculating loop and needle, staining tray, immersion oil, lens paper, bibulous paper, microscope, and glassware marking pencil.

## Procedure Lab One

### Separation of the Bacteria in Mixed Unknown Culture

1. Prepare a Trypticase soy agar broth subculture of the unknown and refrigerate following incubation. You will use this culture if contamination of the test culture is suspected during the identification procedure.



NG= No growth; G= Growth; A/G= Acid and gas; A= Acid only


**Figure 70.1 (continued) Schema for the identification of gram-positive bacteria**

2. Prepare a Gram-stained smear of the original unknown culture. Examine the smear and record your observations in the Lab Report.
3. Prepare four-way streak inoculations (see Experiment 3) on the following media for the separation of the microorganisms in the mixed cultures:
  - a. Trypticase soy agar for observation of colonial characteristics.
  - b. Phenylethyl alcohol agar for isolation of gram-positive bacteria.
  - c. MacConkey agar for isolation of gram-negative bacteria.
4. Incubate all the plates in an inverted position and then subculture for 24 to 48 hours at 37°C.

## Procedure Lab Two

### Preparation of Pure Cultures

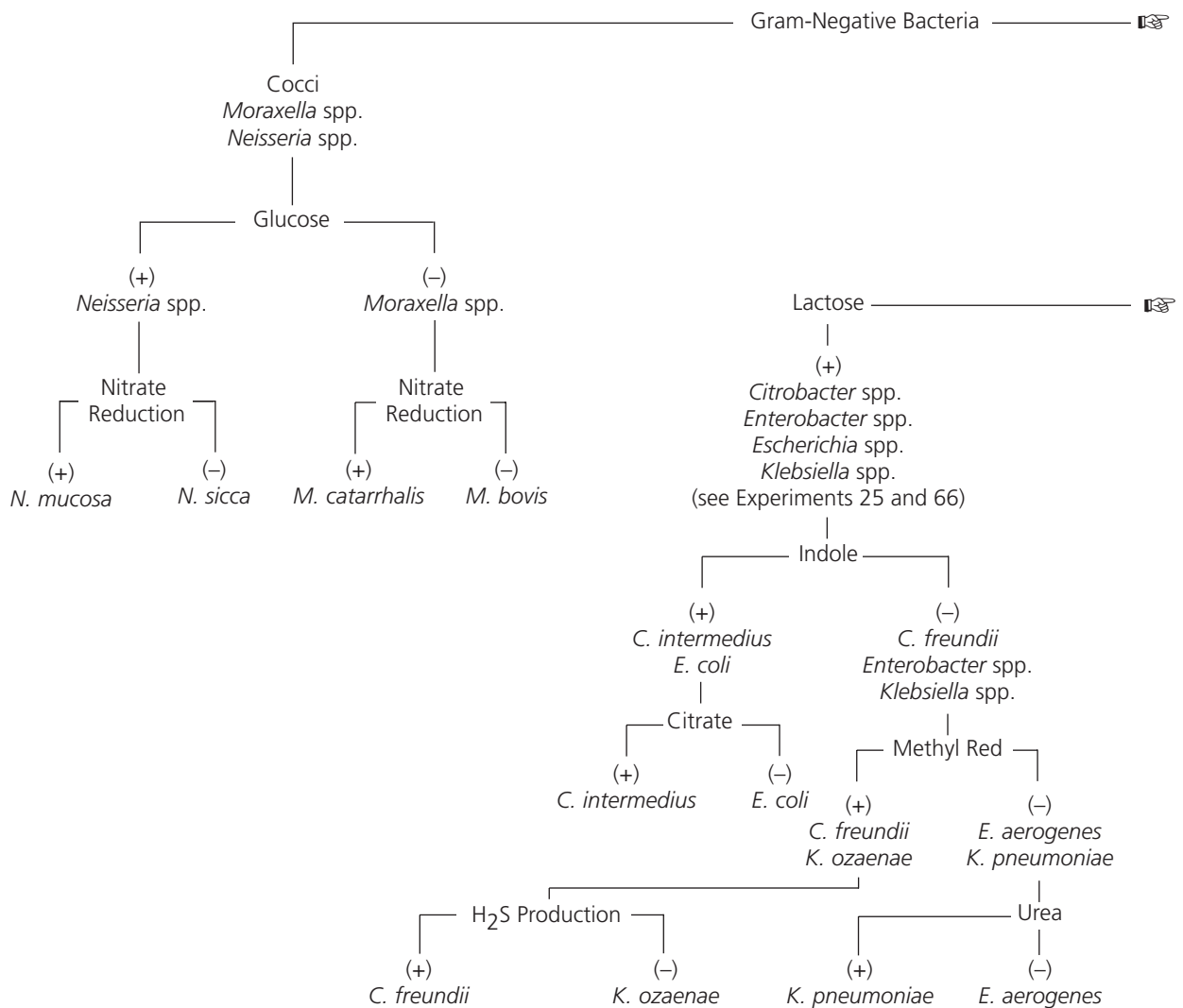1. Isolate a discrete colony on both the phenylethyl alcohol agar plate and the MacConkey agar plate and aseptically transfer each onto a Trypticase soy agar slant (see Experiment 3).
2. Incubate the Trypticase soy agar slants for 24 to 48 hours at 37°C.



NG=No growth; G=Growth; A/G=Acid and gas; A=Acid only

**Figure 70.2** Schema for the identification of gram-negative bacteria


## Procedure Lab Three
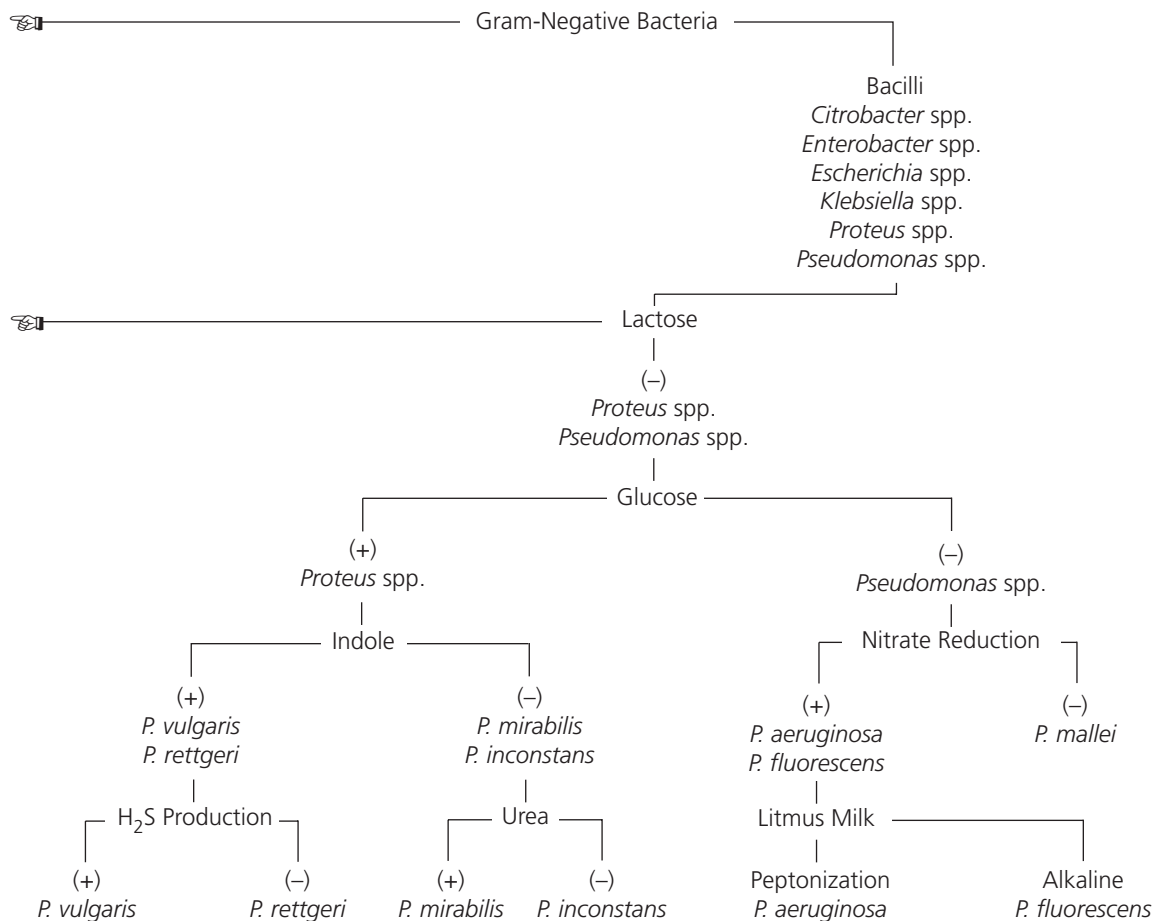
### Preparation of Pure Cultures

1. Examine the Trypticase soy agar plate for the appearance of discrete colonies. Select two colonies that differ in appearance and record their colonial morphologies in the Lab Report.
2. Examine the phenylethyl alcohol and MacConkey agar plates. Record your observations in the Lab Report.

### Identification of Unknown Bacterial Species

1. Prepare a Gram-stained smear from each of the Trypticase soy agar slant cultures to verify their purity by means of the Gram reaction and cellular morphology. Examine the smears and record your observations in the Lab Report.
2. If each Gram-stained preparation is not solely gram-positive or gram-negative, repeat the steps in Labs One and Two, using the



NG= No growth; G= Growth; A/G= Acid and gas; A= Acid only

**Figure 70.2 (continued) Schema for the identification of gram-negative bacteria**

refrigerated Trypticase soy agar subculture as the test culture.

3. If the isolates are deemed to be pure on the basis of their cultural and cellular morphologies, continue with the identification procedure. During this period and in subsequent sessions, use the dichotomous keys in Figures 70.1 and 70.2 to select and perform the necessary biochemical tests on each of your isolates for identification of their species. Incubate all cultures for 24 to 48 hours at 37°C.

## Procedure Lab Four

### Identification of Unknown Bacterial Species

1. Examine all the biochemical test cultures and record your observation and results in the Lab Report.



*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

### Separation of the Bacteria in Mixed Unknown Culture

Record your observations of the Gram-stained smear of your mixed unknown culture in the chart below.

Organism	Cellular Morphology	Gram Reaction
1		
2		

### Preparation of Pure Cultures

- Select from the Trypticase soy agar plates two colonies that differ in appearance, and record their colonial morphologies.

Isolate 1:

Isolate 2:

- Record your observations of the phenylethyl alcohol and MacConkey agar plates in the chart below.

Medium	Growth (+) or (-)	Colonial Morphology	Coloration of Medium
Phenylethyl alcohol agar			
MacConkey agar			

### Identification of Unknown Bacterial Species

- Record your observations of the Gram-stained smears of the Trypticase soy agar cultures obtained from the phenylethyl alcohol and MacConkey agar plates in the chart below.

Agar Slant From	Cellular Morphology	Gram Reaction
Phenylethyl alcohol agar plate culture		
MacConkey agar plate culture		

2. Record your observations and results of all the biochemical tests in the charts below.

GRAM-POSITIVE ISOLATE		
Biochemical Test	Observation	Result

Unknown gram-positive organism: \_\_\_\_\_

GRAM-NEGATIVE ISOLATE		
Biochemical Test	Observation	Result

Unknown gram-negative organism: \_\_\_\_\_

# Immunology

## LEARNING OBJECTIVES

Once you have completed the experiments in this section, you should be familiar with

1. The basic principles of nonspecific (innate) and specific (acquired) immunity.
2. Serological procedures that demonstrate immunological reactions of agglutination and precipitin formation.
3. Rapid immunodiagnostic screening procedures.

## Introduction

**Immunity, or resistance**, is a state in which a person, either naturally or by some acquired mechanism, is protected from contracting certain diseases or infections. The ability to resist disease may be innate (nonspecific), or it may be adaptive (also called acquired, or specific), when the disease state is stimulated in the host.

Innate immunity is native or natural. It is inborn and provides the basic mechanisms that defend the host against intrusion of foreign substances or agents of disease. This defense is not restricted to a single or specific foreign agent, but it provides the body with the ability to resist many pathological conditions. The mechanisms responsible for this native immunity include the **mechanical barriers**, such as the skin and mucous membranes; **biochemical factors**, such as antimicrobial substances present in the body fluids; and the more sophisticated process of **phagocytosis** and action of the **reticuloendothelial** system.

Adaptive immunity, **cell-mediated** and **humoral**, is acquired by the host in response to the presence of a single or particular foreign substance, usually protein, called an **antigen**

(immunogen). In humoral immunity, antigens that penetrate the mechanical barriers of the host, namely the skin and mucous membranes, stimulate formation of **antibodies**. The function of the antibodies is to bind to the specific antigens that are responsible for their production and to inactivate or destroy them. Antibodies are a group of homologous proteins called **immunoglobulins**, which are found in serum and represent five distinct classes: immunoglobulin G (IgG), immunoglobulin A (IgA), immunoglobulin M (IgM), immunoglobulin D (IgD), and immunoglobulin E (IgE).

The primary immunological complexes (antigen + antibody) are as follows:

1. **Agglutination:** This type of reaction uses specific antibodies, **agglutinins**, that are formed in response to the introduction of particulate antigens into host tissues. When these particulate antigens combine with a homologous antiserum, a three-dimensional mosaic complex occurs. This is called an agglutination reaction and can be visualized microscopically and in some cases macroscopically.

**2. Precipitin formation:** This reaction requires specific antibodies, **precipitins**, that are formed in response to the introduction of soluble, nonparticulate antigens into host tissues. These antibodies, when present in serum, form a complex with the specific homologous nonparticulate antigen and result in a visible precipitate.

Advances in chemistry, especially immunology, have enabled us to study the interaction of antigens and immunoglobulins outside the body, in a laboratory setting. These advances have provided an immunological discipline known as **serology**, which studies these *in vitro* reactions that have diagnostic, therapeutic, and epidemiological implications.

In the experiments to follow, you will study several serological procedures based on the principles of agglutination and precipitin formation for the detection of serum antibodies or antigens. The techniques presented in these experiments span a spectrum of methods, ranging from basic reactions to more sophisticated forms of antigen and antibody interactions.



**Note that some of the experimental protocols use positive and negative controls provided in the test kits to demonstrate the desired immunological reactions. These controls do not represent the source of potential pathogens capable of inducing infection in students and instructional staff. The rationale for this design is that body fluids, particularly blood of unknown origin, may serve as a major vehicle for the transmission of infectious viral agents. Thus, our concern with the spread of AIDS and hepatitis precludes the use of blood as a test specimen in a college laboratory.**

It is further suggested that your instructor present experiments that use positive and negative control test kits as laboratory demonstrations. This will reduce the cost of the required materials, which may otherwise be prohibitive at many academic institutions, but will still allow you and your fellow students to observe the advances in immunological serology.

## LEARNING OBJECTIVE


Once you have completed this experiment, you should be able to

1. Demonstrate a precipitin reaction by means of the ring test.

## Principle

The **ring** or **interfacial test** is a simple serological technique that illustrates the precipitin reaction in solution. This antigen-antibody reaction can be demonstrated by the formation of a visible precipitate, a flocculent or granular turbidity, in the test fluid. Antiserum is introduced into a small-diameter test tube, and the antigen is then carefully added to form a distinct upper layer. Following a period of incubation of up to 4 hours, a ring of precipitate forms at the point of contact (interface) in the presence of the antigen-antibody reaction. The rate at which the visible ring forms depends on the concentration of antibodies in the serum and the concentration of the antigen.

To detect the precipitin reaction, a series of dilutions of the antigen is used because both insufficient (zone of antibody excess) and excessive (zone of antigen excess) amounts of antigen will prevent the formation of a visible precipitate (zone of equivalence), as shown in **Figure 71.1**. In addition, you will be able to determine the optimal antibody:antigen ratio by the presence of a pronounced layer of granulation at the interface of the antiserum and antigen solution. This immunological reaction is illustrated in **Figure 71.2**.

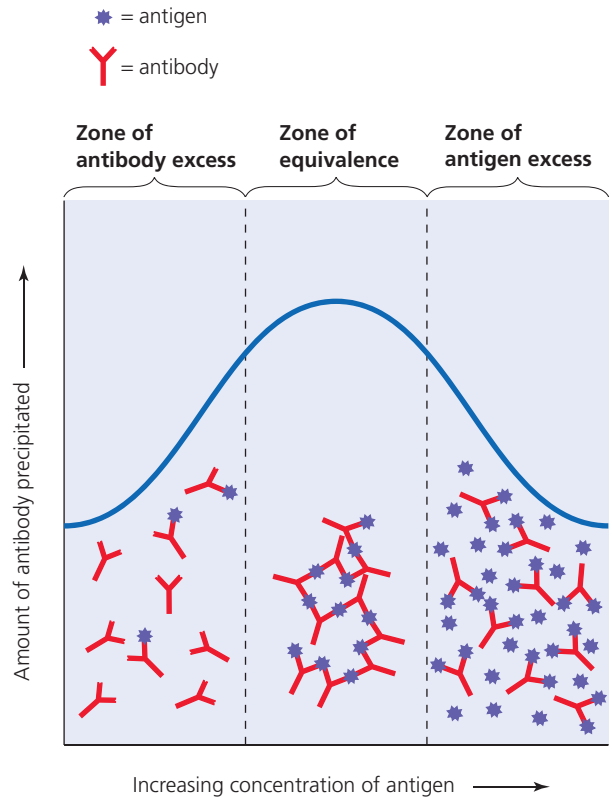


**Figure 71.1** The precipitin reaction

## CLINICAL APPLICATION

### Criminology

The precipitin reaction is a serological test in which an antibody reacts with a specific soluble antigen to form a visible precipitate ring in the tube. This test is mainly used today in criminology for the identification of human blood or other bloodstains, in cases of disputed parentage, and for the determination of the cause of death.

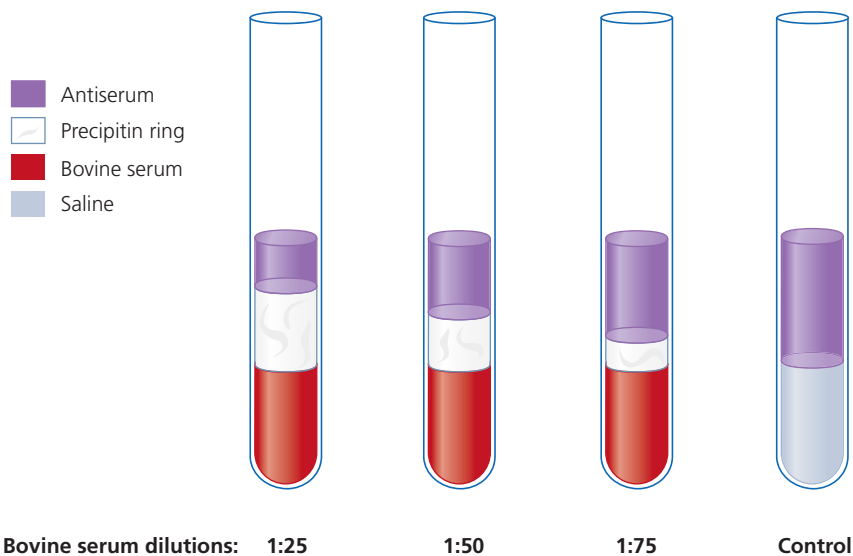

**Figure 71.2** Ring test: Precipitin reactions

## AT THE BENCH



## Materials

### Reagents

Physiological saline (0.85% NaCl), and commercially available bovine globulin antiserum and normal bovine serum diluted to 1:25, 1:50, and 1:75 with physiological saline. The normal bovine serum contains the antigen (bovine globulin), to which antibodies were made commercially in another animal species and provided as antiserum to bovine globulin.

### Equipment


Serological test tubes (8 × 75 mm), 0.5-ml pipettes, serological test tube rack, mechanical pipetting device, glassware marking pencil, and 37°C incubator.


## Procedure

1. Label three serological test tubes according to the antigen dilution to be used (1:25, 1:50, and 1:75) and the fourth test tube as a saline control.
2. Using a different 0.5-ml pipette, transfer 0.3 ml of each of the normal bovine serum dilutions into its appropriately labeled test tube.
3. Using a clean 0.5-ml pipette, transfer 0.3 ml of saline into the test tube labeled as control.
4. Carefully overlay all four test tubes with 0.3 ml of bovine globulin antiserum. To prevent mixing of the sera, tilt the test tube and allow the antiserum to run down the side of the test tube.
5. Incubate all test tubes for 30 minutes at 37°C.
6. Examine all test tubes for the development of a ring of precipitation at the interface. Indicate the presence or absence of a ring in the Lab Report.
7. Determine and record the antigen dilution that produced the greatest degree of precipitation; this is indicative of the optimal antibody:antigen ratio.





3.  Why is it essential to use a series of antigen dilutions in this procedure?

4.  How would you explain the absence of visible precipitate?

## LEARNING OBJECTIVE

Once you have completed this experiment, you should be able to

1. Demonstrate the agglutination reaction by means of the febrile antibody test and an antibody titer determination.

## Principle


The **febrile antibody test** is used in the diagnosis of diseases that produce febrile (fever) symptoms. Some of the microorganisms responsible for febrile conditions are salmonellae, brucellae, and rickettsiae. **Febrile antigens**, such as endotoxins, enzymes, and other toxic end products, elaborated by these organisms are used specifically to detect or exclude the homologous antibodies that develop in response to these antigens during infection.

In this procedure, the antigen is mixed on a slide with the serum being observed. Cellular clumping is indicative of the presence of homologous antibodies in the serum; the absence of homologous antibodies is indicated when there is no visible clumping. Only the febrile antigens and antibodies of *Salmonella* spp. will be used.

**Figure 72.1** shows a positive and a negative agglutination reaction.

The second part of this experiment is designed to illustrate that agglutination reactions such as the febrile antibody test can be used to identify an unknown microorganism through serotyping. A specific antiserum prepared in a susceptible, immunologically competent laboratory animal is mixed with a variety of unknown bacterial antigen preparations on slides. The bacterial antigen that is agglutinated by the antiserum is identified and confirmed to be the agent of infection.

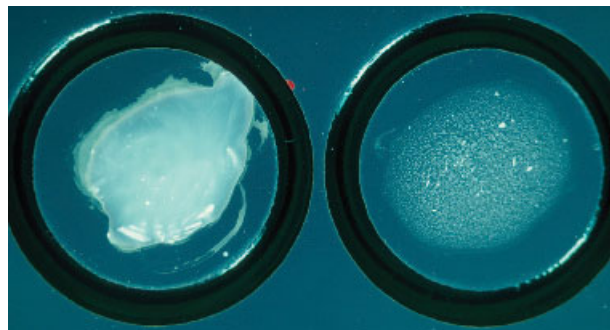These tests are strictly qualitative. A quantitative result can be obtained by performing the **antibody titer test**, which measures the



(a)

(b)

**Figure 72.1 Agglutination reaction.** (a) Visible clumping indicates the presence of homologous antibodies in the serum, and a positive reaction. (b) The lack of visible clumping indicates the absence of homologous antibodies, and a negative reaction.

concentration of an antibody in the serum and thus allows the physician to follow the course of an infection. The patient's serum is titrated (diluted), and the decreasing concentrations of the antiserum are mixed with a constant concentration of homologous antigen. The end point of the test will occur in the test tube containing the serum having the highest dilution showing agglutination.

## CLINICAL APPLICATION

### Febrile Disease Diagnosis

Febrile antigens are standardized suspensions of bacteria or bacterial antigens used to qualify or quantify specific serum antibodies that develop during some febrile infections. Based on the Widal Agglutination test for the diagnosis of typhoid fever, serum from patients can be tested for the presence of antibodies correlating to infectious diseases such as brucellosis, salmonellosis, and some rickettsial infections.



## Materials

### Cultures

Number-coded, washed saline suspensions of *Escherichia coli*, *Proteus vulgaris*, *Salmonella typhimurium*, and *Shigella dysenteriae*.

### Reagents

Physiological saline (0.85% NaCl), commercial preparations of *Salmonella typhimurium* H antigen, and *Salmonella typhimurium* H antiserum (Abcam, Inc.).

### Equipment

Bunsen burner, inoculating loop, glass microscope slides, 13 × 100-mm test tubes, sterile 1-ml pipettes, mechanical pipetting device, applicator sticks, glassware marking pencil, microscope, and waterbath.

## Procedure

### Febrile Antibody Test

1. With a glassware marking pencil, make two circular areas about ½ inch in diameter on a microscope slide. Label the circles A and B.
2. To Area A, add 1 drop of *S. typhimurium* H antigen and 1 drop of 0.85% saline. Mix the two with an applicator stick.
3. To Area B, add 1 drop of *S. typhimurium* H antigen and 1 drop of *S. typhimurium* H antiserum. Mix the two with a clean applicator stick.
4. Pick up the slide, and with two fingers of one hand, rock the slide back and forth.
5. Observe the slide both macroscopically and microscopically, under low power, for cellular clumping (agglutination).
6. Indicate the presence or absence of macroscopic and microscopic agglutination, and draw a representative field of Areas A and B in the Lab Report.

## Serological Identification of an Unknown Organism


1. Prepare two microscope slides as in the previous procedure. Label the four areas on the slides with the numbers of your four unknown cultures.
2. Into each area on both slides, place 1 drop of *S. typhimurium* H antiserum.
3. With a sterile inoculating loop, suspend a loopful of each number-coded unknown culture in the drop of antiserum in its appropriately labeled area on the slides.
4. Pick up the slides and slowly rock them back and forth.
5. Observe both slides macroscopically and microscopically, under low power, for agglutination.
6. In the Lab Report, indicate the presence or absence of macroscopic and microscopic agglutination in each of the suspensions. Also, indicate the suspension that is indicative of a homologous antigen-antibody reaction.

### Determination of Antibody Titer

Refer to **Figure 72.2** when reading the following instructions.

1. Place a row of 10 test tubes (13 × 100 mm) in a rack and number the tubes 1 through 10.
2. Pipette 1.8 ml of 0.85% saline into the first tube and 1 ml into each of the remaining nine tubes.
3. Into Tube 1, pipette 0.2 ml of *Salmonella typhimurium* H antiserum. Mix thoroughly by pulling the fluid up and down in the pipette. *Note: Avoid vigorous washing.* The antiserum has now been diluted 10 times (1:10).
4. Using a clean pipette, transfer 1 ml from Tube 1 to Tube 2 and mix thoroughly as described. Using the same pipette, transfer 1 ml from Tube 2 to Tube 3. Continue this procedure through Tube 9.
5. Discard 1 ml from Tube 9. Tube 10 will serve as the antigen control and therefore will not contain antiserum.

## PROCEDURE

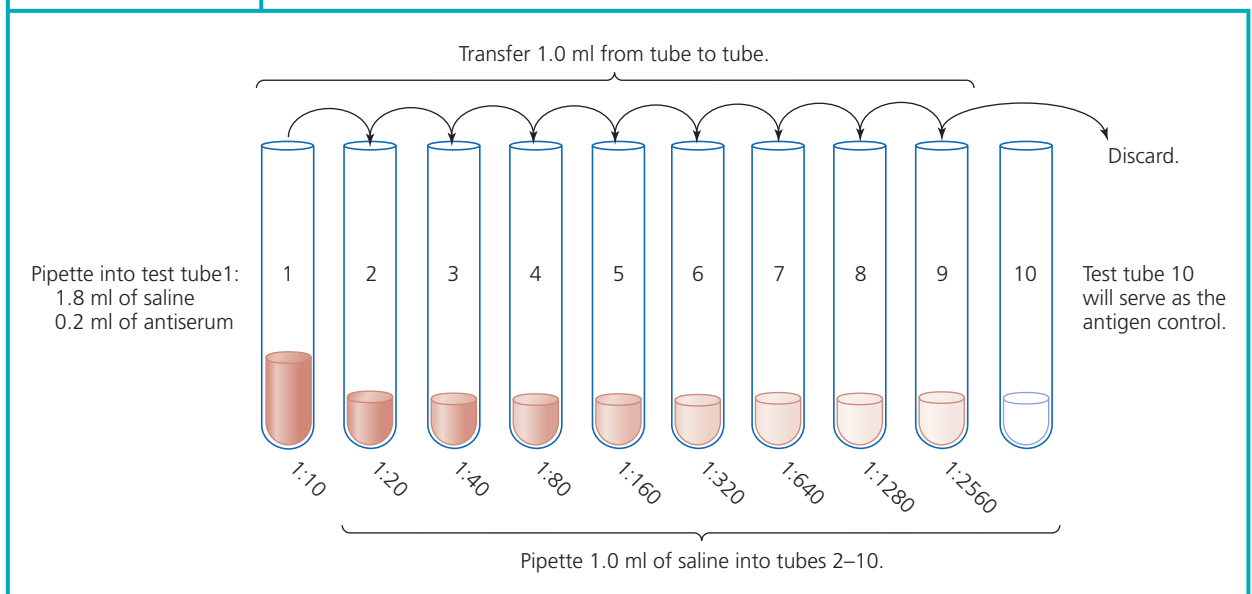

**Figure 72.2 Antibody titer test.** Serial dilution of *Salmonella typhimurium* H antibody

6. The antiserum has been diluted during this twofold dilution to give final dilutions of 1:10, 1:20, 1:40, 1:80, 1:160, 1:320, 1:640, 1:1280, and 1:2560.
7. Add 1 ml of the *Salmonella typhimurium* H antigen suspension adjusted to an absorbance of 0.5 at 600 nm to all tubes.
8. Mix the contents of the test tubes by shaking the rack vigorously.
9. Incubate the test tubes in a 55°C waterbath for 2 to 3 hours.
10. In the Lab Report, indicate the presence or absence of agglutination in each of the antiserum dilutions. Also, indicate the end point of the reaction.



*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results
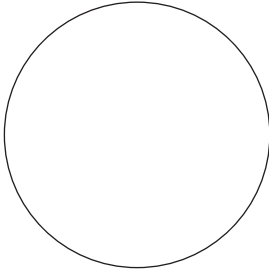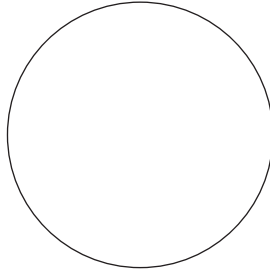
### Febrile Antibody Test

Draw the appearance of the mixture and of the control.	<b>Area A</b>	<b>Area B</b>
	 Saline <i>S. typhimurium</i> H antigen	 <i>S. typhimurium</i> H antiserum <i>S. typhimurium</i> H antigen
Macroscopic agglutination (+) or (-)		
Microscopic agglutination (+) or (-)		

### Serological Identification of an Unknown Organism

AGGLUTINATION				
Cell Antigen	Antiserum	Macroscopic (+) or (-)	Microscopic (+) or (-)	Homologous Antigen-Antibody Reaction
Unknown No: ____	<i>S. typhimurium</i> H			
Unknown No: ____	<i>S. typhimurium</i> H			
Unknown No: ____	<i>S. typhimurium</i> H			
Unknown No: ____	<i>S. typhimurium</i> H			

## Determination of Antibody Titer

Tube	Dilution	Agglutination	Titer
1	1:10		
2	1:20		
3	1:40		
4	1:80		
5	1:160		
6	1:320		
7	1:640		
8	1:1280		
9	1:2560		
10	Antigen control		

## Review Questions

1. What are febrile antibodies?

What is their clinical significance?

2. What is the purpose of determining an antibody titer?

3.  Why does the antibody titer determination use twofold dilutions of the antiserum rather than 10-fold dilutions?

## LEARNING OBJECTIVES

Once you have completed this experiment, you should be able to

1. Perform an antigen-antibody reaction by using fluorescein-labeled antibodies.

## Principle

The **fluorescent antibody technique**, introduced by Coombs in the mid-1950s, is a rapid and reliable procedure to demonstrate agglutination reactions. It has significant merit for use in the identification of microorganisms or their resultant antibodies. For example, the **fluorescent treponemal antibody-absorption (FTA-ABS) test** is used to diagnose syphilis using *Treponema pallidum* as the antigen for detection of syphilitic antibodies in the patient's serum.

This technique requires the use of a specific antibody that has been tagged at the Fc region with a fluorescent dye, such as fluorescein isocyanate or fluorescein isothiocyanate. When the tagged antibody is applied to the antigen preparation, as in the direct method, or to an antigen-antibody complex, as in the indirect method, a microprecipitate forms at the site of the antigen and exhibits a yellow-green fluorescence when viewed under a fluorescence microscope (Figure 73.1).

In this experiment, you will use the direct immunofluorescence method to demonstrate the antigen-antibody reaction and to identify an unknown antigen. The fluorescein-tagged antibody is added to a heat-fixed bacterial smear. A physiologically buffered saline wash removes any uncombined fluorescent antibody. The resultant agglutination reaction, if present, will be demonstrated by a green fluorescence when viewed under a fluorescence microscope.

## CLINICAL APPLICATION

## Diagnosis of Immunobullous Skin Diseases

Immunofluorescence testing is a necessary part of the accurate diagnosis of immunobullous skin diseases and other immune diseases of the skin resulting from autoantibodies to epidermal cells or tissues. Direct immunofluorescence is the diagnostic gold standard in the assessment of patients with these diseases. A fixed and sectioned skin biopsy from a patient's lesion is treated with fluorescent marker conjugated antibodies to human immunoglobulins to test for the presence of autoantibodies.

## AT THE BENCH



## Materials

## Cultures

24-hour brain heart infusion broth cultures of Group A *Streptococcus pyogenes* and Group D *Enterococcus faecalis*; numbered, unknown mixed broth cultures of Group A *S. pyogenes*/*Escherichia coli* and Group D *E. faecalis*/*Escherichia coli*.

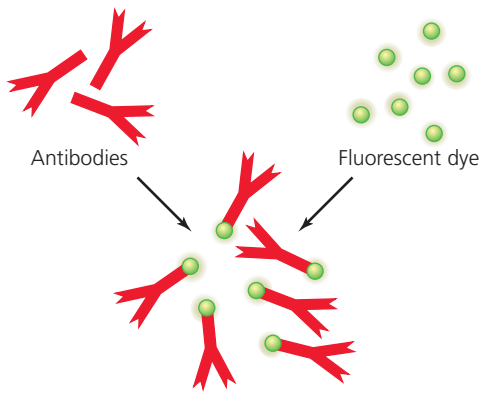
## Reagents

Fluorescent antibody *Streptococcus* Group A (Difco 2318-56-6), fluorescent antibody *Enterococcus* Group D (Difco 2319-56-5), phosphate-buffered saline, and buffered glycerol. *Note: These fluorescent antibodies can also be obtained from Abcam, Inc., Cambridge, Mass.*


## Equipment

Microscope slides, coverslips, Petri dishes, U-shaped glass rods to fit into Petri dishes, filter paper, Coplin jar, glassware marking






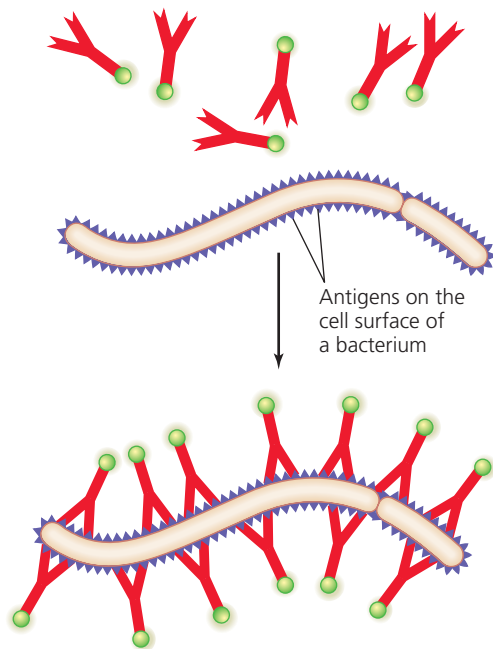(a) Antibodies are tagged with a fluorescent dye.



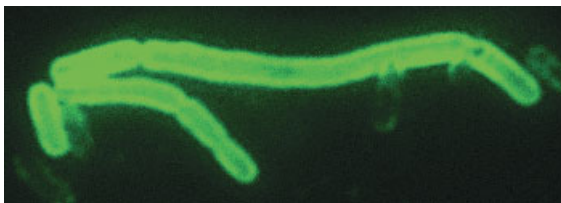(b) Tagged antibodies are introduced to a bacterium and attach to antigens on the cell surface.



(c) Photomicrograph of bacterial immunofluorescence. Only bacteria that bind with the antibodies will fluoresce.

**Figure 73.1** Fluorescent antibody technique

pencil, bibulous paper, sterile pipettes, and fluorescence microscope.

## Procedure

1. With a glassware marking pencil, label two slides *S. pyogenes* and *E. faecalis*, respectively. Divide the third slide in half and label it as “mixed unknown.”
2. Prepare a heat-fixed smear of each known test organism on its appropriately labeled slide. On the slide labeled “mixed unknown,” make a smear on each half of the slide using the unknown mixed culture.
3. On the slides labeled *S. pyogenes* and *E. faecalis*, add 1 drop of each respective fluorescent antibody and spread gently over the surface of the smear. On the slide of the mixed unknown smears, label one side FA-A and the other side FA-D, add 1 drop of each fluorescent antibody to its respectively labeled smear, and allow to spread evenly over the smears.
4. Place moistened filter paper in the Petri dishes, insert the U-shaped glass rod (for slide support), and place the prepared slides on the slide supports. Cover the Petri dishes and incubate for 35 minutes at 25°C.
5. Remove the slides from the Petri dishes and wash away excess antibody with 1% phosphate-buffered saline.
6. Immerse the slide in a Coplin jar containing 1% phosphate-buffered saline for 10 minutes at 25°C.
7. Blot the slides dry with bibulous paper.
8. To each slide add 1 drop of buffered glycerol, cover with a coverslip, and examine under a fluorescence microscope as directed by your instructor.
9. Record your observations in the Lab Report.

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

- Based on your microscopic observations, indicate in the chart the color of each smear and the presence (+) or absence (-) of fluorescence.

Culture	Color	Fluorescence (+) or (-)
Group A: <i>S. pyogenes</i>		
Group D: <i>E. faecalis</i>		
Mixed unknown with Group A antibodies		
Mixed unknown with Group D antibodies		

- Based on your observations, identify your mixed unknown culture.

Culture Number:

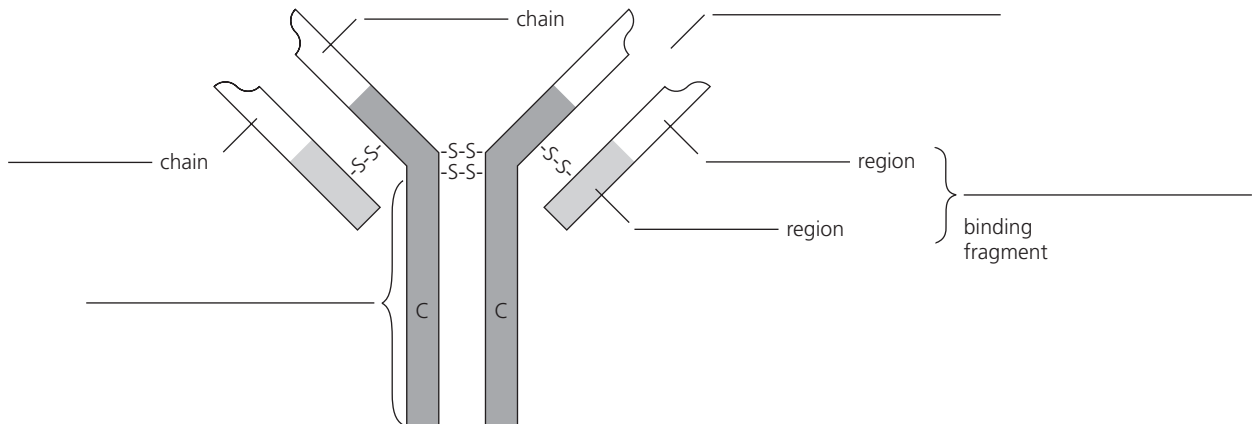
Organisms:

## Review Questions

- What does the presence of fluorescence indicate?
- Why is it necessary to wash away excess labeled antibody before viewing the preparation microscopically?

3. Briefly explain the direct fluorescent antibody procedure.

4. Label all parts of the antibody below, and show where it is tagged with fluorescein isothiocyanate.



5. Aside from syphilis, what diseases can be diagnosed through immunofluorescent techniques?

6. How would you use immunofluorescence to identify Group A *Streptococcus pyogenes* organisms isolated from a patient and differentiate them from other morphologically similar streptococci and staphylococci?

## LEARNING OBJECTIVE

Once you have completed this experiment, you should be able to

1. Demonstrate a method for the identification of either an antigen or an antibody by use of an enzyme-labeled antibody test procedure.

## Principle

The **enzyme-linked immunosorbent assay (ELISA)** procedure is a widely accepted method that is used for the detection of specific antigens or antibodies. The procedure is predicated on the use of an enzyme-linked (labeled) specific antibody to demonstrate the agglutination reaction for the interpretation of the test result. This test can be performed as a double-antibody technique or as an indirect immunosorbent assay. The former method is used for the detection of test antigens; the latter is used for the detection of the test antibodies. In both methods the reactions are carried out in a well of a plastic microtiter plate.

The double-antibody system requires that the unlabeled antibody be allowed to adsorb to the inner surface of the plastic well in the microtiter plate. Any unbound antibody is then washed away, and a specific test antigen is added to the well. If the antigen binds with the antibody adhering to the walls of the well, this immunocomplex will not be removed by the subsequent washing for the removal of any unbound antigen. An enzyme-linked antibody, specific for the antigen, is now added. If the antigen is present in the well, this labeled antibody binds to the antigen, forming an antibody-antigen-antibody complex. Any unbound enzyme-linked antibody is again removed by

washing. This is followed by the addition of a substrate that is capable of producing a colored end product upon its reaction with the enzyme. The resultant enzymatically produced color change may be observed by eye or spectrophotometrically.

The indirect immunosorbent test procedure is similar to the double-antibody technique in that it requires the use of an enzyme-linked antibody. However, an antigen, rather than an antibody, is adsorbed onto the inner surface of the well.

Enzyme-linked immunosorbent assays are used extensively for the diagnosis of human infectious diseases. Included among these are viral infections such as AIDS, influenza, respiratory syncytial viral infection, and rubella. Bacterial infections such as syphilis, brucellosis, salmonellosis, and cholera can also be ascertained by means of this technique. This procedure also can be used for the detection of drugs in blood or tissues.

In this experiment you will use the **Directigen™ Flu A Test** (Becton, Dickinson and Company) to demonstrate the application of an *in vitro* enzyme immunoassay. This rapid, qualitative test employs an enzyme immunomembrane filter assay to detect influenza A antigen extracted from nasopharyngeal or pharyngeal specimens of symptomatic patients. These specimens are added to a ColorPAC™ test device, and any influenza A antigen present is nonspecifically bound to the membrane surface. Detector enzyme conjugated to monoclonal antibodies specific for the influenza A nucleoprotein antigen is bound to the trapped antigen following its addition to the ColorPAC membrane. Two substrates are then added sequentially and allowed to incubate for 5 to 30 minutes prior to determination of the result.

In the experimental procedure to be followed, the positive control will simulate the nasopharyngeal specimen of a symptomatic patient and will be indicative of a positive result. A pharyngeal swab sample of an asymptomatic individual will be used to illustrate a negative result.

## CLINICAL APPLICATION

### Lyme Disease

The ELISA test is commonly used in the diagnosis of Lyme disease for the detection of antibodies to *Borrelia burgdorferi*. Because of the test's sensitivity, it can sometimes produce false positive results and it is not used as the sole basis for diagnosis of Lyme disease. It is generally followed up by a Western Blot test to confirm the diagnosis.

## AT THE BENCH



## Materials

### Cultures

Directigen Flu A positive control and pharyngeal swab specimen from an asymptomatic individual.

### Media

Per designated student group or demonstration: one test tube with 2-ml of sterile saline.

### Equipment

Directigen Flu A Test kit, sterile cotton swabs, sterile 0.2-ml (200- $\mu$ l) pipette, mechanical pipetting device, and disposable gloves.

## Procedure

*Note: This test may be performed as a demonstration for economic reasons or conservation of laboratory time.*



**Wear disposable gloves during the procedure.**

### Preparation of Negative Result by Use of a Pharyngeal Specimen

1. Using a sterile cotton swab, obtain a specimen from the palatine tonsil by rotating the swab vigorously over its surface.
2. Immerse the cotton swab into a test tube containing 2 ml of sterile saline. Mix well. Remove as much liquid from the swab as possible by pressing the swab against the inner surface of the tube. Discard the swab into a container of disinfectant.

3. Using a 0.2-ml (200- $\mu$ l) pipette and a mechanical pipetting device, transfer 124  $\mu$ l of the pharyngeal specimen into a DispensTube™ provided in the kit.
4. Gently mix and add 8 drops of Reagent A into the DispensTube. Mix well.
5. Insert a tip, provided in the kit, into the DispensTube. Dispense all of the extracted specimen into the ColorPAC test well in drops with the sealed flow controller in position. Allow for complete adsorption.
6. Gently mix and rapidly add drops of Reagent 1 until the test well is filled. Allow sufficient time for complete adsorption.
7. Remove the flow controller from the ColorPAC well and discard it into a container of disinfectant.
8. Gently mix and add 4 drops of Reagent 2 onto the ColorPAC membrane. Allow sufficient time for complete adsorption.
9. Gently mix and add 4 drops of Reagent 3 onto the ColorPAC membrane. Allow sufficient time for complete adsorption. Let stand for 2 minutes.
10. Rapidly add enough drops of Reagent 4 to fill the ColorPAC well. Allow sufficient time for complete adsorption.
11. Gently mix and add 4 drops of Reagent 5 onto the ColorPAC membrane. Allow sufficient time for complete adsorption.
12. Gently mix and add 4 drops of Reagent 6 onto the ColorPAC membrane. Allow sufficient time for complete adsorption. *Note: The membrane will turn yellow.*
13. Gently mix and add 4 drops of Reagent 7 onto the ColorPAC membrane. Allow sufficient time for complete adsorption.
14. Wait at least 5 minutes, but no longer than 30 minutes, and then read the results in a well-lighted area.
15. Observe the appearance of the inner surface of the ColorPAC test wells and record your results in the Lab Report.

### Preparation of Positive Result by Use of Positive Control

1. Dispense 4 drops of the positive control, provided in the test kit, into a DispensTube.
2. Repeat Steps 4 through 15 as outlined previously for the preparation of the negative pharyngeal specimen.

Name: \_\_\_\_\_

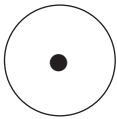
Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results


Record your results below based on the following interpretations of your observations:



**Positive test** (antigen present): The appearance of a purple triangle (of any intensity) on the ColorPAC membrane indicates the presence of the influenza A antigen in the specimen. A purple dot may be evident in the center of the triangle. The background area should be grayish white.




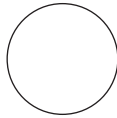**Negative test** (no antigen detected): The appearance of a purple dot on the ColorPAC membrane indicates the absence of the influenza A antigen in the specimen. The background area should be grayish white.




**Uninterpretable test:** The absence of a purple dot, a purple triangle, or an incomplete purple triangle indicates an uninterpretable test.

**Negative pharyngeal specimen:**

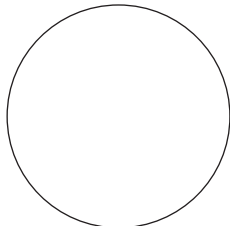

Result

**Positive control specimen:**

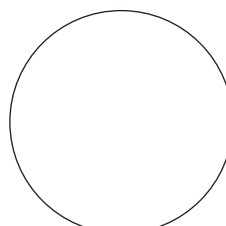

Result

## Review Question

-  Why is the ELISA test used to screen human serum for the AIDS virus, while the Western blot procedure is used only as the confirmation test?

*This page intentionally left blank*

Sexually transmitted diseases (STDs) represent a diverse group of infectious syndromes that share the same mode of transmission, direct sexual contact. Their etiological agents represent a broad spectrum of pathogenic microorganisms that include bacteria, viruses, yeasts, and protozoa. The bacterial STDs include **gonorrhea**, **syphilis**, **nongonococcal urethritis**, and **lymphogranuloma venereum**. The representative viral infections are **genital herpes**, **genital warts**, **hepatitis B**, and the latest member of this group, **AIDS**. The protozoal and fungal infections, namely **trichomoniasis** and **candidiasis**, are diseases of lesser magnitude in the spectrum of STDs.

The experimental procedures that follow were chosen to demonstrate some of the rapid tests that are currently available for the diagnosis of selected STDs, specifically syphilis, genital herpes, and the chlamydial infections. In the methods that follow, you will perform modified procedures in the absence of clinical specimens. Commercially available positive and negative controls will be used to simulate clinical materials. *It is suggested that any of these tests, if performed, should be done as demonstrations.*

## PART A Rapid Plasma Reagin Test for Syphilis

### LEARNING OBJECTIVE

*Once you have completed this experiment, you should be able to*

1. Perform a rapid screening procedure for diagnosis of syphilis.

### Principle

*Treponema pallidum*, the causative agent of **syphilis**, is a tightly coiled, highly motile, delicate spirochete that can be cultivated only in rabbit

tissue cultures or rabbit testes. The organisms are resistant to common staining procedures and are best observed under darkfield microscopy.


Syphilis is a systemic infection that, if untreated, progresses through three clinical stages. The first stage, primary syphilis, is characterized by the formation of a painless papule, called a **chancre**, at the site of infection. Secondary syphilis represents the systemic extension of the infection and presents itself in the form of a **maculopapular rash**, malaise, and lymphadenopathy. Following this stage, the disease becomes self-limiting, and the patient appears asymptomatic until the development of tertiary syphilis. In this final stage, life-threatening complications may develop as a result of the extensive cardiovascular and nervous tissue damage that has ensued.

The **rapid plasma reagin (RPR) test**, which has to a large extent replaced the VDRL (Venereal Disease Research Laboratory) agglutination test, determines the presence of **reagin**, the nonspecific antibody present in the plasma of individuals with a syphilitic infection. The reagin appears in the plasma within 2 weeks of infection and will remain at high concentrations until the disease is eradicated. In the RPR test, if the reagin is present in the blood, it will react with a soluble antigen bound to carbon particles to produce a macroscopically visible antigen, or carbon-antibody complex. This procedure has several advantages over the VDRL test:

1. The serum does not require inactivation by heat for 30 minutes.
2. The serum may be obtained from a finger puncture, unlike the VDRL test, which requires a venous blood sample.
3. The required materials, which include the antigen suspension with a dispensing bottle, diagnostic cards, and capillary pipettes, are all contained in individual kits that do not require additional equipment and are disposable.

In the qualitative form of the RPR test, the patient's blood serum and the carbon-bound antigen suspension are mixed within a circle on the diagnostic card. In the presence of a positive



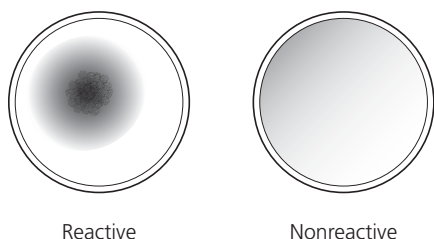

**Figure 75.1** Test card showing results of the rapid plasma reagin test

(reactive) serum, the antigen-antibody complex will produce a macroscopically visible black agglutination reaction. The macroscopic appearance of a light-gray suspension, devoid of any form of agglutination, is indicative of a negative (nonreactive) serum (**Figure 75.1**).

Since this is a nonspecific test, false-positive results may be obtained. It is believed that the reagin is an antibody against tissue lipids in general. Therefore, it may be present in uninfected individuals due to the release of lipids resulting from normally occurring wear and tear of body tissues. It has also been found that serum levels of reagin are elevated during the course of other infectious diseases such as viral pneumonia, lupus erythematosus, infectious mononucleosis, yaws, and pinta. The serum of patients with a reactive RPR result is subjected to additional serological testing, such as the FTA-ABS (fluorescent treponemal antibody-absorption) test, or the TPI (*Treponema pallidum* immobilization) test, using the *Treponema pallidum* bacterium as an antigen to detect specific antibodies that are also present in the serum during syphilitic infection.

## CLINICAL APPLICATION

### Spirochetes

The rapid plasma regain (RPR) test detects non-specific antibodies in the blood of patients that may indicate the spirochete *Treponema pallidum* that causes syphilis. This test is used to screen asymptomatic patients, diagnose symptomatic infection, and track the progress of disease over the treatment period. High incidence of false positives due to cross reactivity and false negatives due to low antibody titers requires further testing using the Venereal Disease Research Lab (VDRL) test in many clinical labs.

## AT THE BENCH



## Materials

### Reagents

Commercially prepared syphilitic serum 4+ and nonsyphilitic serum.

### Equipment

RPR test kit (Inverness Medical Professional Diagnostics), disposable gloves, and rotating machine (optional).

## Procedure

**⚠ Wear disposable gloves.**

1. Label circles on the diagnostic plastic card as reactive and nonreactive.
2. Using a capillary pipette with an attached rubber bulb, draw the reactive serum up to the indicated mark (0.05 ml).
3. Expel the serum directly onto the card in the circle labeled reactive serum. With a clean applicator stick, spread the serum to fill the entire circle.
4. Repeat Steps 2 and 3 for the nonreactive serum.
5. Shake the dispensing bottle to mix the suspension. Hold the bottle with attached 20-gauge needle in a vertical position and dispense 1 drop onto each circle containing the test serum.
6. If a mechanical rotator is available, place the card on the rotator set at 100 rpm or rotate the card back and forth manually for 8 minutes.
7. In the presence of direct light, while tilting the card back and forth, determine the presence or absence of black clumping in each of the serum-antigen mixtures. Record your observations and the reaction as (+) or (-) in the Lab Report.

## PART B Genital Herpes: Isolation and Identification of Herpes Simplex Virus

### LEARNING OBJECTIVE

*Once you have completed this experiment, you should be able to*

1. Perform a tissue culture procedure for the growth and identification of the herpes simplex virus.

### Principle

The double-stranded DNA herpes simplex virus (HSV) is the etiological agent of a variety of human infections. Included among these are **herpes labialis**, fever blisters around the lips; **keratoconjunctivitis**, infection of the eyes; **herpes genitalis**, eruptions on the genitalia; **herpes encephalitis**, a severe infection of the brain; and **neonatal herpes**. The herpes simplex virus is divided into two antigenically distinct groups, **HSV-1** and **HSV-2**. The former is most frequently implicated with infections above the waist, whereas the latter is predominantly responsible for genital infections.

Primary infection with HSV-2 manifests itself with the appearance of vesicular lesions, characterized by itching, tingling, or burning sensations on or within the male and female genitalia. These vesicles regress spontaneously within 2 weeks. Following this symptomatic phase, however, the virus reverts to a latent state in the sacral nerves and remains quiescent until exacerbated by some environmental factor. With no chemotherapeutic cure presently available, recurrent genital herpes with subclinical symptoms is common.

Detection of the herpes simplex virus requires the use of tissue culture techniques. The presence of the virus is then determined by the development of cytopathogenic effects in these cultures, such as the detection of intranuclear inclusion bodies. In recent years, these time-consuming, specialized procedures have been greatly facilitated by the availability of immunoenzymatic reagents for the identification of this clinically significant virus.

The **Cellmatics™ HSV Detection System** is a self-contained system providing for both the growth and identification of the virus from clinical

specimens. In this procedure, the provided tissue culture tubes are inoculated with the clinical sample. Following a 24-hour incubation period and fixation, the presence of HSV antigens is determined by the addition of anti-HSV antibodies, which specifically bind to the HSV antigens. To demonstrate this antigen-antibody complex, a secondary antibody, substrate, and chromogen are added. Following this staining process, HSV-positive cultures viewed microscopically will exhibit brown-black areas of viral infection on a clear background of unstained cells.

In this exercise, you will perform a modified procedure. In the absence of a clinical specimen, the actual culturing and fixation process will not be performed. Instead, the positive and negative commercially available controls will be used to simulate the clinical samples.

### CLINICAL APPLICATION

#### Culturing HSV

Genital herpes is caused by herpes simplex-2 virus (HSV-2). It is spread from person to person during sexual contact. The infection is transmitted by means of viral shedding, which may occur even when no signs or symptoms appear. A swab sample of a vesicular lesion or from the site of a previous lesion is taken from the patient and sent to the clinical or infectious disease lab for identification. This virus can only be grown in tissue culture and not cultivated on or in other laboratory media.

### AT THE BENCH



### Materials

#### Cultures

Cellmatics HSV Positive and Negative Control Tubes.

#### Reagents

Cellmatics Immunodiagnostic Reagents Kit, distilled water (Difco Labs).

#### Equipment

5-ml pipettes, mechanical pipetting device, and microscope.

## Procedure

1. Warm immunodiagnostic reagents to room temperature.
2. Drain all fluid from the positive and negative control tubes.
3. Using a 5-ml pipette, wash the culture tubes twice with 5 ml of distilled water and drain.  
*Note: When washing, exercise care to prevent disruption of the monolayer.*
4. Add 10 drops of primary antiserum (Vial 1).  
*Note: When adding reagents, hold the vials vertically to ensure proper delivery.*
5. Incubate the *tightly capped* tubes in a *horizontal* position for 15 minutes at 37°C. To ensure complete coverage of the monolayer, occasionally rock the tubes gently during incubation.
6. Wash three times with 5 ml of distilled water and drain.
7. Add 10 drops of secondary antibody (Vial 2).
8. Incubate for 15 minutes at 37°C as described in Step 5.
9. Wash three times with 5 ml of distilled water and drain.
10. Add 10 drops of substrate (Vial 3) and 2 drops of chromogen (Vial 4). Mix gently.
11. Incubate for 15 minutes at 37°C as described in Step 5.
12. Wash three times with 5 ml of distilled water and drain.
13. Examine microscopically for the presence of stained cells at 40× and 100× magnifications. Scan the entire stained monolayer of both culture tubes for the presence of brown to blackish-brown stained cells. HSV infection is indicated by the presence of dark-colored cells when viewed against an unstained background of normal cells.
14. Record your observations in the Lab Report.

## PART C Detection of Sexually Transmitted Chlamydial Diseases

### LEARNING OBJECTIVE

*Once you have completed this experiment, you should be able to*

1. Perform an immunofluorescent procedure for diagnosis of *Chlamydia* infections.

## Principle

Members of the genus *Chlamydia* are a group of obligate intracellular parasites. Although they were once believed to be viruses, their morphological and physiological characteristics more closely resemble bacteria, so they are now considered small bacteria. Chlamydiae are gram-negative, nonmotile, thick-walled, spherical organisms possessing both DNA and RNA that reproduce by means of binary fission. Their dependence on living tissues for cultivation and their lack of an ATP-generating system emulate the characteristics of viruses, but their bacterial nature is further affirmed by their sensitivity to antibiotic therapy. *Chlamydia trachomatis*, the human pathogen, is now recognized to be responsible for two sexually transmitted diseases, **nongonococcal urethritis** (NGU) and **lymphogranuloma venereum** (LGV). The incidence of both diseases in contemporary society is increasing dramatically.

NGU is a urethritis (inflammation of the urethra) with symptoms similar to, but less severe than, those of gonorrhea. Undiagnosed and untreated infections may lead to **epididymitis** and **proctitis** in men and **cervicitis, salpingitis, and pelvic inflammatory disease** in women. Nongonococcal urethritis is also caused by other bacteria such as *Ureaplasma urealyticum*, and

*Mycoplasma hominis*, as well as the protozoan *Trichomonas vaginalis*. LGV, the most severe of the genital chlamydial infections, initially develops with a painless lesion at the portal of entry, the genitalia. Systemic involvement is evidenced by swelling of the regional lymph nodes, which become tender and suppurative before disseminating the organisms to other tissues. In the absence of chemotherapeutic intervention, scarring of the lymphatic vessels can cause their obstruction, leading to **elephantiasis**, enlargement of the external genitalia in men, and narrowing of the rectum in women.

The **MicroTrak® Direct Specimen Test** is a rapid, immunofluorescent procedure for the detection of *C. trachomatis*. The procedure circumvents the need to culture these organisms in susceptible tissues prior to their identification. This slide test is designed to detect elementary bodies, the infectious particles produced during the life cycle of this organism, by the use of a staining reagent, a fluorescein-labeled monoclonal antibody specific for the principal protein of the *C. trachomatis* outer-membrane. In this procedure, a slide smear is prepared from the clinical specimen. Following fixation, when the slide is exposed to the Direct Specimen Reagent, the antibody binds to the organisms. Their presence is then determined by the appearance of apple-green chlamydiae against a red background of counterstained cells when viewed under a fluorescent microscope.

## CLINICAL APPLICATION

### Treating Chlamydia

Chlamydial infections are the most commonly reported sexually transmitted diseases. More than fifty million infections occur worldwide with three million cases occurring in the United States annually. Any sexually active person can contract Chlamydia, but it most frequently occurs in teenagers and young adults. The incidence appears higher in females than in males. Chlamydia may be transmitted by an infected mother to her newborn during birth. If a mother's medical or sexual history indicates possible exposure, a urogenital swab will be used to collect a sample for testing.

## AT THE BENCH



## Materials

### Cultures

Commercially prepared positive and negative control slides.

### Reagents

MicroTrak Direct Specimen Test for *Chlamydia trachomatis* (VWR Scientific).

### Equipment

Fluorescent microscope.

## Procedure

1. Stain the positive and negative control slides with the MicroTrak reagent for 15 minutes.
2. Incubate slides for 15 minutes.
3. Rinse the slides with distilled water.
4. Air-dry the slides.
5. Examine the slides under a fluorescent microscope for the presence of apple-green particles indicative of chlamydiae. The particles are evident against a reddish background of counterstained cells.
6. Record your results in the Lab Report.

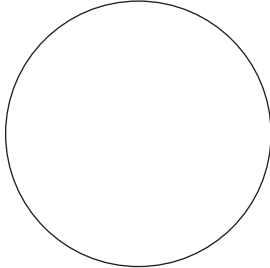
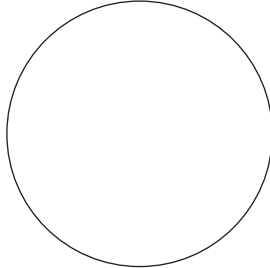
*This page intentionally left blank*

Name: \_\_\_\_\_

Date: \_\_\_\_\_ Section: \_\_\_\_\_

## Observations and Results

### PART A: Rapid Plasma Reagin Test for Syphilis

	Reactive Serum	Nonreactive Serum
Appearance of serum-antigen mixture		
Reaction (+) or (-)		
Draw the observed reaction.		

### PART B: Genital Herpes: Isolation and Identification of Herpes Simplex Virus

Indicate in the chart below the presence (+) or absence (-) of dark-stained patches.

CONTROL CULTURES			
Negative		Positive	
40×	100×	40×	100×

### PART C: Detection of Sexually Transmitted Chlamydial Diseases


Record your results below, indicating the presence (+) or absence (-) of the apple-green chlamydiae on each of the slides.

Positive control slide: \_\_\_\_\_

Negative control slide: \_\_\_\_\_

## Review Questions

1.  Why is an adult who has a high antibody titer to herpes simplex virus 2 (HSV-2) subject to recurrent genital herpes infections?

2.  A 20-year-old college student was informed following a physical examination that her blood test for syphilis was reactive. She indicated that she was a virgin and had never received a blood transfusion. A repeat RPR test was positive, but the TPI test was negative. How would you explain these bizarre results, and what is the clinical status of this patient?

Microbiologists are required to perform a variety of laboratory techniques, including preparing and diluting solutions; expressing concentrations of chemicals, antibiotics, and antiseptics in solution; making quantitative determinations of cell populations based on the standard method for plate counting; and making serial dilutions to accommodate the latter procedure. These techniques commonly involve the use of very large or very small numbers (e.g., 9,000,000,000 or 0.000000009), which can be so cumbersome to manipulate that errors may result. Therefore, it is essential for microbiologists to have a good command of scientific exponential notation known as **scientific notation**.

The basis for this system is predicated on the fact that all numbers can be expressed as the product of two numbers, one of which is the power of the number 10. In scientific notation, the small superscript number next to the 10 is called the **exponent**. Positive exponents tell us how many times the number must be multiplied by 10, while negative exponents indicate how many times a number must be divided by 10 (that is, multiplied by one-tenth).

For example, a number written using the exponential form designated as scientific notation would appear as  $7.5 \times 10^3$ , meaning that  $7.5 \times 10 \times 10 \times 10 = 7500$ . Appendix Table 1.1 shows both large and small numbers written in the exponential form.

Appendix Table 1.1		Scientific (Exponential) Notation	
NUMBERS GREATER THAN ONE		NUMBERS LESS THAN ONE	
1,000,000,000	$= 1 \times 10^9$	0.000 000 001	$= 1 \times 10^{-9}$
100,000,000	$= 1 \times 10^8$	0.000 000 01	$= 1 \times 10^{-8}$
10,000,000	$= 1 \times 10^7$	0.000 000 1	$= 1 \times 10^{-7}$
1,000,000	$= 1 \times 10^6$	0.000 001	$= 1 \times 10^{-6}$
100,000	$= 1 \times 10^5$	0.000 01	$= 1 \times 10^{-5}$
10,000	$= 1 \times 10^4$	0.000 1	$= 1 \times 10^{-4}$
1000	$= 1 \times 10^3$	0.001	$= 1 \times 10^{-3}$
100	$= 1 \times 10^2$	0.01	$= 1 \times 10^{-2}$
10	$= 1 \times 10^1$	0.1	$= 1 \times 10^{-1}$
1	$= 1 \times 10^0$	1	$= 1 \times 10^0$

*Note:* The exponent to which the power of 10 is raised is equal to the number of zeros to the right of 1.

*Note:* The exponent to which the power of 10 is raised is equal to the number of zeros to the left of 1 plus 1.



## Multiplication

**Rule:** To multiply two numbers that are written in scientific notation (exponential form) you must **add** the exponents.

Using numbers larger than 1:

$$75 \times 1200 = 90,000$$

*Scientific notation:*  $(7.5 \times 10^1) \times (1.2 \times 10^3) = 9 \times 10^4$

*Addition of exponents:*  $1 + 3 = 4$

Using numbers less than 1:

$$0.75 \times 1200 = 900$$

*Scientific notation:*  $(7.5 \times 10^{-1}) \times (1.200 \times 10^3) = 9 \times 10^2$

*Addition of exponents:*  $(-1 + 3 = 2)$

$$0.75 \times 0.12 = 0.09$$

*Scientific notation:*  $(7.5 \times 10^{-1}) \times (1.2 \times 10^{-1}) = 9 \times 10^{-2}$

*Addition of exponents:*  $(-1) + (-1) = -2$

## Division

**Rule:** To divide two numbers in scientific notation, you must **subtract** the exponents.

$$75,000 \div 1,200,000 = 0.0625$$

*Scientific notation:*  $(7.5 \times 10^4) \div (1.2 \times 10^6) = 6.25 \times 10^{-2}$

*Subtraction of exponents:*  $(4 - 6 = -2)$

$$7,500 \div .012 = 625,000$$

*Scientific notation:*  $(7.5 \times 10^3) \div (1.2 \times 10^{-2}) = 6.25 \times 10^5$

*Subtraction of exponents:*  $3 - (-2) = 5$

As you practice the use of scientific notation with large and small numbers, you will become more proficient and more comfortable with this system of scientific calculation.

# Methods for the Preparation of Dilutions

In microbiology laboratories as in other science laboratories, solutions must be diluted to achieve a desired final concentration of the active material contained in that solution. A **solution** may be defined as a mixture of two or more substances (**solute**) in which the molecules of the solute are evenly distributed and will not separate on standing, or precipitate from the solution. Solutes are dissolved in a solvent or diluent, such as water, alcohol, or some other vehicle in which the solute is soluble. Solutions are usually referred to as stock solutions and may be diluted by a variety of methods, depending upon the experimental requirements. Some of these methods are listed as follows:

1. A **dilution factor** must be determined first in order to dilute a solution. This dilution factor tells us how many times a solution must be diluted and is calculated by dividing the **initial concentration (IC)** of the solution by the **final concentration (FC)** desired.

$$\text{IC} \div \text{FC} = \text{dilution factor}$$

**Example:** You wish to dilute a 10% stock solution to a final concentration of 2%.

$$10\% \div 2\% = 5 (\text{dilution factor})$$

Take 1.0 ml of the 10% stock solution plus 4.0 ml of diluent (solvent), which equals a total of 5.0 ml. Thus each ml of the final solution will contain 2% solute.

2. Another method is used when a specific volume composed of a specific concentration is required.
  - a.  $\frac{\text{IC}}{\text{FC}} = 10 (\text{dilution factor})$
  - b.  $\frac{\text{volume needed}}{\text{concentration required}} = \text{amount of initial solution needed}$
  - c.  $\text{volume needed} - \text{amount of initial solution} = \text{amount of diluent}$

**Example:** You have a 50% concentrated solution and you need 200 ml of a 5% solution.

- a.  $\frac{\text{IC}}{\text{FC}} = \frac{50\%}{5\%} = 10 (\text{dilution factor})$
  - b.  $\frac{\text{volume needed}}{\text{concentration required}} = \frac{200 \text{ ml}}{5\%} = 40 \text{ ml}$
  - c. 40 ml of 50% IC + 160 ml of diluent = 200 ml of a solution such that each ml will contain 5% solute rather than the original 50% in the stock solution
3. The ability to prepare large dilutions is absolutely essential for work in the microbiology laboratory. This method requires that large dilutions be prepared in two steps.

**Example:** A solution contains 1.0 g per ml of an active material and needs to be diluted to a final concentration of 1.0  $\mu\text{g}$  per ml. A 1,000,000 ( $1 \times 10^6$ )-fold dilution must be made. It is not practical to make such a dilution in one step since 999,999 ml of diluent would be required. This type of dilution is made as follows:

- a. Dilute 1 ml of the stock solution 1000 times:  
1.0 ml + 999 ml of diluent = 1000  $\mu\text{g/ml}$
  - b. Dilute the solution containing 1000  $\mu\text{g/ml}$  another 1000 times:  
1 ml of 1000  $\mu\text{g/ml}$  + 999 ml diluent = 1.0  $\mu\text{g/ml}$
4. When working with large molecules such as proteins, there will be times when you will be required to make large dilutions of the sample to be contained in a specific volume.

**Example:** You need to make 50 ml of a 1/20,000 dilution of albumin.

$$\frac{\text{final dilution}}{\text{volume needed}} = \frac{20,000}{50} = 400 (\text{dilution factor})$$

- a. 1.0 ml of albumin + 399 ml of diluent = 1/400 dilution

b.  $1.0 \text{ ml of a } 1/400 \text{ dilution} + 49 \text{ ml of diluent} = 50 \text{ ml of a solution; each ml contains } 1/20,000 \text{ of albumin.}$

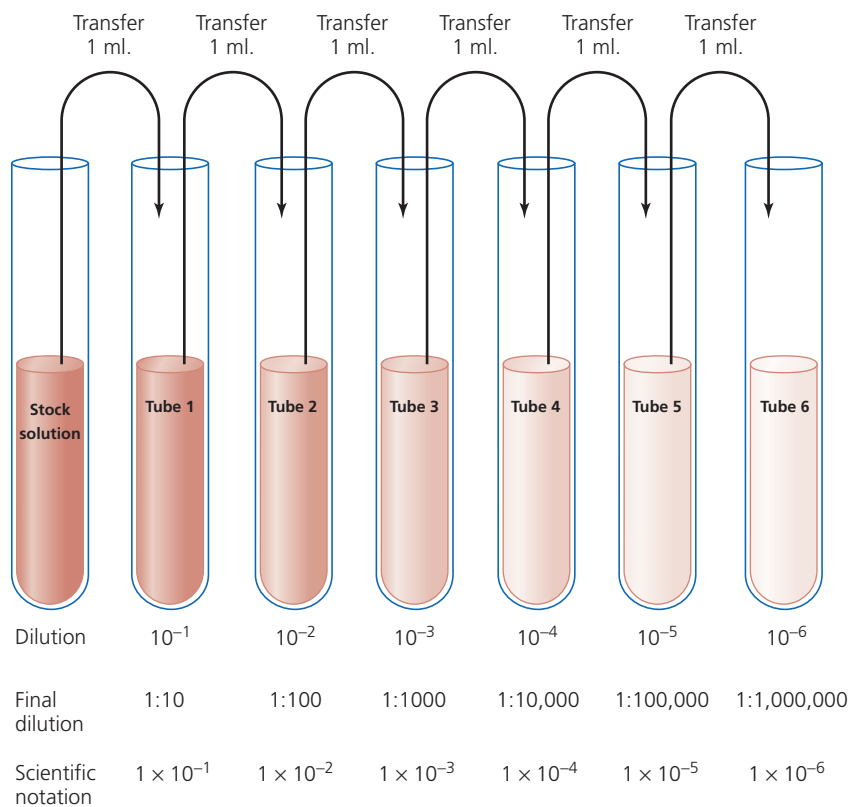
$$50 (\text{volume}) \times 400 (\text{dilution factor}) = 20,000$$

5. Perhaps the most useful type of dilution used in microbiology and immunology is the **serial dilution**. This is essential when small volumes of material are needed. This type of dilution procedure has many uses in the microbiology laboratory, especially for the determination of the total number of cells in culture (Experiments 20 and 21), the number of viral plaques found in suspensions of viruses (Experiment 39), the antibody titer (Experiment 72), and in other immunological studies. The procedure requires the use of dilution blanks containing a known volume of diluent (distilled water, saline, etc.) to which a specific volume of the sample is added. To facilitate the ease of calculations, dilutions are usually made in multiples of 10. For example: 1.0 ml of a sample is added to a 9.0-ml dilution blank (1.0 ml + 9.0 ml = 10) and is recorded as a 1:10 dilution.

It has been statistically determined that greater accuracy is achieved with very large dilutions made from a series of smaller dilutions. The procedure for the performance of a serial dilution has been explained and illustrated in Experiment 20. For the convenience of the student, it is illustrated in **Appendix Figure 2.1**:

1. All dilution blanks contain 9.0 ml of diluent.
2. A fresh pipette is used for each dilution, and the used pipettes are placed in a beaker of disinfectant.
3. After delivery of the sample, the tubes are mixed thoroughly before the next dilution is made.
4. Pippetting by mouth is not allowed. Only mechanical pipette aspirators may be used.

The stock solution in Appendix Figure 2.1 has been diluted 1 million times. In other words, 1.0 ml from Tube 6 will contain 1/1,000,000 of the sample contained in the stock solution.



**Appendix Figure 2.1** Serial dilution

The formulas of the media used in the exercises in this manual are listed alphabetically in grams per liter of distilled water unless otherwise specified. Sterilization of the media is accomplished by autoclaving at 15 lb pressure for 15 minutes unless otherwise specified. Most of the media are available commercially in powdered form, with specific instructions for their preparation and sterilization.

### Ammonium sulfate broth (pH 7.3)

Ammonium sulfate	2.0
Magnesium sulfate · 7H <sub>2</sub> O	0.5
Ferric sulfate · 7H <sub>2</sub> O	0.03
Sodium chloride	0.3
Magnesium carbonate	10.0
Dipotassium hydrogen phosphate	1.0

### Bacteriophage broth 10× (pH 7.6)

Peptone	100.0
Beef extract	30.0
Yeast extract	50.0
Sodium chloride	25.0
Potassium dihydrogen phosphate	80.0

### Basal salts agar\* and broth (pH 7.0)

0.5 M sodium diphosphate	100.0 ml
1.0 M potassium dihydrogen phosphate	100.0 ml
Distilled water	800.0 ml
0.1 M calcium chloride	1.0 ml
1.0 M magnesium sulfate	1.0 ml
Ammonium sulfate	2.0
*Agar	15.0

*Note: Swirl until completely dissolved, autoclave, and cool. Aseptically add 10.0 ml of 1% sterile glucose.*

### Bile esculin (pH 6.6)

Beef extract	3.0
Peptone	5.0
Esculin	1.0
Oxgall	40.0
Ferric citrate	0.5
Agar	15.0

### Blood agar (pH 7.3)

Infusion from beef heart	500.0
Tryptose	10.0
Sodium chloride	5.0
Agar	15.0

*Note: Dissolve the above ingredients and autoclave. Cool the sterile blood agar base to 45°C to 50°C. Aseptically add 50 ml of sterile defibrinated blood. Mix thoroughly, avoiding accumulation of air bubbles. Dispense into sterile tubes or plates while liquid.*

### Brain heart infusion (pH 7.4)

Infusion from calf brain	200.0
Infusion from beef heart	250.0
Peptone	10.0
Dextrose	2.0
Sodium chloride	5.0
Disodium phosphate	2.5
Agar	1.0

### Bromcresol purple dextrose fermentation broth (pH 7.2)

Bacto <sup>®</sup> casitone	10
Dextrose	5
Bromcresol purple (0.2%)	0.01

Bromcresol purple (0.2%) is made separately and filter sterilized. 5 ml is aseptically added to the medium.

*Note: Autoclave at 12 lb pressure for 15 minutes.*

### **Bromcresol purple lactose fermentation broth (pH 7.2)**

Bacto casitone	10
Lactose	5
Bromcresol purple (0.2%)	0.01
Bromcresol purple (0.2%) is made separately and filter sterilized. 5 ml is aseptically added to the above medium.	

*Note: Autoclave at 12 lb pressure for 15 minutes.*

### **Bromcresol purple maltose fermentation broth (pH 7.2)**

Bacto casitone	10
Maltose	5
Bromcresol purple (0.2%)	0.01
Bromcresol purple (0.2%) is made separately and filter sterilized. 5 ml is aseptically added to the medium.	

*Note: Autoclave at 12 lb pressure for 15 minutes.*

### **Bromcresol purple sucrose fermentation broth (pH 7.2)**

Bacto casitone	10
Sucrose	5
Bromcresol purple (0.2%)	0.01
Bromcresol purple (0.2%) is made separately and filter sterilized. 5 ml is added to the medium aseptically.	

*Note: Autoclave at 12 lb pressure for 15 minutes.*

### **Campy-BAP agar (pH 7.0)**

Trypticase peptone	10.0
Thiotone™	10.0
Dextrose	1.0
Yeast extract	2.0
Sodium chloride	5.0
Sodium bisulfide	0.1
Agar	15.0
Vancomycin	10.0 mg
Trimethoprim lactate	5.0 mg
Polymyxin B sulfate	2500.0 IU
Amphotericin B	2.0 mg
Cephalothin	15.0 mg
Defibrinated sheep blood	10.0%

*Note: Aseptically add the antibiotics and defibrinated sheep blood to the sterile, molten, and cooled agar.*

### **Chocolate agar (pH 7.0)**

Proteose peptone	20.0
Dextrose	0.5
Sodium chloride	5.0
Disodium phosphate	5.0
Agar	15.0

*Note: Aseptically add 5.0% defibrinated sheep blood to the sterile and molten agar. Heat at 80°C until a chocolate color develops.*

### **Crystal violet agar (pH 7.0)**

Bacto beef extract	3
Bacto peptone	5
Bacto crystal violet	0.00014
Bacto agar	15

*Note: 1.0 ml of a crystal violet stock solution may be added to the base medium. Stock solution: 14 mg of crystal violet dye dissolved in 100 ml of distilled water.*

### **Decarboxylase broth (Moeller) (pH 6.0)**

Peptone	5.0
Beef extract	5.0
Dextrose	0.5
Bromcresol purple	0.01
Cresol red	0.005
Pyridoxal	0.005
Distilled water	1000.0 ml
To make amino acid-specific medium, add one of the amino acids below; dispense in 3- to 4-ml amounts and autoclave at 121°C for 10 minutes.	
L-lysine dihydrochloride or L-arginine monohydrochloride or L-ornithine dihydrochloride	10 g/l

### **Deoxyribonuclease (DNase) agar (pH 7.3)**

Deoxyribonucleic acid	2.0
Phytane	5.0
Sodium chloride	5.0
Trypticase	15.0
Agar	15.0

### **Endo agar (pH 7.5)**

Peptone	10.0
Lactose	10.0
Dipotassium phosphate	3.5
Sodium sulfite	2.5
Basic fuchsin	0.4
Agar	15.0

### **Eosin-methylene blue agar (Levine) (pH 7.2)**

Peptone	10.0
Lactose	5.0
Dipotassium phosphate	2.0
Agar	13.5
Eosin Y	0.4
Methylene blue	0.065

### **Gel diffusion agar**

Sodium barbital buffer	100.0 ml
Noble agar	0.8

### Glucose acetate yeast sporulation agar (pH 5.5)

Glucose	1
Yeast extract	2
Sodium acetate (with 3H <sub>2</sub> O)	5
Bacto agar	15

### Glucose salts broth (pH 7.2)

Dextrose	5.0
Sodium chloride	5.0
Magnesium sulfate	0.2
Ammonium dihydrogen phosphate	1.0
Dipotassium hydrogen phosphate	1.0

### Glycerol yeast extract agar supplemented with aureomycin (pH 7.0)

Glycerol	5.0 ml
Yeast extract	2.0
Dipotassium phosphate	1.0
Agar	15.0

*Note: Aseptically add aureomycin, 10 µg per ml, to the sterile, molten, and cooled agar.*

### Grape juice broth

Commercial grape or apple juice	
Ammonium biphosphate	0.25%

*Note: Sterilization not required when using a large yeast inoculum.*

### Hay infusion broth

Hay infusion broth preparations are prepared 1 week ahead of the laboratory session in which they will be used. Into a 2000-ml beaker place about 800 ml of water and two to three handfuls of dry grass or hay (obtained from a farm or storage barn). During the incubation period, the infusion should be aerated by passing air through a rubber tube attached to an air supply. This preparation is sufficient for a class and can be dispensed in 50-ml beakers.

### Hektoen enteric agar (pH 7.1)

Peptic digest of animal tissue	12.0
Yeast extract	3.0
Bile salt	9.0
Lactose	12.0
Sucrose	12.0
Salicin	2.0
Sodium chloride	5.0
Sodium thiosulfate	5.0
Ferric ammonium citrate	1.5
Bromthymol blue	0.064
Acid fuchsin	0.5
Agar	13.5

### Inorganic synthetic broth (pH 7.2)

Sodium chloride	5.0
Magnesium sulfate	0.2
Ammonium dihydrogen phosphate	1.0
Dipotassium hydrogen phosphate	1.0

### KF broth (pH 7.2)

Polypeptone	10.0
Yeast extract	10.0
Sodium chloride	5.0
Sodium glycerophosphate	10.0
Sodium carbonate	0.636
Maltose	20.0
Lactose	1.0
Sodium azide	0.4
Phenol red	0.018

### Lactose fermentation broth 1× and 2×\* (pH 6.9)

Beef extract	3.0
Peptone	5.0
Lactose	5.0

\*For 2× broth use twice the concentration of the ingredients.

### Litmus milk (pH 6.8)

Skim milk powder	100.0
Litmus	0.075

*Note: Autoclave at 12 lb pressure for 15 minutes.*

### Luria-Bertani (Miller) agar base (pH 7.0)

Pancreatic digest of casein	10.0
Yeast extract	5.0
Sodium chloride	0.5
Agar	15.0

### Luria-Bertani (Miller) broth (pH 7.0)

Tryptone	10.0
Yeast extract	5.0
Sodium chloride	10.0

### MacConkey agar (pH 7.1)

Bacto peptone	17.0
Proteose peptone	3.0
Lactose	10.0
Bile salts mixture	1.5
Sodium chloride	5.0
Agar	13.5
Neutral red	0.03
Crystal violet	0.001

**Mannitol salt agar (pH 7.4)**

Beef extract	1.0
Peptone	10.0
Sodium chloride	75.0
<i>d</i> -Mannitol	10.0
Agar	15.0
Phenol red	0.025

**m-Endo broth (pH 7.5)**

Yeast extract	6.0
Thiotone peptone	20.0
Lactose	25.0
Dipotassium phosphate	7.0
Sodium sulfite	2.5
Basic fuchsin	1.0

*Note: Heat until boiling; do not autoclave.*

**m-FC broth (pH 7.4)**

Biosate™ peptone	10.0
Polypeptone peptone	5.0
Yeast extract	3.0
Sodium chloride	5.0
Lactose	12.5
Bile salts	1.5
Aniline blue	0.1

*Note: Add 10 ml of rosolic acid (1% in 0.2N sodium hydroxide). Heat to boiling with agitation; do not autoclave.*

**Milk agar (pH 7.2)**

Skim-milk powder	100.0
Peptone	5.0
Agar	15.0

*Note: Autoclave at 12 lb pressure for 15 minutes.*

**Minimal agar (pH 7.0)****Minimal agar, supplemented with streptomycin and thiamine\*****Solution A (pH 7.0)**

Potassium dihydrogen phosphate	3.0
Disodium hydrogen phosphate	6.0
Ammonium chloride	2.0
Sodium chloride	5.0
Distilled water	800.0 ml

**Solution B (pH 7.0)**

Glucose	8.0
Magnesium sulfate • 7H <sub>2</sub> O	0.1
Agar	15.0
Distilled water	200.0 ml

*Note: Autoclave Solutions A and B separately and combine.*

\*To Solution B, add 0.001 g of thiamine prior to autoclaving. To the combined sterile and molten medium, add 50 mg (1 ml of 50 mg per ml) sterile streptomycin solution before pouring agar plates.

**MR-VP broth (pH 6.9)**

Peptone	7.0
Dextrose	5.0
Potassium phosphate	5.0

**Mueller-Hinton agar (pH 7.4)**

Beef, infusion	300.0
Casamino acids	17.5
Starch	1.5
Agar	17.0

**Mueller-Hinton tellurite agar (pH 7.4)**

Casamino acids	20.0
Casein	5.0
L-tryptophan	0.05
Potassium dihydrogen phosphate	0.3
Magnesium sulfate	0.1
Agar	20.0

*Note: Aseptically add 12.5 ml of tellurite serum to the sterile, 50°C molten agar.*

**Nitrate broth (pH 7.2)**

Peptone	5.0
Beef extract	3.0
Potassium nitrate	5.0

**Nitrite broth (pH 7.3)**

Sodium nitrite	2.0
Magnesium sulfate • 7H <sub>2</sub> O	0.5
Ferric sulfate • 7H <sub>2</sub> O	0.03
Sodium chloride	0.3
Sodium carbonate	1.0
Dipotassium hydrogen sulfate	1.0

**Nitrogen-free mannitol agar\* and broth (pH 7.3)**

Mannitol	15.0
Dipotassium hydrogen phosphate	0.5
Magnesium sulfate	0.2
Calcium sulfate	0.1
Sodium chloride	0.2
Calcium carbonate	5.0
*Agar	15.0

**Nutrient agar\* and broth (pH 7.0)**

Peptone	5.0
Beef extract	3.0
*Agar	15.0

**Nutrient gelatin (pH 6.8)**

Peptone	5.0
Beef extract	3.0
Gelatin	120.0

**Peptone broth (pH 7.2)**

Peptone	4.0
---------	-----

**Phenol red dextrose broth (pH 7.3)**

Trypticase	10.0
Dextrose	5.0
Sodium chloride	5.0
Phenol red	0.018

*Note: Autoclave at 12 lb pressure for 15 minutes.*

**Phenol red inulin broth (pH 7.3)**

Trypticase	10.0
Inulin	5.0
Sodium chloride	5.0
Phenol red	0.018

*Note: Autoclave at 12 lb pressure for 15 minutes.*

**Phenol red lactose broth (pH 7.3)**

Trypticase	10.0
Lactose	5.0
Sodium chloride	5.0
Phenol red	0.018

*Note: Autoclave at 12 lb pressure for 15 minutes.*

**Phenol red sucrose broth (pH 7.3)**

Trypticase	10.0
Sucrose	5.0
Sodium chloride	5.0
Phenol red	0.018

*Note: Autoclave at 12 lb pressure for 15 minutes.*

**Phenylalanine agar (pH 7.3)**

Yeast extract	3.0
Dipotassium phosphate	1.0
Sodium chloride	5.0
DL-phenylalanine	2.0
Bacto agar	12.0
Distilled water	1000.0 ml

*Note: Completely dissolve ingredients in boiling water. Dispense in tubes, autoclave, and cool in slanted position.*

**Phenylethyl alcohol agar (pH 7.3)**

Trypticase	15.0
Phytane	5.0
Sodium chloride	5.0
β-Phenylethyl alcohol	2.0
Agar	15.0

**Potato dextrose agar (pH 5.6)**

Infusion from potatoes	200.0
Bacto dextrose	20.0
Bacto agar	15.0

**Sabouraud agar (pH 5.6)****Sabouraud agar supplemented with chlortetracycline (Aureomycin)\***

Peptone	10.0
Dextrose	40.0
Agar	15.0

\*Aseptically add Aureomycin, 10 μg per ml, to the sterile, molten, and cooled medium.

**Salt medium — Halobacterium**

Sodium chloride	250.0
Magnesium sulfate • 7H <sub>2</sub> O	10.0
Potassium chloride	5.0
Calcium chloride • 6H <sub>2</sub> O	0.2
Yeast extract	10.0
Tryptone	2.5
Agar	20.0

*Note: The quantities given are for preparation of 1-liter final volume of the medium. In preparation, make up two solutions, one involving the yeast extract and tryptone and the other the salts. Adjust the pH of the nutrient solution to 7. Sterilize separately. Mix and dispense aseptically.*

**SIM agar (pH 7.3)**

Peptone	30.0
Beef extract	3.0
Ferrous ammonium sulfate	0.2
Sodium thiosulfate	0.025
Agar	3.0

**Simmons citrate agar (pH 6.9)**

Ammonium dihydrogen phosphate	1.0
Dipotassium phosphate	1.0
Sodium chloride	5.0
Sodium citrate	2.0
Magnesium sulfate	0.2
Agar	15.0
Bromthymol blue	0.08

**Snyder test agar (pH 4.8)**

Tryptone	20.0
Dextrose	20.0
Sodium chloride	5.0
Bromcresol green	0.02
Agar	20.0

**Sodium chloride agar, 7.5% (pH 7.0 )**

Bacto beef extract	3.0
Bacto peptone	5.0
Sodium chloride	7.5
Bacto agar	15.0



**Sodium chloride broth, 6.5% (pH 7.0)**

Brain heart infusion broth	100.0 ml
Sodium chloride	6.5

**Starch agar (pH 7.0)**

Peptone	5.0
Beef extract	3.0
Starch (soluble)	2.0
Agar	15.0

**Thioglycollate, fluid (pH 7.1)**

Peptone	15.0
Yeast extract	5.0
Dextrose	5.0
L-cystine	0.75
Thioglycollic acid	0.3 ml
Agar	0.75
Sodium chloride	2.5
Resazurin	0.001

**Tinsdale agar (pH 7.4)**

Proteose peptone, No. 3	20.0
Sodium chloride	5.0
Agar	20.0

*Note: Following boiling, distribute in 100-ml flasks. Autoclave, cool to 55°C, add 15 ml of rehydrated Tinsdale enrichment to each 100 ml, and mix thoroughly before dispensing.*

**Top agar (for Ames test)**

Sodium chloride	5.0
Agar	6.0

**Tributylin agar (pH 7.2)**

Peptone	5.0
Beef extract	3.0
Agar	15.0
Tributylin	10.0

*Note: Dissolve peptone, beef extract, and agar while heating. Cool to 90°C, add the tributyrin, and emulsify in a blender.*

**Triple sugar-iron agar (pH 7.4)**

Beef extract	3.0
Yeast extract	3.0
Peptone	15.0
Proteose peptone	5.0
Lactose	10.0
Saccharose	10.0
Dextrose	1.0
Ferrous sulfate	0.2
Sodium chloride	5.0
Sodium thiosulfate	0.3
Phenol red	0.024
Agar	12.0

**Trypticase nitrate broth (pH 7.2)**

Trypticase	20.0
Disodium phosphate	2.0
Dextrose	1.0
Agar	1.0
Potassium nitrate	1.0

**Trypticase soy agar (pH 7.3)**

Trypticase	15.0
Phytane	5.0
Sodium chloride	5.0
Agar	15.0

**Tryptone agar\* and broth**

Tryptone	10.0
Calcium chloride (reagent)	0.01–0.03 M
Sodium chloride	5.0
*Agar	11.0

**Tryptone soft agar**

Tryptone	10.0
Potassium chloride (reagent)	5.0 ml
Agar	9.0

**Urea broth**

Urea broth concentrate (filter-sterilized solution)	10.0 ml
Sterile distilled water	90.0 ml

*Note: Aseptically add the urea broth concentrate to the sterilized and cooled distilled water. Under aseptic conditions, dispense 3-ml amounts into sterile tubes.*

**Yeast extract broth (pH 7.0)**

Peptone	5.0
Beef extract	3.0
Sodium chloride	5.0
Yeast extract	5.0

## Agarose gel, for Electrophoresis

Agarose	1 g
Tris-borate buffer (1×)	125 ml

*Note: Melt agarose, being careful not to overboil. Cover loosely with foil and hold at room temperature, or place in 60°C waterbath until ready for use.*

## Barritt's reagent, for detection of acetylmethylcarbinol

### Solution A

Alpha-naphthol	5.0 g
Ethanol, absolute	95.0 ml

*Note: Dissolve the alpha-naphthol in the ethanol with constant stirring.*

### Solution B

Potassium hydroxide	40.0 g
Creatine	0.3 g
Distilled water	100.0 ml

*Note: Dissolve the potassium hydroxide in 75 ml of distilled water. The solution will become warm. Allow to cool to room temperature. Add the creatine and stir to dissolve. Add the remaining water. Store in a refrigerator.*

## Biotin-histidine solution, for Ames test

L-Histidine HCl	0.5 mM
Biotin	0.5 mM
Distilled water	10.0 ml

## Buffered glycerol (pH 7.2), for immunofluorescence

Glycerin	90.0 ml
Phosphate buffered saline	10.0 ml

## Diphenylamine reagent, for detection of nitrates

Dissolve 0.7 g diphenylamine in a mixture of 60 ml concentrated sulfuric acid and 28.8 ml of distilled water. Cool and slowly add 11.3 ml of concentrated hydrochloric acid. Allow to stand for 12 hours. Sedimentation indicates that the reagent is saturated.

## Endonuclease buffers

### Buffer 1: EcoRI buffer

Tris-HCl (pH 7.5)	50 mM
MgCl <sub>2</sub>	10 mM
NaCl	100 mM
Triton® X-100	0.02%
BSA	0.1 mg/ml

### Buffer 2: HindIII buffer

Tris-HCl (pH 8.5)	10 mM
MgCl <sub>2</sub>	10 mM
KCl	100 mM
BSA	0.1 mg/ml

### Buffer 3: BamHI buffer

Tris-HCl (pH 8.0)	10 mM
MgCl <sub>2</sub>	5 mM
KCl	100 mM
2-Mercaptoethanol	1 mM
Triton X-100	0.02%
BSA	0.1 mg/ml

## Ferric chloride reagent

Ferric chloride	10.0 g
Distilled water	100.0 ml

## Gram's iodine, for detection of starch

As in Gram's stain

## Hydrogen peroxide, 3%, for detection of catalase activity

*Note: Refrigerate when not in use.*

## Kovac's reagent, for detection of indole

p-Dimethylaminobenzaldehyde	5.0 g
Amyl alcohol	75.0 ml
Hydrochloric acid (concentrated)	25.0 ml

*Note: Dissolve the p-dimethylaminobenzaldehyde in the amyl alcohol. Add the hydrochloric acid.*

**Loading Dye 6×, for gel electrophoresis**

Glycerol (50%)	6 ml
Bromphenol blue (2%)	1 ml
Xylene cyanol (2%)	1 ml
Distilled water	1000 ml

*Note: This can be stored in the refrigerator indefinitely.*

**McFarland Barium Sulfate Standards, for API Staph-Ident procedure**

Prepare 1% aqueous barium chloride and 1% aqueous sulfuric acid solutions. Using the following table, add the amounts of barium chloride and sulfuric acid to clean 15- × 150-mm screw-capped test tubes. Label the tubes 1 through 10.

**Preparation of McFarland Standards**

Tube	Barium Chloride 1% (ml)	Sulfuric Acid 1% (ml)	Corresponding Approximate Density of Bacteria (million/ml)
1	0.1	9.9	300
2	0.2	9.8	600
3	0.3	9.7	900
4	0.4	9.6	1200
5	0.5	9.5	1500
6	0.6	9.4	1800
7	0.7	9.3	2100
8	0.8	9.2	2400
9	0.9	9.1	2700
10	1.0	9.0	3000

**Methyl cellulose, for microscopic observation of protozoa**

Methyl cellulose	10.0 g
Distilled water	90.0 ml

**Methylene blue stain (0.025%)**

Methylene blue 1% stock solution (1 g + 99 ml distilled H <sub>2</sub> O)	10 ml
Distilled water	390 ml

**Methyl red solution, for detection of acid**

Methyl red	0.1 g
Ethyl alcohol	300.0 ml
Distilled water	200.0 ml

*Note: Dissolve the methyl red in the 95% ethyl alcohol. Dilute to 500 ml with distilled water.*

**Nessler's reagent, for detection of ammonia**

Potassium iodide	50.0 g
Distilled water (ammonia-free)	35.0 ml

Add saturated aqueous solution of mercuric chloride until a slight precipitate persists.

Potassium hydroxide (50% aqueous)	400.0 ml
-----------------------------------	----------

*Note: Dilute to 1000 ml with ammonia-free distilled water. Let stand for 1 week, decant supernatant liquid, and store in a tightly capped amber bottle.*

**Nitrate test solution, for detection of nitrites**

<i>Solution A, Sulfanilic acid</i>	
Sulfanilic acid	8.0 g
Acetic acid, 5 N: 1 part glacial acetic acid to 2.5 parts distilled water	1000.0 ml
<i>Solution B, Alpha-naphthylamine</i>	
Alpha-naphthylamine	5.0 g
Acetic acid, 5 N	1000.0 ml

**Oxidase Reagent**

Tetramethyl- <i>p</i> -phenylene diamine dihydrochloride	1.0 g
Distilled water	100 ml

**Orthonitrophenyl-β-D-galactoside (ONPG), for enzyme induction**

0.1 M sodium phosphate buffer (pH 7.0)	50.0 ml
ONPG (8 × 10 <sup>-4</sup> M)	12.5 mg

***p*-Aminodimethylaniline oxalate, for detection of oxidase activity**

<i>p</i> -Aminodimethylaniline oxalate	0.5 g
Distilled water	50.0 ml

*Note: To dissolve fully, gently warm the solution.*

**Phosphate-buffered saline, 1% (pH 7.2–7.4), for immunofluorescence**

<i>Solution A</i>	
Disodium phosphate	1.4 g
Distilled water	100.0 ml
<i>Solution B</i>	
Sodium dihydrogen phosphate	1.4 g
Distilled water	100.0 ml

*Note: Add 84.1 ml of Solution A to 15.9 ml of Solution B. Add 8.5 g of sodium chloride and enough distilled water to make 1 liter.*

**Rabbit plasma, for detection of coagulase activity**

*Note: Store vials at 2°C to 8°C. Reconstitute by the addition of 7.5 ml of sterile water.*

**Sodium barbital buffer, for immunofluorescence**

Sodium barbital	6.98 g
Sodium chloride	6.0 g
1 N hydrochloric acid	27.0 ml
Distilled water, to fill to 1000 ml	

**Toluidine blue solution, 0.1%, for detection of DNase activity**

1% toluidine blue solution	0.1 ml
Distilled water	99.9 ml

**Tris-acetate buffer 10×**

Tris base	48.4 g
Glacial acetic acid	11 g
EDTA (0.5 M)	20 ml
Distilled water	1000 ml

*Note: Add ingredients to 1 liter volumetric flask and then add water to volume.*

**Tris-acetate buffer 1×**

Tris-acetate buffer 10× (see previous entry)	100 ml
Distilled water	900 ml

*Note: Buffer can be stored indefinitely at room temperature.*

**Tris-borate buffer 5×**

Tris base	54 g
Boric acid	27.5 g
EDTA (0.5M, pH 8.0)	20 ml
Distilled water	1000 ml

**Tris-borate buffer 1× working solution**

Tris-borate buffer 5× (see previous entry)	200 ml
Distilled water	800 ml

**Acid-Fast Stain****Carbol fuchsin (Ziehl's)***Solution A*

Basic fuchsin (90% dye content)	0.3 g
Ethyl alcohol (95%)	10.0 ml

*Solution B*

Phenol	5.0 g
Distilled water	95.0 ml

*Note: Mix Solutions A and B. Add 2 drops of Triton X per 100 ml of stain for use in heatless method.*

**Acid Alcohol**

Ethyl alcohol (95%)	97.0 ml
Hydrochloric acid	3.0 ml

**Methylene blue**

Methylene blue	0.3 g
Distilled water	100.0 ml

**Capsule Stain****Crystal violet (1%)**

Crystal violet (85% dye content)	1.0 g
Distilled water	100.0 ml

**Copper sulfate solution (20%)**

Copper sulfate ( $\text{CuSO}_4 \cdot 5\text{H}_2\text{O}$ )	20.0 g
Distilled water	80.0 ml

**Fungal Stains****Lactophenol-cotton-blue solution**

Lactic acid	20.0 ml
Phenol	20.0 g
Glycerol	40.0 ml
Distilled water	20.0 ml
Aniline blue	0.05 g

*Note: Heat gently in hot water (double boiler) to dissolve; then add aniline blue dye.*

**Water-iodine solution**

Gram's iodine (as in Gram's stain)	10.0 ml
Distilled water	30.0 ml

**Gram Stain****Crystal violet (Hucker's)***Solution A*

Crystal violet (90% dye content)	2.0 g
Ethyl alcohol (95%)	20.0 ml

*Solution B*

Ammonium oxalate	0.8 g
Distilled water	80.0 ml

*Note: Mix Solutions A and B.*

**Gram's iodine**

Iodine	1.0 g
Potassium iodide	2.0 g
Distilled water	300.0 ml

**Ethyl alcohol (95%)**

Ethyl alcohol (100%)	95.0 ml
Distilled water	5.0 ml

**Safranin**

Safranin O	0.25 ml
Ethyl alcohol (95%)	10.0 ml
Distilled water	100.0 ml

**Negative Stain****Nigrosin**

Nigrosin, water-soluble	10.0 g
Distilled water	100.0 ml

*Note: Immerse in boiling waterbath for 30 minutes.*

Formalin	0.5 ml
----------	--------

*Note: Filter twice through double filter paper.*

**Spore Stain****Malachite green**

Malachite green	5.0 g
Distilled water	100.0 ml

**Safranin**

Same as in Gram stain

## Cultures

### Bacteria

*Alcaligenes faecalis*  
*Alcaligenes viscolactis*  
*Aquaspirillum itersonii*  
*Bacillus cereus*  
*Bacillus megaterium*  
*Bacillus stearothermophilus*  
*Branhamella catarrhalis* (*Moraxella catarrhalis*)  
*Citrobacter freundii*  
*Citrobacter intermedius*  
*Clostridium butylicum* (*Clostridium beijerinckii*)  
*Clostridium sporogenes*  
*Corynebacterium kutscheri*  
*Corynebacterium xerosis*  
*Enterobacter aerogenes*  
*Enterococcus faecalis*  
*Escherichia coli*  
*Escherichia coli* ATCC e 23724  
*Escherichia coli* ATCC e 23725  
*Escherichia coli* ATCC e 23735  
*Escherichia coli* ATCC e 23740  
*Escherichia coli* ATCC 39991  
*Escherichia coli* ATCC 53100  
*Escherichia coli* B  
*Halobacterium salinarium*  
*Klebsiella ozaenae*  
*Klebsiella pneumoniae* ATCC e 15574  
*Lactobacillus casei*  
*Lactobacillus delbrueckii*  
*Lactobacillus fermenti*  
*Lactococcus lactis*  
*Leuconostoc mesenteroides*  
*Micrococcus luteus*  
*Micrococcus varians*  
*Moraxella bovis*  
*Moraxella catarrhalis* (*Branhamella catarrhalis*)  
*Mycobacterium smegmatis*  
*Neisseria mucosa*  
*Neisseria sicca*  
*Proteus inconstans*  
*Proteus mirabilis*  
*Proteus rettgeri*  
*Proteus vulgaris*  
*Pseudomonas aeruginosa*  
*Pseudomonas fluorescens*  
*Pseudomonas mallei*  
*Pseudomonas savastanoi*  
*Salmonella typhimurium*  
*Salmonella typhimurium* ATCC e 29631  
*Serratia marcescens*  
*Shigella dysenteriae*  
*Staphylococcus aureus* ATCC e 25923  
*Staphylococcus aureus* ATCC e 27659  
*Staphylococcus aureus* ATCC e 27660  
*Staphylococcus aureus* ATCC e 27661  
*Staphylococcus aureus* ATCC e 27691  
*Staphylococcus aureus* ATCC e 27693  
*Staphylococcus aureus* ATCC e 27697  
*Staphylococcus epidermidis*  
*Staphylococcus saprophyticus* ATCC e 15305  
*Streptococcus agalactiae*  
*Streptococcus bovis*  
*Streptococcus mitis*  
*Streptococcus pneumoniae*  
*Streptococcus pyogenes* ATCC 12385  
*Streptococcus* var. Lancefield Group E  
*Streptomyces griseus*  
*Treponema denticola*

### Fungi

*Alternaria* sp.  
*Aspergillus niger*  
*Candida albicans*  
*Cephalosporium* sp.  
*Cladosporium* sp.  
*Fusarium* sp.  
*Mucor mucedo*  
*Penicillium notatum* (*penicillium chrysogenum*)  
*Rhizopus stolonifer*  
*Rhodotorula rubra*  
*Saccharomyces cerevisiae*  
*Saccharomyces cerevisiae* var. *ellipsoideus*  
*Schizosaccharomyces octosporus*  
*Selenotila intestinalis*

## Viruses

T<sub>2</sub> coliphage

## Prepared Slides

### Bacteria

*Aquaspirillum itersonii*

*Bacillus subtilis*

*Spirillum itersonii*

*Staphylococcus aureus*

### Fungi

*Saccharomyces cerevisiae*

## Protozoa

*Balantidium coli*

*Entamoeba histolytica*

*Giardia lamblia (Giardia intestinalis)*

*Plasmodium vivax*

*Trypanosoma gambiense*

## Archaea

*Halobacterium salinarium*

## Other

Blood smear

## Illustration Credits

All illustrations by Precision Graphics unless otherwise noted.

- 35.1: From Harold S. Brown, *Basic Clinical Parasitology*, 4th ed. New York: Appleton-Century-Crofts, 1975.
- 58.2: Carolina Biological Supply Company/Precision Graphics.
- 68.2: Wampole Laboratories Division, Carter-Wallace, Inc., Cranbury, NJ.

## Photo Credits

- Icon.1: Elnur Shutterstock  
Icon.2: Shutterstock
- 0.1: James Cappuccino  
3.2: James Cappuccino  
5.1: Charles D. Winters Photo Reseachers, Inc  
8.1: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education  
9.3a: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education  
9.3b: Jennifer M. Warner, Biology Department, University of North Carolina at Charlotte  
9.3c: Michael Abbey Photo Reseachers, Inc.  
10.1: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education  
11.1a: James Cappuccino  
11.1b: Centers for Disease Control and Prevention  
11.2a-d: David Alexander, University of Portland  
12.2: James Cappuccino  
13.2: Steven R. Spilatro, Department of Biology, Marietta College, Marietta, Ohio  
13.4b: James Cappuccino  
14.2: L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education  
15.1a-b: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education  
15.2a: James Cappuccino  
15.2b: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education  
15.2c: James Cappuccino  
15.3a: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education  
15.3b-c: James Cappuccino  
16.2: David Alexander, University of Portland  
19.3a-e: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education  
20.1a: Hausser Scientific  
20.3: L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education  
20.4: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education  
20.6: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education  
21.4: L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education  
22.1: David Alexander, University of Portland  
22.3: James Cappuccino  
22.5: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education  
23.5a-d: James Cappuccino  
24.2a-e: James Cappuccino  
25.3a-c: James Cappuccino  
25.5a-c: James Cappuccino  
25.8a-c: James Cappuccino  
25.10: James Cappuccino  
26.2a-c: James Cappuccino  
27.2a-b: James Cappuccino  
28.1a-e: James Cappuccino  
29.2a-d: James Cappuccino  
30.1a-b: David Alexander, University of Portland  
30.1c: Brenda Grafton Wellmeyer Department of Biology, Lone Star College North Harris  
31.1: David Alexander, University of Portland  
32.2a-c: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education  
32.4a-c: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education



- 34.1: Eric V. Grave Photo Researchers, Inc.
- 34.2a: Biophoto Associates Photo Researchers, Inc.
- 34.2b: M.I. Walker Photo Researchers, Inc.
- 35.2a-b: Centers for Disease Control and Prevention
- 35.3a: Centers for Disease Control and Prevention
- 35.3b: Biophoto Associates Photo Researchers, Inc.
- 35.4: Eric V. Grave Photo Researchers, Inc.
- 35.5: Marilise Brittes Rott
- 35.6: Centers for Disease Control and Prevention
- 36.2: Jared Martin, Center for Indoor Air Research, Texas Tech University Health Science Center, Lubbock, Texas
- 36.3: Leonard Lessun FBP, Photo Researchers, Inc.
- 36.4: James Cappuccino
- 36.5: Perennou Nuridsany Photo Researchers, Inc.
- 36.6: Biophoto Associates Photo Researchers, Inc.
- 37.1a: James Cappuccino
- 37.1b: Biophoto Associates Photo Researchers, Inc.
- 37.1c: John Durham Science Photo Library, Photo Researchers, Inc.
- 37.2: Nhu Nguyen
- 37.4: L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education
- 39.1: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education
- 41.1a: Courtesy of STERIS Corporation
- 44.2: James Cappuccino
- 44.4: James Cappuccino
- 45.2: James Cappuccino
- 45.3: Ashley Rex, Antimicrobial Test Laboratories
- 46.1: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education
- 49.2: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education
- 49.3a-c: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education
- 50.2: L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education
- 59.4: L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education
- 59.6: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education
- 61.2: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education
- 62.1a-b: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education
- 62.2: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education
- 62.3: James Cappuccino
- 62.4: James Cappuccino
- 62.5a: James Cappuccino
- 62.5b: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education
- 63.1: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education
- 63.2a-b: James Cappuccino
- 63.3: James Cappuccino
- 63.4: James Cappuccino
- 64.1a-b: James Cappuccino
- 64.1c: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education
- 64.2: James Cappuccino
- 64.3: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education
- 64.4: James Cappuccino
- 64.5: James Cappuccino
- 65.1: James Cappuccino
- 66.2: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education
- 66.3a-b: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education
- 67.1: Experiment performed and photographed by L. Brent Selinger, Department of Biological Sciences, University of Lethbridge, Pearson Education
- 69.2: Becton dickinson and Company
- 72.1: LeBeau Custom Medical Stock Photo, Newscom
- 73.1c: Centers for Disease Control

- A**
- Abbé condenser, 37, 38f, 46
- Abscesses, 21, 77, 431
- Absidia* spp., soil populations of, 345
- Absorbance measurements, 100–102
- Abundance of growth (cultural characteristic), 29
- Acetic acid bacteria, 162f
- Acetobacter* spp., 162f
- Acetone, 162f, 344
- Acetylmethylcarbinol, 174, 175f
- Acid curd, 192, 216t
- Acid-alcohol, 81, 82f, 83
- Acid-fast stain, 81–83
- clinical application, 83
  - counterstain, 81–82
  - decolorizing agent, 81
  - of mycobacteria, 83f
  - primary stain, 81
  - procedure, 82f
  - smear preparation, 83
- Acidic stains, 55–56, 57f
- in negative staining, 71f
  - picric acid in, 56f
- Acidophiles, 121, 121f
- Acids/alkylating agents, mechanism of action/uses, 311t
- Acinetobacter* spp., TSI reaction, 168f
- Acne, 431
- Acromonium* spp., 351
- Actinomyces odontolyticus*, dental caries and, 419
- Actinomycetes, soil sample enumeration, 345–348, 347f
- Adaptive enzymes, 365
- Adaptive immunity, 487
- Additive (indifferent) effect of drug combinations, 296, 297, 297f, 299t
- Adenosine triphosphate (ATP), 98, 155
- Adsorption:
  - in animal viruses, 264
  - in bacteriophages, 261
- Aerial mycelium, 239
- Aerobic cellular respiration, 125, 127f, 161f, 164f, 205
- Aerobic microorganisms:
  - identifying via growth distribution, 127f
  - oxygen requirements, 125
  - redox potentials, 131f
- Aerotolerant anaerobes:
  - oxygen requirements, 125
  - distribution of growth, 127f
- Aflatoxin, 237
- African sleeping sickness, 229, 230t, 234f
- Agar, 1, 2f
- blood, 106–107, 107f, 108–109, 126, 423, 424f, 447
  - chocolate, 424, 424f, 425, 425f
  - crystal violet, 105–106, 108–109
  - eosin-methylene blue, 107f
  - Hektoen enteric, 454f
  - MacConkey, 105–106, 107f, 108–109, 111, 126, 454f, 463
  - mannitol salt, 105–106, 107f, 108–109, 108f, 432
  - milk, 137
  - Mueller-Hinton, 294–295
  - phenylethyl alcohol, 105, 106f
  - properties of, 1
  - SIM, 172, 183
  - sodium chloride (7.5%), 105
  - as solid media, 2f
  - starch, 155, 155f
  - tributylin, 156, 156f
- Agar deep tubes, 2, 2f
- Agar plates, 2, 2f
- enrichment culture technique and, 359
  - growth patterns on nutrient, 22, 29–31
  - inoculation of nutrient, 9t
  - Sabouraud, 239, 240, 240f, 242, 248, 425f
- Agar plating method for cell counts, 138–139, 139f
- Agar slants, 2, 2f
- cultural characteristics, 29, 30f
  - cultural characteristics of selected bacteria, 216f
  - growth patterns on nutrient, 29, 30f
- Agarose gel, 400–401, 402f
- Agarose gel electrophoresis, 410, 411, 411f
- Agglutination, 487
- Agglutination reactions, 487
- febrile antibody test, 493–495, 493f
  - fluorescent antibody technique and, 499–500, 500f
- Agglutinins, 487
- Agriculture, biotechnology applications, 387–388
- AIDS (acquired immunodeficiency syndrome), 507, 511
- Alcaligenes faecalis*, 205
- biochemical characteristics, 216f
  - carbohydrate fermentation, 163–164
  - cultivation of, 100
  - cultural characteristics (agar slant), 216f
  - IMViC test series and, 171f
  - litmus milk reaction, 191–194, 216f
  - nitrate reduction test, 197–198, 198f
  - pH requirements, 122
  - TSI agar test, 167–168, 168f
- Alcaligenes* spp.:
  - as food-borne organism, 317
  - TSI reaction, 168f
- Alcaligenes viscolactis*, capsule stain, 91, 92f
- Alcohol production, biochemical pathway, 323f
- Alcohols, mechanism of action/uses, 310t
- Algae, pond water, 45f
- Alkaline reaction in litmus milk reaction, 192, 193f
- Alkalophiles, 121, 121f
- Alpha hemolysis, 107, 107f, 423, 424f
- identifying streptococcal pathogens by, 439, 440f
  - identifying *Streptococcus pneumoniae*, 447, 447f
  - Streptococcus* exhibiting, 107f
- Alternaria* spp., 255, 256t
- soil populations, 345
- Amebic dysentery, 229, 230t, 231t, 233f, 234f
- Ames test for identifying chemical carcinogens, 381–383, 383f
- Amine, 6, 209
- Amino acids, 156, 157
- decarboxylase test, 209–211
  - phenylalanine deaminase test, 209–212, 211f, 212f
  - as zwitterions, 122
- Ammonia, 167
- phenylalanine deamination, 211, 212f
- Ammonification by soil microbes, 343, 344f
- Ammonium hydroxide, 163
- Amoeba*, 221, 224t
- Amoeba proteus*, 221
- Ampicillin, 299t
- Amycolatopsis* spp., 351
- Amylase, 155
- Anaerobes:
  - aerotolerant, 125
  - cultivation of, 131–134
  - facultative, 125, 127f, 131f
  - obligate, 125
  - oxygen's toxicity to, 125
- Anaerobic cellular respiration, 98, 131, 161f
- Analytical profile index (API) system, 455, 456f, 457, 458
- Animal feed/hides as source of food contamination, 317
- Animal viruses, 263–264
- Anionic agents, mechanism of action/uses, 311t
- Anopheles* mosquito, 229, 232f
- Antibiotics, 293
- as antimicrobial agents, 293–299, 293t
  - drug combination synergism, 296–299
  - effective, selection of, 295
  - Kirby-Bauer Antibiotic Sensitivity Test, 294–295, 294–296, 297
  - microbe-produced, 344, 351–353, 352f
  - penicillin activity experiment, 303–305, 304f
  - prototypic, 293t
  - resistance to, 377–378
- Antibodies, 487
- ELISA (enzyme-linked immunosorbent assay), 503–504
- Antibody titer test, 493, 494, 495f
- Antigen-antibody reactions:
  - agglutination and, 487
  - precipitin formation, 488
- Antigens, 487
- ELISA (enzyme-linked immunosorbent assay), 503–504
  - febrile, 493–495, 493f
- Antimetabolites, 294
- Antimicrobial agents, 293–313
- antibiotics, 293–299, 293t
  - antiseptics, 277, 309–313
  - chemotherapeutics, 277, 293–302
  - disinfectants, 277, 309–313

- Antiseptics, 277, 309–313  
 effectiveness of, 309, 312–313  
 mechanism of action/uses, 309–311t  
 susceptibility test, 313f
- Antiserum, 489
- API (analytical profile index) system, 455, 456f, 457, 458
- Apples, fermentation of, 323–324, 323f
- Aquaspirillum itersonii*, 40, 42
- Arborescent (cultural characteristic), 29, 30f
- Arizona* spp., TSI reaction, 168f
- Aseptic inoculation and transfer, 15
- Ascomycetes, 237, 238t
- Ascospores, 247, 248f
- Ascus, 239, 240f, 247, 248f
- Aseptic techniques:  
 defined, 3  
 in subculturing procedures, 14f
- Aspergillus niger*:  
 conidiophore and conidia of, 240f  
 cultivation of, 126, 240f, 242, 242f  
 moist heat to control, 280–281, 281f  
 oxygen requirements, 126  
 radiation to control, 289–290
- Aspergillus* spp.:  
 characteristics of, 238t  
 industrial uses of, 344  
 as soil microbes, 344, 345
- Atmospheric oxygen, *See* Oxygen:
- ATP (adenosine triphosphate), 98, 155
- Autoclaves, 2f, 280, 280f  
 performance testing, 280
- Autotrophs, carbon and, 97
- Auxochrome, 55, 55f, 56f
- Azotobacter* spp., 343
- B**
- Bacillary dysentery, 106, 327, 453
- Bacilli:  
 coliform, *See* Coliform bacteria  
 enteric, *See* Enteric bacteria  
 intestinal, 168f
- Bacilli (rod-shaped) bacteria, 65f, 66, 67f
- Bacillus amyloliquefaciens* H (BamHI), 408, 409t
- Bacillus anthracis*, 88, 91
- Bacillus cereus*, 31, 65, 89  
 biochemical characteristics, 154f, 216f  
 cultivation of, 133  
 cultural characteristics (agar slant), 216f  
 exoenzymatic activities, 157–158  
 IMViC test series, 171f  
 Kirby-Bauer Antibiotic Sensitivity Test, 294–296  
 moist heat to control, 279–281, 281f  
 radiation to control, 289–290
- Bacillus* spp., 75, 87  
 as food-borne organism, 317  
 industrial uses of, 344  
 as soil microbes, 344, 351  
 as vegetative cells or spores, 87
- Bacillus stearothermophilus*, 116, 117
- Bacillus subtilis*, 40, 54
- Bacitracin, 293t  
*Bacillus* spp. and, 344
- Bacitracin sensitivity test, 441, 441f, 442
- Bacteremia, 473
- Bacteria:  
 biochemical characteristics, 216f  
 cell wall of, 277  
 cultivation needs, 99–104  
 cultural characteristics (agar slant), 216f  
 cultural characteristics of, 30f  
 differential staining, 87–91  
 drug susceptibility test, 294–296, 296f  
 genus identification of unknown culture, 215–217, 216f  
 Gram-negative identification, 479, 482–483f  
 Gram-positive identification, 479, 480–481f  
 identification of unknown genus, 215–217, 216t  
 identification of unknown species, 479–483  
 media for routine cultivation of, 99–104  
 morphology, 67f  
 normal flora in human body, *See* Normal flora:  
 oxygen requirements, 125–129  
 pH requirements, 121  
 resistance to moist heat, 279  
 sexually transmitted chlamydial diseases caused by, 510, 511–512  
 shapes/arrangements of, 65f  
 soil sample enumeration, 345–348, 347f  
 specialized media for cultivation of, 105–109  
 spore-forming, 88, 88f  
 temperature requirements, 115–120  
 transduction in, 364, 371, 391  
 vegetative cells, 87  
 viruses infecting, 261–264, *See also* Bacteriophages (phages)
- Bacteria daughter cells, transferring a colony of, 22
- Bacterial cell counts:  
 chemical methods for, 137, 138  
 direct microscopic counts, 137–138  
 electronic cell counters, 138  
 in food products, 319–320  
 serial dilution-agar plate analysis, 138–139, 138f–141f  
 soil sample, 345–348, 347f  
 spectrophotometric analysis, 138
- Bacterial conjugation, 364, 371–373, 391
- Bacterial genetics, 363–383  
 chemical carcinogenicity screening, 381–383  
 enzyme induction, 365–367, 366f  
 isolation of mutants, 377–378  
 recombination via conjugation, 371–373
- Bacterial growth curve, 145–148  
 generation time, 145, 146, 147f  
 stages of typical, 145–146, 145f
- Bacterial morphology, 65–66, 67f
- Bacterial pneumonia, encapsulated, 91
- Bacterial population growth studies, 145
- Bacterial smears, 59–95, 60f  
 from a broth medium, 60–61  
 preparation of, 59–95, 61f  
 proper slide preparation, 60  
 from a solid medium, 60–61
- Bacterial staining, 55–95  
 acid-fast, 81–83  
 differential, 57f, 75  
 gram, 74–80  
 negative, 71–73  
 simple, 57f  
 smear preparation prior to, 60–61, 61f
- Bacteriophages (phages), 261–273, 262f  
 cultivation/enumeration of, 265–267, 265f, 267f  
 infectious process of, 261–262, 263f  
 isolation of coliphages, 271–273, 272f  
 lambda, 262, 407–411, 409t  
 lytic/lysogenic cycles of, 262–263, 263f  
 structural elements, 261, 262f  
 virulent, 262–263, 263f
- Bacteriuria, determination of, 467–469
- Bacteroides* spp., 423  
 as food-borne organism, 317
- Bacturcult culture tube, 469f  
 interpretation of colony counts, 469t
- Baker's yeast, 247
- Balantidium coli*, 221, 229, 231, 231t, 234f, 327
- Barritt's reagent, 174, 175f, 177f
- Basic dyes, mechanism of action/uses, 311t
- Basic stains, 55–56, 57f, 311t  
 methylene blue, 56f  
 use for bacterial staining, 57
- Basidiomycetes, 237, 238t
- Basidium, 239f
- Bausch & Lomb Spectronic 20 spectrophotometer, 100, 101, 101f
- BBL Septi-Chek System, 475, 475f
- Beaded (cultural characteristic), 29
- Beef extract, 99
- Beggiatoa* spp., 344
- Benzalkonium chloride, mechanism of action/uses, 311t
- Benzene, 55, 55f, 56f
- Benzylpenicillin (penicillin G), 303
- Bergey's Manual of Systematic Bacteriology*, 29, 215, 479
- Beta hemolysis, 107, 107f, 423, 424f  
 identifying streptococcal pathogens by, 439, 440f
- Beta-galactosidase, 191f
- Beta-galactosidase induction, 365, 366f
- Beta-lactam ring, 303, 303f
- Beta-propiolactone:  
 as aseptic technique, 3f  
 mechanism of actions/uses, 311t
- Beverages, wine production, 247, 323–324, 344
- Bile, 447
- Bile esculin test, 441
- Bile salts, 447
- Bile solubility test, 447
- Binary fission, 45, 88f, 145, 510
- Biochemical activities of microorganisms, 184f  
 amino acid utilization, 209–212  
 carbohydrate fermentation, 161–164  
 catalase test, 201–204  
 exoenzymatic, 155–159  
 in genus identification of unknown culture, 215–217, 216f  
 hydrogen sulfide test, 154f, 183–184  
 IMViC test series, 154f, 171–181  
 litmus milk reaction, 191–196  
 nitrate reduction test, 197–200  
 oxidase test, 205–206  
 triple sugar-iron agar test, 167–168  
 urease test, 187–190
- Biochemical characteristics of bacteria, 216f
- Biooxidative pathways, 98, 125, 161f
- Biotechnology, 387–416  
 beneficial applications of, 387, 389  
 experiments:  
 bacterial transformation, 389–396, 392f  
 gel electrophoresis, 388, 397, 398, 401, 408  
 isolation of bacterial plasmids, 397–406  
 restriction analysis, 407–416  
 plasmids and, 390–391, 397–398, 401, 407
- Blastomyces*, 238t
- Blood agar, 106–107, 107f, 108–109, 126  
 detecting normal skin flora on, 424  
 detecting *Streptococcus pneumoniae* on, 447, 447f  
 detecting throat flora on, 423, 424f

- types of hemolysis exhibited on, 107*f*, 424*f*, 441*f*
- Blood sample analysis, 473–476
- Blood smear, microscopic examination of, 42, 231
- Blow-out pipette, 5*f*
- Body tube of microscope, 37–38, 38*f*
- Boiling heat, to control microbial growth, 278*f*
- Boils, 431
- Botrytis* spp., as food-borne organism, as food-borne organism, 317
- Breed smears, 138
- Brewer jar, 132*f*
- Brightfield microscope:
  - compound, 36
  - features of, 35
  - limitation of, 35
- Bronchopneumonia, 439, 440*f*
- Broth culture, 13
  - bacterial smears and, 60, 61*f*
  - nutrient, cultural characteristics, 21, 31–32
- Broth media, 1, 2*f*
  - defined, 1
  - glucose salts, 99, 101
  - growth in, 30*f*
  - inorganic synthetic, 99, 101
  - nutrient, 99
- Broth-to-slant transfer, 13, 14*f*
- Brownian movement, 46
- Brucella* spp., 473
- Budding, in yeasts, 247, 248*f*
- Buffering system, 121–122
- Buffers, 121
- Bunsen burners, 2*f*, 3, 15, 21, 22
- Butano, 344
- Butter consistency (cultural characteristic), 29
- ## C
- Cadaverine, 210, 210*f*
- Calcium caseinate, 192
- CAMP test, 441, 441*f*, 443
- Campylobacter jejuni*, 327, 473
  - diseases associated with, 463
  - isolation/presumptive identification, 463–464, 464*f*
- CampyPak jar, 463, 464*f*
- Candida*, 238*t*, 256, 257*t*
- Candida albicans*, 238*t*, 247*f*, 248, 423, 467
- Candidiasis, 507
- Capsular swelling, 448
- Capsule staining, 89–95, 93*f*
  - procedure, 92*f*
- Capsules:
  - defined, 91
  - of *Streptococcus pneumoniae*, 447
- Carbenicillin, 299*t*
- Carbohydrate fermentation, 154*f*, 161–164
  - aerobic respiration pathway vs., 163, 164*f*
  - biooxidative pathways and, 161*f*
  - dental caries and, 419
  - Embden-Meyerhof (glycolytic) pathway, 161, 162*f*
  - gas production in, 161, 163*f*
  - pyruvic acid and, 161, 162*f*
  - by *Saccharomyces cerevisiae*, 116
  - temperature and, 116
  - tests for, 163, 163*f*
  - triple sugar-iron agar test and, 167–168, 168*f*
- Carbol fuchsin, 66, 81, 82*f*, 83
- Carbon:
  - autotrophs and, 97
  - heterotrophs and, 97
- Carbon cycle, 344
- Carbon dioxide:
  - anaerobics and increased growth, 138
  - autotrophs and, 91
  - carbohydrate fermentation and, 163, 163*f*
  - decarboxylation and production of, 209, 210, 210*f*
- Carbon dioxide gas, of lactose fermentation, 191, 193*f*
- Carbuncles, 431
- Carcinogenicity, 381
- Carcinogens, Ames test for identifying, 381–383
- Cardinal (significant) temperature points, 115, 115*f*
- Casein, 156, 192
- Casein hydrolysis, 154*f*, 156, 158
- Catalase, 201
- Catalase test, 154*f*, 201–204
- Cationic agents, mechanism of action/uses, 311*t*
- Cefoxitin, 299*t*
- Cellmatics™ HSV Detection System, 509
- Cell-mediated immunity, 487
- Cell-membrane damage from antimicrobial agents, 277–278
- Cellular respiration, 98
  - aerobic, 161*f*, 162, 164*f*, 205
  - anaerobic, 98, 131, 161*f*
  - defined, 98, 131
  - fermentation vs., 98, 161*f*, 162–164
  - oxygen requirements for, 125–127, 127*f*
- Cellulitis, 439
- Cell-wall injury from antimicrobial agents, 277
- Cephalosporium* spp., 255
  - morphology/microscopic appearance, 257*t*
- Cephalothin, 299*t*
- Cercomonas, 221, 224*f*
- Cervicitis, NGU and, 510
- Chancre, 507
- Chemically defined media, 99
- Chemicals, Ames test for identifying carcinogenic, 381–383, 383*f*
- Chemotherapeutic agents, 277, 293–302
  - antibiotics, 293, 293*t*
  - antiseptics, 277, 293–302
  - determining penicillin activity/inhibition, 303–305
  - disinfectants, 277, 293–302
  - Kirby-Bauer Antibiotic Sensitivity Test, 294–298, 295*f*, 296*f*
  - zone of inhibition, 294–295, 296*f*, 299*t*
  - susceptibility test standards for selected, 299*t*
  - synergistic drug combination test, 296–299, 299*t*
  - synthetic drugs, 293, 294
- Chemotrophs, 98
- Chlamydia, 510
  - treating, 511
- Chlamydia trachomatis*, 510–511
- Chlamydiae, defined, 510
- Chlamydial diseases, detection of, 510–511
- Chlamydomonas, 45*f*
- Chloramine, mechanism of action/uses, 310*t*
- Chloramphenicol, 293*t*
- Chlorine compounds, mechanism of action/uses, 310*t*
- Chocolate agar:
  - detecting skin flora on, 425, 425*f*
  - detecting throat flora on, 423, 424*f*
- Cholera, 327
- Chromium-sulfuric acid method, 132*f*
- Chromogen, 55, 55*f*, 56*f*, 57*f*
- Chromogenic microorganisms, pigmentation of, 29
- Chromophore, 55*f*
- Ciliophora, 221, 224*f*
  - parasitic, 231*t*
- Circular (cultural characteristic), 29, 30*f*
- Circular DNA, 397–398, 398*f*
- Citrase, 177, 177*f*
- Citrate:
  - as carbon source for cellular energy, 175
  - enzymatic degradation, 175–176, 176*f*
- Citrate permease, 175
- Citrate utilization test, 154*f*, 171, 175–176, 176*f*, 177*f*, 216*f*
- Citric acid, 344
- Citrobacter freundii*:
  - decarboxylase test, 209–211, 211*f*
  - TSI reaction, 168*f*
- Citrobacter* spp.:
  - amino acid utilization, 209
  - as food-borne organism, 317
- Cladosporium* spp., 255, 257*t*
  - morphology/microscopic appearance, 257*t*
  - soil populations of, 345
- Clindamycin, 299*t*
- Cloning vectors, plasmids as, 397
- Clonorchis sinensis*, 327
- Closed circular DNA, 397
- Clostridia* spp., 88
- Clostridia, use of pyruvic acid, 162*f*
- Clostridium*, 87
- Clostridium acetobutylicum*, 162*f*, 344
- Clostridium beijerinckii*, 162*f*
- Clostridium sporogenes*, 89
  - cultivation of, 126, 133
  - oxygen requirements, 126
- Clostridium* spp., 344, 423
  - as food-borne organism, 317
  - as vegetative cells or spores, 87
- Coagulase test, to differentiate staphylococcal species, 432, 432*t*, 433*f*
- Coagulase, 431
- Coarse-adjustment knobs, of microscope, 37, 38*f*
- Cocci (spherical-shaped) bacteria, 65*f*, 66, 67*f*
- Coccidioides immitis*, 52
- Coccidioides* spp., 344
- Coenzymes, 98
- Cold temperatures, to control microbial growth, 278*f*
- Coliform bacteria:
  - differential/selective media and, 106
  - in food products, 319–320
  - water sample tests, 329–334, 330*f*, 334*f*, 337–340
- Coliphages, isolation of, 271–273, 272*f*
- Colonial pigmentation test, to differentiate staphylococcal species, 432*f*
- Colony growth characteristics, examining to aid identification, 31
- Competent cell, 389
- Competitive inhibition, 294
- Completed test (water analysis), 330*f*, 331, 333, 336
- Complex media, 99
- Computer-assisted systems for identifying Enterobacteriaceae, 454, 455, 457*f*
- Condenser adjustment knob, microscopes, 37
- Condenser knob (microscope), 38*f*, 39, 40, 87
- Condenser (microscope), 35–39, 38*f*, 39*f*
- Confirmed test (water analysis), 329, 330*f*, 332–333, 335
- Conidia, 239*f*, 240*f*

- Conidiophore, 239f, 240f
- Conjugation, bacterial, 364, 371–373, 391
- Consistency (cultural characteristic), 29
- Constitutive enzymes, 365
- Convex elevation (cultural characteristic), 30f, 31
- Copper sulfate (20%), 89, 91, 92f
- Corynebacterium xerosis*, 425
- biochemical activities and, 154f
- cultural and biochemical characteristics, 216f
- cultural/biochemical characteristics, 216f
- oxygen requirements, 126
- Coulter Counter, 138
- Counterstain, 76–77, 82f
- for acid-fast stain, 81–82
- for gram stain, 75
- for spore stain, 87–88
- Crateriform (cultural characteristic), 30f, 31
- Cresols, mechanism of action/uses, 309t
- Crossing over (exchange of genetic material), 371
- Crowded-plate technique, 351, 352f
- Cryptococcus* spp., 344
- Crystal violet agar, 105–106, 108–109
- Crystal violet (Hucker's) stain, 66, 76, 76f, 105–106, 108–109
- Crystal violet, mechanism of action/uses, 311t
- C-substance hapten, 439
- Cultivation chambers, 2f, 4–6
- Cultivation of microorganisms:
- of anaerobics, 131–134
  - bacteriophages, 265–267
  - counting viable cells, 137–142
  - enumeration of viable cells, 137–142
  - fungi (molds), 239–242
  - media used in, 99–104
  - nutritional requirements, 97–98
  - oxygen requirements, 125–127, 127f
  - oxygen requirements, 98
  - pH requirements, 98, 121–122f
  - temperature requirements, 98, 115–117
- Cultural characteristics of microorganisms, 29–32, 30f
- on nutrient agar plates, 29–31, 30f
  - on nutrient agar slants, 29, 30f, 216f
  - on nutrient broth cultures, 21, 31–32
  - on nutrient gelatin, 31
- Culture media, 1
- Culture tubes, 2f, 3
- Culture vessels, 4f
- Cultures, 65
- medium, 1, *See also* Media
  - mixed, 19
  - pure, 1, 2f, 19–22, 24f, 359
  - stock, 22
  - transfer instruments, 2f, 3–4, 5f
  - transfer techniques, 13–15
  - vessels, 3, 4f
- Curd formation, 192, 193f
- Curvularia* spp., 344
- Cysteine, 183, 184f
- Cysteine desulfurase, 183, 184f
- Cystitis, 431, 467
- Cysts of parasitic protozoa, 230–231t, 234f
- Cytochrome oxidase, 205
- Cytoplasmic damage from antimicrobial agents, 278f
- D**
- Dakin's fluid, mechanism of action/uses, 310t
- Darkfield microscope, 35
- Decarboxylase test, 209–212, 210f, 211f
- Decarboxylation and production, 209, 210f
- Decline (death) phase, in typical bacterial growth curve, 145, 145f
- Decolorizing agent, 75, 76f
- for acid-fast stain, 57f, 81
  - capsule stain, 57f, 89, 90f
  - gram stain, 57f, 75
  - for spore stain, 87
- DeLong shaker flask, 4f
- Denitrification by soil microbes, 343, 344f
- Dental caries, 419
- determining susceptibility to, 419–420, 420f
- Deoxyribonuclease (DNase) test, 432, 432f, 433f
- Desiccation of microbes, 278f
- Desulfotomaculum*, 87
- Detergents/soaps, mechanism of action/uses, 311t
- Deuteromycetes (Fungi Imperfecti), 237, 238t
- soil populations of, 345
- Dextran, dental caries and, 419, 420f
- Dextranucrase, 419
- Dextrins, 155
- Dextrose fermentation, by various bacteria, 216f
- Diamine, 209
- Diaphragm (microscope), 37, 38f, 39–40, 41f
- Diarrhea, protozoa causing, 234f
- Differential staining, 57f, 75
- acid-fast stain, 81–83
  - capsule staining, 89–95, 92f, 93f
  - gram stain, 75–80
  - spore stain, 87–89, 93f
  - visualizing cell structures by, 87–95
- Differential/selective media, 105–106, 107f
- Differentiation of bacteria, and staining techniques, 57
- Dilution:
- defined, 140f
  - serial dilution-agar plate procedure, 138–139, 138f–141f
- Dilution bottle method, 147f
- Dilution factor, defined, 140f
- Dilution-plating procedures for use in bacterial growth curve, 148f
- Dipeptides, 156
- Diphtheroids, *See* *Corynebacterium xerosis*
- Diplobacillus, 65f
- Direct microscopic cell counts, 137–138, 137f
- Direct Specimen Reagent (MicroTrak®), 511
- Directigen Flu A test, 503–505
- Directigen test, 441, 443
- Disinfectants, 277, 309–313
- effectiveness of, 309, 312–313
  - mechanism of action/uses, 309–311t
- DNA (deoxyribonucleic acid), 98, *See also* Plasmids
- bacterial transformation and, 389, 390f
  - circular, 397–398
  - closed and nicked, 397–398, 398f
  - damage:
    - by antimicrobial agents, 278
    - by radiation, 289–290
  - microbial enzyme systems for repairing
    - damaged, 290
    - mutations, 363–364, 377–378
    - recombinant DNA technology, 387–416
    - research, usefulness of bacteria to, 363
    - restriction endonucleases, 407–412, 407f
- DNA viruses, 261
- DNase, 431
- DNase test, to differentiate staphylococcal species, 432, 432f, 433f
- Double-stranded DNA herpes simplex virus (HSV), 509
- Dracunculus medinensis* (guinea worm), 327
- Drug combinations, synergistic effects, 296–298, 297f
- Dry consistency (cultural characteristic), 29
- Dry heat:
- aseptic technique and, 3f
  - to control microbial growth, 278f
  - sterilization time, 279
- Dyes, *See* Acidic stains; Stains:
- Dysentery:
- amebic, 229, 230t, 231t, 233f, 234f
  - bacillary, 106, 327, 453
- E**
- Echinulate (cultural characteristic), 29, 30f
- EcoRI* endonuclease, 407, 408, 409t
- palindrome for, 407f
- Effuse (cultural characteristic), 29, 30f
- Electromagnetic radiation to control microbial growth, 278f, 289–290
- Electromagnetic spectrum, 289–290, 289f
- Electron microscopes, 35, 36
- Electronic cell counters, 138
- Elephantiasis, 511
- Elevation (cultural characteristic), 30f, 31
- ELISA (enzyme-linked immunosorbent assay), 503–505
- Embden-Meyerhof (glycolytic) pathway, 161, 162f
- Emulsifiers, mechanism of action/uses, 311t
- Encapsulated bacterial pneumonia, 91
- ENCISE II computer-assisted system, 444, 445f, 457–458
- Endocarditis, 431, 439
- Endoenzymes (intracellular enzymes), 153–154, 154f
- Endonucleases, 407–408, 407f, 408f
- Endospores, 87, 88, 88f
- Energy sources (metabolic) of microbes, 98
- Enriched media, 106–107
- Enrichment culture technique, isolation of *Pseudomonas* by, 357–362
- Entamoeba histolytica*, 221, 229, 230t, 231, 233f, 327
- Enteric bacteria:
- cultivation of, 106
  - differential/selective media, 106
  - diseases associated with, 453
  - IMViC test series, 154f
  - pathogenic, 171
    - multitest systems for isolating
      - identifying, 453–458
  - pH range and, 121
  - as source of food contamination, 317
  - triple sugar-iron agar test, 167–168
  - uses of pyruvic acid by, 162f
  - waterborne disease from, 327
- Enteritis, 463
- Enterobacter aerogenes*, 91, 92f, 456
- biochemical activities and, 154f
  - cultivation of, 106, 108
  - cultural/biochemical characteristics, 216f
  - glucose fermentation, 174, 175f
  - hydrogen sulfide test, 154f, 183–184, 184f
  - IMViC test series, 171–177
  - methyl red test, 173–174, 177f
  - use of pyruvic acid by, 162f
  - Voges-Proskauer test, 174–175, 177f
- Enterobacter* spp., 209
- cultural characteristics, 216f

- as food-borne organism, 317
  - as normal intestinal flora, 453
  - osmotic pressure to control, 285–286
  - TSI reactions, 168*f*
  - Enterobacteriaceae, 171–181, 417
    - amino acid utilization, 209
    - computer-assisted systems for identifying, 454, 455, 456*f*
    - IMViC test series, 171–182
    - oxidation test, 205
    - pathogenic, 453–458
    - TSI agar test, 167–168, 168*f*
  - Enterococcus faecalis*, 470, 499
    - cultivation of, 108, 126
    - Kirby-Bauer Antibiotic Sensitivity Test, 294–296, 294–298, 297*f*
    - oxygen requirements, 126
  - Enterotoxin, 431
  - Enterotube II Multitest System, 454, 454*f*, 457–458
  - Entire margin (cultural characteristic), 29
  - Enumeration of microorganisms:
    - bacterial growth curve, 145–148
    - bacteriophages, 265–267
    - counting methods, 137–148
    - in food products, 319–320
    - of soil microbes, 345–348, 347*f*
  - Enzyme induction experiment, 365–367, 366*f*
  - Enzyme-linked immunosorbent assay (ELISA), 503–505
  - Enzymes, 115
    - adaptive, 365
    - in anaerobics, 131
    - constitutive, 365
    - defined, 154. *See also* Biochemical activities
    - extracellular, 153–154, 153–158, 154*f*
    - hydrolytic, 153–158, 154*f*
    - inactivation by antimicrobial agents, 278
    - inducible, 365
    - intracellular, 153–154, 154*f*
    - microbe-produced, for industry, 344
    - repair of radiation-damaged DNA by, 290
  - Eosin-methylene blue agar (Levine), 106, 107*f*
  - Epididymitis, NGU and, 510
  - Episome, *See* Plasmids
  - Eremothecium ashbyii*, vitamins and, 344
  - Erysipelas, 439
  - Erythrocytic stage, 230, 232*f*
  - Erythromycin, 299*t*
  - Escherichia coli*, 21, 22, 25, 26, 31, 65, 193*f*
    - antimicrobial drug combinations, 297–298
    - biochemical activities, 154
      - endoenzymatic, 161–212
      - exoenzymatic, 157–158
    - biochemical characteristics, 216*f*
    - blood specimen analysis, 473–476
    - cultivation of, 100, 108
    - decarboxylase test, 209–211, 211*f*
    - enteric pathogen identification, 456
    - enumeration of viable cells, 137–142
    - experiments:
      - antibiotic-producing microbes and, 351–353
      - bacterial growth curve, 145–148
      - bacteriophage cultivation/ enumeration, 265–268
      - blood specimen analysis, 473–476
      - carbohydrate fermentation, 161–164, 163*f*
      - conjugation in strains of, 371–373
      - cultivation of, 100, 105, 106, 116, 122, 133, 139
      - cultural characteristics, 216*f*
      - enzyme induction, 365–367, 366*f*
      - febrile antibody test, 493–495, 493*f*
      - identification of unknown genus, 215–217, 216*f*
      - isolating plasmids from, 397–401
      - isolating streptomycin-resistant mutant of, 377–380
      - Kirby-Bauer Antibiotic Sensitivity Test, 294–298, 295
    - Escherichia coli*:
      - experiments:
        - litmus milk reactions, 191–194
    - Escherichia coli*:
      - experiments:
        - methyl red test, 173–174, 177*f*
        - nitrate reduction test, 197–198, 198*f*
        - oxidase test, 205–206, 206*f*
    - Escherichia coli*:
      - experiments:
        - pH requirements, 122
        - phenylalanine deaminase test, 211–212, 212*f*
        - temperature requirements, 116
    - Escherichia coli*:
      - experiments:
        - urease test, 187–188, 188*f*
    - febrile antibody test, 493–495, 493*f*
    - as food-borne organism, 317
    - genetic map of, 371*f*
    - Gram-negative stain of, 75*f*
    - identification of unknown genus, 215–217, 216*f*
    - IMViC test series, 154*f*, 171–177
    - as indicator of fecal pollution, 327
    - methyl red test, 173–174, 177*f*
    - as normal intestinal flora, 453
    - soil sample enumeration, 345–348, 347*f*
    - spontaneous mutation rate in, 363
    - TSI agar test, 167–168, 168*f*
    - urine specimen analysis, 470
    - Voges-Proskauer test, 174–175, 177*f*
  - Escherichia coli* B, coliphage isolation, 271–273
  - Escherichia coli*, strains of, 365
  - Escherichia* spp., as food-borne organism, 317
  - Essential metabolite, 294
  - Ester bonds, 155
  - Ethyl alcohol (95%), 76, 76*f*
    - mechanism of action/uses, 310*t*
  - Ethylene oxide, mechanism of action/uses, 311*t*
  - Ethylhydrocupreine hydrochloride, 447
  - Eubacteriales, soil populations of, 345
  - Euglena*, 221, 224, 224*f*, 224*t*
  - Excision repair system, 290
  - Exoenzymes, *See* Extracellular enzymes
  - Extracellular enzymes (exoenzymes), 153–158, 154*f*
    - metabolic activities, 155–158
  - Eye conjunctiva, normal flora of, 423
  - Eyepiece lens of microscope, 37
- ## F
- Facultative anaerobes:
    - carbohydrate fermentation by, 161–164
    - distribution of growth, 127*f*
    - oxygen requirements for, 125
    - redox potentials, 131
  - Facultative thermophiles, temperature requirements, 116
  - Fasciolopsis buski*, 327
  - Fastidious microorganisms, 99
    - enriched media to cultivate, 106
    - Streptococcus* spp. as, 106
  - Fatty acids, 155
  - Febrile antibody test, 493–495, 493*f*
  - Febrile antigens, 493–495, 493*f*
  - Fecal contamination of water, 327, 337–340
    - testing, 329–333
  - Fermentation:
    - as biooxidative pathway, 161*f*
    - carbohydrate, 116, 154*f*, 161–164
    - cellular respiration vs., 161–164, 161*f*, 164*f*
    - of insulin, pneumococci and, 448, 449
    - lactose, 103, 106, 191–192, 193*f*
    - wine production and, 323–324
  - Ferrous sulfate, 167
  - Fertility factor (F factor), 371
  - Filamentous margin (cultural characteristic), 30*f*, 31
  - Filiform (cultural characteristic), 29, 30*f*
  - Filtration:
    - aseptic technique, 3*f*
    - to control microbial growth, 278*f*
  - Fine-adjustment knobs, of microscope, 37, 38*f*
  - Fission, in yeasts, 247
  - Flagella, protozoan, 221
  - Flat elevation (cultural characteristic), 30*f*, 31
  - Fleming, Alexander, 303
  - Flocculent growth, 31
  - Fluid thioglycollate medium, 131, 132*f*, 133, 133*f*, 134
  - Flukes, waterborne diseases and, 327
  - Fluorescent antibody technique, 499–500, 500*f*
  - Fluorescent microscopes, 35, 36
  - Fluorescent treponemal antibody absorption (FTA-ABS), 499
  - Food handlers, as contamination source, 317
  - Food microbiology, 153, 317–323
    - bacterial count analysis, 319–320
    - in cheese production, 344
    - role of *Saccharomyces* in, 247
    - in wine production, 247, 323–324, 323*f*
  - Food poisoning, 453
  - Food products, 317–323, 320*f*
  - Food supplements, 247
  - Food utensils, as contamination source, 317
  - Form (cultural characteristic), 29
  - Formaldehyde, mechanism of action/uses, 311*t*
  - Four-way streak-plate technique, 19*f*–20*f*
  - Francisella* spp., 473
  - Free radicals, 125
  - Free spores, 87, 88*f*
  - Free-flowing steam, 278*f*, 279
  - Free-living protozoa, 223–224, 224*t*
  - Freezing, to control microbial growth, 278*f*
  - Fructose, wine production and, 323–324
  - FTA-ABS (fluorescent treponemal antibody absorption) test, 508
  - Fungal infections, ocular micrometers and the diagnosis of, 52
  - Fungi, 237–259, 238*t*, 247–254. *See also* Yeasts
    - cultivation, of molds, 239–242
    - identification of unknown, 255–258
    - industrial importance of, 237
    - molds, 239. *See also* Molds
    - morphology/microscopic appearance, 255–257*t*
    - pathogenic, 239
    - pH requirements, 121
    - resistance to moist heat, 279
    - soil sample enumeration, 345–348, 347*f*
    - toxins of, 237
    - yeasts, 247–254. *See also* Yeasts

Fungi Imperfecti (Deuteromycetes), 237, 238t  
*Fusarium* spp., 255, 256t  
as food-borne organism, 317  
morphology/microscopic appearance, 256t

## G

Galactose, 1  
Gametocytes, *Plasmodium*, 230  
Gamma hemolysis, 106, 107, 107f  
identifying streptococcal pathogens by, 439, 440f  
Gamma radiation to control microbial growth, 278f, 289, 289f  
Gases:  
in carbohydrate fermentation, 163, 163f  
in lactose fermentation, 191, 193f  
GasPak anaerobic system, 131, 132f, 133–134, 134f  
GasPak jar, 463, 464f  
Gastroenteritis, 431, 453  
Gel electrophoresis, 388, 397, 398, 400–401, 401, 408, 409–411  
Gelatin hydrolysis, 154f, 156–157, 157f, 158, 216f  
Gelatinase, 157, 431  
Generation time of bacterial growth, 146, 147f  
Genetic engineering, 387–388, *See also* Biotechnology  
Genetic recombination experiment, 371–373  
Genetic variability mechanisms, 363–364  
Genetics (bacterial), 363–364  
chemical carcinogen screening, 363–364, 381–383  
enzyme induction, 365–367, 366f  
isolation of mutants, 377–378  
recombination via conjugation, 371–373  
Genital herpes, 507  
Genital herpes, isolation of herpes simplex virus (HSV), 509  
Genital warts, 507  
Gentamicin, 299t  
Germ theory of disease, 417  
Germination, 87, 88f  
*Giardia intestinalis*, 221, 229, 231, 231t, 234f, 234t  
waterborne disease and, 327  
*Giardia lamblia*, *See* *Giardia intestinalis*  
Glomerulonephritis, 439, 467  
*Glossina* (tsetse fly), 229  
Glucose:  
in cellular energy production, 155, 173  
in chemically defined media, 99  
degradation in Embden-Meyerhof (glycolytic) pathway, 161, 162f  
Glucose fermentation, 174f, 175f  
carbohydrate fermentation and, 161–164  
methyl red test, 173–174, 177f  
in saliva, dental caries and, 419, 420f  
triple sugar-iron agar test, 167, 168f  
wine production and, 323–324  
Glucose salts broth, 99, 101  
Glycerol, 155  
Glycolytic (Embden-Meyerhof) pathway, 161, 162f  
Glycosidic bonds, 155  
Golden Era of Microbiology, 417  
Gonorrhea, 507  
Gradient-plate technique, 377–378, 377f  
Grain stain, characteristics of selected bacteria, 216f  
Gram, Hans Christian, 75  
Gram stain, 75, 75–80  
counterstain, 75

as diagnostic staining procedure, 77  
procedure, 78f  
reagents, 76f  
staining procedure, 76f  
steps in, 77  
Gram-negative identification, 479, 482–483f  
Gram-positive cells, peptidoglycan layer, 75  
Gram-positive identification, 479, 480–481f  
Gram-positive stain, of streptococci, 75f  
Gram's iodine, 76, 76f  
Grapes, fermentation of, 323–324, 323f  
Griffin, Fred, 389, 390f  
Group A streptococcal pathogens, 439, 441  
Group B streptococcal pathogens, 441  
Group C streptococcal pathogens, 439–440  
Group D streptococcal pathogens, 441–442  
Growth, *See* Microbial growth control  
Guinea worm (*Dracunculus medinensis*), 327

## H

*Haemophilus influenzae*, 91, 389, 473  
*Haemophilus influenzae* D (*HindIII*), 408, 409t  
Hallucinogens, fungi-produced, 237  
*Halobacterium salinarium*, 285–286  
Halogens, mechanism of action/uses, 310t  
Halophiles, 285  
Hand washing, effectiveness of, 7–9  
Hanging-drop preparation to view living microbes, 45–47  
Hay infusions, 50f  
microscopic examination in the diagnosis of, 45–46  
Heat:  
dry (hot air), 278f, 279  
moist, to control microbial growth, 3f, 278f, 279–281  
Heavy metals, mechanism of action/uses, 310t  
Hektoen enteric agar, 454f  
Helminth diseases, 327  
Hemolysins, 431, 440  
Hemolysis, 106, 107, 107f  
Hemolysis tests, to differentiate staphylococcal species, 432t  
Hepatitis B, 507  
Herpes encephalitis, 509  
Herpes genitalis, 509  
Herpes labialis, 509  
Herpes simplex virus (HSV):  
culturing, 509  
detection of, 509  
isolation of, 509  
Herpes simplex-2 virus (HSV-2), 509  
*Heteronema*, 221, 224t  
Heterotrophs:  
carbon and, 97  
protozoa as, 221  
Hexachlorophene, mechanism of action/uses, 309t  
Hexylresorcinol, mechanism of action/uses, 309t  
High-frequency (Hfr) recombinants, 371–373  
Histidine, 381  
*Histoplasma capsulatum*, 52  
Hormones, genetic engineering of, 387–388  
Hospital-acquired infections,  
preventing, 7  
HSV-1, 509  
HSV-2, 509  
Hucker's (crystal violet) stain, 66, 76, 76f, 105–106, 108–109  
Humoral immunity, 487

Hydrogen gas, of lactose fermentation, 191, 193f  
Hydrogen peroxide, 201, 289  
Hydrogen sulfide:  
fermentative pathways, 183, 184f  
triple sugar-iron agar test, 167, 168f  
Hydrogen sulfide test, 154f, 173–174, 174f  
Hydrolytic enzymes, 153  
metabolic activities, 155–158  
Hypertonic solutions, 285, 286f  
Hypertonicity, 278f  
Hyphae, spore and vegetative, 239, 239f  
Hypotonic solutions, 285, 286f  
Hypotonicity, 278f

## I

Identification of microorganisms:  
blood specimen analysis, 473–476  
of *Campylobacter* strains, 463–464, 464f  
enrichment culture technique and, 359  
of enteric bacteria, 154f, 167, 168f, 171–177, 453–458  
of fungi, 247f, 248f, 255–257f  
of molds, 255–257t  
staphylococcal pathogens, 431–438  
of streptococcal pathogens, 439–444  
of *Streptococcus pneumoniae*, 447–449, 447f  
of unknown bacterial genus, 215–216, 216f  
of unknown bacterial species, 479–483  
urine specimen pathogens, 467–470  
Illumination, microscopes and, 37, 39–40  
Immunity, 487  
Immunology, 487–514  
agglutination reactions, 493–495, 493f, 507–508, 507f  
antibody titer test, 493, 494, 495f  
enzyme-linked immunosorbent assay (ELISA), 503–505  
fluorescent antibody technique, 493–495, 493f  
precipitin reaction (ring test), 489–490, 489f, 490f  
rapid plasma reagin test for syphilis, 507–8  
Impetigo, 431  
IMViC test series, 154f, 171–181  
citrate utilization test, 175–176  
indole production test, 171–173  
methyl red test, 173–174  
summary of reactions, 177f  
Voges-Proskauer test, 174–175, 175f  
Incubators, 2f, 4  
Indole production test, 154f, 171–172, 173f, 216f  
reactions summary, 177f  
Induced mutations, 363  
Induction, enzyme, 365–367, 366f  
Industrial microbiology, 153  
enrichment culture technique used in, 357  
soil microbes and, 344  
Infectious diseases, ELISA test to diagnose, 503–505  
Infectious process:  
of animal viruses, 263–264  
of bacteriophages, 261–262, 263f  
Infundibuliform (cultural characteristic), 30f, 31  
Innate (nonspecific) immunity, 487  
Inoculating needle/loop, 15  
Inoculation:  
aseptic, 15  
shake-tube, 125–126, 127f, 132f

- Inorganic salts, 98  
 Inorganic synthetic broth, 99, 101  
 Instruments, for culture transfers, 3, 5f, 13, 14f  
 Insulin:  
   fermentation test, 449  
   genetically engineered, 387  
 Insulin fermentation test, 448  
 Interfacial (ring) test, precipitin reaction, 489–490, 489f, 490f  
 Intestinal flora:  
   normal, 171, 423, 453  
   TSI reactions, 167–168, 168f  
 Intracellular enzymes (endoenzymes), 153–154, 154, 154f  
 Iodine compounds, mechanism of action/uses, 310t  
 Ionizing radiation to control microbial growth, 3f, 289  
 Iris diaphragm (microscope), 37, 39, 40  
 Irregular (cultural characteristic), 21, 29, 30f  
 Isolation of cultures, as diagnostic technique, 21  
 Isolation of herpes simplex virus (HSV), 509  
 Isolation of microorganisms:  
   of antibiotic-producing microbes, 351–353  
   of bacterial plasmids, 397–401  
   blood specimen analysis, 473–476  
   of *Campylobacter* strains, 463–464, 464f  
   of coliphages from raw sewage, 271–276  
   of *E. coli* mutant, 377–380  
   of enteric bacteria, 453–458  
   of *Pseudomonas*, 357–362  
   of streptomycin-resistant mutant, 377–381  
   urine specimen pathogens, 467–470  
 Isopropyl alcohol, mechanism of action/uses, 310t  
 Isotonic solutions, 285, 286f
- K**
- Kanamycins, 299t, 344  
 Keratoconjunctivitis, 509  
 Keto acid phenylpyruvic acid, 211, 212f  
 Kinyoun method, 81  
 Kirby-Bauer antibiotic sensitivity test, 294–296, 294f, 296f, 297  
*Klebsiella pneumoniae*, 91, 456, 470  
   biochemical activities, 154  
   cultural and biochemical characteristics, 216f  
   IMViC test series and, 171  
   methyl red test, 173–174, 177f  
   Voges-Proskauer test, 174–175, 177f  
*Klebsiella* spp.:  
   amino acid utilization, 209  
   identification of, 456f  
   as normal intestinal flora, 454  
   TSI reaction, 168f  
   Voges-Proskauer test, 174–175, 177f  
 Kovac's reagent, 172, 172f  
 Krebs cycle, 162f, 163, 164f
- L**
- Laboratory techniques (overview):  
   aseptic techniques, 3, 3f  
   cultivation chambers, 2f, 4–6  
   cultural characteristic designations, 29–32, 30f  
   culture transfers, 13–15, 22  
   equipment, 2–6, 2f, 4f, 5f  
   isolation of discrete colonies, 19–21  
   isolation of pure cultures, 2f
- media, 1–3, 2f, *See also* Media  
   Petri dishes, 2, 3  
   refrigerator uses, 6  
   test tubes, 2, 3, 4f  
   transfer instruments, 2f, 3–4  
 Lactic acid, 191  
   dental caries and, 419  
 Lactic acid bacteria, use of pyruvic acid by, 162f  
*Lactobacillus acidophilus*, dental caries and, 419  
*Lactobacillus* spp.:  
   as food-borne organism, 317  
   use of pyruvic acid by, 162f  
*Lactococcus lactis*:  
   biochemical activities, 154  
   cultural and biochemical characteristics, 216f  
   IMViC test series and, 171f  
   litmus milk reaction, 191–194  
 Lactose, 365, 366f  
 Lactose fermentation, 106  
   in coliform bacteria testing, 329, 331  
   enteric microbes capable/incapable of, 171f  
   in litmus milk reaction, 191, 193f  
   triple sugar-iron agar test, 167, 168f  
   by various bacterial, 216f  
 Lactose (*lac*) operon, 365, 366f, 391  
 Lag phase, in typical bacterial growth curve, 145, 145f  
 Lambda (temperate) phages, 262, 409–411, 409t  
 Lancefield serogroups, 439, 440f  
 Lancefield streptococcal serogroups, 439, 440f  
 Latex agglutination test, 433–435, 438  
 Leeuwenhoek, Antoni van, 35  
 Leucine (leu), 371  
*Leucostoc mesenteroides*, capsule stain, 91, 92f  
 Leukocidin, 431, 440  
 Light repair system, 290  
 Light sources in microscopes, 37, 39–40  
 Lipases, 155, 156, 156f, 431  
 Lipid hydrolysis, 154f, 155–156, 156f, 158, 158f, 216f  
 Lipids, 155  
 Lipolysis, 156  
 Litmus, 191  
 Litmus milk reactions, 154f, 191–194, 193f  
   enrichment culture technique and, 359  
   of various bacteria, 216f  
 Litmus reduction, 191, 192f, 193f  
 Liver, chemical carcinogens and, 381  
 Lobar pneumonia, 439, 447  
 Lobate margin (cultural characteristic), 29, 30f  
 Logarithmic (log) phase, 145, 145f  
*Lymphogranuloma venereum* (LGV), 507, 510–511  
 Lyophilization, to control microbial growth, 278f  
 Lysine decarboxylase, 209, 210  
 Lysis, 285  
 Lysogenic cells, 263  
 Lysogenic life cycle of bacteriophages, 263f  
 Lytic life cycle of bacteriophages, 262, 263f
- M**
- MacConkey agar, 105–106, 107f, 108–109, 111, 126, 454f, 463  
 Macrogametocytes, 230, 232f  
 Maculopapular rash, 507  
 Magnification, linear, 39f  
 Magnification principles, 37–38  
 Malachite green stain, 87  
 Malaria, 229–230, 232f, 233f  
 Maltase, 155  
 Maltose, 155  
 Mandelic acid, 357–362  
 Mannitol salt agar, 105–106, 107f, 108–109  
   detecting skin flora on, 424, 425f  
   to differentiate staphylococcal species, 432, 432t  
   inoculation procedure to cultivate bacteria, 108f  
 Margins (cultural characteristic), 29–31, 30f  
 Mastigophora, 221, 224t  
   parasitic, 230–231t  
 Maximum growth temperature, 115  
 Mechanical barriers, immunity and, 487  
 Mechanical pipette aspirators, 4, 5f  
 Media, 1–3, 2f  
   agar and, 1  
   chemically defined, 99  
   complex, 99  
   differential/selective, 105–106, 107f  
   enriched, 106–107, 107f  
   measuring turbidity as indicator of growth, 99–100  
   for routine cultivation of bacteria, 99–101  
   selective, 105–106, 106f  
   specialized, 99–100, 105–109  
 Medical microbiology, 417  
   blood specimen analysis, 473–476, 474f  
   *Campylobacter* isolation/presumptive identification, 463–464  
   dental caries susceptibility test, 419–420, 420f  
   enrichment culture technique used in, 357  
   genus identification of unknown culture, 215–217, 216f  
   identifying enteric bacterial pathogens, 453–458  
   identifying normal skin flora, 424–425, 425f  
   identifying normal throat flora, 423–424, 424f  
   identifying staphylococcal pathogens, 431–438, 432t  
   identifying streptococcal pathogens, 439–444  
   identifying *Streptococcus pneumoniae*, 447–449, 447f  
   species identification of unknown culture, 479–483  
   urine specimen analysis, 467–470, 473–476, 474f  
 Meiosis, 371  
 Membrane filters, 337  
   in quantitative water analysis, 337–340, 339f  
 Meningitis, neonatal, 439  
 Merbak, mechanism of action/uses, 310t  
 Mercurochrome, mechanism of action/uses, 310t  
 Mercury compounds, mechanism of action/uses, 310t  
 Merozoites, 229, 230, 232f  
 Merthiolate, mechanism of actions/uses, 310t  
 Mesophiles, optimum growth temperature, 115f, 116  
 Metabolic energy sources for microbes, 98  
 Metabolism (cellular), 153, *See also* Biochemical activities  
 Metabolites, 294  
 Metallic ions needed by microbes, 98  
 Metaphen, mechanism of action/uses, 310t  
 Methyl red test, 154f, 171, 173–174, 177f, 216f



- Methylene blue stain, 56f, 81, 82f, 83
- MIC (minimal inhibitory concentration), 303–304
- Microaerophiles:
  - distribution of growth, 127f
  - oxygen requirements, 125, 127f
- Microbial flora, 1, 7
- Microbial growth:
  - bacterial growth curves and, 145–148
  - generation time and, 145, 145f
  - indicators of, 99–100
  - media used to cultivate, 99–100
  - nutritional needs, 97–98
  - oxygen needs, 125–127
  - pH needs, 98, 121–122
  - physical factors affecting, 98
  - temperature needs, 98, 115–117
- Microbial growth control, 277–313
  - chemical methods, 277, 278f, 293–313, 309–313
    - antiseptics, 277, 309–313
    - chemotherapeutic agents, 293–302
    - disinfectants, 277, 309–313
    - penicillin activity/inhibition, 303–305
  - microbicidal vs. microbistatic, 277
  - moist heat, 279–284
  - penicillin activity/inhibition, 303–305
  - physical methods, 277–284, 277–290, 278f
- Microbicidal effect, 277
- Microbiology, 35, 36
  - defined, 37
- Microbistatic effect, 277
- Micrococcus luteus*, 20, 21, 22, 25, 26, 31
  - biochemical activities and, 154
  - biochemical characteristics, 216f
  - cultivation of, 133
  - cultural characteristics, 216f
  - IMViC test series and, 171f
- Micrococcus* spp., as food-borne organism, 317
- Microgametocytes, *Plasmodium*, 230, 232f
- Microorganisms, 1, 13
  - biochemical activities, 153–219, 154f
  - biochemical activities of, 154
  - cultural characteristics, 29–32
  - fastidious, 99
  - fermentation abilities, 161–164
  - microscopic measurement, 51–54
  - motility of, 183–184
  - nutritional requirements, 99–102
  - oxygen (atmospheric) requirements, 125–127
  - visualization of, in living state, 55
- Micropipette, 5f
  - using, 398–399, 409
- Microscopes/microscopy, 35–54
  - Abbé condenser, 37
  - body tube, 37, 38
  - brightfield, 35, 36
  - compound, 38f
  - darkfield, 35
  - electron, 35, 36
  - examination of living microorganisms using a hanging-drop preparation or a wet mount, 45–50
  - examination of stained cell preparations, 39–42
  - fluorescent, 35, 36
  - illumination, 37, 39–40
  - linear magnification, 39t
  - magnification, 37–38
  - mechanical stage, 37
  - microorganisms, measurement of, 51–54
  - parfocal, 40
  - phase-contrast, 36
  - resolving power/resolution of lenses (microscope), 38–39
  - stage, 37
  - theoretical principles of, 37–40
  - use/care of, 40–42
- MicroTrak® Direct Specimen Test, 511
- Milk:
  - breed smears and, 138
  - litmus milk reactions, 191–194
- Minimal inhibitory concentration (MIC), 303
- Minimum growth temperature, 115
- Mixed culture, isolation of discrete colonies from, 19–21
- Moist heat:
  - as aseptic technique, 3f
  - autoclaving, 279–281, 280f
  - to control microbial growth, 3f, 278f, 279–281
  - pasteurization and, 278f, 279–280
  - sterilization and, 279–280
  - tyndallization and, 278f, 279–280
- Molds:
  - colonial appearance, 241–242, 257t
  - cultivation, 239–242
  - enumeration of soil, 346, 347f
  - identification, 255–256t
  - morphology, 239, 239f, 240f, 255–256t, 257t, 425f
  - pH requirements, 121
  - soil sample enumeration, 345–348, 347f
- Molecular biology, 387
- Monilia, 248
- Mononucleosis, 508
- Moraxella* spp., 423
- Mordant, 75, 76–77, 76f
  - Gram's iodine as, 76, 76f
- Morton test tube closures, 3
- Mosquito (*Anopheles*), 229, 232f
- Most probable number (MPN) test, 329, 332f, 334t
- Motility of microorganisms, detection of, 183–184
- Mouse virulence test, 448
- Mouth, normal flora of, 419–420, 420f, 423
- MPN (most probable number) test, 329, 332f, 334t
- Mucoid consistency (cultural characteristic), 29
- Mucor mucedo* 240f, 241, 242
- Mucor* spp., 241, 255
  - morphology/microscopic appearance, 255t
  - soil populations of, 345
- Mueller-Hinton agar, 294–295
  - detecting native skin flora, 425
  - detecting throat flora, 423, 424f
  - Kirby-Bauer Antibiotic Sensitivity Test, 294–295
- Multitest systems for identifying enteric microbes, 453–458
- Must (wine production), 323
- Mutagenicity, 381
- Mutations, genetic, 363–364, 377–378
- Mycelium, aerial and vegetative, 239
- Mycobacterium leprae*, 81
- Mycobacterium smegmatis*, 31, 351
  - acid-fast stain, 83
  - antiseptic susceptibility test, 312–313
  - Kirby-Bauer Antibiotic Sensitivity Test, 294–296
- Mycobacterium tuberculosis*, 81
- Mycology, 237
- Mycoplasma hominis*, 511
- Mycoses, superficial and systemic, 237
- ## N
- N*-acetylglucosamine, 75–76
- N*-acetylmuramic acid, 75–76
- Napiform (cultural characteristic), 30f, 31
- Needles, transfer, 2, 3, 5f
- Negative staining, 71–73
  - and acidic stains, 71
  - encapsulated invaders, detecting, 71
  - procedure, 71–72, 72f
- Neisseria gonorrhoeae*, 389
- Neisseria meningitidis*, 473
- Neisseria perflava*, 475
- Neisseria* spp.:
  - detecting on chocolate agar, 423, 424f, 425, 425f
  - oxidase test to differentiate, 205
- Neonatal herpes, 509
- Neonatal meningitis, 439
- Neubauer rulings, 137
- Neutrophiles, 121, 121f
- Nicked circular DNA, 397–398, 398f
- Nitrate reductase, 197
- Nitrate reduction test, 154f, 197–198, 198f
- Nitrates, 197, 198f, 343, 344, 344f
- Nitrification, by soil microbes, 343, 344f
- Nitrites, 197, 198f, 343, 344f
- Nitrobacter* spp., 343
- Nitrogen cycle, 343–344, 344f
- Nitrogen fixation, 343–344, 344f
- Nitrogen metabolism, 209
- Nitrogen requirements, microbial, 97
- Nitrosomonas*, 343
- Nocardia*, acid-fast stain, 83
- Nongonococcal urethritis (NGU), 507, 510
- Nonmetallic elements needed by microorganisms, 97–98
- Normal flora:
  - intestinal, 171, 453
  - locations in human body, 423
  - of mouth, 419–420, 420f, 423
  - of skin, 423–427
- Nosepiece of microscope, 37
- Nosocomial infections, preventing, 7
- Novobiocin sensitivity test, 432, 432f, 432t, 433t
- Nuclear stain, 57f
- Nucleases, 440
- Nucleic acids, 55–56, 98
- Nucleoproteins, 261
- Numerical aperture of microscope lens, 37
- Nutrient agar plates:
  - cultural characteristics of, 29–31, 30f
  - E. coli* and *S. aureus* growth on, 106f
  - inoculation of, 9t
- Nutrient agar slants, cultural characteristics of, 29, 30f
- Nutrient broth media, 99, 101
  - cultural characteristics of, 31
- Nutrient gelatin, cultural characteristics, 31
- Nutritional needs, of microbial cells, 97–98
- Nutritional requirements of microorganisms, 97–98
  - cultivation media and, 99–101
- ## O
- Objective lens of microscope, 37–38
  - oil-immersion, 39, 39f, 40
  - overall linear magnification, 39t
  - relationship between working distance, diaphragm opening and, 40, 41f

- Obligate anaerobes, atmospheric oxygen and, 125
- Obligate thermophiles, 116
- Ocular lens of microscope, 37
- Overall linear magnification, 39*t*
- Ocular micrometer, 51, 51*f*  
calibration and use of, 51*f*  
and the diagnosis of fungal infections, 52
- Oidiophore, 239*f*
- Oocyst, *Plasmodium*, 230
- Ookinete, 230, 232*f*
- Opaque (cultural characteristic), 29
- Operon, 365  
lactose (*lac*), 365, 366*f*, 391
- Optical characteristics, 29
- Optimum growth pH range, 121–122, 121*f*
- Optimum growth temperature, 115–117, 115*f*
- Optochin sensitivity test, 447–448, 448*f*, 449
- Organic mercurials, mechanism of action/uses, 310*t*
- Ornithine decarboxylase, 209
- Osmosis, 285–286, 286*f*  
using to control microbial growth, 278*f*, 285–286
- Osmotic pressure of the environment, 285
- Osteomyelitis, 431
- Oven heat, to control microbial growth, 278*f*
- Overall linear magnification, 39*t*
- Oxidase test, 154*f*, 205–206, 206*f*, 216*f*  
cytochrome oxidase, 205
- Oxidation-reduction (redox) reactions, 131, 131*f*
- Oxygen requirements, microbial, 98  
classifying according to need for, 125  
procedure for determining, 125–127, 127*f*
- P**
- PABA (p-Aminobenzoic acid), 294*f*  
p-Aminobenzoic acid (PABA), 294*f*  
p-Aminodimethylaniline oxalate, 205, 424
- Paracasein, 192
- Paraffin plug technique, 132*f*
- Paragonimus westermani*, 327
- Paramecium caudatum*, 221, 224*f*
- Paramecium* spp., 224*t*
- Parasites, viruses as, 261
- Parasitic helminths, 327
- Parasitic protozoa, 229–236, 230–231*t*  
characteristics, 230–231*t*  
life cycle, 229–230, 232*f*  
waterborne disease from, 327
- Paratyphoid bacilli, 106
- Paratyphoid fever, 327, 453
- Parfocal microscopes, 40
- Pasteurization, 278*f*  
moist heat, 279–280  
osmotic pressure, 278*f*, 285–286  
temperature in, 279*f*, 280
- Pathogenic microbes, 417–418, *See also* Medical microbiology  
blood specimen analysis, 473–476, 474*f*  
in enterobacteriaceae family, 171  
food-borne, 317–320, 453  
germ theory of disease and, 417  
IMViC test series to identify, 171–177  
normal flora vs., 171  
occasional pathogens, 171  
soil-borne, 344  
staphylococcal, 431–438  
streptococcal, 439–444, 440*t*  
transformation in, 389  
in urinary trace, 467–470  
in urinary tract, 467–470  
waterborne, 327  
yeasts as, 247
- Pathogens, 171
- Peaches, fermentation of, 323–324
- Pears, fermentation of, 323–324, 323*f*
- Pellicle, 30*f*, 31
- Pelvic inflammatory disease, 510
- Penetration (infection), 261–262  
in animal viruses, 264  
in bacteriophages, 261–262
- Penicillin, 293*t*, 344  
determining activity of, 303–305, 303*t*
- Penicillin G (benzylpenicillin), 299*t*  
molecular structure, 303, 303*t*
- Penicillinase (beta-lactamase), 303–305  
*Penicillium chrysogenum*, 240, 240*f*, 241, 242, 243  
*Penicillium notatum*, *See Penicillium chrysogenum*
- Penicillium* spp., 241, 255, 351  
in cheese production, 344  
as food-borne organism, 317  
morphology/microscopic appearance, 256*t*  
soil populations of, 345
- Peptide bonds, 156, 157*f*
- Peptidoglycan, 75–76
- Peptone, 99, 156, 163, 167
- Peptonization, 192, 193*f*
- Petri dishes, 2, 3  
turntable, 20*f*
- Petroff-Hausser counting chamber, 137, 137*f*
- PFU (plaque-forming unit), 265, 265*f*
- pH:  
dental caries and, 419–420  
requirements of microorganisms, 98, 121–122, 121*f*
- Phages, *See* Bacteriophages
- Phagocytosis, 487
- Phase-contrast microscopes, 36
- Phenol, mechanism of action/uses, 309*t*
- Phenol red, 106, 167
- Phenolic compounds, mechanism of action/uses, 309–310*t*, 310*t*
- Phenylalanine deaminase test, 209–211, 211*f*
- Phenylethyl alcohol agar, 105, 106*f*  
*E. coli* and *S. aureus* growth on, 106*f*
- Phosphorus requirements, microbial, 98
- Phototrophs, 98
- Physical factors affecting microorganisms:  
atmospheric oxygen, 98  
pH of environment, 98  
temperature, 98
- Physical factors affecting microorganisms:  
atmospheric oxygen, 125–130  
pH of environment, 121–124  
temperature, 115–120
- Physical methods to control microbial growth, 277–290, 278*f*  
electromagnetic radiation methods, 289–290  
moist heat, 279–284  
osmotic pressure, 285–286
- Picric acid, 55  
chemical formation of, 56*f*
- Pigmentation (cultural characteristic), 29  
temperature and, 116, 116*f*
- Pipettes, 4  
mechanical aspirators, 4, 5*f*  
micropipette, 398–399, 409
- Plaque, 265
- Plaque-forming unit (PFU), 265, 265*f*
- Plasmids, 371  
isolation of bacterial, 397–401
- Plasmodium* spp., 221
- Plasmodium vivax*, 221, 229–230, 231, 232*f*, 233*f*
- Plasmolysis, 285, 286*f*
- Plums, fermentation of, 323–324, 323*f*
- Pneumonia, 431  
lobar, 439, 447
- Point mutations, 363
- Polylinker, 397
- Polymyxin, 293*t*
- Polypeptides, 156, 192
- Population growth curve stages, 145, 145*f*
- Pour plate-loop dilution technique, 2*f*, 20
- Pour-plate technique, 20*f*  
to count viable cells, 138–139, 138*f*
- Povidone-iodine solution (Betadine), 310*t*
- Precipitin formation, 488
- Precipitin reaction (ring test), 489–490, 489*f*, 490*f*
- Precipitins, 488
- Pre-erythrocytic stage, 229, 232*f*
- Presumptive test (water analysis), 329, 330*f*, 332, 332*f*, 335
- Primary stain:  
acid-fast stain, 81, 82*f*, 83  
capsule stain, 90–91, 90*f*, 91*f*  
defined, 75  
gram stain, 76*f*, 77*f*  
spore stain, 87
- Proctitis, NGU and, 510
- Prophage, 263
- PROS (Profile Recognition System), 455
- Proteases, 156, 157*f*
- Protein hydrolysis, 157*f*
- Proteins:  
and acidic stains, 55–56  
as microbial energy source, 164*f*
- Proteins, as microbial energy source, 192
- Proteolysis, 156  
as litmus milk reaction, 192, 193*f*
- Proteus* spp., 423  
amino acid utilization, 209  
as food-borne organism, 317  
as normal intestinal flora, 453  
TSI reaction, 168*f*
- Proteus vulgaris*, 456, 470, 493  
biochemical activities, 154*f*  
decarboxylase test, 209–212  
hydrogen sulfide test, 154*f*, 183–184, 184*f*  
IMViC test series and, 171*f*, 172  
Kirby-Bauer Antibiotic Sensitivity Test, 294–296  
phenylalanine deaminase test, 209–211, 211*f*  
TSI agar test, 167–168, 168*f*  
urease test, 187–188, 188*f*
- Protozoa, 221–236  
free-living, 223–224, 224*t*  
as heterotrophs, 221  
parasitic, 229–236, 230–231*t*  
waterborne disease from, 327  
pond water, 223  
taxonomy, 221, 224*f*
- Pseudomonadales, soil populations of, 345
- Pseudomonas aeruginosa*, 31, 470, 473  
biochemical activities, 154*f*  
enozymatic activities, 157–158  
experiments:  
antibiotic-producing microbes and, 351–353  
IMViC test series, 171*f*  
litmus milk reaction, 191–194  
Kirby-Bauer Antibiotic Sensitivity Test, 294–298  
nitrate reduction test, 197–198, 198*f*  
oxidase test, 205–206, 206*f*  
TSI reaction, 167–168, 168*f*

- Pseudomonas denitrificans*, vitamins and, 344  
*Pseudomonas savastanoi*, cultivation of, 116  
*Pseudomonas* spp., 317, 344, 358f, 423  
 enrichment culture technique used in, 357–362  
 as food-borne organism, 317  
 oxidase test to differentiate, 205  
 TSI reaction, 168f
- Pseudopods, 221
- Psychrophiles, optimum growth temperature, 115f, 116
- Puerperal fever, 439
- Pure cultures, 1, 2f  
 apparatus used in, 1, 2f  
 defined, 1, 13, 15  
 isolation of, from spread-plate/  
 streak-plate preparation, 22  
 isolation techniques, 19–22, 359  
 procedure for the preparation of, 24f
- Pyelitis, 467
- Pyelonephritis, 431, 467
- Pyrogallic acid technique, 132f
- Pyruvic acid, variation in uses of, 161, 162f
- ## Q
- Quaternary ammonium compounds, mechanism of action/uses, 311t
- Quebec colony counter, 138–139, 139f, 140f
- Quellung (Neufeld) reaction, 448, 449
- Quinolone, 293t
- ## R
- Radiation:  
 to control microbial growth, 3f, 278f, 289–290  
 gamma, 289  
 ionizing forms of, 289  
 x-radiation, 289
- Raifampin, 293t, 299t
- Raised elevation (cultural characteristic), 30f, 31
- Rapid plasma reagin (RPR) test, 507, 507–8, 508f
- Rapid testing methods, 418, 432–433, 511, 515
- Reagent:  
 Barritt's, 174, 175f, 177f  
 Kovac's, 172
- Reagin, 507, 508f
- Recombinant DNA technology, 387–388, *See also* Biotechnology
- Redox (oxidation-reduction) reactions, 131, 131f
- Refractive index (microscope), 39, 39f
- Refrigerators, 2f, 6, 289–290  
 microbistatic effects, 278f
- Rennet curd, 192
- Rennin, 192
- Replication:  
 in bacteriophages, 262, 265  
 viral, 261
- Reproduction, in yeasts, 247, 248f
- Resident flora, 7
- Resistance (immunity), 487
- Resolving power/resolution of lenses (microscope), 38–39
- Resorcinol, mechanism of action/uses, 309t
- Respiratory tract:  
 infections of, 439  
 normal flora of, 423
- Restriction analysis, 407f, 408f, 409–411, 409f
- Restriction endonucleases, 407–410, 407f
- Reticuloendothelial system, 487
- Rheumatic fever, 439
- Rhizobium* spp., as soil microbes, 344  
 soil populations, 345
- Rhizoid (cultural characteristic), 29, 30f
- Rhizopus*, morphology/microscopic appearance, 255t
- Rhizopus* spp., as food-borne organism, 317
- Rhizopus stolonifer*: 240f  
 characteristics of, 238t  
 cultivation of, 240–242  
 morphology/microscopic appearance, 255t  
 as soil microbe, 344
- Rhodotorula rubra*, 246–247, 247f
- Rifampin, 293t, 299t
- Ring test, precipitin reaction, 489–490, 489f, 490f
- Ringworm infection, 237
- RNA (ribonucleic acid), 98
- RNA viruses, 261
- ## S
- Sabouraud agar plate, 239, 240, 240f, 242, 248  
 detecting normal skin flora on, 425, 425f
- Sabouraud broth culture, 116
- Sabouraud broth plate, 255t
- Saccate (cultural characteristic), 30f, 31
- Saccharomyces carlsbergensis*:  
 brewing industry uses, 247  
 use of pyruvic acid by, 162f
- Saccharomyces cerevisiae*, 40, 247, 248f, 248–249, 475  
 brewing industry uses, 247  
 budding, 246f  
 carbohydrate fermentation by, 126  
 cultivation of, 116, 122, 126  
 fermentation study, 249  
 food industry uses, 247  
 moist heat to control, 280–281, 281f  
 morphology, 247f, 248–249  
 reproduction study, 248f, 249  
 temperature requirements, 116, 126  
 use of pyruvic acid by, 162f  
 wine production and, 323–324
- Saccharomyces ellipsoideus*, wine industry uses, 247
- Safranin, 76–77, 94
- Saliva, glucose fermentation in, 419
- Salmonella* spp.:  
 amino acid utilization, 209  
 as food-borne organism, 317  
 as human pathogen, 453, 473  
 TSI reactions, 168f
- Salmonella typhi*, 327, 473  
 use of pyruvic acid by, 162f
- Salmonella typhimurium*, 453, 456, 464  
 Ames test on two strains, 381–383, 383f  
 antibody titer test, 493, 494, 495f  
 biochemical activities, 154f  
 carbohydrate fermentation, 163–164, 163f  
 cultivation of, 108  
 hydrogen sulfide test, 154f, 183–184, 184f  
 IMViC test series and, 171f  
 TSI agar test, 167–168, 168f
- Salmonella typhiphimurium*, 456, 464, 473
- Salpingitis, NGU and, 510
- Sarcodina, 221, 224t, 230t
- Scanning electron microscopes, 36
- Scarlet fever, 439
- Schistosoma haematobium*, 467
- Schistosomiasis, 327
- Schizogony, 229
- Schizonts, 230, 232f
- Schizosaccharomyces*, 247
- Schizosaccharomyces octosporus*, 248f, 248–249
- Sediment, 30f, 31
- Segmenters, 230, 232f
- Selective media, 105–106, 106f
- Selenitla intestinalis*, 248–249
- Semisolid medium, 1, 2f
- Semmelweis, Ignatz, 7
- Sensi-Disc™ dispensers, 295
- Separation of bacteria, and staining techniques, 57
- Septic shock, 473
- Septicemia, 431  
 isolation/identification of pathogen, 473–476, 474f
- Septi-Chek System, 475, 475f
- Serial dilution-agar plate analysis, 138–139, 138f–141f
- Serogroups, 439
- Serologic classification of Lancefield, 439–440, 440t
- Serology, 488  
 ring (interfacial test), 489–490, 490f
- Serrate margin (cultural characteristic), 30f, 31
- Serratia marcescens*, 15, 20, 22, 25, 26, 473  
 cultivation of, 116, 117  
 four-way streak-plate technique, 19f, 20f  
 radiation to control, 289–290  
 temperature and pigment production in, 116
- Serratia* spp., as food-borne organism, 317
- Sewage (raw):  
 coliform bacteria testing, 331–333, 337, 338  
 coliphage isolation, 271–273, 272f
- Sexually transmitted chlamydial diseases, detection of, 510
- Shake-culture technique, 127f, 132f
- Shake-tube inoculation, 125–126, 127f, 132f
- Shaking waterbaths, 4–6
- Shigella dysenteriae*, 168f, 327, 456, 493  
 biochemical activities experiments, 154  
 cultural and biochemical characteristics, 216f  
 hydrogen sulfide test, 154f, 183–184, 184f
- Shigella dysenteriae*, IMViC test series and, 171f
- Shigella dysenteriae*, TSI agar test, 167–168, 168f
- Shigella* spp.:  
 as food-borne organism, 317  
 as human pathogen, 453
- Shigella* spp., TSI reactions, 168f
- Shigellosis, 453
- Signet rings, 230, 232f
- Silver compounds, mechanism of action/uses, 310t
- Silver nitrate, mechanism of action/uses, 310t
- SIM agar, 172, 183
- Simple staining, 57f, 65–67  
 clinical application, 65  
 procedure, 66, 66f  
 purpose of, 65
- Size (cultural characteristic), 29
- Skin:  
 identifying normal flora, 423–426  
 infections, 431, 439  
 infections of, 439  
 sodium chloride (6.5%) test, 442, 442f  
 sodium chloride (7.5%) agar, 105
- Skin lesions, 431
- Slant-to-agar deep transfer, 13, 14f
- Slant-to-broth transfer, 13, 14f

- Snyder test, 419–420, 420f
- Soaps/detergents, mechanism of action/uses, 311t
- Sodium chloride (6.5%) test, 442, 442f
- Sodium chloride (7.5%) agar, 105
- Sodium hypochlorite, mechanism of action/uses, 310t
- Sodium tetradecyl sulfate, mechanism of action/uses, 311t
- Soil microbiology, 343–359
- carbon cycle, 344
  - enumeration methods, 345–348, 347, 347f
  - food contamination and, 317
  - isolation of antibiotic-producing microbes, 351–353
  - nitrogen cycle, 343–344, 344f
  - Pseudomonas* isolation by enrichment culture technique, 357–362
  - sulfur cycle, 344
- Solid medium, 2f
- Solutions, hypertonic/hypotonic/isotonic, 285, 286f
- Sonic vibrations, to control microbial growth, 278f
- Spectrophotometer:
- to determine cell counts, 138
  - to measure microbial growth, 100, 100f, 101–102, 101f
  - to plot bacterial growth curves, 148f
- Spirilla (spiral-shaped) bacteria, 65f, 66, 67f
- Spirochetes, 508
- Spontaneous mutations, 363
- Sporangia (mold), 239f
- Sporangiophore, 239f
- Spore coats, 87
- Spore stain, 93f
- counterstain for, 87–88
  - primary stain for, 87
  - Schaeffer-Fulton method, 87–89
- Spores, 87
- fungal (mold), 239f
- Sporogamy, 230, 232f
- Sporogenesis, 87
- Sporozoa, 221, 229–230, 232f
- Sporozoites, 229, 230, 232f
- Spread-plate technique, 2f, 373f
- Stage micrometer, 51, 51f
- Stained cell preparations, examination of, 39–42
- Staining techniques, 57, 57f
- simple, 65–67, 66f
- Stains:
- acidic, 55, 57f
  - basic, 55–56, 57f
  - chemical composition of, 55f
  - defined, 55
  - Gram, 75–80
  - picric acid, 55
  - staining techniques, 56, 57f
- STAPH-IDENT system, 432–433, 434, 435
- profile register for, 435–436t
- Staphylococcal enteritis, 431
- Staphylococcal pathogens, identifying, 431–436, 431f, 431t, 432t
- Staphylococcus aureus*, 40, 65, 83, 473
- biochemical activities, 154
  - carbohydrate fermentation, 163–164, 163f
  - cultural and biochemical characteristics, 216f
  - detecting on mannitol salt agar, 424, 425f, 432
  - DNase test, 432, 433f
  - enozymatic activities, 157–158
  - experiments, 303–305
  - acid-fast stain, 83
  - antibiotic-producing microbes and, 351–353
  - coagulase test, 432, 433f
  - cultivation of, 106, 106f, 108, 126
  - detecting on mannitol salt agar, 424, 425f, 433
  - drug combinations, 296–298, 297f
  - IMViC test series, 171f
  - moist heat to control growth, 279–281, 281f
  - osmotic pressure to control, 285–286
  - penicillin activity, 303–305, 304f
  - isolation of, 21
  - Kirby-Bauer Antibiotic Sensitivity Test, 294–298
  - laboratory tests to identify, 432t
  - osmotic pressure to control, 285–286
  - oxygen requirements, 126
- Staphylococcus epidermidis*, 423, 431, 433
- cultivation of, 108
  - laboratory tests to identify, 423, 432t
  - Novobiocin sensitivity test and, 432, 433f
- Staphylococcus saprophyticus*, 431, 433
- Staphylococcus* spp.:
- detecting with sodium chloride agar, 105
  - as food-borne organism, 317
  - identifying normal skin flora, 423, 424f
  - identifying pathogenic, 431–436, 431f, 473
  - shape and arrangement of, 67f
- Staphylokinase, 431
- Starch, 155
- Starch agar, 155
- Starch agar plate, 155f
- Starch hydrolysis, 154f, 155, 155f, 158, 216f
- Stationary phase, in typical bacterial growth curve, 145, 145f
- Steam heat, to control microbial growth, 278f, 279–280
- Stemphylium* spp., soil populations, 345
- Stentor* spp., 224t
- Sterilization:
- to control microbial growth, 278f, 279–280
  - defined, 3
  - by dry heat, 279
  - by moist heat, 279–281
  - techniques, 3f
- Steroids, microbe-produced, 344
- Stock cultures, 22
- Stratiform (cultural characteristic), 30f, 31
- Streak-plate technique, 2f
- Streptococcal pneumoniae*, 447–449, 447f
- Streptococci, Gram-positive stain of, 75f
- Streptococcus agalactiae*, 440t, 442
- Streptococcus bovis*, 439, 440t
- Streptococcus equi*, 440t
- Streptococcus mitis*, 440t, 442t
- cultivation of, 100, 108
- Streptococcus mutans*, 419, 420f
- dental caries susceptibility and, 419
- Streptococcus pneumoniae*, 389, 390f, 423, 439, 447–449, 447f
- as cause of lobar pneumonia, 447
- Streptococcus pyogenes*, 21, 439, 440t, 442, 499
- Streptococcus sanguis*, 440t
- Streptococcus* spp., 107, 107f, 417, 423
- hemolytic activities of, 439, 440t
  - identifying pathogenic, 439–444, 440t, 473
  - use of pyruvic acid by, 162f
- Streptococcus* var. Lancefield Group E, 108
- Streptococcus viridans*, 439, 447
- Streptolysin O, 107
- Streptolysin S, 107
- Streptomyces griseus*, 351
- Streptomyces* spp.:
- as antibiotic-producers, 344, 351
  - soil sample enumeration, 345–348, 347f
- Streptomycin, 293t, 299t, 344
- isolation of mutant resistant to, 363, 377–378
- Stretobacillus, 65f
- Subacute endocarditis, 439
- Subculturing, 3, 5f
- defined, 13
  - procedure, 14f
- Sucrose fermentation, 216f
- triple sugar-iron agar test, 167–168, 168f
- Sulfadiazine, 294, 294f
- Sulfanilamides, chemical similarity to PAB, 294f
- Sulfanilic acid, 197, 198f
- Sulfisoxazole-trimethoprim drug combination, 297, 298
- Sulfonamides, 299t
- Sulfur, oxidation by soil microbes, 344
- Sulfur requirements, microbial, 97–98
- Superoxide dismutase, 201
- Surface-active agents, mechanism of action/uses, 311t
- Surgical hand washing, 7
- Synergism, 296
- Synergistic effect of drug combinations, 296–298, 297f
- Synthetic drugs, 293–294
- Syphilis, 499, 507, 511
- defined, 507
  - rapid plasma reagin test for, 507–8
- ## T
- Taxonomy:
- biochemical characteristics and, 153–154, 215–217, 216f
  - cultural characteristics, 216f
  - defined, 215
  - protozoa, 221, 224f
- Teeth:
- dental caries susceptibility experiment, 419–420, 420f
  - normal flora of, 423
- Temperate (lambda) phages, 262, 407–411, 409t
- Temperature:
- agar and, 1
  - in autoclaving, 279–281, 280f
  - cellular enzymes affected by, 115, 279
  - cultivation chambers, 4, 6
  - moist heat and, 278f, 279–281
  - in pasteurization, 279f, 280
  - requirements for cultivating microbes, 2, 98, 115–117
  - in tyndalization, 279–280
- Tergitol, 81
- Test tubes, 2, 3, 4f, 4f
- Tetracycline, 293t, 299t, 344
- drug combination with trimethoprim, 297, 298
- Tetracyclines, 293t
- TFTC (too few to count) designation, 142
- Thermophiles, optimum growth temperature, 115f, 116
- Thiamine (thi), 371
- Thiosulfates, in hydrogen sulfide production, 183, 184f
- Threonine (thr), 371
- Throat, identifying normal flora, 423–427

- Throat infection, 439  
 Thrush infection, 248  
 Thymine dimerization, 278*f*, 289–290  
 Thymol, mechanism of action/uses, 310*t*  
 Tincture of green soap, mechanism of action/uses, 311*t*  
 Tincture of iodine, mechanism of action/uses, 310*t*  
 Tinsdale agar plate, 423, 424*f*, 425  
 TNTC (too numerous to count) designation, 142  
 Tobramycin, 299*t*  
 To-deliver pipette, 5*f*  
 Tonsillitis, 439  
 Too few to count (TFTC) designation, 9, 142  
 Too numerous to count (TNTC) designation, 9, 142  
*Torula* spp., 255, 257*t*  
 Toxins, fungi-produced, 237  
 TPI (*Treponema pallidum* immobilization) test, 508  
 Transduction, genetic, 364, 371, 391  
 Transfer instruments, 3–4, 5*f*  
 Transfer loop, 5*f*  
 Transfer needle, 5*f*  
 Transfer of genetic material, 364  
 Transformation, genetic, 364, 371, 389–394, 389*f*  
   F. Griffin's experiments, 389*f*, 390*f*  
   natural, 389  
 Transient flora, 7  
 Translucent (cultural characteristic), 29  
 Transmission electron microscopes, 36  
 Transparent (cultural characteristic), 29  
*Treponema pallidum*, 499, 507  
 Tributyrin agar plate, 156, 156*f*  
*Trichomonas* spp., 423  
*Trichomonas vaginalis*, 221, 467, 511  
 Trichomoniasis, 507  
*Trichothecium* spp., as food-borne organism, 317  
 Triglycerides, 155, 156*f*  
 Trimethoprim, 299*t*  
 Trimethoprim/sulfamethoxazole, 299*t*  
 Trimethoprim/sulfisoxazole, 297, 298  
 Trimethoprim/tetracycline, 297, 298  
 Triple sugar-iron (TSI) agar test, 154*f*, 167–168  
*Trophozoites*, 230, 230–231*t*, 232*f*  
*Trypanosoma gambiense*, 221, 229, 230*f*, 230*t*, 231, 234*f*  
*Trypanosoma* spp., 221, 229  
 Trypticase soy agar, 100, 108  
 Trypticase™ soy agar plates, 21  
 Tryptophan, indole production test, 171–172, 173*f*, 177*f*  
 Tryptophanase, 171, 172*f*  
 Tsetse fly (*Glossina*), 229  
 TSI (triple sugar-iron) agar test, 154*f*, 167–168  
 Tuberculosis, microscopic examination in the diagnosis of, 42  
 Turbidity of cultures, 100  
   as growth indicator, 99, 138  
   spectrophotometric analysis, 100–101, 100*f*, 101*f*, 138  
 Tween™, 295  
 Tyndalization, 279–280  
 Typhoid bacilli, 106  
 Typhoid fever, 327, 453
- U**  
 Ultraviolet light, 289, 289*f*  
 Umbonate elevation (cultural characteristic), 30*f*, 31  
 Undulate margin (cultural characteristic), 29, 30*f*  
 Uniform fine turbidity, 31  
 Upper respiratory tract:  
   infections of, 439  
   normal flora of, 423  
 Urea:  
   broth, 187*f*  
   enzymatic degradation, 187*f*  
 Urea, enzymatic degradation, 187  
*Ureaplasma urealyticum*, 510  
 Urease, 187, 187*f*  
 Urease test, 154*f*, 187–188, 188*f*  
 Urinary tract infections, 431, 439–440, 467  
   normal flora of, 423  
   pathogenic possibilities, 467  
 Urine specimen analysis, 463*t*, 467–470, 468*f*, 469*t*
- V**  
 Vaginitis, 467  
 Vancomycin, 299*t*  
 VDRL (Venereal Disease Research Laboratory) agglutination test, 507  
 Vegetative cells, 87, 88*f*  
 Vegetative mycelium, 239  
 Venereal Disease Research Laboratory (VDRL) agglutination test, 508  
*Vibrio cholerae*, 327  
*Vibrio* spp., 463  
 Virulent bacteriophage particles, 262–263, 263*f*  
   isolation of, 261–273, 271–276  
 Viruses, 261–273  
   animal, 263–264  
   characteristics, 261  
   coliphage isolation, 271–273  
   electron microscope and, 261  
   replication of, 261  
   sexually transmitted chlamydial diseases caused by, 510, 511–512  
 Visualization of bacteria, and staining techniques, 57
- Vitamins:  
   microbe-produced, 344  
   required by microbes, 98  
 Voges-Proskauer test, 154*f*, 171, 174–175, 175*f*, 177*f*, 216*f*  
*Vorticella*, 224*t*
- W**  
 Water:  
   as decolorizing agent, 87  
   nutritional requirements of microbes, 98  
   as source of food contamination, 317  
 Water microbiology, 327–342  
   membrane filter method analysis, 337–340, 339*f*  
   standard bacterial analysis, 329–334, 330*f*, 334*t*  
 Waterbath for moist heat experiment, 281*f*  
 Waterbaths, 2*f*  
   shaking, 4–6  
 Wetting agents, mechanism of action/uses, 311*t*  
 Wine production, 237, 247, 323–324, 323*f*  
 Wire loops, 2*f*, 3–4  
 World Health Organization (WHO), 327
- X**  
 X-radiation to control microbial growth, 278*f*, 289, 289*f*  
 Xylo, 40
- Y**  
 Yeast extract broth, 99, 101  
 Yeasts, 247–254  
   asexual reproduction, 248*f*  
   cells, 248*f*  
   cultivation of, 247–254  
   fermentation, 248–250, 252  
   importance to food industry, 247  
   morphology, 248–249, 251, 255–257*t*, 425*f*  
   pathogenic, 247  
   reproduction, 247, 248*f*, 249, 252  
   reproductive structures of, 248*f*  
   use of pyruvic acid by, 162*f*  
   in wine production, 323, 323*f*
- Z**  
 Ziehl-Neelsen method, 81  
 Zone of antibody excess, 489, 489*f*  
 Zone of antigen excess, 489, 489*f*  
 Zone of inhibition (Kirby-Bauer test), 294–295, 296*f*, 299*t*  
 Zwitterions, 122  
 Zygomycetes, 237, 238*t*  
   soil populations, 345

# LABORATORY SAFETY RULES AND REGULATIONS

The following basic safety rules should be observed at all times in the laboratory.

1. Upon entering the laboratory, place coats, books, and other paraphernalia in specified locations—never on bench tops.
  2. Keep doors and windows closed during the laboratory session to prevent contamination from air currents.
  3. At the beginning and termination of each laboratory session, wipe bench tops with a disinfectant solution provided by the instructor.
  4. Do not place contaminated instruments, such as inoculating loops, needles, and pipettes, on bench tops. Loops and needles should be sterilized by incineration, and pipettes should be disposed of in designated receptacles.
  5. On completion of the laboratory session, place all cultures and materials in the disposal area as designated by the instructor.
  6. Rapid and efficient manipulation of fungal cultures is required to prevent the dissemination of their reproductive spores in the laboratory environment.
  7. Wash your hands with liquid detergent, rinse with 95% ethyl alcohol, and dry them with paper towels upon entering and prior to leaving the laboratory.
  8. Wear a paper cap or tie back long hair to minimize its exposure to open flames.
  9. Wear a laboratory coat/apron and disposable gloves, as determined by your instructor, to protect your clothing and hands from contamination and from coming into direct contact with staining solutions.
  10. Closed shoes should be worn at all times in the laboratory setting.
  11. Never apply cosmetics or insert contact lenses in the laboratory.
  12. Do not smoke, eat, or drink in the laboratory. These activities are absolutely prohibited.
  13. Carry cultures in a test-tube rack when moving around the laboratory. Likewise, keep cultures in a test-tube rack on the bench tops when not in use. This serves a dual purpose: to prevent accidents and to avoid contamination of yourself and the environment.
  14. Never remove media, equipment, or, especially, microbial cultures from the laboratory. Doing so is absolutely prohibited.
  15. Immediately cover spilled cultures or broken culture tubes with paper towels and then saturate them with disinfectant solution. After 15 minutes of reaction time, remove the towels and dispose of them in a manner indicated by the instructor.
  16. Report accidental cuts or burns to the instructor immediately.
  17. Never pipette by mouth any broth cultures or chemical reagents. Doing so is strictly prohibited. Pipetting is to be carried out with the aid of a mechanical pipetting device only.
  18. Do not lick labels. Use only self-stick labels for the identification of experimental cultures.
  19. Speak quietly and avoid unnecessary movement around the laboratory to prevent distractions that may cause accidents.
- Additional precautions for handling body fluids:
20. Disposable gloves must be worn during the manipulation of these test materials.
  21. Immediate hand washing is required if contact with any of these fluids occurs and also upon removal of gloves.
  22. Masks, safety goggles, and laboratory coats should be worn if an aerosol might be formed or splattering of these fluids is likely to occur.
  23. Spilled body fluids should be decontaminated with a 1:10 dilution of household bleach, covered with paper toweling, and allowed to react for 10 minutes before removal.
  24. Test specimens and supplies in contact with these fluids must be placed into a container of disinfectant prior to autoclaving.

# SERIAL DILUTION

This figure illustrates the basic steps of serial dilution and appears here as a convenient reference for this commonly used technique. Turn back to this figure whenever you need a quick reminder of the procedure. For a more detailed discussion of serial dilutions, refer to either Experiment 20 or Appendix 2.

