

**PHILIPS**

*Eagle Eye Platinum*

RX digital IVUS catheter

Fast, plug-and-play  
simplicity



**Eagle Eye  
Platinum ST  
catheter**

**Eagle Eye  
Platinum  
catheter**


Designed for exceptional deliverability and ease of use

- Soft, tapered tip with lowest available entry profile and choice of two lengths<sup>1</sup>
- GlyDx hydrophilic coating
- Long rapid exchange lumen for pushability
- Radial access appropriate; fits through all 5F guides<sup>2</sup>


Radiopaque markers for quick, convenient length estimation

- Three markers not offered by other IVUS catheters
- 10 mm spacing facilitates length estimation without a pullback device or marker wire.


Provides advanced imaging on the CoreVision platform with SyncVision co-registration

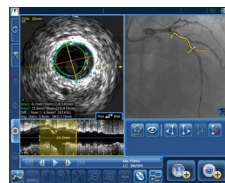


**SyncVision IVUS  
Co-registration**

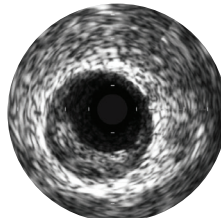

**Grayscale IVUS imaging**

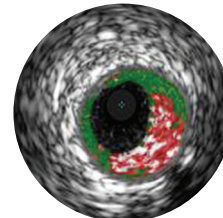


**VH IVUS imaging<sup>3</sup>**

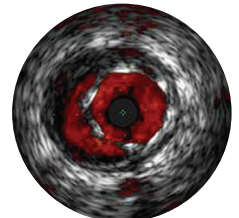


**ChromaFlo imaging**

IVUS is now supported by large randomized and observational studies.<sup>4,5</sup> It is associated with a 74% change in PCI strategy, improved outcomes and quality metrics.<sup>4,5,6,7</sup>

**Eagle Eye Platinum / Platinum ST  
RX digital IVUS catheter**

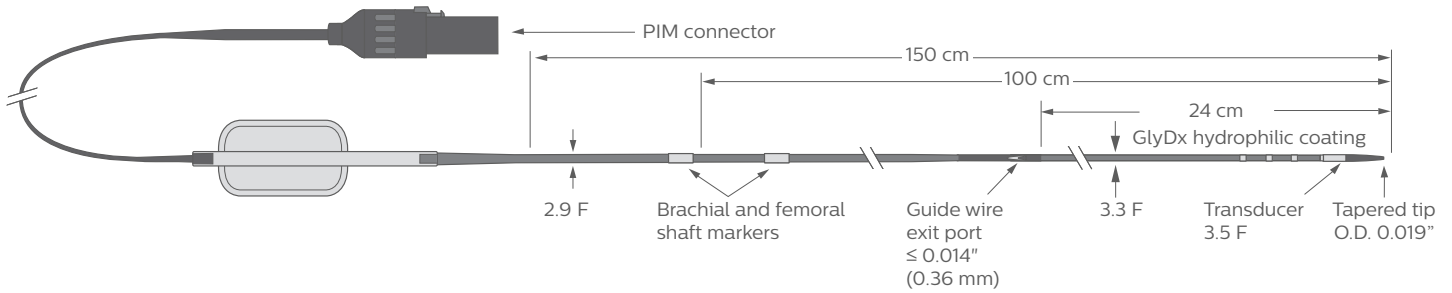

**Distal tip and RO marker detail**

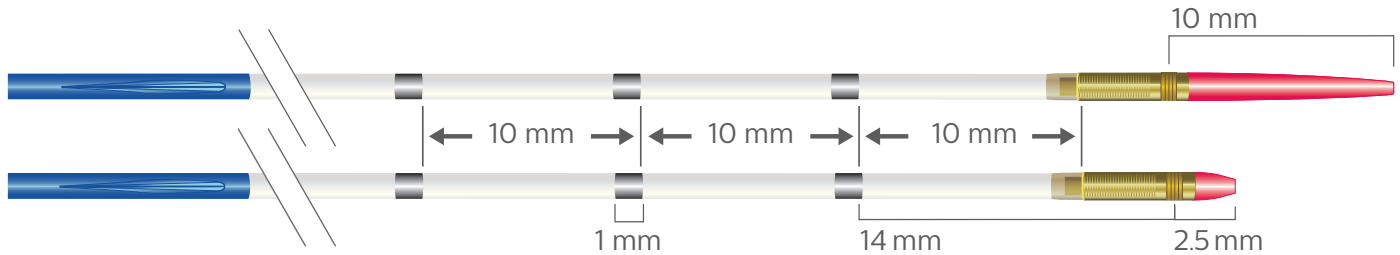

**Technical specifications**

Minimum guide catheter	Maximum guide wire	Maximum field of view	Working length	Frequency	Ordering information
5F (I.D. ≥ .056")	0.014	20 mm	150 cm	20 MHz	85900P Platinum 85900PST Platinum ST

- 0.019" entry profile, data on file; tip to transducer lengths offered include 2.5 mm and 10 mm
- Fits through guide catheters with inner diameters as low as 0.056"; data on file at Philips
- Safety and effectiveness of system for use in the characterization of vascular lesions and tissue types has not been established.
- Elgendy IY et al. Outcomes with intravascular ultrasound-guided stent implantation: a meta-analysis of randomized trials in the era of drug eluting stents. *Circ Cardiovasc Interv.* 2016;9:e003700
- Ahn JM, Kang SJ, Yoon SH, et al. Meta-analysis of outcomes after intravascular ultrasound-guided versus angiography-guided drug-eluting stent implantation in 26,503 patients enrolled in three randomized trials and 14 observational studies. *Am J Cardiol.* 2014;113:1338-1347. Hyperlink "[http://www.ajconline.org/article/S0002-9149\(14\)00549-9/abstract](http://www.ajconline.org/article/S0002-9149(14)00549-9/abstract)"
- Witzenbichler B, et al. Relationship between intravascular ultrasound guidance and clinical outcomes after drug-eluting stents: The ADAPT-DES study. *Circulation.* 2014 Jan; 129,4:463-470.
- Singh V, Badheka AO, Arora S, et al. Comparison of in-hospital mortality, length of hospitalization, costs, and vascular complications of percutaneous coronary interventions guided by ultrasound versus guided by angiography. *Am J Cardiol.* Online 18 Feb 2015.

