

The Custom Hipo Hip Joint

Available in 3 different sizes, each individually custom-made. The Hipo-klein for toddlers progressing from crawling to standing, the Hipo-mittel for children, and the Hipo-gross-Aqua for adult hip disarticulation patients.

The Hipo-gross-Aqua is the first waterproof prosthetic hip joint available for hip disarticulation patients.



Illustration: Hipo-gross-Aqua (CAD Animation)

Special Features

- All Hipo hip joints are **problem free around sand and dirt.**
- The Hipo-gross-Aqua for the **nicest time of year.**
- Individually custom fabricated and suitable for patients with **difficult stump conditions.**
- Ideal as an **alternate prosthesis.**

The System

An individual, distinctive hip joint which can be modified on-site by the attending prosthetist. The hip joint is easily lengthened for patients who are still growing.

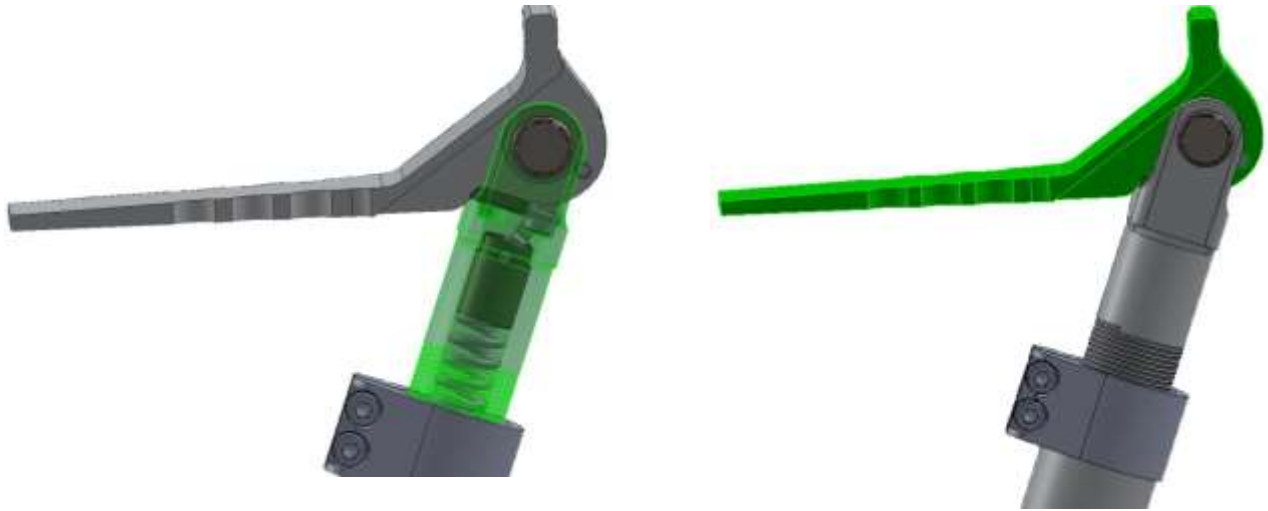


Illustration top left:

Aligned according to the patient with adjustable spring force and rotational.

Illustration top right:

Pelvic connection custom-made according to cast specifications.

Illustration bottom right:

Adapter with neutral alignment angle for your chosen knee joint.



Ordering

The Hipo hip joint is fabricated specifically for your patient. Therefore, it is necessary to send us the cast and knee joint, as well as relevant measurements.

Delivery

Delivery time is approx. 4 -6 weeks after receipt of order. The Hipo hip joint is shipped in the neutral alignment position. Fabrication of the prosthesis and individual adjustments are made on-site by the attending prosthetist.

For further information:

Jeff Matwichuk

and Thomas Böckh

Telephone +49 89 90475750